

Agilent D5000 ScreenTape System Quick Guide

The Agilent 2200 TapeStation system is an automated platform for simpler, faster and more reliable electrophoresis. It is made up of 3 elements:

- 2200 TapeStation system (G2964AA) or 2200 TapeStation Nucleic Acid system (G2965AA)
- D5000 ScreenTape (5067-5588) with D5000 Reagents (5067-5589)
- Agilent Software packages (2200 TapeStation Controller Software, and TapeStation Analysis Software)

Kits

The D5000 ScreenTape system is designed for analysing DNA molecules from 100 – 5000 bp.

Specifications

Analytical Specifications	D5000 ScreenTape assay and reagents
Sizing Range	100 – 5000 bp
Typical Resolution	400 – 5000 bp: 15 %
Sensitivity ¹	0.1 ng/μL
Sizing Precision ²	5 % CV
Sizing Accuracy ^{2,3}	±10 %
Quantification Precision	0.1 – 1 ng/μL: 15 % CV; 1 – 50 ng/μL: 10 % CV
Quantification Accuracy ⁴	±20 %
Quantification Range	0.1 – 50 ng/μL
Maximum Sample Buffer Strength	250 mM KCl, 250 mM Tris-HCl, 125 mM NaCl, 50 mM Acetate, 25 mM MgCl ₂ , 25 mM BSA, 25 mM Guanidine-HCl
Physical Specifications	
Multiple ScreenTape (96 well) compatibility	Not compatible
Analysis Time	16 samples: <25 min
Samples per consumable	1 ladder + 15 samples
Sample volume required	1 μL
Kit stability	4 months
Kit size	112 samples/box

¹ signal-to-noise >3 (single peak)

² Determined using the D5000 Ladder as sample

³ Accuracy for run using single ScreenTape device

⁴ Measured against 2100 Bioanalyzer system



Agilent D5000 ScreenTape System Quick Guide

Storage Conditions

- Reagent vials: 2 – 8 °C (36 – 46 °F)
- The ScreenTape device: 2 – 8 °C (36 – 46 °F) (if you run less than 16 lanes, store used ScreenTape device upright at 2 – 8 °C (36 – 46 °F) for a maximum of 2 weeks.)
- *Never* freeze the ScreenTape device - any ScreenTape device which is accidentally frozen should be discarded.

Kit Components

Part Number	Name	Color	Amount
5067-5588	D5000 ScreenTape		7 ScreenTape devices
5067-5589	D5000 Reagents		2 vials
	• D5000 Ladder	●	10 µL
	• D5000 Sample Buffer	●	1300 µL
5067-5590	D5000 Ladder	●	1 vial
		●	10 µL

Additional Consumables Required for the 2200 TapeStation Instrument

- Loading tips (5067-5152 or 5067-5153)
- Optical Tube 8x Strip (401428) and Optical Cap 8x Strip (401425).
- Vortex mixer (See note below)

Additional Material Required (Not Supplied)

- Volumetric pipette
- Centrifuge

NOTE

2200 TapeStation instruments are supplied with an optional IKA MS3 vortexer which includes a 96-well plate adaptor suitable for 8-way strips.

Safety Information

WARNING

Toxic agents

The handling of solvents, samples and reagents can hold health and safety risks.

- When using/handling the ScreenTape device and working with these substances observe appropriate safety procedures (for example by wearing goggles, safety gloves and protective clothing).
- Always follow good laboratory practices and adhere to the guidelines established in your laboratory.
- Refer to product material safety datasheets for further information.
- The volume of substances should be reduced to the minimum required for the analysis.

CAUTION

Damage to the 2200 TapeStation instrument

→ Use only the recommended consumables and reagents with the 2200 TapeStation system.

General Information on Working with DNA

NOTE

- For best results ensure that all reagents are allowed to equilibrate to room temperature for 30 min prior to use.
- When pipetting Sample Buffer, ensure that excess buffer droplets are removed from the tip before transfer to the sample tubes. Care must be taken due to the viscosity of Sample Buffer.
- When pipetting small volumes ensure that no sample remains within the tip.
- When adding sample buffer to sample, please ensure that they are mixed correctly. To achieve this, gently mix several times with additional pipetting, then cap the tubes, vortex mix using IKA vortexer and adaptor at 2000 rpm for 1 min.
- Briefly centrifuge to collect the contents at the base of the tubes.
- *Improper mixing can lead to quantification errors.*

Essential Measurement Practices

Environmental conditions	<ul style="list-style-type: none"> • Optimal operating temperature: 20 °C (68 °F) • Ambient operating temperature: 15 - 30 °C (59 - 86 °F)
Steps before use on the TapeStation instrument	<ul style="list-style-type: none"> • Equilibrate each vial to room temperature. • Vortex mix each vial and briefly spin. • 'Flick' ScreenTape device to eliminate bubbles in the separation channel, which could interfere with sample loading.
Steps during sample preparation	<ul style="list-style-type: none"> • Keep reagents at room temperature during sample preparation.
Storage after use on the TapeStation instrument	<ul style="list-style-type: none"> • Store all reagent vials and ScreenTape devices at 2 – 8 °C (36 – 46 °F) • Never store reagent vials or ScreenTape devices at room temperature or below 0 °C (32 °F). • If you run less than 16 lanes, store used ScreenTape device upright at 2 – 8 °C (36 – 46 °F) for a maximum of 2 weeks.
Pipette carefully	<ul style="list-style-type: none"> • Always pipette reagents against the side of the sample tube. • If using a standard pipette ensure that no residual material is left on the outside of the tip.
Mix properly after each pipetting step	<ul style="list-style-type: none"> • Mix = Vortex the PCR tubes or 96-well plate using Agilent approved IKA vortexer and adaptor at 2000 rpm for 1 min. • Spin = Move the samples to the bottom of the tubes/wells by pulsing in a centrifuge.

Agilent D5000 ScreenTape System Quick Guide

Prepare TapeStation System D5000

Parts required	p/n	Description
	5067-5588	D5000 ScreenTape

- 1 Launch the 2200 TapeStation Controller Software.
- 2 Load single D5000 ScreenTape device and loading tips into the 2200 TapeStation instrument

Sample Preparation D5000 ScreenTape Assay


Parts required	p/n	Description
	5067-5589	D5000 Reagents

- 1 Allow reagents to equilibrate at room temperature for 30 min
- 2 Vortex mix before use
- 3 Prepare ladder by mixing 10 μ L D5000 Sample Buffer (●) with 1 μ L D5000 Ladder (●) in a tube strip.

NOTE

Use a fresh ladder for each run. No electronic ladder is available for the D5000 assay.

- 4 Prepare sample by mixing 10 μ L D5000 Sample Buffer (●) with 1 μ L DNA sample in the same tube strip.
- 5 Spin down, then vortex using IKA vortexer and adaptor at 2000 rpm for 1 min
- 6 Spin down to position the sample at the bottom of the tube.

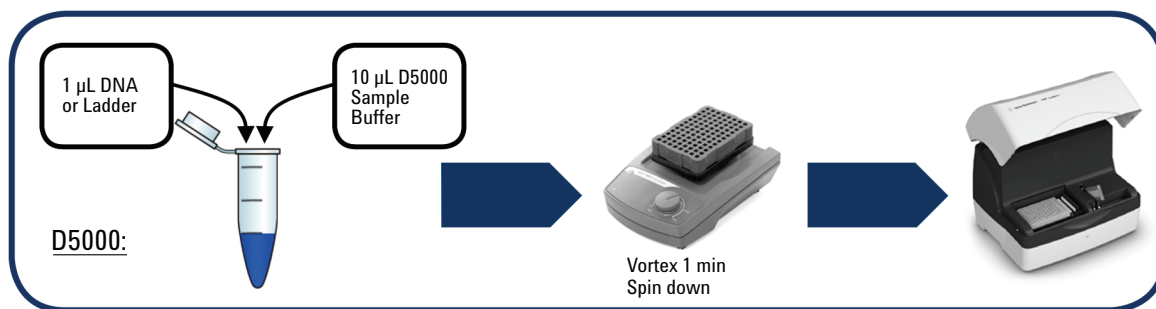

Sample Analysis

- 1 Load samples into the 2200 TapeStation instrument.
- 2 Select the required samples on the 2200 TapeStation Controller Software.
- 3 Click **Start** and specify a filename with which to save your results.

Limited Use Label License

Some products within this system contain SYBR® Green I, which is licensed from Molecular Probes Inc. for use in research and development only. SYBR® is a registered trademark of Molecular Probes, Inc.

For Research Use Only

Not for use in Diagnostic Procedures.

Technical Support

For technical support, please visit www.agilent.com/genomics/contact

Further Information

Visit Agilent Technologies' web site. It offers useful information, support and current developments about the products and technology: www.agilent.com/genomics/tapestation



G2964-90050

Part Number: G2964-90050 Rev. B

Edition 09/2015

Printed in Germany

© Agilent Technologies, Inc. 2014, 2015

Agilent Technologies

Hewlett-Packard-Straße 8

76337 Waldbronn, Germany