




LC AND LC/MS

Your Essential Resource for Columns & Supplies

The Measure of Confidence

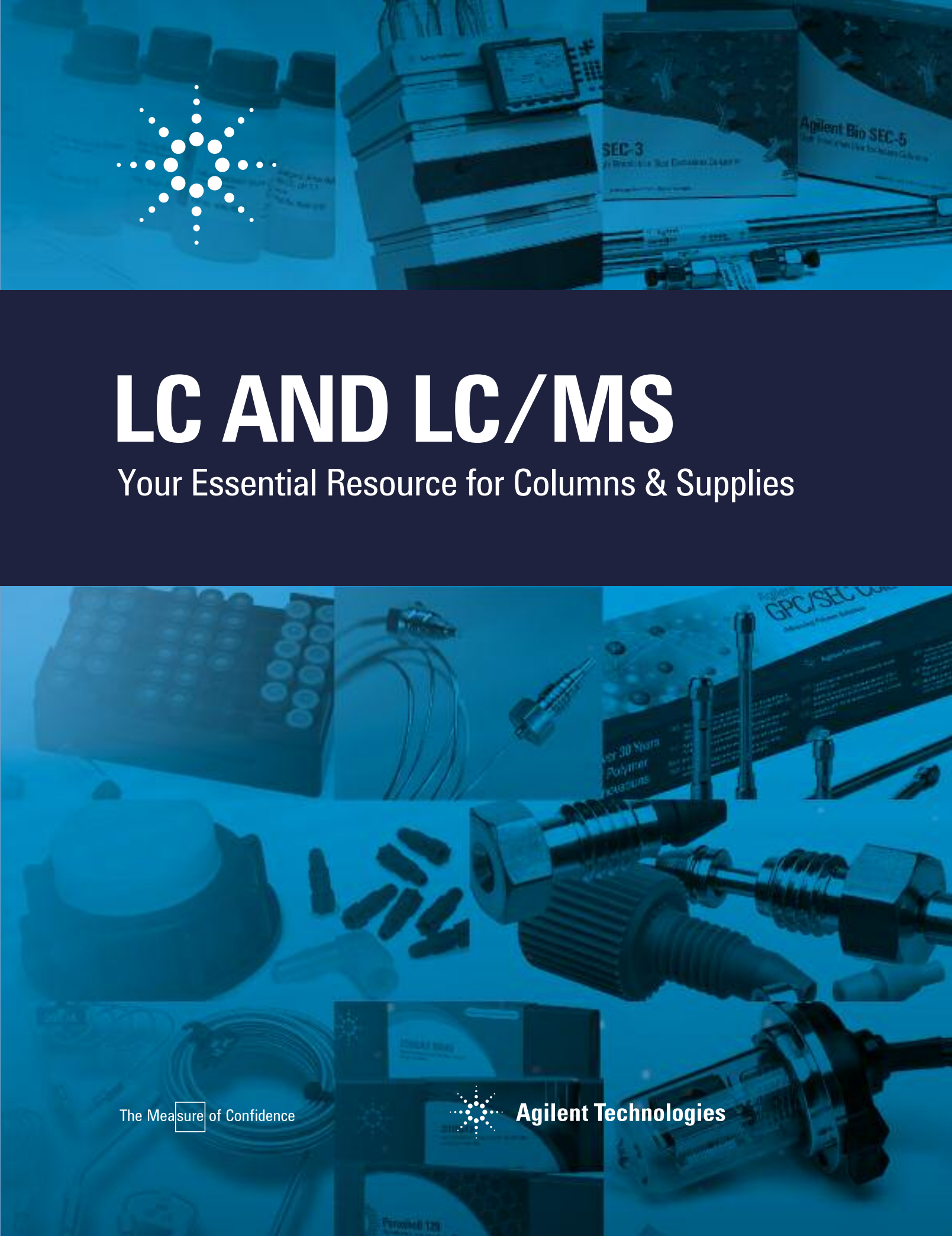

Agilent Technologies

LC and LC/MS Troubleshooting

HPLC Troubleshooting

Symptom Type	Possible Cause	Solution
Baseline disturbance at void time	Positive/negative – Difference in refractive index of injection solvent	Use mobile phase for sample solvent
Detector leaks	Plugged inlet frit	Replace seals/gaskets
Drifting baseline	Positive direction – Contaminant buildup/elution	Flush column, clean up sample, use pure solvents
	Positive/negative – Difference in refractive index of injection solvent	Use mobile phase for sample solvent
	Negative direction (gradient) – Absorbance of "A" mobile phase solvent	Use non-absorbing or HPLC-grade or better solvent
	Positive direction (gradient) – Absorbance of "B" mobile phase solvent	Use non-absorbing or HPLC-grade or better solvent
	Random – Temperature changes	Insulate column and tubing
	Random – Temperature changes	Thermostat column and tubing
	Wavy or undulating – Temperature changes in room	Monitor room temperature and control
Ghost peaks	Peaks from previous injection	Flush column to remove contaminants
	Contamination	Sample cleanup or pre-fractionation
	Unknown interferences in samples	Sample cleanup or pre-fractionation
	Ion-pair – Upset equilibrium	Prepare sample in actual mobile phase to minimize disturbance
	Peptide mapping – Oxidation of TFA	Prepare fresh daily; use anti-oxidant
	Reversed-phase – Contaminated water	Check suitability of water by running different amount through reversed-phase column and measure peak height with elution; use HPLC grade solvents
	Spikes – Bubbles in solvent	De-gas solvents
High column backpressure	Column blockage, adsorbed sample	Better sample cleanup; use guard column
	Mobile phase viscosity too high	Use lower viscosity solvents or higher temperature
	Particle size too small	Use larger d_p packing
	Plugged inlet frit	Replace column
	Plugged inlet frit	Reverse solvent flow
Leaks	Subtle – White powder at fitting/loose fitting	Tighten fittings, cut tubing, or replace ferrules
Leaks, injection valve	Catastrophic – Worn valve rotor	Replace rotor in valve
Leaks, column or other fittings	Catastrophic – Loose fittings	Tighten or replace fittings
Leak, pump	Catastrophic – Pump seal failure	Replace pump seal

(Continued)

HPLC Troubleshooting

Symptom Type	Possible Cause	Solution
Negative peaks	RI detector – solute refractive index less than solvent	No problem; reverse polarity to make positive
	UV detector – solute absorbance less than mobile phase	Use mobile phase with lower UV absorbance; do not recycle solvent too long
Noisy baseline	Random – Contaminant buildup	Flush column; clean up sample; use HPLC-grade solvent
	Continuous – Detector lamp problem	Replace detector lamp
	Occasional – External electrical interference	Use voltage stabilizer for LC system
Peak doubling	Sample volume too large	Reduce the volume e.g. by half and re-inject
	Injection solvent too strong	Use weaker injection solvent or mobile phase
	Blocked frit	Replace and use 0.5 µm porosity in-line filter
	Column void or channeling	Replace column; for some columns, fill in void with packing
	Unswep injector flowpath	Replace injector rotor
	Void at head of column	Replace column, top off column with packing
	Column overloaded with sample	Use higher capacity stationary phase
		Increase column diameter Decrease sample size
Peak tailing	Single peak – interfering components	Sample cleanup; pre-fractionation
	Beginning of peak doubling	See "peak doubling"
	Unswep dead volumes	Minimize number of connections
		Ensure injector seal is tight
		Ensure fittings are properly seated
	Basic compounds – Silanol interactions	Choose endcapped bonded phase Switch to polymeric phase
	Basic substances – Silanol interactions	Use stronger mobile phase or add competing base (e.g. TMA)
	Silica-based – Column degradation	Use specialty column; polymeric column or sterically protected

(Continued)

HPLC Troubleshooting

Symptom Type	Possible Cause	Solution
Peaks are broad	Injection volume too large	Decrease solvent strength of injection solvent to focus solute
	Peak dispersion in injector valve	Introduce air bubble in front/back of sample to decrease dispersion
	Sampling rate of data system too slow	Increase frequency of sampling
	Slow detector time constant	Adjust time constant to match peak width
	Mobile phase viscosity too high	Increase column temperature
	Detector cell volume too large	Use smallest possible cell volume with no heat exchanger in system
	Injector volume too large	Decrease injection volume
	Long retention times	Use gradient elution or stronger mobile phase
Pressure fluctuation	Leaky check valve	Replace check valve
	Pump seal leaks	Replace pump seals
	Buildup of particulates	Filter sample; in-line filter; filter mobile phase
Pressure increasing	Buildup of particulates	Filter sample; in-line filter; filter mobile phase
	Water/organic systems – buffer precipitation	Test buffer-organic mixtures; ensure compatibility
Retention beyond total permeation volume	Size exclusion – Specific interactions	Add mobile phase modifiers or change solvent
Retention times changing	Column temperature varying	Thermostat column; insulate column; ensure lab temperature constant
	Equilibration time insufficient with gradient run or changes in isocratic mobile phase	Make sure at least 10 column volumes pass through column after solvent change or gradient conclusion
	Selective evaporation of mobile phase component	Less vigorous helium sparging; keep solvent reservoirs covered; prepare fresh mobile phase
	Buffer capacity insufficient	Use >20 mM concentration of buffer
	Inconsistent on-line mobile phase mixing	Ensure gradient system delivering constant composition; check vs. manual prep of mobile phase
	Contamination buildup	Occasionally flush column with strong solvent to remove contaminants
	First few injections – Adsorption on active sites	Condition column by initial injection of concentrated sample

(Continued)

HPLC Troubleshooting

Symptom Type	Possible Cause	Solution
Retention times decreasing	Flow rate increasing	Check pump to make sure correct; if not, reset
	Column overloaded with sample	Decrease sample size
	Loss of bonded stationary phase	Keep mobile phase pH between 2 and 8.5
Retention times increasing	Flow rate is slowing	Fix leaks in liquid lines, replace pump seals, check for pump cavitation or air bubbles
	Active sites on silica packing	Use mobile phase modifier
	Loss of bonded stationary phase	Keep mobile phase pH between 2 and 8.5
	Mobile phase composition changing	Make sure mobile phase container is covered
	Active sites on silica packing	Add competing base to mobile phase
	Active sites on silica packing	Use higher coverage packing for stationary phase
Sensitivity problem	Peaks are outside of linear range of detector	Dilute/concentrate to bring into linear region
	First few sample injections – Absorption of sample in loop or column	Condition loop/column with concentrated sample
	Autosampler flow lines blocked	Check flow and make sure there are no blockages
	Injector sample loop underfilled	Make sure that loop is overfilled with sample
	Sample-related losses during preparation	Use internal standard during sample prep; optimize sample prep method
Slow column equilibration times (ion-pairing)	Equilibration time slow for long-chain ion-pairing reagents	Use shorter alkyl chain ion-pair reagent

LC/MS Troubleshooting

Symptom Type	Solution
No peaks	Spray from the nebulizer
	Make sure capillary voltage is set correctly
	Make sure LC/MSD is tuned correctly
	Make sure LC/MSD pressures are within normal ranges
	Check drying gas flow and temperature
	Make sure fragmentor is set correctly
Poor mass accuracy	Recalibrate the mass axis
	Make sure ions used for tuning span mass range of sample ions and show strong stable signals
Low signal	Check the solution chemistry; make sure solvent is appropriate for sample
	Make sure sample is fresh and has been stored correctly
	Make sure LC/MSD is tuned correctly
	Check the nebulizer condition
	Clean the capillary entrance
	Check the capillary for damage and contamination
Unstable signal	Make sure drying gas flow and temperature are correct for the solvent flow
	Make sure solvent is thoroughly degassed
	Make sure LC backpressure is steady; this indicates a steady solvent flow

(Continued)

LC/MS Troubleshooting

Symptom Type	Solution
High spectral noise	Use appropriate mass filter values
	Check spray shape; nebulizer may be damaged or set incorrectly
	Make sure drying gas flow and temperature are correct for the solvent flow
	Make sure solvent is thoroughly degassed
	Make sure LC backpressure is steady; this indicates a steady solvent flow
	If you are using water as part of the mobile phase, make sure it is de-ionized ($> 18 \text{ M}\Omega \text{ cm}$)
Droplets, not spray, exiting the nebulizer	Make sure nebulizing gas pressure is set high enough for the LC flow
	Check position of needle in nebulizer
	Stop solvent flow and remove nebulizer assembly
	Examine end of nebulizer for damage
No flow	Make sure LC is on and there is sufficient solvent in correct bottle
	Check for LC error messages
	Check for blockages
	Repair or replace any blocked components
	Check for leaks
	Make sure MS stream selector valve is set to LC to MSD
Undesired fragmentation	(APCI vs. Electrospray)
	APCI temperature is too high
	Fragmentor voltage is set too high

BioPharmaceutical Applications

NEW!

Protein digest analysis

Column: ZORBAX 300SB-C18
858750-902
2.1 x 100 mm, 1.8 μ m

Mobile Phase: A: 0.1% TFA in water
B: 0.085% TFA in ACN

Flow Rate: 0.5 mL/min

Pressure: 640 bar

Gradient: 2% B 1 min, 2-45% B 8.8 min,
45-95% B 0.2 min, 95% B 2 min,
98-2% B 0.2 min, 2% B 1.8 min


Temperature: 50 °C

Detector: Agilent 1290 Infinity LC

Injection: 5 μ L

Sample: Protein digest

Sample Conc: 1 mg/mL

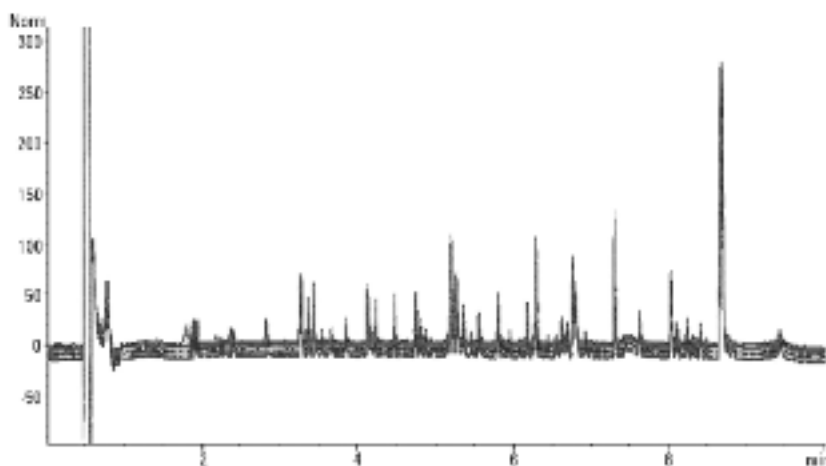

Overlaid chromatograms of 30 runs of a protein digest on an Agilent ZORBAX RRHD 300SB-C18 column.

NEW!**Analysis of oxidized insulin chains**

Column: ZORBAX RRHD 300SB-C18
857750-902

2.1 x 50 mm, 1.8 µm

Mobile Phase: A: 0.1% TFA in water
B: 80% ACN + 0.01% TFA in water

Flow Rate: 1.0 mL/min

Pressure: 650-700 bar


Gradient: 33-50% B, 0-4 min; 33% B, 4-5 min

Detector: UV, 280 nm
Agilent 1290 Infinity LC

Sample: Insulin, oxidized insulin chain A and chain B from bovine pancreas (Sigma Aldrich, St. Louis, MO)

Sample Conc: 1 mg/mL

Injection: 2 µL

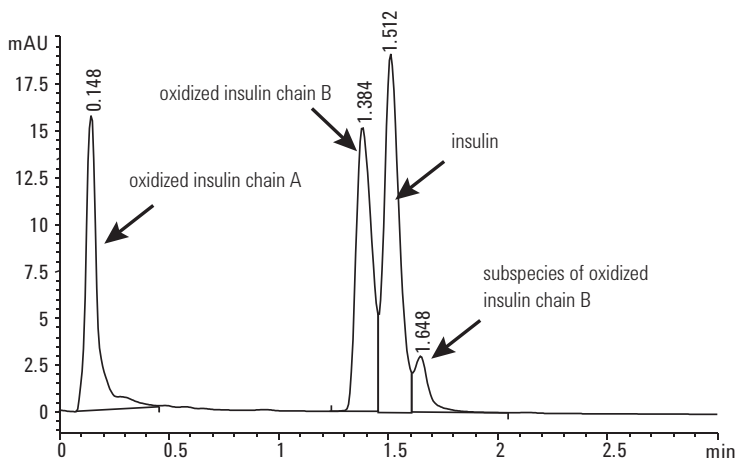

Insulin and oxidized insulin A and B chains are resolved quickly but insulin and oxidized chain B often co-elute.

NEW!**Fast separation of recombinant human erythropoietin**

Column: ZORBAX RRHD 300SB-C18
857750-902

2.1 x 50 mm, 1.8 µm

Mobile Phase: A: 0.1% TFA in water
B: 0.01% TFA in ACN

Flow Rate: 1.0 mL/min

Pressure: 650 bar


Gradient: 5 to 100% B solvent from 0 to 2.5 min

Detector: UV, 280 nm
Agilent 1290 Infinity LC

Sample: Recombinant human EPO protein (rEPO)

Sample Conc: 1.0 mg/mL

Injection: 3 µL

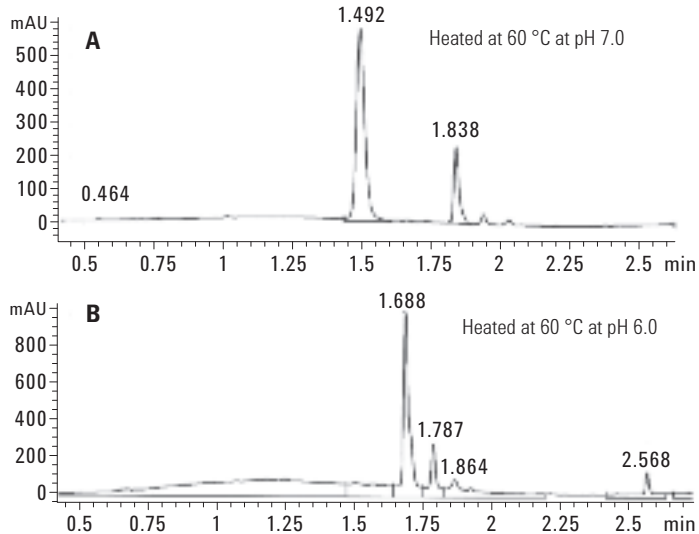

Heat-treated rEPO protein are well resolved by the Agilent ZORBAX RRHD 300SB-C18 column. The column separated these heat-treated rEPO proteins.

NEW!**Separation optimization for ultra fast analysis of reduced and alkylated monoclonal antibody**

Column: ZORBAX RRHD 300SB-C8
858750-906
2.1 x 100 mm, 1.8 μ m

Mobile Phase: (Various)
A: H₂O + 0.1% TFA (v/v)
B: n-propanol:ACN:H₂O (80:10:10) + 0.1% TFA (v/v)

Injection: 1-3 μ L



Flow Rate: 0.5 mL/min

Gradient: Multi-segmented
A (optimized for speed): 0 min-20% B, 3 min-35% B, 4 min-40% B, 5 min-40% B, 5.1 min-90% B, 5.5 min-90% B, 6 min-25% B
B (optimized for resolution): 0 min-25% B, 15 min-32% B, 16 min-32% B, 17 min-90% B, 17.5 min-90% B, 18 min-25% B

Temperature: 75 °C

Detector: UV, 225 nm
Agilent 1290 Infinity LC

For consecutive chromatographic runs, a 2-minute post run was added to re-equilibrate the column.

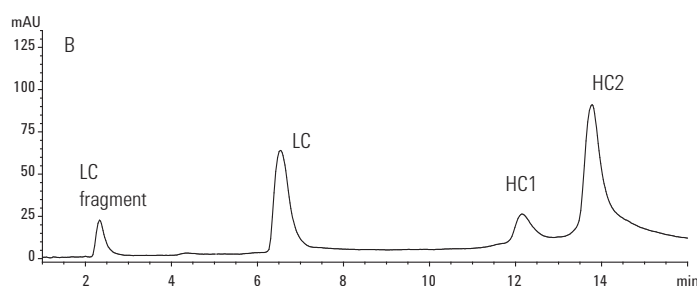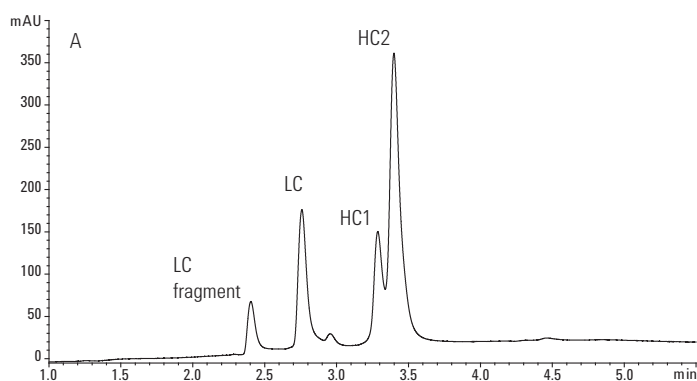

Comparison of two optimized gradients for the ultra fast separation of reduced and alkylated monoclonal antibodies on an Agilent ZORBAX RRHD 300SB-C8 column. The top panel details a rapid separation of the light and heavy chain variants in a shortened run time of less than 4 minutes. The bottom panel displays complete baseline resolution of the two heavy chain variants during a longer runtime using a shallower gradient profile. Both separations were performed at 75 °C and completed with a fast 90% 1-propanol wash step (UV not shown).

NEW!**Column reproducibility – 200 injections of reduced monoclonal antibody using an Agilent ZORBAX RRHD 300SB-C3 column**

Column: Agilent ZORBAX RRHD 300SB-C3
858750-909
2.1 x 100 mm, 1.8 μ m

Temperature: 75 °C
Detector: UV, 280
Agilent 1290 Infinity LC


Mobile Phase: A: 0.1% TFA in water
B: 80% n-propyl alcohol, 10% ACN,
9.9% water and 0.1% TFA

Sample: Reduced monoclonal antibody (IgG1) (1.0 mg/mL) -
Agilent BL05 IgG1

Flow Rate: 0.4 mL/min

Injection: 2 μ L

Gradient: 0 min-1% B, 2 min-20% B, 5 min-50% B,
7 min-50% B, 8.0 min-90% B,
8.3 min-1% hold for 2 min

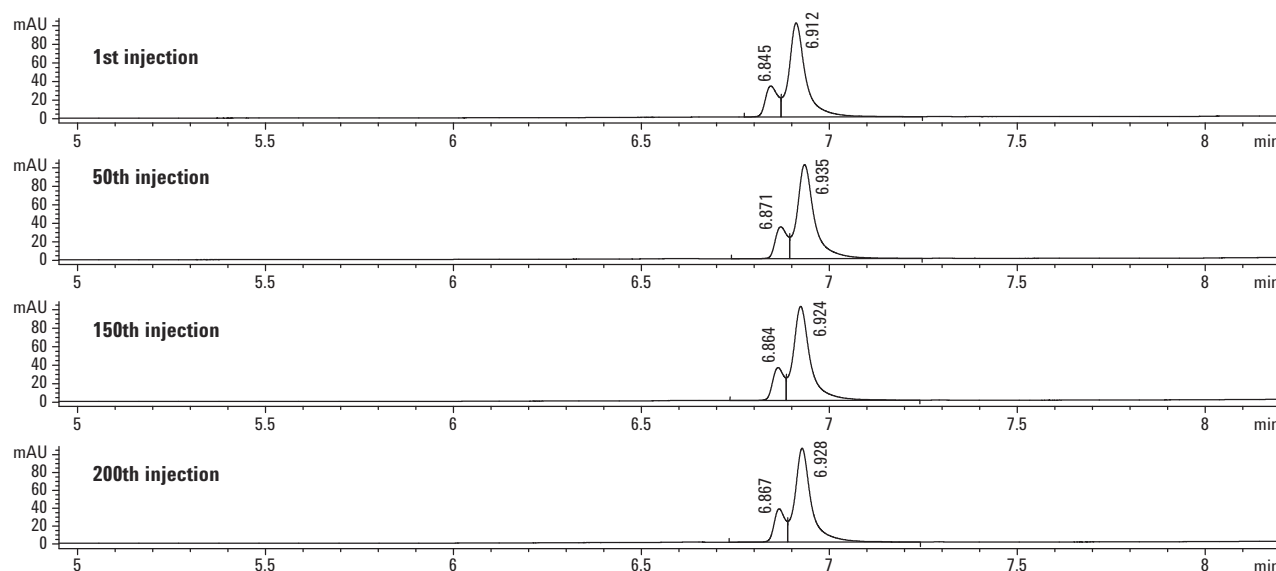

Reduced and alkylated mAb profiling during 200 repeated injections.



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

NEW!**Gradient optimizations for ultra fast analysis of reduced monoclonal antibody**

Column: Agilent ZORBAX RRHD 300SB-Diphenyl 858750-944 2.1 x 100 mm, 1.8 μ m

Mobile Phase: A: 0.1% TFA in water
B: 80% propyl alcohol, 10% ACN, 9.9% water and 0.1% TFA

Flow Rate: 0.5 mL/min


Gradient: 1st condition: 0 min-1% B, 2 min-20% B, 5 min-70% B
2nd condition: 0 min-1% B, 2 min-20% B, 5 min-50% B

Temperature: 74 °C

Detector: UV, 280 nm

Sample: Reduced monoclonal antibody (IgG1) (1.0 mg/mL) - BioCreative IgG1

Injection: 2 μ L

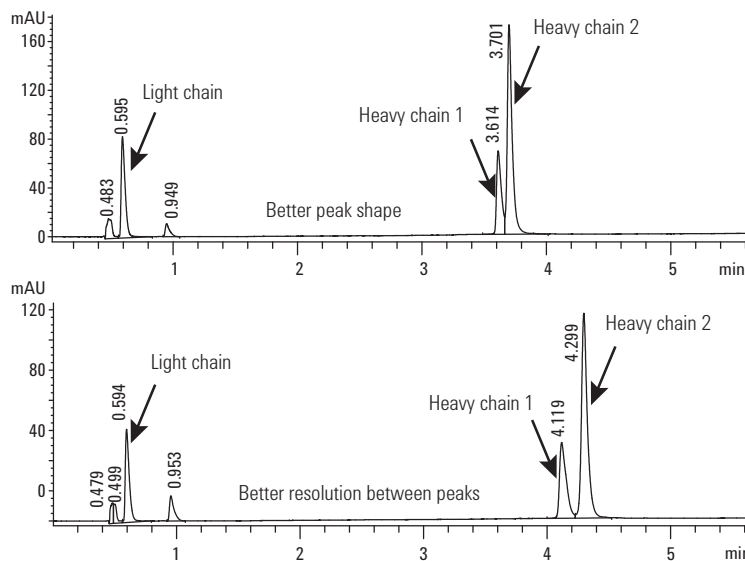

Comparison of two ultra-fast separations of reduced monoclonal antibodies was achieved on a Agilent ZORBAX RRHD 300SB-Diphenyl under different optimized conditions. The top panel separation delivered narrow peak widths with shorter retention times. The bottom panel separation displays higher resolution between the two heavy chain peaks, but with less efficiency.



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

NEW!**Ultra high speed and high resolution
of intact monoclonal antibodies**

Column: Agilent ZORBAX RRHD 300-Diphenyl
858750-944
2.1 x 100 mm, 1.8 μ m

Mobile Phase: A: 0.1% TFA in water
B: 80% n-propyl alcohol,
10% ACN,
9.9% water and 0.1% TFA


Flow Rate: 1.0 mL/min

Temperature: 74 °C

Detector: UV, 280 nm

Sample: Monoclonal antibody (IgG1) (1.0 mg/mL) -
BioCreative IgG1 and Agilent Standard IgG1

Injection: 1 μ L

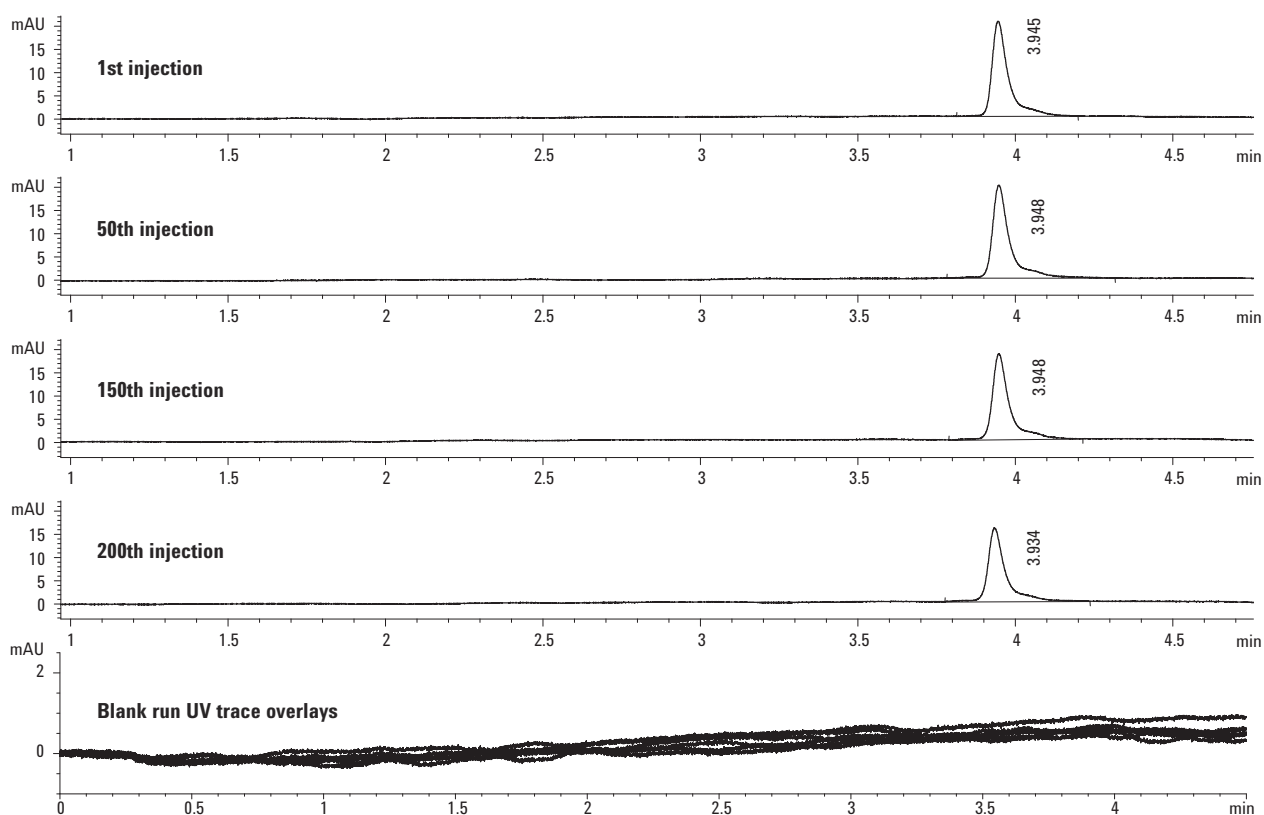

Details of intact mAb profiling during 200 repeated injections. Intact mAb separations shown were collected at 1, 50, 150, and 200th run intervals. The bottom panel displays 5 UV blank run trace overlays collected every 20th run during the column evaluation (**note:** overlay traces are scaled to 2 mAU).

NEW!**Optimizing protein separations with
Agilent weak cation-exchange columns**

Column: Agilent Bio WCX, stainless steel
5190-2453
4.6 x 250 mm, 10 μ m

Column: Agilent Bio WCX, stainless steel
5190-2445
4.6 x 250 mm, 5 μ m

Mobile Phase: A: water
B: 1.6 M NaCl
C: 40.0 mM NaH_2PO_4
D: 40.0 mM Na_2HPO_4
By combining predetermined proportions of C and D, 20 mM buffer solutions at the desired pH range were produced (proportions determined using Buffer Advisor software)

Flow Rate: 1.0 mL/min


Gradient: 0 to 50% B, 0 to 20 min
50% B, 20 to 25 min
0% B, 25 to 35 min

Temperature: Ambient

Detector: UV, 220 nm
Agilent 1260 Infinity Bio-inert Quaternary LC

Sample: Ovalbumin, Ribonuclease A, Cytochrome c, Lysozyme

Sample Conc: 2 mg/mL (in 20 mM sodium phosphate buffer, pH 6.0)

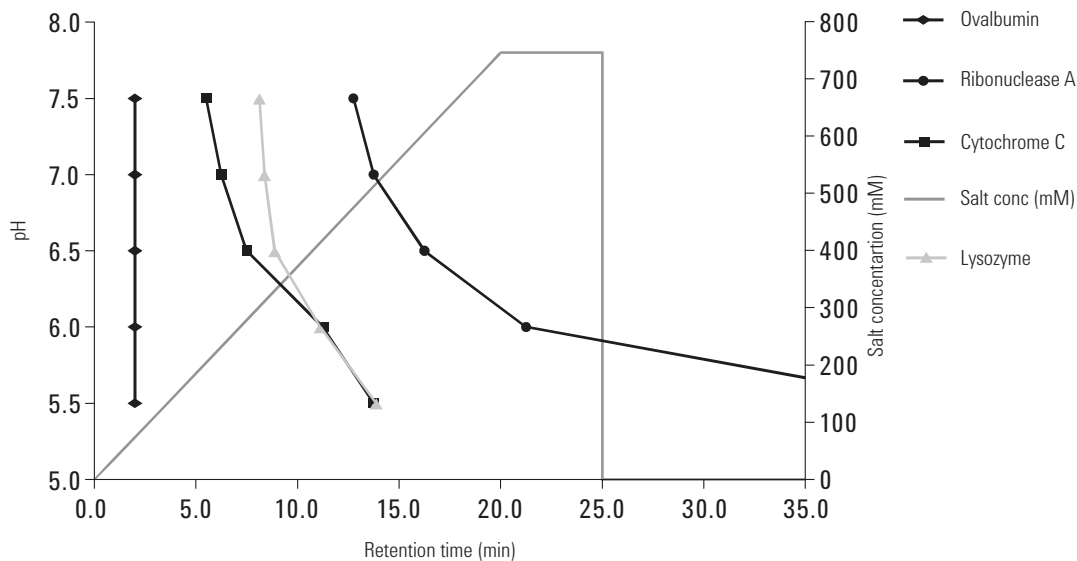

Effect of pH on retention time of protein standards using a Agilent Bio WCX column.

NEW!**Improved resolution with smaller particle size with Agilent weak cation-exchange columns**

Column: Agilent Bio WCX, stainless steel
5190-2453
4.6 x 250 mm, 10 µm

Column: Agilent Bio WCX, stainless steel
5190-2445
4.6 x 250 mm, 5 µm

Mobile Phase: A: water
 B: 1.6 M NaCl
 C: 40.0 mM NaH₂PO₄
 D: 40.0 mM Na₂HPO₄
 By combining predetermined proportions of C and D, 20 mM buffer solutions at the desired pH range were produced (proportions determined using Buffer Advisor software)


Gradient: 0 to 50% B, 0 to 20 min
 50% B, 20 to 25 min
 0% B, 25 to 35 min

Temperature: Ambient

Detector: UV, 220 nm
 Agilent 1260 Infinity Bio-inert Quaternary LC


Sample: Ovalbumin, Ribonuclease A, Cytochrome c, Lysozyme

Sample Conc: 2 mg/mL (in 20 mM sodium phosphate buffer, pH 6.0)

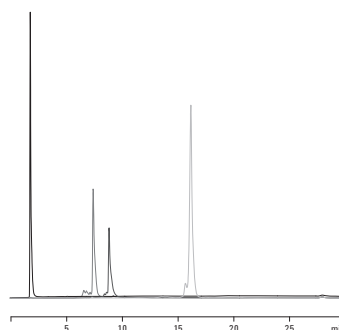

Separation of protein standards at pH 6.5 using an Agilent Bio WCX, NP10 column.


1. Ovalbumin
2. Ribonuclease A
3. Cytochrome c
4. Lysozyme

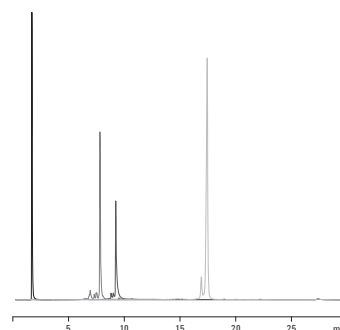

Separation of protein standards at pH 6.5 using an Agilent Bio WCX, NP5 column.



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

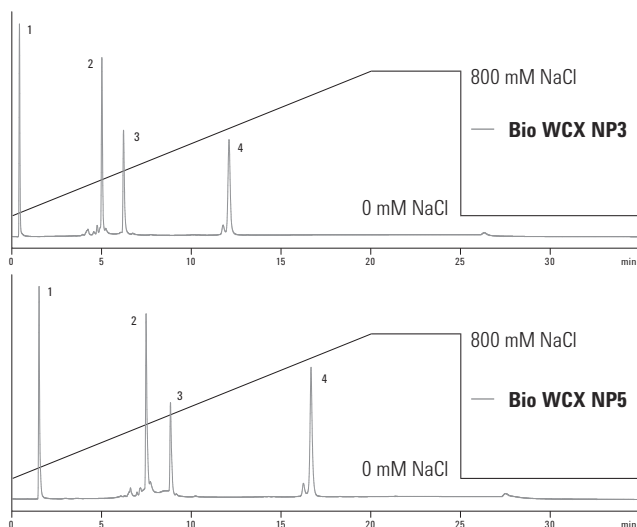NEW!
Faster separations using Agilent weak cation-exchange columns


Protein separation on Agilent Bio WCX NP5 versus Agilent Bio WCX NP3.

Column: Agilent Bio WCX, stainless steel
5190-2445
4.6 x 250 mm, 5 μ m

Column: Agilent Bio WCX, stainless steel
5190-2443
4.6 x 50 mm, 3 μ m

Column: Agilent Bio WCX, stainless steel
5190-2441
4.6 x 50 mm, 1.7 μ m

Mobile Phase: A: 20 mM sodium phosphate, pH 6.5
B: A + 1.6 M NaCl

Flow Rate: 1.0 mL/min

Gradient: 0 to 50% B


Temperature: Ambient

Detector: UV, 220 nm
Agilent 1260 Infinity Bio-inert Quaternary LC



Sample: Ovalbumin, Ribonuclease A, Cytochrome c, Lysozyme

Sample Conc: 0.5 mg/mL

1. Ovalbumin
2. Ribonuclease A
3. Cytochrome c
4. Lysozyme

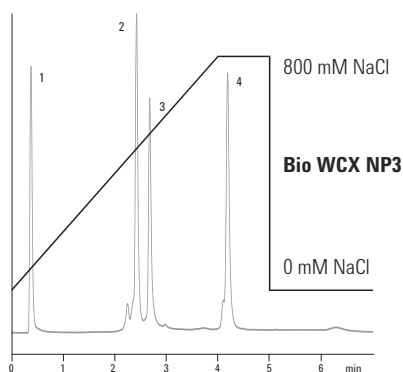

Comparison of Agilent Bio WCX NP3 versus Agilent Bio WCX NP1.7 (flow rate 1.0 mL/min).

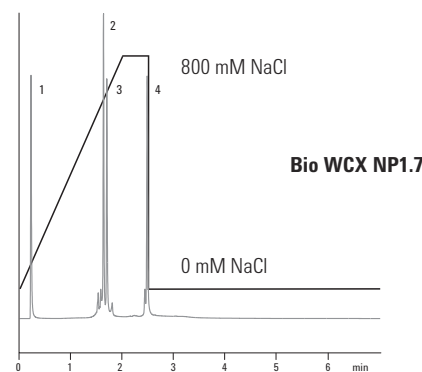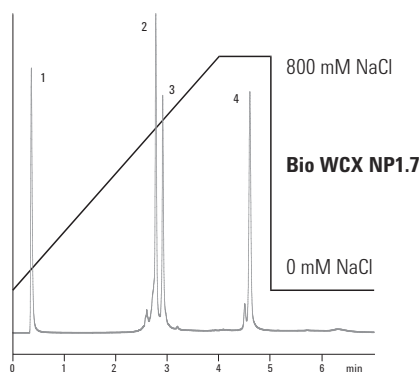

Agilent Bio WCX NP1.7 for protein separations under 3 minutes (flow rate 1.7 mL/min).

NEW!**pH gradient elution for improved separation of monoclonal antibody charged variants**

Column: Bio MAb, stainless steel
5190-2405
4.6 x 250 mm, 5 μ m

Mobile Phase: A: water
B: 1.6 M NaCl
C: 100 mM NaH_2PO_4
D: 100 mM Na_2HPO_4
By combining predetermined proportions of C and D, buffer solutions at the desired pH range were produced at the selected buffer strengths.

Flow Rate: 1.0 mL/min


Gradient: pH 6.0 to 8.0, 0 to 20 minutes
0 to 800 mM NaCl, 20 to 25 minutes
800 mM NaCl, 25 to 30 minutes

Temperature: Ambient

Detector: UV, 220 nm
Agilent 1260 Infinity Bio-inert Quaternary LC

Sample: IgG monoclonal antibody

Sample Conc: 2 mg/mL (in 20 mM sodium phosphate buffer, pH 6.0)

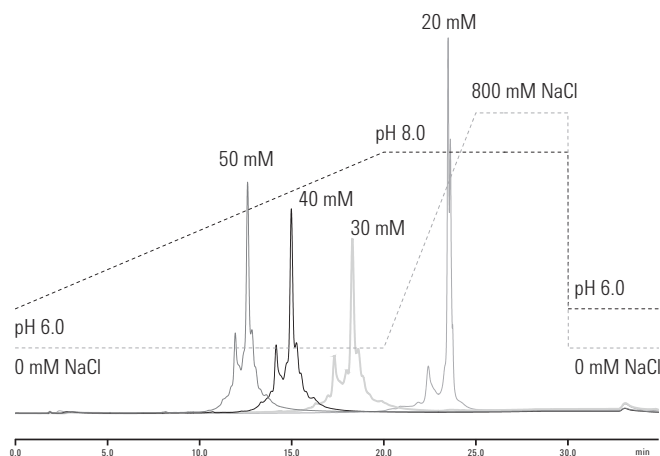

Chromatograms of IgG monoclonal antibody at different ionic strengths.



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

NEW!

Separation of recombinant human erythropoietin (rEPO) using Agilent Bio SEC-3

Column: Bio SEC-3, 100Å
5190-2503
4.6 x 300 mm, 3 µm


Mobile Phase: 150 mM sodium phosphate buffer, pH 7.0

Flow Rate: 0.35 mL/min

Detector: UV, 225 nm
Agilent 1260 Infinity Bio-inert Quaternary LC

Sample: Recombinant human EPO protein (rEPO)

Sample Conc: 1.0 mg/mL

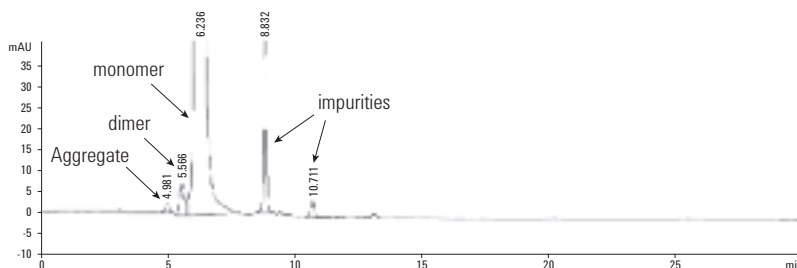

Consistent ion-exchange MAb separation

Column: Bio MAb, PEEK
5190-2411
2.1 x 250 mm, 5 µm

Buffer: A: Sodium phosphate buffer, 20 mM
B: Buffer A + 400 mM NaCl

Gradient: 15-35% Buffer B from 0-30 min


Flow Rate: 0.65 mL/min

Sample: CHO-humanized MAb, 1 mg/mL

Injection: 2.5 µL

Detector: UV, 220 nm

Temperature: Ambient

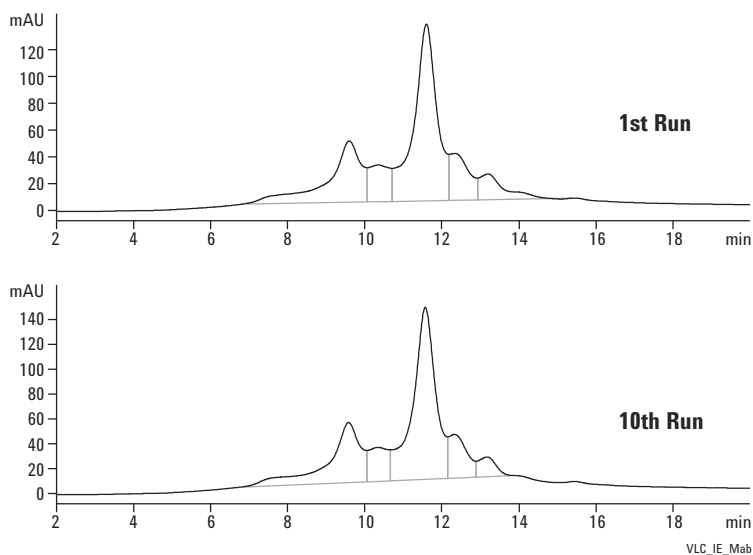

Intact MAb monomer and dimer separation

Column: Bio SEC-3, 300Å
5190-2511
7.8 x 300 mm, 3 µm

Buffer: Sodium phosphate buffer, pH 7.0, 150 mM

Gradient: 0-100% Buffer A from 0-30 min


Flow Rate: 1.0 mL/min

Sample: CHO-humanized MAb, 5 mg/mL – intact

Injection: 5 µL

Detector: UV, 220 nm

Temperature: Ambient

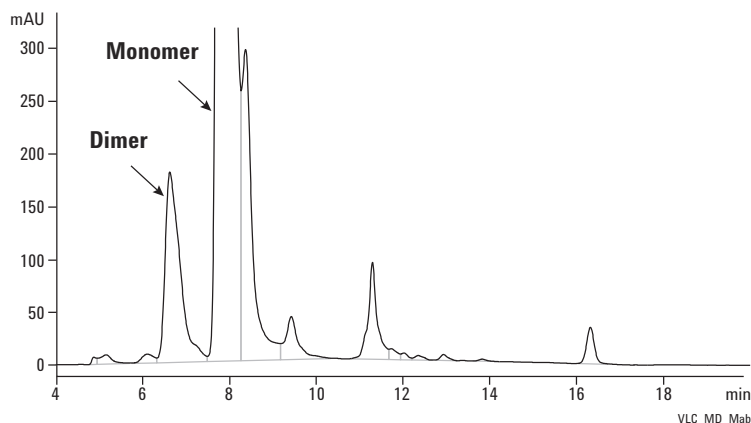**Separation of heated, stressed MAb**

Column: Bio SEC-3, 300Å
5190-2511
7.8 x 300 mm, 3 µm

Buffer: Sodium phosphate buffer, pH 7.0,
150 mM +150 mM sodium sulfate

Gradient: 0-100% Buffer A from 0-30 min


Flow Rate: 1.0 mL/min

Sample: CHO-humanized MAb, 5 mg/mL – stressed at 60 °C

Injection: 5 µL

Detector: UV, 220 nm

Temperature: Ambient

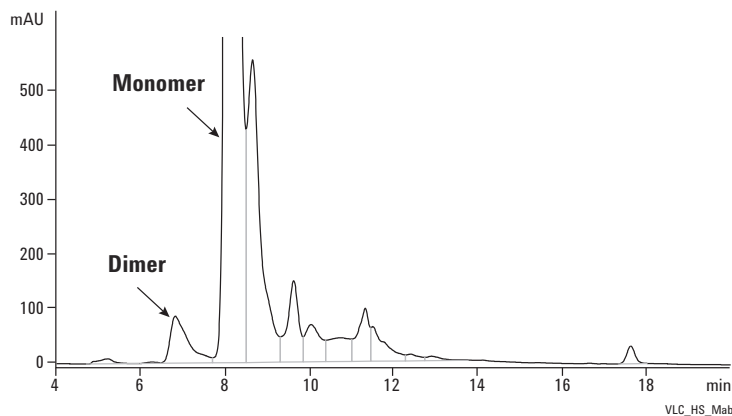

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library


Nucleosides, purines and pyrimidines

Column: Eclipse Plus Phenyl Hexyl
959993-912
4.6 x 150 mm, 5 µm

Mobile Phase: 1% MeOH: 99% 20 mM Ammonium Acetate, pH 4.5

Flow Rate: 1 mL/min

Detector: UV, 254 nm

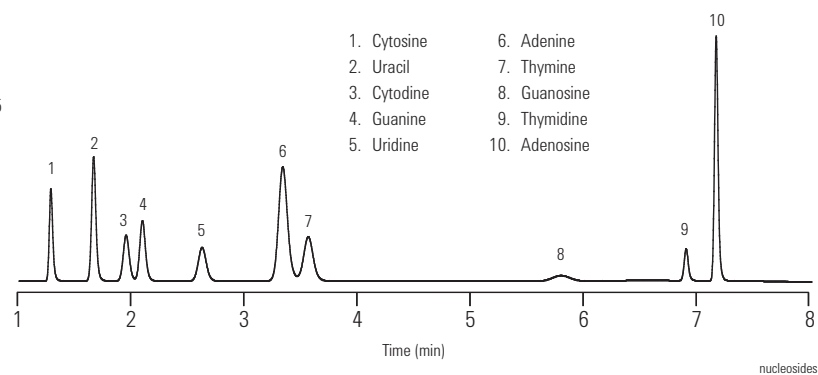**Amino acid standard separation on Eclipse Plus C18**

Column: Eclipse Plus C18
959763-902
2.1 x 150 mm, 3.5 µm

Mobile Phase: A: 10 mM Na₂HPO₄, 10 mM Na₂B₄O₇, 0.5 mM NaN₃, pH 8.2
B: acetonitrile: methanol: water (45:45:10) (v/v/v)

Flow Rate: 0.42 mL/min

Temperature: 40 °C

Detector: UV, 338 nm, then switch to 280 nm at 15.7 min


Sample: 900 pmol Amino Acids with extended Amino Acids and Internal Standards (500 pmol)

Derivatization: Automated, online, OPA / FMOC

- | | | |
|--------|---------|---------|
| 1. ASP | 9. ARG | 17. PHE |
| 2. GLU | 10. ALA | 18. ILE |
| 3. ASN | 11. TYR | 19. LEU |
| 4. SER | 12. CY2 | 20. LYS |
| 5. GLN | 13. VAL | 21. HYP |
| 6. HIS | 14. MET | 22. SAR |
| 7. GLY | 15. NVA | 23. PRO |
| 8. THR | 16. TRP | |

Gradient

Time (min)	% B
0	2
0.5	2
20	57
20.1	100
23.5	100
23.6	2
25	stop

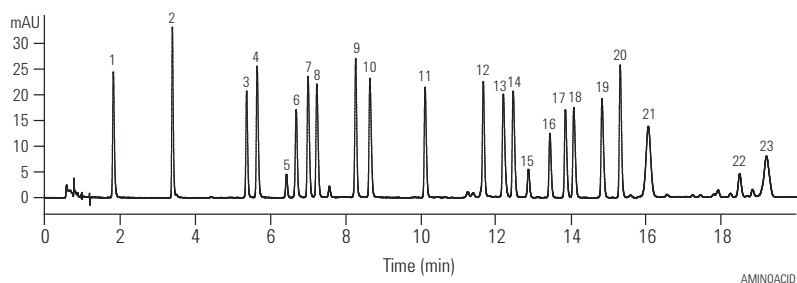

Antibodies: Fast separation of IgM and IgG antibodies

Column: ZORBAX GF-250
884973-701
4.6 x 250 mm, 4 µm


Mobile Phase: 200 mM Sodium Phosphate (pH 7), 0.01% Azide

Flow Rate: 0.94 mL/min

Temperature: Ambient

Detector: UV, 230 nm

Sample: 2.5 µL (1 mg/mL)

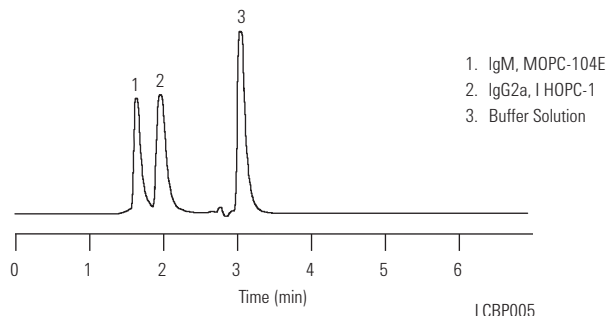**Glycosylated proteins:****Large molecules on Poroshell 300SB-C18 and 300SB-C8**

Column A: Poroshell 300SB-C18
661750-902
1.0 x 75 mm, 5 µm

Column B: Poroshell 300SB-C8
661750-906
1.0 x 75 mm, 5 µm

Column C: ZORBAX 300SB-C18
865630-902
1.0 x 50 mm, 3.5 µm

Mobile Phase: A: 0.1% TFA in H₂O
B: 0.07% TFA in ACN


Flow Rate: A, B: 0.454 mL/min
C: 0.071 mL/min

Gradient: A, B: 0 min 5% B
10 min 100% B
C: 0 min 5% B
50 min 100% B

Temperature: 70 °C

Detector: DAD 212 nm, 1.7 µL flow cell, <0.01 min peak width

Sample: Large glycosylated proteins

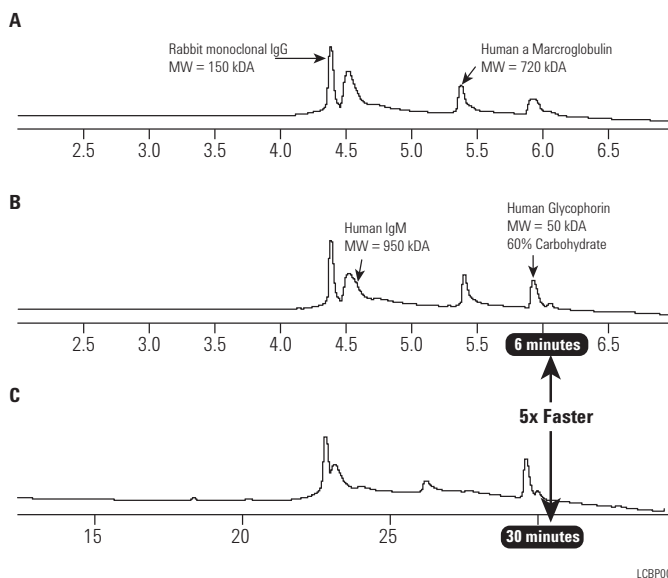

Courtesy of:
Novartis AG, Basel.
Dr. Kurt Forrer
Patrik Roethlisberger



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

HSA tryptic digest on ZORBAX Rapid Resolution HT 1.8 μ m

Column A: ZORBAX SB-C18
883700-922
2.1 x 150 mm, 5 μ m

Column B: ZORBAX SB-C18
822700-902
2.1 x 50 mm, 1.8 μ m

Mobile Phase: A: Water w/0.1% TFA
B: ACN w/0.1% TFA


Flow Rate: A: 0.2 mL/min
B: 0.5 mL/min

Gradient: A: 2 to 50% B in 70min
B: 2 to 50% B in 30min

Temperature: 50 °C

Detector: UV, 214 nm

Sample: HSA tryptic digest, 8 μ L of 15 pmol/ μ L
(120 pmol on column)

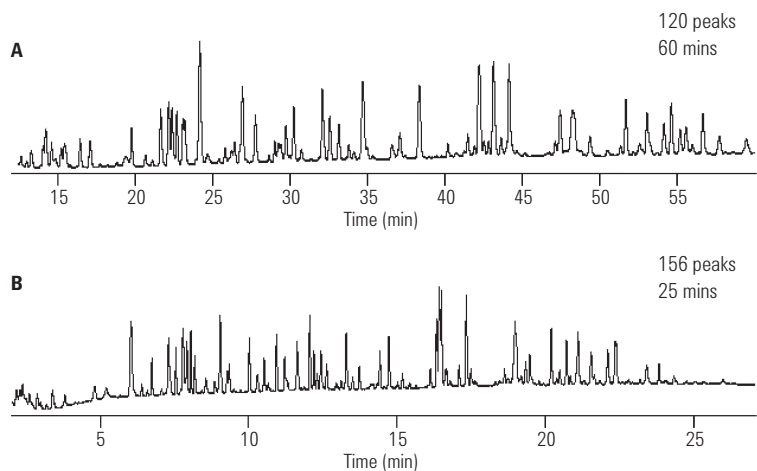

LCBP013

Human serum: Low abundance protein isolation and identification from 1-D gel band by LC/MS

Column: ZORBAX 300SB-C18
Trap: 0.3 x 5 mm, 5 μ m, 5065-9913
Analytical: 0.3 x 150 mm,
5 μ m, 5064-8263


Mobile Phase: A: Water + 0.1% Formic acid
B: Acetonitrile + 0.1% Formic acid

Flow Rate: 6 μ L/min

Gradient: 0 min 3% B
5 min 3% B (loading)
50 min 45% B
52 min 80% B
57 min 80% B
60 min 3% B

Sample: Band from 1-D in gel digest

Base Peak Chromatogram

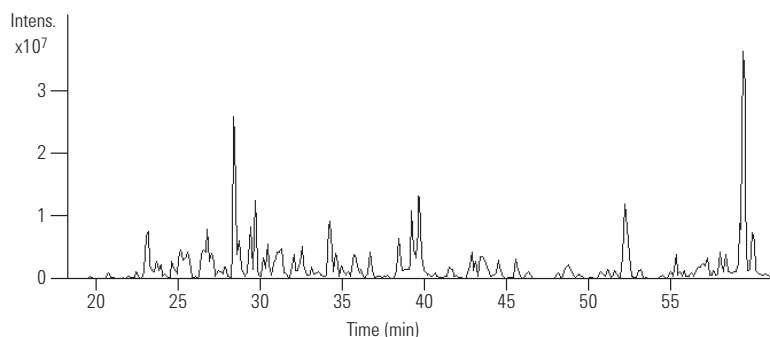

LCBP014

Proteins Identified

1. α -1-Antichymotrypsin
2. Antithrombin-III Precursor
3. Complement Factor B Precursor

Sample Preparation of Human Serum:

Major serum proteins removed using Multiple Affinity Removal

Column: 4.6 x 100 mm, P/N 5185-5985

Followed by 1-D gel digest

Monoclonal IgG1 chains: Separation on Poroshell 300SB-C8

Column: Poroshell 300SB-C8
660750-906
2.1 x 75 mm, 5 µm

Mobile Phase: A: 90% water: 10% ACN + 3 mL/L of MW 300 PEG
B: 10% water: 90% ACN + 3 mL/L of MW 300 PEG


Flow Rate: 1.0 mL/min

Gradient: 0 min 25% B
10 min 40% B
10.1 min 25% B
12 min 25% B

Temperature: 70 °C

Sample: Monoclonal IgG1

*Courtesy of:
Novartis AG, Basel.
Dr. Kurt Forrer
Patrik Roethlisberger*

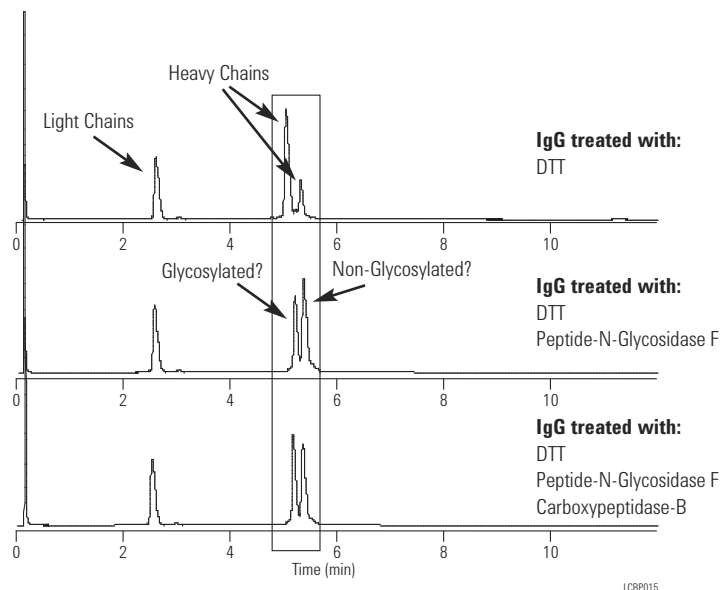

LCBP015

Use ZORBAX Extend-C18 for alternate selectivity at high pH

Column: ZORBAX Extend-C18
773700-902
2.1 x 150 mm, 5 µm

Mobile Phase: A: 0.1% TFA in Water
B: 0.085% TFA in 80% ACN


A: 20 mM NH₄OH in Water
B: 20 mM NH₄OH in 80% ACN

Flow Rate: 0.25 mL/min

Gradient: 5-60% B in 20 min

Temperature: 25 °C

MS Conditions: Pos. Ion ESI-VI 70V, Vcap 4.5 kV
N₂ – 35 psi, 12 L/min, 300 °C
4 µL (50 ng each peptide)

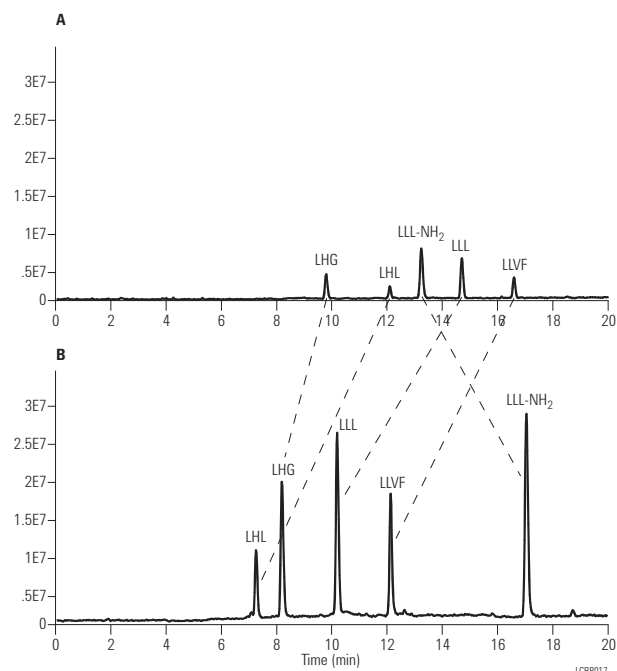

LCBP017

The Extend column can be used for high pH separations of peptides. At high and low pH, very different selectivity can result. Just by changing pH, a complementary method can be developed and it is possible to determine if all peaks are resolved. The Extend column can be used at high and low pH, so the complementary separation can be investigated with one column. Better MS sensitivity for this sample is also achieved at high pH.

Nucleosides: Separation of deoxy and ribonucleosides

Column: ZORBAX SB-C3

883975-909

4.6 x 150 mm, 5 µm


Mobile Phase: 4 mM Ammonium Phosphate (pH 4.0 with Phosphoric Acid)

Flow Rate: 2.0 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

Sample: 2 µL (1.6 µg each)

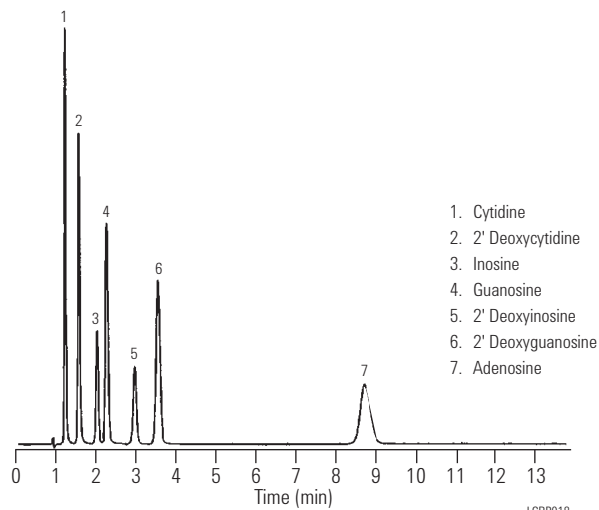

Nucleotides: Separation of mononucleotides

Column: ZORBAX SAX

880952-703

4.6 x 250 mm, 5 µm


Mobile Phase: 0.1 M $\text{NH}_4\text{H}_2\text{PO}_4$

Flow Rate: 2.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Orotic Acid, UMP, GMP, XMP

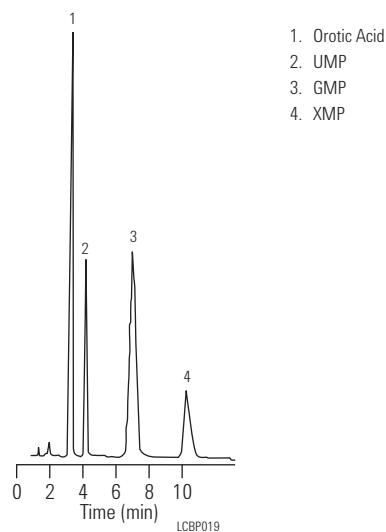

Separation of basic peptides on Bonus-RP versus traditional Alkyl phase

Column A: ZORBAX Bonus-RP
883668-901
4.6 x 150 mm, 5 µm

Column B: Alkyl C8

Mobile Phase: A: 0.010 M ammonium phosphate, pH 7/0.050 M sodium perchlorate
B: 0.010 M ammonium phosphate/0.050 M sodium perchlorate in 50% ACN

Flow Rate: 1.0 mL/min


Gradient: 0-100% B in 50 min

Temperature: 40 °C

Detector: 215 nm

Sample: Basic 11-residue peptides with net +1, +2, +3, +4 positive charges at neutral pH

C1: Ac-Gly-Gly-Gly-Leu-Gly-Gly-Ala-Gly-Gly-Leu-Lys-amide
C2: Ac-Lys-Tyr-Gly-Leu-Gly-Gly-Ala-Gly-Gly-Leu-Lys-amide
C3: Ac-Gly-Gly-Ala-Leu-Lys-Ala-Leu-Lys-Gly-Leu-Lys-amide
C4: Ac-Lys-Tyr-Ala-Leu-Lys-Ala-Leu-Lys-Gly-Leu-Lys-amide

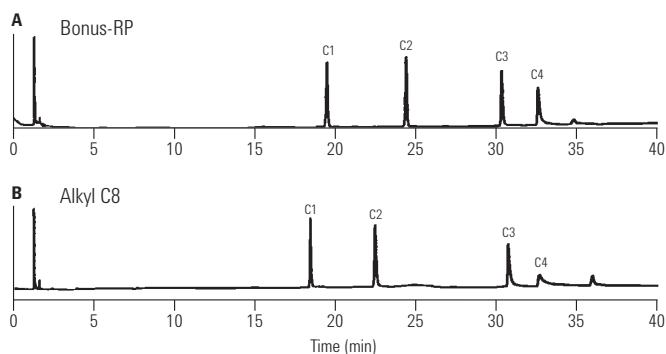

LCBP020

Peptides: Effect of TFA concentration

Column: ZORBAX 300SB-C8
883995-906
4.6 x 150 mm, 5 µm

Mobile Phase: A: Water and TFA
B: ACN and TFA


Flow Rate: 1.0 mL/min

Gradient: 0 min 0% B
30 min 30% B

Temperature: 40 °C

Detector: UV, 254 nm

Sample: Peptide Standards S1-S5, decapeptides differing slightly in hydrophobicity, 6 µL

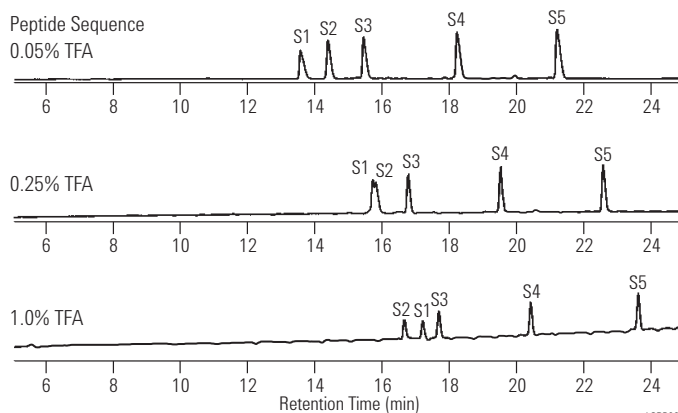

LCBP021

Exploiting chemical stability – TFA concentration


Column: PLRP-S 100Å
PL1512-5500
4.6 x 250 mm, 5 µm

Mobile Phase: A: TFA (various %) in water
B: TFA (various %) in ACN

Gradient: Linear 12-40% B in 15 min

Flow Rate: 1.0 mL/min

Detector: ELS (neb=75 °C, evap=85 °C, gas=1.0 SLM)

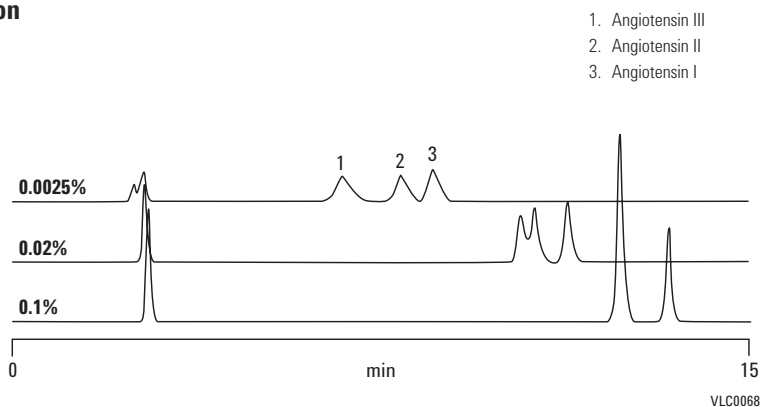

Peptides: Separation of Antiotensins I, II, III with TFA and NH₄OH

Column: ZORBAX Extend-C18
773700-902
2.1 x 150 mm, 5 µm

Mobile Phase: A: Acidic Conditions
A: 0.1% TFA in water
B: 0.085% TFA in 80% ACN
B: Basic Conditions
A: 10 mM NH₄OH in water
B: 10 mM NH₄OH in 80% ACN



Flow Rate: 0.2 mL/min

Gradient: 15-50% B in 15 min

Temperature: 35 °C

MS Conditions: Pos. Ion ESI - Vf 70V, Vcap 4.5 kV
N₂-35 psi, 12 L/min, 325 °C

Sample: 2.5 µL sample (50 pmol each)

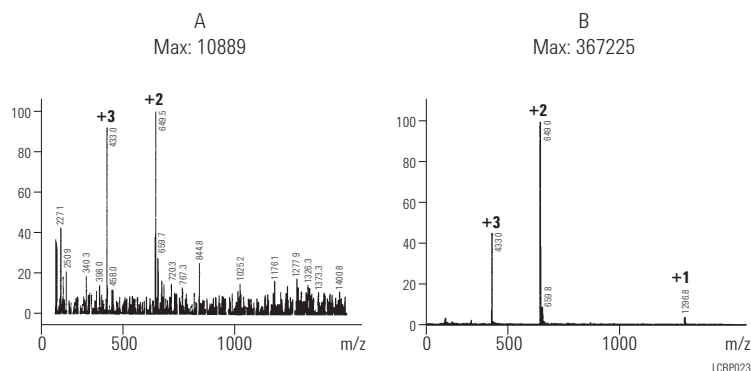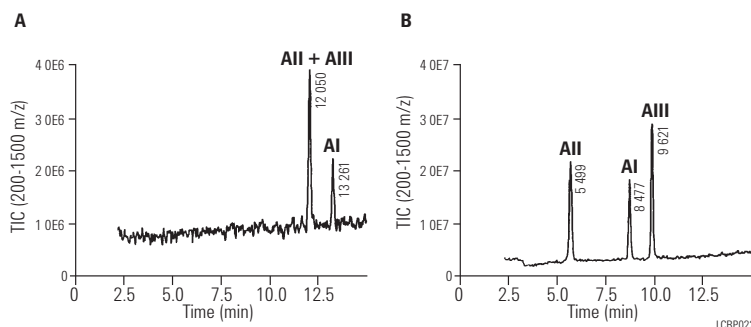

Peptides/proteins:
Equivalent gradient separations

Column: ZORBAX 300SB-C8
 883995-906
 4.6 x 150 mm, 5 µm

Column: ZORBAX 300SB-C8
 883750-906
 2.1 x 150 mm, 5 µm

Mobile Phase: A: 95% Water: 5% ACN with 0.1% TFA
 B: 5% Water: 95% ACN with 0.085% TFA


Flow Rate: A: Analytical
 1 mL/min
 B: Narrow Bore
 0.2 mL/min

Gradient: 10-60% B in 30 min

Temperature: 35 °C

Detector: UV, 215 nm

Sample: 10 µL injection, concentration 2-6 µg

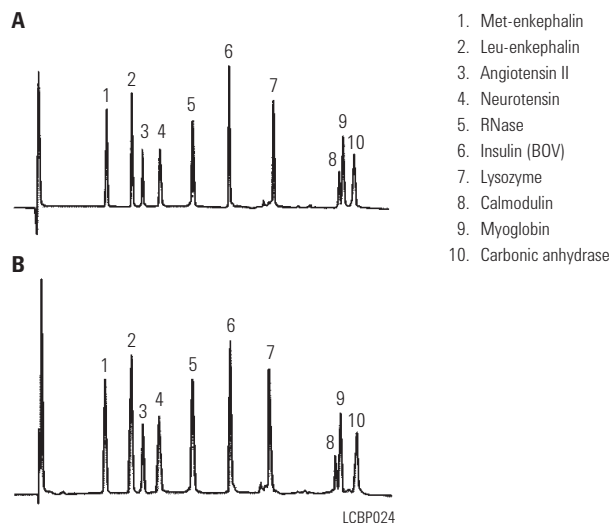

Peptides/proteins:
Effect of elevated temperature

Column: ZORBAX 300SB-C3
 883995-909
 4.6 x 150 mm, 5 µm

Mobile Phase: A: 5:95 ACN:Water with 0.10% TFA (v/v%)
 B: 95:5 ACN:Water with 0.085% TFA (v/v%)


Flow Rate: 1.0 mL/min

Gradient: 15-53% in 20 min, post time 12 min

Temperature: Ambient – 60 °C

Detector: UV, 215 nm

Sample: Polypeptides

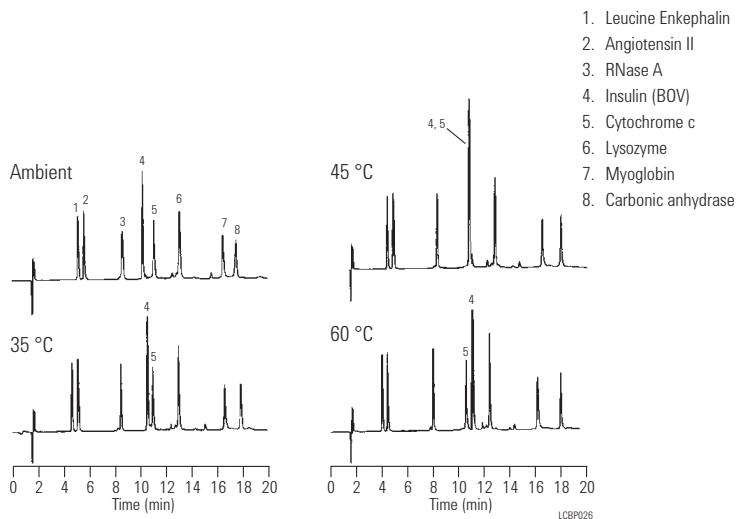

Separation of polypeptides in under 1 minute

Column: Poroshell 300SB-C18
660750-902
2.1 x 75 mm, 5 µm

Mobile Phase: A: 0.1% TFA, H₂O
B: 0.07% TFA, ACN

Flow Rate: 3 mL/min

Gradient: 0-100% B in 1.33 min

Temperature: 70 °C


Detector: DAD 215/16 nm, ref = 310/10 nm

Sample: Peptides/proteins, 0.5 µL

Mixer bypassed with P/N G1312-67301; Loop-bypass program

Sample (peptides/proteins)

1. gly-tyr	0.125 mg/mL
2. Val-tyr-val	0.5 mg/mL
3. Met-enkephalin	0.5 mg/mL
4. Leu-enkephalin	0.5 mg/mL
5. Angiotensin II	0.5 mg/mL
6. RNase A	1 mg/mL
7. Cytochrome C	1 mg/mL
8. Holotransferrin	1 mg/mL
9. Apomyoglobin	1 mg/mL

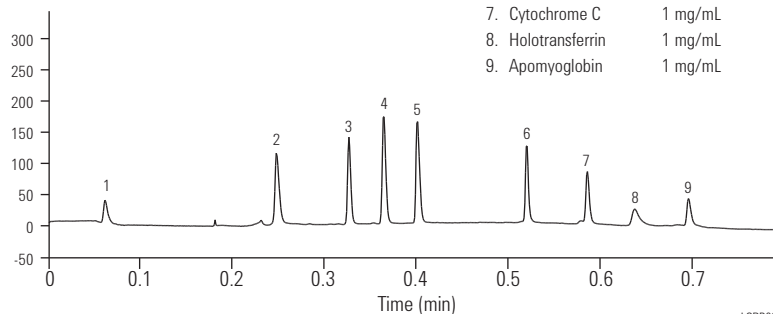

LCBP030

Fast, high-resolution separation of peptides and proteins with Poroshell 300SB-C18

Column: Poroshell 300SB-C18
660750-902
2.1 x 75 mm, 5 µm

Mobile Phase: A: 0.1% TFA
B: 0.07% TFA in ACN

Flow Rate: 3.0 mL/min (360 bar pressure)


Gradient: 5-100% B in 1.0 min

Temperature: 70 °C

Detector: UV, 215 nm

Spaces between solutes indicate good peak capacity for rapidly separating complex samples.

1. Angiotensin II
2. Neurotensin
3. RNase
4. Insulin
5. Lysozyme
6. Myoglobin
7. Carbonic Anhydrase
8. Ovalbumin

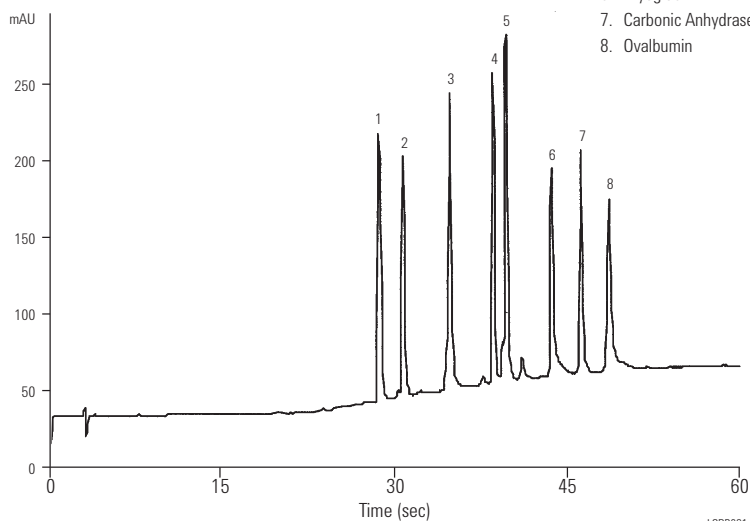

LCBP031

**Peptide RP-HPLC/ESI-MS
using NH_4OH mobile phase
yields both positive and negative ion spectra**

Column: ZORBAX Extend-C18
773700-902
2.1 x 150 mm, 5 μm


Flow Rate: 0.25 mL/min

Gradient: 5-60% B in 20 min

Temperature: 25 °C

MS Conditions: Pos. Ion ESI – Vt 70 V, Vcap 4.5 kV,
 N_2 – 35 psi, 12 L/min, 300 °C
TIC 150-1500 m/z

Sample: 4 μL (50 ng each peptide)

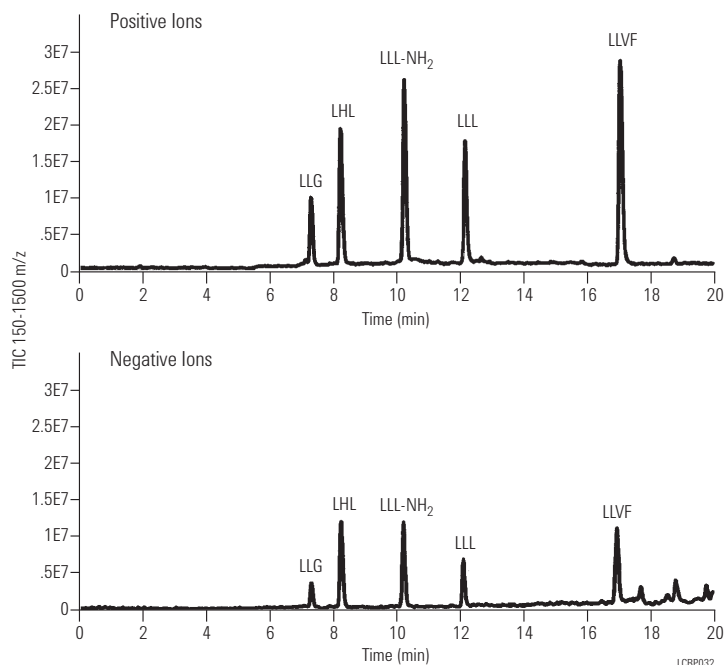

**Comparison of $\text{A}\beta$ peptide RP-HPLC
separations at low and high pH**

Column: ZORBAX Extend-C18
773700-902
2.1 x 150 mm, 5 μm

Mobile Phase: A: 0.1% TFA in water
B: 0.085% TFA in 80% ACN

Flow Rate: 0.25 mL/min

Gradient: 29-41% B in 30 min

Temperature: 80 °C

Detector: UV, 210 nm

Sample: 5 μL sample (100 pmol each)

Mobile Phase: A: 20 mM NH_4OH in water
B: 20 mM NH_4OH in 80% ACN


Flow Rate: 0.25 mL/min

Gradient: 26-38% B in 30 min

Temperature: 25 °C

Detector: UV, 210 nm

Sample: 5 μL sample (100 pmol each)

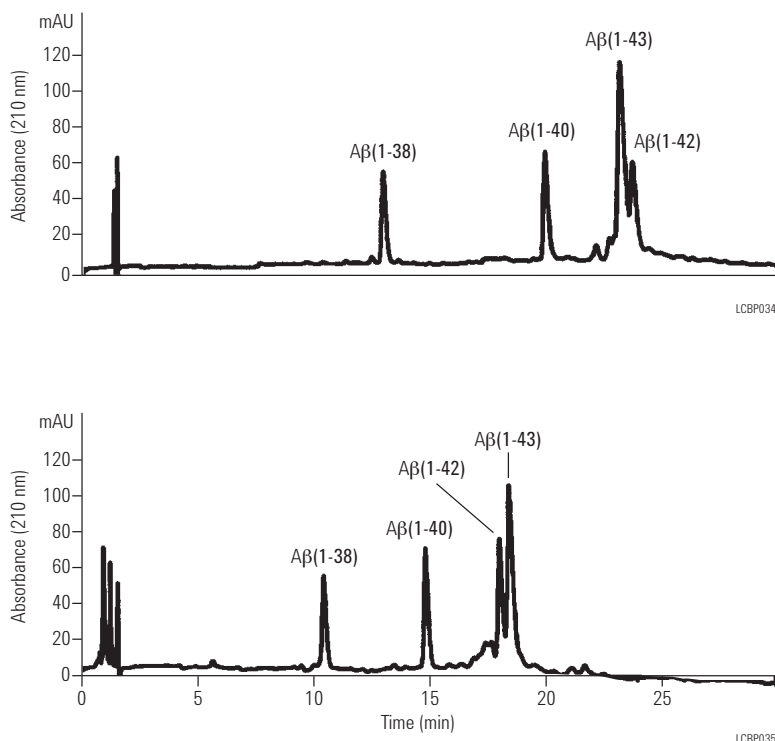

Selectivity comparison of TFA and NH_4OH for peptide RP-HPLC\ESI-MS analysis

Column: ZORBAX Extend-C18
773700-902
2.1 x 150 mm, 5 μm

Mobile Phase: TFA Conditions:
A: 0.1% TFA in water
B: 0.085% TFA in 80% ACN
 NH_4OH Conditions:
A: 20 mM NH_4OH in water
B: 20 mM NH_4OH in 80% ACN


Flow Rate: 0.25 mL/min

Gradient: 5-60% B in 20 min

Temperature: 25 °C

MS Conditions: Pos. Ion ESI – V_f 70V, V_{cap} 4.5 kV,
 N_2 – 35 psi, 12 L/min., 300 °C
TIC 150-1500 m/z

Sample: 4 μL (50 ng each peptide)

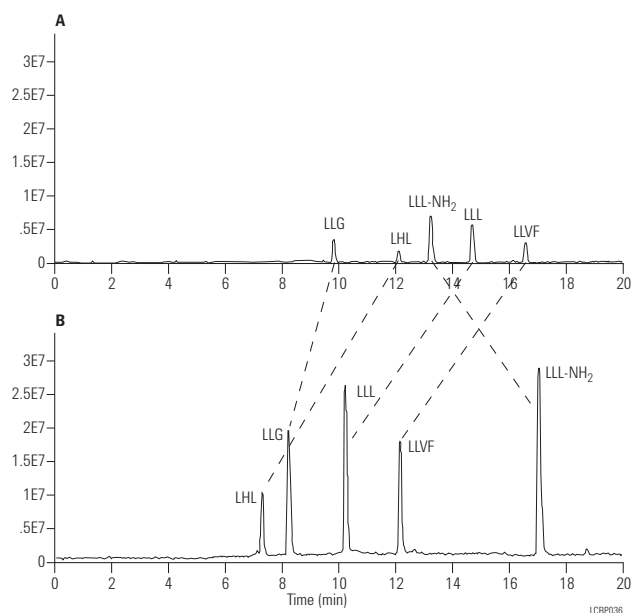

Peptide phosphorylation sites LC and LC/MS using Capillary LC columns

Column: ZORBAX 300SB-C18
5064-8268
0.5 x 150 mm, 3.5 μm

Mobile Phase: A: Water + 0.1% Formic acid
B: Acetonitrile + 0.1% Formic acid

Flow Rate: 5.5 $\mu\text{L}/\text{min}$


Gradient: 5-55% B in 50 min, to
85% B from 55-57 min

Detector: UV, 206 nm


MS Conditions: LC/MS: Pos. Ion ESI with LC/MSD trap
V_{cap}: 4000 V
Drying gas flow: 7 L/min
Drying gas temperature: 250 °C
Nebulizer: 15 psi
Capillary Exit Volt: 50 V Max
Accum Time: 300 ms
Total Averages: 3
Isolation Width: 3 m/z
Frag Amplitude: 1.0 V

Sample: Beta case in digest, 100 nL (4 pmol)


MS

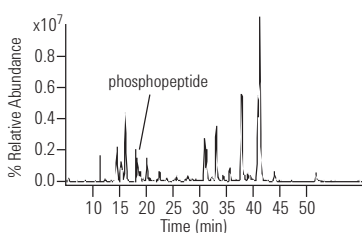


Full Scan MS

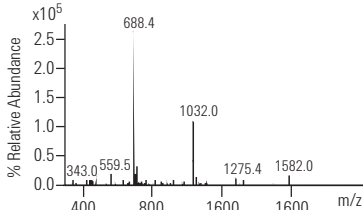

UV

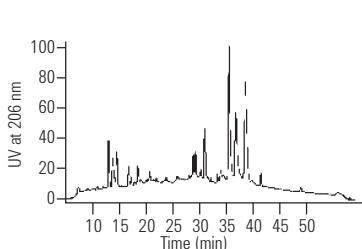

MS/MS of $[\text{M}+2\text{H}]^{2+}$ at m/z 1032

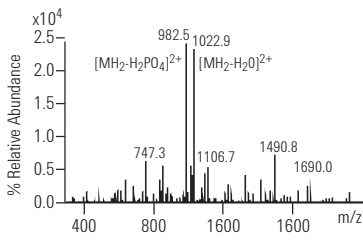

Proteins: Effect of bonded phase, RP

Column A: ZORBAX 300SB-C8
883995-906
4.6 x 150 mm, 5 µm

Column B: ZORBAX 300SB-CN
883995-905
4.6 x 150 mm, 5 µm

Mobile Phase: A: 0.1% TFA in Water,
B: 0.1% TFA in 50/50 ACN/Water


Flow Rate: 1.0 mL/min

Gradient: 1. 46-96% B in 60 min 23-48% ACN
2. 50-86% B in 60 min 25-43% ACN

Temperature: 50 °C

Detector: UV, 210 nm

Sample: Wheat proteins, including w-gliadins

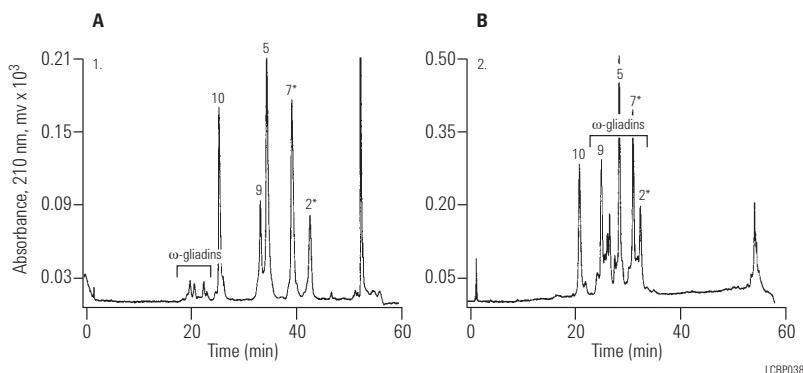**Proteins: Effect of bonded phase**

Column A: ZORBAX RRHD 300SB-C18
883995-902
4.6 x 150 mm, 5 µm

Column B: ZORBAX 300SB-C8
883995-906
4.6 x 150 mm, 5 µm

Column C: ZORBAX 300SB-C3
883995-909
4.6 x 150 mm, 5 µm

Column D: ZORBAX 300SB-CN
883995-905
4.6 x 150 mm, 5 µm

Mobile Phase: A: 0.1% TFA in H₂O
B: 0.09% TFA in 80% ACN/20% Water


Flow Rate: 1.0 mL/min

Gradient: 25-70% B in 40 min

Temperature: 60 °C

Detector: UV, 210 nm

Sample: Polypeptides, 3 µg each

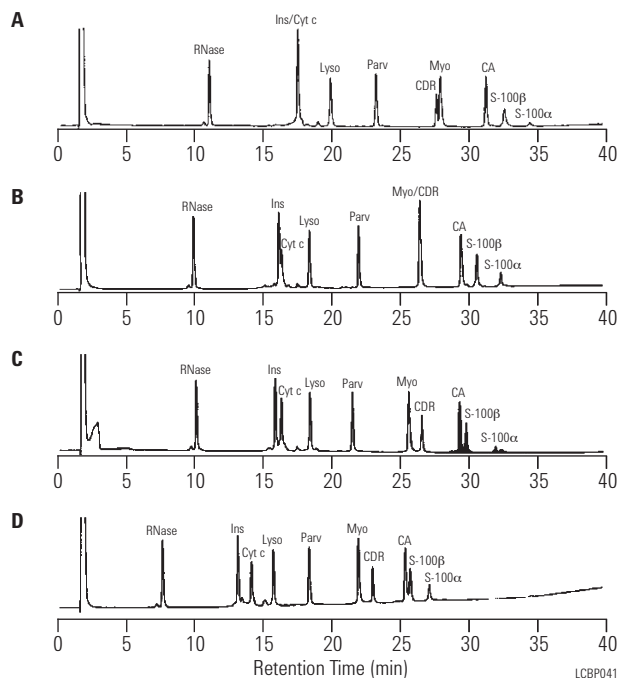

Standard proteins by reversed-phase


Column: PLRP-S 4000Å
PL1512-1803
4.6 x 50 mm, 8 µm

Mobile Phase: A: 0.1% TFA in 95% water:5% ACN
B: 0.1% TFA in 5% water:95% ACN

Gradient: Linear 18-60% B in 1 min

Flow Rate: 4.0 mL/min

Detector: UV, 280 nm

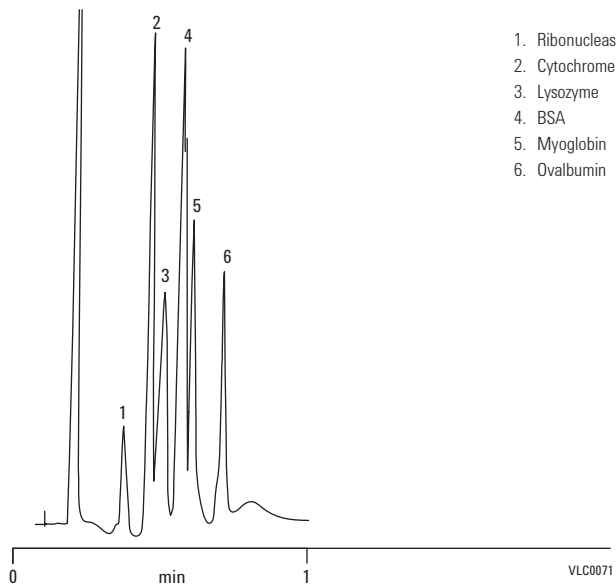

Standard ion-exchange protein separation


Column: PL-SAX 1000Å
PL1551-1502
4.6 x 50 mm, 5 µm

Mobile Phase: A: 10 mM Tris HCl pH 8
B: A+0.35 M NaCl pH 8

Gradient: 0-100% B in 20 min

Flow Rate: 1.0 mL/min

Detector: UV, 220 nm

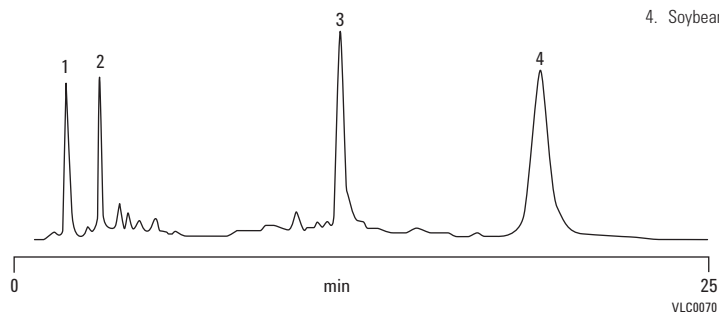

Deoxynucleosides: Using rapid resolution 3.5 μm columns

Column A: ZORBAX SB-CN
883975-905
4.6 x 150 mm, 5 μm


Column B: ZORBAX SB-CN
835975-905
4.6 x 50 mm, 3.5 μm

Mobile Phase: A: 0.1% TFA
B: 90/10 v/v Methanol/Water (0.1% TFA)
Isocratic, 97.5% A, 2.5% B

Flow Rate: 1.0 mL/min

Temperature: 30 °C

Detector: UV, 254 nm

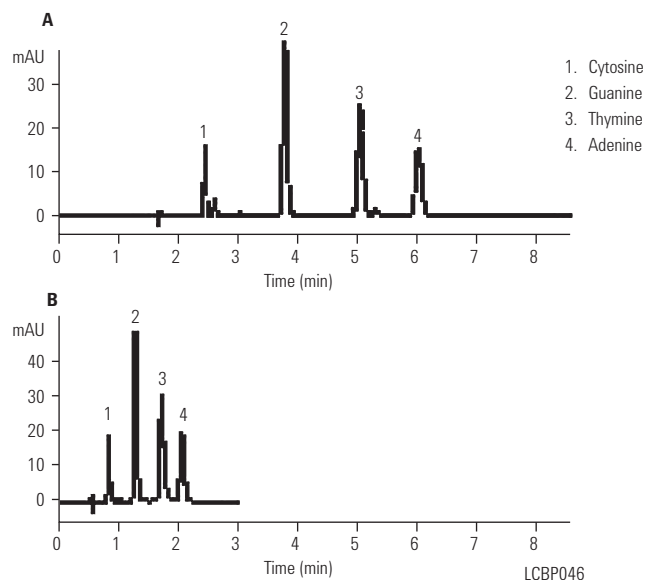

BSA tryptic digest on RRHT

Column: ZORBAX SB-C18
820700-902
2.1 x 150 mm, 1.8 μm

Mobile Phase: A: 0.1% TFA, 5% ACN
B: 0.08% TFA, 95% ACN


Flow Rate: 0.5 mL/min

Gradient: A: Time 0% B 5 min, Time 30% B 60 min
B: Time 0% B 5 min, Time 45% B 60 min
C: Time 0% B 5 min, Time 67.5% B 60 min

Temperature: 80 °C

Detector: UV, 214 nm

Sample: BSA tryptic digest

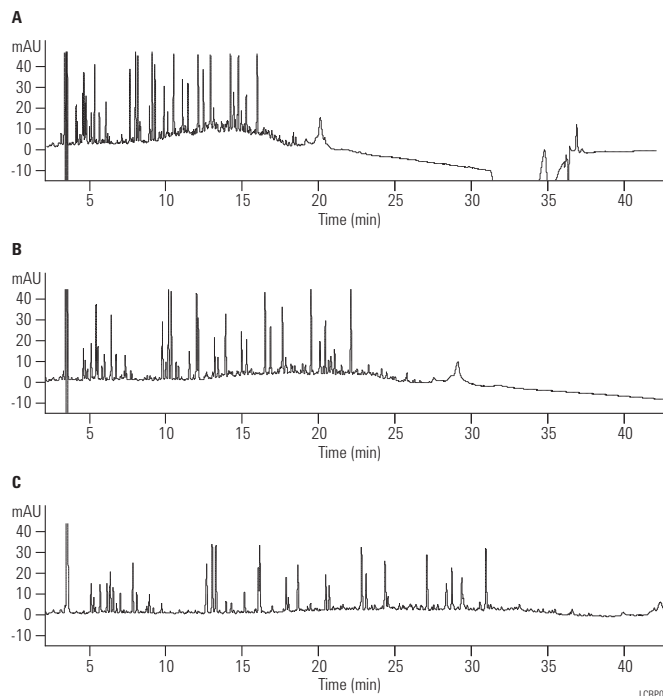


Catecholamines

Column: PLRP-S 100Å
PL1111-3500
4.6 x 150 mm, 5 µm

Mobile Phase: 95% 25 mM citric acid,
25 mM Na₂HPO₄, 1 mM heptane
sulfonic acid:5% ACN, pH 2.85

Flow Rate: 1.0 mL/min

Detector: UV, 280 nm

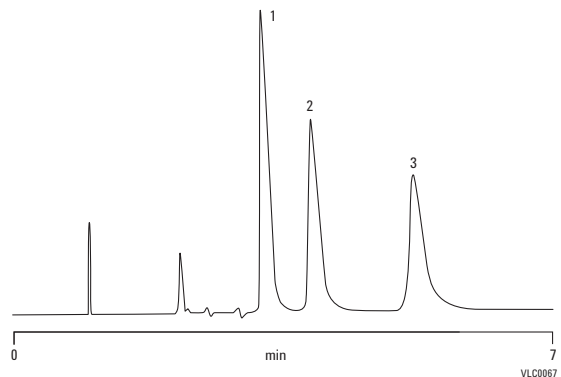

1. Noradrenaline
2. Adrenaline
3. Dopamine

Whey proteins in dairy samples – milk

Column: PLRP-S 300Å
PL1512-3801
4.6 x 150 mm, 8 µm


Mobile Phase: A: 0.1% TFA in 99% water:1% ACN
B: 0.1% TFA in 1% water:99% ACN

Gradient: 36-48% B, 0-24 min, 48-100% B, 24-30 min
100% B, 30-35 min, 100-36% B, 35-40 min

Flow Rate: 1.0 mL/min

Injection Volume: 10 µL

Detector: UV, 220 nm

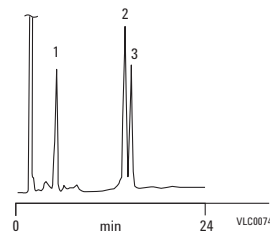

1. α-Lactalbumin
2. β-Lactoglobulin (B chain)
3. β-Lactoglobulin (A chain)

Temperature as a tool to enhance mass transfer and improve resolution of oligonucleotides in ion-pair reversed-phase HPLC

Column: PLRP-S 100Å
PL1512-1300
4.6 x 50 mm, 3 µm


Mobile Phase: A: 100 mM TEAA
B: 100 mM TEAA in 25% ACN

Gradient: 5% change in buffer B over 5 min

Flow Rate: 1.0 mL/min

Temperature: 35 °C, 60 °C, or 80 °C

Detector: UV, 254 nm

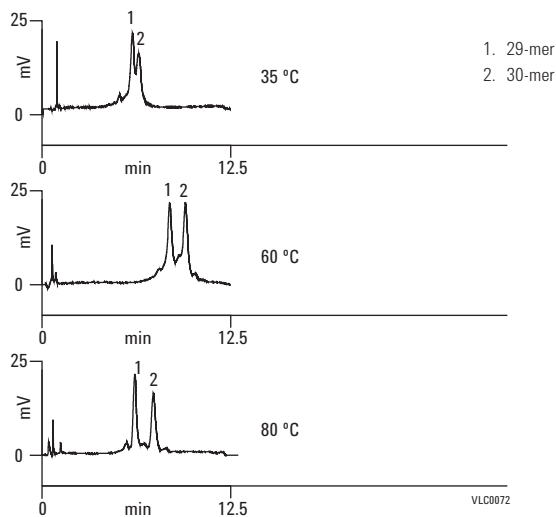

Hydrophilic purine/pyrimidine separation


Column: ZORBAX SB-Aq
883975-914
4.6 x 150 mm, 5 µm

Mobile Phase: 50 mM NaOAc, pH 4.6

Flow Rate: 2.0 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

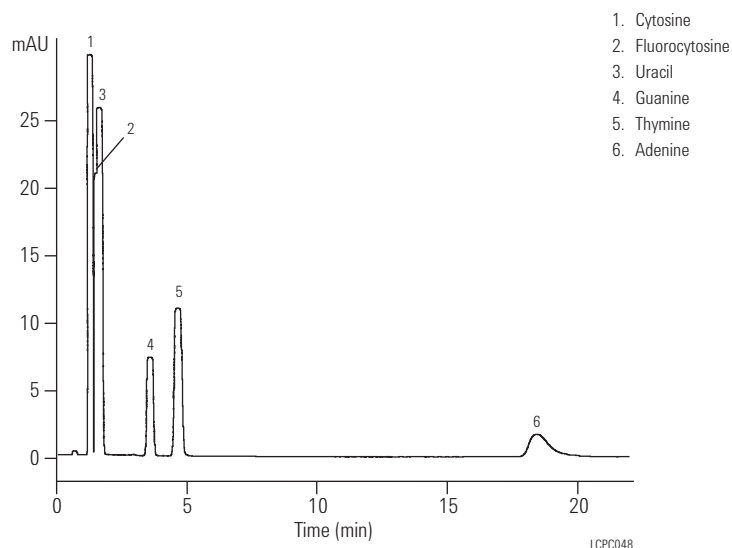

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Chemical/Industrial Applications

Analysis of biocides in hand sanitizer

Column: ZORBAX RRHD Eclipse Plus C18
959757-902
2.1 x 50 mm, 1.8 μ m

Mobile Phase: A: H₂O (0.5% TFA)
B: ACN (0.04% TFA)


Flow Rate: 1.7 mL/min

Gradient: Time 0.0 95/5 A/B DAD: 275 nm (0 min)
Time 1.0 55/45 A/B 225 nm (0.46 min)
Time 1.1 0/100 A/B 255 nm (0.67 min)

Sample: 1 μ L injection of 50 ppm std.


Temperature: 30 $^{\circ}$ C

Biocides Standards

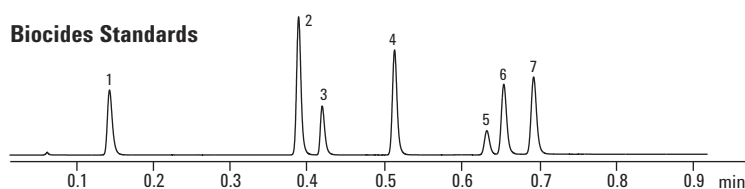

1. Kathon 1A
2. Kathon 1B
3. Carbendazim
4. 1,2-Benzisothiazol-3(2H)-one
5. 2-Phenoxyethanol
6. Benzoic Acid
7. Methyl Paraben

Hand Sanitizer Sample - Injected Neat

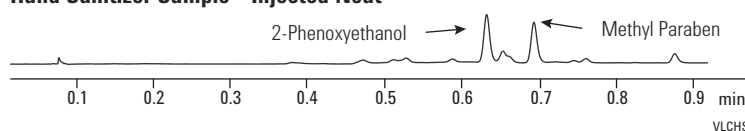

Triton X-114: Decreasing run-time by changing bonded phase

Column A: ZORBAX SB-C3
883975-909
4.6 x 150 mm, 5 μ m

Column B: ZORBAX SB-C18
883975-902
4.6 x 150 mm, 5 μ m


Mobile Phase: MeOH and H₂O (as indicated)

Flow Rate: 1.0 mL/min

Temperature: 50 $^{\circ}$ C

Detector: UV, 225 nm

Sample: Triton X-114

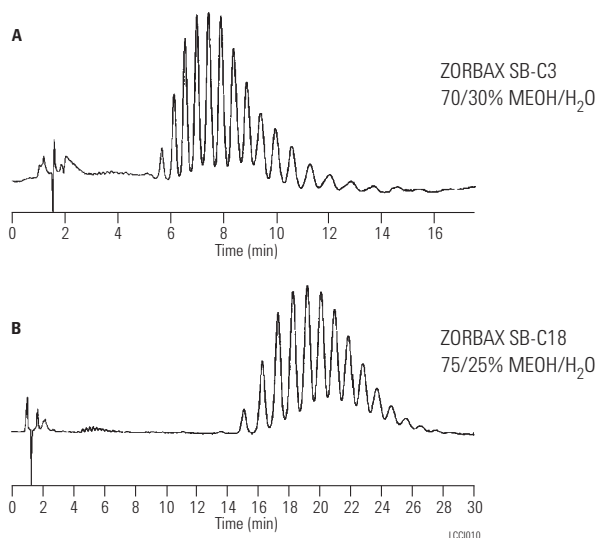

Organic acids separated on ZORBAX SB-Aq


Column: ZORBAX SB-Aq
883975-914
4.6 x 150 mm, 5 µm

Mobile Phase: 99% 20 mM NaH₂PO₄, pH 2, 1% ACN

Flow Rate: 1.0 mL/min

Temperature: 35 °C

Detector: UV, 210 nm

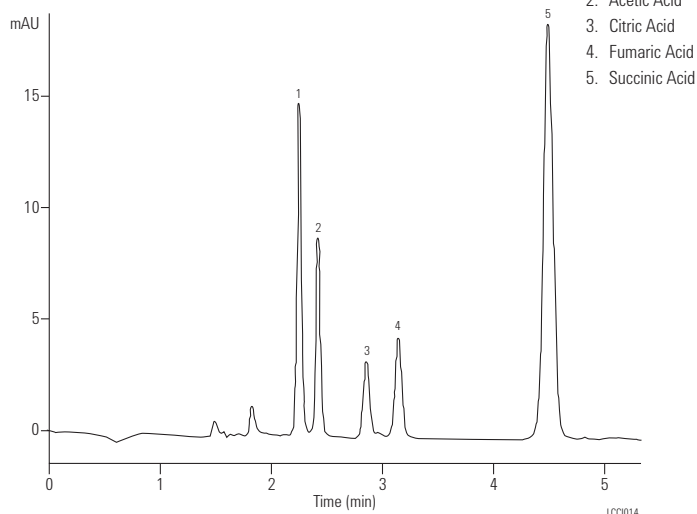**Brij 35**

Column: PLRP-S 100Å
PL1111-3500
4.6 x 150 mm, 5 µm

Mobile Phase: A: Water
B: ACN


Gradient: 0-100% B in 40 min

Flow Rate: 0.8 mL/min

Injection Volume: 10 µL

Sample Conc: 1 mg/mL

Detector: ELS (neb=50 °C, evap=70 °C, gas=1.5 SLM)

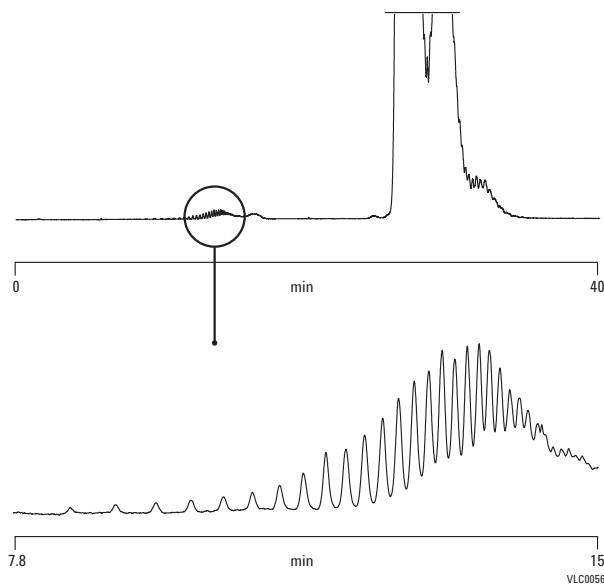

Alcohols and aliphatic compounds

Column: Hi-Plex H
PL1170-6830
7.7 x 300 mm, 8 µm


Mobile Phase: Water

Flow Rate: 0.6 mL/min


Temperature: 40 °C

Detector: 356-LC RI

1. Methanol
2. Ethanol
3. Isopropanol
4. Tert-butanol
5. n-Propanol

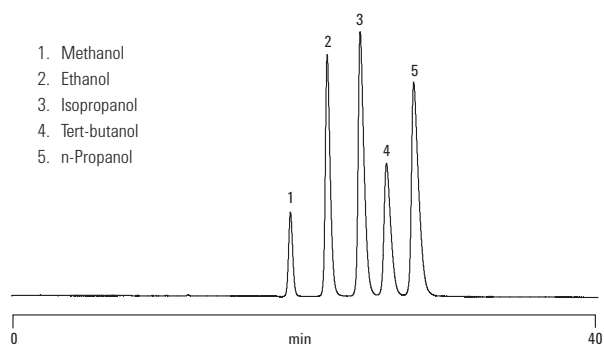

1. Ethanol
2. Acetone
3. Tert-butanol

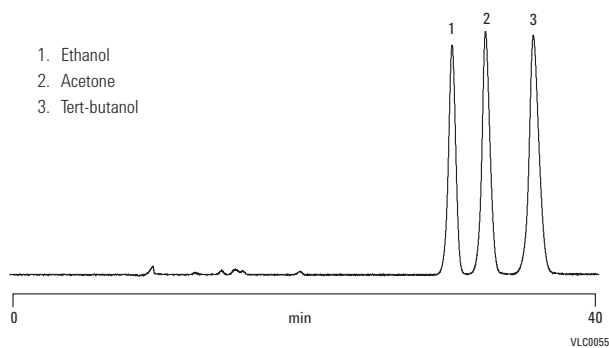

VLC0055



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Environmental Applications

NEW!

Fast LC/MS/MS analysis of group 4 pharmaceuticals from EPA-1694

Column: ZORBAX RRHD HILIC Plus
959758-901
2.1 x 100 mm, 1.8 μ m


Mobile Phase: A: 10 mM ammonium acetate in water, pH 6.7
B: acetonitrile

Flow Rate: 1 mL/min

Detector: Agilent 1290 Infinity LC with an
Agilent 6410 Triple Quadrupole Mass Spectrometer

MS Conditions: TCC: 25 $^{\circ}$ C
dMRM, ESI positive mode, cycle time 35 ms
Drying Gas: 9 L/min, 300 $^{\circ}$ C
Nebulizer Pressure: 40 psig
Capillary Voltage: 4000

Sample: 0.1 μ L injection of 0.1 mg/mL each in
acetonitrile/water (3:1): cimetidine, albuterol,
ranitidine and metformin

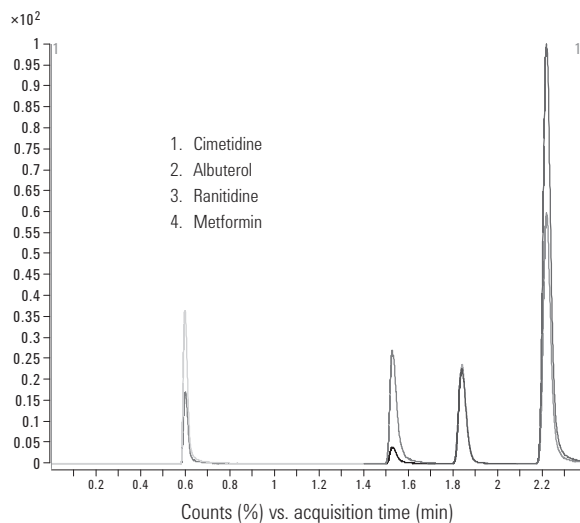
NEW!

Separation of azo dye degradation products

Column A: Poroshell 120 EC-C18
695775-902
2.1 x 100 mm, 2.7 μ m

Column B: Poroshell 120 SB-C18
685775-902
2.1 x 100 mm, 2.7 μ m

Column C: Poroshell 120 Phenyl-Hexyl
695775-912
2.1 x 100 mm, 2.7 μ m


Column D: Poroshell 120 Bonus RP
685775-901
2.1 x 100 mm, 2.7 μ m

Flow Rate: 0.4 mL/min

Gradient: 15 to 100% MeOH over 10 min

Solvent: 10 mM Ammonium acetate, pH 4.8

1. Aniline
2. o-Toluidine
3. Methoxyaniline
4. Chloroaniline
5. Benzidine
6. Dimethylbenzidine
7. 3,3'-Dimethoxybenzidine
8. Naphthylamine
9. Dichlorobenzidine

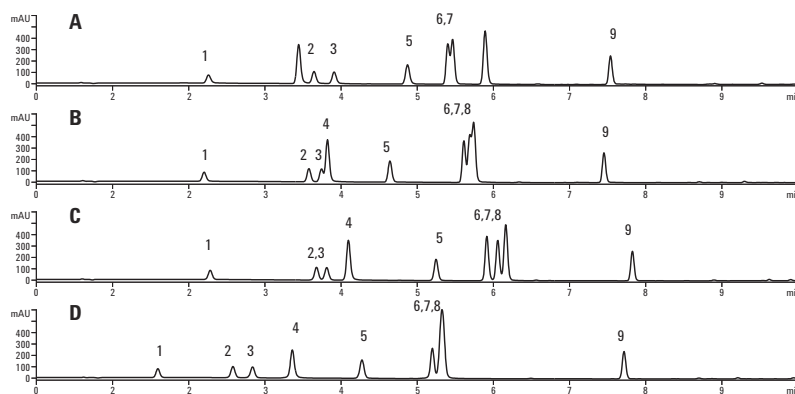

Comparison of phenols separation with Poroshell 120

Column: Poroshell 120 EC-C18
699975-902
4.6 x 50 mm, 2.7 µm

Mobile Phase: A: Water with 0.1% Formic Acid
B: Acetonitrile

Gradient: Time %B
0.8 5%
6.8 60%


1200 SL controlled temperature
at 25 °C 2 mm flow cell

Column: Poroshell 120 EC-C18
695975-902
4.6 x 100 mm, 2.7 µm

Mobile Phase: A: Water with 0.1% Formic Acid
B: Acetonitrile

Gradient: Time %B
2.0 5%
17 60%

1200 RRLLC SL controlled temperature
at 25 °C 2 mm flow cell

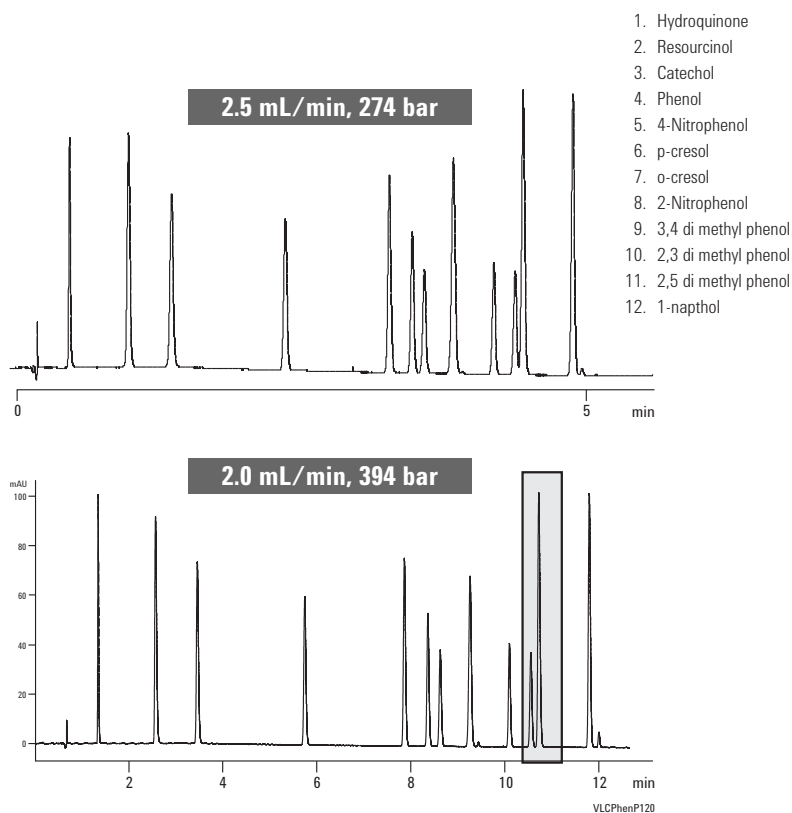

DNPH: Derivatized Aldehydes obtained from air

Column: ZORBAX ODS
884950-543
4.6 x 250 mm, 5 µm

Mobile Phase: A: 100% Water
B: 100% ACN


Flow Rate: 1.0 mL/min

Gradient: 60-75% B in 30 min; Wash: From 75-100% B in 5 min, after 5 min return to 60% B

Temperature: 35 °C

Detector: UV, 230 nm

Sample: DNPH Derivatized Aldehydes

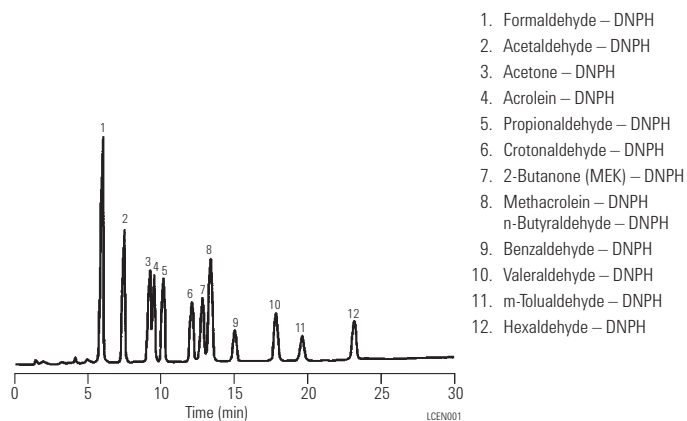

Amitrol in water by LC/MS, 0.05 ppb

Column: ZORBAX SB-C18
863954-302

3.0 x 150 mm, 3.5 μ m

Mobile Phase: A: 10 mM ammonium acetate
B: MeOH


Flow Rate: 0.4 mL/min

Gradient: 0 min, 65% B; 10 min, 65% B;
15 min, 100% B; 20 min, 65% B

Temperature: 30 °C

MS Conditions: Ionization Mode: APCI, positive polarity
SIM parameters: Ion: 213 Amitrol
Ion: 259 IS
Fragmentor: 100 V
SIM Resolution: Low
Vaporizer: 325 °C
Drying Gas (N_2): 5.0 L/min
Gas Temperature: 350 °C
Nebulizer pressure: 60 psig
Vcap: 4000 V
Corona: 4.0 uA

Sample: Amitrol in water, 100 μ L

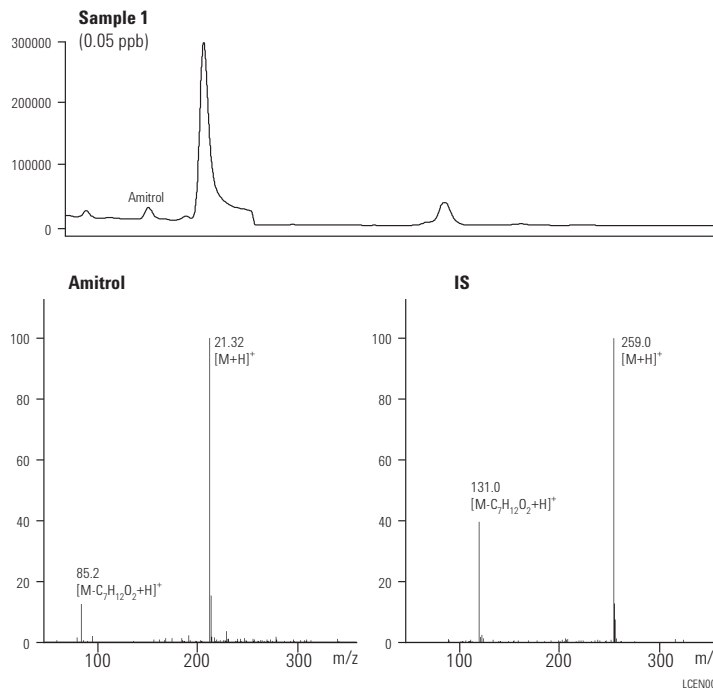**Anilines, substituted: Rapid separation**

Column: ZORBAX Rx/SB-C8
866953-906

4.6 x 75 mm, 3.5 μ m


Mobile Phase: 20% ACN/80% 25 mM phosphate buffer, pH 2.5

Flow Rate: 3.0 mL/min

Temperature: 60 °C

Detector: UV, 254 nm

Sample: Anilines

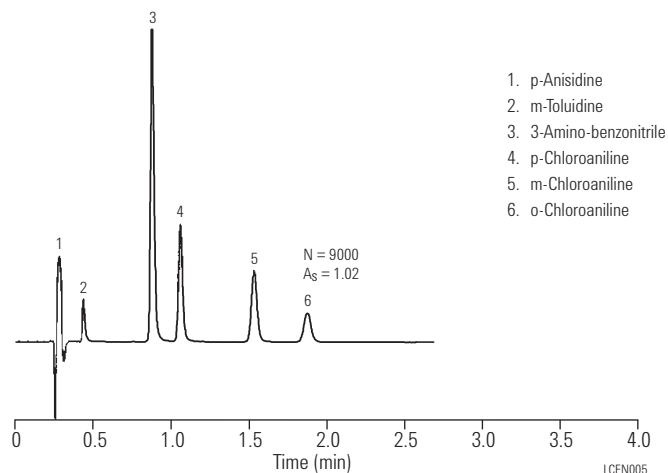

**Explosives and related compounds:
Qualitative and quantitative analysis**

Column A: ZORBAX SB-C18
883700-922
2.1 x 150 mm, 5 µm

Column B: ZORBAX SB-CN
883700-905
2.1 x 150 mm, 5 µm

Mobile Phase: A = ACN + 5% H₂O + 5 mM CF₃COONH₄
B = H₂O + 5% ACN + 5 mM CF₃COONH₄,
pH 2.7 (CF₃COOH)

Flow Rate: 0.23 mL/min

Gradient: A:
0 min 80% B
2 min 80% B
10 min 70% B
20 min 65% B
25 min 60% B
35 min 30% B
40 min 30% B
42 min 80% B


B:
0 min 80% B
1 min 80% B
15 min 70% B
30 min 20% B
35 min 20% B
37 min 80% B

Temperature: 18 °C

Detector: UV, 210, 240, 360 nm, wavelength
switching for each compound

Sample: 10 µL of 19 explosive compounds
in ACN/H₂O (20/80)

- | | |
|-------------------------------|--------------------------------|
| 1. Picric acid | 11. 4-Amino-4,6-dinitrotoluene |
| 2. 4-Amino-2-nitrotoluene | 12. 2-Nitrotoluene |
| 3. 2-Amino-6-nitrotoluene | 13. 2,6-Dinitrotoluene |
| 4. RDX | 14. 4-Nitrotoluene |
| 5. 2-Amino-4-nitrotoluene | 15. 3-Nitrotoluene |
| 6. HMX | 16. 2,4,6-Trinitrotoluene |
| 7. 1,3-Dinitrobenzene | 17. Tetryl |
| 8. 1,3,5-Trinitrobenzene | 18. Diphenylamine |
| 9. 2-Amino-4,6-dinitrotoluene | 19. Hexyl |
| 10. 2,4-Dinitrotoluene | |

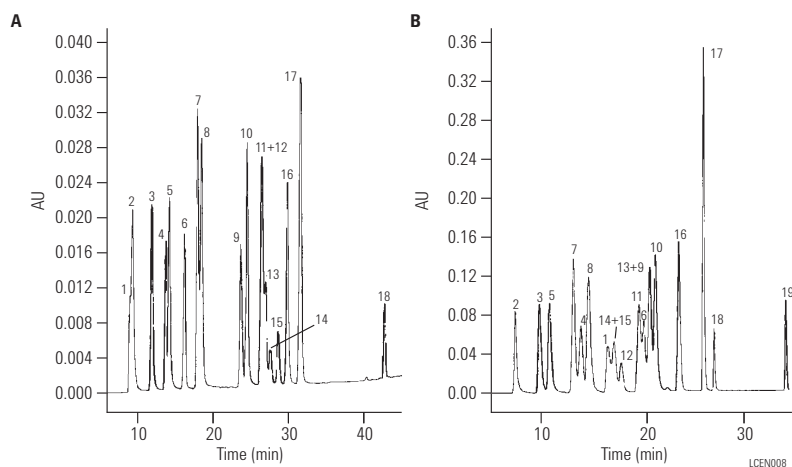

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Explosives from soil extract

Column: ZORBAX SB-C18
880975-302
3.0 x 250 mm, 5 µm


Mobile Phase: Methanol/Water (50/50) (v/v)

Flow Rate: 0.3 mL/min

Temperature: Ambient

Detector: UV, 230 nm

Sample: 10 µL explosives mix

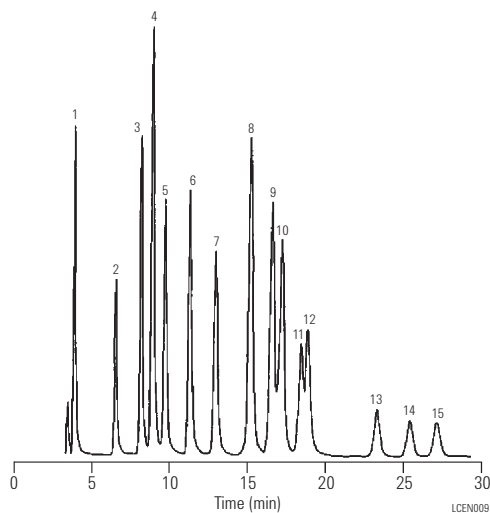

1. Octogen (HMX)
2. Hexogen (RDX)
3. 2-Amino-6-nitrotoluene
4. 1,3,5-Trinitrobenzene
5. 2-Amino-4-nitrotoluene
6. 1,3-Dinitrobenzene
7. Tetryl
8. 2,4,6-Trinitrotoluene
9. 4-Amino-2,6-dinitrotoluene
10. 2-Amino-4,6-dinitrotoluene
11. 2,6-Dinitrotoluene
12. 2,4-Dinitrotoluene
13. 2-Nitrotoluene
14. 4-Nitrotoluene
15. 3-Nitrotoluene

Herbicides on different bonded phases

Column A: ZORBAX SB-CN
883975-905
4.6 x 150 mm, 5 µm

Column B: ZORBAX SB-Phenyl
883975-912
4.6 x 150 mm, 5 µm

Column C: ZORBAX SB-C8
883975-906
4.6 x 150 mm, 5 µm


Mobile Phase: 35% ACN, 65% Water

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Herbicides

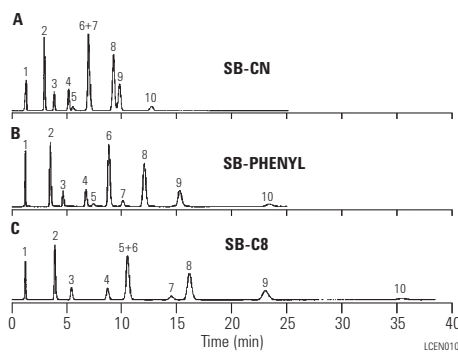

1. Bentazon
2. Tebuthiuron
3. Simazine
4. Atrazine
5. Prometon
6. Diuron
7. Propazine
8. Propanil
9. Prometryne
10. Metolachlor

Herbicide/pesticide standards: Effect of bonded phase

Column: Eclipse XDB-C8
993967-906
4.6 x 150 mm, 5 µm

Mobile Phase: Water/Acetonitrile

Flow Rate: 1.0 mL/min

Gradient: 20-60% in 15 min

Temperature: 50 °C
40 °C
30 °C
20 °C

Detector: DAD 240

Sample: Herbicide & pesticide standards

Column: Eclipse XDB-C18
993967-902
4.6 x 150 mm, 5 µm

Mobile Phase: Water/Acetonitrile



Flow Rate: 1.0 mL/min

Gradient: 20-60% in 15 min

Temperature: 50 °C
40 °C
30 °C
20 °C

Detector: DAD 240

Sample: Herbicide & pesticide standards

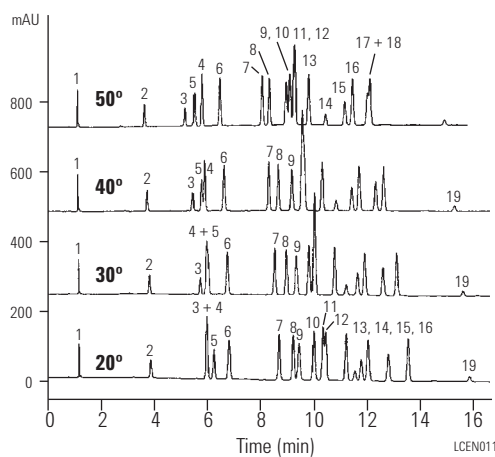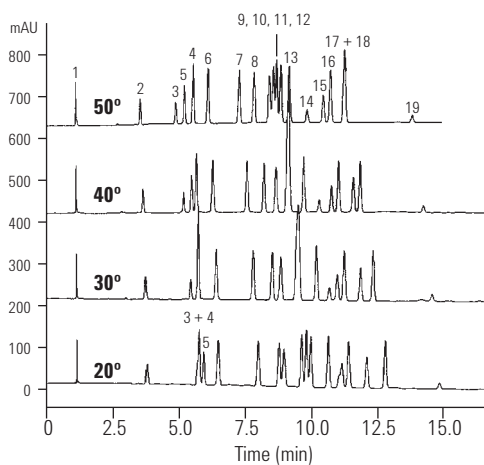

1. Desethyldeisopropylatrazine
2. Desethylatrazine
3. Benzthiazuron
4. Hexazinon
5. Metoxuron
6. Simazine
7. Methabenzthiazuron
8. Simazine
9. Atrazine
10. Isoproturon
11. Diuron
12. Monoluron
13. Metobromuron
14. Metazachlor
15. Propazine
16. Sebutylazine
17. Terbutylazine
18. Linuron
19. Metolachlor



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Separation of EPA 610 PAH Mix

Column: Eclipse PAH
959990-318
3.0 x 250 mm, 5 µm

Mobile Phase: A: Water
B: Acetonitrile
Initial %B = 40


Flow Rate: 0.85 mL/min

Gradient:

Time (Min)	%B
0.00	45
17.5	100
24.0	100
25.5	40
27.5	40
Stop Time = 25.0	

Temperature: 25 °C

Detector: 220, 4 nm No Ref.; Stop time = 26.0 min

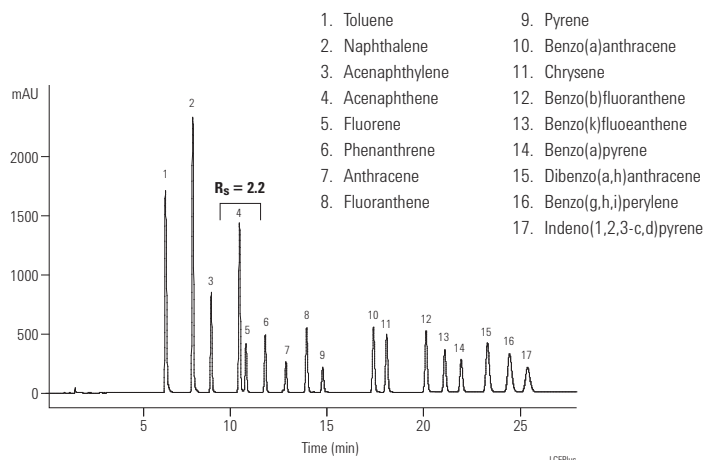

Polycyclic aromatic hydrocarbons according to EPA Method 610


Column: Pursuit PAH
A7001100X046
4.6 x 100 mm, 3 µm

Sample: NIST 16473 Standard

Mobile Phase: A: ACN:water, 25:75
B: ACN

Flow Rate: 2.0 mL/min

Detector: UV, 254 nm

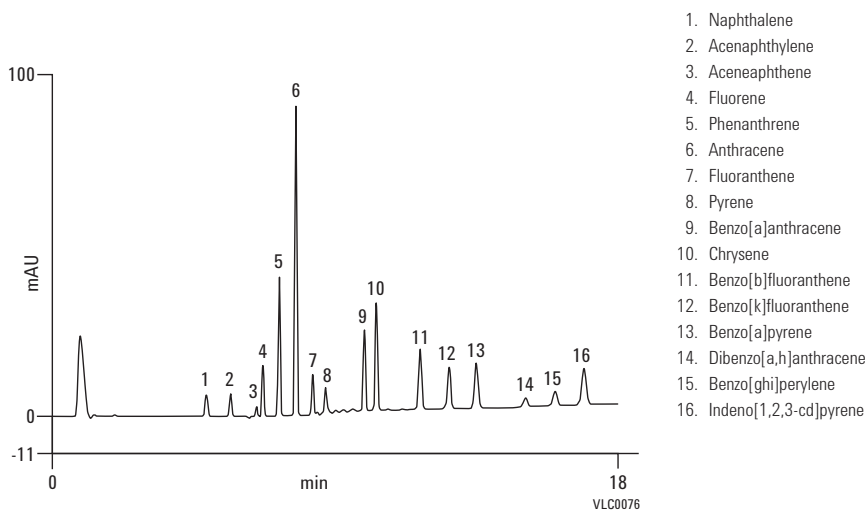

NEW!**Rapid method development for 18 PAH compounds with an Agilent RRHD Eclipse PAH column**

Column: ZORBAX RRHD Eclipse PAH
959758-918
2.1 x 100 mm, 1.8 μ m

Mobile Phase: A: Water
B: Acetonitrile

Flow Rate: 0.84 mL/min

Gradient: 40-100% B, gradient time (t_g) varies from 1 to 20 min;
isocratic hold at 100% B for 2 min,
re-equilibrate column at 40% B for 3 min


Temperature: 25 $^{\circ}$ C

Detector: Agilent 1290 Infinity LC

MS Conditions: Sig = 220, 4 nm; Ref = Off

Sample: 0.5 μ L injection of diluted Agilent PAH Mixture
(P/N 8500-6035) spiked with thiourea as a v_0 marker

- | | |
|-----------------------------|-----------------------------|
| 1. Thiourea (V_0 marker) | 10. Pyrene |
| 2. Toluene | 11. Benzo(a)anthracene |
| 3. Naphthalene | 12. Chrysene |
| 4. Acenaphthylene | 13. Benzo(b)fluoranthene |
| 5. Acenaphthene | 14. Benzo(k)fluoranthene |
| 6. Fluorene | 15. Benzo(a)pyrene |
| 7. Phenanthrene | 16. Dibenzo(a,h)anthracene |
| 8. Anthracene | 17. Benzo(g,h,i)perylene |
| 9. Fluoranthene | 18. Indeno(1,2,3-c,d)pyrene |

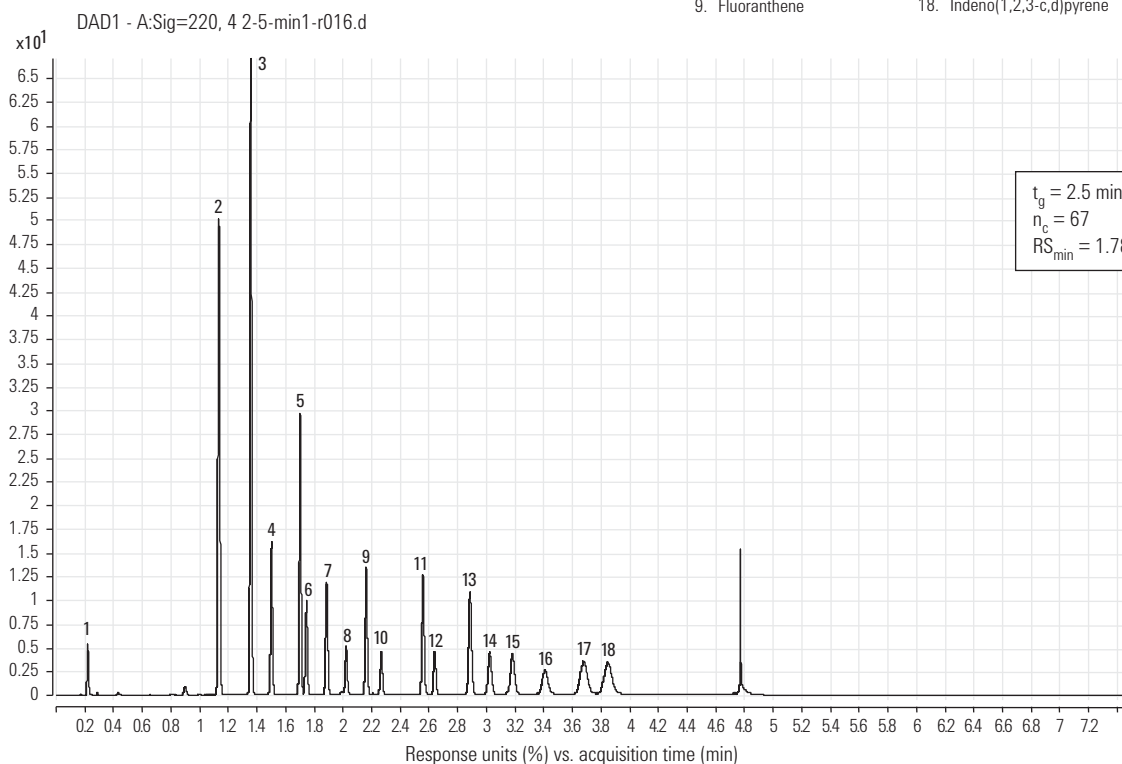

Gradient times are rapidly screened for the separation of 18 compounds.

Separation of 20 PAHs on Eclipse PAH

Column: Eclipse PAH
959964-918
4.6 x 100 mm, 1.8 µm

Mobile Phase: A: Water
B: Acetonitrile

Flow Rate: 1.8 mL/min


Gradient:

Time (Min)	% B
0	40
6	100
9.5	100
10	40

Stop Time = 12

Temperature: 25 °C

Detector: 230, 8 nm No Ref.; Data rate 0.2 s,
micro flow cell

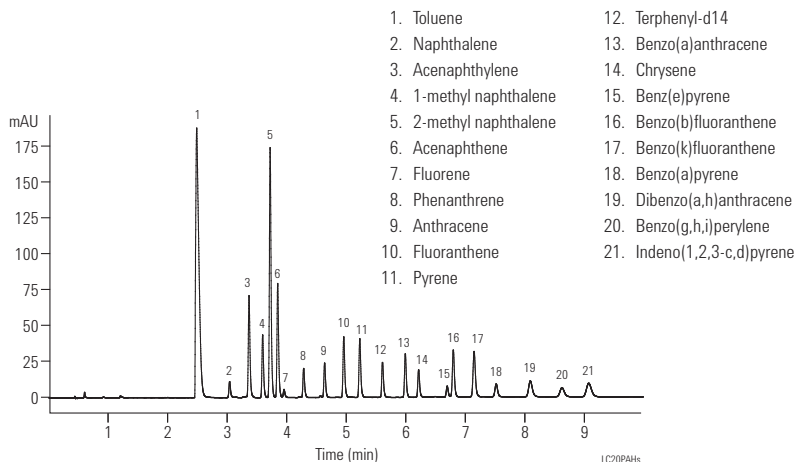
**Herbicides: Rapid separation**

Column: Eclipse XDB-C18
933975-902
4.6 x 30 mm, 3.5 µm

Mobile Phase: MeOH:H₂O (60:40)

Flow Rate: 2 mL/min

Temperature: Ambient

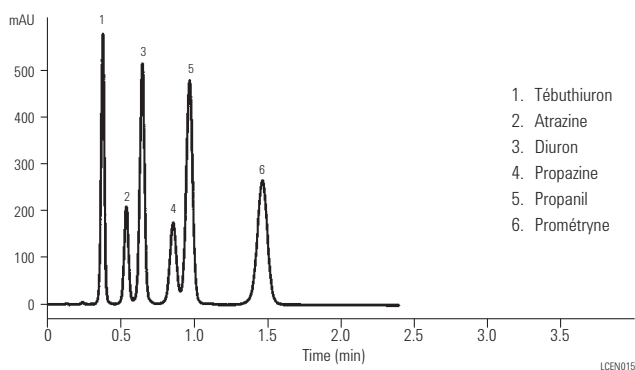**Phenoxyacid herbicides**


Column: Pursuit XRs C8
A6010150X046
4.6 x 150 mm, 5 µm

Mobile Phase: MeCN:water+0.1% HCOOH, 50:50

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 220 nm

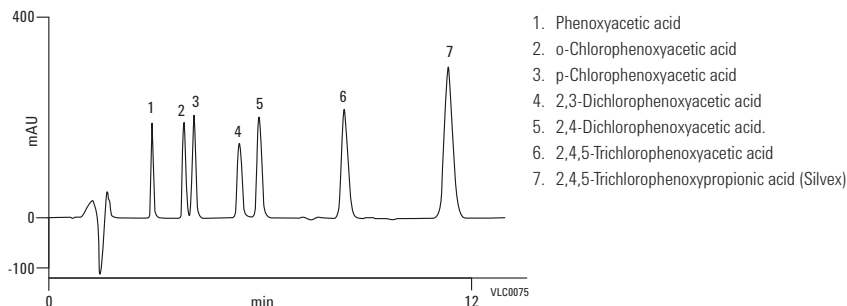

Triazine pesticides on Bonus-RP and Alkyl C8 phase

Column: ZORBAX Bonus-RP
883668-901
4.6 x 150 mm, 5 µm

Mobile Phase: MeOH: 0.1% TFA (70:30)*



Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: 254 nm

Sample: Triazine pesticides, 2 µL

1. Prometryne
2. Tebuthiuron
3. Atrazine
4. Propazine
5. Diuron
6. Propanil
7. Dacthal

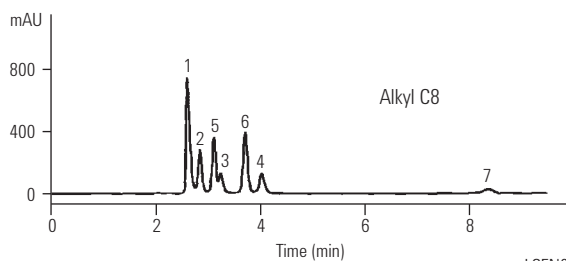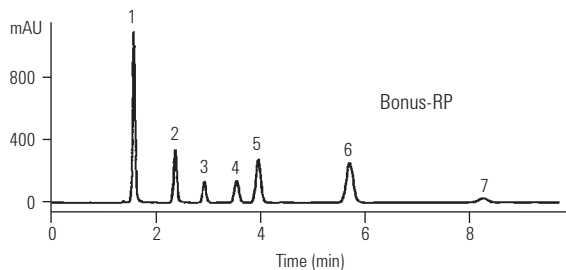

LCEN01

* For low pH work with Bonus-RP, a TFA mobile phase is often preferred over phosphate, and is compatible with LC/MS.

Phenols, substituted

Column: ZORBAX SB-C18
883975-902
4.6 x 150 mm, 5 µm

Mobile Phase: 20% ACN/80% 0.01 M H₃PO₄ to 45% ACN in 7.5 min


Flow Rate: 1.5 mL/min

Gradient: 80% ACN in 2.0 min

Temperature: 35 °C

Detector: UV, 254 nm

Sample: Phenols

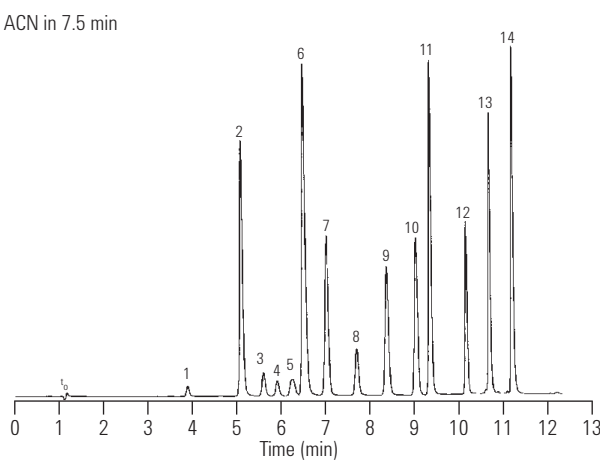

LCEN020

1. Phenol
2. 4-Nitrophenol
3. m-Cresol
4. o-Cresol
5. 2-Chlorophenol
6. 2,4-Dinitrophenol
7. 2-Nitrophenol
8. 2,4-Dimethylphenol
9. 4-Chloro-3-methylphenol
10. 2,4-Dichlorophenol
11. 2-Methyl-4,6-dinitrophenol
12. 2,4,6-Trichlorophenol
13. 2,3,4,6-Tetrachlorophenol
14. Pentachlorophenol

Plant hormones: Rapid gradient elution separation

Column: ZORBAX Rx/SB-C8
866953-906
4.6 x 75 mm, 3.5 µm


Mobile Phase: A: Water with 0.1% TFA
B: Acetonitrile with 0.1% TFA

Flow Rate: 3.0 mL/min

Temperature: 60 °C

Detector: UV, 245 nm

Sample: Plant hormones

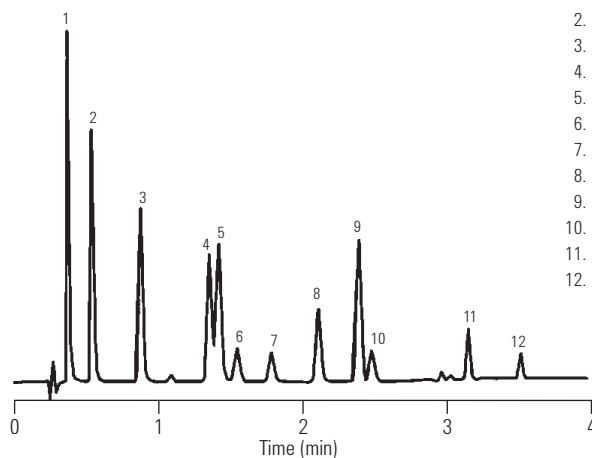

1. Kinetin
2. n-6-Benzyl adenine
3. 3-Indole acetic acid
4. 1-Naphthyl acetamide
5. 3-Indole propionic acid
6. o-Chlorophenoxy acetic acid
7. p-Chlorophenoxy acetic acid
8. 3-Indole butyric acid
9. 1-Naphthyl acetic acid
10. o-Chlorophenoxy propionic acid
11. 3,4,5-Trichlorophenoxy acetic acid
12. 3,4,5-Trichlorophenoxy propionic acid

LCEN022

VX nerve agent metabolites by LC/MS-IS standard (C13 labeled)

Column: ZORBAX NH2
860700-708
2.1 x 50 mm, 5 µm


Mobile Phase: 1:1 (20 mM Ammonium Acetate pH 4.5/Acetonitrile)

Flow Rate: 0.5 mL/min, 1 µL injection (prepared std in ACN)

Temperature: 35 °C

Detector: ESI-Negative Ion, Gas Flow 12 L/min, Nebulizer 60 psi

Sample	MW
1. Cyclohexyl methylphosphonic acid	178
2. Pinacolyl methylphosphonic acid	180
3. Isopropyl methylphosphonic acid	138
4. Ethyl methylphosphonic acid	124

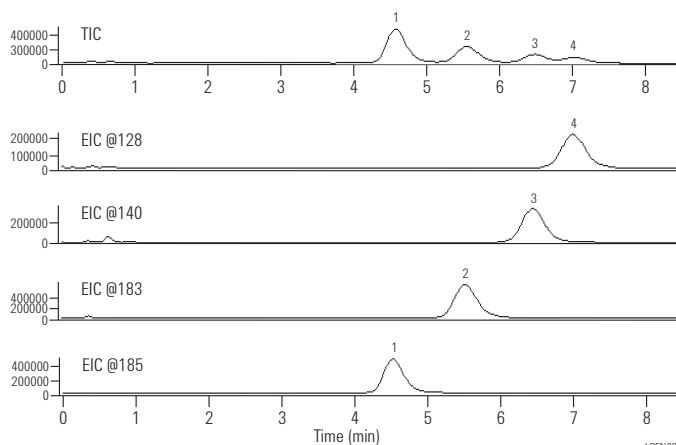

LCEN025



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Food and Consumer Product Applications

NEW!**Blueberry anthocyanin analysis**

Column A: Poroshell 120 SB-C18
687975-902
4.6 x 75 mm, 2.7 μ m


Column B: ZORBAX SB-C18
863953-902
4.6 x 150 mm, 3.5 μ m

Column C: ZORBAX SB-C18
880975-902
4.6 x 250 mm, 5 μ m

Flow Rate: 1 mL/min

Detector: Agilent 1260 Rapid Infinity LC

Blueberry anthocyanin analysis on totally porous and superficially porous StableBond C18 columns. Overlay of anthocyanin method with 250 mm 5 μ m, 150 mm 3.5 μ m, and 75 mm 2.7 μ m at 1 mL/min.

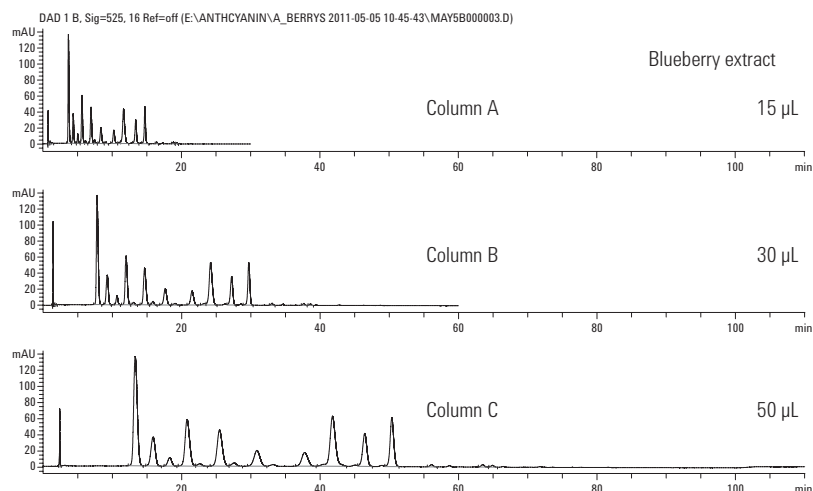**NEW!****Analysis of pesticide residues in green tea**

Column: Poroshell 120 EC-C18
695775-902
2.1 x 100 mm, 2.7 μ m

Mobile Phase: A: 5 mM FA in water
B: 5 mM FA in ACN


Flow Rate: 0.4 mL/min

Gradient: 5% B in 1 min, 50% B in 3 min,
90% B in 7 min, 90% B in 8 min,
5% B in 8.2 min, 5% B in 9 min

Temperature: 30 $^{\circ}$ C


MRM chromatograms of 50 ng/g fortified sample processed by EN method.

- | | |
|------------------|---------------------|
| 1. Acephate | 7. Propoxur |
| 2. Pymetrozine | 8. Carbaryl |
| 3. Carbenfazole | 9. Cyprodinil |
| 4. Thiabendazole | 10. Ethoprophos |
| 5. Imidacloprid | 11. Penconazole |
| 6. Imazalil | 12. Kresoxim-methyl |
| | IS TPP |

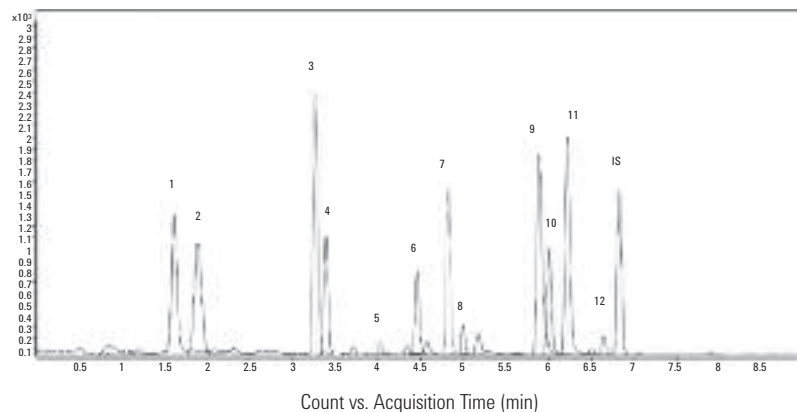

NEW!

**An overlay of the original ZORBAX Eclipse Plus 5 μ m method and Agilent Poroshell 120 method.
All 11 peaks on Poroshell 120 are resolved by the time the first peak elutes on the original
5 μ m ZORBAX Eclipse Plus method**

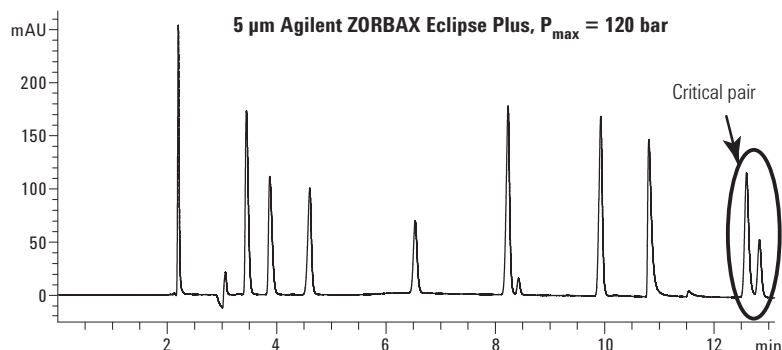


Column: Eclipse Plus C18
959990-902
4.6 x 250 mm, 5 μ m

Mobile Phase: A: 20 mM ammonium acetate, pH 4.80
B: acetonitrile

Flow Rate: 1.000 mL/min

Gradient: 14% B at t_0 , ramp to 52% B in 12.0 min

Temperature: 30 °C

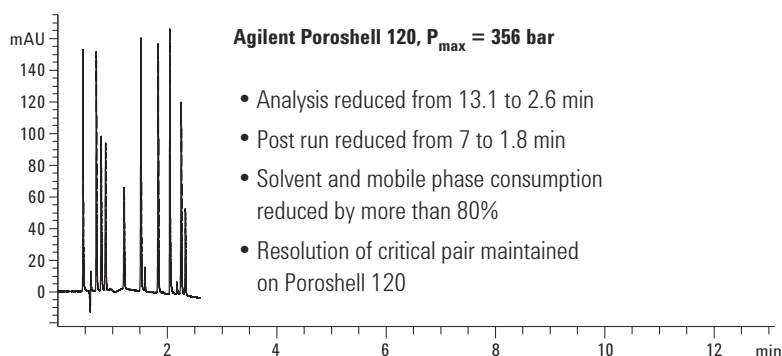

Column: Poroshell 120 EC-C18
695975-302
3.0 x 100 mm, 2.7 μ m

Mobile Phase: A: 20 mM ammonium acetate, pH 4.80
B: acetonitrile

Flow Rate: 0.851 mL/min

Gradient: 14% B at t_0 , ramp to 52% B in 2.1 min

Temperature: 30 °C



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

NEW!**Fast analysis of sulfa drugs**

Column: Eclipse Plus C18
959990-902
4.6 x 250 mm, 5 μ m


Column: Poroshell 120 EC-C18
695975-902
4.6 x 100 mm, 2.7 μ m

Gradient: Formic acid/acetonitrile

Detector: Agilent 1100 Series LC

Sample: Ten sulfa drugs

A separation of ten sulfa drugs scaled from an Agilent ZORBAX Eclipse Plus C18 column to an Agilent Poroshell 120 EC-C18 column showing analysis time decreased from 30 min to 8 min using a formic acid/acetonitrile gradient.

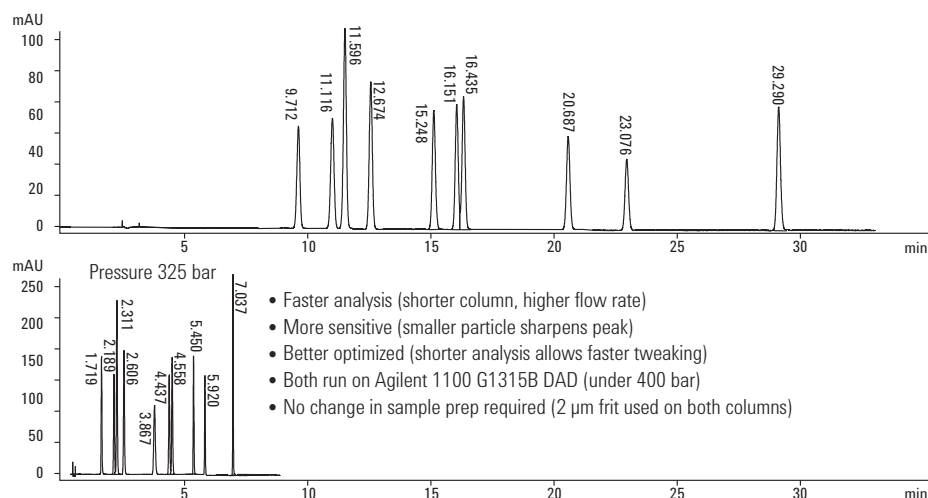

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

NEW!**Determination of anthocyanins in blueberries**

Column: ZORBAX RRHD Eclipse Plus C18
959758-902
2.1 x 100 mm, 1.8 μ m

Column: ZORBAX RRHD Eclipse Plus Phenyl-Hexyl
959758-912
2.1 x 100 mm, 1.8 μ m

Column: ZORBAX RRHD SB-Aq
858700-914
2.1 x 100 mm, 1.8 μ m

Column: ZORBAX RRHD SB-Phenyl
858700-912
2.1 x 100 mm, 1.8 μ m

Mobile Phase: A: 5% HCOOH in H₂O
B: CH₃CN

Flow Rate: 0.65 mL





Gradient: 10-50% B in 15 min

Detector: Agilent 1290 Infinity LC

MS Conditions: DAD: Sig = 525, 8 nm; Ref = Off
MS2 Scan: ESI + 200-1000
Scan time: 100 ms, 0.2 amu step
Fragmentor: 180 V
Drying gas: 10 L/min, 350 °C
Nebulizer Pressure: 50 psig
Capillary Voltage: 3500

Sample: 5 μ L injection of blueberry extract

1. Cyanidin, m/z 286
2. Peonidin, m/z 300
3. Delphinidin, m/z 302
4. Petunidin, m/z 316
5. Malvidin, m/z 330

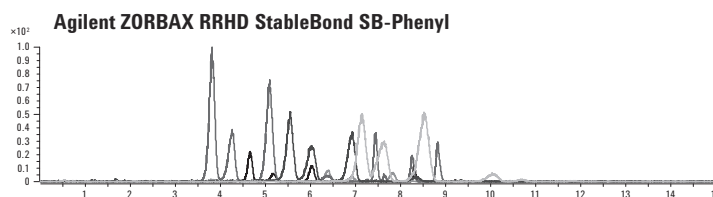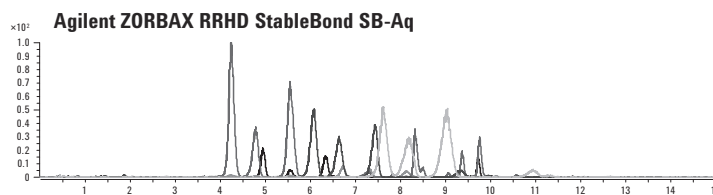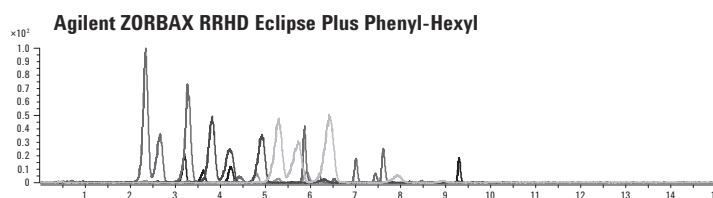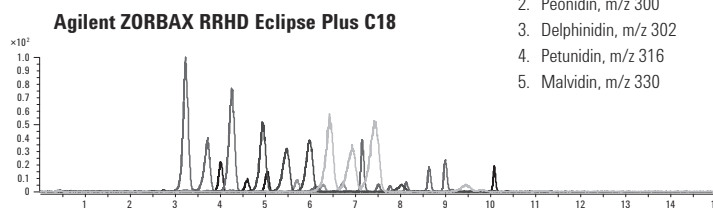

Counts (%) versus Acquisition time (min)

Separation of Azo Dyes


Column: Eclipse Plus Phenyl Hexyl
959996-912
4.6 x 100 mm, 5 µm

Mobile Phase: A: 10 mM Ammonium Acetate, pH 4.7
B: MeOH

Flow Rate: 1.5 mL/min

Gradient: Time (Min): %B:
0 25
5 50

Detector: UV, 254 nm

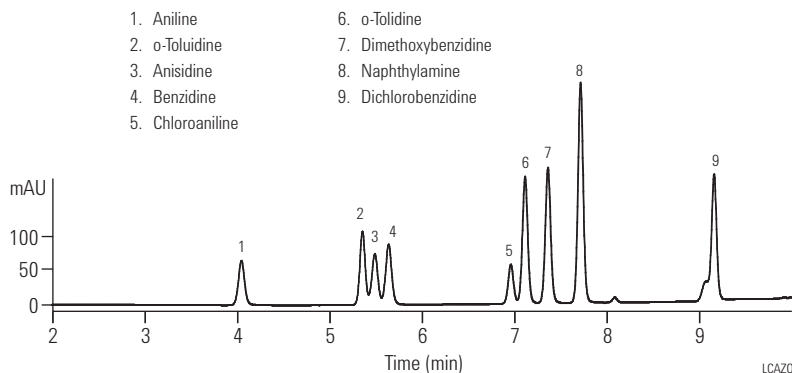
**Anthocyanins from blueberries:
High-efficiency high-speed separation**

Column A: ZORBAX SB-C18
880975-902
4.6 x 250 mm, 5 µm

Column B: ZORBAX SB-C18
863953-902
4.6 x 150 mm, 3.5 µm

Column C: ZORBAX SB-C18
866953-902
4.6 x 75 mm, 3.5 µm

Mobile Phase: A: 3% Phosphoric acid
B: 100% MeOH

Flow Rate: 1.0 mL/min


Gradient: As shown

Temperature: 30 °C


Detector: UV, 525 nm

Sample: Natural anthocyanins


Time	Percent B
0 min	23% B
35 min	26% B
97 min	60% B

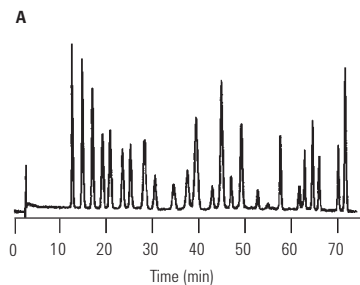

Time	Percent B
0 min	23% B
21 min	26% B
58.2 min	60% B

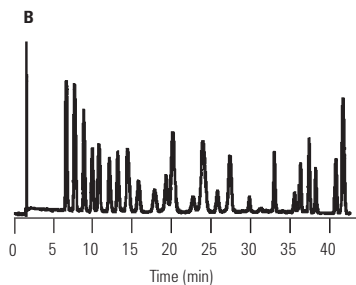

Time	Percent B
0 min	23% B
10.5 min	26% B
29.1 min	60% B

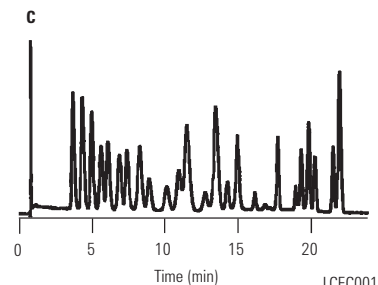

Aromatics II

Column: Eclipse XDB-Phenyl
963967-912
4.6 x 150 mm, 3.5 µm


Mobile Phase: H₂O: MeOH, 40:60

Flow Rate: 1.0 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

Sample: Aromatic Sample

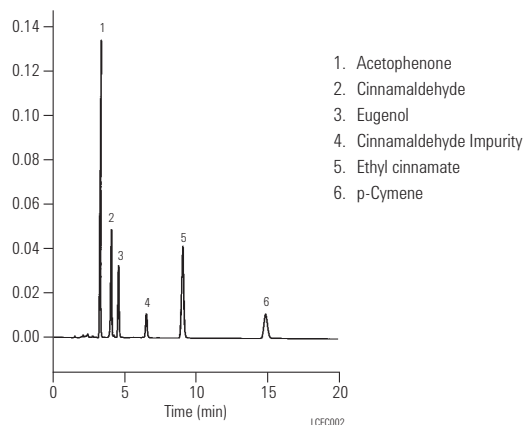**Aspartame: Metabolites and applications**

Column: ZORBAX SB-C18
866953-902
4.6 x 75 mm, 3.5 µm


Mobile Phase: 85/15, 0.1% TFA/ACN

Flow Rate: 1.0 mL/min


Temperature: 35 °C

Detector: UV, 210 nm


Sample: Aspartame

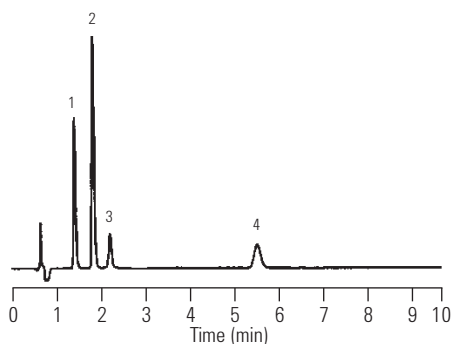**Aspartame and Its Metabolites**

1. Phenylalanine
2. 5-benzyl-3,6-dioxo-2-piperazineacetic acid
3. Aspartic acid-phenylalanine dipeptide
4. Aspartame

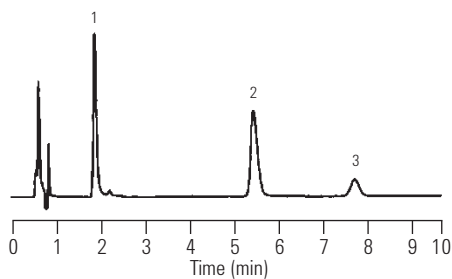**Diet Coke**

1. Caffeine
2. Aspartame
3. Unknown

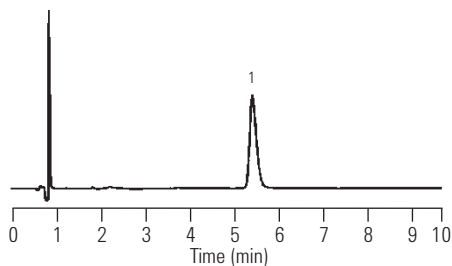**Equal Sweetener**


1. Aspartame

Carbohydrates: Carbohydrate standards


Column: ZORBAX Carbohydrate Analysis
843300-908
4.6 x 150 mm, 5 µm

Mobile Phase: 63% CH₃CN/H₂O
Flow Rate: 0.5 mL/min

Detector: Agilent RID
Sample: Carbohydrate standard:
A: 25 ng/ L, 1 µL injected
B: 500 pg/ L, 50 µL injected

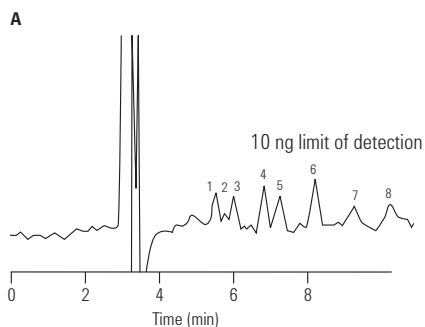Carbohydrates: Separation showing high sensitivity

1. Ribose
2. Rhamnose
3. Xylose
4. Fructose
5. Glucose
6. Sucrose
7. Maltose
8. Lactose
9. Raffinose

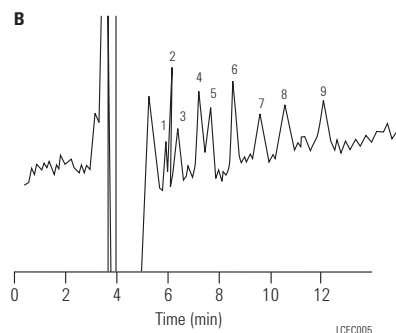Sensitivity of high injection volume (50 µL)**Carbohydrates: Effect of mobile phase strength**

Column: ZORBAX NH2
880952-708
4.6 x 250 mm, 5 µm


Mobile Phase: ACN/Water, as indicated
Flow Rate: 1.0 mL/min

Temperature: Ambient


Detector: RI

Sample: Mono- and Disaccharides

ACN/H₂O: 70/30

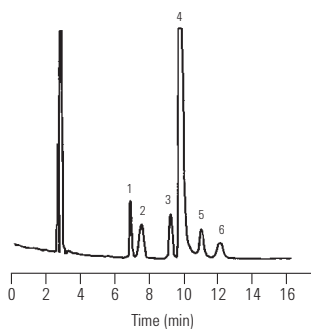

ACN/H₂O: 75/25

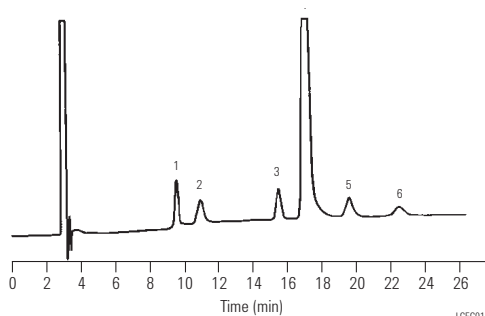

1. Fructose
2. Glucose
3. Saccharose
4. Palatinose
5. Trehalulose
6. Isomaltose

Carbohydrates in colas

Column: ZORBAX Carbohydrate Analysis
843300-908
4.6 x 150 mm, 5 µm


Mobile Phase: 75% ACN:25% H₂O

Flow Rate: 2.0 mL/min

Temperature: 30 °C

Detector: RID

Sample: No dilution
A: COLA, Fountain
B: COLA, Can, Brand A
C: COLA, Brand B
D: COLA, Brand B, diet

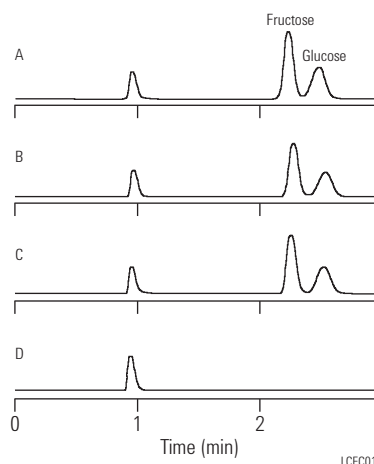

LCFC013

Carbohydrates: Sugar alcohols

Column: ZORBAX Carbohydrate Analysis
843300-908
4.6 x 150 mm, 5 µm


Mobile Phase: 75% ACN:25% H₂O

Flow Rate: 2.0 mL/min

Temperature: 30 °C

Detector: RID

Sample: Chewing gum, sugar-free

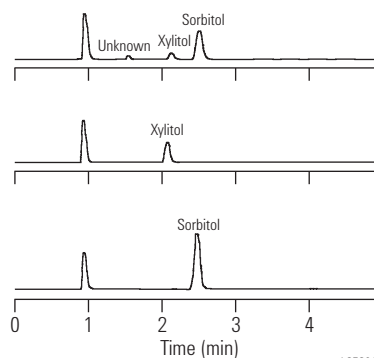

LCFC014

Carbohydrates in juices

Column: ZORBAX Carbohydrate Analysis
843300-908
4.6 x 150 mm, 5 µm


Mobile Phase: 75% ACN/25% H₂O

Flow Rate: 2.0 mL/min

Temperature: 30 °C

Detector: RID

Sample: Diluted to 0.1X in 50:50 ACN:H₂O

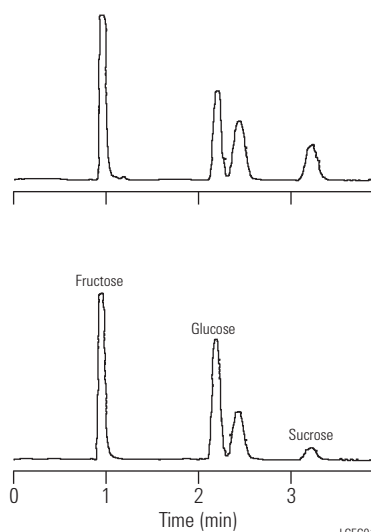

LCFC016

Apple Drink

36.8% Fructose
24.9% Sucrose
38.3% Glucose

Apple Juice

58.7% Fructose
9.9% Sucrose
33.4% Glucose

Carbohydrates in milk

Column: ZORBAX Carbohydrate Analysis
843300-908
4.6 x 150 mm, 5 µm

Mobile Phase: 75% ACN/25% H₂O

Flow Rate: 2.0 mL/min

Temperature: 30 °C


Detector: RID

Sample: Partitioned between CH₂Cl₂: H₂O

Milk (2%)



100% Lactose-Free Milk

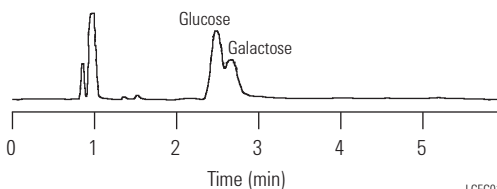

LCFC015

Flavoring agents

Column: ZORBAX SB-Phenyl
860975-912
2.1 x 50 mm, 5 µm


Mobile Phase: 0.3% TFA: ACN, 65:35

Flow Rate: 0.3 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Cool mint Listerine sample

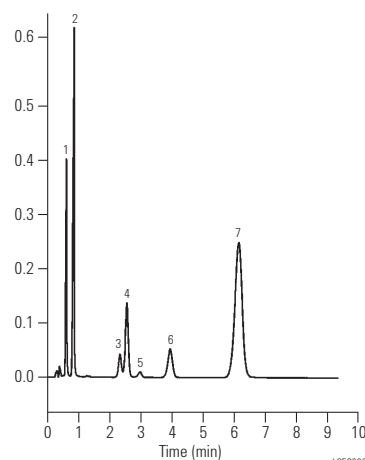

1. Unknown
2. Benzoic acid
3. Methyl salicylate
4. Carvone
5. Unknown
6. Thymol
7. Anethole

LCFC006

Food colors, FD&C

Column: ZORBAX Eclipse XDB-C18
935967-902
4.6 x 50 mm, 3.5 µm

Mobile Phase: A: 0.1% TFA, pH to 4.4 with TEA, B: MeOH


Flow Rate: 1.0 mL/min

Gradient: 17 to 100% B/4 min

Temperature: Ambient

Detector: UV, 254 nm

- | | | |
|------------------|-----------------|--------|
| 1. Yellow #5 | C16H9N4Na3O9S2 | MW=534 |
| 2. Red #40 | C18H14N2Na2O8S2 | MW=496 |
| 3. Blue #1 | C37H34N2Na2O9S3 | MW=760 |
| 4. Propylparaben | C10H12O3 | MW=180 |
| 5. Red #3 | C20H414Na2O5 | MW=878 |

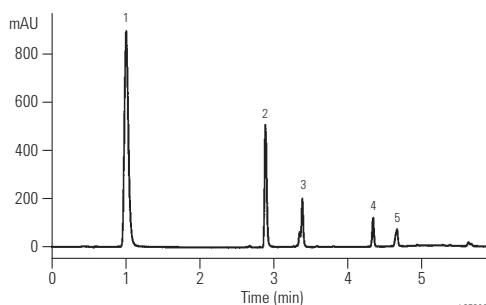

LCFC007

Neutraceuticals: Extract from green tea

Column: ZORBAX SB-C8
863953-906

4.6 x 150 mm, 3.5 µm

Mobile Phase: 75% 0.1% Trifluoroacetic acid: 25% Methanol


Injection: 1 mL/min

Temperature: 40 °C


Detector: UV, 280 nm

Sample: Green tea extract, 5 µL

Catechin Mixture

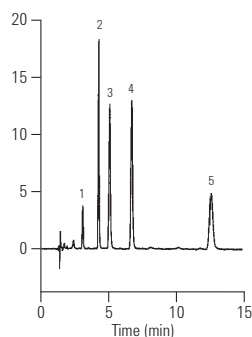

Green Tea Extract

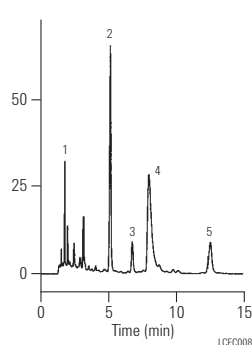

1. Epigallocatechin
2. Epicatechin
3. Epigallocatechin gallate
4. Catechol
5. Epicatechin gallate

Tocopherols by LC/MS with APPI

Column: Eclipse XDB-C18
993967-302

3.0 x 150 mm, 5 µm


Mobile Phase: 97% MeOH: 3% 10 mM CH₃COONH₄

Flow Rate: 0.5 mL/min

Temperature: 40 °C

MS Conditions: MS: Agilent 1100MSD SL
Ionization: APPI (Positive)
Scan range: m/z 100-500
Vcap: 1500 V
SIM ion: base peak
Drying gas: 7 L/min at 350 °C
Nebulizer gas: 60 psi
Vaporizer temp: 350 °C
Fragmentor: 140 V
EM gain: 4

Sample Volume: 10 µL

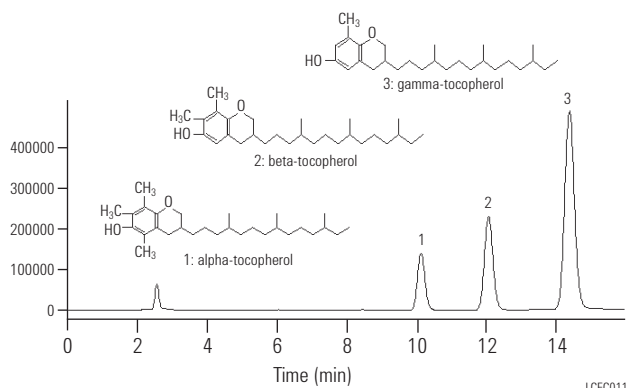

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Sugars in plain and milk chocolate


Column: Hi-Plex Pb
PL1170-6820
7.7 x 300 mm, 8 µm

Mobile Phase: Water

Flow Rate: 0.6 mL/min

Temperature: 80 °C

Detector: RI

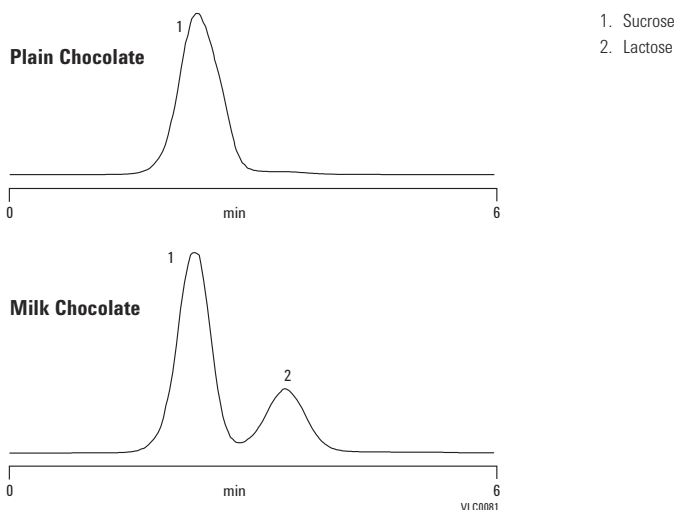

Sugars

Column: Hi-Plex K
PL1170-6860
7.7 x 300 mm, 8 µm


Sample: Sugars mixture (all 10 mg/mL), 20 µL injection

Mobile Phase: Water

Flow Rate: 0.6 mL/min

Temperature: 85 °C

Detector: 356-LC RI

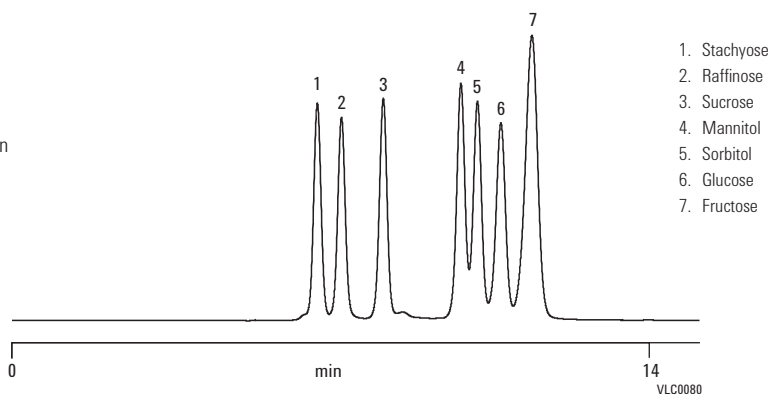

Parabens: High speed separation

Column: ZORBAX SB-C18 Rapid Resolution
Cartridge
833975-902
4.6 x 30 mm, 3.5 µm


Mobile Phase: 0.1% H₃PO₄:ACN, (50:50)

Flow Rate: 2 mL/min

Temperature: Top: ambient, bottom: 60 °C

Detector: UV, 254 nm with standard flow cell (13 µL)

Sample: Parabens, 1 µL

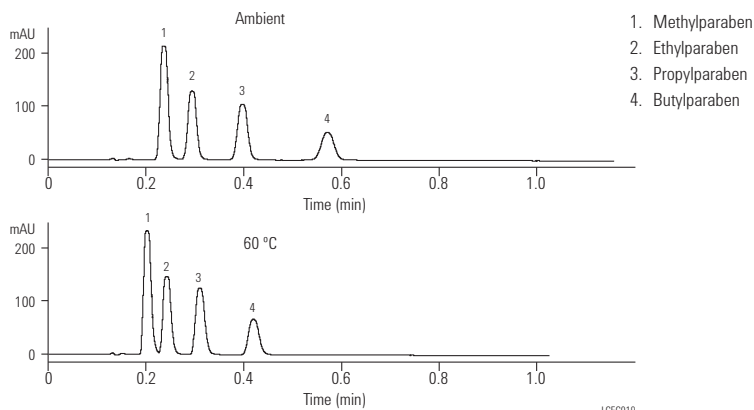

Separation of vitamin D2/D3


Column: Eclipse PAH
959941-918
4.6 x 50 mm, 1.8 µm

Mobile Phase: 92% MeOH, 8% water

Flow Rate: 2 mL/min

Temperature: 40 °C

Detector: 325 nm for VA/280 nm for VD and VE

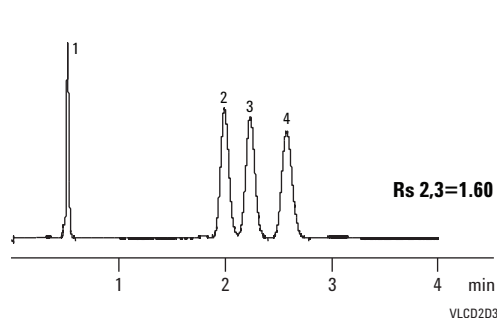

1. Vitamin A
2. Vitamin D2
3. Vitamin D3
4. Vitamin E (a-VE)

Fat-soluble vitamins on ZORBAX Eclipse XDB-C8

Column: Eclipse XDB-C8
993967-906
4.6 x 150 mm, 5 µm


Mobile Phase: 5/95 Water/MeOH

Flow Rate: 1.0 mL/min

Temperature: A: Ambient
B: 50 °C

Detector: UV, 280 nm

Sample: Fat-soluble vitamins

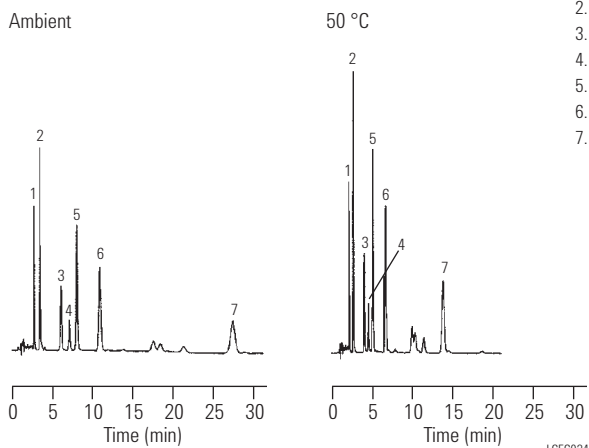

1. Retinol
2. Retinol acetate
3. Vitamin D3
4. γ-Tocopherol
5. α-Tocopherol
6. Tocopherol acetate
7. Retinol palmitate

Water-soluble vitamins

Column: ZORBAX SB-C8
883975-906
4.6 x 150 mm, 5 µm

Mobile Phase: A: 50 mM Sodium Phosphate, pH 2.5/MeOH (90/10)
B: 50 mM Sodium Phosphate, pH 2.5/MeOH (10/90)


Flow Rate: 1.0 mL/min

Gradient: 0-70% B in 18 min

Temperature: Ambient

Detector: UV, 245 nm

Sample: Water-soluble vitamins

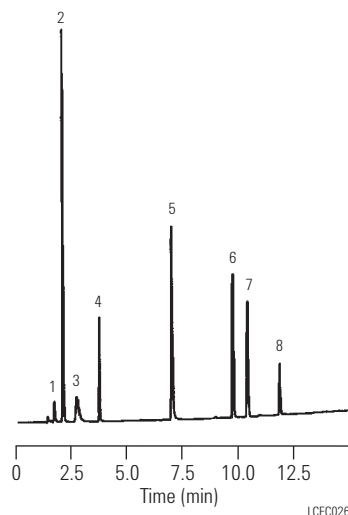

1. B₁-Thiamine
2. Vitamin C
3. B₃-Niacin
4. B₆-Pyridoxine
5. Pantothenic acid
6. Folic acid
7. B₁₂-Cyanocobalamin
8. B₂-Riboflavin

**Water-soluble vitamins:
High speed separation using ion-pairing**

Column: ZORBAX Rx/SB-C8
866953-906
4.6 x 75 mm, 3.5 µm


Mobile Phase: 10 mM Hexane Sulfonate with 0.1%
Phosphoric Acid: MeOH (74:26)

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 245 nm

Sample: Water-soluble vitamins

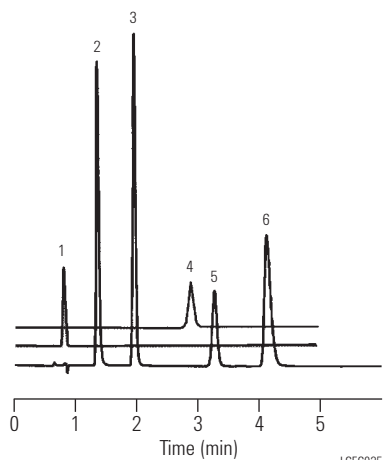

1. Vitamin C
2. B₃-Niacin
3. B₆-Pyridoxine
4. Folic acid
5. B₂-Riboflavin
6. B₁-Thiamine

**Water-soluble vitamins using the USP 23
method**

Column: ZORBAX SB-C18
880975-902
4.6 x 250 mm, 5 µm


Mobile Phase: 7.2 mM Hexane Sulfonate/MeOH/Acetic Acid
(73/27/1) (ratio to 101)

Flow Rate: 1.0 mL/min

Temperature: 30 °C

Detector: UV, 280 nm

Sample: Water-soluble vitamins

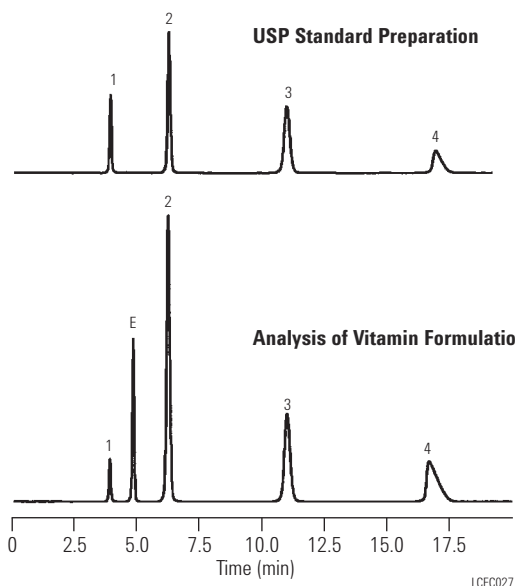

1. B₃-Niacin
2. B₆-Pyridoxine
3. B₂-Riboflavin
4. B₁-Thiamine
- E. Excipient

Water-soluble B vitamins separated on ZORBAX SB-Aq


Column: ZORBAX SB-Aq
883975-914
4.6 x 150 mm, 5 µm

Mobile Phase: 5% MeOH/95% water (0.1% TFA)

Flow Rate: 2.0 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

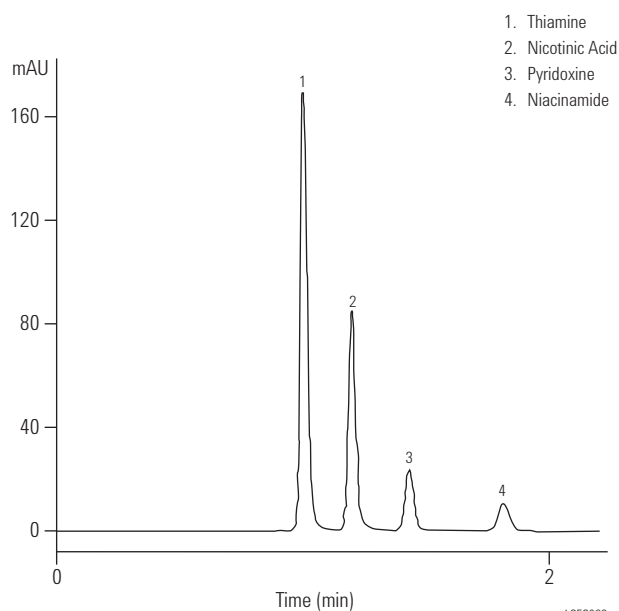

Sunscreen ingredients: Perform conventional, fast and ultra-fast separations on the same column family

Column A: Eclipse XDB-C18
993967-902
4.6 x 150 mm, 5 µm
6 µL inj

Column B: Eclipse XDB-C18
961967-902
4.6 x 100 mm, 3.5 µm
4 µL inj

Column C: Eclipse XDB-C18
927975-902
4.6 x 50 mm, 1.8 µm
2 µL inj


Mobile Phase: A: 15% water
B: 85% MeOH

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Sunscreens

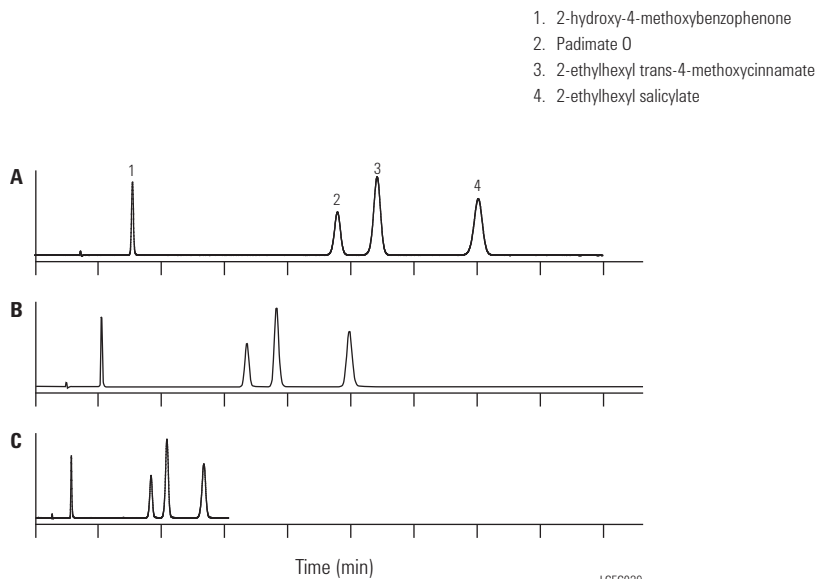

Fast vitamin E analysis on Rapid Resolution HT


Column A: Eclipse XDB-C18
927975-902
4.6 x 50 mm, 1.8 μ m

Column B: Eclipse XDB-C18
993967-902
4.6 x 150 mm, 5 μ m

Mobile Phase: A: 5% water
B: 95% MeOH

Flow Rate: 3 mL/min, 1 mL/min

Temperature: Ambient

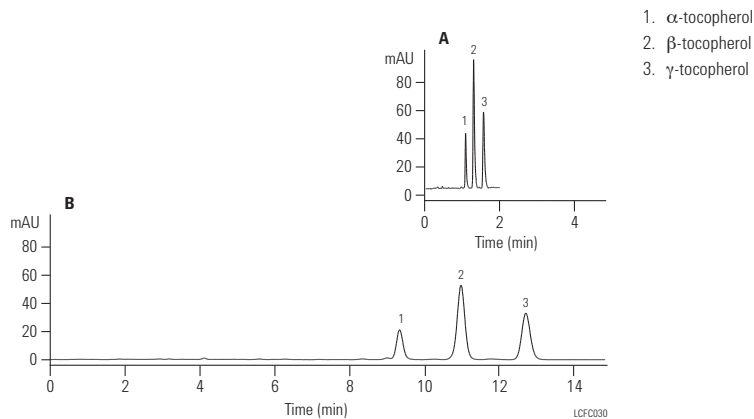

Theobromine in beverages

Column: ZORBAX SB-C18
827975-902
4.6 x 50 mm, 1.8 μ m


Mobile Phase: A: 92% 0.1% formic acid
B: 8% 0.1% formic acid in ACN

Flow Rate: 1.5 mL/min

Temperature: Ambient

Detector: UV, 254 nm, flow cell 2 μ L,
3 mm flow path

Sample: Theobromine

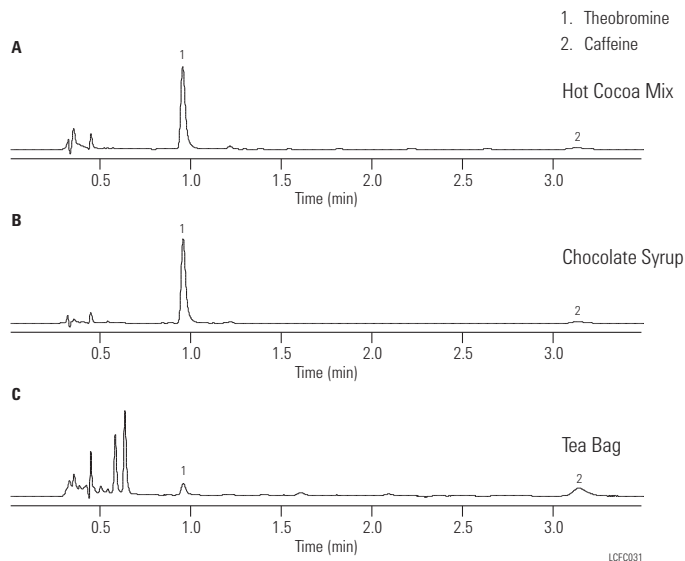


Benzoic acid/sorbic acid

Mobile Phase: 0.1% formic acid in water:
0.1% formic acid in MeCN, 80:20

Flow Rate: 0.7 mL/min

Detector: UV, 254 nm

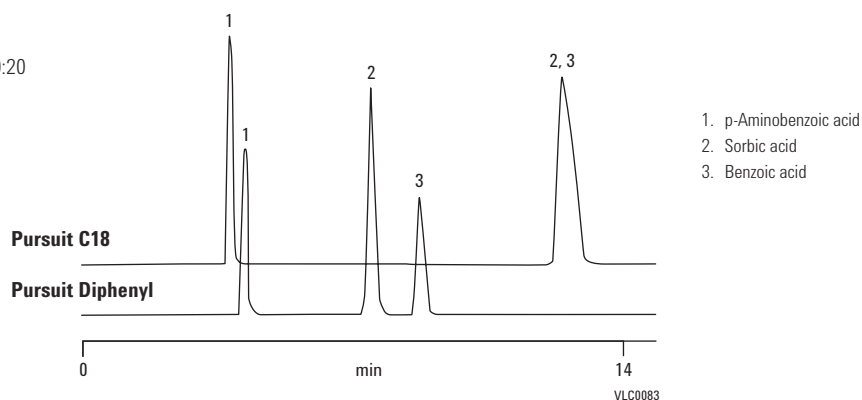**Quantification and qualification of vitamin C and citric acid in fresh grapefruit juice**


Column: **PLRP-S 100Å**
PL1512-5500
4.6 x 250 mm, 5 µm

Sample: Diluted 1:50 in eluent

Mobile Phase: 0.2M NaH₂PO₄, pH 2.14

Flow Rate: 0.5 mL/min

Detector: UV, 220 nm

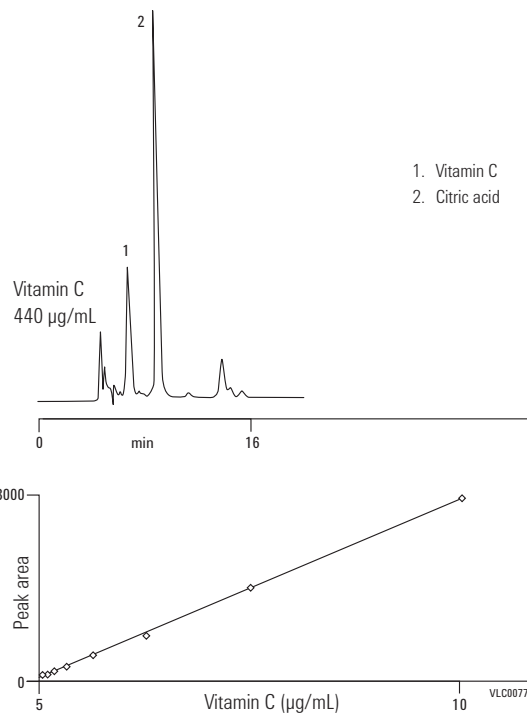

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Rose wine

Column: Hi-Plex H
PL1170-6830
7.7 x 300 mm, 8 µm


Mobile Phase: 0.004M H₂SO₄

Flow Rate: 0.4 mL/min

Pressure: 13 bar

Temperature: 75 °C

Detector: RI

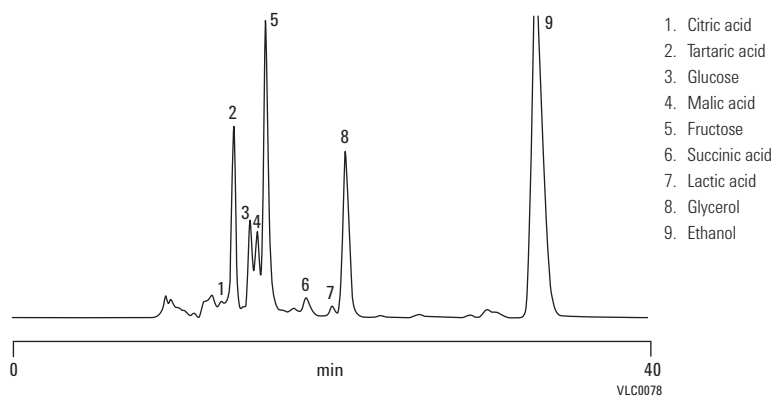

Sports drink

Column: Hi-Plex Na
PL1171-6140
7.7 x 300 mm, 10 µm


Sample: High energy orange flavor non-carbonated sports drink

Mobile Phase: Water

Flow Rate: 0.3 mL/min

Temperature: 80 °C

Detector: RI

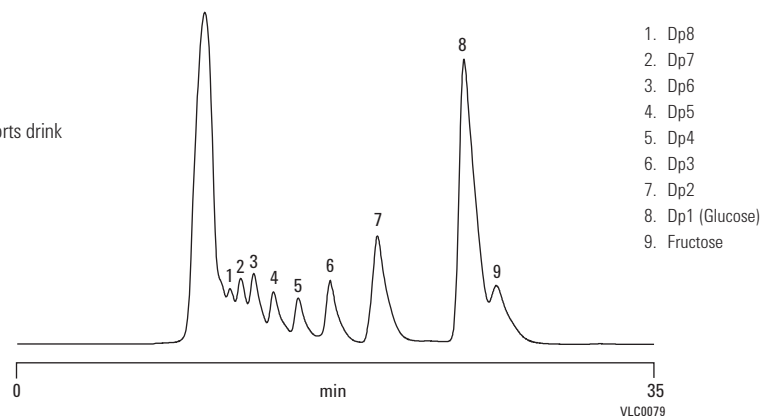

Oligosaccharides


Column: Hi-Plex Ca (Duo)
PL1F70-6850
6.5 x 300 mm, 8 µm

Mobile Phase: DI water

Flow Rate: 0.5 mL/min

Temperature: 90 °C

Detector: RI

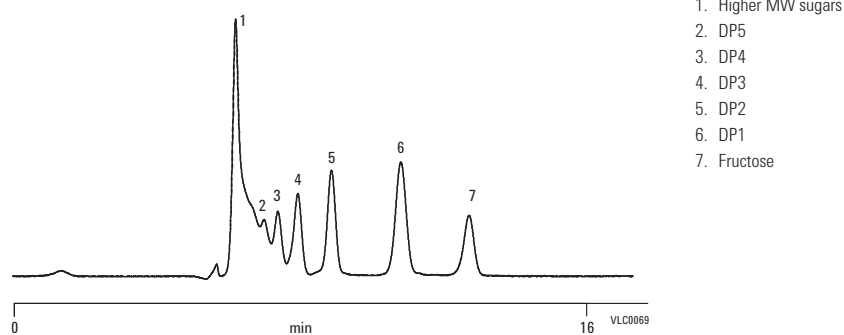

Pharmaceutical Applications


NEW!

Fast analysis of cefepime and related impurities

Column: Poroshell 120 EC-C18
697975-902
4.6 x 75 mm, 2.7 µm

Column: Eclipse Plus C18
959990-902
4.6 x 250 mm, 5 µm

Detector: Agilent 1200 Infinity Series

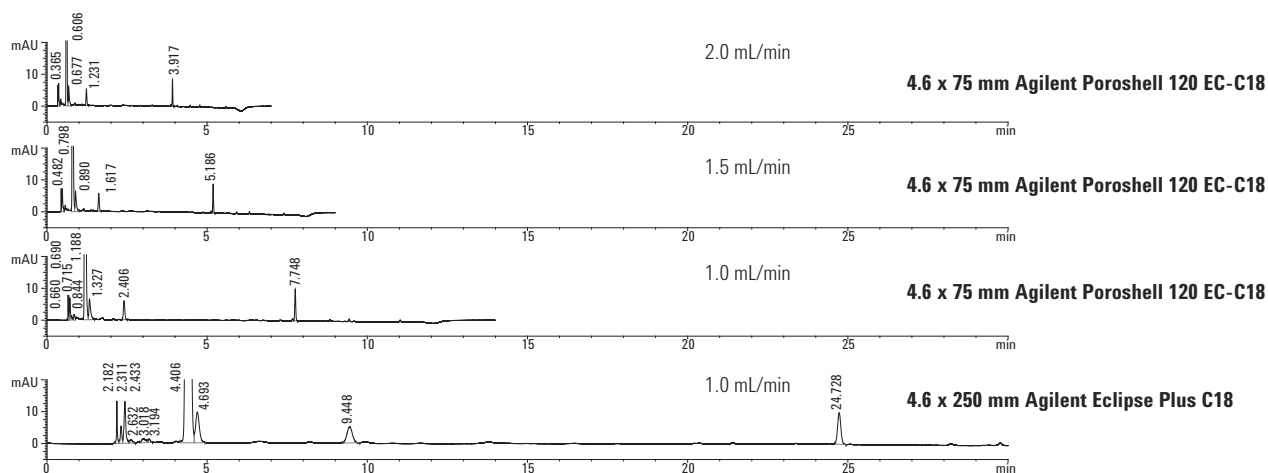

NEW!

Naproxen analysis

Column A: Eclipse Plus C18
959993-902
4.6 x 150 mm, 5 µm


Column B: Poroshell 120 EC-C18
697975-902
4.6 x 50 mm, 2.7 µm

Mobile Phase: 50:49:1 MeCN:H₂O:Glacial acetic acid

Flow Rate: 1.2 mL/min

Injection: Column A: 20 µL
Column B: 6.7 µL

Injection: Naproxen

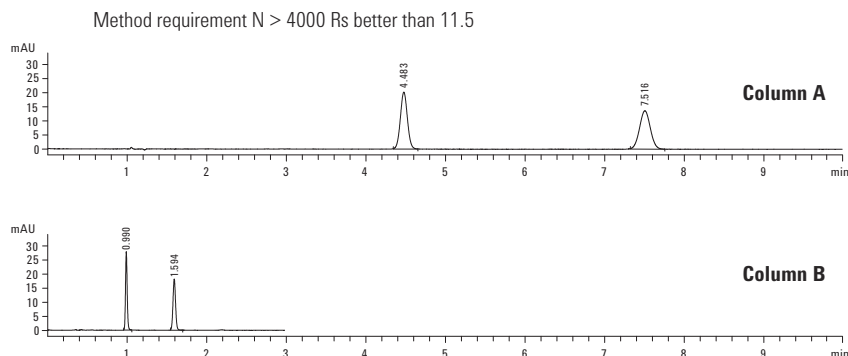

4-fold reduction in analysis time for this method when transferring to Poroshell 120.

NEW!


Analysis of water soluble vitamins in multivitamin tablets

Column: Poroshell 120 EC-C18
697975-902
4.6 x 75 mm, 2.7 µm

Flow Rate: 1.5 mL/min

Gradient: 0 min-1% B, 0.5 min-12% B,
0.52 min-30% B,
3.5 min-30% B, 4.5 min-1% B

Injection: 5 µL

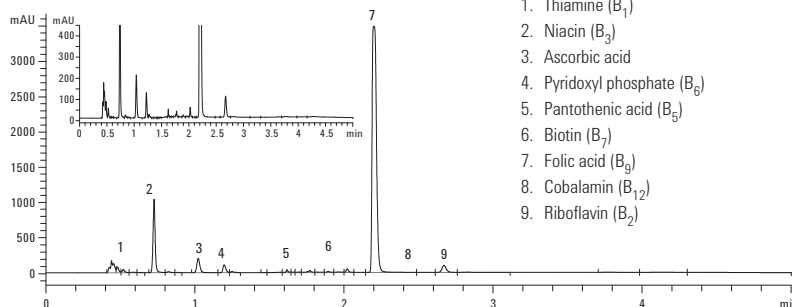

NEW!





Fast method for ginseng analyses scaled from a traditional method

Column: Eclipse Plus C18
959993-902
4.6 x 150 mm, 5 µm

Column: Poroshell 120 EC-C18
695975-302
3.0 x 100 mm, 2.7 µm

Detector: 1200 Infinity Series

Sample: Ginsenoside

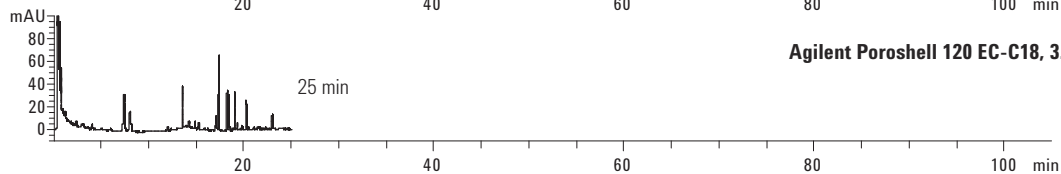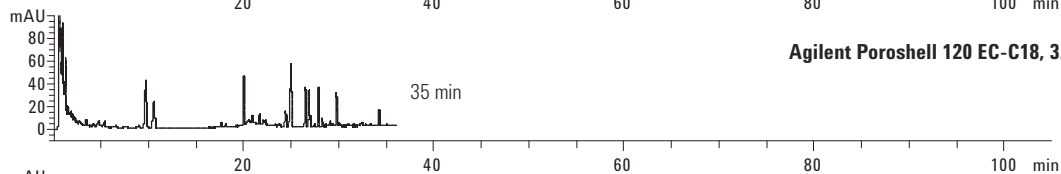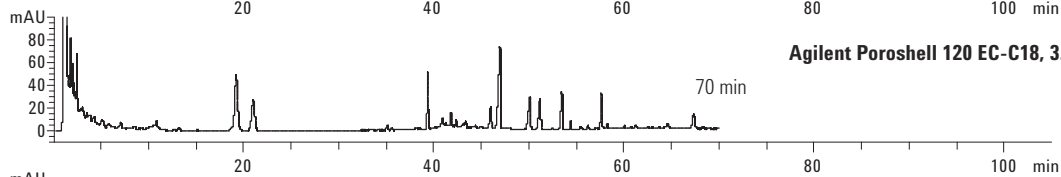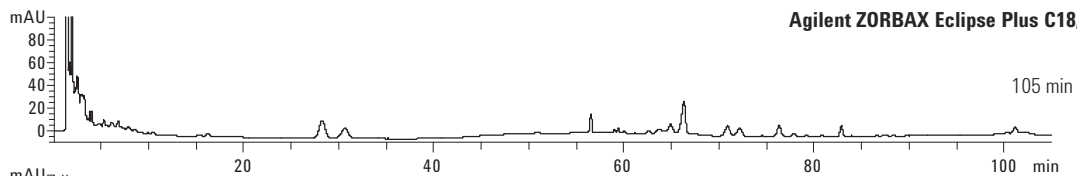

NEW!

Separation of 8 steroids

Column A: Poroshell 120 EC-C18
695775-902
2.1 x 100 mm, 2.7 µm

Column B: Poroshell 120 SB-C18
685775-902
2.1 x 100 mm, 2.7 µm

Column C: Poroshell 120 Phenyl-Hexyl
695775-912
2.1 x 100 mm, 2.7 µm


Column D: Poroshell 120 Bonus RP
685775-901
2.1 x 100 mm, 2.7 µm

Mobile Phase: 0.1% formic acid
in both water and MeOH

Flow Rate: 0.4 mL/min, 25 °C,
2.1 x 100 mm 40 °C

Gradient: 40-80% MeOH in 14 min

1. Hydrocortisone
2. β-Estradiol
3. Androstadiene 3,17 dione
4. Testosterone
5. Ethynylestradiol
6. Estrone
7. Norethindrone acetate
8. Progesterone

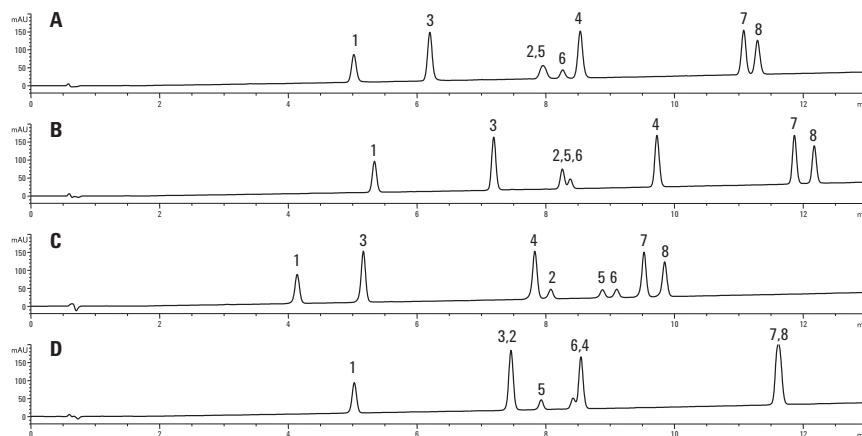

NEW!

Mixture of beta blockers

Column A: Poroshell 120 Bonus RP
685775-901
2.1 x 100 mm, 2.7 µm

Column B: Poroshell 120 Phenyl-Hexyl
695775-912
2.1 x 100 mm, 2.7 µm

Column C: Poroshell 120 EC-C18
695775-902
2.1 x 100 mm, 2.7 µm

Column D: Poroshell 120 SB-C18
685775-902
2.1 x 100 mm, 2.7 µm


Mobile Phase: 10 mM pH 3.8 NH₄HCO₂, methanol

Flow Rate: 0.35 mL/min

Gradient: 90% B to 30% B over 12 min

* Nadolol is isobaric and elutes as two peaks.

1. Atenolol
2. Pindolol
3. Nadolol
4. Metoprolol
5. Acebutolol
6. Propranolol
7. Alprenolol

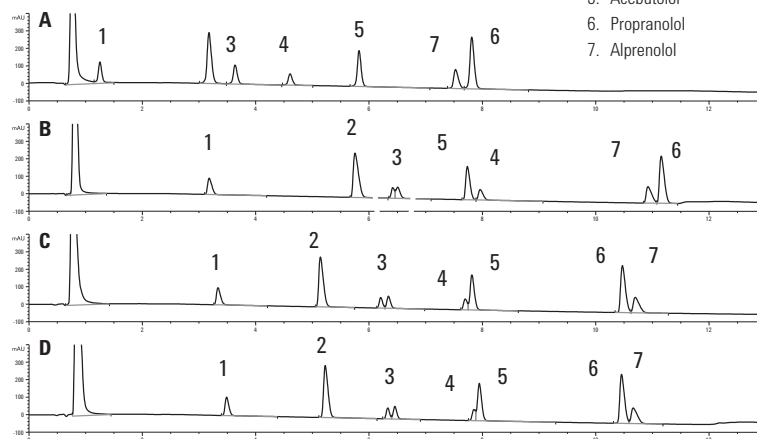

NEW!

Several ZORBAX RRHD 1.8 μ m selectivities facilitate method development

Column: ZORBAX RRHD Eclipse Plus C18
959758-902
2.1 x 100 mm, 1.8 μ m

Column: ZORBAX RRHD Eclipse XDB-C18
981758-902
2.1 x 100 mm, 1.8 μ m

Column: ZORBAX RRHD SB-C18
858700-902
2.1 x 100 mm, 1.8 μ m

Column: ZORBAX RRHD Extend-C18
758700-902
2.1 x 100 mm, 1.8 μ m

Mobile Phase: A: H₂O
B: CH₃CN, each with 0.1% HCOOH


Detector: Agilent 1290 Infinity LC with an
Agilent 6410 Triple Quadrupole Mass Spectrometer

MS Conditions: TCC: 30 °C
MS Source: Electrospray AP-ESI
Drying-gas temperature and flow: 325 °C, 12 L/min
Nebulizer gas pressure: 35 psi
Capillary voltage: 3000 V

Sample: Four endocannabinoid fatty amides:
Arachidonoylglycerol (AEA)
2-Arachidonoylglycerol (2-AG)
Palmitoylethanolamide (PEA)
Oleoylethanolamide (OEA)

* The second black peak is an impurity, believed to be
1,3-arachidonoylglycerol, a rearrangement of 2-AG

The selectivity of four Agilent ZORBAX RRHD C18 columns is compared
using a method for endocannabinoids.

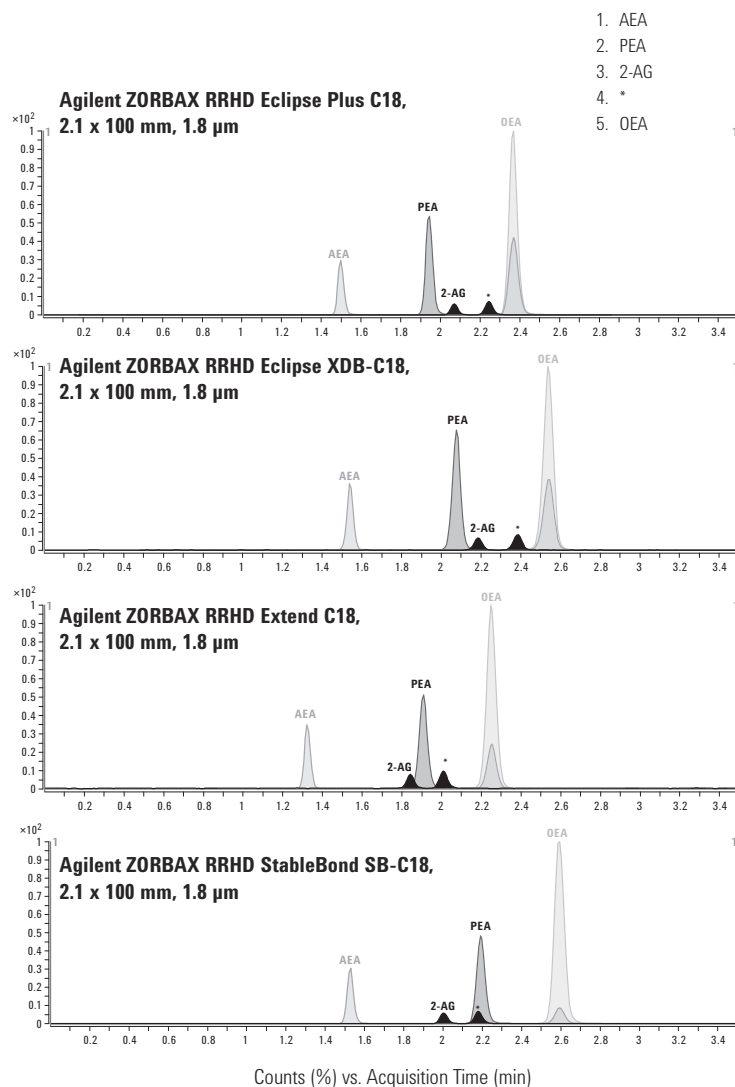

For a comprehensive listing of chromatograms searchable by compound name, visit our online
Chromatogram Library at www.agilent.com/chem/library

Fast analysis 11 common compounds found in analgesics


Column: Poroshell 120 EC-C18
695975-902
4.6 x 100 mm, 2.7 µm

Mobile Phase: A : Water + 0.1% formic acid
B: ACN

Flow Rate: 3.5 mL/min

Temperature: 40 °C

Detector: DAD 254 nm

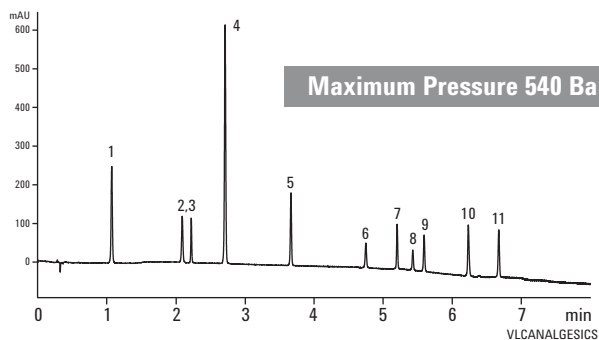

1. Acetaminophen
2. Caffeine
3. 2-Acetamidophenol
4. Acetamide
5. Phenacetin
6. Sulindac
7. Piroxicam
8. Tolmetin
9. Ketoprofen
10. Diflusal
11. Diclofenac

Faster analysis of USP Method for simvastatin tablet

Column A: Eclipse Plus C18
959990-902
4.6 x 250 mm, 5 µm


Column B: Poroshell 120 EC-C18
697975-902
4.6 x 75 mm, 2.7 µm

Mobile Phase: 65% CH₃CN,
35% 3.9 g/L NaH₂PO₄ (pH 4.5)

Flow Rate: 1.5 mL/min for 5 µm column
2.8 mL/min for 2.7 µm Poroshell 120 column

Temperature: 45 °C

Detector: DAD Sig = 238, 8
Ref = 360, 100 nm

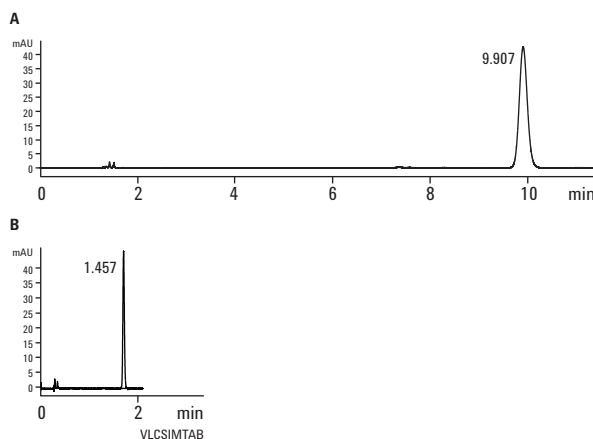

	USP Requirement	5 µm (1.5 mL/min)	2.7 µm (2.8 mL/min)
T_R	N/A	9.907	1.457
k'	> 3.0	5.962	5.122
N	> 4500	16939	14439
T_f	< 2.0	1.09	1.10

Faster separation of sulfa drugs

Column A: Eclipse Plus C18
959990-902
4.6 x 250 mm, 5 µm


Time	%B
0	8
33	33
35	33

Column B: Poroshell 120 EC-C18
695975-902
4.6 x 100 mm, 2.7 µm

Time	%B
0	8
12	33
13.2	33

Mobile Phase: A: 0.1% formic acid in Water
B: 0.1% formic acid in ACN

Flow Rate: 1 mL/min

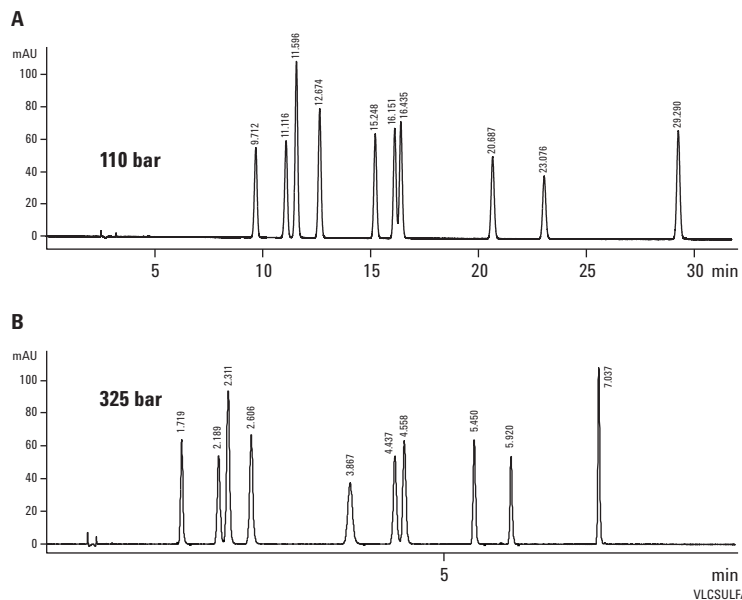

Separation of pharmaceutical cardiac drugs


Column: Eclipse Plus C18
959996-902
4.6 x 100 mm, 5 µm

Mobile Phase: A: 20 mM Ammonium Acetate, pH 4.8
B: ACN

Flow Rate: 1 mL/min

Gradient: 10-90% in 10 min

Detector: UV, 254 nm

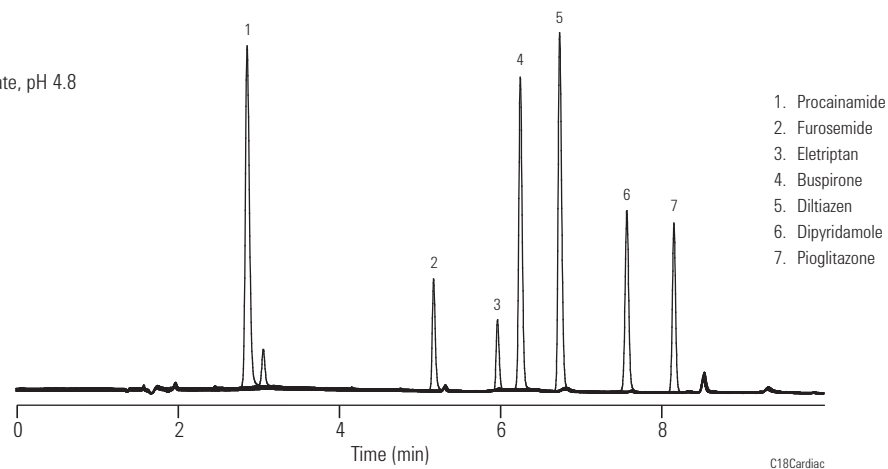

Fast and ultra-fast analysis of basic compounds

Column: Eclipse Plus C18
959941-902
4.6 x 50 mm, 1.8 μ m


Mobile Phase: A: 50% 8 mM K_2HPO_4 , pH 7
B: 50% ACN

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Berberine, 0.4 mg/mL, 2 μ L

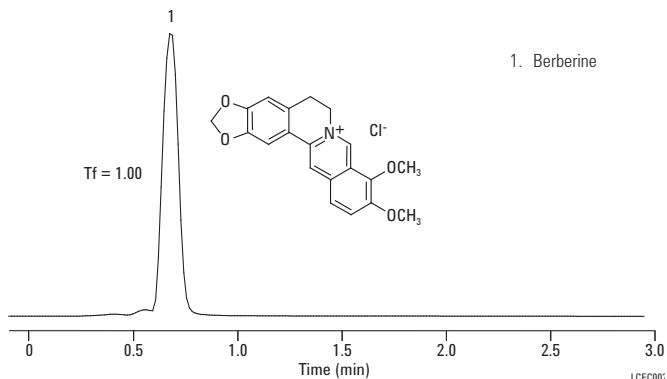

Xanthines: Higher resolution, same selectivity with RRHT

Column A: ZORBAX SB-C18
846975-902
4.6 x 50 mm, 5 μ m

Column B: ZORBAX SB-C18
827975-902
4.6 x 50 mm, 1.8 μ m


Mobile Phase: A: 92% 0.1% formic acid
B: 8% 0.1% formic acid in ACN

Flow Rate: 1.5 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Xanthines

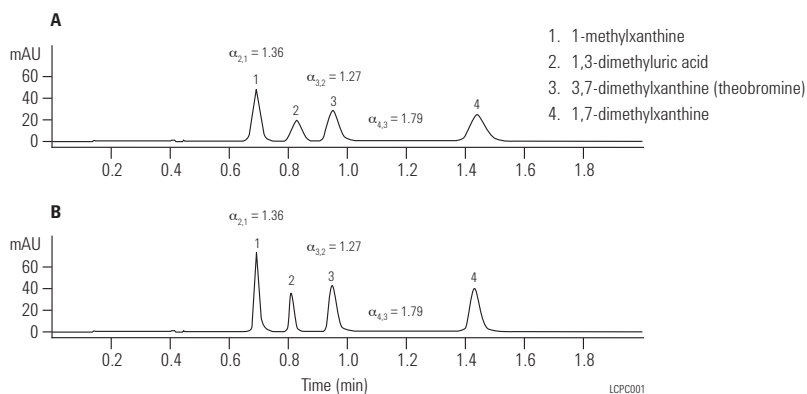

Antihistamines: Fast separations on RRHT Extend-C18

Column A: ZORBAX Extend-C18
773450-902
4.6 x 150 mm, 5 μ m

Column B: ZORBAX Extend-C18
727975-902
4.6 x 50 mm, 1.8 μ m


Mobile Phase: A: 30% 50 mM pyrrolidine buffer
B: 70% MeOH

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 220 nm

Sample: Antihistamines

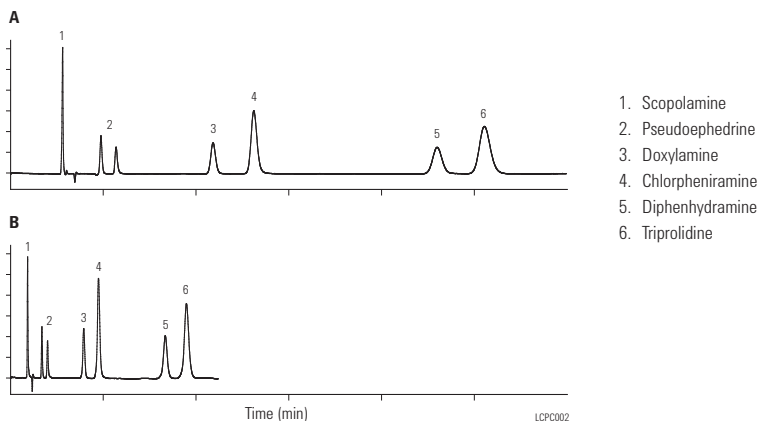

Ibuprofen: Optimizing selectivity with RRHT Columns

Column A: SB-C8
827975-906
4.6 x 50 mm, 1.8 μ m

Column B: Eclipse XDB-C8
927975-906
4.6 x 50 mm, 1.8 μ m


Mobile Phase: A: 63% water
B: 37% acetonitrile + 1.8 mL H_3PO_4

Flow Rate: 2.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Ibuprofen oral suspension

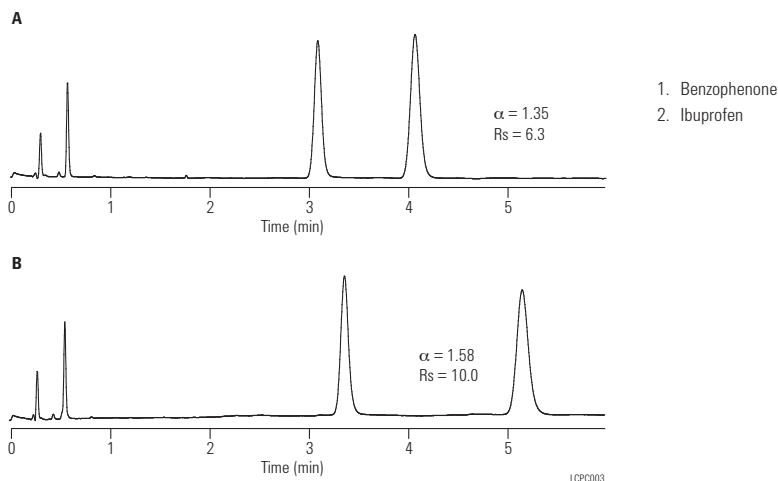

Analgesics

Column: Pursuit XR_s Diphenyl
A6020150X046
4.6 x 150 mm, 5 μ m


Mobile Phase: A: Water+0.1% HCOOH
B: MeCN+0.1% HCOOH

Gradient: 25-80% B in 20 min

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

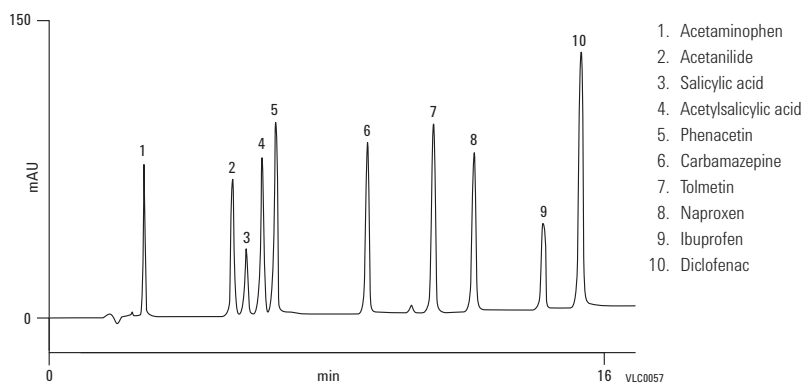

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Anesthetics, local: Bonded phase selectivity

Column A: ZORBAX SB-C18
883975-902
4.6 x 150 mm, 5 µm

Column B: ZORBAX SB-C8
883975-906
4.6 x 150 mm, 5 µm

Column C: ZORBAX SB-C3
883975-909
4.6 x 150 mm, 5 µm

Column D: ZORBAX SB-Phenyl
883975-912
4.6 x 150 mm, 5 µm

Column E: ZORBAX SB-CN
883975-905
4.6 x 150 mm, 5 µm

Mobile Phase: A: 50 mM NaH₂PO₄ pH 2.5 in 95% H₂O/5% ACN
B: 50 mM NaH₂PO₄ pH 2.5 in 47% H₂O/53% ACN


Flow Rate: 1.5 mL/min

Gradient: 0-100% B in 18.8 min





Temperature: 26 °C

Detector: UV, 254 nm

Sample: 10 µL, 10 µg/mL

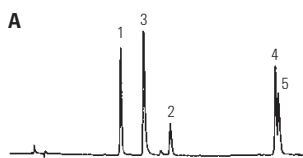

1. Procaine
2. Lidocaine
3. d-Cinchonine
4. Butacaine
5. Tetracaine

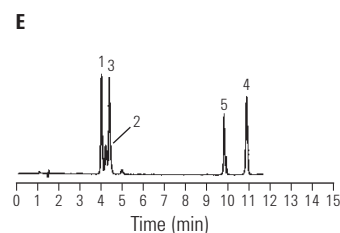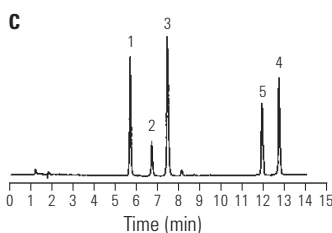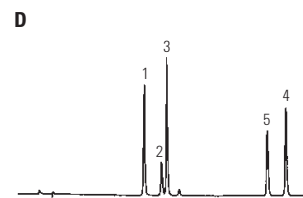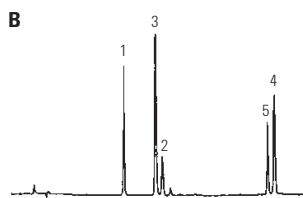

LCPC005

Local anesthetics


Column: Pursuit XRs C8
A6010150X046
4.6 x 150 mm, 5 µm

Mobile Phase: 65:35 MeOH:5 mM NH₄CO₃, pH 10

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 210 nm

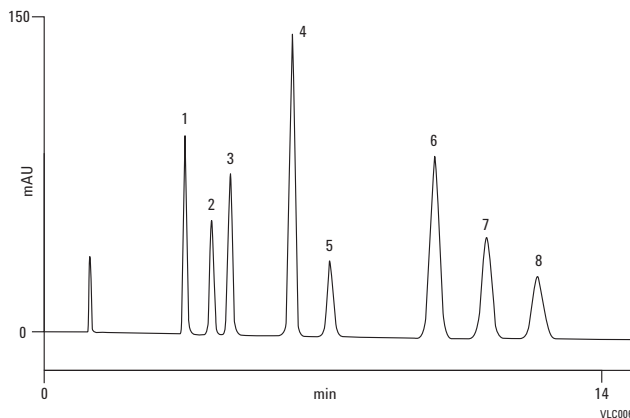

1. Benzocaine
2. Procaine
3. Chlorocaine
4. Mepivacaine
5. 4-Hydroxypropivacaine
6. Cocaine
7. Lidocaine
8. Ropivacaine

VLC0063

Antibiotics: High speed separation

Column: ZORBAX Rx/SB-C8
866953-906
4.6 x 75 mm, 3.5 μ m

Mobile Phase: 8.0% acetonitrile/92% 0.1% aqueous TFA


Flow Rate: 3.0 mL/min

Gradient: 45-70% B in 35 min

Temperature: 60 °C

Detector: UV, 260 nm

Sample: 1 μ L containing 0.40, 0.36, 0.10 and 0.37 μ g each of 1-4 resp.

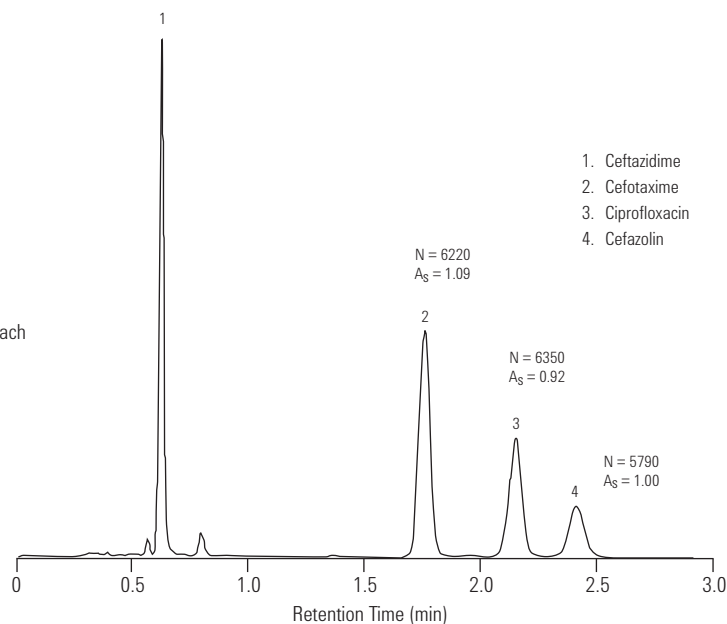

LCPC007

Antibiotics: Lincomycin and Clindamycin by LC-APCI-MS LC-TIC

Column: ZORBAX SB-C18 cartridge
823700-902
2.1 x 30 mm, 1.8 μ m

Mobile Phase: Gradient: 15-50% B in 1 min, hold for 1.5 min,
A: 0.2% formic acid pH 2.8
B: ACN + 0.2% formic acid

Flow Rate: 0.5 mL/min


Gradient: Post time: 1.5 min

Temperature: Ambient

Detector: APCI, Positive ion

MS Conditions: Peak width: 0.10 min
Scan: 150-600 Da, step 0.1
Fragmentor: 70
Gas Temp: 350 °C
Vaporizer: 350 °C
Drying gas: 12 L/min
Nebulizer pres: 50 psi
Vcap: +3000 V
Corona: 4.0 μ A

Sample: Antibiotics, 1 μ L

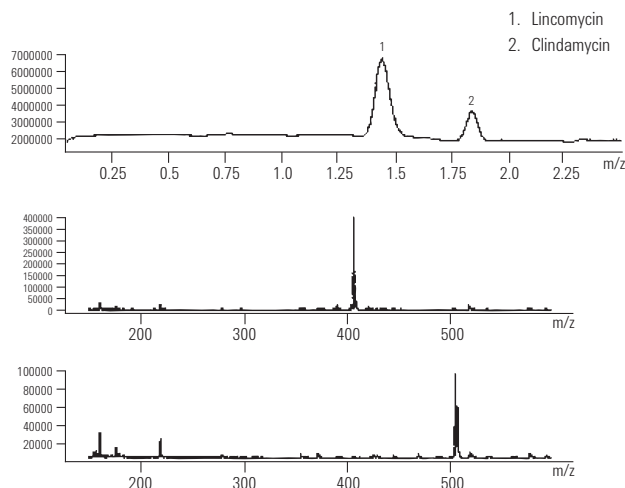

LCPC008

Antifungal medications

Column: ZORBAX Bonus-RP
883668-901
4.6 x 150 mm, 5 µm


Mobile Phase: 35% 25 mM NaH₂PO₄, Dibasic (pH 6.5 with H₃PO₄);
65% ACN

Flow Rate: 1 mL/min

Temperature: Ambient

Detector: UV, 220 nm

Sample: Antifungals, 2 µL

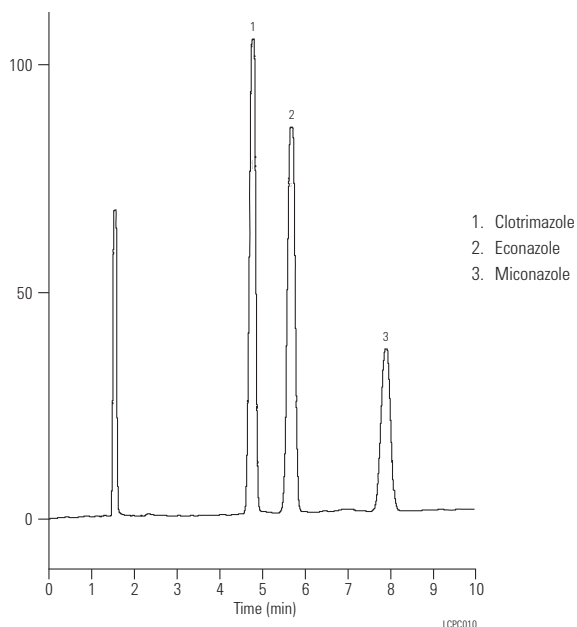

Antifungals


Column: Pursuit XRs Diphenyl
A6020150X046
4.6 x 150 mm, 5 µm

Mobile Phase: Water+0.1% HCOOH:
MeCN+0.1% HCOOH, 80:20

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

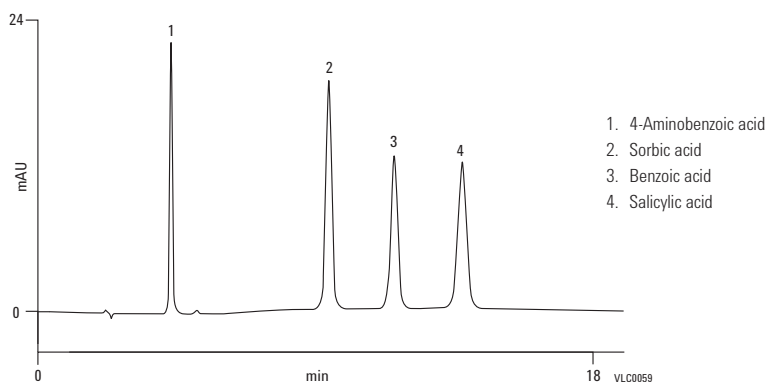

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

**Analgesics: Non-steroidal anti-inflammatory drugs:
Narrow bore separation**

Column: Eclipse XDB-C8
993700-906
2.1 x 150 mm, 5 µm


Mobile Phase: 50/50, 25 mM Sodium Phosphate
(pH 7.0 with Phosphoric Acid), MeOH

Flow Rate: 0.2 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

Sample: 2 µL, 10 µg/mL

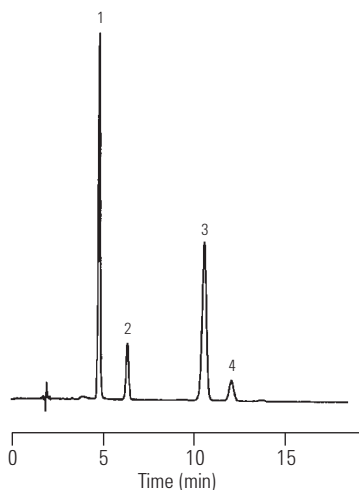

NSAID	pK _a
1. Phenacetin	2.2
2. Tolmetin	3.5
3. Phenylbutazone	4.4
4. Fenoprofen	4.5

LCPC011

Separation of small molecule anorectics

Column A: ZORBAX Bonus-RP
883668-901
4.6 x 150 mm, 5 µm

Column B: Traditional Alkyl C8 Phase


Mobile Phase: 25 mM K₂HPO₄, pH 7.2/MeOH: ACN (50:50), 45/55

Flow Rate: 1 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Anorectics "Fen-phen", 5 µL

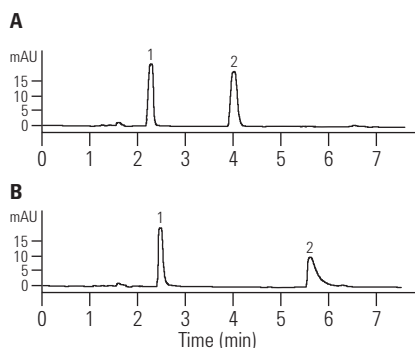

1. Phentermine
2. Fenfluramine

LCBP004



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Aromatic acids/benzoic acids: Selectivity differences

Column A: ZORBAX SB-C8
880975-906
4.6 x 250 mm, 5 µm

Column B: ZORBAX SB-Phenyl
880975-912
4.6 x 250 mm, 5 µm

Column C: ZORBAX SB-CN
880975-905
4.6 x 250 mm, 5 µm

Mobile Phase: 30-45% methanol in 25 mM Na Phosphate, pH 2.5

A: 45% Methanol

B: 40% Methanol


C: 30% Methanol

Flow Rate: 1.0 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

Sample: Benzoic acids

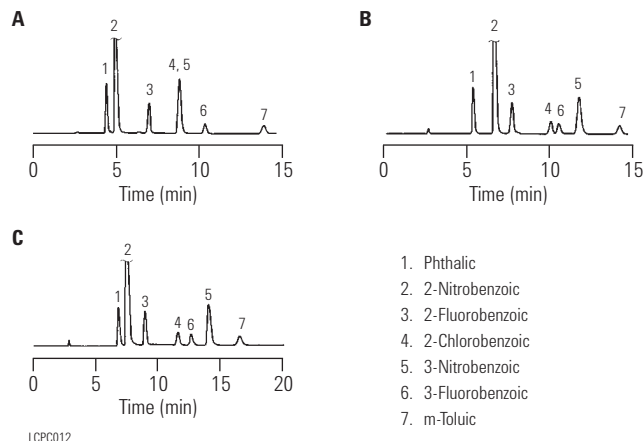

Catecholamines/biogenic amines: Rapid separation using ion-pair reagents

Column: ZORBAX Rx/SB-C8
866953-906
4.6 x 75 mm, 3.5 µm

Mobile Phase: 0.14 M sodium phosphate,
20 mM EDTA,
0.75 mM octyl sulfonate,
9% methanol pH 3.5


Flow Rate: 1.5 mL/min

Temperature: 26 °C

Detector: 0.75 V vs Ag/AgCl with electro-chemical detection

Sample: 10 µg/mL each standard; volume
20 µL (2 g tissue sample)
A: Standards (2pmol; DHBA 5pmol)
B: Mouse Satrium
C: Mouse Neocortex

- | | |
|--------------------------------------|----------------------------------|
| 1. DOPA-Dihydroxyphenylalanine | 6. HIAA-Hydroxyindoleacetic acid |
| 2. DHBA-Dihydroxybenzyl amine | 7. EP-Epinephrine |
| 3. DOPAC-Dihydroxyphenyl acetic acid | 8. HVA-Homovanillic acid |
| 4. NE-Norepinephrine | 9. 5-HT-Hydroxytryptamine |
| 5. DA-Dopamine | 10. 3-MT-Methoxytyrosine |

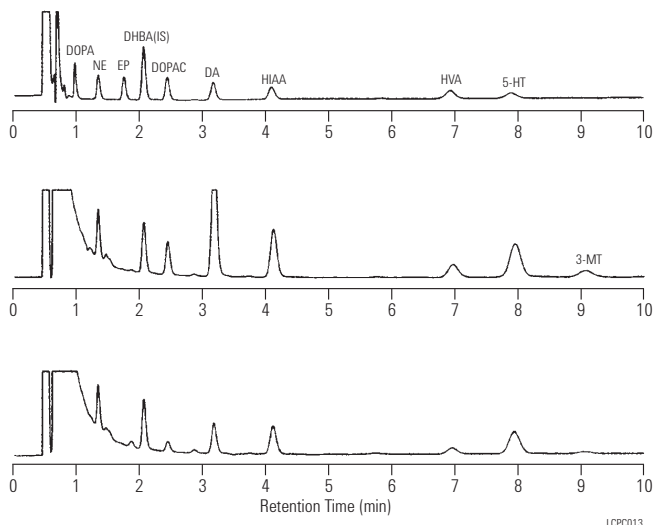

Chiral ethiazide (diuretic drug) separation

Column: Ultron ES-OVM Chiral
702111651
4.6 x 150 mm, 5 µm


Mobile Phase: 20 mM KH₂PO₄ (pH 4.6)

Flow Rate: 1.0 mL/min

Temperature: 25 °C

Detector: UV, 220 nm

Sample: 20 µL containing 0.35 µg Ethiazide

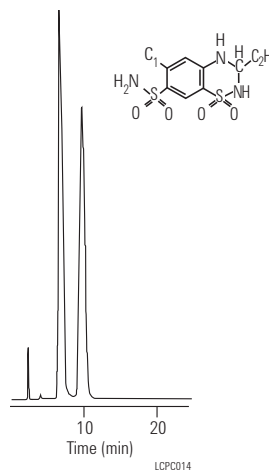

Chiral separation of fluoxetine enantiomers (Prozac)

Column: Ultron ES-OVM Chiral
702111651
4.6 x 150 mm, 5 µm


Mobile Phase: 25/75 (v/v) EtOH / 20 mM KH₂PO₄, pH 5.5
(adjusted with NaOH)

Flow Rate: 0.8 mL/min

Temperature: Ambient

Detector: UV, 225 nm

Sample: Mixture fluoxetine (Prozac) enantiomers

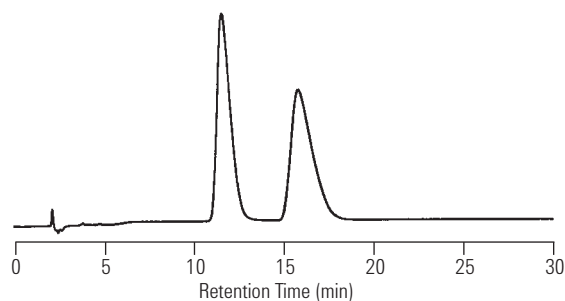

Courtesy of D.S. Ristry and V.S. Sharp, Eli Lilly and Co.



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Goldenseal and related alkaloids on Rapid Resolution Eclipse XDB-C18

Column: Eclipse XDB-C18
963967-902
4.6 x 150 mm, 3.5 µm

Mobile Phase: 68% 30 mM ammonium acetate,
14 mM TEA, pH ~4.85
32% Acetonitrile


Flow Rate: 1.0 mL/min

Temperature: 30 °C

Detector: 230 nm

Sample: Goldenseal and related alkaloids

Alkaloids, such as the active components in Goldenseal and other related plants, are quickly and accurately separated using isocratic conditions on an Eclipse XDB-C18 Rapid Resolution column.

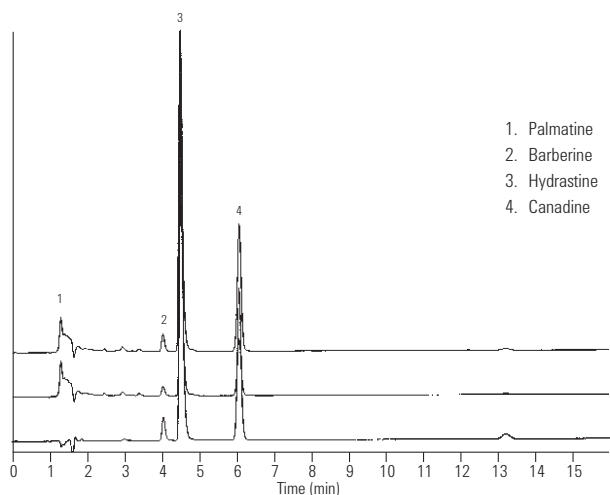

Components of green tea separated on Rapid Resolution StableBond SB-C8

Column: ZORBAX SB-C8
863953-906
4.6 x 150 mm, 3.5 µm

Mobile Phase: 75% 0.1% TFA : 25% MeOH


Flow Rate: 1.0 mL/min

Temperature: 40 °C

Detector: 280 nm

Sample: Green tea

Nutraceuticals, such as the components of green tea, are quickly separated on a StableBond SB-C8 Rapid Resolution column.

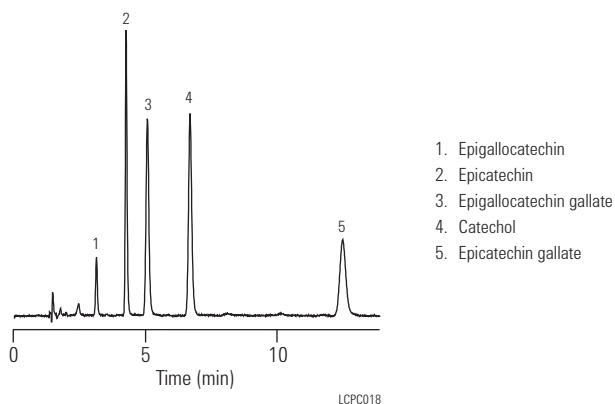

Chiral separation of hexobarbital


Column: Chiradex
79925CB-584
4.0 x 250 mm, 5 µm

Mobile Phase: Methanol/water, 20:80

Flow Rate: 1.0 mL/min

Detector: UV, 220 nm

Sample: Hexobarbital

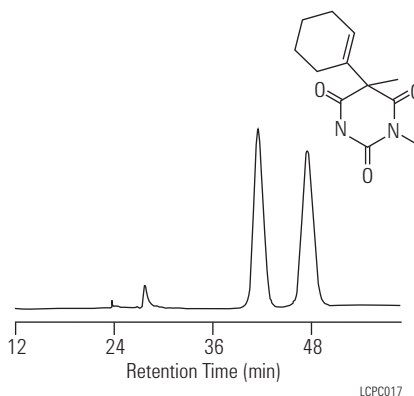

Chiral separation of S- and R-Norfluoxetine

Column: Ultron ES-OVM Chiral
724111653
4.6 x 250 mm, 10 µm


Mobile Phase: 6/94 (v/v) MeOH / 20 mM KH₂PO₄

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 225 nm

Sample: 50 µg/mL of 2:3 mixture R- : S-Norfluoxetine

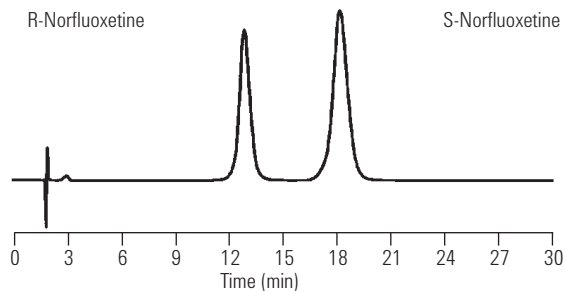

Courtesy of D.S. Ristry and V.S. Sharp, Eli Lilly and Co.

LCPC019

Chiral separation of salbutamol

Column: Ultron ES-Pepsin
822111631A
4.6 x 150 mm, 5 µm


Mobile Phase: 20 mM phosphate buffer, pH 6.0

Flow Rate: 1.0 mL/min

Temperature: 25 °C

Detector: UV, 220 nm

Sample: 20 µL containing 0.35 µg salbutamol mixture

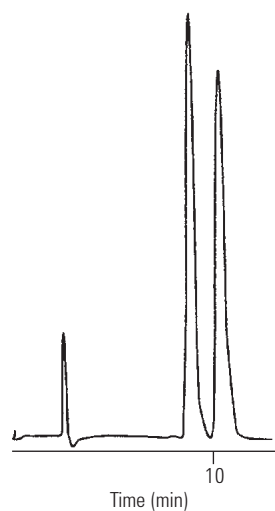

LCPC020



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Chiral separation of tolperison enantiomers

Column: Ultron ES-OVM Chiral
702111651
4.6 x 150 mm, 5 µm


Mobile Phase: 20 mM KH₂PO₄ (pH 5.5), C₂H₅OH (100/4 v/v)

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 220 nm, 0.04 AUFS

Sample: Tolperison, 5 µL

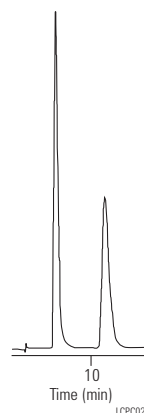

LCPC022

Chiral separation of atenolol

Column: Ultron ES-Pepsin
822111631A
4.6 x 150 mm, 5 µm


Mobile Phase: 20 mM phosphate buffer, pH 6.0/Ethanol (99/1)

Flow Rate: 1.0 mL/min

Temperature: 25 °C

Detector: UV, 220 nm, 0.04 AUFS

Sample: 1.5 µL, 0.25 mg/mL, atenolol racemic mixture

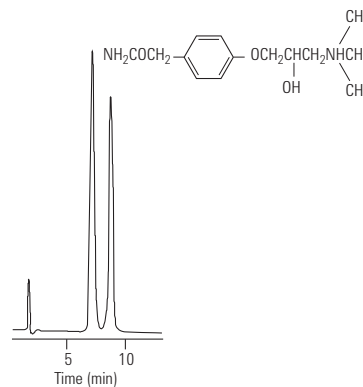

LCPC021

Cocaine and metabolites


Column: ZORBAX Rx-SIL
883975-901
4.6 x 150 mm, 5 µm

Mobile Phase: MeOH: NH₄ Acetate, 25 mM, pH 6 (70:30)

Flow Rate: 1.0 mL/min

Temperature: 40 and 80 °C

Detector: UV, 210 nm

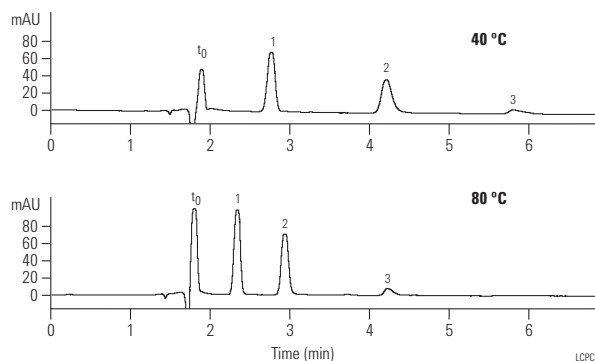

LCPC023

Aspirin and cough remedy

Column: Eclipse XDB-C8
993967-906
4.6 x 150 mm, 5 µm


Mobile Phase: (75:25) 25 mM Na₂HPO₄ (pH 3.0): ACN

Flow Rate: 1.0 mL/min

Temperature: 40 °C

Detector: UV, 254 nm

Sample: 5 µL, 10 µg/mL

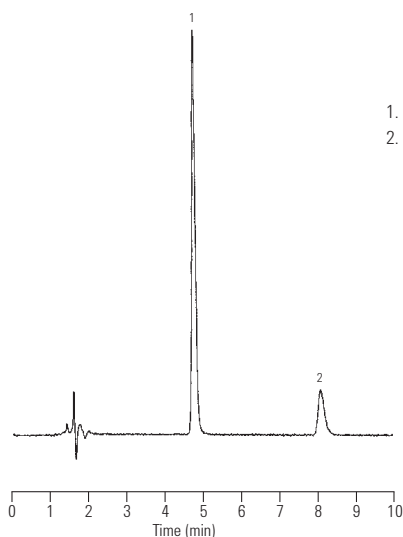

1. Acetylsalicylic acid
2. Dextromethorphan

LCPC024

Cough formula mixture: Fast and efficient separation

Column A: ZORBAX SB-CN
866953-905
4.6 x 75 mm, 3.5 µm

Column B: ZORBAX SB-CN
883975-905
4.6 x 150 mm, 5 µm


Mobile Phase: 20/80, Acetonitrile/150 mM Na Citrate, pH 2.6

Flow Rate: 1.5 mL/min, 1.0 mL/min


Temperature: 35 °C

Detector: UV, 270 nm

Sample: 2 µL, cough formula

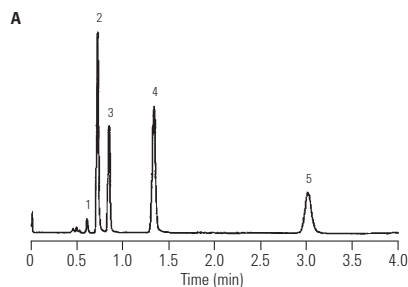

1. Maleic acid
2. Pseudo-ephedrine
3. Acetaminophen
4. Chlorpheniramine
5. Dextromethorphan

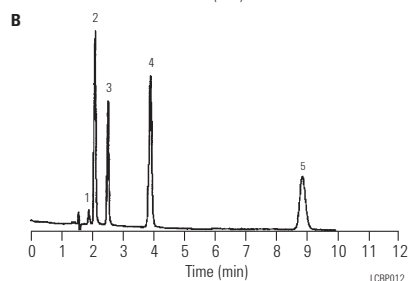

LCBP012

Guaifenesin: USP analysis of guaifenesin

Mobile Phase: 40% Methanol:60% Water:1.5% Glacial Acetic Acid

Flow Rate: 1.0 mL/min

Temperature: 25 °C


Sample: Guaifenesin

A: 8 µL


B: 2 mL

Column: **Eclipse XDB-C18**
990967-902
4.6 x 250 mm, 5 µm

Peak	TR	N	Rs
1	6.63	12,737	0
2	11.19	18,552	15.8

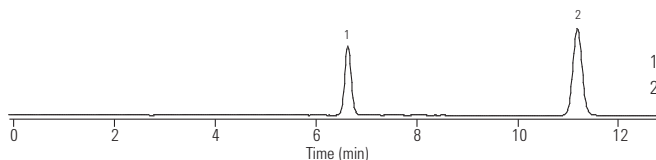

1. Guaifenesin: 0.04 mg/mL
2. Benzoic Acid: 0.10 mg/mL

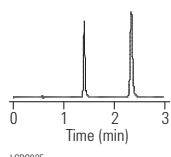

Column: **Eclipse XDB-C18**
922975-902
4.6 x 50 mm, 1.8 µm

Peak	TR	N	Rs
1	1.4	11,421	0
2	2.33	12,909	12.3

Minimum Resolution Required = 3.0

Metronidazole: Updating USP methods

Column A: **ZORBAX C8**
883952-706
4.6 x 150 mm, 5 µm

Column B: **Eclipse XDB-C8**
993967-906
4.6 x 150 mm, 5 µm

Column C: **Eclipse XDB-C8**
963967-906
4.6 x 150 mm, 3.5 µm


Mobile Phase: 80/20, Water/Methanol

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm


Sample: Metronidazole

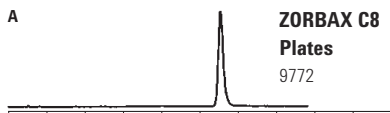

ZORBAX C8
Plates
9772

USP TF
1.30


Particle Size
5 µm

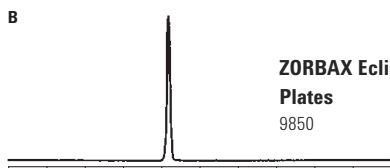

ZORBAX Eclipse XDB-C8
Plates
9850

USP TF
0.98

Particle Size
5 µm

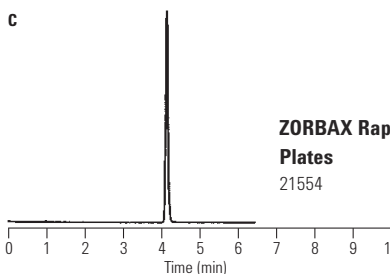

ZORBAX Rapid Resolution Eclipse XDB-C8
Plates
21554

USP TF
1.13

Particle Size
3.5 µm

Morphine and metabolites: Extracted blood plasma sample separation

Column: ZORBAX SB-C18
863953-902
4.6 x 150 mm, 3.5 µm


Mobile Phase: 97/3 70 mM KH₂PO₄ + 1 mM EDTA/ACN, pH 4.5

Flow Rate: 1.5 mL/min

Temperature: Ambient

Detector: A: Electrochemical, 720 mV
B: Fluorescence, Ex = 285 nm, Em = 352 nm

Sample: 50 µL
Morphine-3-glucuronide 125 ng
Morphine-6-glucuronide 62.5 ng
Morphine 12.5 ng in Water

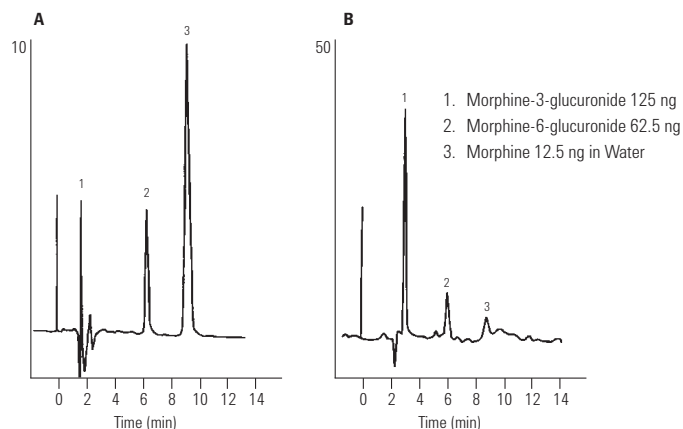

Courtesy of J. Visser, Center for Pharmacy, Univ. Groningen, The Netherlands.

LCPC027

Opiates (drugs of abuse) by LC/MS

Column: ZORBAX SB-AQ
830990-914
2.1 x 150 mm, 3.5 µm

Mobile Phase: A: Acetonitrile with 0.1% formic acid
B: Water with 0.1% formic acid


Flow Rate: 0.25 mL/min

Gradient: 0 min 10% B
5 min 35% B
5.1 min 100% B

MS Conditions: Time of Flight (TOF)
Standard with calibrant delivery system
providing constant low flow of ~ 2 µM purine
and HP-921 calibrant to dual ESI for
continuous auto-calibration

Sample: Opiates

XIC of +TOF MS

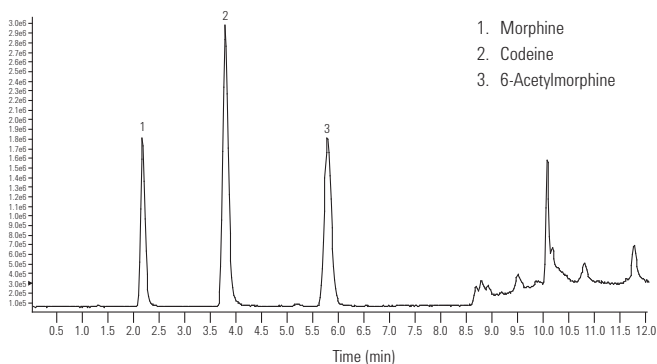

LCPC028



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

NEW!**Comparing HILIC and RPLC of morphine using Agilent ZORBAX RRHD columns with UHPLC/MS**

Column: Agilent ZORBAX Eclipse Plus C18
2.1 x 100mm, 5 μ m
(Custom column)

Column: ZORBAX RRHD HILIC Plus
959758-901
2.1 x 100 mm, 1.8 μ m

Mobile Phase: A: 10 mM NH_4HCO_2 , pH 3.2
B: $\text{CH}_3\text{CN}/100$ mM NH_4HCO_2 , pH 3.2 (9:1)
Column A: 10% B isocratic
Column B: 70% B isocratic

Flow Rate: Column A: 0.4 mL/min
Column B: 1 mL/min

Pressure: Column A: 90 bar
Column B: 810 bar


Temperature: 25 $^\circ\text{C}$

Detector: Agilent 1290 Infinity LC with an
Agilent 6410A Triple Quadrupole Mass Spectrometer

MS Conditions: MS Source: Positive ESI, capillary 4000 V, drying gas temperature, flow rate and nebulizer pressure vary with mobile phase flow rate
MS Acquisition: Selected ion mode (SIM), delta EMV 200 V, MS dwell time varies with mobile phase flow rate
Software: Agilent MassHunter versions B.03.01, B.02.00 AND B.03.01 were used for data acquisition, qualitative, and quantitative analyses, respectively

Sample: 2 μL injection of 1 $\mu\text{g}/\text{mL}$ each of morphine, normorphine, morphine-3- β -D-glucuronide: HILIC sample was prepared in CH_3CN ; RPLC sample was prepared in H_2O

HILIC mode with UHPLC columns cuts analysis time in half, while improving sensitivity by more than a factor of 10, compared to traditional LC columns in RPLC mode with MS detection.

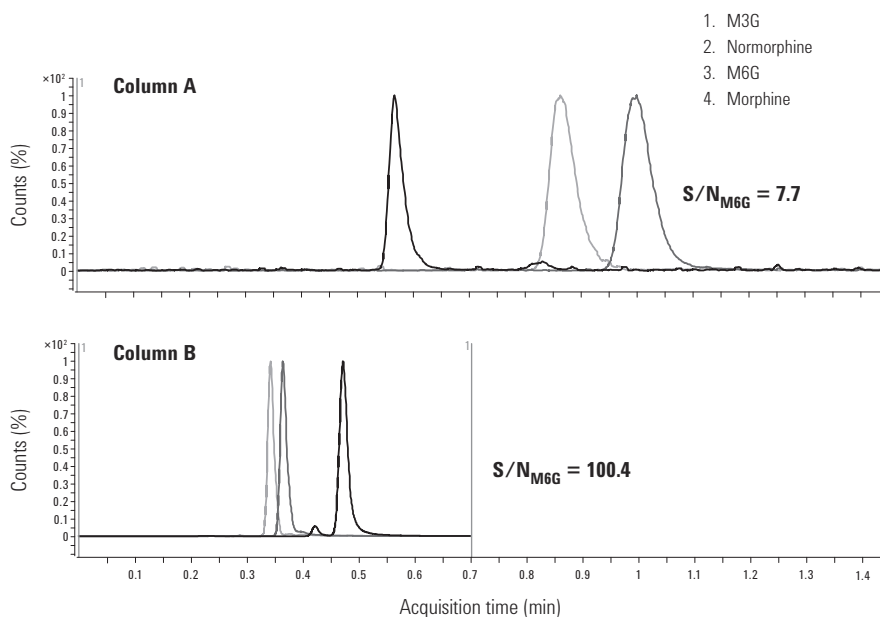

Neutraceuticals: Hypericin separation in St. John's Wort

Column: Eclipse XDB-C8
993967-906
4.6 x 150 mm, 5 µm


Mobile Phase: 23% 25 mM Na₂HPO₄,
Dibasic (pH 7.0 with H₃PO₄): 77% MeOH

Flow Rate: 1.0 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

Sample: Neutraceuticals

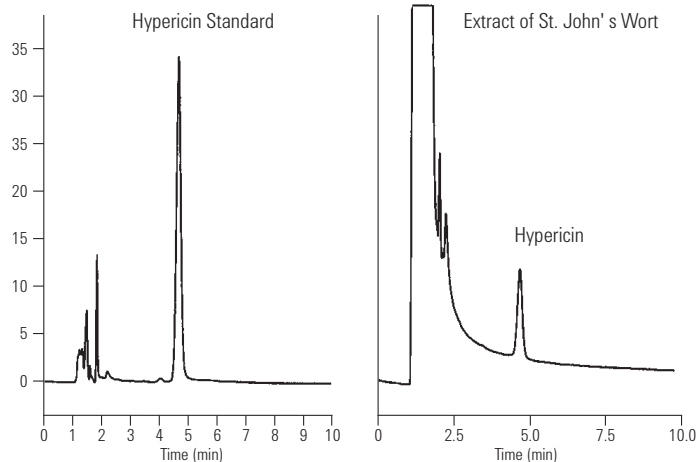

LCPC029

Pharmaceuticals: Rapid, high sensitivity LC and LC/MS of 11 drugs

Column: Eclipse XDB-C18
925700-902
2.1 x 50 mm, 1.8 µm

Mobile Phase: A: 10 mM NH₄ Formate (pH = 3.6)
B: ACN with 10 mM NH₄ Formate


Flow Rate: 0.6 mL/min

Gradient: 5% B to 70% B in 7.5 min, to 95% B in 8.5 min

Temperature: 65 °C

Detector: UV, 230 nm and MSD Trap SL

MS Conditions: Pos. Dry Gas: 345 °C
Neb.: 45 psi
HV Cap: 3500 V
Range: 100-700
Average: 5 Spectra
ICC: 30000
Charge Con: On
Smart Par. Settings: Tar Mas: 250 m/z
Comp. Stab.: 100%
Trap Drive: 100%
Frag. Options: Smart Frag: On
Frag. Width: 10 m/z

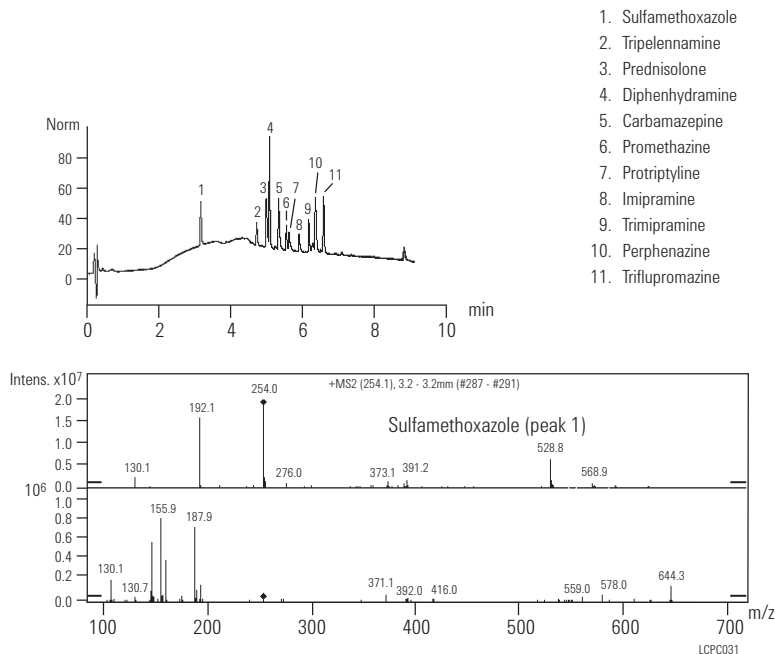

Hormones/steroids

Column: ZORBAX RRHT SB-C18
823975-902
4.6 x 30 mm, 1.8 µm


Mobile Phase: 50% 20 mM NaH₂PO₄, pH 2.8: 50% ACN

Flow Rate: 1.0 mL/min

Temperature: RT

Detector: UV, 230 nm

Sample: Hormones/steroids

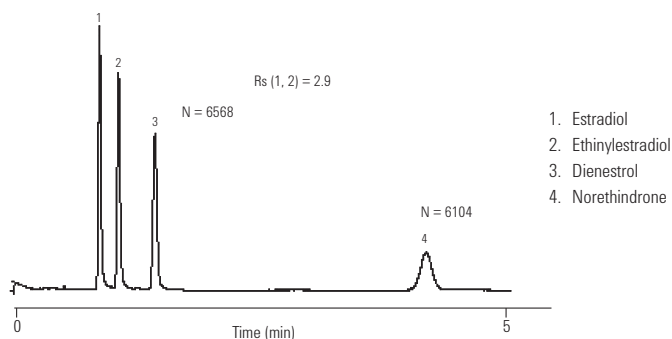

LCPC034

Steroids: Separation

Column: Eclipse XDB-CN
993967-905
4.6 x 150 mm, 5 µm


Mobile Phase: 40:60 ACN:Water

Flow Rate: 1.0 mL/min

Temperature: 25 °C

Detector: UV, 205 nm

Sample: 1. Norethindrone 0.514 mg/mL
2. Progesterone 0.407 mg/mL
3. Mestranol 0.057 mg/mL

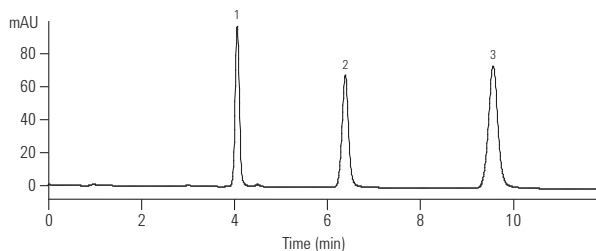

LCPC036



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Steroids

Column A: Eclipse XDB-Phenyl
963967-912
4.6 x 150 mm, 3.5 µm

Column B: Eclipse XDB-C18
993967-902
4.6 x 150 mm, 5 µm

Mobile Phase: H₂O:ACN, 60:40



Flow Rate: 1.0 mL/min

Temperature: 35 °C

Detector: UV, 254 nm

Sample:

1. Prednisolone
2. Corticosterone
3. 11 -hydroxyprogesterone
4. Cortisone acetate
5. Deoxycorticosterone
6. 17 hydroxyprogesterone
7. Progesterone

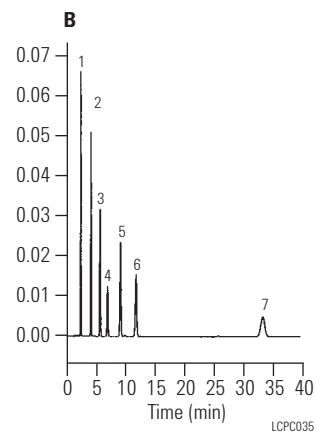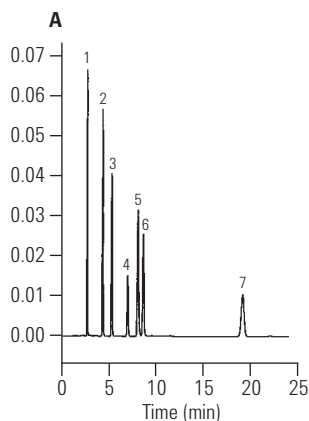

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Triamcinolone – USP analysis of triamcinolone


Column: Eclipse XDB-C18
923975-902
4.6 x 30 mm, 1.8 µm

Mobile Phase: 47% Methanol:53% Water


Flow Rate: 1.5 mL/min

Temperature: 25 °C

Sample: Triamcinolone, 1 µL

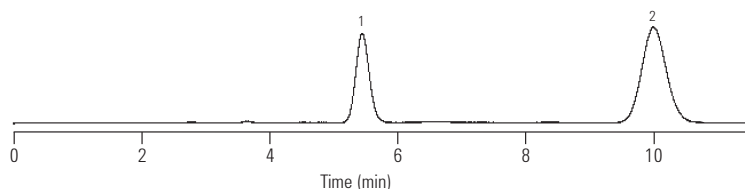

Peak	TR	N	Rs
1	5.45	3199	0
2	9.99	3212	8.1

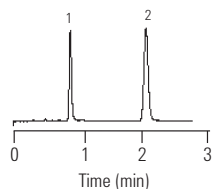

1. Triamcinolone: 0.2 mg/mL
2. Hydrocortisone: 0.3 mg/mL
Minimum Resolution Required = 3.0

Peak	TR	N	Rs
1	0.89	3256	0
2	2.07	4851	11.8

LCPC038

Separation of highly basic antidepressants above their pKa in free base form (pKa 9.5-9.7)


Column: ZORBAX Extend-C18
773450-902
4.6 x 150 mm, 5 µm

Mobile Phase: 75% Methanol / 25% 50 mM Pyrrolidine Buffer, pH 11.5

Flow Rate: 0.5 mL/min

Temperature: 40 °C

Detector: UV, 215 nm

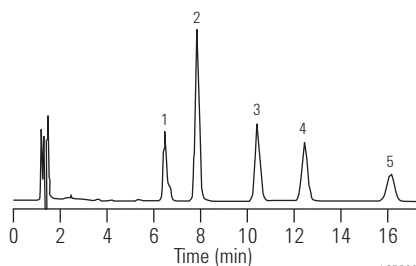

1. Doxepin
2. Imipramine
3. Nortriptyline
4. Amitriptyline
5. Trimipramine

LCBP007

Basic drugs can often be separated in their charged form at low pH with StableBond or at mid-range pH with Eclipse XDB or Bonus -RP columns. With Extend-C18, you can separate at high pH to improve solubility, improve retention, or obtain different selectivity.

Antidepressants, tricyclic: Comparative separation

Column A: ZORBAX Bonus-RP
883668-901
4.6 x 150 mm, 5 µm

Column B: Brand A Polar-linked C8

Column C: Brand B Polar-linked C18




Mobile Phase: ACN: 20 mM Na Citrate, pH 6 (60:40)

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Tricyclic antidepressants (u= uracil)

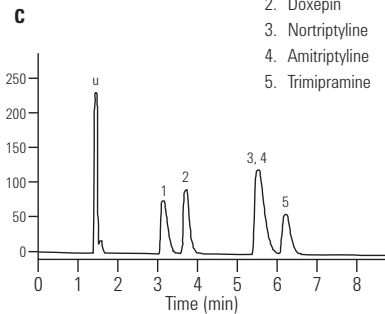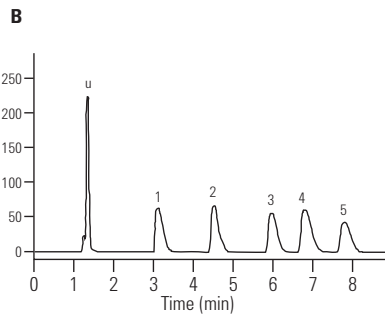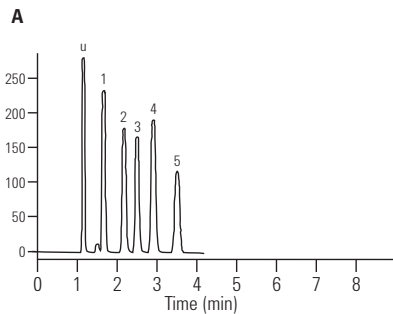

1. Propranolol
2. Doxepin
3. Nortriptyline
4. Amitriptyline
5. Trimipramine

LCBP011

Tricyclic antidepressants

Column: Eclipse XDB-C8
993967-906
4.6 x 150 mm, 5 µm


Mobile Phase: 38/62 THF/25 mM Potassium Phosphate, pH7

Flow Rate: 1.0 mL/min

Temperature: 23 °C

Detector: UV, 254 nm

Sample: 10 µL, Antidepressant mix, 10 µg/mL

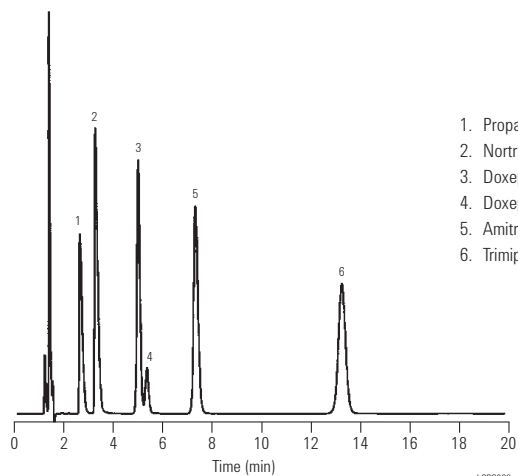

1. Propranolol
2. Nortriptyline
3. Doxepin
4. Doxepin dimer
5. Amitriptyline
6. Trimipramine

LCPC039

Tricyclic antidepressants and metabolites: Effect of pore size

Column A: ZORBAX SB-C18
863953-902
4.6 x 150 mm, 3.5 μ m

Column B: ZORBAX RRHD 300SB-C18
883995-902
4.6 x 150 mm, 5 μ m


Mobile Phase: 40/60, 25 mM Phosphate Buffer,
10 mM Triethylamine, pH 6.2/ACN

Flow Rate: 1.2 mL/min


Temperature: Ambient

Detector: UV, 254 nm

Sample: 10 μ L, Antidepressant mix, 10 μ g/mL

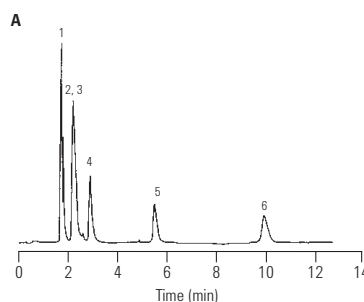

1. trans- 10-OH - Nortriptyline
2. trans- 10-OH - Amitriptyline
3. cis- 10-OH - Nortriptyline
4. cis- 10-OH - Amitriptyline
5. Nortriptyline
6. Amitriptyline

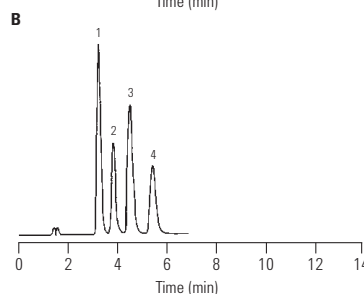

LCPCD40

Ulcer treatment drugs at intermediate pH

Column: ZORBAX Bonus-RP
883668-901
4.6 x 150 mm, 5 μ m


Mobile Phase: Na citrate, 20 mM, pH 6.1: MeOH, (80:20)

Flow Rate: 1.0 mL/min

Temperature: Ambient

Detector: UV, 220 nm

Sample: Ulcer treatment drugs

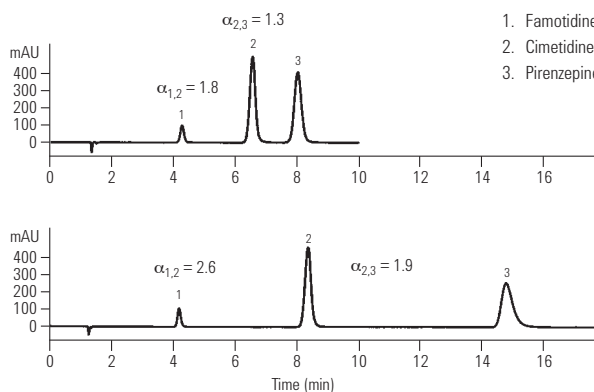

LCPCD42



For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Urine, LSD analysis by LC/MS

Column: Eclipse XDB-C8
960967-906
2.1 x 50 mm, 5 µm

Mobile Phase: 15 : 85, ACN : 10 mM Ammonium Formate, pH 3.7


Flow Rate: 0.3 mL/min

Temperature: 30 °C

Detector: MS

MS Conditions: SIM mode, Ions: 324.2, 223.1, 208.1
Fragmentor (dynamically ramped) 100V at 324.2,
148V at 223.1, 170V at 208.1

Sample: LSD

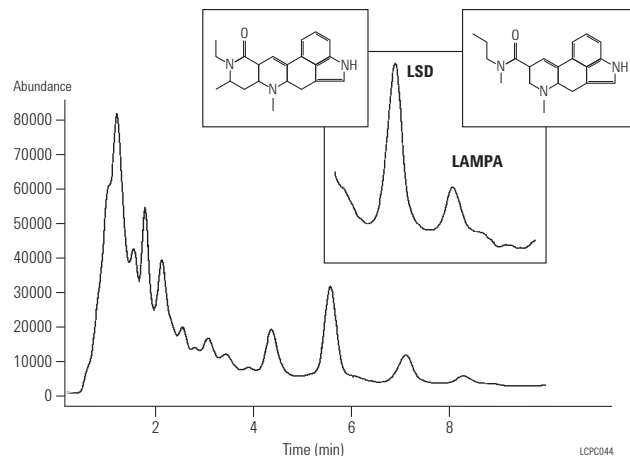

Hughes, J.M., C.A. Miller and S.M. Fischer, "Development of a Method for the Forensic Analysis of LSD in Urine", presented at the ASMS, Palm Springs, June 1997.

USP method: Glyburide and internal standard, progesterone

Column: Eclipse XDB-C8
990967-906
4.6 x 250 mm, 5 µm


Mobile Phase: 45/55, 50 mM Ammonium Phosphate/ACN, Final pH 5.35

Flow Rate: 1.5 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: 5 µL, 10 ug/mL each of standard

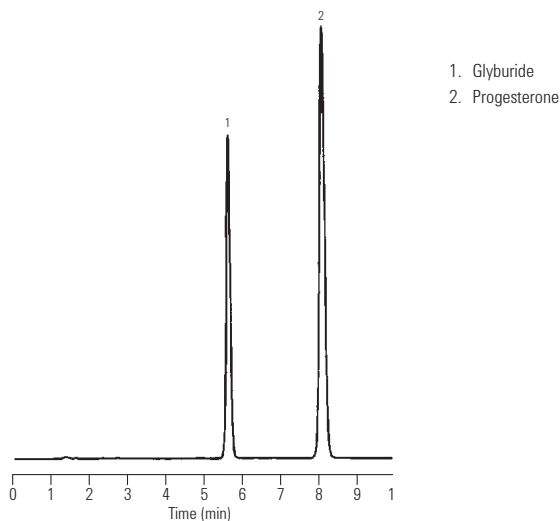

Dexamethasone, USP method: Rapid analysis

Column A: ZORBAX SB-C8
880975-906
4.6 x 250 mm, 5 µm

Column B: ZORBAX Rx/SB-C8
866953-906
4.6 x 75 mm, 3.5 µm


Mobile Phase: A = Water, B = ACN; Isocratic 30% B

Flow Rate: 2.0 mL/min

Temperature: Ambient

Detector: UV, 254 nm

Sample: Dexamethasone
10 µL and 5 µL, 10 µg/mL

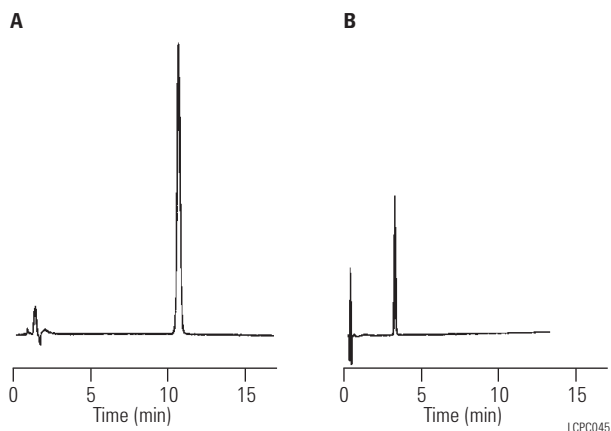

USP analysis of tetracyclines

Column: PLRP-S 100Å
PL1512-5500
4.6 x 250 mm, 5 µm


Sample: 20 mg tetracycline in 25 mL 0.01M HCl

Mobile Phase: 60 g 2-Methyl-2-propanol + 200 mL UHP water +
400 mL 0.2 M K₂HPO₄ at pH 8 + 50 mL 10 g/L
tetrabutylammonium hydrogen sulphate at pH 8 +
10 mL 40 g/L sodium edetate at pH 8, made up to
1000 mL with water (adjust pH with dilute NaOH)

Flow Rate: 1.0 mL/min

Temperature: 60 °C

Detector: UV, 254 nm

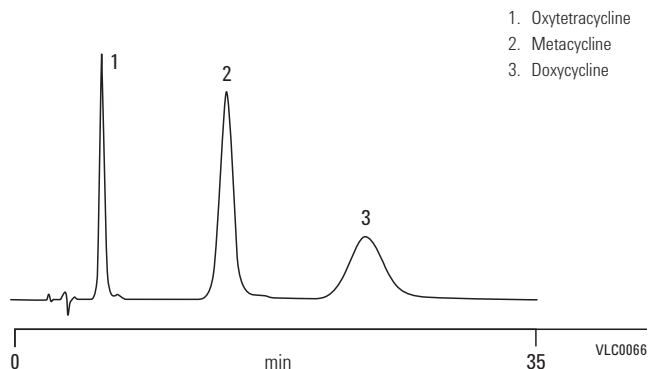

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Warfarin: USP chromatographic purity method using Eclipse XDB-CN

Column: Eclipse XDB-CN
993967-905
4.6 x 150 mm, 5 µm


Mobile Phase: 32:68:1 Acetonitrile:Water:Glacial Acetic Acid

Flow Rate: 1.5 mL/min

Temperature: 25 °C

Detector: UV, 260 nm

Sample: Warfarin, 2 µL

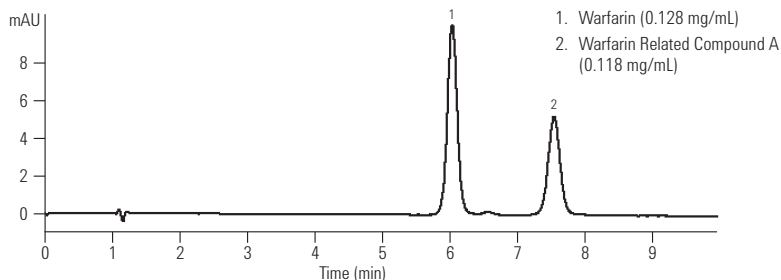

LCPC047

Ten cardiac drugs on Rapid Resolution HT SB-C18

Column: SB-C18
829975-902
4.6 x 150 mm, 1.8 µm

Mobile Phase: A: 0.1% TFA, 5% ACN
B: 0.08% TFA, 95% ACN


Flow Rate: 2 mL/min

Gradient: 0.0 min 12.5% B
10.5 min 60% B
12.0 min 60% B

Temperature: 70 °C

Detector: UV, 230 nm

Sample: Cardiac drugs

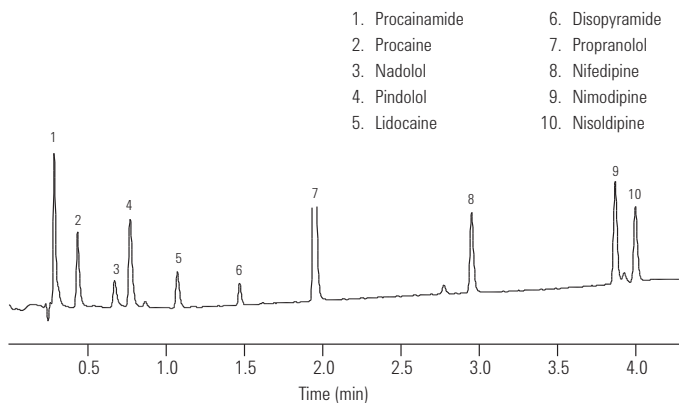

LCPC049

Sulfonamides – Fast analysis with RRHT columns

Column: SB-C18
824700-902
2.1 x 30 mm, 1.8 µm


Mobile Phase: A: 90% 0.1% formic acid
B: 10% 0.1% formic acid in MeOH

Flow Rate: A: 0.2 mL/min
B: 0.8 mL/min

Temperature: 35 °C

Detector: TIC, Single Quad

Sample: Sulfonamides

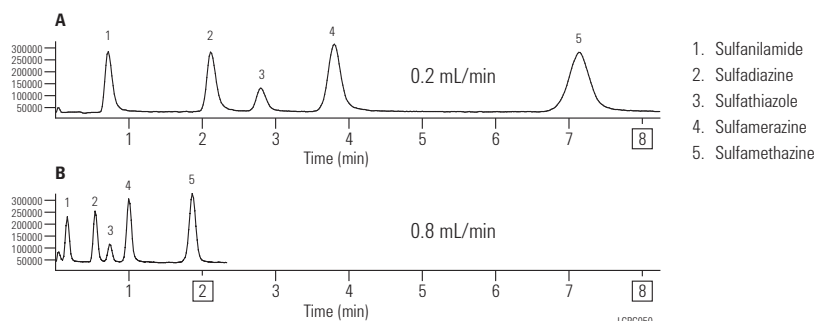

LCPC050

Sulfa drugs

Column: Pursuit XRs Ultra C8
A7511100X020
2.0 x 100 mm, 3.0 µm


Mobile Phase: A: Water+0.1% TFA
B: MeCN+0.1% TFA

Gradient: 10% B for 10 min,
ramp to 45% B in 1 min and hold for 1 min,
return to 10% B in 1 min and hold for 1 min

Flow Rate: 0.65 mL/min

Temperature: Ambient

Detector: UV, 254 nm

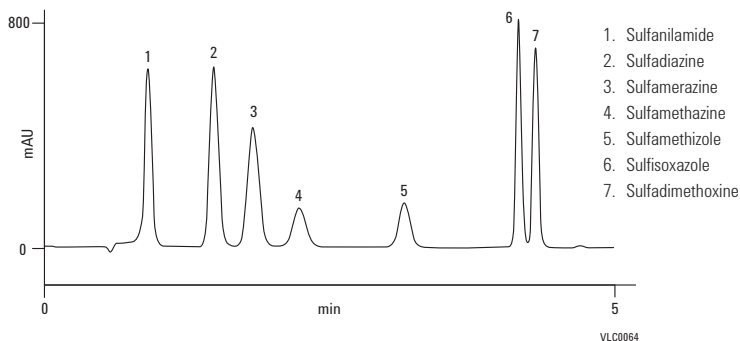


Sulfa drugs

Column: PLRP-S 100Å
PL1111-3500
4.6 x 150 mm, 5 µm

Mobile Phase: Potassium sulfate:
ACN 7:1, pH 2.2

Flow Rate: 1.0 mL/min

Detector: UV, 254 nm

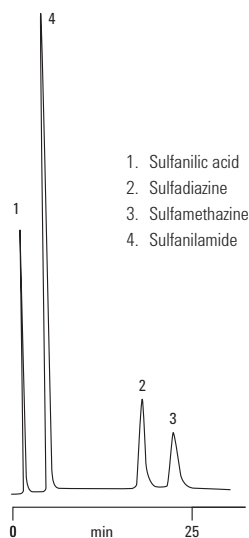


Column: PLRP-S 100Å
PL1111-3500
4.6 x 150 mm, 5 µm

Mobile Phase: Disodium tetraborate: ACN 6:1,
pH 9.3

Flow Rate: 1.0 mL/min

Detector: UV, 254 nm

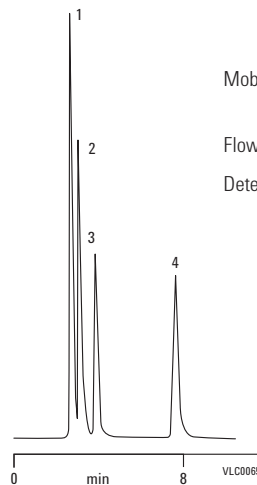

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Fast analysis of Pindolol

Column A: ZORBAX SB-CN
863953-905
4.6 x 150 mm, 3.5 μ m

Column B: ZORBAX SB-CN
827975-905
4.6 x 50 mm, 1.8 μ m


Mobile Phase: A: 70% 50 mM Na Acetate
B: 30% ACN

Flow Rate: 1 mL/min

Temperature: Ambient

Detector: UV, 219 nm

Sample: Pindolol, 2 μ L

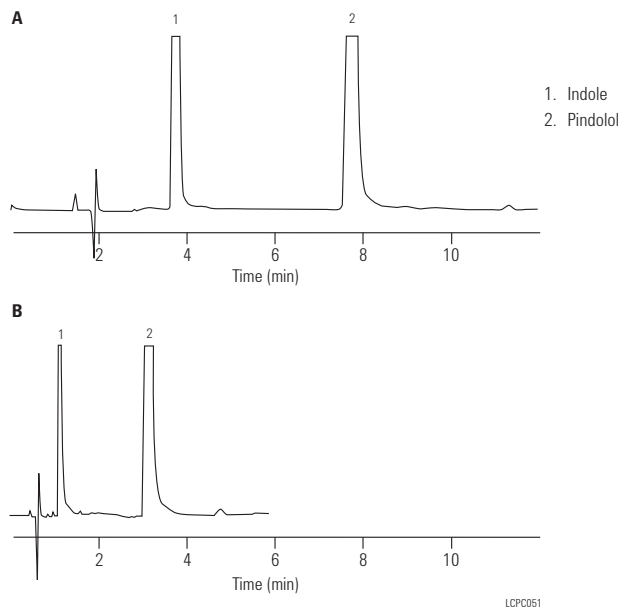

Lamotrigine


Column: Pursuit XRs Ultra C8
A7511100X020
2.0 x 100 mm, 3.0 μ m

Mobile Phase: ACN:water, 25:90 for 1 min

Flow Rate: 0.2 mL/min

Injection Volume: 5 μ L, 50% MeOH

Detector: MS

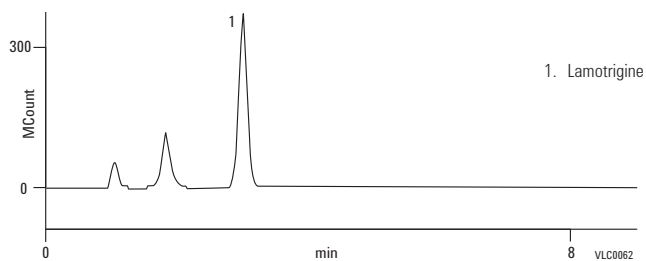

Barbiturates

Column: PLRP-S 100Å
PL1512-5500
4.6 x 250 mm, 5 µm


Mobile Phase: Water

Flow Rate: 1.0 mL/min

Temperature: 200 °C

Detector: UV, 220 nm

Courtesy: Smith, RM, Burgess, RJ, Cheinthaovorn, O and Stuttard, JR (1999) Superheated water: a new look at chromatographic eleunts for reversed-phase liquid chromatography. LCGC Europe, January 1999, 30-36. Used with permission.

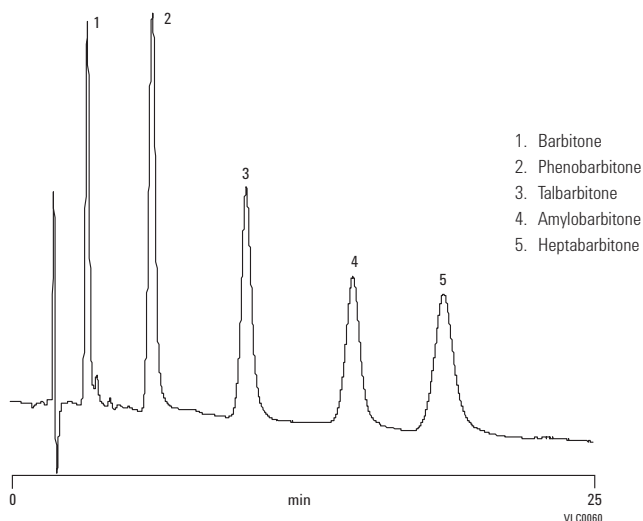

Analysis of ciprofloxacin and ciprofloxacin metabolites

Column: PLRP-S 100Å
PL1111-3500
4.6 x 150 mm, 5 µm


Mobile Phase: 74% 20 mM TCA:22%
ACN:4% MeOH adjusted to pH 3

Flow Rate: 1.0 mL/min

Detector: UV, 277 nm

Krol GJ, Noe, AJ and Beerman, D (1986) Liquid chromatographic analysis of ciprofloxacin and ciprofloxacin metabolites in body fluids. Journal of Liquid Chromatography, 9(13), 2897-2919. Reprinted with permission of the publisher (Taylor & Francis Group, www.informaworld.com).

- Blank urine sample containing known concentrations of internal standard, ciprofloxacin and its metabolites
- Blank urine sample containing only internal standard

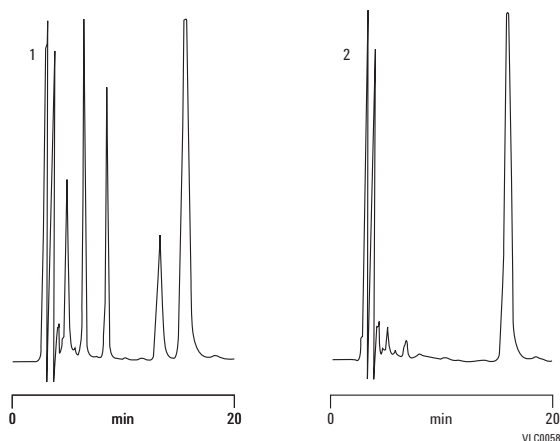

For a comprehensive listing of chromatograms searchable by compound name, visit our online Chromatogram Library at www.agilent.com/chem/library

Product Index

Agilent Supplies

Autosamplers	
Injection Valves	61
Kits	66-67
Loop Capillaries	64
Maintenance Schedule	61
Metering Device	64
Needles and Needle Seats	62-63
Trays	65
Bio-inert Quaternary	
Autosampler	101-102
Column Compartment	102
Detectors	102
Fittings	102
Pumps	101
Valves	102
Capillaries	
1200 and 1100 Prep LC Systems	26
1220/1120 Infinity Series LC Systems	33
1260 Bio-inert LC System	30
1260/1200/1100 Infinity Series LC Systems	23
1290 Infinity Series LC Systems	24
1290 Valve Head	25-27
Loop	32
Miscellaneous	27-29
PEEK Coated Fused Silica Capillaries – 100 µL/min Flow	32
PEEK Coated Fused Silica Capillaries – 20 µL/min Flow	31
PEEK Coated Fused Silica Capillaries for Nano LC	31
Chip LC	100
Detectors	
80 nL and 500 nL Flow Cell	95-96
Diode Array Detector (DAD)/ Multiple Wavelength Detector (MWD)	93
Evaporative Light Scattering Detector (ELSD)	98
Fluorescence Detector (FLD)	99
Lamps	90
Maintenance	89
Maintenance Kits	97
Preparative Flow Cells	96
Refractive Index Detector (RID)	98
Variable Wavelength Detector (VWD)	91-92
Fittings and Unions	
Capillary and Fittings Kits	36-39
Fittings	40-41
How to Align Connection Properly	45
How to Prepare Connection	44
How to Tighten Correctly	43
Fraction Collectors	
Capillary Kits and Needles	70
Collecting Tubes and Trays	69
G1364D Micro Fraction Collector	71
Maintenance Schedule	68
Well Plate Trays	69
In-Line Filters	12
LC/MS	
Chemicals	107
Cleaning Supplies	108
Foreline Pump	107

Gas Purifiers	108
Instrument Supplies	105-106
Maintenance Schedule	103
Preventive Maintenance Kit	104
Quiet Cover	109
Standards	
Analyzer Kit	110
Calibrant Mixes	110
Kits	110
Tools	108
Loop Capillaries	32
Pumps	
1290 Infinity LC Pump Supplies	51
Frits and Adapters	54
Inlet Valves	49
Kits	59-60
Outlet Valves	49
Pistons and Seals	50
Purge Valves	48
Routine Maintenance	47
Safety Caps	55
Seal Wash	52
Solvent Filters	54
Solvent Reservoir	53
Vacuum Degassers	58
Rack for LC Systems	14
Solvent Filters/Degassers	13
Standards	15
Thermostatted Column Compartment	
Capillary Kits	83-88
Tools	11
Troubleshooting	546-549
Tubing	
Accessories	35
PEEK	34
Rigid Capillary Tubing	35
Unions	42
Valves	
Injection Valves	73
Maintenance Notes	72
Manual Injection Valves with Position Sensing Switches	76-77
Sample Loops	78-79
Switching Valves	
External Switching Valve	74
Internal Switching Valve	75
Syringes for Manual Injection	
Fitted Plungers	80
PTFE-Tipped Plungers	81

Applications

BioPharmaceutical	552
Chemical/Industrial	580
Environmental	583
Food and Consumer Products	594
Pharmaceutical	611

CE and CE/MS

Alignment Interface	190
CE/MS Accessories	
Adapter Kit	194
Capillaries	195
Sprayer Kit	195
Capillaries	
Alignment Interface	190
Bulk Fused-Silica	181
CEP Coated	185
Capillary Cassette	190
Capillary Electrochromatography (CEC)	188-189
Cross-linked and Bonded µSIL	186-187
Extended Light Path	
(Bubble Cell) Bare Fused-Silica	179-180
Polyvinyl Alcohol (PVA) Coated	182-184
Standard Bare Fused-Silica	178
Universal Bare Fused-Silica	181
Capillary Cassette	190
High Sensitivity Detection Cell	191-192
Instrument Supplies	199
Solution Kits	
Cation	173
Forensics	175
Inorganic Anion	171-172
Organic Acids	174
µPAGE	176-177
Standards & Reagents	
CZE Buffers for Charged Analytes	196
CZE Buffers for Proteins	197
Capillary Conditioning Solutions	196
MEKC Buffers for Neutral and Charged Analytes	197
Plating Bath Analysis Buffer	197
Ultra Pure CE Water	196
µPAGE Buffer Solutions and Oligo Standards	197
System Start-up and Test Kits	198
Troubleshooting	201-203
Vial Rack	200
Vials and Caps	199
Window Etching Tool	200

Columns for Biomolecule Separations

Affinity Chromatography	
Agilent Bio-Monolith Protein A	434-435
Multiple Affinity Removal System	438-440
Multiple Affinity Removal System Starter Kits	441
Protein Fractionation System	437
Proteomic Reagents	442
mRP-C18 High-Recovery Protein	442-443
Analysis Scale	
Capillary and Nano	452
MicroBore (1.0 mm id)	461
ZORBAX Bio-SCX Series II	458
Citations	
PL-SAX	487
PL-SCX	487
PLRP-S	486-487
Poroshell 300	486
ZORBAX 300	485

Column Selection	
Amino Acid Analysis.....	362
Broad Distribution.....	363
DNA and RNA Oligonucleotide.....	360
Natural and Synthetic Peptides.....	358
Peptide Mapping.....	356
Protein.....	353
Ion-Exchange Chromatography	
Agilent Bio IEX.....	402
Agilent Bio MAB.....	399
Agilent Bio-Monolith.....	412-413
PL-SAX Strong Anion-Exchange.....	406
PL-SCX Strong Cation-Exchange.....	410
Literature.....	477-484
Method Development	
Bio Ion-Exchange Column.....	447
High Sensitivity Capillary Column.....	451
Reversed-Phase LC/MS Column.....	446
SEC Column.....	449-450
ZORBAX Columns.....	444
Purification - Prep HPLC	
PL-SAX and PL-SCX.....	472-473
PLRP-S.....	467-468
Peptide Purification	
VariPure IPE.....	476
VariTide RPC.....	475
ZORBAX PrepHT.....	466
Reversed-Phase HPLC	
PLRP-S.....	387
Poroshell 120.....	385
Poroshell 300.....	380-381
ZORBAX 300Extend-C18.....	376
ZORBAX 300StableBond.....	367
ZORBAX Eclipse	
Amino Acid Analysis (AAA).....	394-395
ZORBAX RRHD 300-Diphenyl.....	374
Size Exclusion Chromatography (SEC)	
Agilent Bio SEC-3.....	418
Agilent Bio SEC-5.....	424
ProSEC 300S.....	428
ZORBAX GF-250 and GF-450.....	431

Columns for Small Molecule Separations

Column Selection	
Column and Mobile Phase Guidelines.....	211-213
Method Development.....	220-221
Quick Guide to Reversed-Phase.....	207-209
Small Molecule Overview.....	218-219
Fast Columns for Reversed-Phase HPLC/UHPLC	
Fast Guards for UHPLC.....	246
Poroshell 120.....	228
ZORBAX Rapid Resolution	
High Definition (RRHD) 1.8 μ m.....	233
ZORBAX Rapid Resolution	
High Throughput (RRHT) 1.8 μ m.....	239
Guard Columns.....	223-224
Kits for Analytical HPLC	
ZORBAX Method Development Kits.....	284-285
ZORBAX Method Validation Kits.....	286-288
Oligo Solutions	
StratoSpheres DNA Cartridges.....	347
TOP Cartridges.....	348-349

Other HPLC Techniques	
Hi-Plex for Carbohydrate Analysis.....	340
ZORBAX HILIC Plus.....	324
ZORBAX Ion-Exchange.....	333
ZORBAX Normal-Phase.....	326-328
Preparative HPLC	
Bulk Materials.....	321
Load & Lock.....	322
Polaris Prep.....	321
Prep LC.....	311-312
Pursuit and Pursuit XRs Prep.....	319
ZORBAX PrepHT.....	314
Reversed-Phase HPLC	
HC-C18(2).....	304-305
PLRP-S.....	306
Polaris.....	298
Pursuit.....	287-288
TC-C18(2).....	304-305
ZORBAX 80Å Extend-C18.....	274
ZORBAX 80Å StableBond.....	264
ZORBAX Bonus-RP.....	278
ZORBAX Eclipse PAH.....	254
ZORBAX Eclipse Plus.....	248
ZORBAX Eclipse XDB.....	256
ZORBAX Original.....	283
ZORBAX Rx.....	272
USP Designations.....	343-346

CrossLab Supplies

Autosampler Maintenance.....	119
Autosampler Syringes.....	116-118
CTC HPLC Autosamplers	
Autosampler Syringes.....	169
Capillaries and Tubing.....	130
Detector Lamps.....	113
Dionex HPLC Systems	
Autosampler Supplies.....	158
Autosampler Syringes.....	157
Buffer Loops.....	160
Capillaries.....	166-168
Detector Lamps.....	157
Fittings, Ferrules and Unions.....	168
In-Line Filters.....	168
Performance Maintenance Kits.....	161-165
Pump Maintenance Procedure.....	126
Pump Supplies.....	158-159
Sample Loops.....	160
Tubing.....	168
Valve Replacement Parts.....	159
Fittings.....	131
In-Line Filters.....	132
Performance Maintenance Kits.....	129
Pump Maintenance.....	124
Pump Supplies.....	122-123
Sealing Mats.....	115
Shimadzu HPLC Systems	
Autosampler Syringes.....	150
Capillaries.....	154
Detector Lamps.....	150
Fittings, Ferrules and Unions.....	155
In-Line Filters.....	156
Performance Maintenance Kits.....	154

Pump Maintenance Procedure.....	125
Pump Supplies.....	151-152
Sample Loops.....	153
Tubing.....	155
Valve Replacement Parts.....	152-153
Troubleshooting.....	133-136
Valve Supplies.....	127-128
Waters HPLC Systems	
Autosampler Syringes.....	137-138
Capillaries.....	148
Detector Lamps.....	137
Detector Supplies.....	142
Fittings, Ferrules and Unions.....	149
Performance Maintenance Kits.....	144-158
Pump Maintenance Procedure.....	124
Pump Supplies.....	139-142
Sample Loops.....	143
Valve Replacement Parts.....	142
Well Plates.....	115

GPC/SEC Columns and Standards

Aqueous	
PL aquagel-OH Analytical.....	525
PL aquagel-OH Preparative.....	528
Column Accessories.....	529
Organic	
PLgel Individual Pore Size.....	505
PLgel MIXED.....	498-501
PLgel MIXED-LS.....	502
PLgel MiniMix.....	504
PLgel Preparative.....	506
Polymer Standards	
EasiCal.....	536
EasiVial.....	532
Polyacrylic Acid.....	545
Polyethylene Glycol/Oxide.....	541-542
Polymethylmethacrylate.....	539-540
Polysaccharides.....	543-544
Polystyrene.....	537-538
Setting Up a GPC/SEC System.....	491-495
Special Application	
EnviroPrep.....	507
MesoPore.....	520
OligoPore.....	521
PL HFIPgel.....	509
PL Rapide.....	510-511
PLgel Olexis.....	508
PlusPore.....	514-515
PolarGel.....	512
PolyPore.....	516
ResiPore.....	518

Part Number Index

0100-0549	107	0101-1250	79	1535-4045	74-75, 77	443905-902	313
0100-0900	42	0101-1251	79	1535-4046	77	443910-901	313
0100-0969	105	0101-1252	79	1535-4860	79	443910-902	313
0100-1259	41	0101-1253	77	1535-4970	104, 107	446905-101	313
0100-1516	40, 92, 96, 99, 106	0101-1254	77	1535-5045	77	446905-102	313
0100-1543	195	0101-1255	77	1535-5082	77	446905-301	313
0100-1597	107	0101-1258	75	160-2644-5	181	446905-302	313
0100-1631	40	0101-1267	61, 73	160-2650-5	181	446905-901	313
0100-1710	11, 58	0101-1268	61, 73	160-2660-5	181	446905-902	313
0100-1847	42	0101-1288	74-75	190-0131	181	446910-102	313
0100-1849	61, 73	0101-1360	75	190-0231	181	449905-101	313
0100-1850	75, 77	0101-1362	75	190-0331	181	449905-102	313
0100-1851	61, 73-75, 77, 101-102	0101-1385	61, 73	190-0431	181	449905-301	313
0100-1852	75	0101-1409	74-75, 77, 102	191-1311	177	449905-302	313
0100-1853	61, 73	0101-1415	74-75	191-3211	177	449905-501	313
0100-1854	75, 108	0101-1416	61, 73, 101	191-5211	177	449905-502	313
0100-1855	75, 104, 108	0101-1417	74-75, 77	192-1311	176	449905-901	313
0100-1859	77	0101-1421	74-75	192-3211	176	449905-902	313
0100-1860	77	0101-1422	61, 73	192-5211	176	449910-902	313
0100-1921	79	0101-2415	61, 73	194-8111	186	5001-3702	82
0100-1922	79	01018-22707	48	196-7203	186	5001-3726	14
0100-1923	79	01018-23702	11	197-7202	186	5001-3743	52
0100-1924	79	01018-60025	54	199-2602	186	5021-1816	28
0100-2051	104, 108	01018-68722	59	204310	442	5021-1817	28
0100-2086	41	01048-87302	28	2140-0585	199	5021-1818	28
0100-2087	75	01078-87302	32, 64	2140-0590	90	5021-1819	28-29
0100-2088	61, 73	01078-87305	28	2140-0600	99	5021-1820	27
0100-2089	75	01080-68702	15	2140-0813	90	5021-1821	27
0100-2175	40	01080-68704	15	2140-0820	90	5021-1822	27
0100-2231	61, 73	01080-83202	40	280959-904	12	5021-1823	27, 92
0100-2233	75	01090-27609	12	280959-907	12	5021-1845	262, 270
0100-2298	42	01090-68702	12	3150-0509	13	5021-1866	11
0100-2304	106	01090-68703	12	3150-0576	13	5022-2133	42
0100-2410	35	01090-87304	28	3150-0577	13	5022-2141	195
0100-2441	42	01090-87306	23, 33	3150-0619	199	5022-2144	42
0101-0376	79	01090-87610	23	3150-0944	54	5022-2145	42
0101-0377	79	01090-87611	23	3162-0178	106	5022-2146	95
0101-0378	79	01100-68700	11	3162-1056	107	5022-2155	58
0101-0379	79	05971-80103	106	3162-1057	107	5022-2159	27
0101-0620	77	05980-60051	108	400510	442	5022-2165	12
0101-0623	77	0890-1727	107	410910-101	313	5022-2166	12
0101-0921	61, 73	0890-1761	34	410910-102	313	5022-2175	64
0101-1050	61, 73	0890-1762	34	410910-301	313	5022-2184	42, 95
0101-1219	79	0890-1763	34, 102	410910-302	313	5022-2185	12
0101-1226	79	0890-1915	34, 106	410910-501	313	5022-2188	51
0101-1227	79	0905-1163	199	410910-502	313	5022-2192	48, 54
0101-1228	79	0905-1175	52	413910-101	313	5022-6503	82
0101-1229	79	0905-1192	48	413910-102	313	5022-6509	29
0101-1230	79	0905-1294	64	413910-301	313	5022-6510	29
0101-1231	77	0905-1420	51	413910-302	313	5022-6531	69
0101-1232	77	0905-1516	54	413910-501	313	5022-6532	69
0101-1233	77	0905-1599	64	413910-502	313	5022-6533	69
0101-1234	79	0905-1717	64	419910-301	313	5022-6534	69
0101-1235	79	0905-1718	52	419910-302	313	5022-6536	41
0101-1236	79	0905-1719	51	419910-501	313	5022-6538	65, 69
0101-1237	79	0905-1731	101	419910-502	313	5022-6539	65, 69
0101-1238	79	12102300	349	420212-901	313	5022-6541	71
0101-1239	79	12102301	349	420212-902	313	5022-6542	71
0101-1240	79	126-1012	187	420420-901	313	5022-6543	71
0101-1241	79	126-1013	187	420910-901	313	5022-6544	71
0101-1242	79	126-1713	187	420910-902	313	5022-6546	71
0101-1243	79	127-1012	187	440905-901	313	5023-0208	71
0101-1244	79	127-1712	187	440905-902	313	5023-0209	71
0101-1245	79	127-1713	187	440910-901	313	5023-0213	71
0101-1246	79	1400-0563	107	440910-902	313	5023-0214	71
0101-1247	79	14251921	349	443905-101	313	5023-0215	71
0101-1248	79	1460-2571	104-105, 108	443905-102	313	5023-0238	71
0101-1249	79	1520-0401	195	443905-901	313	5023-0271	12

5023-0282.....	11	5061-3362.....	35	5064-8264.....	373, 460	5065-9978.....	34, 52
5023-1803.....	42	5062-2418.....	40	5064-8265.....	373, 460	5067-1540.....	41
5041-2168.....	54, 199	5062-2461.....	23	5064-8266.....	373, 460	5067-1547.....	41
5042-1385.....	69	5062-2462.....	23-24, 33-34, 99	5064-8267.....	373, 460	5067-1551.....	12
5042-1386.....	69	5062-2463.....	34, 99	5064-8268.....	373, 460	5067-1553.....	12
5042-1388.....	69	5062-2478.....	395	5064-8269.....	263	5067-1555.....	12
5042-1389.....	69	5062-2479.....	395	5064-8270.....	373, 460	5067-1557.....	41
5042-6454.....	69	5062-2483.....	34, 58	5064-8271.....	263	5067-1558.....	41
5042-6458.....	69	5062-2484.....	52	5064-8273.....	12	5067-1562.....	12
5042-6459.....	69	5062-2486.....	11	5064-8286.....	263	5067-1565.....	54
5042-6461.....	34	5062-8517.....	54, 199	5064-8287.....	263	5067-1581.....	77
5042-6462.....	34	5062-8522.....	23, 92	5064-8288.....	263	5067-1582.....	100
5042-6463.....	34	5062-8524.....	172, 175	5064-8291.....	263	5067-1584.....	100
5042-6470.....	69	5062-8529.....	192	5064-8293.....	64	5067-1585.....	100
5042-6476.....	15	5062-8534.....	58	5064-8294.....	373, 460	5067-1595.....	84
5042-6478.....	199	5062-8535.....	23, 33, 92	5064-8295.....	373, 460	5067-1596.....	85
5042-6491.....	199	5062-8541.....	40	5064-8296.....	263	5067-1597.....	85
5042-6500.....	40	5062-8544.....	199	5064-8297.....	263	5067-4104.....	77
5042-8502.....	69	5062-8562.....	49	5064-8298.....	263	5067-4105.....	77
5042-8507.....	52	5062-8571.....	196	5064-8300.....	373, 460	5067-4107.....	75
5042-8517.....	42, 71	5062-8572.....	196	5065-4402.....	69	5067-4108.....	75
5042-8518.....	42, 71	5062-8573.....	196	5065-4410.....	41	5067-4111.....	75
5042-8519.....	42, 71	5062-8574.....	197	5065-4420.....	15	5067-4112.....	75
5042-8922.....	58	5062-8575.....	172, 196	5065-4421.....	53	5067-4113.....	75
5042-8954.....	52	5062-8576.....	172, 174, 196	5065-4422.....	41, 95	5067-4114.....	61, 73
5042-8957.....	41	5062-8577.....	196	5065-4423.....	41	5067-4117.....	74-75
5042-9954.....	11	5062-8578.....	172-175, 196	5065-4426.....	40	5067-4118.....	74-75
5042-9967.....	11	5062-8587.....	64	5065-4427.....	460	5067-4121.....	75
5043-0221.....	56	5062-8588.....	82	5065-4454.....	40	5067-4124.....	51
5043-0222.....	56	5063-6502.....	77	5065-4459.....	373, 460	5067-4131.....	74, 101
5043-0223.....	56	5063-6506.....	63	5065-4460.....	373, 460	5067-4132.....	74-75
5043-0224.....	56	5063-6510.....	174	5065-4461.....	373, 460	5067-4134.....	74-75
5043-0225.....	56	5063-6511.....	172	5065-4462.....	373, 460	5067-4137.....	75
5043-0226.....	56	5063-6512.....	188	5065-4463.....	373, 460	5067-4141.....	77
5043-0227.....	56	5063-6513.....	188	5065-4464.....	379, 460	5067-4142.....	74-75
5043-0228.....	57	5063-6514.....	198	5065-4465.....	379, 460	5067-4144.....	74-75
5043-0229.....	57	5063-6515.....	198	5065-4466.....	379, 460	5067-4146.....	74-75
5043-0230.....	57	5063-6520.....	198	5065-4467.....	379, 460	5067-4148.....	74-75
5043-0231.....	57	5063-6526.....	82	5065-4468.....	384, 460	5067-4158.....	77
5043-0232.....	56-57	5063-6531.....	53	5065-4498.....	66	5067-4159.....	74-75, 102
5043-0233.....	57	5063-6535.....	188	5065-4499.....	60	5067-4170.....	75
5043-0234.....	56	5063-6536.....	188	5065-4500.....	12	5067-4174.....	48
5043-0235.....	57	5063-6540.....	188	5065-9901.....	12	5067-4202.....	77
5043-0236.....	57	5063-6541.....	188	5065-9908.....	110	5067-4601.....	83-85
5043-0237.....	57	5063-6544.....	188	5065-9910.....	373, 460	5067-4607.....	29
5043-0238.....	57	5063-6586.....	50, 64	5065-9911.....	373, 460	5067-4608.....	28-29
5043-0239.....	57	5063-6589.....	51, 64	5065-9912.....	458, 460	5067-4609.....	29
5043-0242.....	56	5063-6591.....	40	5065-9913.....	373, 460	5067-4633.....	37
5043-0243.....	56	5063-6592.....	95	5065-9914.....	373, 460	5067-4638.....	12
5043-0255.....	57	5063-6593.....	41, 95	5065-9915.....	373, 460	5067-4646.....	38, 84, 86
5043-0272.....	56	5063-6597.....	15, 99	5065-9922.....	60	5067-4647.....	25
5043-0300.....	56	5063-6598.....	58	5065-9923.....	373, 460	5067-4648.....	28
5043-0828.....	56	5063-6599.....	58	5065-9924.....	373, 460	5067-4649.....	25
5043-0829.....	56	5064-8203.....	173	5065-9926.....	29	5067-4650.....	25
5043-0830.....	56	5064-8205.....	173	5065-9927.....	29	5067-4651.....	25
5043-0831.....	56	5064-8206.....	173	5065-9931.....	28	5067-4653.....	26
5043-0832.....	56	5064-8208.....	175	5065-9932.....	28	5067-4657.....	24
5061-3303.....	40, 99	5064-8209.....	175	5065-9933.....	28-29	5067-4658.....	24
5061-3304.....	39	5064-8211.....	11	5065-9935.....	27	5067-4659.....	24
5061-3315.....	39	5064-8220.....	15	5065-9937.....	36	5067-4660.....	24
5061-3327.....	99	5064-8236.....	197	5065-9938.....	36	5067-4661.....	51
5061-3328.....	99	5064-8253.....	271	5065-9939.....	37	5067-4662.....	67
5061-3329.....	99	5064-8254.....	271	5065-9942.....	458, 460	5067-4669.....	24
5061-3330.....	395	5064-8255.....	271	5065-9947.....	37, 82	5067-4670.....	24
5061-3331.....	395	5064-8256.....	271	5065-9948.....	35, 41	5067-4678.....	50
5061-3332.....	395	5064-8257.....	271	5065-9950.....	35, 41	5067-4682.....	38, 86
5061-3333.....	395	5064-8258.....	271, 460	5065-9952.....	52	5067-4684.....	25
5061-3334.....	395	5064-8259.....	373, 460	5065-9963.....	28	5067-4685.....	28
5061-3335.....	395	5064-8260.....	271	5065-9964.....	27	5067-4686.....	25
5061-3337.....	395	5064-8261.....	271	5065-9967.....	40	5067-4687.....	25
5061-3339.....	395	5064-8262.....	271	5065-9971.....	34, 42	5067-4688.....	28
5061-3361.....	35	5064-8263.....	373, 460	5065-9976.....	34	5067-4689.....	26

PART NUMBER INDEX

5067-4695.....	50, 64, 101-102	5069-3637.....	412	5188-6557.....	439	5190-2433.....	405
5067-4697.....	51	5069-3639.....	435	5188-6558.....	439	5190-2434.....	405
5067-4699.....	59	5080-5400.....	108	5188-6559.....	439	5190-2435.....	405
5067-4703.....	32, 64	5133001.....	349	5188-6560.....	439	5190-2436.....	405
5067-4710.....	32, 64	5133005.....	349	5188-6562.....	439	5190-2439.....	405
5067-4716.....	51	5180-4108.....	40, 99	5188-8283.....	441	5190-2440.....	405
5067-4717.....	51	5180-4114.....	40, 99	5188-8825.....	439	5190-2441.....	405
5067-4728.....	48-49	5181-1507.....	199	5188-8826.....	439	5190-2442.....	405
5067-4729.....	39, 87	5181-1512.....	199	5190-0443.....	110	5190-2443.....	405
5067-4730.....	87	5181-1513.....	199	5190-0469.....	110	5190-2444.....	405
5067-4733.....	40	5181-1541.....	71	5190-0488.....	15	5190-2445.....	405
5067-4735.....	25	5181-8836.....	200	5190-0551.....	110	5190-2446.....	405
5067-4737.....	28	5182-0567.....	199	5190-0554.....	110	5190-2447.....	405
5067-4738.....	40	5182-1530.....	90	5190-0555.....	110	5190-2448.....	405
5067-4739.....	40	5182-9697.....	199	5190-0556.....	110	5190-2451.....	405
5067-4741.....	42, 102	5182-9710.....	106	5190-0917.....	90, 199	5190-2452.....	405
5067-4744.....	25	5183-2021.....	285	5190-0924.....	11	5190-2453.....	405
5067-4745.....	25	5183-2022.....	285	5190-1401.....	108	5190-2454.....	405
5067-4746.....	26	5183-4619.....	199	5190-1431.....	284	5190-2455.....	405
5067-4767.....	88	5183-4623.....	199	5190-1432.....	284	5190-2456.....	405
5067-4769.....	88	5183-4624.....	285	5190-1433.....	284	5190-2459.....	405
5067-4777.....	29	5183-4625.....	285	5190-1434.....	284	5190-2460.....	405
5067-4778.....	29	5183-4626.....	285	5190-1435.....	284	5190-2461.....	405
5067-4779.....	29	5183-4627.....	285	5190-1436.....	284	5190-2462.....	405
5067-4780.....	29	5183-4669.....	200	5190-1443.....	104	5190-2463.....	405
5067-4781.....	29	5183-4670.....	200	5190-1480.....	80	5190-2464.....	405
5067-4782.....	29	5185-5807.....	285	5190-1484.....	80	5190-2465.....	405
5067-4798.....	58	5185-5808.....	285	5190-1485.....	80	5190-2466.....	405
5067-4800.....	87	5185-5809.....	285	5190-1486.....	80-81	5190-2467.....	405
5067-5103.....	38, 88	5185-5810.....	285	5190-1492.....	81	5190-2468.....	405
5067-5104.....	25	5185-5920.....	269, 372, 463	5190-1494.....	80	5190-2471.....	405
5067-5106.....	25	5185-5921.....	261	5190-1499.....	81	5190-2472.....	405
5067-5107.....	25	5185-5922.....	282	5190-1501.....	80	5190-2473.....	405
5067-5109.....	26	5185-5923.....	277	5190-1505.....	81	5190-2474.....	405
5067-5110.....	26	5185-5968.....	384, 463	5190-1508.....	80	5190-2475.....	405
5067-5111.....	26	5185-5984.....	439	5190-1512.....	81	5190-2476.....	405
5067-5112.....	26	5185-5985.....	439	5190-1515.....	80	5190-2479.....	405
5067-5113.....	26	5185-5986.....	441	5190-1520.....	81	5190-2480.....	405
5067-5120.....	28	5185-5987.....	441	5190-1522.....	80	5190-2481.....	405
5067-5189.....	39	5185-5988.....	441	5190-1526.....	81	5190-2482.....	405
5067-5378.....	58	5185-5990.....	441	5190-1558.....	81	5190-2483.....	405
5067-5380.....	58	5185-5991.....	441	5190-1560.....	81	5190-2484.....	405
5067-5383.....	58	5188-2743.....	13	5190-1561.....	81	5190-2485.....	405
5068-0001.....	75	5188-2744.....	13	5190-1562.....	81	5190-2486.....	405
5068-0002.....	75	5188-2745.....	13	5190-1564.....	81	5190-2487.....	405
5068-0004.....	48	5188-2746.....	13	5190-1571.....	81	5190-2488.....	405
5068-0005.....	48	5188-5217.....	439	5190-2401.....	401	5190-2491.....	405
5068-0006.....	74-75	5188-5218.....	439	5190-2402.....	401	5190-2492.....	405
5068-0007.....	61, 73	5188-5230.....	439	5190-2403.....	401	5190-2493.....	405
5068-0008.....	74-75	5188-5231.....	442	5190-2404.....	401	5190-2494.....	405
5068-0011.....	74-75	5188-5249.....	441	5190-2405.....	401	5190-2495.....	405
5068-0012.....	74-75	5188-5250.....	441	5190-2406.....	401	5190-2496.....	405
5068-0040.....	74-75, 102	5188-5251.....	441	5190-2407.....	401	5190-2499.....	405
5068-0041.....	74-75, 102	5188-5252.....	441	5190-2408.....	401	5190-2500.....	405
5068-0044.....	74-75	5188-5253.....	441	5190-2411.....	401	5190-2501.....	423
5068-0045.....	74-75, 102	5188-5254.....	441	5190-2412.....	401	5190-2502.....	423
5068-0052.....	77	5188-5289.....	439	5190-2413.....	401	5190-2503.....	423
5068-0053.....	77	5188-5321.....	69	5190-2414.....	401	5190-2504.....	423
5068-0060.....	74-75, 77, 101-102	5188-5322.....	69	5190-2415.....	401	5190-2505.....	423
5068-0067.....	74-75	5188-5332.....	439	5190-2416.....	401	5190-2506.....	423
5068-0076.....	74-75	5188-5333.....	439	5190-2419.....	401	5190-2507.....	423
5068-0077.....	74-75	5188-5334.....	439	5190-2420.....	401	5190-2508.....	423
5068-0082.....	77	5188-5336.....	439	5190-2421.....	405	5190-2509.....	423
5068-0093.....	74-75	5188-5341.....	439	5190-2422.....	405	5190-2510.....	423
5068-0095.....	74-75	5188-6408.....	439	5190-2423.....	405	5190-2511.....	423
5068-0097.....	74-75	5188-6409.....	439	5190-2424.....	405	5190-2512.....	423
5068-0115.....	75	5188-6410.....	439	5190-2425.....	405	5190-2513.....	423
5068-0116.....	75	5188-6411.....	439	5190-2426.....	405	5190-2514.....	423
5068-0122.....	48	5188-6510.....	442	5190-2427.....	405	5190-2515.....	423
5068-0123.....	48	5188-6511.....	442	5190-2428.....	405	5190-2516.....	427
5069-3635.....	412	5188-6523.....	110	5190-2431.....	405	5190-2517.....	427
5069-3636.....	412	5188-6529.....	15	5190-2432.....	405	5190-2518.....	427

5190-2519	427	687975-302	229	726975-302	243, 277	8001-0405	150
5190-2520	427	687975-902	229	726975-902	243, 276	8001-0406	150
5190-2521	427	689775-902	229	727700-902	243, 277	8001-0501	151
5190-2522	427	689775-906	229	727975-302	243, 277	8001-0502	151
5190-2523	427	689775-914	229	727975-902	243, 276	8001-0503	151
5190-2524	427	689975-302	229	728700-902	243, 277	8001-0504	151
5190-2525	427	689975-306	229	728975-302	243, 277	8001-0506	152
5190-2526	427	689975-314	229	728975-902	243, 276	8001-0509	154
5190-2527	427	689975-902	229	735700-902	277	8001-0510	152
5190-2528	427	689975-906	229	735953-902	276	8001-0511	152
5190-2529	427	689975-914	229	735954-302	276	8001-0512	152
5190-2530	427	691775-902	229	746450-902	276	8001-0513	152
5190-2531	427	691775-906	229	75400001	349	8001-0514	151
5190-2532	427	691975-302	229	75700001	349	8001-0515	151
5190-2533	427	691975-306	229	7571901C	349	8001-0516	152
5190-2534	427	691975-902	229	75719025	349	8001-0517	154
5190-2535	427	691975-906	229	75719050	349	8001-0519	152
5190-2536	427	693768-901	229	7572915B	349	8001-0520	152
5190-2537	427	693775-901	325	7572915C	349	8001-0521	152
5190-2538	427	693775-902	229, 386	7573915B	349	8001-0522	151
5190-2539	427	693775-906	229	7573915C	349	8001-0527	152
5190-2540	427	693775-912	229	757700-302	238, 277	8001-0528	152
5190-2541	427	693968-301	229	757700-902	238, 277	8001-0529	152
5190-2542	427	693968-901	229	758700-302	238, 277	8001-0530	152
5190-2543	427	693975-301	325	758700-902	238, 277	8001-0531	152
5190-2544	427	693975-302	229, 386	759700-302	238	8001-0532	152
5190-2545	427	693975-306	229	759700-902	238, 277	8001-0533	151
520518-904	305	693975-312	229	760450-902	277	8001-0534	151
520518-905	305	693975-901	325	761600-902	277	8001-0535	152
588905-902	305	693975-902	229, 386	761753-902	277	8001-0601	153
588915-902	305	693975-906	229	761775-902	379	8001-0603	153
588925-902	305	693975-912	229	761973-902	379	8001-0604	153
588935-902	305	695768-901	229	763600-902	277	8001-0607	153
590-3003	200	695775-901	325	763750-902	379	8001-0608	153
590-4000	177, 197	695775-902	229, 386	763953-902	276	8001-0609	153
590-4001	177, 197	695775-906	229	763954-302	276	8001-0610	154
590-4005	177, 197	695775-912	229	763973-902	379	8001-0612	153
59980-20134	107	695968-301	229	764953-302	276	8001-0613	153
59987-20033	107	695968-901	229	764953-902	276	8001-0614	153
59987-20040	105	695975-302	229, 386	765600-902	277	8001-0615	152
6040-0798	107	695975-306	229	765750-902	379	8001-0701	150
6040-0834	104, 107	695975-312	229	765973-902	379	8001-0702	150
660750-902	210, 384	695975-901	325	766953-902	276	8001-0703	150
660750-906	384	695975-902	210, 229, 386	770050-902	277, 318	8001-0704	150
660750-909	384	695975-906	229	770100-902	277, 318	8001-0705	150
661750-902	384, 463	695975-912	229	770150-902	277, 318	8001-0801	153
661750-906	384, 463	697775-902	229	770450-302	276	8001-0802	153
661750-909	384, 463	697775-906	229	770450-902	276	8001-0803	155
670750-902	384	697975-302	229	770995-902	379	8001-0805	155
671750-902	384, 463	697975-306	229	773450-302	276	8001-0806	155
681775-902	229	697975-902	210, 229	773450-902	276	8001-0807	155
681975-302	229	697975-906	229	773700-902	277	8001-0808	156
681975-902	229	699768-901	229	773995-902	379	8001-0809	153
683775-902	229, 386	699775-901	325	79835-87638	29	8001-0810	154
683775-906	229	699775-902	229	79841-87610	27	8001-0812	153
683775-914	229	699775-906	229	7995108-344	262	8001-0813	155
683975-302	229, 386	699775-912	229	7995108-585	262	8001-0814	153
683975-306	229	699968-301	229	7995108-595	262	8001-0816	155
683975-314	229	699968-901	229	7995118-344	262	8001-0817	155
683975-902	229, 386	699975-301	325	7995118-504	262, 270	8001-0818	154
683975-906	229	699975-302	229	7995118-585	262	8001-0819	155
683975-914	229	699975-306	229	7995118-595	262	8001-0821	154
685775-902	229, 386	699975-312	229	7995208-344	270	8001-0822	154
685775-906	229	699975-901	325	7995208-585	270	8001-0823	155
685775-914	229	699975-902	229	7995208-595	270	8001-0824	155
685975-302	229, 386	699975-906	229	7995218-344	270	8001-8020	154
685975-306	229	699975-912	229	7995218-585	270-271	8002-0401	157
685975-314	229	722975-902	245	7995218-595	270	8002-0402	157
685975-902	229, 386	724700-902	243, 277	7995230-344	262	8002-0403	157
685975-906	229	724975-302	243, 277	8001-0401	150	8002-0404	157
685975-914	229	724975-902	243, 276	8001-0402	150	8002-0405	157
687775-902	229	726700-902	243, 277	8001-0403	150	8002-0406	157

PART NUMBER INDEX

8002-0407	157	8002-0917	161	8005-0911	145	820950-908	330
8002-0408	157	8002-0918	162	8005-0912	145	820950-911	433
8002-0412	157	8002-0919	162	8005-0913	144	820950-912	283
8002-0413	157	8002-0921	168	8005-0914	146	820950-913	273
8002-0414	157	8002-0923	164	8005-0915	144	820950-914	273
8002-0415	157	8002-0924	164	8005-0916	144	820950-915	269
8002-0501	158	8005-0414	138	8005-0925	145	820950-916	269
8002-0502	158	8005-0416	138	8005-0926	145	820950-917	269
8002-0515	158	8005-0417	138	8005-0927	146	820950-918	372
8002-0516	159	8005-0418	138	8005-0928	144	820950-919	329
8002-0517	159	8005-0419	138	8005-0929	145	820950-920	269
8002-0601	158	8005-0420	138	8005-0930	144	820950-921	372
8002-0602	159	8005-0422	138	8010-0440	169	820950-922	269
8002-0603	159	8005-0423	138	8010-0441	137, 169	820950-923	372
8002-0604	159	8005-0508	141	8010-0442	137, 169	820950-924	283, 372
8002-0605	159	8005-0512	142	8010-0443	137, 169	820950-925	261
8002-0607	159	8005-0513	141	8010-0444	137, 169	820950-926	261
8002-0608	159	8005-0514	139	8010-0445	137, 169	820950-927	261
8002-0610	165	8005-0515	140	8010-0446	137, 169	820950-928	282
8002-0611	165	8005-0516	140	8010-0448	137, 169	820950-930	277
8002-0701	157	8005-0523	140	8010-0449	169	820950-931	394
8002-0702	157	8005-0524	141	8010-0450	137, 169	820950-932	379
8002-0703	157	8005-0525	141	8010-0455	137, 169	820950-933	269
8002-0704	157	8005-0526	141	8010-0456	137, 169	820950-935	261, 329
8002-0705	157	8005-0527	140	8010-0457	137, 169	820950-936	253
8002-0706	157	8005-0528	142, 145	8010-0458	137, 169	820950-937	253
8002-0802	168	8005-0529	140	8010-0459	137, 169	820950-938	253
8002-0803	168	8005-0530	140	8010-0460	137, 169	820950-939	253, 255
8002-0805	168	8005-0531	140	8010-0467	137, 169	820999-901	225, 253, 255, 261, 269, 273, 277, 282-283, 305, 329-330, 334, 372, 379, 384, 394, 433
8002-0806	168	8005-0532	142	8010-0468	137, 169		
8002-0808	158	8005-0533	140	820212-911	433		
8002-0809	158	8005-0535	139	820212-914	273, 318		
8002-0810	158	8005-0536	139	820212-915	269, 273, 317	821075-918	384
8002-0811	160	8005-0537	141	820212-918	317, 373, 466	821075-920	384
8002-0815	167	8005-0538	139	820212-919	318, 329	821075-924	384
8002-0816	167	8005-0539	140	820212-920	269, 317	821125-915	269, 273
8002-0817	167	8005-0540	140	820212-921	317, 373, 466	821125-918	372
8002-0818	166	8005-0541	139	820212-924	317, 373, 466	821125-919	329
8002-0819	166	8005-0601	142	820212-925	261, 318	821125-924	269, 372
8002-0820	166	8005-0602	142	820212-926	261, 318	821125-926	261
8002-0821	166	8005-0603	142	820212-928	282, 318	821125-928	282
8002-0822	166	8005-0604	142	820212-930	277, 318	821125-930	277
8002-0823	166	8005-0605	142	820212-933	269, 317	821125-932	379
8002-0824	166	8005-0702	137	820385-901	318	821125-933	269
8002-0825	166	8005-0704	137	820400-901	261, 269, 273, 277, 282-283, 313, 317-318, 329-330, 373, 433, 466	821125-935	261, 329
8002-0826	166	8005-0705	137			821125-936	253
8002-0831	166	8005-0812	148	820444-901	261, 269, 273, 277, 282, 313, 317-318, 329-330, 373, 433, 466	821125-937	253
8002-0832	166	8005-0822	148			821125-938	253
8002-0833	166	8005-0823	148			821125-939	253, 255
8002-0834	166	8005-0824	148			821700-902	245, 270
8002-0835	167	8005-0825	148	820555-901	245, 263, 270	821700-932	245, 270
8002-0837	166	8005-0826	148	820675-111	330, 433	821725-901	238, 243, 246, 253
8002-0856	160	8005-0835	148-149	820675-112	261	821725-902	238, 244, 246, 269
8002-0857	160	8005-0836	149	820675-115	269, 273, 283	821725-903	238, 243, 246, 261
8002-0858	160	8005-0837	149	820675-119	329-330	821725-904	238, 244, 246, 269
8002-0859	160	8005-0838	143, 149, 156	820675-124	269, 283, 372	821725-911	229, 246
8002-0860	160	8005-0839	143	820700-902	244, 269	821725-912	229, 246
8002-0901	164	8005-0840	143	820700-905	244, 269	821725-913	229, 246
8002-0902	164	8005-0841	143	820700-906	244, 269	821725-914	229, 246
8002-0903	164	8005-0842	143	820700-912	244, 269	821975-902	245, 270
8002-0904	163	8005-0843	143	820750-901	243, 246, 252	821975-932	245, 270
8002-0905	164	8005-0844	143	820750-902	244, 246, 267	822700-902	270
8002-0906	163	8005-0845	143	820750-903	243, 246, 260	822700-932	270
8002-0907	161	8005-0846	143	820750-904	244, 246, 267	822975-902	244-245, 270
8002-0908	161	8005-0901	146	820750-911	229, 246	822975-906	245, 270
8002-0909	162	8005-0902	146	820750-912	229, 246	822975-932	270
8002-0910	162	8005-0903	146	820750-913	229, 246	823700-902	245, 270
8002-0911	161	8005-0904	146	820750-914	229, 246	823700-932	245, 270
8002-0912	162	8005-0905	147	820950-901	330	823750-901	238, 243, 246, 253
8002-0913	163	8005-0906	147	820950-902	283	823750-902	238, 244, 246, 268
8002-0915	165	8005-0907	146	820950-905	283, 330	823750-903	238, 243, 246, 260
8002-0916	165	8005-0908	147	820950-906	283	823750-904	238, 244, 246, 268

823750-911	229, 246	829975-302	244, 268	858700-314	238	863750-906	372
823750-912	229, 246	829975-305	244, 268	858700-902	238, 268	863953-902	267
823750-913	229, 246	829975-306	244, 268	858700-905	238, 268	863953-905	267
823750-914	229, 246	829975-312	244, 268	858700-906	238, 268	863953-906	267
823975-902	245, 270	829975-902	244, 267	858700-912	238, 268	863953-912	267
823975-932	245, 270	829975-905	244, 267	858700-914	238	863953-914	267
824700-902	244, 269	829975-906	244, 267	858750-902	372	863954-302	267
824700-905	244, 269	829975-912	244, 267	858750-906	372	863954-305	267
824700-906	244, 269	829975-914	244, 267	858750-909	372	863954-306	267
824700-912	244, 269	830668-901	244, 281	858750-944	372, 374-375	863954-312	267
824700-914	244, 269	830975-906	244, 267	858768-901	238, 281	863954-314	267
824975-302	244, 268	830990-902	268	859700-302	238, 268	863967-302	273
824975-305	244, 268	830990-906	268	859700-306	238, 268	863967-902	273
824975-306	244, 268	830990-914	268	859700-902	238, 268	863973-902	210, 372
824975-902	244, 267	831975-902	270	859700-905	238, 268	863973-905	372
824975-905	244, 267	831975-906	270	859700-906	238, 268	863973-906	372
824975-906	244, 267	831975-932	270	859700-912	238, 268	863973-909	372
824975-912	244, 267	831975-936	270	859700-914	238	863974-302	372
824975-914	244, 267	832975-902	267	859768-901	238, 281	863974-306	372
825700-902	245, 270	832975-906	267	860700-304	334	863974-309	372
825700-932	245, 270	833975-902	270	860700-704	334	864668-301	281
825975-902	245, 270	833975-906	270	860700-708	330	864668-901	281
825975-932	245, 270	833975-912	270	860950-902	372	865600-902	269
826700-902	244, 269	833975-932	270	860950-905	372	865600-906	269
826700-906	244, 269	833975-936	270	860950-906	372	865608-901	282
826975-302	244, 268	834975-902	267	860950-909	372	865630-902	372, 463
826975-306	244, 268	834975-906	267	860975-902	268	865630-906	372, 463
826975-902	244, 267	835668-901	281	860975-905	268	865750-902	372
826975-906	244, 267	835975-902	267	860975-906	268	865750-906	372
827668-301	244, 281	835975-905	267	860975-909	268	865973-902	210, 372
827668-901	244, 281	835975-906	267	860975-912	268	865973-905	372
827700-901	244, 329	835975-912	267	860975-914	268	865973-906	372
827700-902	244, 269	835975-914	267	861600-902	269	865973-909	372
827700-905	244, 269	840140-901	261, 269, 273, 283, 329-330, 372, 433	861600-906	269	8660-0827	108
827700-906	244, 269			861608-901	282	866668-901	281
827700-912	244, 269	840300-908	330	861630-902	372, 463	866735-902	268
827700-914	244, 269	843300-908	330	861700-901	282	866953-302	267
827768-901	244, 282	846952-704	334	861753-902	268	866953-902	267
827975-301	244, 329	846975-202	267	861753-905	268	866953-905	267
827975-302	244, 268	846975-902	267	861753-906	268	866953-906	267
827975-305	244, 268	846975-906	267	861753-912	268	866953-912	267
827975-306	244, 268	846975-914	267	861753-914	268	866953-914	267
827975-312	244, 268	8500-1867	107	861767-902	273	866967-902	273
827975-314	244, 268	8500-2236	107	861768-901	282	868050-901	282, 318
827975-901	244	8500-4410	195	861775-902	372	868100-901	282, 318
827975-902	244, 267, 329	8500-6762	15	861775-906	372	868150-901	282, 318
827975-905	244, 267	8500-6782	196	861953-902	267	870050-902	269, 317
827975-906	244, 267	8500-6785	174	861953-905	267	870050-906	269, 317
827975-912	244, 267	8500-6786	197	861953-906	267	870050-914	269, 317
827975-914	244, 267	8500-6787	197	861953-912	267	870100-902	269, 317
828668-301	244, 281	8500-6797	172	861953-914	267	870100-906	269, 317
828668-901	244, 281	8500-6900	174	861954-302	267	870100-914	269, 317
828700-901	244, 329	8500-6917	110	861954-305	267	870150-902	269, 317
828700-902	244, 269	857700-302	238, 268	861954-306	267	870150-906	269, 317
828700-905	244, 269	857700-305	238, 268	861954-309	267	870150-914	269, 317
828700-906	244, 269	857700-306	238, 268	861954-312	267	8710-0004	108
828700-912	244, 269	857700-312	238, 268	861954-314	267	8710-0510	99, 108
828700-914	244, 269	857700-314	238	861967-302	273	8710-0806	108
828768-901	244, 282	857700-902	238, 268	861967-902	273	8710-1534	95
828975-301	244, 329	857700-905	238, 268	861971-901	281	8710-1615	108
828975-302	244, 268	857700-906	238, 268	861973-306	372	8710-1622	108
828975-305	244, 268	857700-912	238, 268	861973-902	372	8710-1924	11, 48
828975-306	244, 268	857700-914	238	861973-906	372	8710-1930	11, 35, 58, 108
828975-309	268	857750-902	372	863600-902	269	8710-1931	11, 35
828975-312	244, 268	857750-906	372	863600-905	269	8710-2699	108
828975-314	244, 268	857750-909	372	863600-906	269	871700-902	268
828975-901	244, 329	857750-944	372, 375	863608-901	282	871700-906	268
828975-902	244, 267	857768-901	238, 281	863630-902	372, 463	871700-914	268
828975-905	244, 267	858700-302	238, 268	863630-906	372, 463	872700-902	268
828975-906	244, 267	858700-305	238, 268	863668-301	281	872700-906	268
828975-912	244, 267	858700-306	238, 268	863668-901	281	873700-902	270
828975-914	244, 267	858700-312	268	863700-901	282	873700-906	270

PART NUMBER INDEX

873700-932	270	880995-206	372	895150-909	317, 373, 466	933975-902	262
873700-936	270	880995-209	372	897150-102	317, 373, 466	933975-906	262
874700-902	268	880995-902	372	897150-106	317, 373, 466	933975-932	262
874700-906	268	880995-905	372	897150-109	317, 373, 466	933975-936	262
875700-902	270	880995-906	372	897250-102	317, 373, 466	934967-902	260
875700-906	270	880995-909	372	897250-105	317, 373, 466	934967-906	260
875700-932	270	881750-902	372	897250-106	317, 373, 466	935967-902	260
875700-936	270	883668-301	281	897250-109	317, 373, 466	935967-906	260
877150-102	269, 317	883668-901	281	921700-902	245, 263	935967-912	260
877150-106	269, 317	883700-704	334	921700-932	245, 263	946975-902	260
877150-114	269, 317	883700-714	334	921975-902	245, 263	946975-906	260
877250-101	318, 329	883700-901	329	921975-932	245, 263	959701-902	253
877250-102	269, 317	883700-902	273	922700-902	245, 262	959701-906	253
877250-105	269, 317	883700-905	268	922700-932	245, 262	959701-912	253
877250-106	269, 317	883700-906	268	922975-902	243, 245, 262	959701-918	253, 255
877250-112	269, 317	883700-909	268	922975-906	245, 262	959731-902	243, 253
877250-114	269, 317	883700-912	268	922975-932	245, 262	959731-906	243, 253
877952-101	317, 330	883700-922	268	923700-902	245, 263	959731-912	243, 253
877952-102	283, 317	883725-901	281	923700-932	245, 263	959733-902	253
877952-105	283, 317, 330	883750-902	372	923975-902	245, 262	959733-906	253
877952-106	283, 317	883750-905	372	923975-932	245, 262	959733-912	253
877952-108	317, 330	883750-906	372	924700-902	243, 261	959741-902	243, 253
877967-102	273, 318	883750-909	372	924700-906	243, 261	959741-906	243, 253
877974-901	433	883952-302	283	924975-302	243, 260	959741-912	243, 253
877974-910	433	883952-701	330	924975-306	243, 260	959741-918	243, 253, 255
878150-101	282, 318	883952-702	283	924975-902	243, 260	959743-901	325
878250-101	282, 318	883952-703	334	924975-906	243, 260	959743-902	253
880668-301	281	883952-704	334	925700-902	245, 263	959743-906	253
880668-901	281	883952-705	330	925700-932	245, 263	959743-912	253
880952-201	330	883952-706	283	925975-902	245, 262	959746-902	253
880952-202	283	883952-708	330	925975-932	245, 262	959746-906	253
880952-203	334	883952-710	283	926700-902	243, 261	959757-301	238, 325
880952-204	334	883952-712	283	926700-906	243, 261	959757-302	238, 252
880952-205	330	883952-714	334	926975-302	243, 260	959757-306	238, 252
880952-206	283	883967-302	273	926975-306	243, 260	959757-312	238
880952-208	330	883967-901	273	926975-902	243, 260	959757-901	238, 325
880952-302	283	883967-902	273	926975-906	243, 260	959757-902	210, 238, 253
880952-701	330	883975-202	267	927700-902	243, 261	959757-906	238, 253
880952-702	283	883975-302	267	927700-906	243, 261	959757-912	238
880952-703	334	883975-305	267	927975-302	243, 260	959758-301	238, 325
880952-704	334	883975-306	267	927975-306	243, 260	959758-302	238, 252
880952-705	330	883975-309	267	927975-902	243, 260	959758-306	238, 252
880952-706	283	883975-312	267	927975-906	243, 260	959758-312	238
880952-708	330	883975-314	267	928700-902	243, 261	959758-901	238, 325
880952-710	283	883975-901	329	928700-906	243, 261	959758-902	238, 253
880952-712	283	883975-902	267	928975-302	243, 260	959758-906	238, 253
880952-714	334	883975-905	267	928975-306	260	959758-912	238
880967-201	267, 273	883975-906	267	928975-902	243, 260	959759-302	238, 252
880967-202	273	883975-909	267	928975-906	260	959759-306	238, 252
880967-302	273	883975-912	267	9300-1747	199	959759-901	238, 325
880967-901	273	883975-914	267	9300-1748	199	959759-902	238, 253
880967-902	273	883995-902	210, 372	9301-0407	99	959759-906	238, 253
880975-201	329	883995-905	372	9301-0656	53	959759-912	238
880975-202	267	883995-906	372	9301-0722	200	959763-902	253
880975-205	267	883995-909	372	9301-0895	13	959763-906	253
880975-209	267	884950-507	283	9301-0978	199	959763-912	253
880975-212	267	884950-526	283	9301-1291	106	959764-902	243, 253
880975-302	267	884950-543	283	9301-1337	58	959764-906	243, 253
880975-305	267	884950-567	267	9301-1420	53	959764-912	243, 253
880975-306	267	884950-577	281	9301-1421	53	959764-918	243, 253, 255
880975-309	267	884973-701	433	9301-1446	99	959790-918	253, 255
880975-312	267	884973-901	433	9301-1450	53	959793-901	325
880975-314	267	884973-902	433	9301-6341	53	959793-902	253
880975-901	329	884975-202	267	9301-6342	53	959793-906	253
880975-902	267	895050-902	317, 373, 466	930990-902	260	959793-912	253
880975-905	267	895050-906	317, 373, 466	930990-906	260	959793-918	253, 255
880975-906	267	895050-909	317, 373, 466	931975-902	262	959794-902	243
880975-909	267	895100-902	317, 373, 466	931975-906	262	959931-902	243, 252
880975-912	267	895100-906	317, 373, 466	931975-932	262	959931-906	243, 252
880975-914	267	895100-909	317, 373, 466	931975-936	262	959931-912	243, 252
880995-202	372	895150-902	317, 373, 466	932967-902	260	959931-918	243, 252, 255
880995-205	372	895150-906	317, 373, 466	932967-906	260	959933-902	252

959933-906	252	961967-905	260	993967-905	260, 329	A2001150X046	300
959933-912	252	961967-906	260	993967-906	260	A2001200X030	301
959936-902	252	963400-902	394	993967-912	260	A2001250C046	302
959936-906	252	963600-902	261	A2000020X020	301	A2001250X020	301
959936-912	252	963600-906	261	A2000030X020	301	A2001250X030	301
959941-302	243, 252	963954-302	260	A2000030X046	300	A2001250X046	300
959941-306	243, 252	963954-305	260	A2000050X020	301	A2001MG	303
959941-312	243, 252	963954-306	260	A2000050X030	301	A2001MG1	303
959941-902	243, 252	963954-312	260	A2000050X046	300	A2002250X046	300
959941-906	243, 252	963967-902	260	A2000100C020	302	A2002250X212	300, 321
959941-912	243, 252	963967-905	260	A2000100C030	302	A2002250X500	300
959941-918	243, 252, 255	963967-906	260	A2000100C046	302	A2002MG	303
959943-901	325	963967-912	260	A2000100R030	302	A2003020X020	301, 332
959943-902	252	965600-902	261	A2000100R046	302	A2003030X020	301, 332
959943-906	252	965600-906	261	A2000100T030	302	A2003050X020	301, 332
959943-912	252	966400-902	394	A2000100T046	302	A2003050X030	301, 332
959943-918	252, 255	966735-902	261	A2000100X020	301	A2003050X046	300, 332
959946-902	252	966954-302	260	A2000100X030	301	A2003050X212	300, 331
959946-906	252	966967-902	260	A2000100X046	300	A2003100X020	301, 332
959951-902	243, 252	966967-905	260	A2000100X212	300	A2003100X030	301, 332
959961-302	252	966967-906	260	A2000100X300	300	A2003100X046	300, 332
959961-306	252	966967-912	260	A2000125X040	301	A2003125X040	301, 332
959961-312	252	970050-902	261, 318	A2000150C046	302	A2003150X020	301, 332
959961-901	325	970050-906	261, 318	A2000150R030	302	A2003150X030	301, 332
959961-902	252	970100-902	261, 318	A2000150R046	302	A2003150X040	301, 332
959961-906	252	970100-906	261, 318	A2000150T030	302	A2003150X046	300, 331-332
959961-912	252	970150-902	261, 318	A2000150T046	302	A2003250X020	301, 332
959961-918	252, 255	970150-906	261, 318	A2000150X020	301	A2003250X030	301, 332
959963-302	252	971700-902	261	A2000150X030	301	A2003250X040	301, 332
959963-306	252	971700-906	261	A2000150X040	301	A2003250X046	300, 332
959963-312	252	972700-902	261	A2000150X046	300	A2003250X212	300, 321, 331
959963-902	210, 252	972700-906	261	A2000150X212	300	A2003MG	303, 332
959963-906	252	973700-902	262	A2000200X046	300	A2003MG2	303, 332
959963-912	252	973700-906	262	A2000250C030	302	A2004250X212	300, 321, 331
959963-918	252, 255	973700-932	262	A2000250C046	302	A2004250X500	300, 331
959964-302	243, 252	973700-936	262	A2000250R046	302	A2004MG	303, 332
959964-306	243, 252	974700-902	261	A2000250T046	302	A2004MG2	303, 332
959964-312	243, 252	974700-906	261	A2000250X020	301	A2005020X020	301, 332
959964-902	243, 252	975700-902	262	A2000250X030	301	A2005030X020	301, 332
959964-906	243, 252	975700-906	262	A2000250X040	301	A2005050X020	301, 332
959964-912	243, 252	975700-932	262	A2000250X046	300	A2005050X030	301, 332
959964-918	243, 252, 255	975700-936	262	A2000250X100	300, 321	A2005050X046	301, 332
959990-318	252, 255	977150-102	261, 318	A2000250X212	300, 321	A2005100X020	301, 332
959990-902	252	977150-106	261, 318	A2000MG	303	A2005100X030	301, 332
959990-906	252	977250-102	261, 318	A2000MG2	303	A2005100X046	300, 332
959990-912	252	977250-106	261, 318	A2001020X020	301	A2005150X020	301, 332
959990-918	252, 255	981757-302	238, 260	A2001030X020	301	A2005150X030	301, 332
959993-302	252	981757-902	238, 261	A2001030X030	301	A2005150X046	300, 332
959993-306	252	981758-302	238, 260	A2001030X046	301	A2005250X020	301, 332
959993-902	252	981758-902	238, 261	A2001050C020	302	A2005250X046	300-301, 332
959993-906	252	981759-302	238, 260	A2001050R020	302	A2005MG	303, 332
959993-912	252	981759-902	238, 261	A2001050T020	302	A2005MG2	303, 332
959993-918	252, 255	990967-202	260	A2001050X020	301	A2006030X020	301
959994-902	243	990967-206	260	A2001050X030	301	A2006050X020	301
959996-902	252	990967-302	260	A2001050X046	301	A2006050X046	300
959996-906	252	990967-305	260	A2001075X046	300	A2006100X020	301
959996-912	252	990967-306	260	A2001100C020	302	A2006100X030	301
959996-918	252, 255	990967-312	260	A2001100R030	302	A2006100X046	300
960967-902	260	990967-902	260	A2001100R046	302	A2006150X020	301
960967-905	260	990967-905	260, 329	A2001100T030	302	A2006150X030	301
960967-906	260	990967-906	260	A2001100T046	302	A2006150X046	300
960967-912	260	990967-912	260	A2001100X020	301	A2006250X020	301
961400-302	394	993400-902	394	A2001100X030	301	A2006250X030	301
961600-902	261	993700-902	260	A2001100X046	300	A2006250X046	300
961600-906	261	993700-905	260, 329	A2001150C020	302	A2006250X100	300, 321
961753-902	261	993700-906	260	A2001150C046	302	A2006250X212	300
961753-905	261	993700-912	260	A2001150R020	302	A2006MG	303
961753-906	261	993967-302	260	A2001150R046	302	A2006MG2	303
961967-302	260	993967-305	260	A2001150T020	302	A2007030X020	301
961967-306	260	993967-306	260	A2001150T046	302	A2007030X030	300
961967-312	260	993967-312	260	A2001150X020	301	A2007050X020	301
961967-902	260	993967-902	260	A2001150X030	301	A2007050X030	301

PART NUMBER INDEX

A2007050X046.....	301	A2014250X046.....	300, 332	A3000150T046.....	294	A3030250X046.....	291
A2007100X020.....	301	A2014MG.....	303, 332	A3000150X020.....	292	A3030250X100.....	291, 320
A2007100X030.....	301	A2014MG2.....	303, 332	A3000150X030.....	291	A3030MG.....	294
A2007100X046.....	300	A2020050X020.....	301	A3000150X039.....	291	A3030MG2.....	294
A2007150X020.....	301	A2020050X046.....	300	A3000150X046.....	291	A3031030X020.....	292
A2007150X030.....	301	A2020100X030.....	301	A3000150X100.....	291	A3031050R020.....	294
A2007150X046.....	300	A2020125X040.....	301	A3000150X212.....	292	A3031050T020.....	294
A2007250X020.....	301	A2020150X020.....	301	A3000250C020.....	293	A3031050X020.....	292
A2007250X030.....	301	A2020150X030.....	301	A3000250C030.....	293	A3031100C046.....	293
A2007250X046.....	300	A2020150X040.....	301	A3000250C046.....	293	A3031100X020.....	292
A2007MG.....	303	A2020150X046.....	300	A3000250X020.....	292	A3031100X046.....	291
A2007MG2.....	303	A2020250X020.....	301	A3000250X030.....	291	A3031150C046.....	293
A2008250X100.....	300	A2020250X030.....	301	A3000250X040.....	291	A3031150R046.....	294
A2008250X212.....	300	A2020250X040.....	301	A3000250X046.....	291	A3031150T046.....	294
A2008MG.....	303	A2020250X046.....	300	A3000250X100.....	291, 320	A3031150X020.....	292
A2010050X020.....	301	A2020250X100.....	300, 321	A3000250X212.....	292	A3031150X046.....	291
A2010100X030.....	301	A2020250X212.....	300	A3000300X039.....	291	A3031250C046.....	293
A2010100X046.....	300	A2020MG.....	303	A3000MG.....	294	A3031250X046.....	291
A2010125X040.....	301	A2020MG2.....	303	A3000MG1.....	294	A3031MG.....	294
A2010150X020.....	301	A2021050X020.....	301	A3000MG2.....	294	A3031MG2.....	294
A2010150X030.....	301	A2021050X030.....	301	A3001020X020.....	292	A3032100X046.....	291
A2010150X040.....	301	A2021050X046.....	301	A3001030X020.....	292	A3032150X046.....	291
A2010150X046.....	300	A2021050X100.....	300	A3001030X046.....	291	A3032250X046.....	291
A2010250X030.....	301	A2021075X020.....	301	A3001050C020.....	293	A3032250X100.....	291, 320
A2010250X040.....	301	A2021100X020.....	301	A3001050C046.....	293	A3032250X212.....	292, 320
A2010250X046.....	300	A2021150X020.....	301	A3001050R046.....	294	A3032250X500.....	292
A2010250X100.....	321	A2021150X030.....	301	A3001050T046.....	294	A3040030X020.....	292
A2010250X212.....	300, 321	A2021250X020.....	301	A3001050X020.....	292	A3040050X020.....	292
A2010MG.....	303	A2021250X046.....	300	A3001050X030.....	292	A3040050X046.....	291
A2010MG2.....	303	A2021MG.....	303	A3001050X046.....	291	A3040100X020.....	292
A2011030X030.....	301	A2021MG2.....	303	A3001100C020.....	293	A3040100X046.....	291
A2011050X020.....	301	A2030050X020.....	301	A3001100C030.....	293	A3040100X212.....	292
A2011075X046.....	300	A2030100X030.....	301	A3001100C046.....	293	A3040150X020.....	292
A2011100X046.....	300	A2030125X040.....	301	A3001100R030.....	294	A3040150X030.....	291
A2011150X020.....	301	A2030150X020.....	301	A3001100T030.....	294	A3040150X046.....	291
A2011150X046.....	300	A2030150X030.....	301	A3001100X020.....	292	A3040250X030.....	291
A2011250X020.....	301	A2030150X040.....	301	A3001100X030.....	291	A3040250X046.....	291
A2011MG.....	303	A2030150X046.....	300	A3001100X046.....	291	A3040250X100.....	320
A2011MG2.....	303	A2030250X020.....	301	A3001150C020.....	293	A3040MG.....	294
A2013020X020.....	301, 332	A2030250X030.....	301	A3001150C030.....	293	A3040MG1.....	294
A2013030X020.....	301, 332	A2030250X040.....	301	A3001150C046.....	293	A3040MG2.....	294
A2013050X020.....	301, 332	A2030250X046.....	300	A3001150R030.....	294	A3041020X020.....	292
A2013050X046.....	300, 332	A2030250X100.....	300, 321	A3001150T030.....	294	A3041030X020.....	292
A2013100X020.....	301, 332	A2030250X212.....	300, 321	A3001150X020.....	292	A3041050X020.....	292
A2013100X030.....	301, 332	A2030MG.....	303	A3001150X030.....	291	A3041050X030.....	292
A2013100X046.....	300, 332	A2031050X020.....	301	A3001150X046.....	291	A3041050X046.....	291
A2013125X040.....	301, 332	A2031050X030.....	301	A3001250X020.....	292	A3041100X020.....	292
A2013150X020.....	301, 332	A2031050X046.....	301	A3001250X030.....	291	A3041100X030.....	291
A2013150X030.....	301, 332	A2031100X020.....	301	A3001250X046.....	291	A3041100X046.....	291
A2013150X040.....	301, 332	A2031150X020.....	301	A3001MG.....	294	A3041150X020.....	292
A2013150X046.....	300, 332	A2031250X020.....	301	A3001MG2.....	294	A3041150X030.....	291
A2013250X020.....	301, 332	A2031250X046.....	300	A3002100X046.....	291	A3041150X046.....	291
A2013250X030.....	301, 332	A2031MG2.....	303	A3002150X046.....	291	A3041200X020.....	292
A2013250X040.....	301, 332	A3000020X020.....	292	A3002150X212.....	292	A3041250X020.....	292
A2013250X046.....	300, 332	A3000030X020.....	292	A3002250X046.....	291	A3041250X046.....	291
A2013250X100.....	300, 321, 331	A3000050X020.....	292	A3002250X100.....	291	A3041MG.....	294
A2013250X212.....	300, 321, 331	A3000050X046.....	291	A3002250X212.....	292	A3041MG1.....	294
A2013MG.....	303, 332	A3000100C020.....	293	A3002250X500.....	292	A3041MG2.....	294
A2013MG2.....	303, 332	A3000100C030.....	293	A3002300X039.....	291	A3050020X020.....	292
A2014020X020.....	301, 332	A3000100C046.....	293	A3002MG.....	294	A3050030X020.....	292
A2014030X020.....	301, 332	A3000100R030.....	294	A3002MG2.....	294	A3050050X020.....	292
A2014050X020.....	301, 332	A3000100T030.....	294	A3030050X020.....	292	A3050050X046.....	291
A2014050X030.....	301, 332	A3000100X020.....	292	A3030100C046.....	293	A3050100X020.....	292
A2014050X046.....	301, 332	A3000100X030.....	291	A3030100X020.....	292	A3050100X030.....	291
A2014100X020.....	301, 332	A3000100X046.....	291	A3030100X046.....	291	A3050100X046.....	291
A2014100X030.....	301, 332	A3000125X040.....	291	A3030150C020.....	293	A3050150X030.....	291
A2014100X046.....	300, 332	A3000150C020.....	293	A3030150C046.....	293	A3050150X046.....	291
A2014150X020.....	301, 332	A3000150C030.....	293	A3030150R046.....	294	A3050150X100.....	291
A2014150X030.....	301, 332	A3000150C046.....	293	A3030150T046.....	294	A3050150X212.....	292
A2014150X046.....	300, 332	A3000150R030.....	294	A3030150X020.....	292	A3050250X046.....	291
A2014250X020.....	301, 332	A3000150R046.....	294	A3030150X046.....	291	A3050250X100.....	291, 320
A2014250X030.....	301, 332	A3000150T030.....	294	A3030250C046.....	293	A3050250X212.....	292

A3050MG	294	A6004250X212	296, 320, 331	A7000100R030	294	G1313-27302	65
A3050MG2	294	A6004250X300	296, 320, 331	A7000100T030	294	G1313-43204	63
A3051020X020	292	A6004250X500	296, 320, 331	A7000150C046	293	G1313-43216	35
A3051030X020	292	A6004MG	297	A7000150R046	294	G1313-44510	65
A3051050X020	292	A6005050X020	296, 331	A7000150T046	294	G1313-44512	65
A3051050X030	292	A6005050X046	295, 331	A7000150X046	291	G1313-44513	65
A3051050X046	291	A6005100X046	295, 331	A7000250C046	293	G1313-60004	65
A3051100X020	292	A6006050X046	295, 331	A7000250R046	294	G1313-68709	66
A3051100X030	291	A6006100X021	296, 331	A7000250T046	294	G1313-68711	67
A3051100X046	291	A6006100X046	295, 331	A7000250X046	291	G1313-68719	66
A3051150X020	292	A6010050X020	296	A7000MG3	297	G1313-68730	66
A3051150X030	291	A6010100X020	296	A7001100C046	293	G1313-87102	62
A3051150X046	291	A6010100X030	295	A7001100R030	294	G1313-87201	62
A3051250X046	291	A6010100X046	295	A7001100R046	294	G1313-87202	62
A3051MG	294	A6010100X212	320	A7001100T030	294	G1313-87203	62
A3051MG2	294	A6010150X020	296	A7001100T046	294	G1313-87300	34
A6000030X020	296	A6010150X030	295	A7001100X020	292	G1313-87303	32, 64
A6000030X212	296	A6010150X040	295	A7001100X030	291	G1313-87304	23, 82
A6000050X020	296	A6010150X046	295	A7001100X046	291	G1313-87305	23, 33, 82
A6000050X046	295	A6010150X212	320	A7001150C046	293	G1314-60081	91-92
A6000050X100	295	A6010250X030	295	A7001150R046	294	G1314-60082	91-92
A6000050X212	296, 320	A6010250X040	295	A7001150T046	294	G1314-60083	91-92
A6000050X300	296	A6010250X046	295	A7001MG3	297	G1314-60086	91-92
A6000100X020	296	A6010250X300	320	A7501030X020	297	G1314-60087	91-92
A6000100X030	295	A6010MG	297	A7501050X020	297	G1314-60100	90
A6000100X046	295	A6010MG2	297	A7501100X020	297	G1314-60101	90
A6000100X212	296, 320	A6011050X020	296	A7501100X030	297	G1314-60182	91-92
A6000100X300	296, 320	A6011050X030	295	A7501150X020	297	G1314-60183	91-92
A6000150X020	296	A6011050X046	295	A7501150X030	297	G1314-60186	91-92
A6000150X030	295	A6011100X020	296	A7511030X020	297	G1314-60187	91-92
A6000150X040	295	A6011100X030	295	A7511050X020	297	G1314-65052	91, 97
A6000150X046	295	A6011100X046	295	A7511100X020	297	G1314-65054	91, 97
A6000150X100	295	A6011150X020	296	A7511150X030	297	G1314-65056	97
A6000150X212	296, 320	A6011150X030	295	A7521030X020	297	G1314-65061	91, 97
A6000150X300	296, 320	A6011MG	297	A7521050X020	297	G1314-87301	92
A6000250X020	296	A6011MG2	297	A7521100X020	297	G1314-87302	92
A6000250X030	295	A6012250X212	296	BHT-4	104, 108	G1315-27705	96
A6000250X040	295	A6020050X020	296	BMT-4	108	G1315-45003	11, 95
A6000250X046	295	A6020050X046	295	G1103-60001	90	G1315-60011	94-95
A6000250X100	295	A6020100X030	295	G1156-68711	83	G1315-60012	94-95
A6000250X212	296, 320	A6020100X046	295	G1156-68712	83	G1315-60015	94-95
A6000250X300	296, 320	A6020100X212	296, 320	G1156-68713	83	G1315-60016	94, 96
A6000MG	297	A6020150X020	296	G1156-68714	83	G1315-60017	94, 96
A6000MG2	297	A6020150X030	295	G1160-68706	83	G1315-60018	94, 96
A6001020X020	296	A6020150X046	295	G1310-68730	59	G1315-60022	94-95
A6001030X030	295	A6020150X300	296	G1310-68731	59	G1315-60024	94-95
A6001030X046	295	A6020250X020	296	G1310-68741	59	G1315-60025	94-95
A6001050X020	296	A6020250X030	295	G1310-68742	59	G1315-67301	96
A6001050X030	295	A6020250X046	295	G1311-60003	53	G1315-67302	96
A6001050X046	295	A6020250X100	295	G1311-60006	12	G1315-68703	31
A6001100X010	296	A6020250X212	296, 320	G1311-60009	48	G1315-68708	31-32
A6001100X020	296	A6020MG	297	G1311-68705	60	G1315-68712	97
A6001100X030	295	A6020MG2	297	G1311-68710	59	G1315-68713	91, 97
A6001100X046	295	A6021030X020	296	G1311-68711	59	G1315-68715	96-97
A6001150X010	296	A6021030X046	295	G1312-60020	49	G1315-68716	94, 96
A6001150X020	296	A6021050X020	296	G1312-60025	49	G1315-68724	94, 96
A6001150X030	295	A6021050X030	295	G1312-60061	48	G1315-68725	96-97
A6001150X046	295	A6021050X046	295	G1312-60066	49	G1315-80001	96
A6001250X020	296	A6021100X010	296	G1312-60067	49	G1315-80002	96
A6001250X046	295	A6021100X020	296	G1312-67305	23	G1315-80003	96
A6001MG	297	A6021100X030	295	G1312-67500	23	G1315-80004	96
A6001MG2	297	A6021100X046	295	G1312-68711	59	G1315-87101	96
A6002050X046S	295	A6021150X020	296	G1312-68716	53	G1315-87302	95
A6002150X300	320	A6021150X030	295	G1312-68726	60	G1315-87303	23
A6002250X046	295	A6021150X046	295	G1312-68730	60	G1315-87305	96
A6002250X100	295, 320	A6021150X100	295	G1312-68741	59	G1315-87306	95
A6002250X212	296, 320	A6021250X020	296	G1312-68755	60	G1315-87311	23, 33, 92, 99
A6002250X300	320	A6021250X046	295	G1312-87303	23	G1315-87312	23
A6002250X500	296, 320	A6021MG	297	G1312-87304	23	G1315-87313	96
A6002MG	297	A6021MG2	297	G1312-87305	28	G1315-87318	96
A6004250X046	295, 331	A6022250X500	320	G1312-87306	28	G1315-87319	95
A6004250X100	295, 331	A7000100C030	293	G1312-87330	54	G1315-87321	95

PART NUMBER INDEX

G1315-87323.....	96	G1364-68723.....	70	G1600-60211.....	178	G1947-20029.....	105
G1315-87325.....	95	G1364-81701.....	71	G1600-60230.....	190	G1947-60103.....	105
G1315-87328.....	96	G1364-83205.....	71	G1600-60232.....	180	G1956-20302.....	106
G1315-87333.....	96	G1364-84516.....	69	G1600-60233.....	180	G1956-80000.....	106
G1315-87338.....	96	G1364-84521.....	69	G1600-60310.....	190	G1958-60098.....	105
G1315-87339.....	95	G1364-84522.....	69	G1600-60311.....	178	G1958-60136.....	105
G1316-27301.....	28	G1364-84523.....	69	G1600-60330.....	190	G1960-80039.....	107
G1316-60001.....	82	G1364-84524.....	69	G1600-60332.....	180	G1960-80060.....	105
G1316-67005.....	75	G1364-84525.....	69	G1600-60400.....	194	G1969-20302.....	106
G1316-67006.....	75	G1364-84531.....	69	G1600-60411.....	178	G1969-60086.....	108
G1316-67007.....	75	G1364-84532.....	69	G1600-60419.....	184	G1969-85000.....	110
G1316-67009.....	75	G1364-86711.....	71	G1600-61132.....	180	G1969-85001.....	107
G1316-68708.....	83	G1364-87201.....	70	G1600-61211.....	178	G1969-85003.....	15, 107
G1316-68710.....	83	G1364-87202.....	70	G1600-61219.....	184	G1969-85010.....	110
G1316-68711.....	83	G1364-87304.....	70-71	G1600-61232.....	173, 180	G1969-85020.....	110
G1316-68716.....	38, 82	G1364-87305.....	70-71	G1600-61239.....	184	G1969-85026.....	107
G1316-68721.....	83	G1364-87306.....	70-71	G1600-61311.....	178	G1972-60025.....	71
G1316-68744.....	82	G1367-60001.....	69	G1600-61332.....	180	G1978-85000.....	15, 107
G1316-80002.....	82	G1367-68730.....	66	G1600-61411.....	178	G1982-85001.....	15
G1316-80003.....	82	G1367-68734.....	66	G1600-61419.....	184	G1982-85002.....	15
G1316-80004.....	82	G1367-68741.....	66	G1600-62132.....	180	G1982-85003.....	15
G1316-83200.....	82	G1367-87012.....	63	G1600-62211.....	172, 178	G2228-68700.....	59
G1316-87300.....	23, 33, 82	G1367-87017.....	63	G1600-62232.....	180	G2250-04500.....	65
G1316-87303.....	27, 82	G1367-87101.....	62	G1600-62311.....	174, 178	G2250-04501.....	65
G1316-87305.....	27	G1367-87102.....	62	G1600-62318.....	185	G2250-04502.....	65
G1316-87306.....	28	G1367-87200.....	62, 70	G1600-62332.....	180	G2250-04503.....	65
G1316-87309.....	23	G1367-87201.....	62	G1600-62402.....	199	G2250-04504.....	65
G1316-87312.....	27	G1367-87202.....	63	G1600-62411.....	178	G2255-68700.....	65, 69
G1316-87313.....	27	G1367-87300.....	32, 64	G1600-62700.....	190	G2255-68709.....	65
G1316-87314.....	27	G1367-87304.....	28	G1600-63200.....	192	G2255-68710.....	65
G1316-87316.....	27	G1375-87301.....	31	G1600-63211.....	178	G2255-68720.....	65
G1316-87317.....	27	G1375-87302.....	31	G1600-63311.....	178	G2255-68730.....	65
G1316-87318.....	27	G1375-87303.....	32, 64	G1600-63411.....	178	G2258-23201.....	54
G1316-87319.....	27	G1375-87304.....	31	G1600-64211.....	175, 178	G2258-60003.....	64
G1316-87321.....	28	G1375-87305.....	32	G1600-64232.....	180	G2258-60011.....	65
G1316-87323.....	28	G1375-87306.....	32	G1600-64311.....	178	G2258-68710.....	63
G1321-60005.....	99	G1375-87308.....	32	G1600-64332.....	180	G2258-87102.....	63
G1321-60007.....	99	G1375-87309.....	31-32	G1600-64411.....	178	G2258-87307.....	34
G1321-60015.....	99	G1375-87310.....	31	G1600-67201.....	199	G2258-87310.....	34
G1322-67300.....	58	G1375-87311.....	32	G1600-67219.....	184, 195	G2258-87311.....	34
G1322-68705.....	58, 60	G1375-87312.....	32	G1600-67220.....	195	G2258-87312.....	34
G1328-87600.....	23	G1375-87315.....	32, 64	G1600-67311.....	194-195	G2258-87313.....	34
G1329-60011.....	65	G1375-87320.....	31	G1600-67312.....	195	G2258-87314.....	34
G1329-68718.....	67	G1375-87321.....	31	G1600-67319.....	184, 195	G2258-87315.....	34
G1329-68727.....	67	G1375-87322.....	31	G1600-68319.....	184, 192	G2258-87316.....	34
G1329-68736.....	67	G1375-87323.....	31	G1600-68714.....	192	G2260-68711.....	32, 64
G1329-68737.....	67	G1375-87324.....	31	G1600-68715.....	192	G2260-87101.....	63
G1329-80001.....	62	G1375-87325.....	31	G1600-68716.....	192	G2260-87201.....	63
G1329-87012.....	62	G1375-87326.....	34	G1600-68723.....	192	G2260-87300.....	26
G1329-87017.....	62	G1375-87327.....	31	G1603A.....	194	G2260-87301.....	26
G1329-87101.....	62	G1376-60003.....	53	G1607-20030.....	195	G2421-60001.....	110
G1329-87103.....	62	G1376-60005.....	59	G1607-60000.....	195	G2423A.....	107
G1329-87300.....	23	G1376-68705.....	60	G1607-60001.....	195	G2424A.....	107
G1329-87302.....	32, 64	G1376-68707.....	60	G1607-60041.....	195	G2425A.....	107
G1353-68750.....	75	G1376-68710.....	60	G1607A.....	195	G2426A.....	107
G1361-22402.....	50	G1377-44900.....	63	G160U-60419.....	184	G2427A.....	105
G1361-23204.....	54	G1377-87000.....	63	G160U-61219.....	184	G2428A.....	105
G1361-23205.....	54	G1377-87001.....	63	G160U-61239.....	184	G2431A.....	110
G1361-60012.....	49	G1377-87002.....	63	G160U-61419.....	184	G2432A.....	110
G1361-60022.....	53	G1377-87201.....	63	G1946-00034.....	107	G2441-80010.....	106
G1361-67302.....	26	G1377-87300.....	32, 64	G1946-20215.....	108	G2453-85050.....	107
G1361-68707.....	60	G1377-87310.....	32, 64	G1946-20301.....	105	G2453-85060.....	107
G1361-68710.....	60	G1379-67310.....	58	G1946-60037.....	105	G2455-85001.....	15
G1362-68706.....	98	G1600-23223.....	199	G1946-60098.....	105	G2571-80103.....	106
G1362-68709.....	98	G1600-60002.....	190	G1946-60157.....	108	G3199B.....	109
G1362-87300.....	98	G1600-60007.....	199	G1946-60180.....	108	G4203-68708.....	11
G1362-87301.....	98	G1600-60013.....	194	G1946-80009.....	105	G4204-40000.....	51
G1364-27107.....	71	G1600-60027.....	192	G1946-80019.....	106	G4204-40005.....	51
G1364-60021.....	71	G1600-60033.....	199	G1946-80049.....	108	G4204-60004.....	51
G1364-68706.....	71	G1600-60132.....	180	G1946-80054.....	108	G4204-60022.....	49
G1364-68711.....	70	G1600-60150.....	190	G1946-85004.....	107	G4208-68700.....	11
G1364-68712.....	70	G1600-60210.....	190	G1946-85021.....	107	G4212-60007.....	94

G4212-60008.....	94	G5611-67300.....	30	PL1110-6320.....	505	PL1210-6120.....	506
G4212-60011.....	94	G5611-67301.....	30	PL1110-6400.....	508	PL1210-6120EPA.....	507
G4212-60022.....	51	G5611-68710.....	37	PL1110-6500.....	500	PL1210-6125.....	506
G4212-60032.....	94	G5611-68741.....	59, 101	PL1110-6504.....	500	PL1210-6130.....	506
G4212-60038.....	94	G5615-60005.....	102	PL1110-6515.....	505	PL1210-6140.....	506
G4212-68001.....	51	G5615-60017.....	102	PL1110-6520.....	505	PL1210-6150.....	506
G4216-68711.....	67	G5615-60018.....	102	PL1110-6525.....	505	PL1210-6160.....	506
G4218-20000.....	98	G5615-60022.....	102	PL1110-6530.....	505	PL1212-1102.....	470
G4218-20001.....	98	G5616-60050.....	30, 102	PL1110-6540.....	505	PL1212-1103.....	470
G4218-20002.....	98	G5664-86703.....	102	PL1110-6550.....	505	PL1212-3100.....	470
G4218-20003.....	98	G5664-86706.....	102	PL1111-3500.....	310, 393	PL1212-3101.....	470
G4218-20004.....	98	G5667-40500.....	102	PL1113-1300.....	518	PL1212-3702.....	470
G4218-40000.....	98	G5667-60310.....	32	PL1113-1320.....	521	PL1212-3703.....	470
G4218-40010.....	98	G5667-60320.....	64, 102	PL1113-1325.....	520	PL1212-3800.....	470
G4218-40011.....	98	G5667-60500.....	30	PL1113-1500.....	516	PL1212-3801.....	470
G4218-40100.....	98	G5667-60501.....	30	PL1113-3100.....	511	PL1212-6100.....	470
G4218-40110.....	98	G5667-60502.....	30	PL1113-3120.....	511	PL1212-6101.....	470
G4218-40130.....	98	G5667-60503.....	30	PL1113-3300.....	511	PL1212-6200.....	470
G4218-40150.....	98	G5667-60504.....	30	PL1113-3500.....	511	PL1212-6201.....	470
G4218-40220.....	98	G5667-60505.....	30	PL1113-6300.....	518	PL1212-6400.....	470
G4218-60100.....	98	G5667-87017.....	63, 102	PL1113-6325.....	520	PL1212-6401.....	470
G4218-68010.....	98	G5667-87200.....	63, 101	PL1113-6500.....	516	PL1212-6800.....	470
G4218-85000.....	15, 98	G6011A.....	109	PL1113-6520.....	521	PL1212-6801.....	470
G4220-20012.....	11	G6012A.....	109	PL1114-1900HFIP.....	509	PL1213-6520.....	521
G4220-24013.....	52	G6013A.....	109	PL1114-6900HFIP.....	509	PL1220-6130.....	528
G4220-26210.....	52	G6014A.....	109	PL1117-1800.....	512	PL1245-1102.....	474
G4220-60006.....	51	G7100-60002.....	190	PL1117-1830.....	512	PL1245-1103.....	411, 474
G4220-60007.....	53	G7100-60007.....	199	PL1117-6800.....	512	PL1245-3102.....	474
G4220-60012.....	51	G7100-60033.....	199	PL1117-6830.....	512	PL1245-3103.....	411, 474
G4220-60015.....	51	G7100-60150.....	190	PL1120-3830.....	511	PL1245-3702.....	411, 474
G4220-60016.....	51	G7100-60210.....	190	PL1120-6520.....	525	PL1245-3703.....	411, 474
G4220-60022.....	49	G7100-60230.....	190	PL1120-6830.....	525	PL1249-1120.....	525, 528
G4220-60028.....	49	G7100-60310.....	190	PL1145-1802.....	474	PL1249-6100.....	528
G4220-60035.....	51	G7100-60330.....	190	PL1145-1803.....	474	PL1249-6140.....	528
G4220-63010.....	52	G7100-60400.....	190, 194	PL1147-1501.....	428	PL1249-6150.....	528
G4220-63015.....	52	G7100-62700.....	190	PL1147-6501.....	428	PL1251-1102.....	409, 474
G4226-60021.....	65	G7100-68705.....	200	PL1149-1530.....	525	PL1251-1103.....	409, 474
G4226-60310.....	32, 64	G7100-68723.....	192	PL1149-1840.....	525	PL1251-3102.....	409, 474
G4226-67001.....	67	PCG931AAKIT.....	322	PL1149-3800.....	511	PL1251-3103.....	409, 474
G4226-68735.....	66	PCG932AAKIT.....	322	PL1149-6240.....	525	PL1251-3702.....	409, 474
G4226-87012.....	63	PCG933AAKIT.....	322	PL1149-6250.....	525	PL1251-3703.....	409, 474
G4226-87020.....	63	PCG93LL500X25.....	322	PL1149-6260.....	525	PL1310-0001.....	529
G4226-87201.....	63	PCG93LL500X25WJ.....	322	PL1149-6800.....	525	PL1310-0002.....	529
G4240-23705.....	100	PCG93LL500X50.....	322	PL1149-6801.....	525	PL1310-0005.....	529
G4240-25206.....	100	PCG93LL500X50WJ.....	322	PL1149-6820.....	525	PL1310-0007.....	529
G4240-43200.....	100	PCG93LL500X75.....	322	PL1149-6840.....	525	PL1310-0008.....	529
G4240-87300.....	100	PCG93LL500X75WJ.....	322	PL1149-6850.....	525	PL1310-0012.....	529
G4240-87301.....	100	PCG93LLSTAND123.....	322	PL1149-6860.....	525	PL1310-0016.....	310, 342, 393
G4240-87302.....	100	PL1012-5A05.....	475	PL1151-1802.....	474	PL1310-0036.....	529
G4240-87303.....	100	PL1013-2100.....	511	PL1151-1803.....	474	PL1310-0048.....	529
G4240-87304.....	100	PL1013-2120.....	511	PL1151-3802.....	474	PL1312-1300.....	393, 463
G4240-87309.....	100	PL1013-2300.....	511	PL1151-3803.....	474	PL1312-1301.....	463
G4240-87310.....	100	PL1013-2500.....	511	PL1170-1810.....	342	PL1312-1500.....	393, 463
G4280-60031.....	48	PL1020-2830.....	511	PL1170-1820.....	342	PL1312-1501.....	463
G4280-60033.....	48	PL1049-2800.....	511	PL1170-1830.....	342	PL1312-1502.....	393, 463
G4280-60061.....	48	PL1110-1120.....	500, 502	PL1170-1840.....	342	PL1312-1503.....	463
G4280-68710.....	60, 66	PL1110-1220.....	500, 502	PL1170-1850.....	342	PL1312-1802.....	393, 463
G4280-68730.....	60, 66	PL1110-1320.....	500	PL1170-1860.....	342	PL1312-1803.....	463
G4280-68750.....	60, 66	PL1110-1400.....	508	PL1170-2820.....	342	PL1312-3300.....	393, 463
G4280-68770.....	60, 66	PL1110-1520.....	500	PL1170-2823.....	342	PL1345-1502.....	411, 463
G4296-68715.....	11	PL1110-6100.....	500	PL1170-6810.....	342	PL1345-1503.....	411, 463
G5611-21503.....	51, 64, 101-102	PL1110-6100LS.....	502	PL1170-6820.....	342	PL1351-1502.....	409, 463
G5611-26210.....	52, 101	PL1110-6115.....	505	PL1170-6830.....	342	PL1351-1503.....	409, 463
G5611-60020.....	49, 101	PL1110-6120.....	505	PL1170-6840.....	342	PL1410-0101.....	529
G5611-60025.....	49, 101	PL1110-6125.....	505	PL1170-6860.....	342	PL1410-0200.....	529
G5611-60061.....	48, 101	PL1110-6130.....	505	PL1171-1140.....	342	PL1410-0301.....	529
G5611-60067.....	49, 101	PL1110-6140.....	505	PL1171-6140.....	342	PL1410-0501.....	529
G5611-60500.....	30	PL1110-6150.....	505	PL1210-1120.....	506	PL1412-4100.....	471
G5611-60501.....	30	PL1110-6160.....	505	PL1210-3120EPA.....	507	PL1412-4101.....	471
G5611-60502.....	30	PL1110-6200.....	500	PL1210-6100.....	506	PL1412-4102.....	471
G5611-60503.....	30	PL1110-6200LS.....	502	PL1210-6104.....	506	PL1412-4103.....	471
G5611-63010.....	101	PL1110-6300.....	500	PL1210-6115.....	506	PL1412-4200.....	471

PART NUMBER INDEX

PL1412-4201	471	PL1512-3702	470	PL1712-3703	470	PL2010-0104	537
PL1412-4400	471	PL1512-3703	470	PL1712-3800	470	PL2010-0105	537
PL1412-4401	471	PL1512-3800	310, 393, 470	PL1712-3801	470	PL2010-0200	535
PL1412-4702	471	PL1512-3801	310, 393, 470	PL1712-6800	470	PL2010-0201	535
PL1412-4703	471	PL1512-3802	310, 393	PL1712-6801	470	PL2010-0202	535
PL1412-4A05	475	PL1512-3803	310, 393	PL1745-3102	474	PL2010-0203	535
PL1412-4K00	471	PL1512-5100	470	PL1745-3103	411, 474	PL2010-0300	535
PL1412-4K01	471	PL1512-5101	470	PL1745-3702	474	PL2010-0301	535
PL1412-4K02	471	PL1512-5102	470	PL1745-3703	411, 474	PL2010-0400	535
PL1412-6100	471	PL1512-5103	470	PL1751-3102	409, 474	PL2010-0401	535
PL1412-6101	471	PL1512-5200	470	PL1751-3103	409, 474	PL2010-0402	535
PL1412-6102	471	PL1512-5201	470	PL1751-3702	409, 474	PL2010-0403	535
PL1412-6103	471	PL1512-5400	470	PL1751-3703	409, 474	PL2010-0501	536
PL1412-6200	471	PL1512-5401	470	PL1812-3102	470	PL2010-0505	536
PL1412-6201	471	PL1512-5500	310, 393	PL1812-3103	470	PL2010-0601	536
PL1412-6400	471	PL1512-5501	310, 393	PL1812-6100	470	PL2010-0605	536
PL1412-6401	471	PL1512-5702	470	PL1812-6101	470	PL2012-0001	538
PL1412-6702	471	PL1512-5703	470	PL1812-6200	470	PL2012-0005	538
PL1412-6703	471	PL1512-5800	310, 393, 470	PL1812-6201	470	PL2012-0010	538
PL1412-6800	471	PL1512-5801	310, 393, 470	PL1812-6400	470	PL2012-1001	538
PL1412-6801	471	PL1512-5802	310, 393	PL1812-6401	470	PL2012-1005	538
PL1412-6A05	475	PL1512-5A05	475	PL1812-6800	470	PL2012-1010	538
PL1412-6K00	471	PL1514-1900HFIP	509	PL1812-6801	470	PL2012-2001	538
PL1412-6K01	471	PL1514-5900HFIP	509	PL1845-2102	474	PL2012-2005	538
PL1412-6K02	471	PL1545-1502	411	PL1845-2103	411, 474	PL2012-2010	538
PL1417-0800	512	PL1545-1503	411	PL1845-3102	411, 474	PL2012-3001	538
PL1417-0830	512	PL1545-1802	411	PL1845-3103	411, 474	PL2012-3005	538
PL1445-4102	411, 474	PL1545-1803	411	PL1851-2102	409, 474	PL2012-3010	538
PL1445-4103	474	PL1545-3102	411, 474	PL1851-2103	409, 474	PL2012-4001	538
PL1445-4702	411, 474	PL1545-3103	411, 474	PL1851-3102	409, 474	PL2012-4005	538
PL1445-4703	411, 474	PL1545-3702	411, 474	PL1851-3103	409, 474	PL2012-4010	538
PL1445-6102	411, 474	PL1545-3703	411, 474	PL1912-1300	310, 393	PL2012-5001	538
PL1445-6103	411, 474	PL1545-3802	411	PL1912-1301	310, 393	PL2012-5005	538
PL1445-6702	411, 474	PL1545-3803	411	PL1912-1500	310, 393	PL2012-5010	538
PL1445-6703	411, 474	PL1545-5102	411, 474	PL1912-1501	310, 393	PL2012-6001	538
PL1451-4102	409, 474	PL1545-5103	411, 474	PL1912-1502	310, 393	PL2012-6005	538
PL1451-4103	409, 474	PL1545-5702	474	PL1912-1503	310, 393	PL2012-6010	538
PL1451-4702	409, 474	PL1545-5703	411, 474	PL1912-1801	310, 393	PL2012-7001	538
PL1451-4703	409, 474	PL1547-1501	428	PL1912-1802	310, 393	PL2012-7005	538
PL1451-6102	409, 474	PL1547-5501	428	PL1912-1803	310, 393	PL2012-7010	538
PL1451-6103	409, 474	PL1551-1502	409	PL1912-3300	310, 393	PL2012-8001	538
PL1451-6702	409, 474	PL1551-1503	409	PL1912-3301	310, 393	PL2012-8005	538
PL1451-6703	409, 474	PL1551-1802	409	PL1912-3500	310, 393	PL2012-8010	538
PL1510-1100	504	PL1551-1803	409	PL1912-3501	310, 393	PL2012-9001	538
PL1510-1200	504	PL1551-3102	409, 474	PL1912-3801	310, 393	PL2012-9005	538
PL1510-1300	504	PL1551-3103	409, 474	PL1912-3802	310, 393	PL2012-9010	538
PL1510-1500	504	PL1551-3702	409, 474	PL1912-3803	310, 393	PL2013-1001	538
PL1510-1504	504	PL1551-3703	409, 474	PL1912-5500	310, 393	PL2013-1005	538
PL1510-5100	504	PL1551-3802	409	PL1912-5501	310, 393	PL2013-1010	538
PL1510-5200	504	PL1551-3803	409	PL1912-5801	310, 393	PL2013-2001	538
PL1510-5300	504	PL1551-5102	409, 474	PL1945-1502	411	PL2013-2005	538
PL1510-5500	504	PL1551-5103	409, 474	PL1945-1503	411	PL2013-2010	538
PL1510-5504	504	PL1551-5702	409, 474	PL1945-1802	411	PL2013-3001	538
PL1512-1300	310, 393	PL1551-5703	409, 474	PL1945-1803	411	PL2013-3005	538
PL1512-1301	310, 393	PL1570-5810	342	PL1945-3802	411	PL2013-3010	538
PL1512-1500	310, 393	PL1612-1801	310, 393	PL1945-3803	411	PL2013-4001	538
PL1512-1501	310, 393	PL1670-0810	342	PL1951-1502	409	PL2013-4005	538
PL1512-1502	310, 393	PL1670-0820	342	PL1951-1503	409	PL2013-4010	538
PL1512-1503	310, 393	PL1670-0830	342	PL1951-1802	409	PL2013-5001	538
PL1512-1801	310, 393	PL1670-0840	342	PL1951-1803	409	PL2013-5005	538
PL1512-1802	310, 393	PL1670-0850	342	PL1951-3802	409	PL2013-5010	538
PL1512-1803	310, 393	PL1670-0860	342	PL1951-3803	409	PL2013-6001	538
PL1512-3100	470	PL1671-0140	342	PL1C12-2502	393	PL2013-6005	538
PL1512-3101	470	PL1712-3100	470	PL1E10-3120EPA	507	PL2013-6010	538
PL1512-3102	470	PL1712-3101	470	PL1E10-6120EPA	507	PL2013-7001	538
PL1512-3103	470	PL1712-3102	470	PL1E12-5A05	475	PL2013-7005	538
PL1512-3200	470	PL1712-3103	470	PL1F70-6830	342	PL2013-7010	538
PL1512-3201	470	PL1712-3200	470	PL1F70-6850	342	PL2013-8001	538
PL1512-3300	310, 393	PL1712-3201	470	PL2010-0100	537	PL2013-8005	538
PL1512-3301	310, 393	PL1712-3400	470	PL2010-0101	537	PL2013-8010	538
PL1512-3401	470	PL1712-3401	470	PL2010-0102	537	PL2013-9001	538
PL1512-3501	310, 393	PL1712-3702	470	PL2010-0103	537	PL2013-9005	538

PL2013-9010.....	538	PL2023-6005.....	540	PL2083-2010.....	542	PL2144-2000.....	545
PL2014-0001.....	538	PL2023-6010.....	540	PL2083-3001.....	542	PL2144-2001.....	545
PL2014-0005.....	538	PL2023-7001.....	540	PL2083-3005.....	542	PL2144-3000.....	545
PL2014-0010.....	538	PL2023-7005.....	540	PL2083-3010.....	542	PL2144-3001.....	545
PL2014-1001.....	538	PL2023-7010.....	540	PL2083-4001.....	542	PL3540-C603VP.....	476
PL2014-1005.....	538	PL2023-8001.....	540	PL2083-4005.....	542	PL3540-D603VP.....	476
PL2014-1010.....	538	PL2023-8005.....	540	PL2083-4010.....	542	PL3540-P603VP.....	476
PL2014-2001.....	538	PL2023-8010.....	540	PL2083-5001.....	542	PL3549-3603VP.....	476
PL2014-2005.....	538	PL2023-9001.....	540	PL2083-5005.....	542	PL3554-1602dAbz.....	347
PL2014-2010.....	538	PL2023-9005.....	540	PL2083-5010.....	542	PL3554-1602dCac.....	347
PL2014-3001.....	538	PL2023-9010.....	540	PL2083-6001.....	542	PL3554-1602dCbz.....	347
PL2014-3005.....	538	PL2024-0001.....	540	PL2083-6005.....	542	PL3554-1602dGdmf.....	347
PL2014-3010.....	538	PL2024-0005.....	540	PL2083-6010.....	542	PL3554-1602dGibu.....	347
PL2014-4001.....	538	PL2024-0010.....	540	PL2083-7001.....	542	PL3554-1602dT.....	347
PL2014-4005.....	538	PL2024-1001.....	540	PL2083-7005.....	542	PL3554-4602dAbz.....	347
PL2014-4010.....	538	PL2024-1005.....	540	PL2083-7010.....	542	PL3554-4602dCac.....	347
PL2014-6001.....	538	PL2024-1010.....	540	PL2083-8001.....	542	PL3554-4602dCbz.....	347
PL2014-6005.....	538	PL2024-2001.....	540	PL2083-8005.....	542	PL3554-4602dGdmf.....	347
PL2014-6010.....	538	PL2024-2005.....	540	PL2083-8010.....	542	PL3554-4602dGibu.....	347
PL2014-7001.....	538	PL2024-2010.....	540	PL2083-9001.....	542	PL3554-4602dT.....	347
PL2014-7005.....	538	PL2070-0200.....	535	PL2083-9005.....	542	RMSN-2.....	108
PL2014-7010.....	538	PL2070-0201.....	535	PL2083-9010.....	542	RMSN-4.....	108
PL2014-8001.....	538	PL2070-0202.....	535	PL2084-0001.....	542		
PL2014-8005.....	538	PL2070-0203.....	535	PL2084-0005.....	542		
PL2014-8010.....	538	PL2070-1001.....	542	PL2084-0010.....	542		
PL2014-9001.....	538	PL2070-1005.....	542	PL2084-1001.....	542		
PL2014-9005.....	538	PL2070-1010.....	542	PL2084-1005.....	542		
PL2014-9010.....	538	PL2070-2001.....	542	PL2084-1010.....	542		
PL2020-0200.....	535	PL2070-2005.....	542	PL2084-2001.....	542		
PL2020-0201.....	535	PL2070-2010.....	542	PL2084-2005.....	542		
PL2020-0202.....	535	PL2070-3001.....	542	PL2084-2010.....	542		
PL2020-0203.....	535	PL2070-3005.....	542	PL2090-1000.....	544		
PL2022-2001.....	540	PL2070-3010.....	542	PL2090-3000.....	544		
PL2022-2005.....	540	PL2070-4001.....	542	PL2090-4000.....	544		
PL2022-2010.....	540	PL2070-4005.....	542	PL2090-5000.....	544		
PL2022-3001.....	540	PL2070-4010.....	542	PL2090-6000.....	544		
PL2022-3005.....	540	PL2070-5001.....	542	PL2090-8000.....	544		
PL2022-3010.....	540	PL2070-5005.....	542	PL2091-1000.....	544		
PL2022-5001.....	540	PL2070-5010.....	542	PL2091-2000.....	544		
PL2022-5005.....	540	PL2070-6001.....	542	PL2091-3000.....	544		
PL2022-5010.....	540	PL2070-6005.....	542	PL2091-4000.....	544		
PL2022-6001.....	540	PL2070-6010.....	542	PL2142-3000.....	545		
PL2022-6005.....	540	PL2070-7001.....	542	PL2142-3001.....	545		
PL2022-6010.....	540	PL2070-7005.....	542	PL2142-5000.....	545		
PL2022-7001.....	540	PL2070-7010.....	542	PL2142-6000.....	545		
PL2022-7005.....	540	PL2070-8001.....	542	PL2142-6001.....	545		
PL2022-7010.....	540	PL2070-8005.....	542	PL2142-7000.....	545		
PL2022-8001.....	540	PL2070-8010.....	542	PL2142-7001.....	545		
PL2022-8005.....	540	PL2070-9001.....	542	PL2142-8000.....	545		
PL2022-8010.....	540	PL2070-9005.....	542	PL2142-8001.....	545		
PL2022-9001.....	540	PL2070-9010.....	542	PL2143-0000.....	545		
PL2022-9005.....	540	PL2071-0001.....	542	PL2143-0101.....	545		
PL2022-9010.....	540	PL2071-0005.....	542	PL2143-2000.....	545		
PL2023-0001.....	540	PL2071-0010.....	542	PL2143-2001.....	545		
PL2023-0005.....	540	PL2071-1001.....	542	PL2143-3000.....	545		
PL2023-0010.....	540	PL2071-1005.....	542	PL2143-3001.....	545		
PL2023-1001.....	540	PL2071-1010.....	542	PL2143-4000.....	545		
PL2023-1005.....	540	PL2071-2001.....	542	PL2143-4001.....	545		
PL2023-1010.....	540	PL2071-2005.....	542	PL2143-5000.....	545		
PL2023-2001.....	540	PL2071-2010.....	542	PL2143-5001.....	545		
PL2023-2005.....	540	PL2071-3001.....	542	PL2143-6000.....	545		
PL2023-2010.....	540	PL2071-3005.....	542	PL2143-6001.....	545		
PL2023-3001.....	540	PL2071-3010.....	542	PL2143-7000.....	545		
PL2023-3005.....	540	PL2080-0200.....	535	PL2143-7001.....	545		
PL2023-3010.....	540	PL2080-0201.....	535	PL2143-8000.....	545		
PL2023-4001.....	540	PL2080-0202.....	535	PL2143-8001.....	545		
PL2023-4005.....	540	PL2080-0203.....	535	PL2143-9000.....	545		
PL2023-4010.....	540	PL2083-1001.....	542	PL2143-9001.....	545		
PL2023-5001.....	540	PL2083-1005.....	542	PL2144-0000.....	545		
PL2023-5005.....	540	PL2083-1010.....	542	PL2144-0101.....	545		
PL2023-5010.....	540	PL2083-2001.....	542	PL2144-1000.....	545		
PL2023-6001.....	540	PL2083-2005.....	542	PL2144-1001.....	545		

Application Title Index

12 phenols analyzed
using a longer (4.6 x 100 mm)
Agilent Poroshell 120 EC-C18 column.....231

A

Adrenocorticosteroids
on Pursuit PFP and C18290

Agilent Bio SEC-3
column length comparison, 150 mm.....422

Agilent Bio SEC-3
column length comparison, 300 mm.....422

Alberta Peptide Institute test mix.....308, 390

Alcohols and aliphatic compounds.....582

Alkyd resin.....519

Amino acid standard separation
on Eclipse Plus C18.....564

Amitrol in water by LC/MS, 0.05 ppb585

An overlay of the original ZORBAX
Eclipse Plus 5 μ m method and Agilent
Poroshell 120 method. All 11 peaks on
Poroshell 120 are resolved by the time
the first peak elutes on the original 5 μ m
ZORBAX Eclipse Plus method595

Analgesics.....618

Analgesics: Non-steroidal anti-inflammatory
drugs: Narrow bore separation622

Analysis of Biocides in Hand Sanitizer.....580

Analysis of Bovine Serum Albumin by light
scattering using ProSEC 300S columns.....429

Analysis of carbohydrates
on Hi-Plex H columns340

Analysis of choline kinase on PL-SAX 4000Å.....407

Analysis of ciprofloxacin
and ciprofloxacin metabolites.....643

Analysis of diazepam on Rx-C18.....272

Analysis of fruit juice.....338

Analysis of low molecular weight polystyrene
and oligomer fractions collected from
OligoPore preparative columns522

Analysis of oxidized insulin chains.....553

Analysis of pesticide residues in green tea.....594

Analysis of representative whey proteins408

Analysis of sugars with high sodium matrix341

Analysis of sweeteners on Hi-Plex Ca columns.....340

Analysis of water soluble vitamins
in multivitamin tablets612

Analytical separation of
low molecular weight polystyrene522

Anesthetics, local: Bonded phase selectivity.....619

Anilines, substituted: Rapid separation585

Anthocyanins from blueberries:
High-efficiency high-speed separation.....598

Antibiotics: High speed separation620

Antibodies: Fast separation of IgM
and IgG antibodies565

Antidepressants, tricyclic:
Comparative separation636

Antifungal medications621

Antifungals290, 621

Antihistamines: Fast separations
on RRHT Extend-C18617

Antibiotics: Lincomycin and Clindamycin
by LC-APCI-MS LC-TIC620

Aromatic acids/benzoic acids:
Selectivity differences623

Aromatics II.....599

Aspartame: Metabolites and applications.....599

Aspirin and cough remedy628

B

BSA tryptic digest on RRHT.....577

Barbiturates.....643

Baseline expansion of a separation
of protein standards.....414

Basic antihistamines
on Extend-C18 at high pH.....275

Benzoic acid/sorbic acid609

Bio-Monolith DEAE column monitors
phage production during fermentation.....415

Blueberry anthocyanin analysis594

Brij 35581

C

Calibration curves.....505

Calibration of the ProSEC 300S
column with globular proteins429

Capillary columns for HPLC
analyses with UV and MS detection455

Carbohydrates in colas.....601

Carbohydrates in juices.....601

Carbohydrates in milk.....602

Carbohydrates: Carbohydrate standards.....600

Carbohydrates:
Effect of mobile phase strength.....600

Carbohydrates: Sugar alcohols.....601

Catecholamines.....578

Catecholamines/biogenic amines:
Rapid separation using ion-pair reagents.....623

Charge isoform analysis
of monoclonal antibodies.....401

Chiral ethiazide (diuretic drug) separation624

Chiral separation of S- and R-Norfluoxetine626

Chiral separation of atenolol.....627

Chiral separation of
fluoxetine enantiomers (Prozac)624

Chiral separation of hexobarbital.....625

Chiral separation of salbutamol626

Chiral separation of tolperison enantiomers627

Cocaine and metabolites627

Column reproducibility – 200 injections of
reduced monoclonal antibody using an
Agilent ZORBAX RRHD 300SB-C3 column.....555

Column stability testing at pH 3 and 60 °C.....257

Column stability testing at pH 7.0.....258

Columns for sample clean-up.....507

Comparing HILIC and RPLC of morphine
using Agilent ZORBAX RRHD columns
with UHPLC/MS631

Comparison of A β peptide RP-HPLC
separations at low and high pH.....573

Comparison of Agilent Bio SEC-3
and competitor column in the analysis
of a monoclonal antibody.....420

Comparison of PolyPore
with conventional individual
pore size GPC columns517

Comparison of phenols separation
with Poroshell 120.....584

Components of green tea separated
on Rapid Resolution StableBond SB-C8625

Consistent ion-exchange
MAb separation400, 562

Corn syrup, Hi-Plex.....339

Cough formula mixture:
Fast and efficient separation628

Cough/cold remedies on ZORBAX 300SCX334

Crude bradykinin prep load469

Crude peptide screen.....476

D

DNPH: Derivatized Aldehydes
obtained from air.....584

Deoxynucleosides: Using rapid resolution
3.5 μ m columns.....577

Determination of anthocyanins in blueberries597

Dexamethasone, USP method: Rapid analysis.....639

Differences in composition of
two alkyl naphthalene sulfonates527

Dimethyl-C18/amide, Bonus-RP280

E

EasiVial PS-H533

Eclipse Plus C18 vs. C8.....250

Eclipse Plus C8 is less retentive
than Eclipse Plus C18.....251

Eliminate tailing and maximize resolution
with Eclipse Plus Columns.....250

Environmental phenols on Poroshell 120.....230

Epoxy resin501

Excellent separation of two phenol
formaldehyde resins with PolarGel-M.....513

Exceptional Lot-to-lot Reproducibility.....426

Exceptional separating power403

Exploiting chemical stability –
NH₄OH concentration.....308, 390

Explosives and related compounds:
Qualitative and quantitative analysis.....586

Explosives from soil extract.....587

Extend-C18 provides good peak shape
at low pH276

F

Fast analysis 11 common compounds
found in analgesics615

Fast analysis of cefepime
and related impurities217, 611

Fast analysis of Pindolol642

Fast analysis of sulfa drugs.....596

Fast and ultra-fast analysis
of basic compounds251, 617

Fast LC/MS/MS analysis of group 4 pharmaceuticals from EPA 1694	583
Fast method for ginseng analyses scaled from a traditional method	612
Fast separation of recombinant human erythropoietin	553
Fast separation of reduced monoclonal antibody	375
Fast vitamin E analysis on Rapid Resolution HT	608
Fast, high-resolution separation of peptides and proteins with Poroshell 300SB-C18	572
Faster analysis of USP Method for simvastatin tablet	615
Faster separation of sulfa drugs	616
Faster separations using Agilent weak cation-exchange columns	560
Fat-soluble vitamins on ZORBAX Eclipse XDB-C8	605
Five different bonded phases provide selectivity options	266
Flavoring agents	602
Food colors, FD&C	602
Four different 300SB bonded phases optimize separation of large polypeptides	371
Fraction analysis – the concentration overload purification	469

G

Glycosylated proteins: Large Molecules on Poroshell 300SB-C18 and 300SB-	565
Goldenseal and related alkaloids on Rapid Resolution Eclipse XDB-C18	625
Good peak shape over a wide pH range with ZORBAX Eclipse XDB	257
Gradient optimizations for ultra-fast analysis of reduced monoclonal antibody	556
Guaifenesin: USP analysis of guaifenesin	629

H

HPLC of 25 bp DNA ladder	307, 388
HPLC separation of 12 phenols performed in just 5 minutes – and under 400 bar – using an Agilent Poroshell 120 EC-C18 column	231
HSA tryptic digest on ZORBAX Rapid Resolution HT 1.8 µm	566
Heparin	526
Herbicide/pesticide standards: Effect of bonded phase	588
Herbicides on different bonded phases	587
Herbicides: Rapid separation	591
High Resolution of 24 Amino Acids Using ZORBAX Eclipse AAA Protocol	396
High purity and high recovery with ZORBAX PrepHT columns	315
High resolution and fast analysis on RRHT Eclipse PAH column	254
High resolution normal-phase separation of octylphenoxy ethanol surfactant on ZORBAX CN	327
High resolution separation of a Poly-T-Oligonucleotide size standard spiked with 10-mer, 15-mer, 30-mer and 50-mer (main peaks)	408

High sensitivity with capillary columns	453
Higher resolution of intact monoclonal antibody	368
Higher resolution of oxidation study	354, 368
Hormones/steroids	633
Human serum: Low abundance protein isolation and identification from 1-D gel band by LC/MS	454, 566
Hyaluronic acid	527
Hydrophilic purine/pyrimidine separation	579

I

Ibuprofen: Optimizing selectivity with RRHT Columns	618
Improved peak shape of basic compounds using Bonus-RP	279
Improved reproducibility of monoclonal antibodies	369
Improved resolution with smaller particle size with Agilent weak cation-exchange columns	559
Increase peak capacity with RRHT columns	240
Increased resolution for peptide mapping	357, 369
Intact MAb monomer and dimer separation	355, 419, 563

L

LC/MS analysis of angiotensin on Extend-C18	377
LC/MS performance test mix for Polaris C8-A	299
Lamotrigine	642
Large fibrous proteins	309, 392
Liquid chromatography phase test mixture (LPTM) on Pursuit C8	290
Local anesthetics	619
Long life at high pH with 300Extend-C18	378
Long life at high pH with Extend-C18	275
Long lifetime of RRHT columns at elevated temperatures	241

M

Mechanical stability of Pursuit XRs	289
Metronidazole: Updating USP methods	629
Microbore HPLC for sensitive peptide analysis	462
MicroBore Poroshell 300 columns provide maximum sensitivity for LC/MS	382
Mixture of beta blockers	613
Monoclonal IgG1 chains: Separation on Poroshell 300SB-C8	383, 567
Morphine and metabolites: Extracted blood plasma sample separation	630

N

Naproxen analysis	611
Natural products – capsaicin and dihydrocapsaicin on Pursuit XRs C18	319
Neutraceuticals: Extract from green tea	603
Neutraceuticals: Hypericin separation in St. John's Wort	632
New levels of sensitivity and resolution	235
Nucleosides, purines and pyrimidines	564
Nucleosides: Separation of deoxy and ribonucleosides	568
Nucleotides: Separation of mononucleotides	568

O

Oligosaccharides	610
Opiates (drugs of abuse) by LC/MS	630
Optimize separations with Eclipse XDB selectivity options	259
Optimizing protein separations with Agilent weak cation-exchange columns	558
Organic acid analysis	338
Organic acids separated on ZORBAX SB-Aq	581
Overlay of UV and Light Scattering 90° for a Sample of γ-globulins, Illustrating Monomer, Dimer, Trimer, and Aggregate Peaks	430
Overlay of UV and light scattering 90° for a sample of γ-globulins, illustrating monomer, dimer, and trimer peaks	430

P

pH gradient elution for improved separation of monoclonal antibody charged variants	561
PLgel LS column	503
PLgel Olexis reveals true modalities across the range of polyolefins	508
Parabens: High speed separation	604
Peak shape and efficiency are better with ZORBAX Eclipse Plus	249
Peptide RP-HPLC/ESI-MS using NH4OH mobile phase yields both positive and negative ion spectra	573
Peptides/proteins: Effect of elevated temperature	370, 571
Peptides/proteins: Equivalent gradient separations	571
Peptides: Effect of TFA concentration	370, 569
Peptides: Separation of antitensins I, II, III with TFA and NH4OH	570
Peptides: Separation of antitensins I, II, III with TFA and NH4OH	632
Phenols, substituted	592
Phenoxyacid herbicides	591
Plant hormones: Rapid gradient elution separation	593
Plasticized PVC	501
Polyamides	509
Polycyclic aromatic hydrocarbons according to EPA Method 610	589
Polyester	519
Polyesterimide	520
Polyethylene Glycol/Oxide standards	542
Polyethylene glycols	307, 388
Polymethylmethacrylate in DMF	517
Polymethylmethacrylate standards	540
Polyol	501
Polyphenylene Sulfides	500
Polystyrene standards	538
Polyurethanes	520
Polyvinyl alcohol	526, 528
Pore size choice: Mouse IgG	421
Pore size choice: Proteins	421
Poroshell 120 EC-C18 for fast UHPLC separations	232
Poroshell 300 columns separate proteins and peptides in seconds	381

Preparative fractionation of a culture filtrate containing amyloglucosidases on Agilent PL-SAX 4000A	473
Preparative scale purification of Leuprolide by concentration overload	359
Preparative separation of low molecular weight polystyrene	522
Protein digest analysis	552
Protein elution pattern on ZORBAX Poroshell 300SB-C8	384
Proteins in a complex sample by 2-D HPLC with Nano HPLC columns	457
Proteins: Effect of bonded phase	575
Proteins: Effect of bonded phase, RP	575
Pullulan polysaccharide standards	544
Purification of a 25-mer trityl-off oligonucleotide and analytical quantitation of the fraction using PLRP-S 100Å, 4.6 x 50 mm	469
Purification of a large oligonucleotide	473

Q

Quantification and qualification of vitamin C and citric acid in fresh grapefruit juice	609
---	-----

R

Rapid analysis of an analgesic tablet, selectivity differences at pH 2.7 and pH 7	250
Rapid method development for 18 compounds with an Agilent RRHD Eclipse PAH column	590
Rapid Resolution HT (RRHT) provides double the efficiency of Rapid Resolution columns	240
Reduce analysis time dramatically with Rapid Resolution HT columns	241
Reduce peptide map analysis time by 90% with Poroshell 300SB	382
Resin analysis by rapid GPC	511
Rose wine	610

S

SB-CN optimizes retention and resolution	266
Scale-up from analytical to prep	321
Selectivity changes for basic compounds with Eclipse XDB and StableBond	258
Selectivity comparison of TFA and NH ₄ OH for peptide RP-HPLC\ESI-MS analysis	574
Selectivity comparison of TFA and NH ₄ OH for peptide RP-HPLC\ESI-MS analysis	574
Selectivity comparison: C18 columns	236
Selectivity comparison: Phenyl columns	237
Selectivity for urea pesticides	259
Selectivity in peptide RP-LC	389
Selectivity test mix for Polaris columns	299
Separation of 20 PAHs on Eclipse PAH	591
Separation of 8 steroids	613
Separation of azo dye degradation products	583
Separation of Azo Dyes	598
Separation of EPA 610 PAH Mix	589
Separation of a tryptic digest on ZORBAX MicroBore 300SB-C18	461
Separation of basic peptides on Bonus-RP versus traditional Alkyl phase	569
Separation of charge variants of human IgG1 with pH gradient	355

Separation of group 4 analytes in EPA 1694 on ZORBAX HILIC Plus column	325
Separation of heated, stressed MAb	563
Separation of highly basic antidepressants above their pKa in free base form (pKa 9.5-9.7)	635
Separation of licorice root on RRHD columns	234
Separation of pharmaceutical cardiac drugs	616
Separation of polypeptides in under 1 minute	572
Separation of protein standards on Agilent 3 µm ion-exchange columns by cation-exchange chromatography	403
Separation of recombinant human erythropoietin (rEPO)	562
Separation of small molecule anorectics	622
Separation of vitamin D2/D3	605
Separation optimization for ultra-fast analysis of reduced and alkylated monoclonal antibody	554
Several ZORBAX RRHD 1.8 µm selectivities facilitate method development	614
Short-chain ZORBAX 300SB-C3 is stable at low pH, high temperature	371
Shorter chain ZORBAX SB-CN is also stable at low pH (pH 2.0, 50 °C)	265
Space	265
Sports drink	610
StableBond SB-C18 shows excellent stability at low pH and high temperature	265
Standard ion-exchange protein separation	407, 576
Standard protein separation	410
Standard proteins by reversed-phase	576
Starches	500
Steroids	634
Steroids: Easy scalability using Agilent Prep columns	312
Steroids: Separation	633
Sub 1 minute separations with RRHD columns	235
Sugars	604
Sugars in plain and milk chocolate	604
Sulfa drugs	641
Sulfonamides – Fast analysis with RRHT columns	640
Sunscreen ingredients: Perform conventional, fast and ultra-fast separations on the same column family	607
Superior loadability on Agilent Prep C18 with basic compounds	312

T

Temperature as a tool to enhance mass transfer and improve resolution of oligonucleotides in ion-pair reversed-phase HPLC	391, 579
Ten cardiac drugs on Rapid Resolution HT SB-C18	640
Theobromine in beverages	608
Tocopherols by LC/MS with APPI	603
Triamcinolone – USP analysis of triamcinolone	635
Triazine pesticides on Bonus-RP and Alkyl C8 phase	592

Tricyclic antidepressants	636
Tricyclic antidepressants and benzodiazepines	289
Tricyclic antidepressants and metabolites: Effect of pore size	637
Triton X-114: Decreasing run-time by changing bonded phase	580
Two samples of melamine resin analyzed by PolarGel-L	513

U

UHPLC efficiency at HPLC pressures	230
Ultra-high speed and high resolution of intact monoclonal antibodies	557
USP analysis of tetracyclines	639
USP method: Glyburide and internal standard, progesterone	638
USP method for sorbitol	341
USP methods for sugar alcohols	339
Ulcer treatment drugs at intermediate pH	637
Urine, LSD analysis by LC/MS	638
Use ZORBAX Extend-C18 for alternate selectivity at high pH	378, 567

V

VX nerve agent metabolites by LC/MS-IS standard (C13 labeled)	593
Virtually eliminate retention time variations	400

W

Warfarin: USP chromatographic purity method using Eclipse XDB-CN	640
Water-soluble B vitamins separated on ZORBAX SB-Aq	607
Water-soluble vitamins	605
Water-soluble vitamins using the USP 23 method	606
Water-soluble vitamins: High speed separation using ion-pairing	606
Weak cation-exchange chromatography for P128 therapeutic protein sample on the Agilent 1260 Bio-inert Quaternary LC system using different cation-exchange columns	404
Whey proteins in dairy samples – milk	391, 578

X

Xanthines: Higher resolution, same selectivity with RRHT	617
--	-----

Z

ZORBAX Bio-SCX Series II provides more retention of small peptides	459
ZORBAX Bonus-RP is stable at low and mid pH	279
ZORBAX Bonus-RP provides unique selectivity	280
ZORBAX Eclipse Plus: Best peak shape in the industry without tailing	249
ZORBAX Nano columns for high sensitivity protein digest analysis by LC/MS	453

Compound Index

A

Acebutolol	613
Acenaphthene	589-591
Acenaphthylene	589-591
Acephate	594
Acetaldehyde – DNPH	584
Acetamide	615
2-Acetamidophenol	235, 615
Acetaminophen	235, 250, 615, 618, 628
Acetanilide	235, 618
Acetate	174-175
Acetic Acid	581
Acetone	582
Acetone – DNPH	584
Acetophenone	599
6-Acetylmorphine	628
Acetylsalicylic acid	235, 250, 618, 628
Acrolein – DNPH	584
Adenine	564, 579
Adenosine	564, 568
Adonitol	339
Adrenaline	578
Alanine	396, 564
Alanine-3-Glycine-4	308, 390
Albuterol	325, 583
Alprenolol	613
2-Amino-5-azotoluene	183
4-Aminobenzoic acid	290, 621
p-Aminobenzoic acid	609
3-Amino-benzonitrile	585
7-Aminoclonazepam	289
2-Amino-4,6-dinitrotoluene	586-587
4-Amino-2,6-dinitrotoluene	587
4-Amino-4,6-dinitrotoluene	586
7-Aminoflunitrazepam	289
2-Aminonaphthalene	183
2-Amino-4-nitrotoluene	586-587
2-Amino-6-nitrotoluene	586-587
4-Amino-2-nitrotoluene	586
Amitriptyline	258, 289, 635-637
cis-10-OH -Amitriptyline	637
trans-10-OH-Amitriptyline	637
Ammonium	173
Amylbenzene	299
Amylobarbitone	643
Anadamine	236
Androstadiene 3,17 dione	613
Anethole	602
Angiotensin I	308, 377, 389-390, 462, 570
Angiotensin II	194, 308, 370, 377, 381, 389-390, 462, 570-572
Angiotensin III	308, 377, 389-390, 570
Aniline	183, 583, 598
p-Cl-Aniline	183
Anisidine	598
p-Anisidine	585

Anthracene	189, 589-591
α -1-Antichymotrypsin	454, 566
Antithrombin-III	454, 566
Apomyoglobin	572
Aprotinin	403
Arabinose	341
Arabitol	339
2-Arachinoylglycerol	381
Arginine	396, 564
Arsenate	175
Arsenite	175
Ascorbic acid	612
Asparagine	396, 596
Aspartame	232, 290, 299, 599
Aspartic acid	396, 596
Aspartic acid-phenylalanine dipeptide	599
Atenolol	613
Atrazine	587-588, 591
Azide	172, 175, 421

B

Barberine	625
Barbital	258
Barbitone	643
Beclomethasone	290
Bendroflumethiazide	189
Bentazon	587
Benzaldehyde – DNPH	584
Benz(e)pyrene	591
Benzidine	183, 583, 598
1,2-Benzisothiazol-3(2H)-one	580
Benzisothiazol-3(2H)-one	235
Benzo(a)anthracene	589-591
Benzo(a)pyrene	589-591
Benzo(b)fluoranthene	589-591
Benzocaine	619
Benzo(g,h,i)perylene	589-591
Benzoic acid	232, 290, 389, 580, 602, 609, 621, 629
Benzo(k)fluorene	589-591
Benzophenone	618
Benzoylcegonine	627
Benzthiazuron	588
n-6-Benzyl adenine	593
5-Benzyl-3,6-dioxo-2-piperazineacetic acid	599
Berberine	251, 617
Biosynthetic human insulin	384
Biotin (B7)	612
Biphenyl	189
Bovine carbonic anhydrase	407, 576
Bromide	172
BSA	382, 425-426, 432, 576-577
Bumetanide	189
Buspirone	616
Butacaine	266, 619
tert-butanol	582

2-Butanone (MEK) – DNPH	584
Butylbenzene	299
Butylparaben	189, 241, 604
n-Butyraldehyde – DNPH	584

C

Caffeine	232, 235, 250, 258, 276, 299, 316, 334, 599, 608, 615
Calcium	173
Calmodulin	571
Canadine	625
Carbamazepine	618, 632
Carbaryl	594
Carbendazim	235, 580, 594
Carbonate	175
Carbonic Anhydrase	370-371, 381-382, 408, 571-572
Carvone	602
Catalase	432
Catechol	230-231, 584, 603, 625
Cefazolin	620
Cefepime	217
Cefotaxime	620
Ceftazidime	620
Celecoxib	235
Cephaclo	280
Cephalexin	280
Cephoxitin	280
Cephuroxime	280
Chlorate	172
Chloride	172, 174-175
Chloroaniline	583, 598
m-Chloroaniline	585
o-Chloroaniline	585
p-Chloroaniline	585
2-Chlorobenzoic	623
Chlorocaine	619
4-Chloro-3-methylphenol	592
5-Chloro-2-methyl-4-isothiazolin-3-one	235
2-Chlorophenol	592
o-Chlorophenoxyacetic acid	591, 593
p-Chlorophenoxyacetic acid	591, 593
o-Chlorophenoxy propionic acid	593
Chlorothiazide	189
Chlorpheniramine	275, 279, 617, 628
Chlorthalidone	189
Chrysene	589-591
Chymotrypsinogen A	410
Cimetidine	325, 583, 637
d-Cinchonine	266, 619
Cinnamaldehyde	599
Ciprofloxacin	620
Citrate	174
Citric Acid	338, 396, 581, 609-610
Clindamycin	620
Clomipramine	289

Clonazepam	289
Clotrimazole	621
Cobalamin (B12)	612
Cocaine	619, 627
Codeine	630
Corticosterone	290
Cortisone	290, 299
Cortisone acetate	290
m-Cresol	590
o-Cresol	230-231, 584, 592
p-Cresol	230-231, 584
Crotonaldehyde – DNPH	584
Cyanide	175
Cyanidin	237, 597
Cyanocobalamin (B12)	605
p-Cymene	599
Cyprodinil	594
Cysteine	396
Cytidine	568
Cytochrome C	183, 382, 370-371, 400, 403, 410, 558-560, 571, 572, 576
Cytidine	564
Cytosine	564, 579

D

Daidzen	266
Dehydroacetic acid	232
Delphinidin	237, 597
2' Deoxycytidine	568
2' Deoxyguanosine	568
2' Deoxyinosine	568
Desethylatrazine	588
Desethyldeisopropylatrazine	588
Desipramine	289
Dextromethorphan	628
Diazepam	272, 289
Dibenzo(a,h)anthracene	589-591
Dichlorobenzidine	583, 598
3,3-Dichlorobenzidine	183
2,4-Dichlorophenol	592
2,3-Dichlorophenoxyacetic acid	591
2,4-Dichlorophenoxyacetic acid	591
Diclofenac	235, 615, 618
Dienestrol	633
Diethylstilbestrol	265
Diflusal	235, 615
Dihydroxy benzylamine	623
Dihydroxyphenyl acetic acid	623
Dihydroxyphenyl alanine	623
Diltiazem	616
Dimethoxybenzidine	598
3,3'-Dimethoxybenzidine	583
Dimethylbenzidine	583
2,3-Dimethyl phenol	230-231, 530
2,4-Dimethyl phenol	592
2,5-Dimethyl phenol	230-231, 530
3,4-Dimethyl phenol	231, 584
1,3-Dimethyluric acid	617
Dimethylxanthine	276

1,7-Dimethylxanthine	617
3,7-Dimethylxanthine (theobromine)	617
1,3-Dinitrobenzene	586-587
2,4-Dinitrophenol	592
2,4-Dinitrotoluene	586-587
2,6-Dinitrotoluene	586-587
Dioctyl phthalate	290, 299
Diphenhydramine	275, 617, 632
Diphenylamine	586
Dipropyl phthalate	241
Dipropylthallate	250
Dipyridamole	616
Disopyramide	640
Diuron	259, 587-588, 591
Dopamine	578, 623
Doxepin	258, 289, 635-636
Doxycycline	639
Doxylamine	275, 279, 617
Dulcitol	339

E

Epinephrine	623
Ecgonine methylester	627
Econazole	621
Eletriptan	616
Epagallo catechin	603
Epicatechin	603, 625
Epicatechin gallate	603, 625
Epigallocatechin	625
Epigallocatechin gallate	603, 625
Estradiol	633
β -Estradiol	613
Estriol	266
Estrone	613
Ethanol	582, 610
Ethinylestradiol	633
Ethoprophos	594
Ethyl cinnamate	599
2-Ethylhexyl <i>trans</i> -4-methoxycinnamate	607
<i>bis</i> -(2-Ethylhexyl) phthalate	507
Ethylhexyl salicylate	259
2-Ethylhexyl salicylate	607
Ethylparaben	189, 272, 604
Eugenol	599
Excipient	606

F

Famotidine	637
Fenfluramine	622
Fenoprofen	622
Fenuron	259
Fibrinogen	309, 392
Flunitrazepam	289
Fluocinolone acetate	290
Fluoranthene	189, 589-591
Fluorene	189, 589-591
Fluoride	172, 175
2-Fluorobenzoic	623

3-Fluorobenzoic	623
Fluorocytosine	579
Folic acid	605-606, 612
Formaldehyde – DNPH	584
Formate	174
Fructose	338, 340, 600, 604, 610
Fumaric Acid	581
Furosemide	616

G

Galactose	340
Genistein	266
γ -Globulin	421
Glucagon	384
Glucose	338, 340-341, 600, 604, 610
Glutamine	396, 564
Glutamate	175
Gly3-Gly4 (N α -acetylated)	308, 390
Glyburide	638
Glycerol	610
Glyceryl Guaicolate	334
Gluconate	396, 564
Glycine	389, 396, 564
Guaifenesin	629
Guanine	564, 579
Guanosine	564, 568

H

Heptabarbital	643
Hexaldehyde – DNPH	584
Hexazinon	588
Hexogen (RDX)	586-587
Hexyl	586
Histidine	396, 564
Homocyclonite	586
Homovanillic acid	623
Holotransferrin	572
Hydrastine	625
Hydrochlorothiazide	189
Hydrocortisone	312, 613, 635
Hydroflumethiazide	189
Hydroquinone	230-231, 584
Hydroxyindoleacetic acid	623
2-Hydroxy-4-methoxybenzophenone	607
4-Hydroxypropivacaine	619
Hydroxyproline	396, 564
5-Hydroxytryptamine	623

I

Ibuprofen	235, 396, 618
IgA	421
IgG2a, I HOPC-1	565
IgM, MOPC-104E	565
Imazalil	594
Imidacloprid	594
Imipramine	289, 632, 635
Indeno(1,2,3-c,d)pyrene	589-591

Indole.....	642
3-Indole acetic acid.....	593
3-Indole butyric acid.....	593
3-Indole propionic acid.....	593
Inosine.....	568
Insulin.....	370-371, 381-382, 462, 571-572
Iodide.....	172
<i>iso</i> -erythritol.....	339
Isoleucine.....	564
Isomaltose.....	600
Isopropanol.....	582
Isoproturon.....	588

K

Kathon 1A.....	580
Kathon 1B.....	580
Ketoprofen.....	235, 615
Kinetin.....	593
Kresoxim-methyl.....	594

L

α -Lactalbumin.....	391, 408, 578
Lactate.....	174-175
Lactic Acid.....	581, 338, 610
Lactoglobulin A.....	183
Lactoglobulin B.....	183
β -Lactoglobulin (A chain).....	391, 408, 578
β -Lactoglobulin (B chain).....	391, 408, 432, 578
Lactose.....	340, 600, 602
Lamotrigine.....	642
Leucine.....	396, 564
Leucine Enkephalin.....	194, 370, 571-572
Lidocaine.....	266, 619, 640
Lincomycin.....	620
Linuron.....	259, 588
Lysine.....	396, 564
Lysozyme.....	183, 370-371, 381-382, 400, 403, 410, 558-560, 571-572, 576

M

Magnesium.....	173
Malate.....	174
Maleate.....	275
Maleic acid.....	628
Malic acid.....	610
Maltose.....	600
Maltotriose.....	340
Malvidin.....	237, 597
Mannitol.....	339-341, 604
Mefanamic acid.....	241
Mepivacaine.....	619
Metacycline.....	639
Metazachlor.....	588
Met-Enkephalin.....	194, 571-572
Metformin.....	325, 583
Methabenzthiazuron.....	588
Methacrolein -DHCP.....	584

Methanol.....	582
3-Methoxytyrosine.....	623
Methionine.....	396, 564
Methoxyaniline.....	583
4-Methoxybenzenesulfonamide.....	289
Methoxychlor, 200 mg/L.....	507
Methyl-3-aminothiophene-2-carboxylate.....	289
2-Methyl-4,6-dinitrophenol.....	592
4,4-Methylene-bis-2-chloroaniline.....	183
2-Methyl-4-isothiazolin-3-one.....	235
1-Methyl naphthalene.....	591
2-Methyl naphthalene.....	591
2-Methyl-5-nitroaniline.....	183
Methyl paraben.....	189, 235, 580, 604
Methyl prednisolone.....	290
Methyl salicylate.....	602
1-Methylxanthine.....	617
Metobromuron.....	588
Metolachlor.....	587-588
Metoprolol.....	613
Metoxuron.....	588
Miconazole.....	621
Molybdate (VI).....	172
Monolinuron.....	259, 588
Monuron.....	259
Morphine.....	630-631
Morphine-6-glucuronide.....	630
Morphine-3-glucuronide.....	630
Myoglobin.....	370-371, 381-382, 407, 410, 421, 432, 571-572, 576

N

Nadolol.....	613, 640
Naphthalene.....	189, 241, 589-591
1-Naphthol.....	230-231, 584
1-Naphthyl acetamide.....	593
1-Naphthyl acetic acid.....	593
Naphthylamine.....	583, 598
Naproxen.....	241, 618
Neurotensin.....	381, 571-572
Niacin (B3).....	605-606, 612
Niacinamide.....	607
Nicotinic Acid.....	607
Nifedipine.....	315, 640
Nimodipine.....	315, 640
Nisoldipine.....	315, 640
Nitrate.....	172, 175
Nitrite.....	172, 175
2-Nitrobenzoic.....	623
3-Nitrobenzoic.....	623
2-Nitrophenol.....	230-231, 584, 592
4-Nitrophenol.....	230-231, 584, 592
2-Nitrotoluene.....	586-587
3-Nitrotoluene.....	586-587
4-Nitrotoluene.....	586-587
Noradrenaline.....	578
Nordiazepam.....	289
Nordoxepin.....	289
Norepinephrine.....	623

Norethindrone.....	633
Norethindrone acetate.....	613
Normorphine.....	631
Nortriptyline.....	250, 258, 289, 635-637
<i>cis</i> -10-OH-Nortriptyline.....	637
<i>trans</i> -10-OH-Nortriptyline.....	637
Norvaline.....	396, 564

O

Octogen (HMX).....	587
Octylmethoxycinnamate.....	250
Oleylethanolamide (OEA).....	328, 614
Orotic Acid.....	568
Ovalbumin.....	381, 407, 412, 421, 558-560, 572, 576
Oxalate.....	174
Oxalic acid.....	338, 389
Oxybenzone.....	250, 259
Oxytetracycline.....	639
Oxytocin.....	462

P

Padimate-O.....	259, 607
Palatinose.....	600
Palmatine.....	625
Palmitoylethanolamide.....	236, 614
Pantothenic acid.....	605, 612
Parvalbumin.....	371
Penconazole.....	594
Pencycuron.....	259
Pentachlorophenol.....	592
<i>iso</i> -Pentane.....	565
Pentylparaben.....	189
Peonidin.....	237, 597
Perphenazine.....	632
Perylene.....	507
Petunidin.....	237, 597
Phenacetin.....	235, 615, 618, 622
Phenanthrene.....	189, 589-591
Phenobarbitone.....	643
Phenol.....	231, 240, 249, 299, 584, 592
Phenoxyacetic acid.....	591
2-Phenoxyethanol.....	235, 580
Phentermine.....	622
Phenylalanine (PHE).....	396, 564, 599
Phenylbutazone.....	622
Phenylephrine.....	334
Phosphate.....	174-175
Phthalic acid.....	623
<i>p</i> -hydroxybenzoic acid.....	232
Picric acid.....	586
Pindolol.....	613, 640, 642
Pioglitazone.....	616
Pirenzepine.....	637
Piroxicam.....	235, 615
Poly-DL-alanine.....	425
Potassium.....	173

Prednisolone	290, 632
Prednisolone acetate	290
Procainamide	616, 640
Procaine	266, 619, 640
Progesterone	613, 638
Promethazine	632
Prometon	587
Prométryne	587, 591
Propanil	587, 591
n-Propanol	582
Propanolol	636
Propazine	587-588, 591
Propionaldehyde – DNPH	584
Propoxur	594
Propranolol	241, 613, 636, 640
Propylparaben	189, 604
Protriptyline	396, 564, 632
Pseudoephedrine	275, 617, 628
Pymetrozine	594
Pyrene	589-591
Pyridine	249, 389
Pyridoxine	605-607
Pyridoxyl phosphate (B6)	612
Pyrilamine	334
Pyroglutamate	174
Pyruvate	174

Q

Quercetin	266
-----------	-----

R

Raffinose	340, 600, 604
Ranitidine	325, 583
Reserpine	290, 299
Resorcinol	230-231, 584
Retinol	605
Retinol acetate	605
Retinol palmitate	605
Rhamnose	600
Riboflavin (B2)	605-606, 612
Ribonuclease (RNase)	371, 381, 571-572
Ribonuclease A	370, 400, 403, 425-426, 558-560, 571-572, 576
Ribose	600
Ropivacaine	619

S

Saccharin	232
Saccharose	600
Salicylic acid	235, 290, 618, 621
Sarcosine	396, 564
Scopolamine	275, 617
Sebutylazine	588
Selenate	175
Serine	396, 564
Simazine	587-588
Sodium	173

Sorbic acid	290, 609, 621
Sorbitol	338-341, 604
Stachyose	340, 604
Succinate	174
Succinic Acid	338, 581, 610
Sucrose	338, 340-341, 600, 604
Sulfachloropyridazine	616
Sulfadiazine	249, 616, 640-641
Sulfadimethoxine	616, 641
Sulfamerazine	249, 616, 640-641
Sulfamethazine	616, 640-641
Sulfamethizole	616, 641
Sulfamethoxazole	249, 258, 616, 632
Sulfamethoxypyridazine	616
Sulfanilamide	249, 640-641
Sulfanilic acid	641
Sulfapyridine	616
Sulfate	172, 174-175
Sulfathiazole	249, 616, 640
Sulfisoxazole	641
Sulfmethazine	249
Sulfur	507
Sulindac	235, 615

T

Talbarbitone	643
Tartarate	175
Tartaric acid	338, 610
Tebuthiuron	587
Terbutylazine	588
Terphenyl-d14	591
o-Terphenyl	299
Testosterone	312, 613
Tetracaine	266, 619
2,3,4,6-Tetrachlorophenol	592
Tetryl	586-587
Theobromine	276, 316, 608
Theophylline	276, 316, 334
Thiabendazole	594
Thiamine pyrophosphate (TPP)	594
Thiamine (B1)	607, 612
Thiocyanate	172
Thiosulfate	172
Thiourea	189, 590
Threonine (THR)	396, 564
Thymidine	564
Thymine	564, 579
Thymol	602
Thyroglobulin	421, 425-426, 432
α-Tocopherol	605, 608
β-Tocopherol	608
γ-Tocopherol	605, 608
Tocopherol acetate	605
o-Tolidine	598
Tolmetin	235, 615, 618, 622
m-Tolualdehyde – DNPH	584
Toluene	240, 299, 589-591
m-Toluic	623
m-Toluidine	585

o-Toluidine	183, 583, 598
Trehalulose	600
Triamcinolone	290, 635
Triamcinolone acetonide	290
2,4,6-Trichlorophenol	592
2,4,5-Trichlorophenoxyacetic acid	591
3,4,5-Trichlorophenoxyacetic acid	593
2,4,5-Trichlorophenoxypropionic acid (Silvex)	591
3,4,5-Trichlorophenoxypropionic acid	593
Triflupromazine	632
Trimipramine	258, 289, 632, 635-636
1,3,5-Trinitrobenzene	586-587
2,4,6-Trinitrotoluene	586-587
Tripelennamine	632
Triphenylene	299
Triprolidine	275, 279, 617
Tryptophan (TRP)	396, 459, 564
Tyrosine (TYR)	396, 459, 564
Tébuthiuron	591

U

Uracil	240-241, 258, 299, 425-426, 564, 579
Uridine	564, 588

V

Valeraldehyde – DNPH	584
Valine	396, 564
Valine3-Glycine4 (Nα-acetylated)	308, 390
Valine3-Valine4 (Nα-acetylated)	308, 390
Valine-tyrosine-valine	572
Vitamin A	605
Vitamin B12	421
Vitamin C	605-606, 609
Vitamin D2	605
Vitamin D3	605
Vitamin E (a-VE)	605

W

Warfarin	640
----------	-----

X

Xanthosine-5'-monophosphate (XMP)	568
Xylitol	339
Xylose	600

Ordering Information

Easy Ordering Terms and Conditions

Discounts and Delivery

Agilent Technologies specializes in fast delivery. In the US, if you call before 2 PM EST, we will ship your order that day. You may also request overnight express delivery before 6 PM EST and you will have your order the next day. Volume discounts on a variety of individual products are offered when the entire quantity is shipped to one address at one time.

A shipping and handling fee will be added to your order unless the purchase is over \$2000 US for orders place online or over \$4000 for orders place via phone. Special shipping (i.e., overnight in the US) is available in most regions at an additional cost.

Agilent is required to collect all state and local sales taxes unless the buyer's tax-exempt certificate is on file with Agilent Technologies. Please be prepared to provide a copy if it is not on file, when placing your order.

Please check with your Agilent Customer Service Representative, local Authorized Distributor, or the Agilent website for current prices, special offers, promotions and discounts when placing your order.

Satisfaction Guaranteed

If you are not satisfied with your Agilent product within the first 60 days, you may return your purchase in its original condition for a full refund or credit. A return policy statement is included in every Agilent shipment and posted under Product Information on the website. In the US and Canada, please call for a Return Authorization form and return instructions at **1-800-227-9770**. If your Agilent product was purchased from a distributor, please contact the distributor.

Shipping Damages

If items are damages in transit, please follow the instructions below:

- If a shipment is visibly damaged upon arrival, do not accept it until the person making the delivery has endorsed the bill of lading with statement for the extent of the damage.
- If any damage is found after unpacking, retain all cartons and inner packaging and immediately request an inspection from the carrier.
- Notify the Agilent Customer Contact Center at **1-800-227-9770** about the damaged shipment so that we can make the appropriate sales adjustment and/or provide you with return instructions (Sales order number, product number and quantity damaged will be needed).

Easy Ways To Order

- Phone: **1-800-227-9770** (option 1, 1) in the US and Canada — Mon-Fri, 8AM to 8PM EST
- Fax: **1-302-633-8901** in the US
- Email: **cag_sales-na@agilent.com** in the US and Canada
- Online: **www.agilent.com/chem** in the US and Canada

Payment Options

- In the US, Visa, MasterCard, Discover and American Express are accepted with a minimum order of \$20 (not applicable in all countries).
- Email ePay@agilent.com to make an electronic payment using the ACH/EFT (Automated Clearing House/Electronic Funds Transfer) method.
- Establish a charge account through your Agilent Customer Service Representative or Your Local Agilent sales office. An account number will be assigned to you for charging your purchases. Payment terms are net 30 days from the invoice date. All orders are subject to credit approval.

We will be happy to supply a price quote via, phone, email or fax if you need it in writing.

Warranties

All Agilent Technologies products in this catalog are designed and manufactured to stringent standards under the Agilent quality system registered to ISO 9001. At Agilent, we back every product with a 90-day warranty and a money-back guarantee. If Agilent receives notice of defects during the warranty period. Agilent shall, at its option, either repair or replace products which prove to be defective. If Agilent is unable, within a reasonable time, to repair or replace any product to a condition as warranted, the buyer shall be entitled to a refund of the purchase price upon return of the product to Agilent. The warranty period for each product begins on the day of shipment.

This warranty shall not apply to any defect, failure, or damage caused by improper use or improper or inadequate maintenance or care. This warranty is exclusive and no other warranty, whether written or oral, is expressed or implied. Agilent specifically disclaims the implied warranties of merchantability and fitness for particular purposes. The remedies provided herein are the buyer's sole and exclusive remedies. In no event shall Agilent be liable for direct, indirect, special, incidental, or consequential damages (including loss of profits) whether based on contract, tort, or any other legal theory.

Agilent Technologies Order Form

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.



Agilent Technologies

Order Date	Purchase Order or Credit Card Number & Expiration	Taxable Y or N?	Yes	No
Name:		If No, please provide Certificate #		
Title:				
Phone:		Fax:		
Company:		Email:		
	Shipping Address	Billing Address		
Company:				
Street:				
Room/Bldg/Dept:				
City:				
State/Province/Country:				
Zip/Postal Code:				
Deliver to:				
Part Number	Description	Quantity	Price	Total Cost
Special Instructions:			Subtotal:	
			Tax:	
			Total:	
For Assistance: Please call the Agilent Technologies Customer Contact Center at 1-800-227-9770 (U.S. and Canada) Email Node: cag_sales-na@agilent.com Fax Number: 302-633-8901 U.S. and Canada Note: 1. All pricing, tax, discount, and availability information is subject to verification by Agilent Technologies. 2. Shipping and handling is free for orders over \$2000 US.				

Agilent Technologies Order Form

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.



Agilent Technologies

[illegible]

Agilent Technologies Order Form

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.



Agilent Technologies

Order Date		Purchase Order or Credit Card Number & Expiration		Taxable Y or N?	Yes	No
Name:				If No, please provide Certificate #		
Title:						
Phone:				Fax:		
Company:				Email:		
		Shipping Address		Billing Address		
Company:						
Street:						
Room/Bldg/Dept:						
City:						
State/Province/Country:						
Zip/Postal Code:						
Deliver to:						
Part Number	Description	Quantity	Price	Total Cost		
Special Instructions:				Subtotal:		
				Tax:		
				Total:		
For Assistance: Please call the Agilent Technologies Customer Contact Center at 1-800-227-9770 (U.S. and Canada) Email Node: cag_sales-na@agilent.com Fax Number: 302-633-8901 U.S. and Canada Note: 1. All pricing, tax, discount, and availability information is subject to verification by Agilent Technologies. 2. Shipping and handling is free for orders over \$2000 US.						



For more information

Buy online:

www.agilent.com/chem/store

Contact us:

www.agilent.com/chem/contactus

This information is subject to change without notice.

© Agilent Technologies, Inc. 2012

Printed in Canada October 31, 2012

5991-1059EN



Scan the QR code
with your smartphone
for more information



The Measure of Confidence



Agilent Technologies