



GC AND GC/MS

Your Essential Resource for Columns & Supplies


20
15/16



GC AND GC/MS

Achieve excellent, reproducible performance for difficult samples

For over 40 years, Agilent has broken new ground with innovations in Gas Chromatography. We continue our leadership tradition by offering the industry's broadest selection of GC and GC/MS columns and supplies. All are manufactured to Agilent's exact specifications to minimize downtime and ensure consistent, high-quality results that you can rely on.

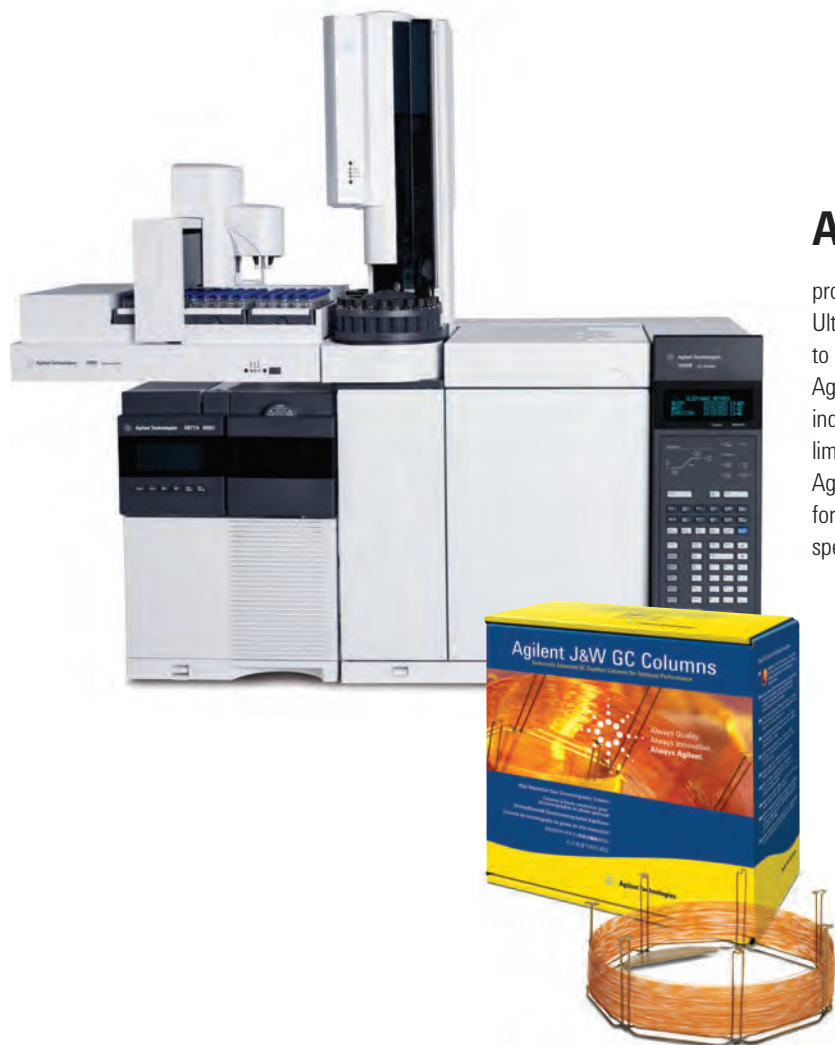

Agilent Ultra Inert solutions

provide the flow path inertness vital to analytical success. Ultra Inert split and splitless liners are manufactured and tested to our highest level of scrutiny to ensure quality and consistency. Agilent J&W Ultra Inert GC columns are tested with the industry's most demanding test probe to reduce detection limits and produce more accurate data for difficult analytes. Agilent GC and GC/MS instruments bring together all elements for trace-level analysis, dramatically improving MS resolution, spectral integrity, and detection limits.




GC and GC/MS supplies

More samples, lower detection levels, with fewer analysts. These demands challenge laboratories to maximize the productivity and performance of their instrumentation. To help you stay ahead, Agilent is continuously improving our extensive portfolio of innovative, award winning GC columns and supplies, designed to help you resolve many of the day to day setbacks encountered in your lab. You can avoid downtime and your time can be better spent on meeting your analytical and business challenges.

For labs pushing the detection limits of trace level analysis on very active compounds, **Agilent Inert Flow Path solutions** ensure a reliably inert flow path for higher sensitivity, accuracy, and reproducibility. Install industry leading Agilent J&W GC columns with new proprietary design GC column nuts and ferrules to simplify your day yet maximize your GC and GC/MS systems output.

- Inert Flow Path components – Ultra Inert GC columns, Ultra Inert liners, Ultra Inert gold seals, UltiMetal Plus Capillary Flow Technology devices with Flexible Metal ferrules – have Agilent proprietary deactivation chemistries to ensure sample integrity.
- "Better Connectivity" with products such as Self Tightening column nuts, UltiMetal Plus Flexible Metal ferrules, and Ultra Inert liners in Touchless packaging improves productivity with ease of use and convenience.
- Full portfolio of premium GC products to support your lab needs – including Agilent CrossLab brand and Agilent Bulk supplies packaging.

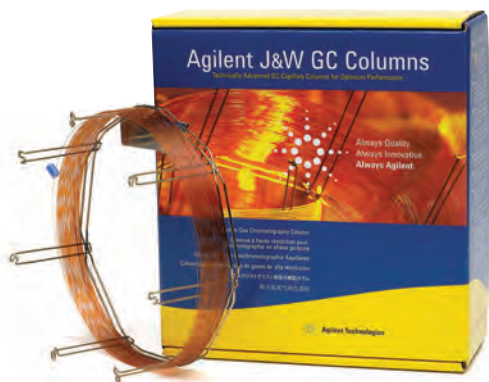

Agilent J&W GC columns

deliver the best inertness for acids, bases, and mixed functional compounds, the lowest bleed levels, and the tightest column-to-column reproducibility. Mass Spec Grade GC columns (VF-ms, DB-ms and HP-ms) give you robust performance, low column bleed, and a wide range of selectivity. LTM column modules combine a fused silica capillary GC column with heating and temperature-sensing components for efficient column heating and cooling. What's more, integrated guard columns protect your analytical columns from non-volatile compounds in the sample matrix.

Table of Contents

Agilent Solutions, Services and Support	4	Detector Systems.....	92
Featured Products	8	Flame Ionization Detector (FID).....	92
Agilent Parts and Supplies	14	Electron Capture Detector (ECD)	100
GC and GC/MS Maintenance Schedule	14	Thermal Conductivity Detector (TCD)	103
Bulk GC Supplies.....	16	Flame Photometric Detector (FPD).....	108
Inlet Septa	17	Nitrogen Phosphorus Detector (NPD)	116
Inlet Liners.....	24	Nitrogen and Sulfur Chemiluminescence Detectors	122
Agilent Ultra Inert Liners	26	GC Standards	124
Capillary Column Ferrules and Nuts	34	7820A GC System.....	125
GC Column Connection Supplies	40	GC/MS Parts and Supplies	141
Capillary Flow Technology Supplies.....	42	MSD Contamination	142
Press-fit Capillary Column Connectors	44	Ion Source	149
Graphpak Capillary Connectors	45	MSD Filaments.....	162
Large Valve Oven.....	46	Vent Valve Supplies	163
Valves and Loops	47	Gas Clean Filters	164
Sample Introduction Systems.....	50	Quadrupole Mass Filter.....	164
Agilent Vials and Closures for GC, GC/MS and GC/HS	56	MSD Electron Multipliers and Replacement Horn	165
Headspace Vials and Closures	56	Vacuum Systems and Pumps	166
High Performance Septa.....	57	Diffusion Pump.....	168
CombiPAL Headspace Vials and Closures	59	Quiet Cover	169
Crimping and Decapping Tools	60	Foreline Pump.....	170
Teledyne Tekmar Purge and Trap Supplies	62	7000 Triple Quadrupole GC/MS	171
Markes Thermal Desorption	65	7200 Q-TOF for GC/MS.....	175
Inlet Systems	67	240-MS Ion Trap Parts and Supplies	178
Split/Splitless Inlets.....	69	220-MS Parts and Supplies.....	180
Multimode Inlet.....	76	GC/MS Standards.....	181
Cool On-Column Inlets.....	78	Agilent Syringes	183
Programmable Temperature Vaporizer (PTV) Inlets ...	81		
Purged Packed Inlets	88		

Agilent CrossLab GC Parts and Supplies	192	Applications	500
Product Introductions.....	193	Environmental	501
Supplies for Bruker, Varian GC Systems.....	206	Hydrocarbons.....	501
Supplies for PerkinElmer GC Systems.....	221	Pesticides and Herbicides.....	506
Supplies for Shimadzu GC Systems.....	227	Semivolatiles.....	532
Supplies for Thermo Scientific GC Systems.....	234	Volatiles	544
Supplies for CTC GC Autosamplers	239	Air Analysis.....	549
Agilent J&W GC Columns	240	Food, Flavor and Fragrance	554
Column Selection	247	Energy and Fuels.....	576
Column Selection Principles	247	Industrial Chemical	602
GC Column Application and Method Guides	262	Forensic Toxicology and Pharma.....	635
Agilent J&W Ultra Inert GC Columns.....	286	Indices	652
Agilent J&W High Efficiency GC Capillary Columns	293	Ordering Information	684
Low-bleed GC/MS Columns	294		
Premium Polysiloxane Columns.....	318		
Polyethylene Glycol (PEG) Columns	351		
Specialty Columns	362		
PLOT Columns	420		
Columns with Non-Bonded Stationary Phases.....	441		
Guard Columns	445		
LTM Column Modules	447		
Fused Silica Tubing	464		
Stainless Steel Tubing.....	469		
Packed GC Columns	470		
Custom GC Column Ordering.....	481		
GC Column Test Standards.....	482		
Column Installation and Troubleshooting.....	484		



PUT MORE THAN 40 YEARS OF RELENTLESS INNOVATION BEHIND YOUR EVERY RESULT

By continually raising the standards for technologies that support your routine analyses, Agilent's R&D efforts have led to breakthroughs such as:

- **New GC columns** that help you achieve higher levels of inertness and column-to-column reproducibility
- **LC column choices** that deliver the sensitivity and reliability you need for demanding applications
- **Cutting-edge sample preparation products** that promote reliable extraction and concentration
- **Fresh atomic and molecular spectroscopy ideas** for identifying and confirming targets and unknowns

Longtime Agilent customers have experienced our commitment firsthand. And now, we look forward to demonstrating how Agilent's approach to relentless innovation can work to your advantage, too.



CHEMICAL ANALYSIS SOLUTIONS

Food

From high-volume pesticide screening in food products to rapid identification of pathogens, Agilent understands the analytical needs of food producers, shippers, and regulators. Utilizing our easy-to-use analyzers and updated screening libraries, customers can quickly develop robust and reliable methods. Agilent's leading gas chromatography and mass spectrometry systems are widely regarded as valuable food testing techniques for an array of different analyses.

Environmental

Agilent offers more than 40 years of environmental testing and regulatory expertise. We help government and private labs with the full range of assays, from routine testing of soils for heavy metals to detection of pharmaceuticals in groundwater, in concentrations down to parts per trillion.

Energy & Chemicals

Agilent collaborates closely with process industry customers to offer analytical systems that meet their needs for separation, detection, throughput, and support. We'll even preconfigure custom or standard analyzers so they arrive at the lab ready-to-go. From crude oil, natural gas, and refining, to specialty chemicals and alternative fuels, Agilent provides the latest technologies and solutions to increase quality, safety, and profitability for energy and chemical labs, while meeting the industry's stringent quality requirements. Agilent leads the way in ASTM collaborations that have evolved – and will continue to evolve – into industry standards.

Forensics


Whether testing for poisons in a forensics investigation, screening athletes for performance enhancing drugs, analyzing samples for recreational drugs, or checking a crime scene for explosive residue – lives and professions may be dependent on the accuracy of your equipment. Agilent Technologies leads the industry with a comprehensive portfolio of workflow solutions that provide the ability to identify, confirm and quantify thousands of substances.

Lab Informatics

The ways labs capture, analyze and share data profoundly affect their efficiency. Agilent offers a rich, integrated suite of software products built on customer-driven architectural values with the Agilent OpenLAB Software Suite. OpenLAB delivers superior performance and connection across multiple systems, providing open systems integration and investment protection. Our commitment is to deliver more value across each step in the life cycle of scientific data – from data collection and analysis to interpretation and management.

Materials Science

Agilent offers a newly expanded portfolio of instruments used for the research, manufacturing and testing of advanced materials, from precision optics to pulp and paper. Tools for atomic spectroscopy, molecular spectroscopy, chromatography, and X-ray crystallography all support continuous progress in materials science.

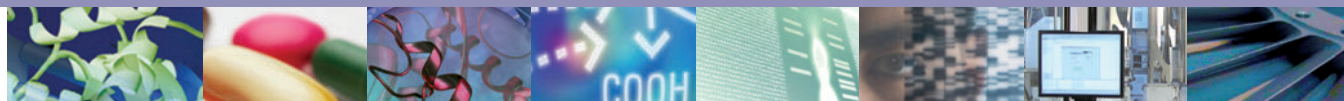

LIFE SCIENCE SOLUTIONS

Biopharmaceutical

Biotherapeutics have enormous potential to improve human health, with growing numbers of protein and antibody therapeutics to address unmet medical needs. At every development stage, from disease research to QA/QC and manufacturing, Agilent can help you make the right choices for moving therapeutics to market. We understand the biopharmaceutical workflow so our product families work together seamlessly, as engines of research, discovery, and development. Agilent columns deliver complete characterization of biomolecules using reversed-phase, size exclusion, ion exchange, and affinity chromatography. Our bio-inert supplies ensure that every part of your workflow delivers the performance you need to optimize your bio-separation.

Pharmaceutical

You need the most efficient processes to evaluate drug candidates, determine efficacy, and ensure safety and compliance during development and manufacture. Agilent has worked with pharma companies for many years to ensure reliability and reproducibility for regulatory compliance, from lab-to-lab and around the world. Our pharma solutions provide high-throughput capability at every stage of the product lifecycle, with automated sample prep, industry-leading U/HPLC systems, the largest family of Fast LC columns, open access LC/MS, spectroscopy, and automated dissolution. A complete family of LC supplies and lamps help optimize every analysis and take day-to-day lab efficiency one step further.

Proteomics

Research into how large sets of proteins affect the health of an organism requires special sets of analytical tools. Agilent has built a formidable arsenal of liquid chromatograph/mass spectrometers, bioinformatics systems, multiple affinity protein removal columns, and OFFGEL electrophoresis for protein identification and protein biomarker discovery. Accurate-Mass mass spectrometry and the microfluidic HPLC-Chip/MS are two Agilent innovations speeding the work of proteomics researchers around the globe.

Metabolomics

Collections of small molecules are increasingly being seen as rich sources of biomarkers, but studying metabolites presents many challenges. The need for speed, accuracy, and powerful interpretation capabilities in looking at chemical profile snapshots is underscored because molecules are constantly entering, leaving or changing within the metabolome. Agilent's GC, LC, and MS portfolios, along with our excellent bioinformatics offerings, user-customizable METLIN metabolite database for LC/MS, and the industry's first commercial GC/MS retention time locked metabolite library align well with the needs of metabolomics researchers.

Genomics

Agilent is a global leader in microarrays, scanners, and NGS reagents used in a wide variety of genomic-based disease research experiments. Our SureSelect and HaloPlex Target Enrichment Systems dominate the category, streamlining next generation sequencing studies. Agilent offers a wide range of catalog CGH and gene expression microarrays and a highly-developed capability to produce custom arrays using our free online design tool, SureDesign. All Agilent microarrays feature highly sensitive, selective 60-mer probes, and, with as many as eight arrays printed on a slide, the cost per sample is cost-efficient.

Life Science Informatics

Mirroring its extensive instrument portfolio, Agilent offers the industry's most extensive suite of bioinformatics software, helping users derive knowledge from complex genomic, proteomic, metabolomic and other biological data. SureCall and CytoGenomics software analyzes NGS and aCGH data and the GeneSpring suite provides multi-omic analysis and visualization capabilities to help compare complex datasets to explore biological questions from multiple perspectives. The GeneSpring suite includes the GX module for microarray-based gene expression and genotyping data, the PA module for Pathway Analysis and multi-omic analysis and the MPP software, which analyzes mass spec data from proteomics and metabolomics experiments.

Lab Automation

To meet the skyrocketing demand for more throughput and automation, Agilent has substantially expanded its lab automation offerings. The Agilent line of liquid handlers and microplate processors are designed to streamline high-volume life science workflows. Agilent is also continually upgrading its advanced autosamplers for LC, GC, LC/MS and GC/MS, adding functionality and speed to reflect the performance of its advanced instruments.

Vacuum Technology

Agilent works with customers to solve vacuum challenges from experiments in high-energy physics to developing systems for nanotechnology. Agilent manufactures vacuum systems used in its own mass spectrometry instruments as well as those of other manufacturers. Agilent's vacuum technology has been proven by the most powerful physics experiment ever built, CERN's Big Bang machine, which was used in the discovery of the Higgs boson.



Get the Agilent Service Guarantee

Should your instrument require service while covered by an Agilent Advantage service agreement, we guarantee repair or we will replace your instrument for free.

No other company offers this level of commitment to keep your lab up and running at peak efficiency.

Agilent Service and Support for Instrument Systems

Focus on what you do best

For over 40 years, Agilent has been building and maintaining the instruments you count on to stay competitive and successful. Trust us to protect your investment with a broad portfolio of services, backed by a global network of experienced service professionals dedicated to the productivity of your lab.

Agilent Advantage Service Plans

The best service available for your Agilent instruments

Agilent offers a flexible range of service plans so that you can choose the level of coverage that is best for your lab.

- **Agilent Advantage Gold** – Priority-one coverage for ultimate uptime and productivity
- **Agilent Advantage Silver** – Comprehensive coverage for dependable laboratory operations
- **Agilent Advantage Bronze** – Total repair coverage at a fixed annual price
- **Agilent Repair Service** – Basic coverage for reliable instrument repair

Agilent Advantage service plans include Agilent Remote Advisor for real-time remote monitoring and diagnostics. Through secure internet connections, you can interact with Agilent service professionals, receive detailed asset reports, and configure text or email alerts to notify you before problems occur – helping you to maximize instrument uptime and optimize laboratory workflows.

And for Agilent-quality service on analytical instruments from other leading manufacturers, Agilent CrossLab services offer the same quality coverage you have come to expect from the expert Agilent engineers you know and trust.



Agilent Compliance Services

Equipment qualification that meets the most stringent requirements

Enterprise Edition Compliance was developed to streamline qualification delivery compliance across your entire lab. Used worldwide in regulated labs, including standards organizations and regulatory agencies, Enterprise Edition enables you to:

- Improve qualification efficiency by harmonizing protocols across platforms to ensure greater efficiency and minimize regulatory risk
- Standardize your entire compliance operation with robust test designs that work with all your instruments
- Add, remove or reconfigure tests based upon your unique user requirements
- Reduce staff review time significantly with consistently formatted, computer generated, tamper-proof reports

Agilent Education and Consulting Services

Our best minds, working for you

Make the most of your instrument with training and consulting from the same experts who designed the instruments, software and processes you use every day.

- Classroom, online, and on-site training in instrument operation, troubleshooting and maintenance
- Customized consulting services to meet your lab's unique needs

The Agilent Value Promise – 10 Years of Guaranteed Value

In addition to continually evolving products, we offer something else unique to the industry – our 10-year value promise guarantee. The Agilent Value Promise guarantees you at least 10 years of instrument use from your date of purchase, or we will credit you with the residual value of the system toward an upgraded model. Not only does Agilent ensure a reliable purchase now, but we also ensure that your investment is just as valuable in the future.


For more detailed information, please go to www.agilent.com/chem/services or contact your local Agilent Services and Support representative.



Technical Support at work for you

Have a hardware, software, application, instrument repair or troubleshooting question? Agilent's technical experts are available to answer your questions. With years of laboratory experience, our technical support specialists can provide in-depth knowledge and experience.

For questions pertaining to supplies found in this catalog, contact your local Agilent office or Authorized Agilent Distributor or visit www.agilent.com/chem/techsupport

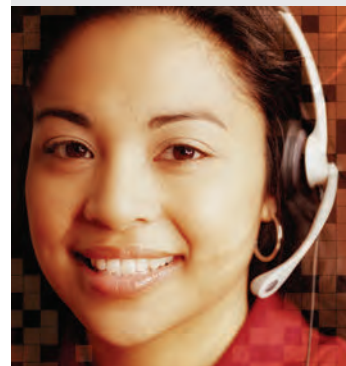

Need more information?

Visit www.agilent.com/chem/contactus to:

- Locate your nearest Agilent office or distributor for expert technical support.
- Get fast sales and product assistance by phone. Simply use the scroll-down menu to select your country.
- Receive email assistance using our convenient online forms.



GC and GC/MS Applications

Industry-specific applications from your partner in chromatography

With over 40 years of chromatography expertise, Agilent is a great resource for all types of applications. In fact, we're developing new ones every day.

Simply turn to the pages listed below for the most current applications based on your area of specialization.

Environmental — you'll learn how to perform critical analyses — such as measuring the levels of atmospheric halocarbons and identifying organochlorine pesticides in soil — while meeting your increasing demands for speed and accuracy. **Turn to page 501.**

Food, Flavors, and Fragrances — we'll discuss how to ensure quality, safety, and regulatory compliance for fragrances, perfumes, and essential oils. Applications focus on chiral compounds, menthol, and FAMES. **Turn to page 554.**


Energy and Fuels — here you'll find applications — such as the analysis of sulfur compounds in propylene — that you can use right away to meet regulatory requirements, improve efficiency, and maintain good environmental stewardship. **Turn to page 576.**

Industrial Chemical — we'll help you maintain product quality — and production efficiency — by sharing the latest applications for alcohols, halogenated hydrocarbons, aromatic solvents, phenols, and inorganic gases. **Turn to page 602.**

Forensic Toxicology and Pharma — we'll bring you fully up-to-date on the newest screening methods for controlled substances such as amphetamines, narcotics, and alcohol. We'll also review the latest techniques for monitoring residual solvents. **Turn to page 635.**

TIPS & TOOLS

Search the ScanView database to find almost 2000 GC applications and standard methods of all types, old and new. Get your free copy of ScanView at www.agilent.com/chem/scanview

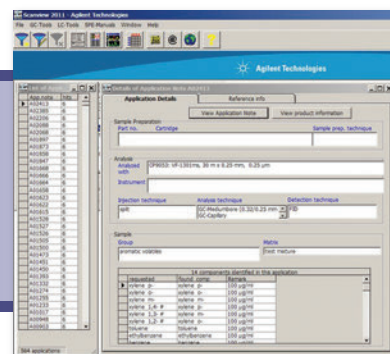

Environmental Applications, Hydrocarbons

Unleaded Gasoline

Column: DB-VRX
124-1534
30 m x 0.45 mm, 2.55 μ m

Carrier: Helium at 109 cm/s (10.4 mL/min), measured at 40 °C

Oven: 40 °C for 2 min,
40-200 °C at 12 °C/min,
200 °C for 5 min

Sampler: Purge and Trap (O.I.A. 4560)

Trap: BTEX (Supelco) at 50 °C during purge

Desorb: 270 °C for 1 min


Injection: LVI (Low Volume Injector)

Detector: A: FID, 250 °C


B: PID (O.I.A. 4430), 200 °C

Sample: 115 ppb gasoline in 5 mL water

A: FID

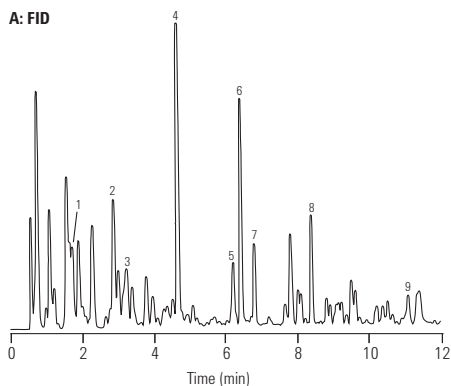

B: PID

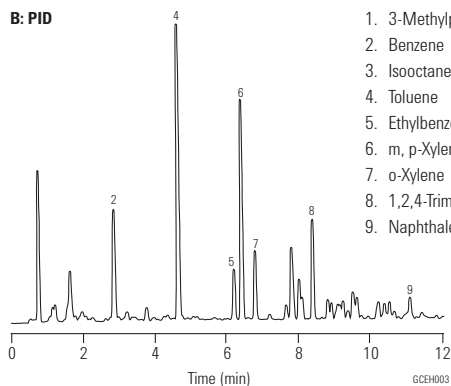

1. 3-Methylpentane
2. Benzene
3. Isooctane
4. Toluene
5. Ethylbenzene
6. m, p-Xylene
7. o-Xylene
8. 1,2,4-Trimethylbenzene
9. Naphthalene

GCEH003

Determination of Chlorophenols in Water and Soil

Column: VF-5ms
CP8961
60 m x 0.32 mm, 0.25 µm

Oven: 60 °C, 30 °C/min to 300 °C

Carrier: He 80 kPa, 0.8 bar, 5.7 psi

Injection: Splitless, initial time: 1 min; Splitflow: 50 mL/min
250 °C
2 µL


Detector: MS
280 °C

Sample: Isohexane

Sample Conc: Standard, 1 µg/mL, derivatization with acetic acid anhydride

Dr. Weßling, Laboratorien GmbH

- | | |
|---------------------------|-------------------------------|
| 1. Phenol | 10. 2,4,6-Trichlorophenol |
| 2. 2-Chlorophenol | 11. 2,3,6-Trichlorophenol |
| 3. 3-Chlorophenol | 12. 2,3,5-Trichlorophenol |
| 4. 4-Chlorophenol | 13. 2,4,5-Trichlorophenol |
| 5. 2,6-Dichlorophenol | 14. 2,3,4-Trichlorophenol |
| 6. 2,4+2,5-Dichlorophenol | 15. 3,4,5-Trichlorophenol |
| 7. 3,5-Dichlorophenol | 16. 2,3,5,6-Tetrachlorophenol |
| 8. 2,3-Dichlorophenol | 17. 2,3,4,6-Tetrachlorophenol |
| 9. 3,4-Dichlorophenol | 18. 2,3,4,5-Tetrachlorophenol |

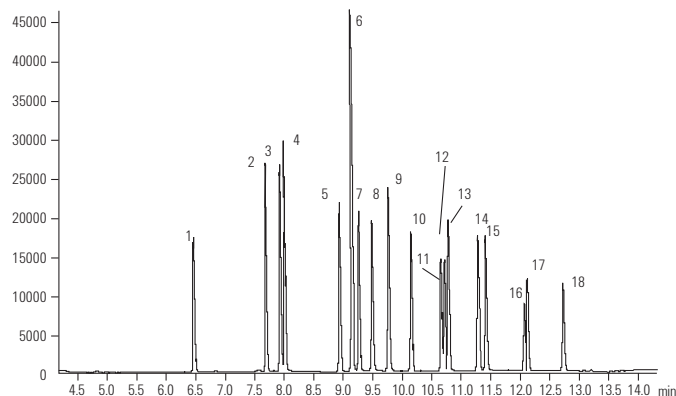**PBDEs by ECD**

Column: DB-XLB
15 m x 0.18 mm, 0.07 µm
Agilent Technologies custom column

Carrier: Hydrogen at 72 cm/s at 100 °C (4.0 mL/min),
constant flow mode

Oven: 100 °C for 0.5 min
100 °C to 300 °C at 30 °C/min
300 °C for 5 min


Injection: Split, 250 °C
Split ratio 20:1

Detector: ECD, 300 °C
Peak, Congener (2.5 mg/mL)

Sample: 1 µL

*Special thanks to AccuStandard, Inc. of New Haven, CT,
for PBDE standards.*

- | | |
|-----------------------------------|--|
| 1. 2,2',4-TriBDE (BDE-17) | 8. 2,2',3,4,4'-PentaBDE (BDE-85) |
| 2. 2,4,4'-TriBDE (BDE-28) | 9. 2,2',4,4',5,6'-HexaBDE (BDE-154) |
| 3. 2,3',4',6-Tetra-BDE (BDE-71) | 10. 2,2',4,4',5,5'-HexaBDE (BDE-153) |
| 4. 2,2',4,4'-Tetra-BDE (BDE-47) | 11. 2,2',3,4,4',5'-HexaBDE (BDE-138) |
| 5. 2,3',4,4'-TetraBDE (BDE-66) | 12. 2,2',3,4,4',5',6'-HeptaBDE (BDE-183) |
| 6. 2,2',4,4',6-PentaBDE (BDE-100) | 13. 2,3,3',4,4',5,6'-HeptaBDE (BDE-190) |
| 7. 2,2',4,4',5-PentaBDE (BDE-99) | 14. DecaBDE (BDE-209) (12.5 mg/mL) |

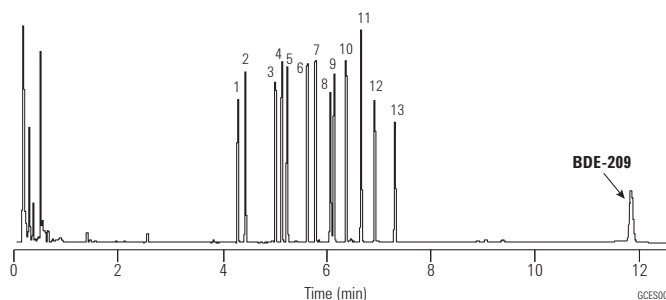

Diesel Fuel

Column: DB-5ms
125-5532
30 m x 0.53 mm, 1.50 μ m

Carrier: Helium at 48.5 cm/s, measured at 60 °C

Oven: 60 °C for 2 min
60-300 °C at 12 °C/min
300 °C for 10 min

Injection: Direct, 280 °C

Detector: FID, 250 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 μ L injection in hexane
A: Standard, 50 ng/component
B: Sample, 0.6 mg/mL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759


Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Syringe: 10 μ L tapered, FN 23-26s/42/HP, 5181-1267

**Diesel fuel standard
50 ng/component**

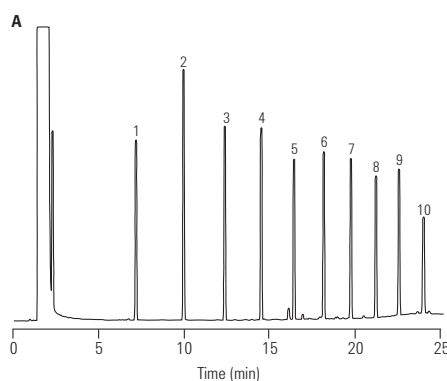

**Diesel fuel
0.6 mg/mL**

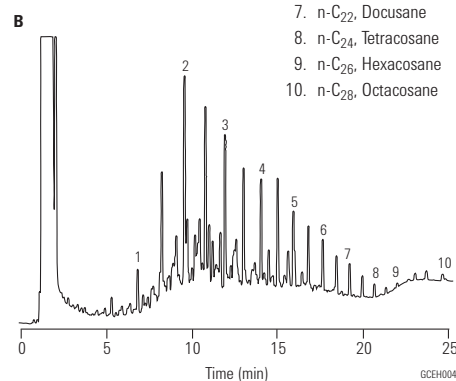

1. n-C₁₀, Decane
2. n-C₁₂, Dodecane
3. n-C₁₄, Tetradecane
4. n-C₁₆, Hexadecane
5. n-C₁₈, Octadecane
6. n-C₂₀, Eicosane
7. n-C₂₂, Docosane
8. n-C₂₄, Tetracosane
9. n-C₂₆, Hexacosane
10. n-C₂₈, Octacosane

Analysis of Polycyclic Aromatic Hydrocarbons


Column: VF-Xms
CP8805
30 m x 0.25 mm, 0.10 μ m

Sample: 1 μ L ca. 3 ng per component on-column

Carrier: Helium, 60 kPa

Injection: Split, T=275 °C

Detector: Agilent Ion Trap MS

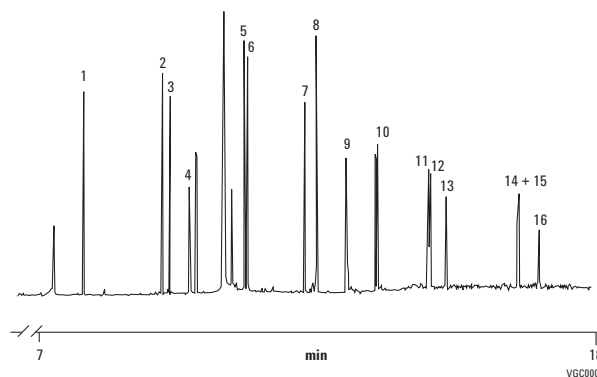

1. Naphthalene
2. Acenaphthylene
3. Acenaphthene
4. Fluorene
5. Phenanthrene
6. Anthracene
7. Fluoranthene
8. Pyrene
9. Chrysene
10. Benzo[a]anthracene
11. Benzo[k]fluoranthene
12. Benzo[b]fluoranthene
13. Benzo[a]pyrene
14. Indeno[1,2,3-cd]pyrene
15. Dibenzo[a,h]anthracene
16. Benzo[g,h,i]perylene

Dioxins and Dibenzofurans

Column: CP-Sil 88
CP6173
50 m x 0.25 mm, 0.20 µm

Sample: 1.0 µL Toluene


Sample Conc: 100-400 pg/µL

Carrier: Helium, 170 kPa (1.7 bar, 24 psi)

Oven: 100 °C to 180 °C to 230 °C, 3 °C/min

Injection: Splitless

Detector: MSD

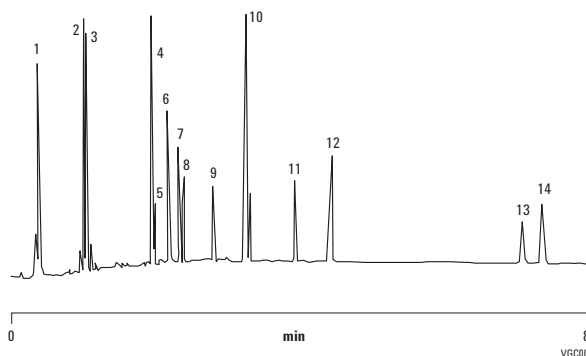

1. 2,3,7,8-TCDD
2. 2,3,7,8-TCDF
3. 1,2,3,7,8-PeCDF
4. 1,2,3,4,7,8-HxCDF
5. 1,2,3,6,7,8-HxCDF
6. 2,3,4,7,8-PeCDF
7. 1,2,3,4,7,8-HxCDD + 1,2,3,7,8-PeCDD
8. 1,2,3,6,7,8-HxCDD
9. 1,2,3,7,8,9-HxCDD
10. 1,2,3,4,6,7,8-HxCDF
11. 2,3,4,6,7,8-HpCDD
12. 1,2,3,4,6,7,8-HpCDD
13. 1,2,3,4,6,7,8,9-OCDF
14. 1,2,3,4,6,7,8,9-OCDD

**78 Semi-volatile Components
on an Agilent J&W DB-UI 8270D**

Column: DB-UI 8270D
122-9732
30 m x 0.25 mm, 0.25 µm

Instrument: Agilent 7890 Series GC

Carrier: Helium, 1.2 mL/min constant flow,
septum, purge 3 mL/min,
purge time on 0.7 min 50 mL/min, gas saver off


Oven: 30 °C (1.0 min), 15 °C/min to 100 °C,
20 °C/min to 240 °C (0.5 min),
15 °C to 325 °C (6.7 min)

Inlet: MMI in nonpulsed splitless mode, 1 µL at 275 °C

Inlet liner: Dual taper direct connect liner

Sampler: Agilent 7693, 10.0 µL syringe (p/n G4513-80216)

Detector: MSD: 325 °C Transfer line,
280 °C source,
150 °C quad, 35-500 amu range

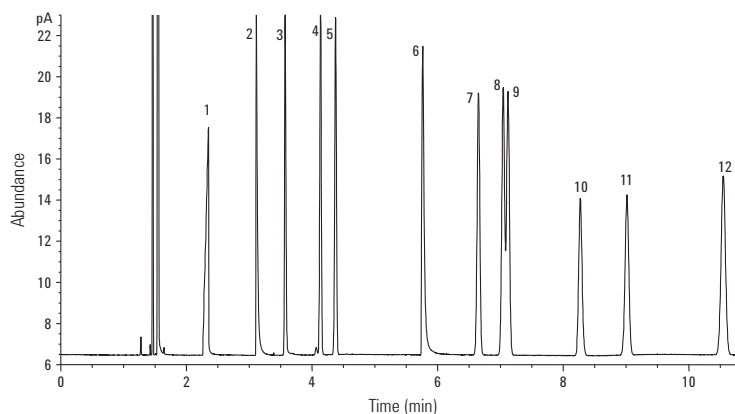

Example total ion chromatogram of a 78 component semi-volatile standard injection with a 10 ng on-column loading for each component.

Polybrominated Diphenyl Ethers (PBDEs)

Column: DB-5ms Ultra Inert
122-5512UI
15 m x 0.25 mm, 0.25 µm

Instrument: Agilent 6890N/5973B MSD


Sampler: Agilent 7683B, 5.0 µL syringe
(p/n 5188-5246),
1.0 µL splitless injection,
5 ng each component on-column

Carrier: Helium 72 cm/s, constant flow

Inlet: Pulsed splitless; 325 °C, 20 psi
until 1.5 min,
purge flow 50 mL/min at 2.0 min

Oven: 150 to 325 °C
(17 °C/min),
hold 5 min

Detector: MSD source at 300 °C,
quadrupole at 150 °C,
transfer line at 300 °C,
scan range 200-1000 amu

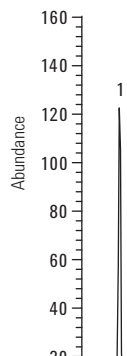

1. BDE-47
2. BDE-100
3. BDE-99
4. BDE-154
5. BDE-153
6. BDE-183
7. BDE-205
8. BDE-209

Suggested Supplies

Liner: Direct connect, dual taper, deactivated, 4 mm id, G1544-80700

Syringe: Autosampler syringe, 0.5 µL, 23 g, cone, 5188-5246

15+1 EU Priority PAHs

**Resolution of Critical Pairs
on an Agilent J&W DB-EUPAH Column**

Column: DB-EUPAH
121-9627
20 m x 0.18 mm, 0.14 µm

Instrument: Agilent 6890N/5975B MSD

Sampler: Agilent 7683B, 5.0 µL syringe, 0.5 µL splitless
injection, injection speed 75 µL/min

Carrier: Helium, ramped flow 1.0 mL/min (0.2 min),
5 mL/min² to 1.7 mL/min

Inlet: 325 °C splitless, purge flow 60 mL/min at 0.8 min

Oven: 45 °C (0.8 min) to 200 °C (45 °C/min),
2.5 °C/min to 225 °C, 3 °C/min to 266 °C,
5 °C/min to 300 °C, 10 °C/min to 320 °C (4.5 min)

Detector: MSD source at 300 °C, quadrupole at 180 °C,
transfer line at 330 °C, scan range 50-550 amu

All 15+1 EU regulated priority PAHs are well resolved with the DB-EUPAH column. Challenging benzo[b,k,j]fluoranthene isomers are baseline resolved, allowing for accurate quantitation of each isomer. In addition, baseline resolution is achieved for critical pairs benz[a]anthracene and cyclopenta[c,d]pyrene, cyclopenta[c,d]pyrene and chrysene, and indeno[1,2,3-cd]pyrene and dibenz[a,h]anthracene. This application demonstrates that the DB-EUPAH column can provide excellent sensitivity and selectivity for the analysis of EU regulated PAHs.


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, dual taper, deactivated, 4 mm id, G1544-80700

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

- | | |
|--------------------------|----------------------------|
| 1. Benzo[c]fluorene | 9. Benz[a]pyrene |
| 2. Benz[a]anthracene | 10. Indeno[1,2,3-cd]pyrene |
| 3. Cyclopenta[c,d]pyrene | 11. Dibenz[a,h]anthracene |
| 4. Chrysene | 12. Benzo[g,h,i]perylene |
| 5. 5-Methylchrysene | 13. Dibenz[a,i]pyrene |
| 6. Benzo[b]fluoranthene | 14. Dibenz[a,e]pyrene |
| 7. Benzo[k]fluoranthene | 15. Dibenz[a,i]pyrene |
| 8. Benzo[j]fluoranthene | 16. Dibenz[a,h]pyrene |

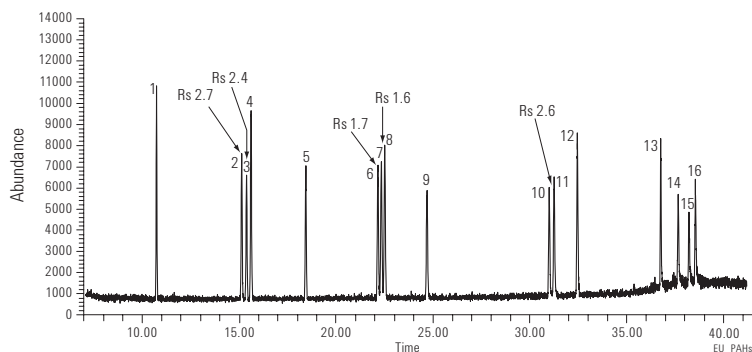

Environmental Applications, Pesticides and Herbicides

Fast CLP Pesticides

Column: DB-CLP1
123-8232
30 m x 0.32 mm, 0.25 µm

Column: DB-CLP2
123-8336
30 m x 0.32 mm, 0.50 µm

Instrument: Agilent 7890 GC with dual µECD

Carrier: Helium, constant flow 3.5 mL/min

Oven: 150 °C (hold 0.2 min), 45 °C/min to 250 °C,
18 °C/min to 300 °C, 30 °C/min to 330 °C, hold 2.5 min


Sampler: Agilent 7693

Injection: 1 µL splitless

Detector: µECD at 340 °C

Sample: 50 ng/mL CLP Pesticides

1. Tetrachloro-m-xylene
2. α-BHC
3. γ-BHC
4. β-BHC
5. Heptachlor
6. δ-BHC
7. Aldrin
8. Heptachlor epoxide
9. γ-Chlordane
10. α-Chlordane
11. Endosulfan I
12. 4,4'-DDE
13. Dieldrin
14. Endrin
15. 4,4'-DDD
16. Endosulfan II
17. 4,4'-DDT
18. Endrin aldehyde
19. Endosulfan sulfate
20. Methoxychlor
21. Endrin ketone
22. Decachlorobiphenyl

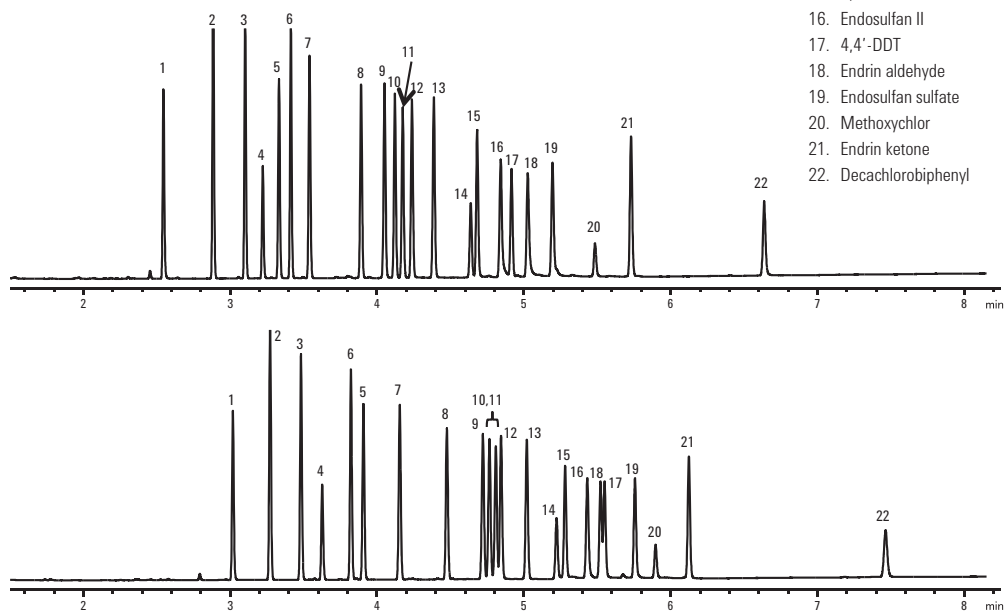

**EPA Method 504.1 – 1,2 dibromoethane (EDB),
1,2-dibromo-3-chloropropane (DBCP),
and 1,2,3-trichloropropane (123TCP)**

Column: DB-CLP1
123-8232
30 m x 0.32 mm, 0.25 µm

Column: DB-CLP2
123-8336
30 m x 0.32 mm, 0.50 µm

Carrier: Helium, constant flow, 3.75 mL/min


Oven: 50 °C, hold 1.5 min, 20 °C/min to 95 °C,
40 °C/min to 175 °C, hold 1.25 min

Injection: 2 µL, splitless, 200 °C


Detector: µECD, 300 °C

Sample: 100 ng/mL EPA 504.1 analytes, 100 ng/mL
chlorinated solvents + trihalomethanes

- | | |
|--------------------------|--|
| 1. Chloroform | 7. 1,1,2-Trichloroethane |
| 2. 1,1,1-Trichloroethane | 8. Dibromochloromethane |
| 3. Carbon tetrachloride | 9. 1,2-Dibromoethane (EDB) |
| 4. Trichloroethane | 10. Bromoform |
| 5. Bromodichloromethane | 11. 1,2,3-Trichloropropane (123TCP) |
| 6. Tetrachloroethane | 12. 1,2-Dibromo-3-chloropropane (DBCP) |

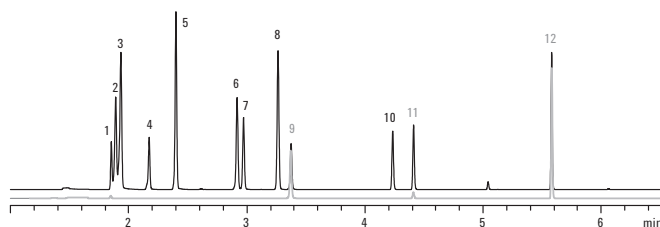

100 ng/mL chlorinated solvents + THMs
100 ng/mL EPA 504.1 analytes

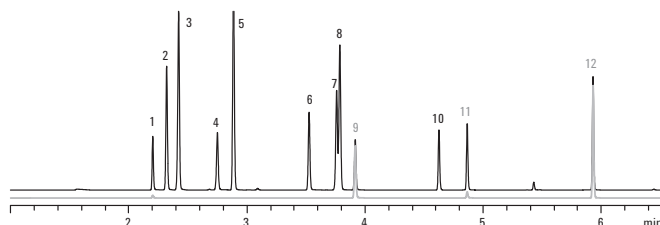

100 ng/mL chlorinated solvents + THMs
100 ng/mL EPA 504.1 analytes

Agilent J&W DB-CLP1/DB-CLP2 columns analyze 1,2-dibromoethane (EDB), 1,2-dibromo-3-chloropropane (DBCP), and 1,2,3-trichloropropane (123TCP) according to EPA Method 504.1 with cooler analysis temperatures allowing a faster GC cycle time.

Organochlorine Pesticides, EPA Method 8081B

Column: DB-CLP1
123-8232
30 m x 0.32 mm, 0.25 µm

Column: DB-CLP2
123-8336
30 m x 0.32 mm, 0.25 µm

Instrument: Agilent 7890 GC with dual µECD

Carrier: Helium at 43.5 cm/s (constant flow)

Oven: 80 °C (hold 0.5 min) to 150 °C at 20 °C/min,
5 °C/min to 235 °C, 15 °C/min to 300 °C, hold 5 min


Sampler: Agilent 7693

Injection: 2 µL, splitless

Detector: µECD at 325 °C

Sample: 50 ng/mL 8081B analytes

- | | |
|----------------------------------|---------------------------------|
| 1. 1,2-Dibromo-3-chloropropane | 24. Heptachlor epoxide |
| 2. Hexachlorocyclopentadiene | 25. γ-Chlordane |
| 3. 1-Bromo-2-nitrobenzene | 26. trans-Nonachlor |
| 4. Etridiazole | 27. α-Chlordane |
| 5. Chloroneb | 28. Endosulfan I |
| 6. Trifluralin | 29. 4,4'-DDE |
| 7. TCMX | 30. Dieldrin |
| 8. Propachlor | 31. Chlorobenzilate (250 ng/mL) |
| 9. Di-allate isomers (250 ng/mL) | 32. Perthane (250 ng/mL) |
| 10. Hexachlorobenzene | 33. Chloropropylate (250 ng/mL) |
| 11. α-BHC | 34. Endrin |
| 12. Pentachloronitrobenzene | 35. Nitrofen |
| 13. γ-BHC | 36. 4,4'-DDD |
| 14. β-BHC | 37. Endosulfan II |
| 15. Heptachlor | 38. 4,4'-DDT |
| 16. Dichloro | 39. Endrin aldehyde |
| 17. Alachlor | 40. Endosulfan sulfate |
| 18. δ-BHC | 41. Captafol |
| 19. Chlorothalonil | 42. Methoxychlor |
| 20. Aldrin | 43. Endrin ketone |
| 21. DCPA | 44. Mirex |
| 22. Isodrin | 45. cis-Permethrin |
| 23. Kelthane | 46. trans-Permethrin |
| | 47. Decachlorobiphenyl |

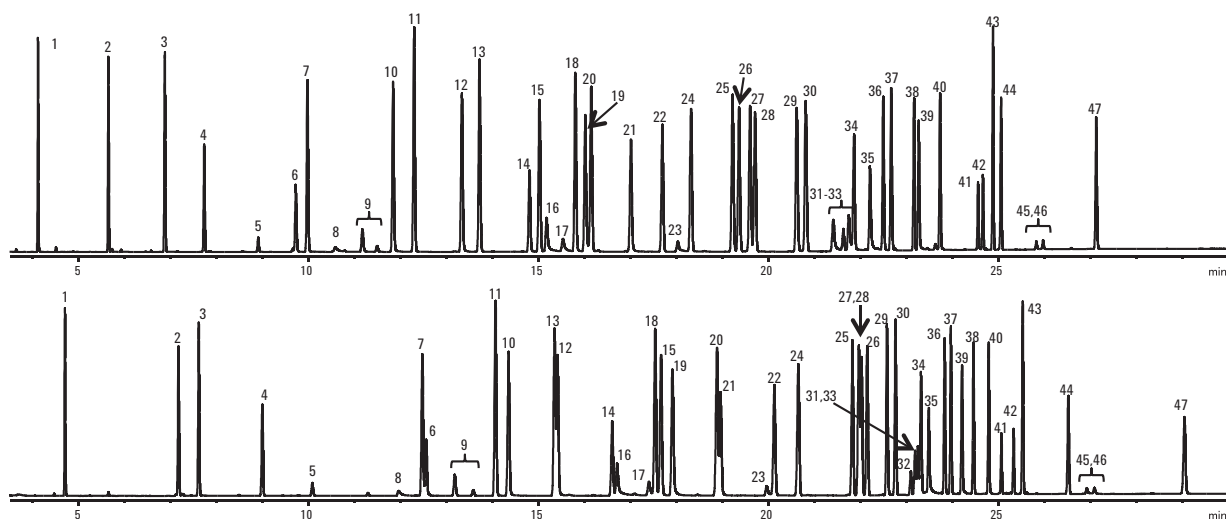

DB-624UI Organic Acid Performance

Column: DB-624 Ultra Inert
123-1334UI
30 m x 0.32 mm, 1.80 µm

Column: Non-Agilent 624, 30 m x 0.32 mm, 1.8 µm

Carrier: Hydrogen, 4 mL/min constant flow


Oven: 70 °C (1 min), then 20 °C/min to 260 °C

Inlet: 250 °C, 1 µL, split 1:200

Inlet liner: 4 mm, glass wool

Detector: FID at 260 °C

Organic acid mix C₁-C₁₀ (6 to 17 ng) on a DB-624UI column (A) and a traditional non-Agilent 624 column (R) after conditioning at 260 °C for 1 h.

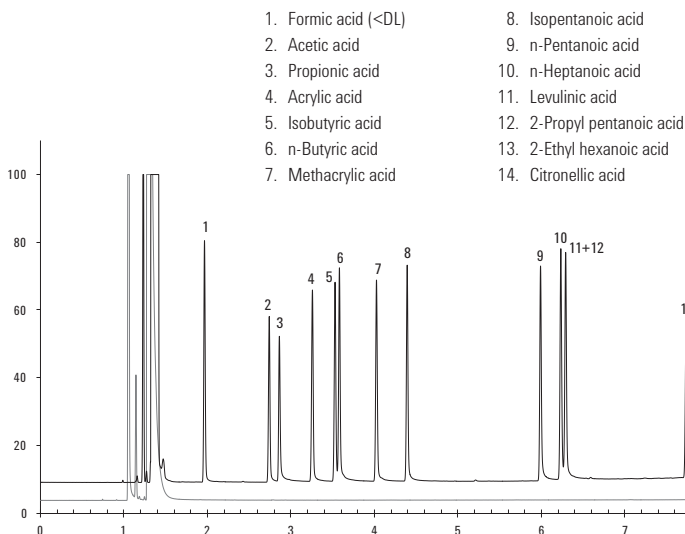
**EPA Method 551 – Chlorinated Solvents,
Trihalomethanes (THMs),
and Disinfection Byproducts (DBPs)**

Column: DB-CLP1
123-8232
30 m x 0.32 mm, 0.25 µm

Carrier: Helium, constant flow, 45 cm/s


Oven: 35 °C, hold 5.75 min, 20 °C/min to 95 °C,
40 °C/min to 200 °C, hold 1.25 min

Column: DB-CLP2
123-8336
30 m x 0.32 mm, 0.50 µm


Injection: 2 µL splitless, 200 °C

Detector: µECD, 300 °C

DB-CLP1

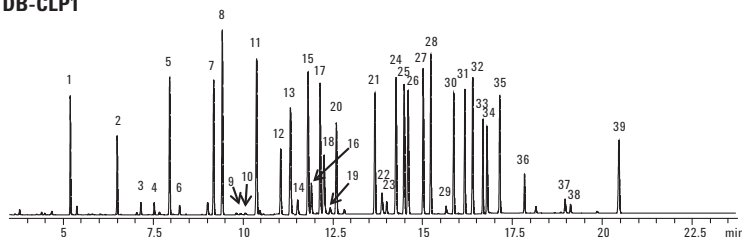

DB-CLP2

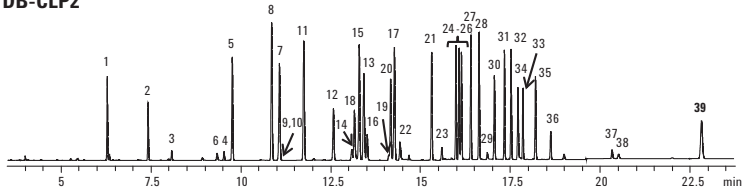

1. Chloroform
2. 1,1,1-Trichloroethane
3. Carbon tetrachloride
4. Trichloroacetonitrile
5. Trichloroethane
6. Chloral hydrate
7. Bromodichloromethane
8. 1,1-Dichloro-2-propanone
9. Dichloroacetonitrile
10. Chloropicrin
11. Tetrachloroethane
12. 1,1,2-Trichloroethane
13. Dibromochloromethane
14. 1,2-Dibromoethane
15. 1,1,1-Trichloro-2-propanone
16. Bromochloroacetonitrile
17. Bromoform
18. 1,2,3-Trichloropropane
19. Dibromoacetonitrile
20. 1,2-Dibromo-3-chloropropane

Analysis of Semivolatiles

Column A: DB-5.625
122-5632
30 m x 0.25 mm, 0.50 μ m

Column B: DB-5.625
121-5622
20 m x 0.18 mm, 0.36 μ m


Carrier: He constant flow mode, 1.1 mL/min

Oven: 40 °C (1 min), 25 °C/min to 320 °C
4.80 min hold


Injection: Splitless 0.5 μ L injected at 300 °C,
QuickSwap pressure 5.0 psi during acquisition,
80.0 psi during backflush with inlet set to
1.0 psi during backflush

Detector: Agilent 5975C Performance Turbo MSD
equipped with 6 mm large-aperture drawout lens,
p/n G2589-20045

Translating 0.25 mm id column method to 0.18 mm id format
results in 32% reduction in analysis time. Resolution of 77 peaks
of interest is also maintained for the faster 0.18 mm id separation.

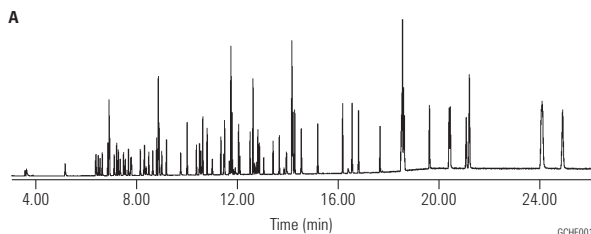

US EPA Method 8270, 5 ng/mL System Performance Check Compounds
Chromatogram using a DB-5.625, 30 m x 0.25 mm, 0.5 μ m

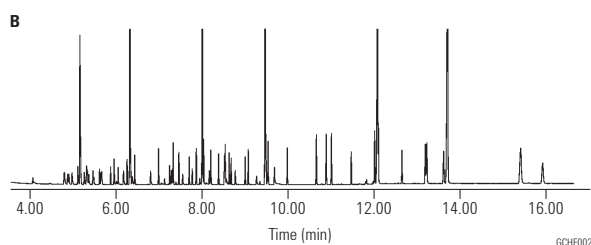

US EPA Method 8270, 5 ng/mL System Performance Check Compounds
Chromatogram using a DB-5.625, 20 m x 0.18 mm, 0.36 μ m

**TIPS & TOOLS**

Learn more about the Agilent 7890B GC System at www.agilent.com/chem/7890BGC



Pesticides, EPA 508.1

Column: DB-35ms
123-3832
30 m x 0.32 mm, 0.25 µm

Column: DB-XLB
123-1236
30 m x 0.32 mm, 0.50 µm

Carrier: Helium at 45 cm/s (EPC in constant flow mode)

Oven: 75 °C for 0.5 min
75-300 °C at 10 °C/min
300 °C for 2 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: µECD, 350 °C
Nitrogen makeup gas
(column + makeup flow = 30 mL/min constant flow)



Sample: 50 pg per component

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

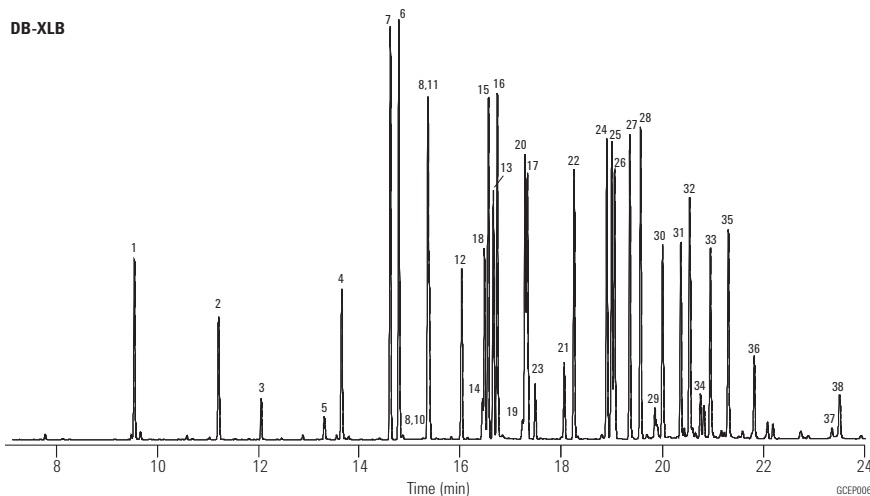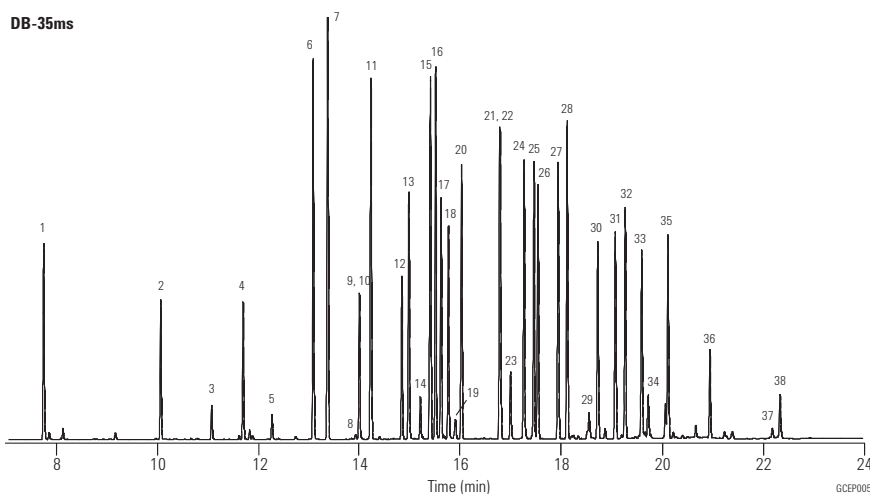

1. Hexachlorocyclopentadiene
2. Etridiazole
3. Chloroneb
4. Trifluralin
5. Propachlor
6. Hexachlorobenzene
7. α-BHC
8. Atrazine
9. Pentachloronitrobenzene
10. Simazine
11. γ-BHC
12. β-BHC
13. Heptachlor
14. Alachlor
15. δ-BHC
16. Chlorothalonil
17. Aldrin
18. Metribuzin
19. Metolachlor
20. DCPA
21. 4,4'-Dibromobiphenyl
22. Heptachlor epoxide
23. Cyanazine
24. γ-Chlordane
25. α-Chlordane
26. Endosulfan I
27. 4,4'-DDE
28. Dieldrin
29. Chlorobenzilate
30. Endrin
31. 4,4'-DDD
32. Endosulfan II
33. 4,4'-DDT
34. Endrin aldehyde
35. Endosulfan sulfate
36. Methoxychlor
37. cis-Permethrin
38. trans-Permethrin

Phenoxy Acid Herbicides – Methyl Derivatives, EPA 8151A

Column: DB-35ms
123-3832
30 m x 0.32 mm, 0.25 µm

Carrier: Helium at 45 cm/s (EPC in constant flow mode)

Oven: 50 °C for 0.5 min
50-100 °C at 25 °C/min
100-320 °C at 12 °C/min
320 °C for 2 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: µECD, 350 °C
Nitrogen makeup gas
(column + makeup flow = 30 mL/min constant flow)


Sample: 50 pg per component

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

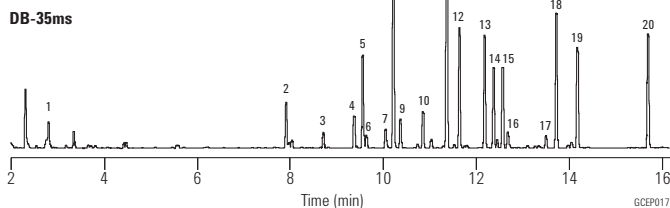

1. Dalapon
2. 3,5-Dichlorobenzoic acid
3. 4-Nitrophenol
4. Methyl-2,4-dichlorophenylacetate (SS)
5. Dicamba
6. MCPP
7. MCPA
8. 4,4'-Dibromooctafluorobiphenyl (IS)
9. Dichloroprop
10. 2,4-D
11. Pentachlorophenol
12. 2,4,5-T,P
13. 2,4,5-T
14. Chloramben
15. Dinoseb
16. 2,4-DB
17. Bentazone
18. DCPA
19. Picloram
20. Acifluorfen

Direct Comparison for Rapid CLP (Contract Laboratory Program) Pesticide Analysis

Column: DB-17ms
121-4722
20 m x 0.18 mm, 0.18 µm

Column: DB-XLB
121-1222
20 m x 0.18 mm, 0.18 µm


Carrier: Hydrogen (69 cm/s at 120 °C,
ramped at 99 mL/min to
106 cm/s at 4.4 min)

Oven: 120 °C (0.32 min); 120 °C/min to 160 °C;
30 °C/min to 258 °C (0.18 min);
38.81 °C/min to 300 °C (1.5 min)


Injection: Split/splitless, 220 °C, pulsed splitless
(35 psi for 0.5 min, purge flow of 40 mL/min
on at 1 min, gas saver flow
20 mL/min on 3 min)

Detector: µECD 320 °C; nitrogen makeup;
constant column + makeup flow 60 mL/min



**DB-17ms primary column
DB-XLB confirmatory column**

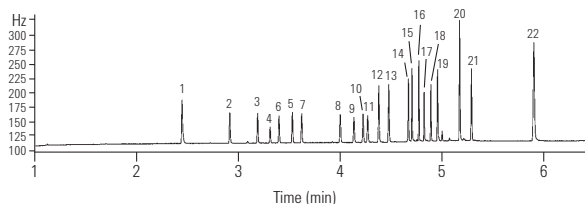

1. Tetrachloro-m-xylene
2. α-BHC
3. γ-BHC
4. β-BHC
5. δ-BHC
6. Heptachlor
7. Aldrin
8. Heptachlor epoxide
9. γ-Chlordane
10. α-Chlordane
11. Endosulfan I
12. 4,4' DDE
13. Dieldrin
14. Endrin
15. 4,4' DDD
16. Endosulfan II
17. 4,4' DDT
18. Endrin aldehyde
19. Endosulfan sulfate
20. Methoxychlor
21. Endrin ketone
22. Decachlorobiphenyl

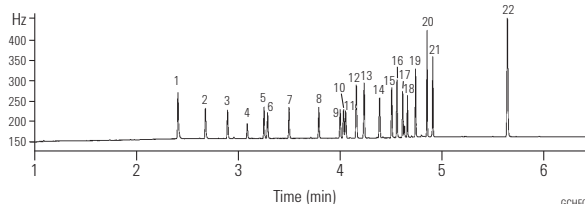

**Vendor R primary column, 20 m x 0.18 mm, 0.18 µm
Vendor R confirmatory column, 20 m x 0.18 mm, 0.14 µm**

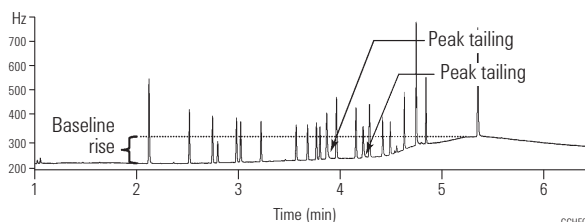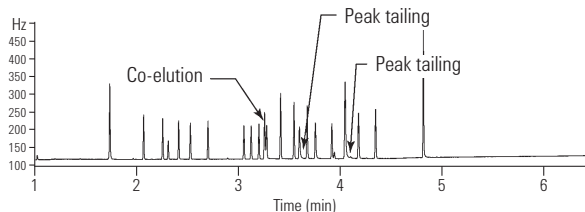

The DB-17ms primary column and DB-XLB confirmatory column sufficiently resolved all the peaks of interest in less than six minutes with sharp, symmetrical peaks and minimal baseline drift. In contrast, vendor R's primary analysis column resolved only 20 of 22 peaks with visible peak tailing. Vendor R's confirmatory column resolved all 22 peaks of interest but with peak tailing and an unacceptable level of temperature dependent baseline drift.

Aroclors 1016-1268 (without 1221)

Column: DB-XLB
121-1232
30 m x 0.18 mm, 0.18 µm


Carrier: Helium at 37 cm/s, measured at 150 °C

Oven: 100 °C for 1 min
100-265 °C at 1.2 °C/min

Injection: Hot on-column, 250 °C

Detector: MSD, 340 °C transfer line, SIM

Sample: 1 µL in isooctane, 12.5 ppm

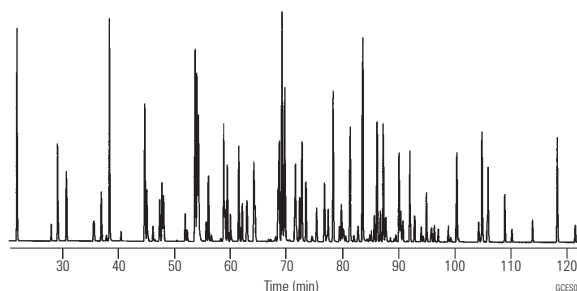**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

CLP Pesticides

Column: DB-35ms
123-3832
30 m x 0.32 mm, 0.25 µm

Column: DB-XLB
123-1236
30 m x 0.32 mm, 0.50 µm


Carrier: Helium at 45 cm/s
(EPC in constant flow mode)

Oven: 110 °C for 0.5 min
110-320 °C at 15 °C/min
320 °C for 2 min


Injection: Splitless, 250 °C
30 s purge activation time

Detector: µECD, 350 °C
Nitrogen makeup gas
(column + makeup flow =
30 mL/min constant flow)

Sample: 50 pg per component

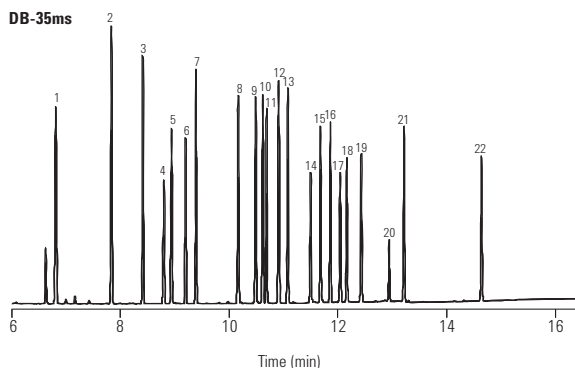

**Complete resolution and confirmation of
22 CLP Pesticides in under 16 minutes!**

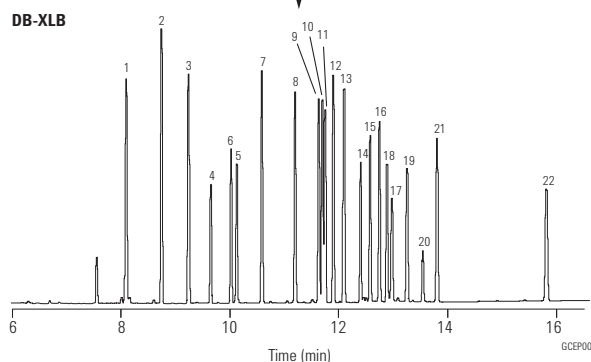

1. Tetrachloro m-xylene (SS)
2. α-BHC
3. γ-BHC
4. β-BHC
5. Heptachlor
6. δ-BHC
7. Aldrin
8. Heptachlor epoxide
9. γ-Chlordane
10. α-Chlordane
11. Endosulfan I
12. 4,4'-DDE
13. Dieldrin
14. Endrin
15. 4,4'-DDD
16. Endosulfan II
17. 4,4'-DDT
18. Endrin aldehyde
19. Endosulfan sulfate
20. Methoxychlor
21. Endrin ketone
22. Decachlorobiphenyl (SS)

SS - Surrogate Standard

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759


Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

High Speed VOC, EPA Method 8260**Column:** DB-VRX
121-1524**20 m x 0.18 mm, 1.00 µm****Carrier:** Helium at 55 cm/s (1.5 mL/min)**Injection:** Split, 150 °C
Split ratio 60:1**Oven:** 45 °C for 3.0 min
45-190 °C at 36 °C/min
190-225 °C at 20 °C/min
225 °C for 0.5 min**Detector:** Agilent 5975 MSD
Scan range: 35-260 amu
Scan rate: 3.25 scans/s
Quad temp: 150 °C
Source temp: 200 °C
Transfer line temp: 200 °C**Sampler:** Purge and trap (Tekmar 3100)
Purge: 11 min
Trap: VoCarb 3000
Preheat: 245 °C
Desorb: 250 °C for 1 min
Bake: 260 °C for 10 min
Line & valve: 100 °C**Sample:** 5 mL

- Halogenated and aromatic analytes at 40 ppb
- Internal standards at 20 ppb
- Polar analytes (i.e., ethers, alcohols and ketones at 100-800 ppb)

- | | | |
|-------------------------------|-------------------------------|----------------------------------|
| 1. Dichlorodifluoromethane | 47. Carbon tetrachloride | 93. Propylbenzene |
| 2. Chloromethane | 48. Chloroacetonitrile | 94. 2-Chlorotoluene |
| 3. Hydroxypropionitrile | 49. Benzene | 95. 4-Chlorotoluene |
| 4. Vinyl chloride | 50. tert-Amylmethyl ether | 96. 1,3,5-Trimethylbenzene |
| 5. Bromomethane | 51. Fluorobenzene (IS) | 97. Pentachloroethane |
| 6. Chloroethane | 52. 2-Pentanone | 98. tert-Butylbenzene |
| 7. Ethanol | 53. Dibromomethane | 99. 1,2,4-Trimethylbenzene |
| 8. Acetonitrile | 54. 1,2-Dichloropropane | 100. sec-Butylbenzene |
| 9. Acrolein | 55. Trichloroethene | 101. 1,3-Dichlorobenzene |
| 10. Trichlorofluoromethane | 56. Bromodichloromethane | 102. Benzyl chloride |
| 11. Isopropyl alcohol | 57. 2-Nitropropane | 103. 1,4-Dichlorobenzene-d4 (IS) |
| 12. Acetone | 58. 1,4-Dioxane | 104. 1,4-Dichlorobenzene |
| 13. Ethyl ether | 59. Epichlorohydrin | 105. Isopropyltoluene |
| 14. 1,1-Dichloroethene | 60. Methyl methacrylate | 106. 1,2-Dichlorobenzene |
| 15. tert-Butyl alcohol | 61. cis-1,3-Dichloropropene | 107. Butylbenzene |
| 16. Acrylonitrile | 62. Propiolactone | 108. 1,2-Dibromo-3-chloropropane |
| 17. Methylene chloride | 63. Bromoacetone | 109. Hexachloroethane |
| 18. Allyl chloride | 64. Pyridine | 110. Nitrobenzene |
| 19. Allyl alcohol | 65. trans-1,3-Dichloropropene | 111. 1,2,4-Trichlorobenzene |
| 20. 1-Propanol | 66. 1,1,2-Trichloroethane | 112. Naphthalene |
| 21. Propargyl alcohol | 67. Toluene-d8 (IS) | 113. Hexachlorobutadiene |
| 22. trans-1,2-Dichloroethene | 68. Toluene | 114. 1,2,3-Trichlorobenzene |
| 23. MTBE | 69. 1,3-Dichloropropane | |
| 24. 1,1-Dichloroethane | 70. Paraldehyde | |
| 25. Propionitrile | 71. Ethyl methacrylate | |
| 26. 2-Butanone | 72. Dibromochloromethane | |
| 27. Diisopropyl ether | 73. 3-Chloropropionitrile | |
| 28. cis-1,2-Dichloroethene | 74. 1,2-Dibromoethane | |
| 29. Methacrylonitrile | 75. Tetrachloroethene | |
| 30. Bromochloromethane | 76. 1,1,1,2-Tetrachloroethane | |
| 31. Chloroform | 77. 1-Chlorohexane | |
| 32. 2,2-Dichloropropane | 78. Chlorobenzene | |
| 33. Ethyl acetate | 79. Ethylbenzene | |
| 34. Ethyl-tert-butyl ether | 80. Bromoform | |
| 35. Methyl acrylate | 81. m-Xylene | |
| 36. Dibromofluoromethane (IS) | 82. p-Xylene | |
| 37. Isobutanol | 83. trans-Dichlorobutene | |
| 38. Dichloroethane-d4 (IS) | 84. 1,3-Dichloro-2-propanol | |
| 39. Pentafluorobenzene | 85. Styrene | |
| 40. 1,2-Dichloroethane | 86. 1,1,2,2-Tetrachloroethane | |
| 41. 1,1,1-Trichloroethane | 87. o-Xylene | |
| 42. 1-Chlorobutane | 88. 1,2,3-Trichloropropane | |
| 43. Crotonaldehyde | 89. cis-Dichlorobutene | |
| 44. 2-Chloroethanol | 90. 4-Bromofluorobenzene (IS) | |
| 45. 1,1-Dichloropropene | 91. Isopropylbenzene | |
| 46. 1-Butanol | 92. Bromobenzene | |

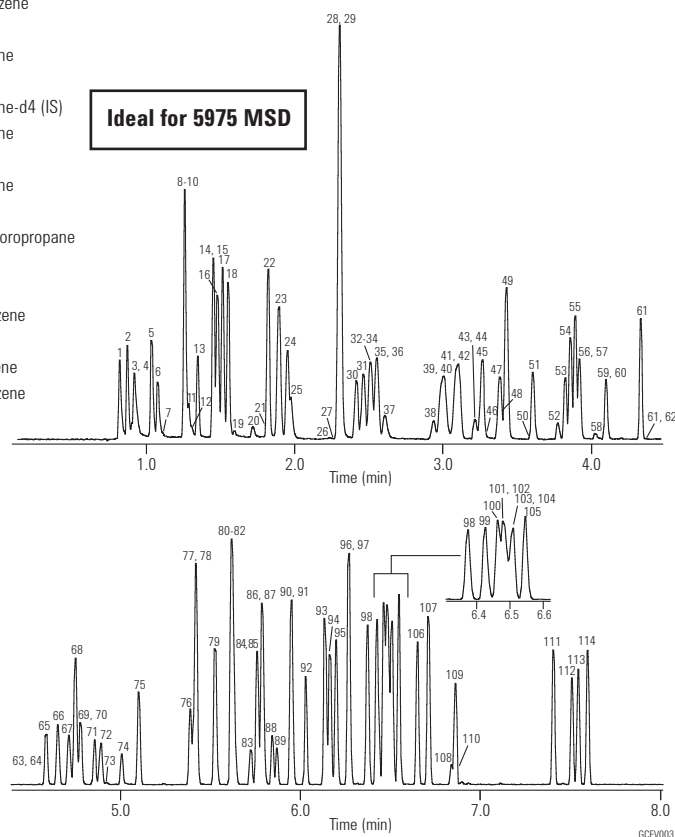**Suggested Supplies****Septum:** 11 mm Advanced Green septa, 5183-4759**Liner:** Direct, 1.5 mm id, 18740-80200**Seal:** Gold plated seal, 18740-20885

PBDEs

Column: DB-XLB
122-1231
30 m x 0.25 mm, 0.10 µm


Carrier: Helium at 38 cm/s at 100 °C (1.2 mL/min),
constant flow mode

Oven: 100 °C for 1 min; 100 °C to 340 °C at 20 °C/min,
340 °C for 12 min

Injection: Cool on-column, oven-track mode

Detector: Agilent 5973 MSD, 325 °C transfer line, EI SIM
(ions monitored: 231.8, 248.0, 327.9, 398.6, 400.5,
405.8, 845.7, 563.6, 643.5, 721.4, 799.3)

Sample: 0.5 µL

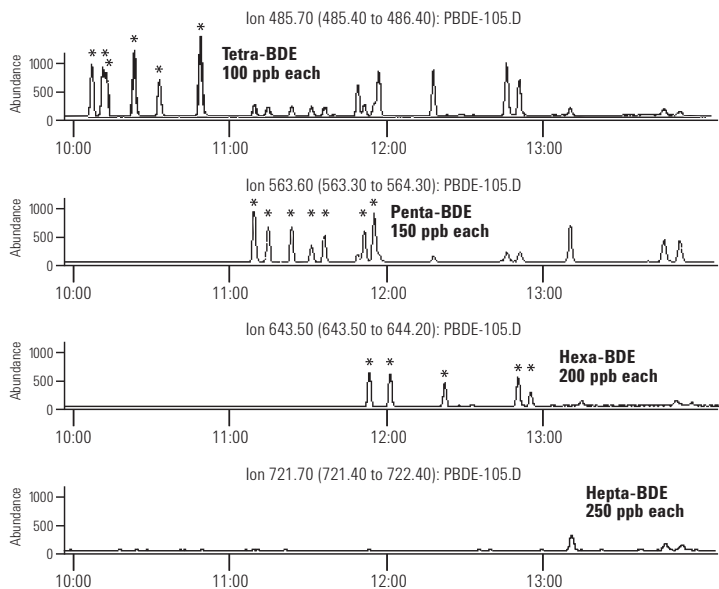

For a complete Application Note, visit www.agilent.com/chem, select "Literature" from the Library and type 5989-0094EN into the "Keyword" field.

EPA Volatiles by GC/MS (Split Injector)

Column: DB-VRX
122-1564
60 m x 0.25 mm, 1.40 µm


Carrier: Helium at 30 cm/s, measured at 45 °C

Oven: 45 °C for 10 min
45-190 °C at 12 °C/min
190 °C for 2 min
190-225 °C at 6 °C/min
225 °C for 1 min

Sampler: Purge and trap (O.I.A. 4560)
Purge: Helium for 11 min at 40 mL/min
Trap: Tenax/Silica Gel/Carbosieve
Preheat: 175 °C
Desorb: 220 °C for 0.6 min

Injection: Split, 110 °C
Split flow 30 mL/min

Detector: MSD, 235 °C transfer line
Full scan 35-260 amu (m/z 44 subtracted)

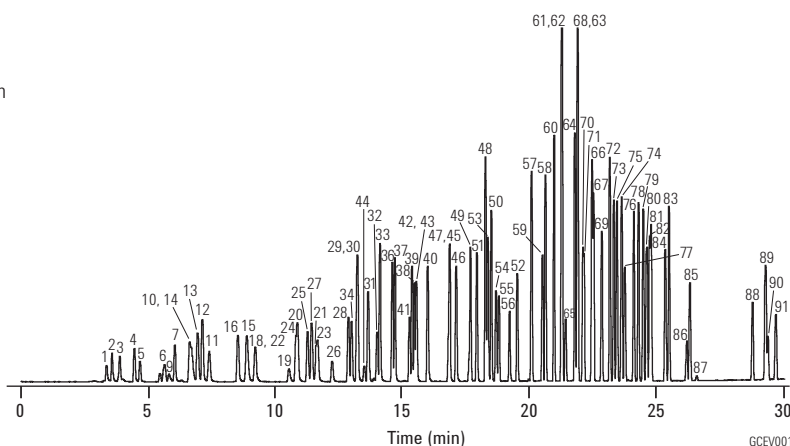Suggested Supplies**Septum:** 11 mm Advanced Green septa, 5183-4759**Liner:** Direct, 1.5 mm id, 18740-80200**Seal:** Gold plated seal kit, 5188-5367

GCEV001

- | | | |
|------------------------------|-----------------------------------|---------------------------------|
| 1. Dichlorodifluoromethane | 32. Carbon tetrachloride | 63. o-Xylene |
| 2. Chloromethane | 33. Benzene | 64. Styrene |
| 3. Vinyl chloride | 34. 1,2-Dichloroethane | 65. Bromoform |
| 4. Bromomethane | 35. 2,2-Dimethylhexane | 66. Isopropylbenzene |
| 5. Chloroethane | 36. Fluorobenzene (IS) | 67. 4-Bromofluorobenzene (SS) |
| 6. Trichlorofluoromethane | 37. 1,4-Difluorobenzene (IS) | 68. 1,1,2,2-Tetrachloroethane |
| 7. Diethyl ether | 38. Trichloroethene | 69. Bromobenzene |
| 8. 1,1-Dichloroethene | 39. 1,2-Dichloropropane | 70. 1,2,3-Trichloropropane |
| 9. Acetone | 40. Methyl methacrylate | 71. trans-1,4-Dichloro-2-butene |
| 10. Iodomethane | 41. Dibromomethane | 72. n-Propylbenzene |
| 11. Carbon disulfide | 42. Bromodichloromethane | 73. 2-Chlorotoluene |
| 12. Allyl chloride | 43. 2-Nitropropane | 74. 1,3,5-Trimethylbenzene |
| 13. Methylene chloride | 44. Chloroacetonitrile | 75. 4-Chlorotoluene |
| 14. Acrylonitrile | 45. cis-1,3-Dichloropropene | 76. tert-Butylbenzene |
| 15. Methyl-tert-butyl ether | 46. 4-Methyl-2-pentanone | 77. Pentachloroethane |
| 16. trans-1,2-Dichloroethene | 47. 1,1-Dichloro-2-propanone | 78. 1,2,4-Trimethylbenzene |
| 17. Hexane | 48. Toluene | 79. sec-Butylbenzene |
| 18. 1,1-Dichloroethane | 49. trans-1,3-Dichloropropene | 80. 1,3-Dichlorobenzene |
| 19. 2-Butanone | 50. Ethyl methacrylate | 81. p-Isopropyltoluene |
| 20. cis-1,2-Dichloroethene | 51. 1,1,2-Trichloroethane | 82. 1,4-Dichlorobenzene |
| 21. 2,2-Dichloropropane | 52. Tetrachloroethene | 83. n-Butylbenzene |
| 22. Propionitrile | 53. 1,3-Dichloropropane | 84. 1,2-Dichlorobenzene |
| 23. Methyl acrylate | 54. 2-Hexanone | 85. Hexachloroethane |
| 24. Methacrylonitrile | 55. Dibromochloromethane | 86. 1,2-Dibromo-3-chloropropane |
| 25. Bromochloromethane | 56. 1,2-Dibromoethane | 87. Nitrobenzene |
| 26. Tetrahydrofuran | 57. 1-Chloro-3-fluorobenzene (IS) | 88. 1,2,4-Trichlorobenzene |
| 27. Chloroform | 58. Chlorobenzene | 89. Hexachlorobutadiene |
| 28. Pentafluorobenzene (IS) | 59. 1,1,1,2-Tetrachloroethane | 90. Naphthalene |
| 29. 1,1,1-Trichloroethane | 60. Ethylbenzene | 91. 1,2,3-Trichlorobenzene |
| 30. 1-Chlorobutane | 61. m-Xylene | |
| 31. 1,1-Dichloropropene | 62. p-Xylene | |

EPA Method 525.2

Column: DB-5ms
122-5532
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 32 cm/s, measured at 45 °C, constant flow mode

Oven: 45 °C for 1 min
45-130 °C at 30 °C/min
130 °C for 3 min
130-180 °C at 12 °C/min
180-240 °C at 7 °C/min
240-325 °C at 12 °C/min
325 °C for 5 min

Injection: Splitless, 300 °C
1.0 min purge activation time
Focus liner

Detector: MSD, 325 °C transfer line
Full scan m/z 45-450

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Sample: Composite mixture of AccuStandard
Method 525.2 standards (M-525.2-SV-ASL,
M-525.2-FS-ASL, M-525.2-CP-ASL,
M-525.2-NP1-ASL, M-525.2-NP2-ASL):
target compounds at 2 ng/µL, IS/SS at 5 ng/µL

Compound	RT	m/z	Compound	RT	m/z	Compound	RT	m/z
1. Isophorone	5.85	82	49. 2,4,5-Trichlorobiphenyl	15.59	256	84. DEF	19.84	57/169
2. 1,3-Dimethyl-2-nitrobenzene (SS)	6.65	134	50. Metribuzin	15.95	198	85. 2,2',4,4',5,6'-Hexachlorobiphenyl	19.90	360
3. Dichlorvos	7.41	109	51. Alachlor	16.14	160	86. Dieldrin	19.92	79
4. Hexachlorocyclopentadiene	8.87	237	52. Simetryn	16.23	213	87. Carboxin	19.97	143
5. EPTC	9.17	128	53. Ametryn	16.33	227/170	88. Endrin	20.43	67/81
6. Mevinphos	10.09	127	54. Heptachlor	16.36	100	89. Chlorobenzilate	20.56	139
7. Butylate	10.18	57/146	55. Prometryn	16.40	241/184	90. Endosulfan II	20.68	195
8. Vernolate	10.42	128	56. Prebane (terbutryn)	16.72	226/185	91. p,p'-DDD	20.77	235/165
9. Dimethyl phthalate	10.45	163	57. Bromacil	16.79	205	92. Endrin aldehyde	21.01	67
10. Terrazole (etridiazole)	10.47	211/183	58. Di-n-butyl phthalate	16.90	149	93. Norflurazon	21.36	145
11. 2,6-Dinitrotoluene	10.56	165	59. 2,2',4,4'-Tetrachlorobiphenyl	17.02	292	94. Benzyl butyl phthalate	21.49	149
12. Tillam (pebulate)	10.61	128	60. Metolachlor	17.11	162	95. Endosulfan sulfate	21.53	272
13. Acenaphthylene	10.65	152	61. Dursban (chlorpyrifos)	17.15	197/97	96. p,p'-DDT	21.61	235/165
14. Acenaphthene-d10 (IS)	11	164	62. Cyanazine	17.23	225/68	97. Hexazinone	21.68	171
15. Chloroneb	11.17	191	63. Dacthal (DCPA methyl ester)	17.27	301	98. Bis(2-ethylhexyl) adipate	21.87	129
16. 2-Chlorobiphenyl	11.19	188	64. Aldrin	17.29	66	99. Triphenylphosphate (SS)	21.98	326/325
17. Tebuthiuron	11.37	156	65. Triadimefon	17.43	57	100. Endrin ketone (breakdown product)	22.52	67/317
18. 2,4-Dinitrotoluene	11.51	165	66. Diphenamid	17.73	72/167	101. 2,2',3,3',4,4',6-Heptachlorobiphenyl	22.59	394/396
19. Molinate	11.68	126	67. MGK-264 (isomer A)	17.78	164/66	102. Benz[a]anthracene	22.66	228
20. Diethyl phthalate	12.21	149	68. MGK-264 (isomer B)	18.11	164	103. Chrysene-d12 (IS)	22.68	240
21. Fluorene	12.35	166	69. Heptachlor epoxide	18.28	81	104. 2,2',3,3',4,4',5,6'-Octachlorobiphenyl	22.70	430/428
22. Propachlor	12.46	120	70. 2,2',3',4,6-Pentachlorobiphenyl	18.34	326	105. Methoxychlor	22.73	227
23. Ethoprop	12.82	158	71. Merphos	18.36	209/153	106. Chrysene	22.74	228
24. Cycloate	12.86	83/154	72. γ-Chlordane	18.88	373	107. Bis(2-ethylhexyl) phthalate	23.10	149
25. Chlorpropham	13.08	127	73. Tetrachlorvinphos (stirifos)	18.95	109	108. Fenarimol	23.80	139
26. Trifluralin	13.14	306	74. Butachlor	19.03	176/160	109. cis-Permethrin	24.38	183
27. α-BHC	13.69	181	75. Pyrene-d10 (SS)	19.13	212	110. trans-Permethrin	24.50	183
28. 2,3-Dichlorobiphenyl	13.74	222/152	76. Pyrene	19.18	202	111. Benzo[b]fluoranthene	25.06	252
29. Hexachlorobenzene	13.77	284	77. α-Chlordane	19.21	375/373	112. Benzo[k]fluoranthene	25.12	252
30. Gesatamine (atraton)	13.99	196/169	78. Endosulfan I	19.22	195	113. Fluridone	25.66	328
31. Prometon	14.14	225/168	79. trans-Nonachlor	19.28	409	114. Benzo[a]pyrene	25.67	252
32. Atrazine	14.26	200/215	80. Fenamiphos	19.33	303/154	115. Perylene-d12 (SS)	25.78	264
33. Simazine	14.27	201/186	81. Napropamide	19.39	72	116. Indeno[1,2,3-c,d]pyrene	27.63	276
34. β-BHC	14.28	181	82. Tricyclazole	19.61	189	117. Dibenzo[a,h]anthracene	27.69	278
35. Pentachlorophenol	14.35	266	83. p,p'-DDE	19.76	246	118. Benzo[g,h,i]perylene	28.11	276
36. Propazine	14.35	214/172						
37. γ-BHC	14.52	181						
38. Terbufos	14.62	57						
39. Pronamide	14.69	173						
40. Diazinon	14.76	137/179						
41. Phenanthrene-d10 (IS)	14.85	188						
42. Chlorothalonil	14.89	266						
43. Phenanthrene	14.92	178						
44. Terbacil	15.02	161						
45. Methyl paraoxon	15.04	109						
46. Disulfoton	15.05	88						
47. Anthracene	15.06	178						
48. δ-BHC	15.20	181						

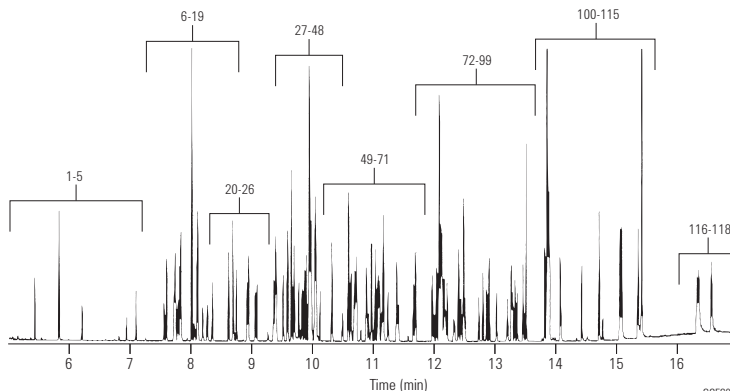

Pesticides and Fire Retardants (US EPA 527)

Column: DB-5ms Ultra Inert
122-5532UI
30 m x 0.25 mm, 0.25 µm


Carrier: Helium, 52 cm/s, constant flow

Oven: 60 °C (1 min) to 210 °C (25 °C/min), 20 °C/min to 310 °C (3 min)

Injection: Splitless, 250 °C, purge flow 50 mL/min at 1 min,
gas saver 80 mL/min on at 3 min

Detector: Transfer line 290 °C, source 300 °C, quad 180 °C

Sample: Pesticide/PBDE standards, 1 ng with 5 ng IS/SS on-column

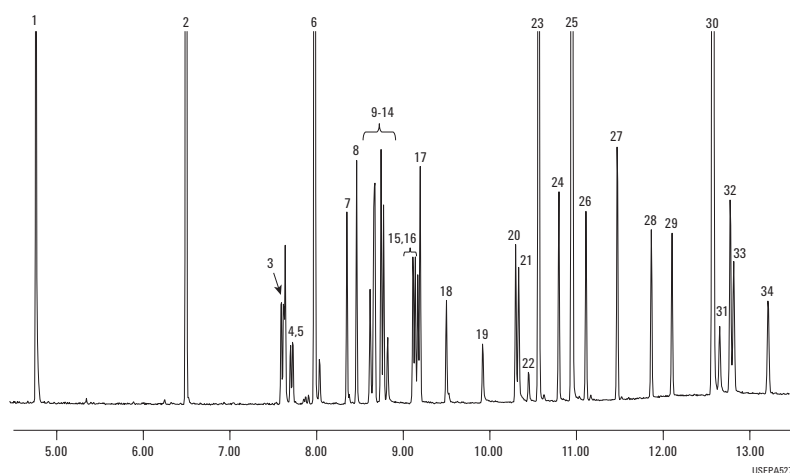

- | | |
|--------------------------------|-------------------------|
| 1. 1,2-Dimethyl-2-nitrobenzene | 18. Fenamiphos |
| 2. Acenaphthalene-D10 | 19. Nitrophen |
| 3. Dimethoate | 20. Norflurazon |
| 4. Atrazine | 21. Kepone |
| 5. Propazine | 22. Hexazinone |
| 6. Anthracene-D10 | 23. Triphenyl phosphate |
| 7. Vinclozoline | 24. Bifenthrin |
| 8. Prometryn | 25. Chrysene-D12 |
| 9. Bromacil | 26. BDE-47 |
| 10. Malathion | 27. Mirex |
| 11. Thiazopyr | 28. BDE-100 |
| 12. Dursban | 29. BDE-99 |
| 13. Benthicarb | 30. Perylene-D12 |
| 14. Parathion | 31. Fenvalerate |
| 15. Terbufos sulfone | 32. Esfenvalerate |
| 16. Bioallethrin | 33. Hexabromobiphenyl |
| 17. Oxychlorane | 34. BDE-153 |

**EPA Method 508.1 –
Chlorinated Pesticides and Herbicides**

Column: DB-CLP1
123-8232
30 m x 0.32 mm, 0.25 µm

Column: DB-CLP2
123-8336
30 m x 0.32 mm, 0.50 µm

Carrier: Helium, constant flow, 35 cm/s

Oven: 80 °C, hold 0.5 min, 26 °C/min to 175 °C, 6.5 °C/min to 235 °C, 15 °C/min to 300 °C, hold 6 min


Injection: 2 µL, splitless, 250 °C

Detector: µCED, 340 °C


Sample: 100 ng/mL EPA 508.1 analytes,
100 ng/mL pesticide surrogate mix

- | | |
|--|---|
| 1. Hexachlorocyclopentadiene | 20. DCPA |
| 2. Etridiazole | 21. Heptachlor epoxide |
| 3. Chloroneb | 22. Cyanazine |
| 4. Trifluralin | 23. Butachlor |
| 5. Tetrachloro-m-xylene (surrogate standard) | 24. γ-Chlordane |
| 6. Propachlor | 25. α-Chlordane |
| 7. Hexachlorobenzene | 26. Endosulfan I |
| 8. α-BHC | 27. 4,4'-DDE |
| 9. Atrazine | 28. Dieldrin |
| 10. Simazine | 29. Chlorobenzilate |
| 11. γ-BHC | 30. Endrin |
| 12. β-BHC | 31. 4,4'-DDD |
| 13. Heptachlor | 32. Endosulfan II |
| 14. Alachlor | 33. 4,4'-DDT |
| 15. δ-BHC | 34. Endrin aldehyde |
| 16. Chlorothalonil | 35. Endosulfan sulfate |
| 17. Aldrin | 36. Methoxychlor |
| 18. Metribuzin | 37. cis-Permethrin |
| 19. Metolachlor | 38. trans-Permethrin |
| | 39. Decachlorobiphenyl (surrogate standard) |

DB-CLP1

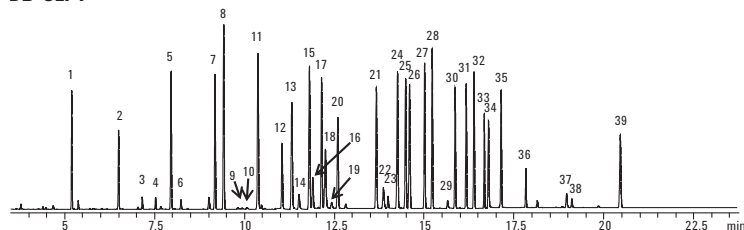

DB-CLP2

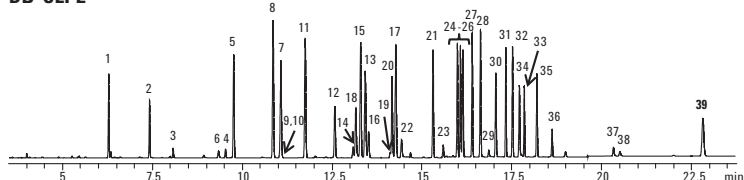

The DB-CLP1 column separates all chlorinated pesticide and herbicide analytes according to EPA Method 505.

Chlorinated Pesticides, EPA Method 508

Column: HP-5ms
19091S-433
30 m x 0.25 mm, 0.25 µm

Carrier: Helium, 24 psi, 45 cm/s (80 °C) constant flow


Oven: 80 °C for 1 min
80-180 °C at 30 °C/min
180-205 °C at 3 °C/min
205 °C for 4 min
205-290 °C at 2 °C/min
290 °C for 2 min

Injection: Splitless
1 min purge delay

Detector: ECD, 320 °C
Nitrogen makeup gas at 60 mL/min
Anode purge 3 mL/min

Sample: 1 µL

- | | | |
|----------------------|------------------------|------------------------|
| 1. Etridiazole | 11. Heptachlor | 21. Endosulfan II |
| 2. Chloroneb | 12. Aldrin | 22. Chlorobenzilate |
| 3. Propachlor | 13. DCPA | 23. 4,4'-DDD |
| 4. Trifluralin | 14. Heptachlor epoxide | 24. Endrin aldehyde |
| 5. α-BHC | 15. γ-Chlordane | 25. Endosulfan sulfate |
| 6. Hexachlorobenzene | 16. Endosulfan I | 26. 4,4'-DDT |
| 7. β-BHC | 17. α-Chlordane | 27. Methoxychlor |
| 8. δ-BHC | 18. Dieldrin | 28. cis-Permethrin |
| 9. γ-BHC | 19. 4,4'-DDE | 29. trans-Permethrin |
| 10. Chlorothalonil | 20. Endrin | |

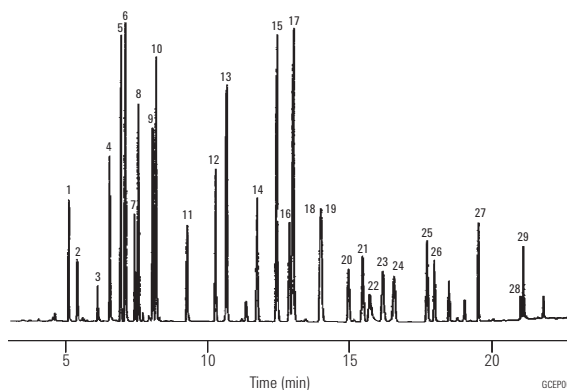**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Organochlorine Pesticides

Column: DB-5
125-5037
30 m x 0.53 mm, 0.50 µm

Carrier: Helium at 30 cm/s (4.0 mL/min)


Oven: 150-275 °C at 4 °C/min
275 °C for 30 min

Injection: Splitless, 250 °C

Detector: ECD, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.7 µL of 100 pg/µL standard in isooctane

- | | |
|--------------------------------------|-----------------------------|
| 1. 2,4,5,6-Tetrachloro-m-xylene (IS) | 12. Dieldrin |
| 2. α-BHC | 13. p,p'-DDE |
| 3. β-BHC | 14. Endrin |
| 4. γ-BHC | 15. Endosulfan II |
| 5. δ-BHC | 16. p,p'-DDD |
| 6. Heptachlor | 17. Endrin aldehyde |
| 7. Aldrin | 18. Endosulfan sulfate |
| 8. Heptachlor epoxide | 19. p,p'-DDT |
| 9. γ-Chlordane | 20. Endrin ketone |
| 10. Endosulfan I | 21. Methoxychlor |
| 11. α-Chlordane | 22. Decachlorobiphenyl (IS) |

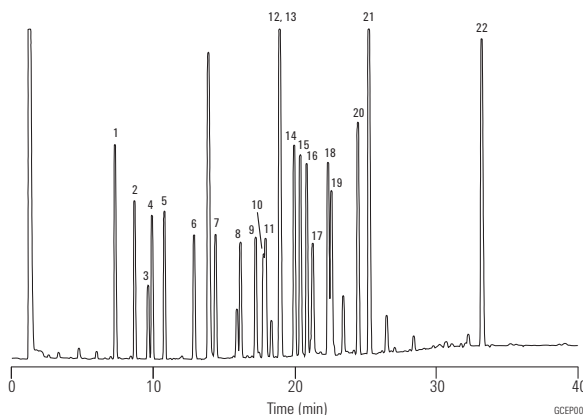**Suggested Supplies**

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Septum: 11 mm Advanced Green septa, 5183-4759

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Organochlorine Pesticides III

Column: DB-1701
125-0737
30 m x 0.53 mm, 0.50 µm

Carrier: Helium at 30 cm/s (4.0 mL/min)


Oven: 150-275 °C at 4 °C/min
275 °C for 30 min

Injection: Splitless, 250 °C

Detector: ECD, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.7 µL of 100 pg/µL standard in isooctane

- | | |
|--------------------------------------|-----------------------------|
| 1. 2,4,5,6-Tetrachloro-m-xylene (IS) | 12. Dieldrin |
| 2. α-BHC | 13. p,p'-DDE |
| 3. β-BHC | 14. Endrin |
| 4. γ-BHC | 15. Endosulfan II |
| 5. δ-BHC | 16. p,p'-DDD |
| 6. Heptachlor | 17. Endrin aldehyde |
| 7. Aldrin | 18. Endosulfan sulfate |
| 8. Heptachlor epoxide | 19. p,p'-DDT |
| 9. γ-Chlordane | 20. Endrin ketone |
| 10. Endosulfan I | 21. Methoxychlor |
| 11. α-Chlordane | 22. Decachlorobiphenyl (IS) |

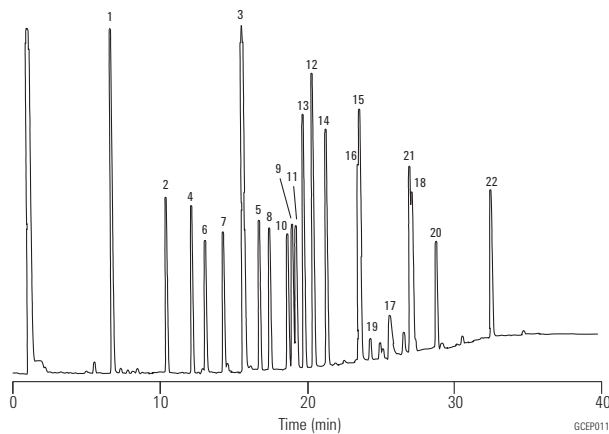**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Organochlorine Pesticides IV

Column: DB-35
125-1937
30 m x 0.53 mm, 0.50 µm

Carrier: Helium at 30 cm/s (4.0 mL/min)

Oven: 150-275 °C at 4 °C/min
275 °C for 30 min

Injection: Splitless, 250 °C

Detector: ECD, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.7 µL of 100 pg/µL standard in isooctane


- | | |
|--------------------------------------|-----------------------------|
| 1. 2,4,5,6-Tetrachloro-m-xylene (IS) | 12. Dieldrin |
| 2. α-BHC | 13. p,p'-DDE |
| 3. β-BHC | 14. Endrin |
| 4. γ-BHC | 15. Endosulfan II |
| 5. δ-BHC | 16. p,p'-DDD |
| 6. Heptachlor | 17. Endrin aldehyde |
| 7. Aldrin | 18. Endosulfan sulfate |
| 8. Heptachlor epoxide | 19. p,p'-DDT |
| 9. γ-Chlordane | 20. Endrin ketone |
| 10. Endosulfan I | 21. Methoxychlor |
| 11. α-Chlordane | 22. Decachlorobiphenyl (IS) |

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

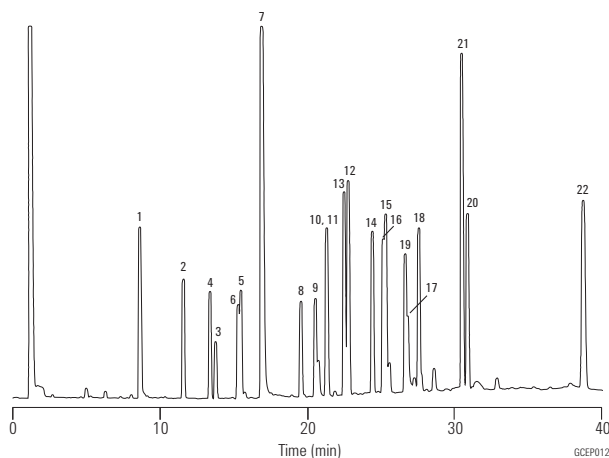

Organochlorine Pesticides, DB-5/DB-1701P

Column: DB-5
123-5032
30 m x 0.32 mm, 0.25 µm

Carrier: Helium at 29.2 cm/s, measured at 150 °C

Oven: 60 °C for 0.5 min
60-140 °C at 20 °C/min
140-280 °C at 11 °C/min
280 °C for 23 min



Column: DB-1701P
123-7732
30 m x 0.32 mm, 0.25 µm

Injection: Splitless, 200 °C

Column: Guard Column
160-2535-10
30 m x 0.32 mm, 0.25 µm

Detector: ECD, 325 °C
Nitrogen makeup gas at 30 mL/min

Sample: 2.0 µL, 20-200 pg/µL

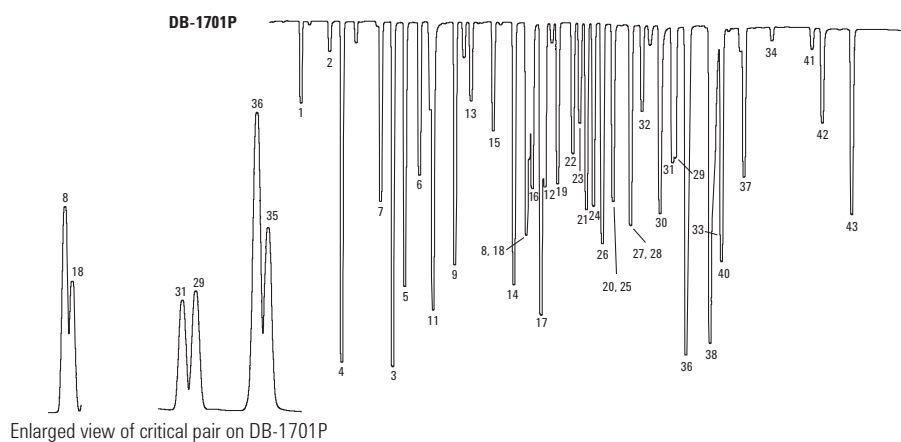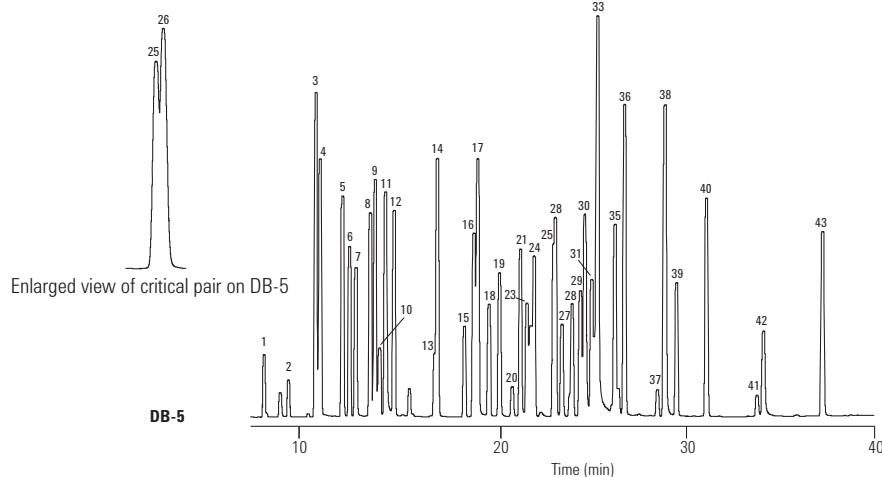


1. Etridiazole
2. Chloroneb
3. Propachlor
4. Tetrachloro-m-xylene (IS)
5. Trifluralin
6. α-BHC
7. Hexachlorobenzene
8. β-BHC
9. γ-BHC
10. Pentachloronitrobenzene
11. p,p'-Dichlorobiphenyl
12. δ-BHC
13. Heptachlor
14. Alachlor
15. Aldrin
16. Chlorpyrifos
17. DCPA
18. Isodrin
19. Heptachlor epoxide
20. Captan
21. γ-Chlordane
22. o,p'-DDE
23. Endosulfan I
24. α-Chlordane
25. Dieldrin
26. p,p'-DDE
27. o,p'-DDD
28. Endrin
29. Endosulfan II
30. Chlorobenzilate
31. p,p'-DDD
32. o,p'-DDT
33. Endrin aldehyde
34. Endrin ketone
35. Carbophenothion
36. p,p'-DDT
37. Endosulfan sulfate
38. Hexabromobenzene (HBB)
39. Methoxychlor
40. Mirex
41. cis-Permethrin
42. trans-Permethrin
43. Decachlorobiphenyl (IS)

Organochlorine Pesticides

Column: VF-17ms
CP8982
30 m x 0.25 mm, 0.25 µm

Sample: 1.0 µL
Sample Conc: 200 µg/mL
Carrier: Helium, 70 kPa
Injection: Splitter, 1:100
Detector: MS, Ion Trap, TIC

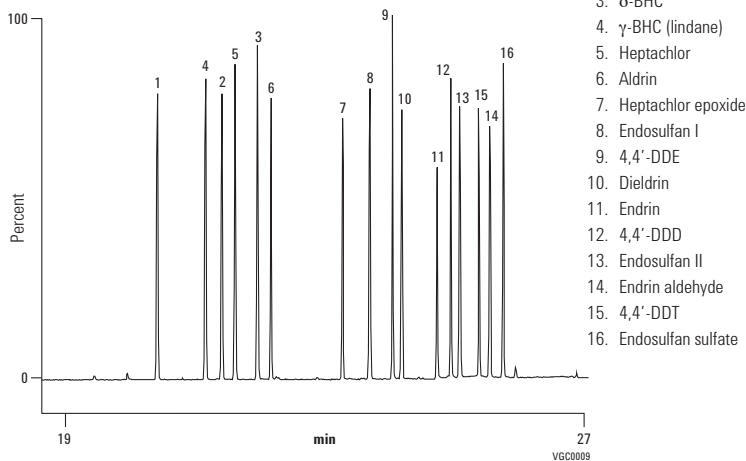**Nitrogen/Phosphorus Containing Pesticides,
EPA Method 507**

Column: HP-5ms
19091S-433
30 m x 0.25 mm, 0.25 µm


Carrier: Helium, 30 cm/s (13.6 psi) pressure program
Oven: 80-178 °C at 30 °C/min
178 °C for 4 min
178-205 °C at 2 °C/min
205-310 °C at 30 °C/min
310 °C for 4 min

Injection: Splitless, 260 °C
1 min purge delay

Detector: NPD, 290 °C
Helium makeup gas at 30 mL/min

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

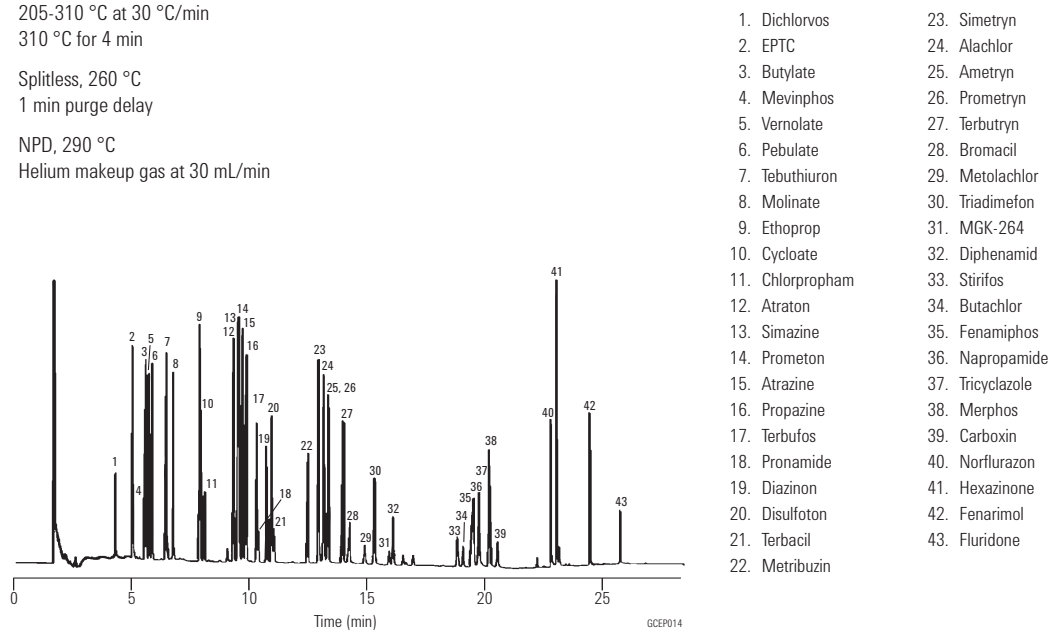

Herbicides I

Column: DB-XLB
122-1232
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 32 cm/s, measured at 50 °C

Oven: 50 °C for 1 min
50-180 °C at 10 °C/min
180-230 °C at 5 °C/min
230-320 °C at 10 °C/min
320 °C for 2 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: MSD, 300 °C transfer line
Full scan 50-400


Sample: 2 µL x 10-50 ng/µL solution
in acetone

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

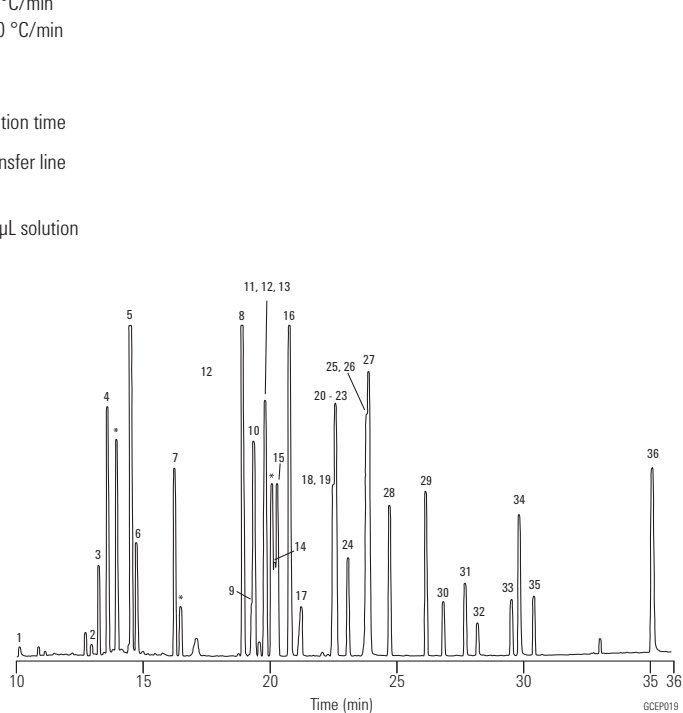

- | | |
|-------------------|------------------|
| 1. Monuron | 19. Propanil |
| 2. Diuron | 20. Ametryn |
| 3. EPTC | 21. Prometryn |
| 4. Dichlobenil | 22. Simetryn |
| 5. Vernolate | 23. Metribuzin |
| 6. Pebulate | 24. Terbutryn |
| 7. Molinate | 25. Metolachlor |
| 8. Sulfallate | 26. Bromacil |
| 9. Atraton | 27. Dacthal |
| 10. Prometon | 28. Diphenamid |
| 11. Atrazine | 29. Butachlor |
| 12. Propazine | 30. Napropamide |
| 13. Simazine | 31. Carboxin |
| 14. Terbutylazine | 32. Tricyclazole |
| 15. Pronamide | 33. Norflurazon |
| 16. Secbumeton | 34. Hexazinone |
| 17. Terbacil | 35. Difolatan |
| 18. Alachlor | 36. Fluridone |

* Impurity

Herbicides II


Column: DB-210
122-0232
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 35 cm/s

Oven: 140-215 °C at 3 °C/min

Injection: Split 1:50, 1 µL

Detector: ECD, 300 °C
Nitrogen makeup gas at 30 mL/min

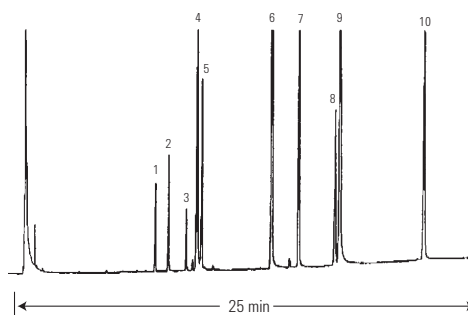

- | |
|-----------------|
| 1. Phorate |
| 2. Ethoprop |
| 3. Terbufos |
| 4. Atrazine |
| 5. Fonofos |
| 6. Propachlor |
| 7. Chlorpyrifos |
| 8. Alachlor |
| 9. Metolachlor |
| 10. Cyanazine |

GCHERB01

C₁ and C₂ Halocarbons (Freons)

Column: GS-GasPro
113-4362
60 m x 0.32 mm


Carrier: Helium at 35 cm/s, constant velocity

Oven: 40 °C for 2 min,
40-120 °C at 10 °C/min
120 °C for 3 min
120-200 °C at 10 °C/min

Injection: Splitless, 250 °C
0.20 min purge activation time

Detector: MSD, 280 °C,
Full scan 45-180 amu

Sample: 1.0 µL of 100 ppm mixture
of AccuStandard M-REF &
M-REF-X in methanol

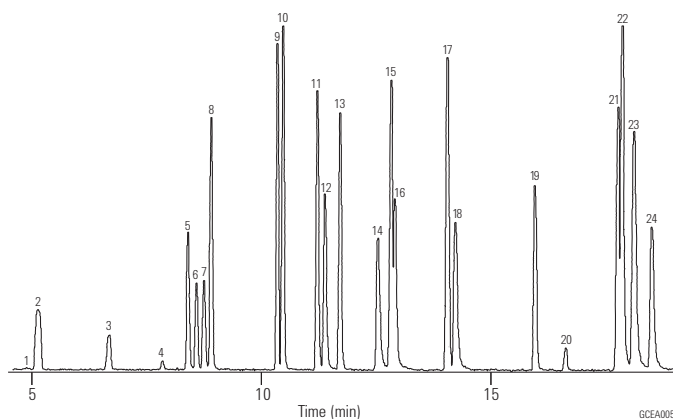**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316
Seal: Gold plated seal, 18740-20885
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Freon #

1. Chlorotrifluoromethane*	13
2. Trifluoromethane	23
3. Bromotrifluoromethane	13B1
4. Chloropentafluoroethane	115
5. Pentafluoroethane	125
6. 1,1,1-Trifluoroethane	143a
7. Dichlorodifluoromethane	12
8. Chlorodifluoromethane	22
9. 1,1,1,2-Tetrafluoroethane	134a
10. Chloromethane	40
11. 1,1,2,2-Tetrafluoroethane	134
12. Bromochlorodifluoromethane	12B1
13. 1,1-Difluoroethane	152a
14. 1,2-Dichloro-1,1,2,2-tetrafluoroethane	114
15. 2-Chloro-1,1,2-tetrafluoroethane	124
16. 1-Chloro-1,1-difluoroethane	142b
17. Dichlorofluoromethane	21
18. Trichlorofluoromethane	11
19. Chloroethane	160
20. Dichloromethane	30
21. 1,1-Dichloro-1-fluoroethane	141b
22. 2,2-Dichloro-1,1,1-trifluoroethane	123
23. 1,1,2-Trichloro-1,2,2-trifluoroethane	113
24. 1,2-Dibromo-1,1,2,2-tetrafluoroethane	114B2

*Peak not shown

Nitrogen Containing Herbicides (EPA Method 507)


Column: DB-35
125-1937
30 m x 0.53 mm, 0.50 µm

Carrier: Helium at 38 cm/s (5 mL/min),
measured at 150 °C

Oven: 60 °C for 1 min
60-290 °C at 15 °C/min
290 °C for 5 min

Injection: Megabore direct, 290 °C, 1 µL of 3 ng/µL standard

Detector: NPD, 290 °C

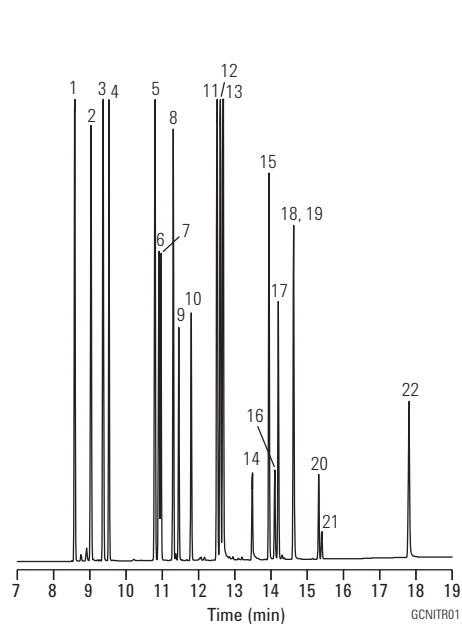

1. Eptam
2. Sutan
3. Vernam
4. Tillam
5. Ordram
6. Treflan
7. Balan
8. Ro-Neet
9. Propachlor
10. Tolban
11. Propazine
12. Atrazine
13. Simazine
14. Terbacil
15. Sencor
16. Dual
17. Paarlan
18. Prowl
19. Bromacil
20. Oxadiazon
21. GOAL
22. Hexazinone


Free Phenols

Column: HP-50+
19091L-433
30 m x 0.25 mm, 0.25 μ m

Carrier: Hydrogen, constant flow 45 cm/s

Injection: Split, 100:1

Detector: FID, 300 °C

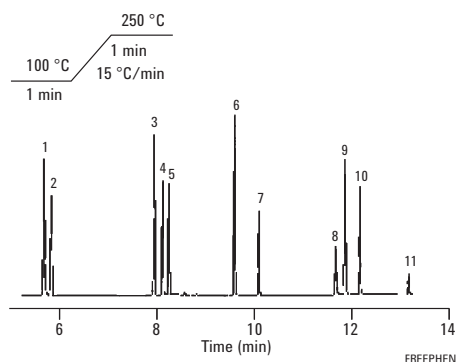


1. Phenol
2. 2-Chlorophenol
3. 2,4-Dimethylphenol
4. 2-Nitrophenol
5. 2,4-Dichlorophenol
6. 4-Chloro-3-methylphenol
7. 2,4,6-Trichlorophenol
8. 2,4-Dinitrophenol
9. 4-Nitrophenol
10. 2-Methyl-4,6-dinitrophenol
11. Pentachlorophenol

EPA 625 Halogenated Pesticides on "1701" Type Phases

Column: VF-1701 Pesticides
CP9070
30 m x 0.25 mm, 0.25 μ m

Oven: 150 °C, 5 °C/min to 275 °C

Injection: Split: T=275 °C
ECD: T=275 °C, 2 pg

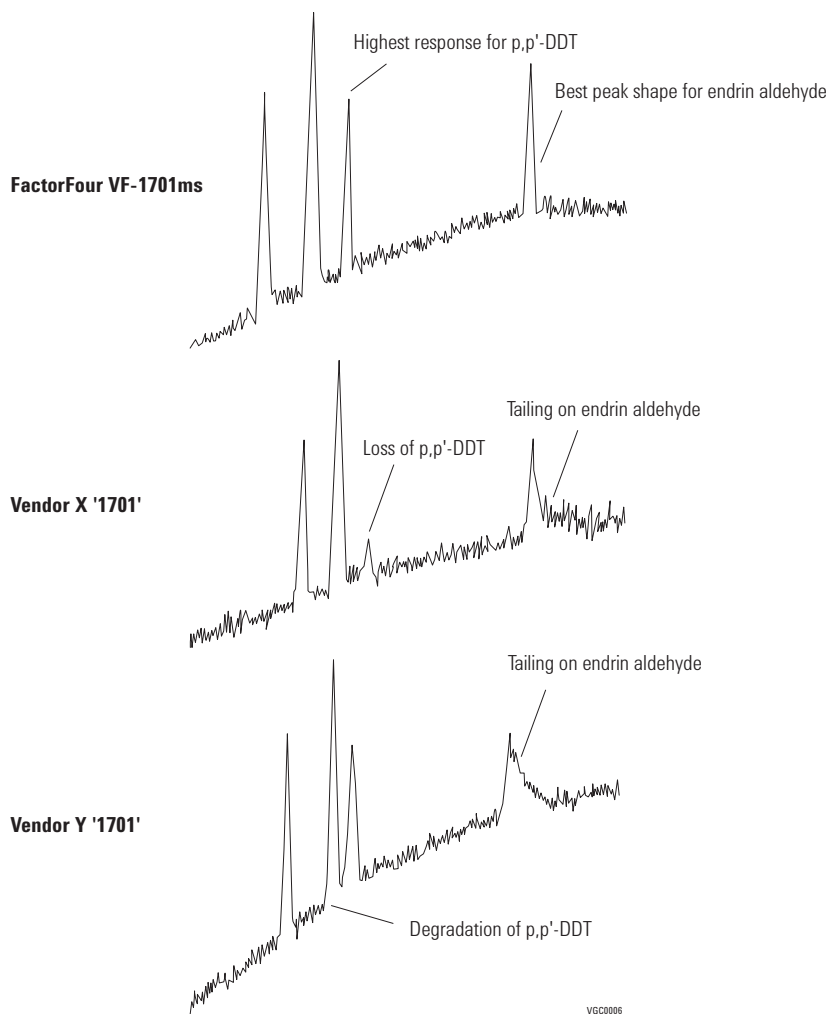

Organochlorine Pesticides to EPA 625 via GC/MS


Column: VF-35ms
CP8877
30 m x 0.25 mm, 0.25 µm

Carrier: Helium, approx. 1.0 mL/min, 60 kPa

Oven: 45 °C + 10 °C/min to 325 °C

Injection: Split/splitless, in split mode, 1:100

Detector: Ion Trap MS

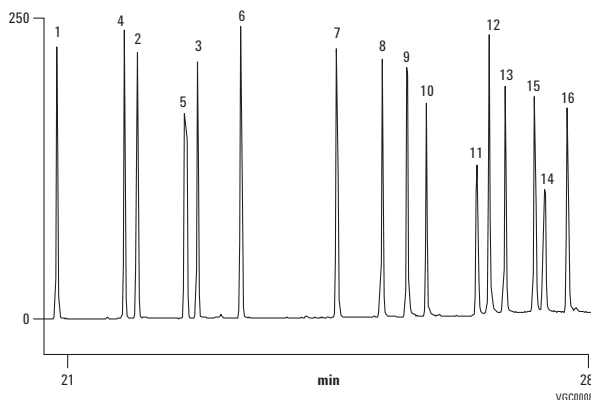

1. α-BHC
2. β-BHC
3. δ-BHC
4. γ-BHC (lindane)
5. Heptachlor
6. Aldrin
7. Heptachlor epoxide
8. Endosulfan I
9. 4,4'-DDE
10. Dieldrin
11. Endrin
12. 4,4'-DDD
13. Endosulfan II
14. Endrin aldehyde
15. 4,4'-DDT
16. Endosulfan sulfate

Organochlorine Pesticides I EPA Method 8081A

Column: DB-35ms
122-3832
30 m x 0.25 mm, 0.25 µm


Carrier: Helium at 35 cm/s, measured at 50 °C

Oven: 50 °C for 1 min
50-100 °C at 25 °C/min
100-300 °C at 5 °C/min
300 °C for 5 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: MSD, 300 °C transfer line
Full scan at m/z 50-500

Sample: 1 µL of 35 µg/mL composite 8081A
standards, AccuStandard Inc.

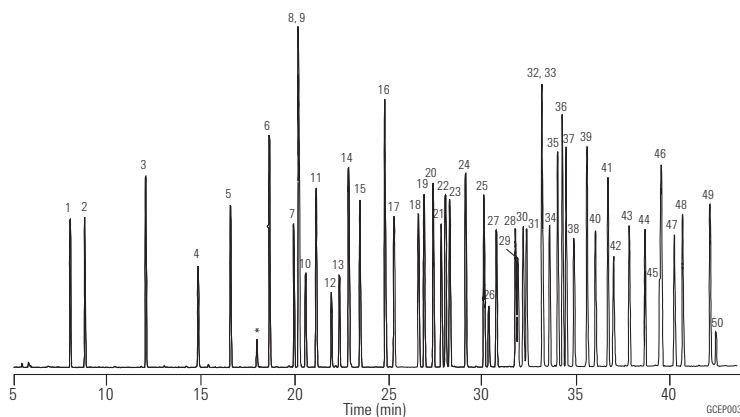

1. 1,2-Dibromo-3-chloropropane
2. 4-Chloro-3-nitrobenzotrifluoride (SS)
3. Hexachloropentadiene
4. 1-Bromo-2-nitrobenzene (IS)
5. Terrazole
6. Chloroneb
7. Trifluralin
8. 2-Bromobiphenyl (SS)
9. Tetrachloro m-xylene (SS)
10. α, α-Dibromo-m-xylene
11. Propachlor
12. Di-allate A
13. Di-allate B
14. Hexachlorobenzene
15. α-BHC
16. Pentachloronitrobenzene (IS)
17. γ-BHC
18. β-BHC
19. Heptachlor
20. Alachlor
21. δ-BHC
22. Chlorothalonil
23. Aldrin
24. Dacthal
25. Isodrin
26. Kelthane
27. Heptachlor epoxide
28. γ-Chlordane
29. trans-Nonachlor
30. α-Chlordane
31. Endosulfan I
32. Captan
33. p,p'-DDE
34. Dieldrin
35. Chlorobenzilate
36. Perthane
37. Chloropropylate
38. Endrin
39. p,p'-DDD
40. Endosulfan II
41. p,p'-DDT
42. Endrin aldehyde
43. Endosulfan sulfate
44. Dibutyl chlorendate (SS)
45. Captafol
46. Methoxychlor
47. Endrin ketone
48. Mirex
49. cis-Permethrin
50. trans-Permethrin

* Breakdown Products
SS - Surrogate Standard
IS - Internal Standard

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Standards used were a composite of individual solutions supplied courtesy of AccuStandard Inc., 25 Science Park, New Haven, CT 06511, 800-442-5290.

Organochlorine Pesticides II EPA Method 8081A

Column: DB-5ms
122-5532
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 35 cm/s, measured at 50 °C

Oven: 50 °C for 1 min
50-100 °C at 25 °C/min
100-300 °C at 5 °C/min
300 °C for 5 min

Injection: Splitless, 250 °C
30 s purge activation time


Detector: MSD, 300 °C transfer line
Full scan at m/z 50-500

Sample: 1 µL of 35 µg/mL composite 8081A standards, AccuStandard Inc.

Standards used were a composite of individual solutions supplied courtesy of AccuStandard Inc., 25 Science Park, New Haven, CT 06511, 800-442-5290.

- | | |
|--|------------------------------|
| 1. 1,2-Dibromo-3-chloropropane | 26. Kelthane |
| 2. 4-Chloro-3-nitrobenzotrifluoride (SS) | 27. Heptachlor epoxide |
| 3. Hexachloropentadiene | 28. γ-Chlordane |
| 4. 1-Bromo-2-nitrobenzene (IS) | 29. trans-Nonachlor |
| 5. Terrazole | 30. α-Chlordane |
| 6. Chloroneb | 31. Endosulfan I |
| 7. Trifluralin | 32. Captan |
| 8. 2-Bromobiphenyl (SS) | 33. p,p'-DDE |
| 9. Tetrachloro m-xylene (SS) | 34. Dieldrin |
| 10. α, α-Dibromo-m-xylene | 35. Chlorobenzilate |
| 11. Propachlor | 36. Perthane |
| 12. Di-allate A | 37. Chloropropylate |
| 13. Di-allate B | 38. Endrin |
| 14. Hexachlorobenzene | 39. p,p'-DDD |
| 15. α-BHC | 40. Endosulfan II |
| 16. Pentachloronitrobenzene (IS) | 41. p,p'-DDT |
| 17. γ-BHC | 42. Endrin aldehyde |
| 18. β-BHC | 43. Endosulfan sulfate |
| 19. Heptachlor | 44. Dibutyl chlorendate (SS) |
| 20. Alachlor | 45. Captafol |
| 21. δ-BHC | 46. Methoxychlor |
| 22. Chlorothalonil | 47. Endrin ketone |
| 23. Aldrin | 48. Mirex |
| 24. Dacthal | 49. cis-Permethrin |
| 25. Isodrin | 50. trans-Permethrin |

* Breakdown Products
SS - Surrogate Standard
IS - Internal Standard

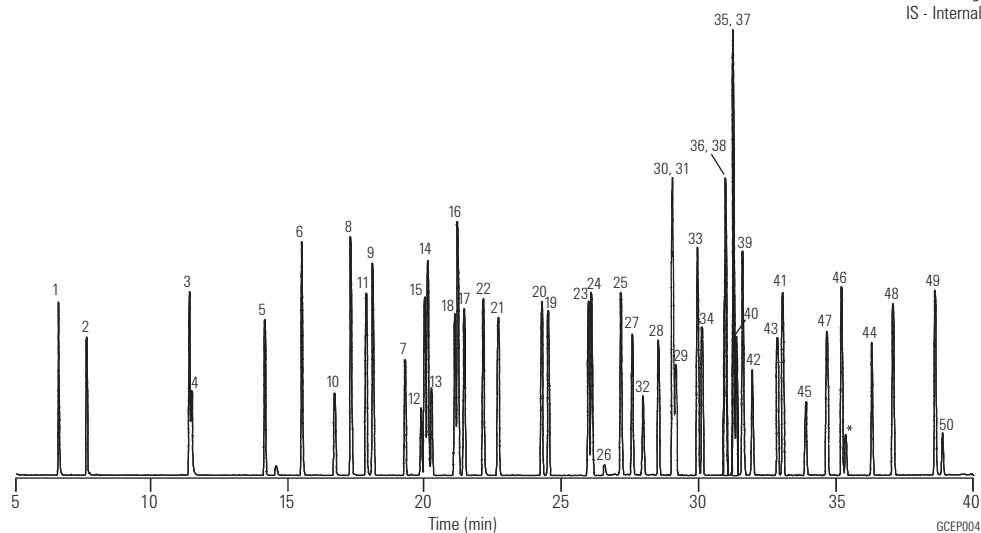

Organophosphorus Pesticides in Apple Matrix

Column: DB-35ms Ultra Inert
121-3822UI
20 m x 0.18 mm, 0.18 µm

Instrument: Agilent 7890 GC/Agilent 5975C Series GC/MSD

Sampler: Agilent 7683B automatic liquid sampler,
5.0 µL syringe (p/n 5181-1273)

CFT Device: Purged 2-way splitter (p/n G3180B)
Split Ratio MSD:FPD = 3:1

MSD Restrictor: 1.2 m x 0.15 mm id deactivated fused silica tubing

FPD Restrictor: 1.4 m x 0.15 mm id deactivated fused silica tubing

PCM 1: 3.8 psi constant pressure

Inlet: 1 µL splitless; 250 °C, purge flow 60 mL/min
at 0.25 min, gas saver on at 2 min 20 mL/min


Carrier: Helium, constant pressure 43.5 psi at 95 °C

Oven: 95 °C (1.3 min), 15 °C/min to 125 °C, 5 °C/min to
165 °C, 2.5 °C/min to 195 °C, 20 °C/min to 280 °C
(3.75 min)


Postrun: 5 min at 280 °C, PCM 1 pressure 70 psi during
backflush, 2 psi inlet pressure during backflush

Detector: 310 °C transfer line, 310 °C source, 150 °C quad

1. Oxydemeton-methyl
2. Methamidophos
3. Mevinphos
4. Acephate
5. Naled
6. Diazinon
7. Dimethoate
8. Chlorpyrifos
9. Malathion
10. Methidathion
11. TPP (surrogate std)
12. Phosmet
- 13.

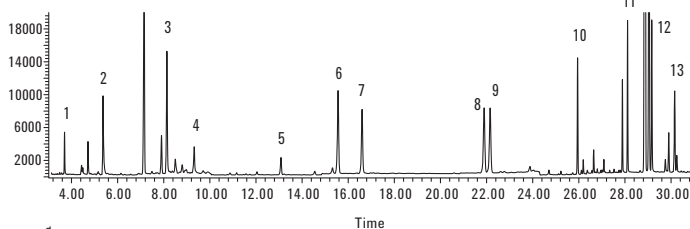

MSD (SIM): 600 ng/mL

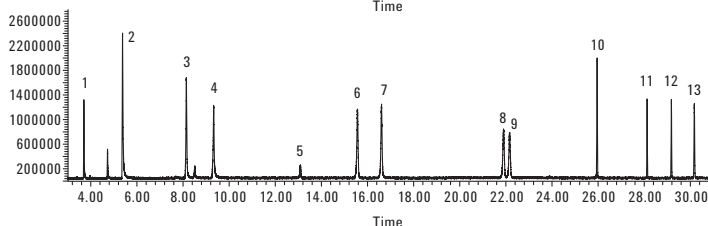

FPD (P): 200 ng/mL

GC/MS-SIM and FPD chromatograms of a matrix matched organophosphorus pesticides standard analyzed on an Agilent J&W DB-35ms UI column. The effluent split ratio is MSD:FPD = 3:1.

Environmental Applications, Semivolatiles

Agilent's Ultra Inert Test Probe Mixture

Column: DB-5ms Ultra Inert
122-5532UI
30 m x 0.25 mm, 0.25 μ m

Carrier: Hydrogen, constant pressure, 38 cm/s


Oven: 65 °C isothermal

Sampler: Agilent 7683B, 0.5 μ L syringe
(p/n 5188-5246), 0.02 μ L split injection

Injection: Split/splitless, 250 °C, 1.4 mL/min; split column flow
900 mL/min; gas saver flow 75 mL/min at 2.0 min

Detector: FID at 325 °C; 450 mL/min air, 40 mL/min hydrogen,
45 mL/min nitrogen makeup

A properly deactivated DB-5ms Ultra Inert column delivers symmetrical peak shapes, along with increased peak heights, which allow for accurate integration and detection of trace analytes.

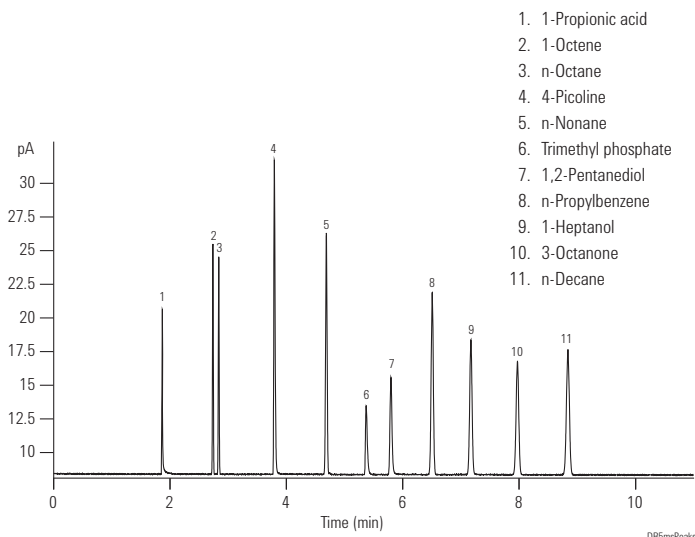

Trace Level Polycyclic Aromatic Hydrocarbon (PAH) Analyses

Column: DB-5ms Ultra Inert
122-5532UI
30 m x 0.25 mm, 0.25 μ m


Carrier: Helium constant flow 30 cm/s

Oven: 40 °C (1 min) to 100 °C (15 °C/min)
10 °C to 210 °C (1 min)
5 °C/min to 310 °C (8 min)

Injection: Split/splitless, 260 °C, 53.7 mL/min total flow,
purge flow 50 mL/min on at 0.5 min,
gas saver flow 80 mL/min on at 3.0 min

Detector: MSD source at 300 °C
Quadrupole at 180 °C
Transfer line at 290 °C
Scan range 50-550 amu

- | | |
|-------------------|----------------------------|
| 1. Naphthalene | 9. Benz[a]anthracene |
| 2. Acenaphthylene | 10. Chrysene |
| 3. Acenaphthene | 11. Benzo[b]fluoranthene |
| 4. Fluorene | 12. Benzo[k]fluoranthene |
| 5. Phenanthrene | 13. Benzo[a]pyrene |
| 6. Anthracene | 14. Indeno[1,2,3-cd]pyrene |
| 7. Fluoranthene | 15. Dibenzo[a,h]anthracene |
| 8. Pyrene | 16. Benzo[g,h,i]perylene |

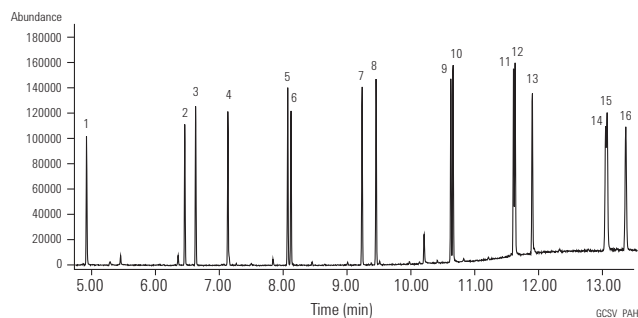

Tetrachlorodibenzo-p-furans


Column A: DB-225
122-2232
30 m x 0.25 mm, 0.25 µm

Column B: DB-225ms
122-2932
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 12 mL/min

Oven: 160-250 °C at 7 °C/min
250 °C until compounds elute

Injection: Splitless, 240 °C

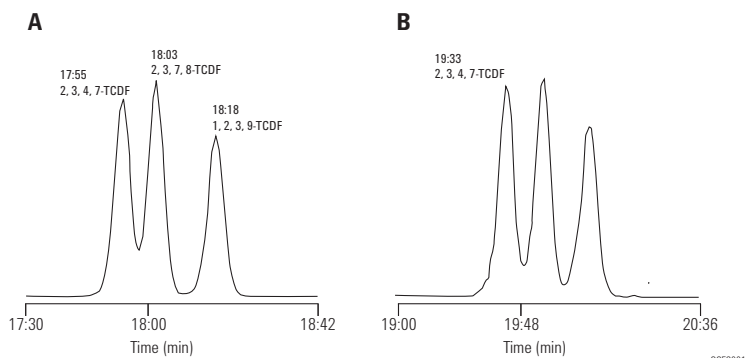

GCE5001

Note the separation between 2,3,7,8-TCDF and 2,3,4,7-TCDF on DB-225 is also easily achievable (and actually a little better) on Agilent J&W DB-225ms.

Congeners in DIN Method PCBs

Column: DB-XLB
122-1236
30 m x 0.25 mm, 0.50 µm

Carrier: Helium at 34.2 cm/s, measured at 150 °C

Oven: 100 °C for 1 min
100-320 °C at 5.6 °C/min

Injection: Hot on-column, 250 °C
Split flow 100 mL/min

Detector: MSD, 300 °C transfer line
SIM of 221.9, 255.9,
291.9, 325.8, 359.8,
395.8, 429.7, 463.7


Sample: 2 µL dilute Aroclor mixture

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

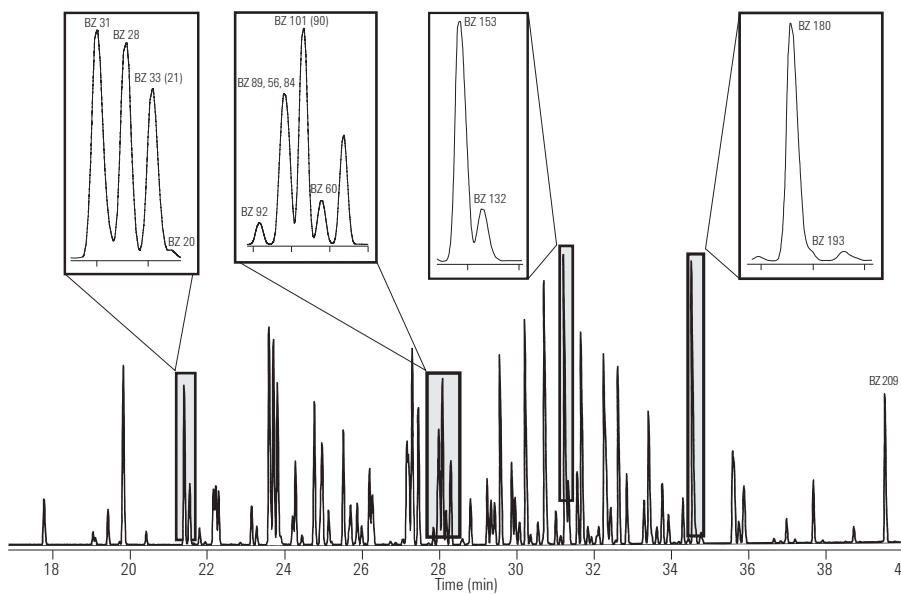

GCE5005

Extended Temperature Program Resolving Congeners 52 and 138

Column: DB-XLB
122-1236
30 m x 0.25 mm, 0.50 µm


Carrier: Helium at 34.2 cm/s, measured at 150 °C

Oven: 100 °C for 1 min
100-275 °C at 1.6 °C/min

Injection: Hot on-column, 250 °C
Split flow 100 mL/min

Detector: MSD, 300 °C transfer line
SIM of 221.9, 255.9, 291.9, 325.8,
359.8, 395.8, 429.7, 463.7

Sample: 2 µL dilute Aroclor mixture

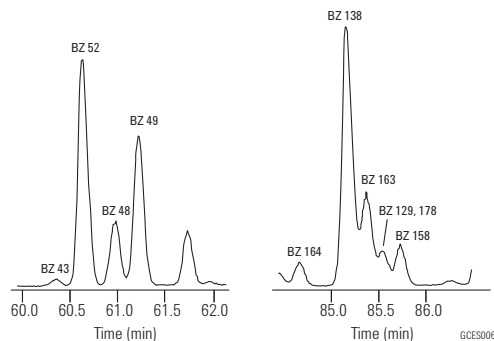

PCBs by EPA Method 8082

Column: DB-35ms
123-3832
30 m x 0.32 mm, 0.25 µm

Column: DB-XLB
123-1236
30 m x 0.32 mm, 0.50 µm



Carrier: Helium at 45 cm/s
(EPC in constant flow mode)

Oven: 110 °C for 0.5 min
110-320 °C at 15 °C/min
320 °C for 5 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: µECD, 350 °C
Nitrogen makeup gas
(column + makeup flow =
30 mL/min constant flow)

Sample: 50 pg per component

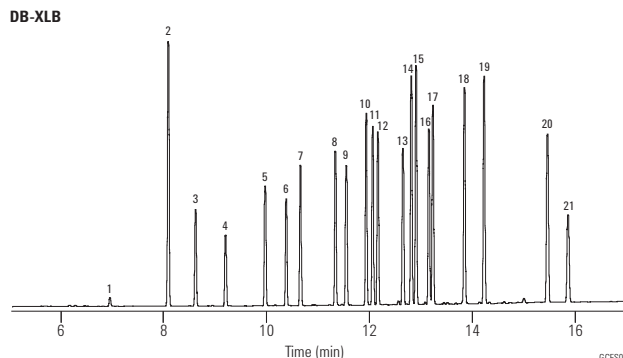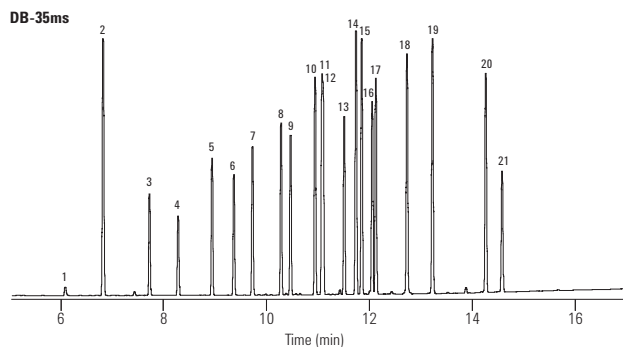

1. IUPAC 1
 2. Tetrachloro-m-xylene (IS/SS)
 3. IUPAC 5
 4. IUPAC 18
 5. IUPAC 31
 6. IUPAC 52
 7. IUPAC 44
 8. IUPAC 66
 9. IUPAC 101
 10. IUPAC 87
 11. IUPAC 110
 12. IUPAC 151
 13. IUPAC 153
 14. IUPAC 141
 15. IUPAC 137
 16. IUPAC 187
 17. IUPAC 183
 18. IUPAC 180
 19. IUPAC 170
 20. IUPAC 206
 21. Decachlorobiphenyl (IS/SS)
- IS/SS - Internal Standard/
Surrogate Standard

Suggested Supplies

Septum: 11 mm Advanced Green septa,
5183-4759

Liner: Splitless, single taper, deactivated,
4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP,
5181-1267


Automated Cleanup of PCB extracts from Waste Oil Using 7696A Sample Prep Workbench

Column: DB-5ms
122-5532
30 m x 0.25 mm, 0.25 µm

Instrument: Agilent 7000 Triple Quadrupole GC/MS system
Carrier: Helium, 1 mL/min constant flow
During backflush: 2 mL/min
Oven: 80 °C (1 min), 10 °C/min to 305 °C, 7.5 min hold
Injection: 1 µL, pulsed splitless
QuickSwap: 28 kPa constant pressure
Backflush: Start at 23.5 min

Detector: MRM mode
CE 25 V, dwell time 100 ms per transition
Trichloro-biphenyls: 256.0 > 186.0; 258.0 > 186.0
Tetrachloro-biphenyls: 293.8 > 222.0; 291.8 > 222.0
Pentachloro-biphenyls: 325.8 > 256.0; 327.8 > 256.0
Hexachloro-biphenyls: 359.9 > 289.9; 361.9 > 289.9
Heptachloro-biphenyls: 393.8 > 323.8; 395.8 > 323.8
Octachloronaphthalene (IS): 404.0 > 404.0 (CE OV)

Sample: Reference sample BCR-449, five aliquots

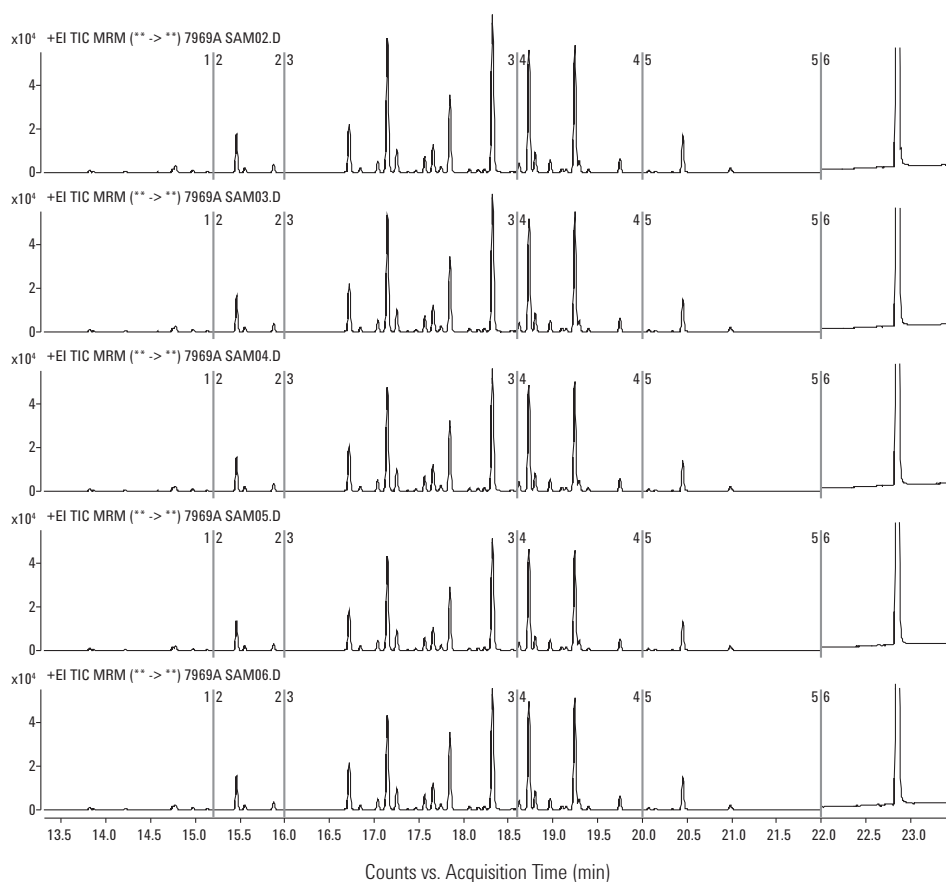

Pyrethrins

Column: DB-1
123-1032
30 m x 0.32 mm, 0.25 µm

Carrier: Helium at 39 cm/s, measured at 150 °C


Oven: 180 °C for 11 min
180-200 °C at 10 °C/min
200 °C for 8 min
200-210 °C at 10 °C/min
210 °C for 18 min
210-245 °C at 30 °C/min
245 °C for 4 min

Injection: Split, 250 °C
Split ratio 1:20

Detector: FID, 300 °C
Helium makeup gas at 30 mL/min

Sample: 1 µL


*Chromatogram courtesy of Khan Nguyen and Richard Moorman
of Sandoz Agro Inc.*

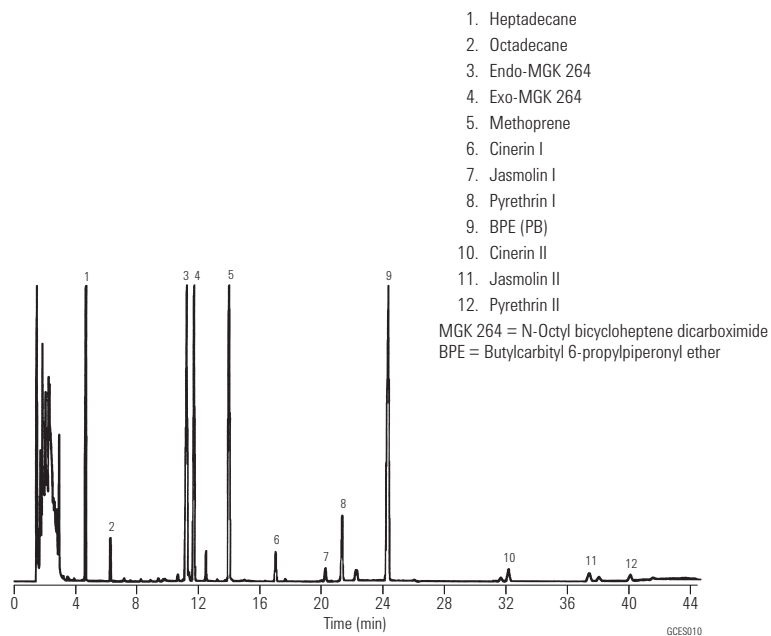

Organotin Compounds I

Column: HP-1
19091Z-012
25 m x 0.32 mm, 0.17 µm

Carrier: Helium, 100 kPa
Oven: 50 °C for 1 min
50-260 °C at 15 °C/min
Injection: Splitless
Detector: AED, 330 °C
Sample: 1 µL

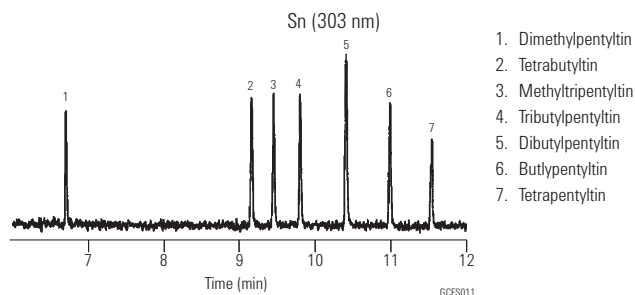
**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Organotin Compounds II

Column: HP-5
19091J-002
25 m x 0.20 mm, 0.11 µm

Carrier: Helium, 0.75 mL/min constant flow
Oven: 60-360 °C at 5 °C/min
Injection: Splitless, 300 °C
Detector: AED, 300 °C
Hg selective at 254 nm
Sample: 1 µL

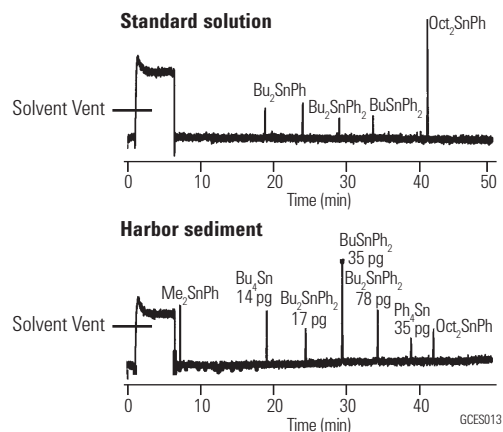**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Semivolatile Compounds, US EPA Method 8270

Column: HP-5ms
19091S-133
30 m x 0.25 mm, 0.50 µm

Carrier: Ramped flow 1.2 mL/min for 0.0 min
Ramp at 99 mL/min to 2.0 mL/min
2.0 mL/min for 0.35 min
Ramp at 10 mL/min to 1.2 mL/min

Oven: 40 °C for 1.0 min
40-100 °C at 15 °C/min
100-240 °C at 20 °C/min
240-310 °C at 10 °C/min

Injection: Splitless, 250 °C
30 mL/min purge flow
at 0.35 min

Detector: 5973 MSD, 310 °C transfer line
Scan range 35-500 amu,
3.25 scans/s

Sample: 1 µL of 50 ng standard


Suggested Supplies

Septum: 11 mm Advanced Green septa,
5183-4759

Liner: Splitless, single taper, deactivated,
4 mm id, 5181-3316

Syringe: 10 µL tapered,
FN 23-26s/42/HP, 5181-1267

- | | | | |
|----------------------------------|---------------------------------|--------------------------------|---------------------------------|
| 1. n-Nitrosodimethylamine | 36. 2,4,5-Trichlorophenol | 52. Fluorene | 68. Terphenyl-d14 |
| 2. Pyridine | 37. 2-Fluorobiphenyl | 53. 4-Nitroaniline | 69. Benzidine |
| 3. 2-Fluorophenol | 38. 2-Chloronaphthalene | 54. 4,6-Dinitro-2-methylphenol | 70. Butylbenzylphthalate |
| 4. Phenol-d5 | 39. 2-Nitroaniline | 55. n-Nitrosodiphenylamine | 71. 3,3'-Dichlorobenzidine |
| 5. Phenol | 40. Dimethyl phthalate | 56. Azobenzene | 72. Benzo[a]anthracene |
| 6. Aniline | 41. 2,6-Dinitrotoluene | 57. 2,4,6-Tribromophenol | 73. Chrysene-d12 |
| 7. Bis(2-chloroethyl) ether | 42. Acenaphthylene | 58. 4-Bromophenyl-phenylether | 74. Chrysene |
| 8. 2-Chlorophenol | 43. 3-Nitroaniline | 59. Hexachlorobenzene | 75. Bis(2-ethylhexyl) phthalate |
| 9. 1,3-Dichlorobenzene | 44. Acenaphthene-d10 | 60. Pentachlorophenol | 76. Di-n-octylphthalate |
| 10. 1,4-Dichlorobenzene-d4 | 45. Acenaphthene | 61. Phenanthrene-d10 | 77. Benzo[b]fluoranthene |
| 11. 1,4-Dichlorobenzene | 46. 2,4-Dinitrophenol | 62. Phenanthrene | 78. Benzo[k]fluoranthene |
| 12. Benzyl alcohol | 47. 4-Nitrophenol | 63. Anthracene | 79. Benzo[a]pyrene |
| 13. 1,2-Dichlorobenzene | 48. Dibenzofuran | 64. Carbazole | 80. Perylene-d12 |
| 14. 2-Methylphenol | 49. 2,4-Dinitrotoluene | 65. Di-n-butyl phthalate | 81. Indeno[1,2,3-cd]pyrene |
| 15. Bis(2-chloroisopropyl) ether | 50. Diethyl phthalate | 66. Fluoranthene | 82. Dibenz[a,h]anthracene |
| 16. 4-Methylphenol | 51. 4-Chlorophenyl-phenyl ether | 67. Pyrene | 83. Benzo[g,h,i]perylene |
| 17. n-Nitroso-di-n-propylamine | | | |
| 18. Hexachloroethane | | | |
| 19. Nitrobenzene-d5 | | | |
| 20. Nitrobenzene | | | |
| 21. Isophorone | | | |
| 22. 2-Nitrophenol | | | |
| 23. 2,4-Dimethylphenol | | | |
| 24. Benzoic acid | | | |
| 25. Bis(2-chloroethoxy) methane | | | |
| 26. 2,4-Dichlorophenol | | | |
| 27. 1,2,4-Trichlorobenzene | | | |
| 28. Naphthalene-d8 | | | |
| 29. Naphthalene | | | |
| 30. 4-Chloroaniline | | | |
| 31. Hexachlorobutadiene | | | |
| 32. 4-Chloro-3-methylphenol | | | |
| 33. 2-Methylnaphthalene | | | |
| 34. Hexachlorocyclopentadiene | | | |
| 35. 2,4,6-Trichlorophenol | | | |

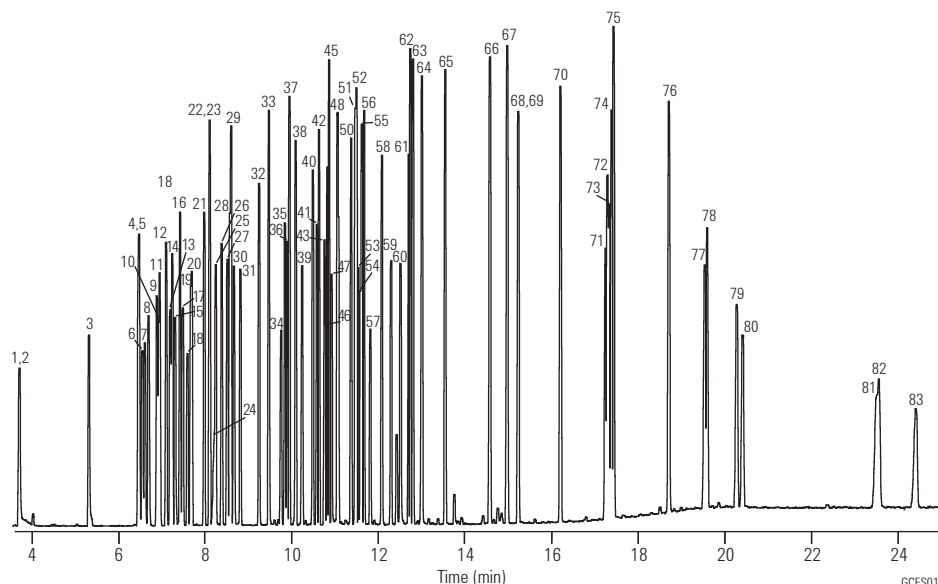

GCES015

A variety of HP-5ms and DB-5ms columns can be used for 8270 and similar semivolatiles applications. The column shown above was chosen to maximize inertness and robustness to residues with a thicker 0.5 µm film, but the price paid is a slightly longer run time.

An HP-5ms, 30 m x 0.25 mm id, 0.25 µm, p/n 19091S-433 would give shorter run times, with slightly less inertness and robustness.

A DB-5ms, 30 m x 0.25 mm id, 0.25 µm, p/n 122-5532, would give slightly less inertness, but offer better resolution of PAHs such as benzo[b]fluoranthene and benzo[k]fluoranthene.

A DB-5ms, 20 m x 0.18 mm x 0.18 µm, p/n 121-5522, can offer significantly reduced run times with a modest loss of inertness.

US EPA Method 8061 (Phthalate Esters)

Column: DB-5ms
121-5522
20 m x 0.18 mm, 0.18 µm

Carrier: Helium at 49 cm/s, measured at 80 °C
constant flow program

Oven: 80 °C for 0.5 min
80-160 °C at 30 °C/min
160-320 °C at 15 °C/min

Injection: Splitless, 300 °C
30 s purge activation time

Detector: MSD, 325 °C transfer line
Full scan m/z 50-400


Sample: 1 µL of 20 ng/µL
Method 8061 mixture (AccuStandard) in hexane

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

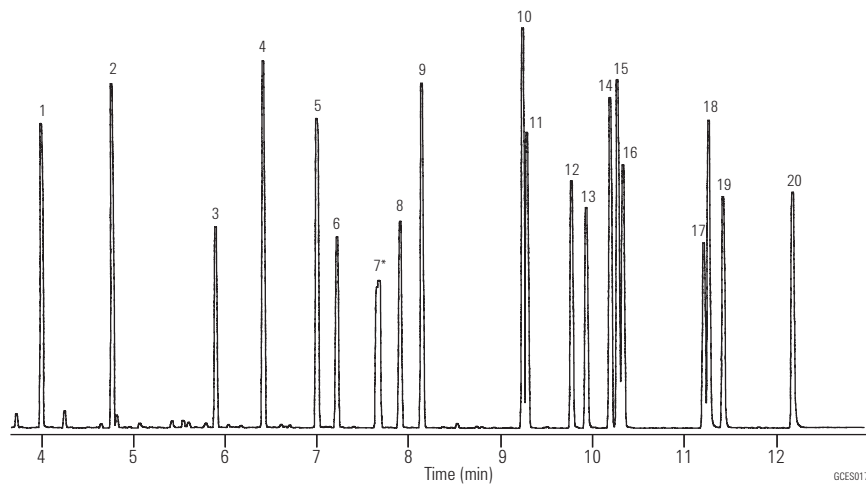

1. Dimethyl phthalate
2. Diethyl phthalate
3. Benzyl benzoate (IS)
4. Diisobutyl phthalate
5. Di-n-butyl phthalate
6. Bis(4-methoxyethyl) phthalate
7. Bis(4-methyl-2-pentyl) phthalate *
8. Bis(2-ethoxyethyl) phthalate
9. Diamyl phthalate
10. Dihexyl phthalate
11. Butyl benzyl phthalate
12. Hexyl 2-ethylhexyl phthalate
13. Bis(2-n-butoxyethyl) phthalate
14. Dicyclohexyl phthalate
15. Bis(2-ethylhexyl) phthalate
16. Diphenyl phthalate (SS)
17. Diphenyl isophthalate (SS)
18. Di-n-octyl phthalate
19. Dibenzyl phthalate (SS)
20. Dinonyl phthalate

* Two isomers

IS - Internal Standard

SS - Surrogate Standard

PAHs

Column: DB-17ms
122-4732
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at: 34.1 cm/s, measured at 150 °C

Oven: 95 °C for 0.5 min
95-340 °C at 5 °C/min
340 °C for 5 min

Injection: Split, 300 °C
Split ratio 1:40

Detector: MSD, 340 °C transfer line
Scan 80-330 amu

Sample: 2 µL, PAH standard

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Ions		Ions	
1. Naphthalene	128	27. 3,6-Dimethylphenanthrene	206, 191
2. 2-Methylnaphthalene	142, 141	28. 1,3-Dinitronaphthalene	126, 218
3. 1-Methylnaphthalene	142, 141	29. 1,5-Dinitronaphthalene	218, 114
4. Azulene	128	30. Fluoranthene	202
5. Acenaphthene	154	31. 2,2'-Dinitrobiphenyl	198, 139
6. Biphenyl	154	32. Pyrene	202
7. 2,6-Dimethylnaphthalene	156, 155	33. 2-Methylfluoranthene	216, 215
8. Acenaphthalene	152	34. 2,3-Benzofluorene	216, 215
9. Dibenzofuran	168, 139	35. Dodecahydrotriphenylene	240, 198
10. Dibenzo-p-dioxin	184	36. 1-Amino-4-nitronaphthalene	188, 115
11. Fluorene	166, 165	37. 9-Phenylanthracene	254, 253
12. 1-Nitronaphthalene	127, 173	38. 1,2-Benzanthracene	228
13. 9,10-Dihydroanthracene	179, 180	39. Chrysene	240
14. 2-Nitronaphthalene	127, 173	40. Benz[a]anthracene-7,12-dione	258, 202
15. 2-Nitrobiphenyl	152, 115	41. 2,7-Dinitrofluorene	256, 163
16. Dibenzothiophene	184	42. Benzo[b]fluoranthene	252
17. Phenanthrene	178	43. Benzo[k]fluoranthene	252
18. Anthracene	178	44. 7,12-Dimethylbenz[a]anthracene	256, 241
19. 3-Nitrobiphenyl	199, 152	45. Benzo[e]pyrene	252
20. 4-Nitrobiphenyl	199, 152	46. Benzo[a]pyrene	252
21. 5,6-Benzoquinoline	179	47. Perylene	252
22. Carbazole	167	48. 3-Methylcholanthrene	268
23. 2-Methylanthracene	192, 191	49. 9,10-Diphenylanthracene	330
24. 1,2,3,4-Tetrahydrofluoranthene	178, 206	50. 1,2,3,4-Dibenzanthracene	278
25. 2-Phenylnaphthalene	204	51. 1,2,5,6-Dibenzanthracene	278
26. 9-Methylanthracene	192, 191	52. Benzo[g,h,i]perylene	276

The DB-35ms is also an excellent column for PAH analysis.

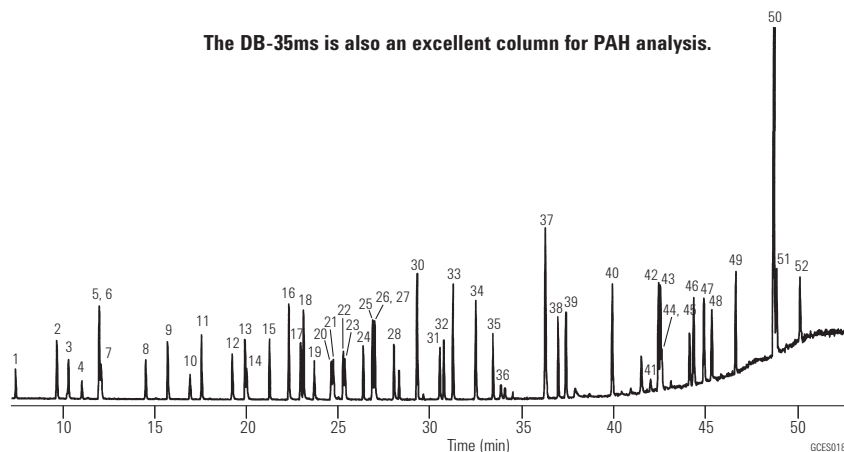

Phenols

Column: DB-5ms
122-5532
30 m x 0.25 mm, 0.25 µm

Column: DB-XLB
122-1232
30 m x 0.25 mm, 0.25 µm

Carrier: He at 1.2 mL/min constant flow



Oven: 40 °C for 2 min
40-100 °C at 40 °C/min
100 °C for 0.50 min
100-140 °C at 2 °C/min
140-340 °C at 30 °C/min

Injection: Pulsed splitless, 200 °C
Pulse pressure & time: 25 psi for 1 min
Purge flow & time: 50 mL/min for 0.25 min
Gas saver flow & time: 20 mL/min for 3 min

Detector: MSD, 320 °C transfer line
Quadrupole at 150 °C
Source at 230 °C

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

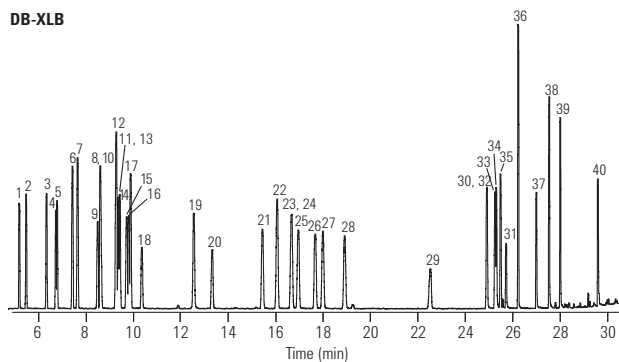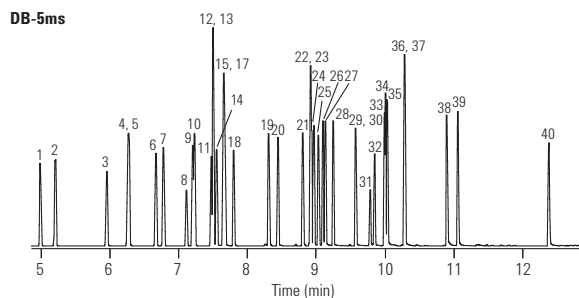

GCES019

1. Phenol
2. 2-Chlorophenol
3. 2-Methylphenol
4. 4-Methylphenol
5. 3-Methylphenol
6. 2-Chloro-5-methylphenol
7. 2,6-Dimethylphenol
8. 2-Nitrophenol
9. 2,4-Dimethylphenol
10. 2,5-Dimethylphenol
11. 2,4-Dichlorophenol
12. 2,3-Dimethylphenol
13. 2,5-Dichlorophenol
14. 2,3-Dichlorophenol
15. 2-Chlorophenol
16. 4-Chlorophenol
17. 3,4-Dimethylphenol
18. 2,6-Dichlorophenol
19. 4-Chloro-2-methylphenol
20. 4-Chloro-3-methylphenol
21. 2,3,5-Trichlorophenol
22. 2,4-Dibromophenol
23. 2,4,6-Trichlorophenol
24. 2,4,5-Trichlorophenol
25. 2,3,4-Trichlorophenol
26. 3,5-Dichlorophenol
27. 2,3,6-Trichlorophenol
28. 3,4-Dichlorophenol
29. 3-Nitrophenol
30. 2,5-Dinitrophenol
31. 2,4-Dinitrophenol
32. 4-Nitrophenol
33. 2,3,5,6-Tetrachlorophenol
34. 2,3,4,5-Tetrachlorophenol
35. 2,3,4,6-Tetrachlorophenol
36. 3,4,5-Trichlorophenol
37. 2-Methyl-4,6-dinitrophenol
38. Pentachlorophenol
39. Dinoseb
40. 2-Cyclohexyl-4,6-dinitrophenol

10 ng/ μ L Semivolatile Checkout Standard on a 20 m x 0.18 mm, 0.36 μ m Agilent J&W DB-UI 8270D Capillary GC Column using an Ultra Inert Liner with Wool

Column: DB-UI 8270D
121-9723
20 m x 0.18 mm, 0.36 μ m

Inlet: S/SL 1 μ L pulsed splitless, 300 °C 44 psi pulse to 1.4 min, purge flow 50 mL/min at 1.42 min, gas saver off

Inlet liner: Agilent Ultra Inert single taper with wool (p/n 5190-2293)

Oven: 40 °C (2.5 min), 25 °C/min to 320 °C (4.8 min)

Carrier: Helium, constant flow 1.58 mL/min set at 40 °C

MSD: 325 °C transfer line, 300 °C source, 150 °C quad, 30-550 amu range



GC/MSD: Agilent 7890 Series GC/5975C Series GC/MSD

Aux EPC: 2 psi with 5 mL/min bleed during run

Sampler: Agilent 7683B, 5.0 μ L syringe (p/n G4513-80206)

Backflush: Post run 3.5 min at 75 psi Aux EPC, 2 psi inlet pressure

- | | |
|---------------------------------|----------------------------|
| 1. N-Nitrosodimethylamine | 16. Simazine |
| 2. Aniline | 17. Atrazine |
| 3. 1,4-Dichlorobenzene-d4 | 18. Pentachlorophenol |
| 4. Isophorone | 19. Terbufos |
| 5. 1,3-Dimethyl-2-nitrobenzene | 20. Chlorothalonil |
| 6. Naphthalene | 21. Phenanthrene-d10 |
| 7. Hexachlorocyclopentadiene | 22. Aldrin |
| 8. Mevinphos | 23. Heptachlor epoxide |
| 9. Acenaphthene-d10 | 24. Endrin |
| 10. 2,4-Dinitrophenol | 25. 4,4'-DDT |
| 11. 4-Nitrophenol | 26. 3,3'-Dichlorobenzidine |
| 12. 2,4-Dinitrotoluene | 27. Chrysene d-12 |
| 13. Fluorene | 28. Benzo[b]fluoranthene |
| 14. 4,6-Dinitro-2-methyl phenol | 29. Perylene-d12 |
| 15. Trifluralin | |

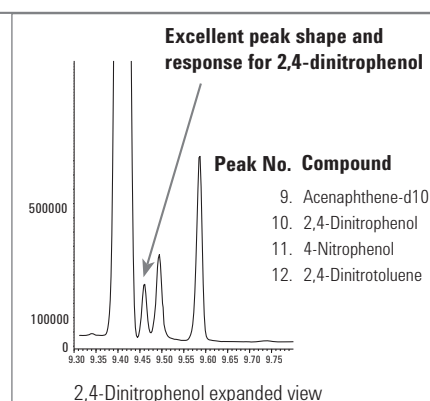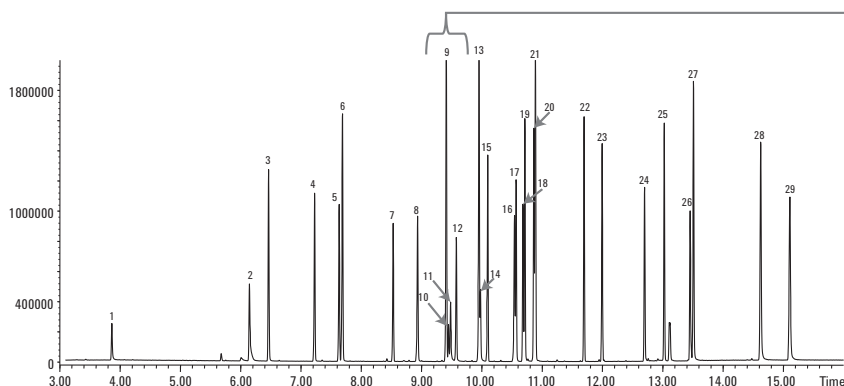

High Resolution Phenol Analysis by GC/MS

Column: VF-5ms
CP8944
30 m x 0.25 mm, 0.25 μ m


Sample Conc: Approx. 5-10 ng per component on-column

Carrier: Helium, 70 kPa

Injection: Split, 1:200, T=275 °C

Detector: Agilent Ion Trap MS

- | | |
|-----------------------------|---|
| 1. Phenol | 11. 2,3,5-Trichlorophenol |
| 2. 2-Chlorophenol | 12. 2,4,6-Trichlorophenol |
| 3. o-Cresol | 13. 2,4,5-Trichlorophenol |
| 4. m-Cresol | 14. 2,3,4-Trichlorophenol |
| 5. p-Cresol | 15. 2,3,6-Trichlorophenol |
| 6. 2-Nitrophenol | 16. 4-Nitrophenol |
| 7. 2,4-Dimethylphenol | 17. 2,4-Dinitrophenol |
| 8. 2,4-Dichlorophenol | 18. 2,3,5,6-Tetrachlorophenol |
| 9. 2,6-Dichlorophenol | 19. 2-Methyl-4,6-dinitrophenol |
| 10. 4-Chloro-3-methylphenol | 20. Pentachlorophenol |
| | 21. 2-sec-Butyl-4,6-dinitrophenol (dionseb) |

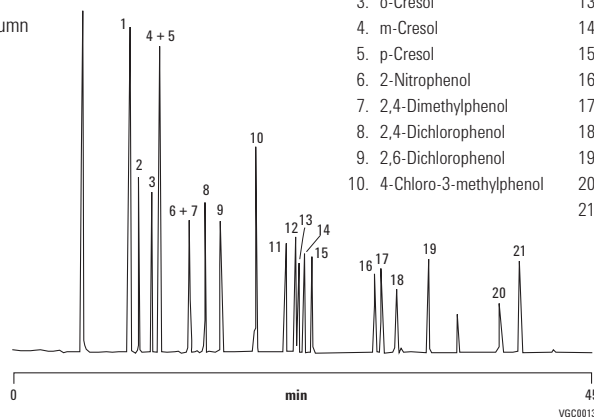

Phenols According to EPA Method 8040

Column: CP-Sil 8 CB
CP7454
50 m x 0.32 mm, 0.25 µm


Sample Conc: 1 ppm

Oven: 80 °C to 200 °C, 8 °C/min

Carrier: H₂, 150 kPa (1.5 bar, 21 psi)

Injection: Split, 100 mL/min

Detector: FID

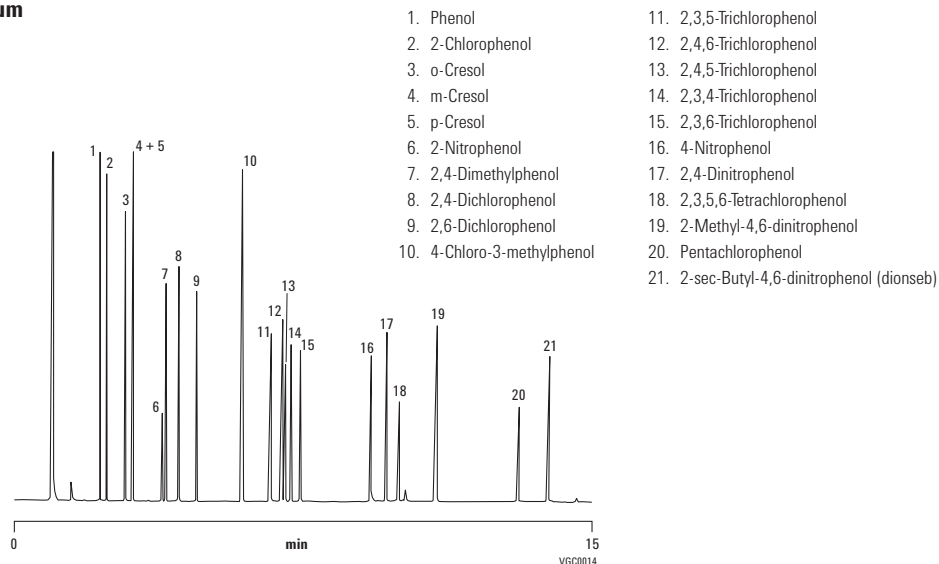**EPA Method 552.2**

Column: DB-35ms
123-3832
30 m x 0.32 mm, 0.25 µm

Column: DB-XLB
123-1236
30 m x 0.32 mm, 0.50 µm


Carrier: Helium at 45 cm/s
(EPC in constant flow mode)

Oven: 40 °C for 0.5 min
40-200 °C at 15 °C/min
200 °C for 2 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: µECD, 350 °C
Nitrogen makeup gas
(column + makeup flow =
30 mL/min constant flow)

Sample: 50 pg per component

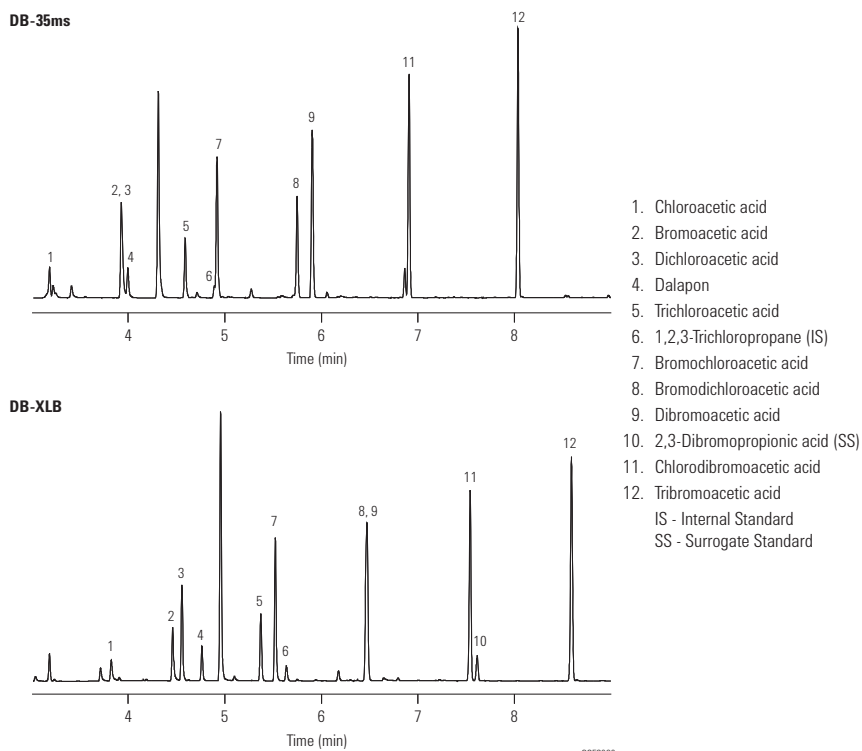**Suggested Supplies**

Septum: 11 mm Advanced Green septa,
5183-4759

Liner: Direct connect, dual taper,
deactivated, 4 mm id,
G1544-80700

Syringe: 10 µL tapered, FN 23-26s/42/HP,
5181-1267

Environmental Applications, Volatiles

Extended Analyte List for EPA Method 8021 (ELCD)

Column: DB-624
124-1374
75 m x 0.45 mm, 2.55 µm

Column: DB-VRX
124-1574
75 m x 0.45 mm, 2.55 µm

Carrier: Helium at 9 mL/min, measured at 35 °C

Oven: 35 °C for 12 min
35-60 °C at 5 °C/min
60 °C for 1 min
60-200 °C at 17 °C/min
200 °C for 5 min

Sampler: Purge and Trap (O.I.A. 4560)

Trap: VoCarb 3000

Preheat: 175 °C

Desorb: 260 °C for 1 min

Injection: J&W LVI (Low Volume Injector), 150 °C

Detector: A: PID (O.I.A. 4430), 200 °C Helium
makeup gas at 20 mL/min
B: ELCD (O.I.A. 4420), with NiCat reaction tube
in the halogen mode, 950 °C reactor temperature

Sample: 20 ppb per component in 5 mL water

Suggested Supplies


Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885


Septum: 11 mm Advanced Green septa, 5183-4759

- | | |
|------------------------------|-----------------------------------|
| 1. Dichlorodifluoromethane | 26. 1,2-Dichloropropane |
| 2. Chloromethane | 27. Dibromomethane |
| 3. Vinyl chloride | 28. Trifluorotoluene (IS) |
| 4. Bromomethane | 29. Bromodichloromethane |
| 5. Chloroethane | 30. 2-Chloroethyl vinyl ether |
| 6. Trichlorofluoromethane | 31. cis-1,3-Dichloropropene |
| 7. 2-Chloropropane (IS) | 32. Toluene |
| 8. 1,1-Dichloroethene | 33. trans-1,3-Dichloropropene |
| 9. Iodomethane | 34. 1,1,2-Trichloroethane |
| 10. Allyl chloride | 35. Tetrachloroethene |
| 11. Methylene chloride | 36. 1,3-Dichloropropane |
| 12. trans-1,2-Dichloroethene | 37. Dibromochloromethane |
| 13. 1,1-Dichloroethane | 38. 1,2-Dibromoethane |
| 14. Chloroprene | 39. 1-Chloro-3-fluorobenzene (IS) |
| 15. cis-1,2-Dichloroethene | 40. Chlorobenzene |
| 16. 2,2-Dichloropropane | 41. 1,1,1,2-Tetrachloroethane |
| 17. Bromochloromethane | 42. Ethylbenzene |
| 18. Chloroform | 43. m-Xylene |
| 19. 1,1,1-Trichloroethane | 44. p-Xylene |
| 20. Carbon tetrachloride | 45. Styrene |
| 21. 1,1-Dichloropropene | 46. o-Xylene |
| 22. Benzene | 47. Bromoform |
| 23. 1,2-Dichloroethane | 48. Isopropylbenzene |
| 24. Fluorobenzene (IS) | 49. cis-1,4-Dichlorobutene |
| 25. Trichloroethene | 50. 1,1,2,2-Tetrachloroethane |
| | 51. Bromobenzene |
| | 52. 1,2,3-Trichloropropane |
| | 53. n-Propylbenzene |
| | 54. 2-Chlorotoluene |
| | 55. 1,3,5-Trimethylbenzene |
| | 56. 4-Chlorotoluene |
| | 57. tert-Butylbenzene |
| | 58. 1,2,4-Trimethylbenzene |
| | 59. sec-Butylbenzene |
| | 60. 1,3-Dichlorobenzene |
| | 61. p-Isopropyltoluene |
| | 62. 1,4-Dichlorobenzene |
| | 63. Benzyl chloride |
| | 64. n-Butylbenzene |
| | 65. 1,2-Dichlorobenzene |
| | 66. Bis(2-chloroisopropyl) ether |
| | 67. 1,2-Dibromo-3-chloropropane |
| | 68. 1,2,4-Trichlorobenzene |
| | 69. Hexachlorobutadiene |
| | 70. Naphthalene |
| | 71. 1,2,3-Trichlorobenzene |

DB-624

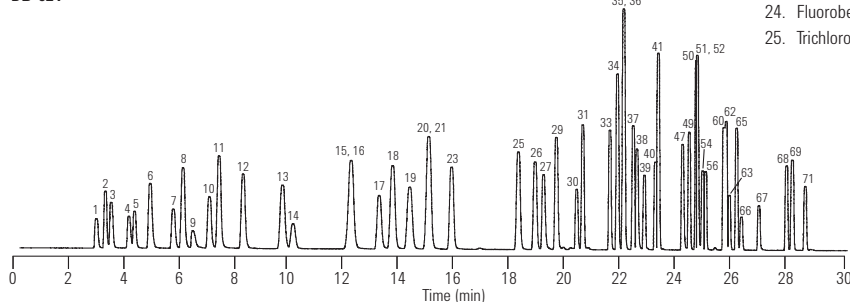

DB-VRX

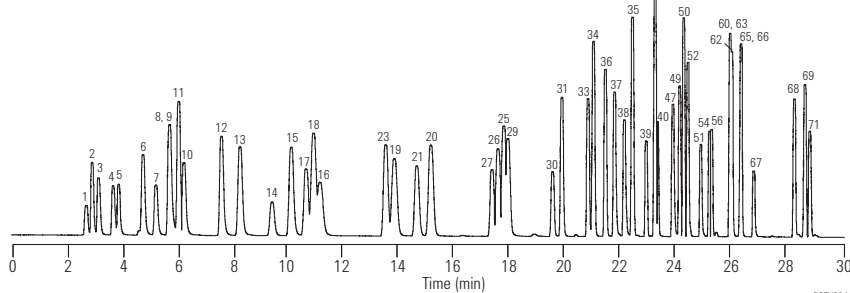

GCEV004

Fast VOC Analysis

Column: DB-624
121-1324
20 m x 0.18 mm, 1.00 µm

Carrier: Helium at 37 cm/s, (constant flow mode)

Oven: 35 °C for 4 min
35-200 °C at 15 °C/min
200 °C for 0.1 min
60-200 °C at 17 °C/min

Sampler: Purge and trap (Tekmar LSC 3000)
Purge: Helium for 11 min at 50 mL/min
Preheat: 250 °C
Desorb: 260 °C for 2 min
Line & valve: 100 °C

Detector: MSD, 250 °C transfer line
Full scan 35-260 amu
3.25 scans per s

Sample: 10 ppb per component in 25 mL water

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

- | | |
|-------------------------------|---------------------------------|
| 1. Dichlorofluoromethane | 27. Dibromochloromethane |
| 2. Chloromethane | 28. 1,2-Dibromomethane |
| 3. Vinyl chloride | 29. Chlorobenzene |
| 4. Bromomethane | 30. 1,1,1,2-Tetrachloroethane |
| 5. Chloroethane | 31. Ethylbenzene |
| 6. Trichlorofluoromethane | 32. m-Xylene |
| 7. 1,1-Dichloroethene | 33. p-Xylene |
| 8. Methylene chloride | 34. o-Xylene |
| 9. trans-1,2-Dichloroethene | 35. Bromoform |
| 10. 1,1-Dichloroethane | 36. Isopropylbenzene |
| 11. 2,2-Dichloropropane | 37. Bromofluorobenzene |
| 12. Bromochloromethane | 38. Bromobenzene |
| 13. Chloroform | 39. n-Propylbenzene |
| 14. 1,1,1-Trichloroethane | 40. 2-Chlorotoluene |
| 15. Carbon tetrachloride | 41. 1,3,5-Trimethylbenzene |
| 16. Benzene | 42. tert-Butylbenzene |
| 17. Fluorobenzene | 43. 1,2,4-Trimethylbenzene |
| 18. Trichloroethene | 44. sec-Butylbenzene |
| 19. 1,2-Dichloropropane | 45. 1,3-Dichlorobenzene |
| 20. Dibromomethane | 46. 4-Isopropyltoluene |
| 21. Bromodichloromethane | 47. 1,4-Dichlorobenzene |
| 22. cis-1,3-Dichloropropene | 48. 1,2-Dichlorobenzene |
| 23. Toluene | 49. 1,2-Dibromo-3-chloropropane |
| 24. trans-1,3-Dichloropropene | 50. 1,2,4-Trichlorobenzene |
| 25. 1,1,2-Trichloroethane | 51. Hexachlorobutadiene |
| 26. Tetrachloroethene | 52. Naphthalene |
| | 53. 1,2,3-Trichlorobenzene |

Ideal for 5975 MSD

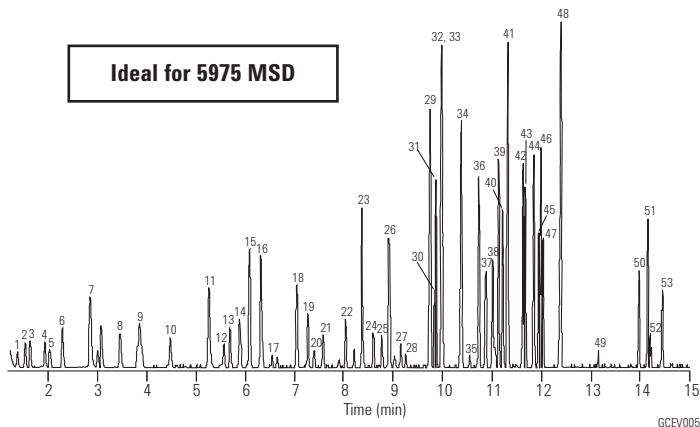

Analysis of Volatile Organic Compounds in Environmental Waters Using the Agilent 7697A Headspace and 7890B/5977A GC/MS

Column: VF-624ms
CP9103
60 m x 0.25 mm, 1.40 µm

Instrument: Agilent 7697A Headspace and 7890B/5977A GC/MS

Carrier: Helium, 11 mL/min, 160 °C

Oven: 32 °C for 2 min, then 10 °C/min to 220 °C for 5 min


Injection: Split, 4:1, 160 °C for 5 min, purge 100 mL/min for 1 min

Detector: 5977A MSD, simultaneous Scan/SIM mode


Sample: Standard VOC mix

Sample Conc: 10 µg/L

	RT, min	CAS Number		RT, min	CAS Number		RT, min	CAS Number
1. Dichlorodifluoromethane	3.387	75-71-8	11. trans-1,2-Dichloroethene	7.069	156-60-5	21. Benzene	9.440	71-43-2
2. Chloromethane	3.734	74-87-3	12. 1,1-Dichloroethane	7.644	75-34-3	22. 1,2-Dichloroethane	9.497	107-06-2
3. Vinyl chloride	3.980	75-01-4	13. Ethyl tert-butyl ether	8.091	637-92-3	23. tert-Amyl methyl ether	9.540	994-05-8
4. Bromomethane	4.390	74-83-9	14. cis-1,2-Dichloroethene	8.353	156-59-2	24. Trichloroethene	10.232	79-01-6
5. Chloroethane	4.788	75-00-3	15. 2,2-Dichloropropane	8.370	594-20-7	25. 1,2-Dichloropropane	10.576	78-87-5
6. Trichlorofluoromethane	5.202	75-69-4	16. Bromochloromethane	8.656	74-97-5	26. Dibromomethane	10.699	74-95-3
7. 1,1-Dichloroethene	5.998	75-34-4	17. Chloroform	8.756	67-66-3	27. Bromodichloromethane	10.884	75-27-4
8. Carbon disulfide	6.338	75-15-0	18. 1,1,1-Trichloroethane	8.995	71-55-6	28. cis-1,3-Dichloropropene	11.437	10061-01-5
9. Dichloromethane	6.701	75-09-2	19. 1,1-Dichloro-1-propene	9.177	563-58-6	29. Toluene	11.890	108-88-3
10. Methyl tert-butyl ether	7.046	1634-04-4	20. Carbon tetrachloride	9.189	56-23-5	30. trans-1,3-Dichloropropene	12.165	10061-02-6
						31. 1,1,2-Trichloroethane	12.443	79-00-5
						32. Tetrachloroethene	12.580	127-18-4
						33. 1,3-Dichloropropane	12.673	142-28-9
						34. Dibromochloromethane	12.981	124-48-1
						35. 1,2-Dibromoethane	13.175	106-93-4
						36. Chlorobenzene	13.830	108-90-7
						37. 1,1,1,2-Tetrachloroethane	13.939	630-20-6
						38. Ethylbenzene	13.934	100-41-4
						39. m and p-Xylene	14.115	108-38-3 & 106-42-3
						40. o-Xylene	14.669	95-47-6
						41. Styrene	14.699	100-42-5
						42. Bromoform	14.994	75-25-2
						43. Isopropylbenzene	15.183	98-82-8
						44. 1,1,2,2-Tetrachloroethane	15.612	79-34-5
						45. Bromobenzene	15.697	108-86-1
						46. 1,2,3-Trichloropropane	15.731	96-18-4
						47. n-Propylbenzene	15.793	103-65-1
						48. 2-Chlorotoluene	15.952	95-49-8
						49. 3-Chlorotoluene	16.042	108-41-8
						50. 1,3,5-Trimethylbenzene	16.048	108-67-8
						51. 4-Chlorotoluene	16.133	106-43-4
						52. tert-Butylbenzene	16.526	98-06-6
						53. 1,2,4-Trimethylbenzene	16.608	95-63-6
						54. sec-Butylbenzene	16.856	135-98-8
						55. 1,3-Dichlorobenzene	17.071	541-73-1
						56. 4-Isopropyltoluene	17.077	99-87-6
						57. 1,4-Dichlorobenzene	17.220	106-46-7
						58. 1,2,3-Trimethylbenzene	17.231	526-73-8
						59. n-Butylbenzene	17.689	104-51-8
						60. 1,2-Dichlorobenzene	17.761	95-50-1
						61. 1,2-Dibromo-3-chloropropane	18.949	96-12-8
						62. 1,3,5-Trichlorobenzene	19.215	108-70-3
						63. 1,2,4-Trichlorobenzene	20.179	120-82-1
						64. Hexachlorobutadiene	20.370	87-68-3
						65. Naphthalene	20.604	91-20-3
						66. 1,2,3-Trichlorobenzene	20.922	87-61-6

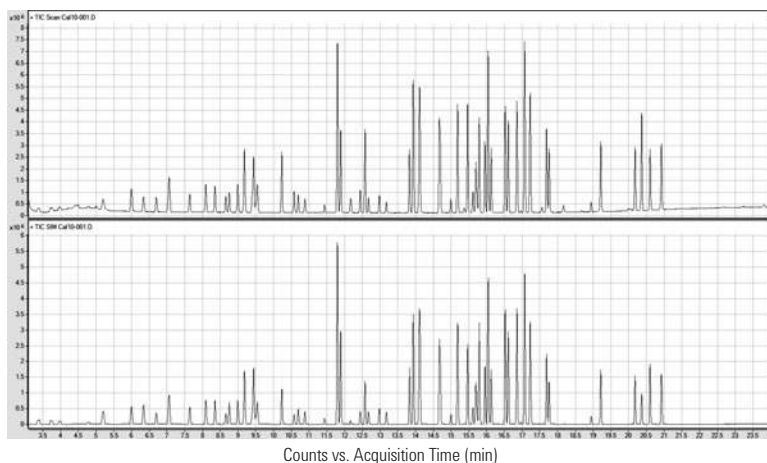

10 µg/L VOC Standard Scan and SIM Traces

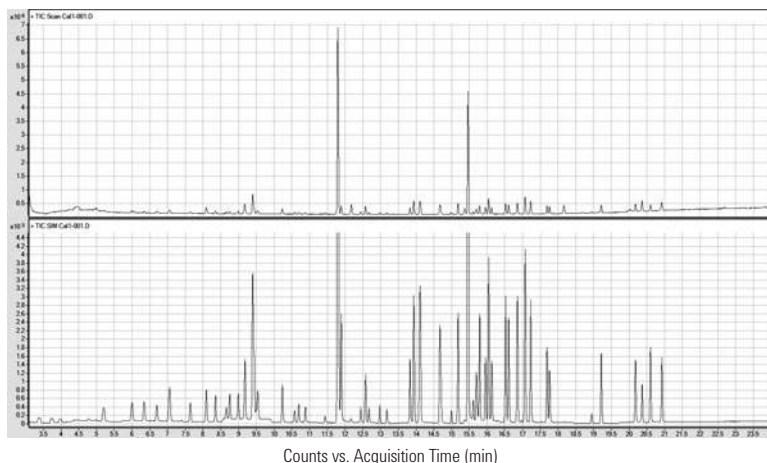

1 µg/L VOC Standard Scan and SIM Traces

EPA Method 551

Column: DB-1
122-1033
30 m x 0.25 mm, 1.00 µm


Carrier: Helium at 24.8 cm/s,
measured at 150 °C

Injection: Splitless, 200 °C
15 s purge activation time

Oven: 35 °C for 9 min
35-40 °C at 10 °C/min
40 °C for 3 min
40-150 °C at 6 °C/min
150 °C for 1 min

Detector: ECD, 300 °C
Sample: 1 µL of 50 pg/µL,
AccuStandard

- | | |
|-----------------------------|---------------------------------|
| 1. Chloroform | 10. Chloropicrin |
| 2. 1,1,1-Trichloroethane | 11. Dibromochloromethane |
| 3. Carbon tetrachloride | 12. Bromochloroacetonitrile |
| 4. Trichloroacetonitrile | 13. 1,2-Dibromoethane |
| 5. Dichloroacetonitrile | 14. Tetrachloroethene |
| 6. Bromodichloromethane | 15. 1,1,1-Trichloropropanone |
| 7. Trichloroethene | 16. 1,2-Dibromopropane (IS) |
| 8. Chloral hydrate | 17. Bromoform |
| 9. 1,1-Dichloro-2-propanone | 18. Dibromoacetonitrile |
| | 19. 1,2-Dibromo-3-chloropropane |

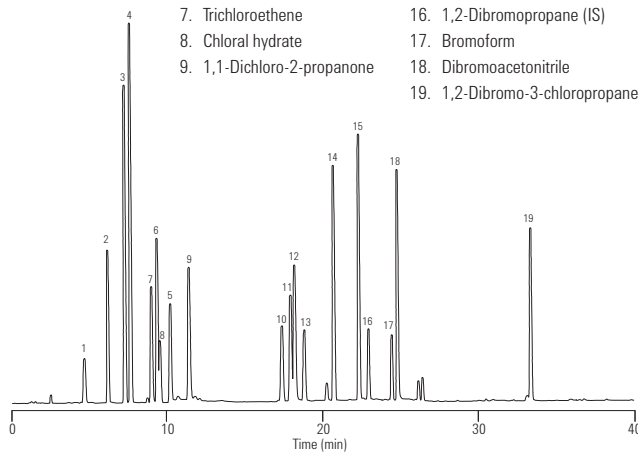

GCEV008

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316
Seal: Gold plated seal, 18740-20885
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

European Red List Volatiles

Column: DB-5.625
122-5632
30 m x 0.25 mm, 0.50 µm

Column: DB-624
122-1334
30 m x 0.25 mm, 1.40 µm

Carrier: Helium at 35 cm/s, measured at 40 °C

Injection: Split, 250 °C
Split ratio 1:50

Oven: 40 °C for 2 min
40-140 °C at 12 °C/min


Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of headspace of neat mixture


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct, 1.5 mm id, 18740-80200
Seal: Gold plated seal, 18740-20885

DB-5.625

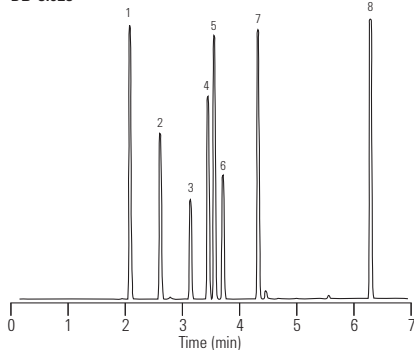

DB-624

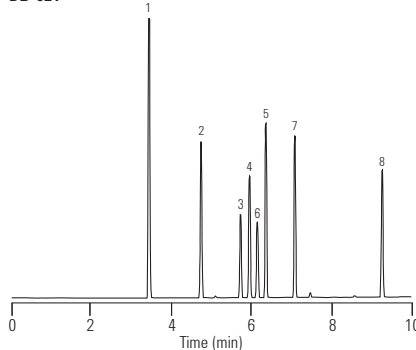

GCEV010

1. 1,1-Dichloroethylene
2. 1,1-Dichloroethane
3. Chloroform
4. 1,1,1-Trichloroethane
5. 1,2-Dichloroethane
6. Carbon tetrachloride
7. Trichloroethylene
8. Tetrachloroethylene

EPA Volatiles by GC/MS (Split Injector)

Column: DB-VRX
122-1564
60 m x 0.25 mm, 1.40 µm

Carrier: Helium at 30 cm/s, measured at 45 °C

Oven: 45 °C for 10 min
45-190 °C at 12 °C/min
190 °C for 2 min
190-225 °C at 6 °C/min
225 °C for 1 min

Sampler: Purge and trap (O.I.A. 4560)
Purge: Helium for 11 min at 40 mL/min
Trap: Tenax/Silica Gel/Carbosieve
Preheat: 175 °C
Desorb: 220 °C for 0.6 min

Injection: Split, 110 °C
Split flow 30 mL/min


Detector: MSD, 235 °C transfer line
Full scan 35-260 amu (m/z 44 subtracted)

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal kit, 5188-5367

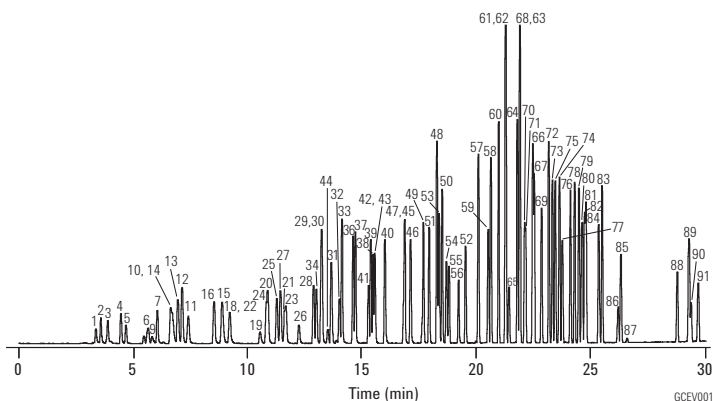

Column: DB-624
122-1364
60 m x 0.25 mm, 1.40 µm

Carrier: Helium at 31 cm/s, measured at 40 °C

Oven: 45 °C for 3 min
45-90 °C at 8 °C/min
90 °C for 4 min
90-200 °C at 6 °C/min
200 °C for 5 min

Sampler: Purge and trap (O.I.A. 4560)
Purge: Helium for 11 min at 40 mL/min
Trap: Tenax/Silica Gel/Carbosieve
Preheat: 175 °C
Desorb: 220 °C for 0.6 min

Injection: Split, 110 °C
Split flow 30 mL/min


Detector: MSD, 235 °C transfer line
Full scan 35-260 amu (m/z 44 subtracted)

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal kit, 5188-5367

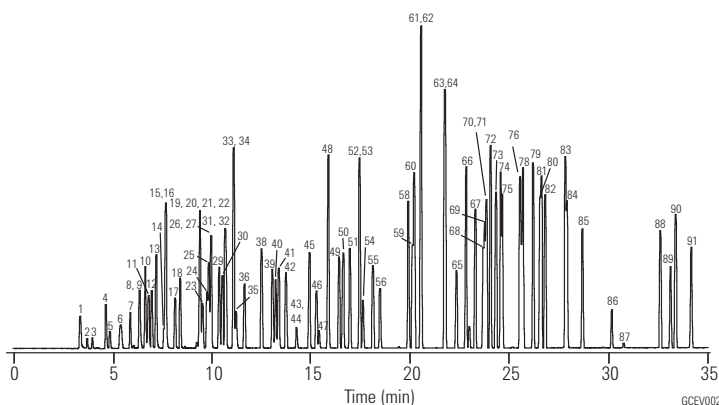

- | | | | | |
|------------------------------|------------------------------|-----------------------------------|---------------------------------|---------------------------------|
| 1. Dichlorodifluoromethane | 20. cis-1,2-Dichloroethene | 39. 1,2-Dichloropropane | 58. Chlorobenzene | 77. Pentachloroethane |
| 2. Chloromethane | 21. 2,2-Dichloropropane | 40. Methyl methacrylate | 59. 1,1,1,2-Tetrachloroethane | 78. 1,2,4-Trimethylbenzene |
| 3. Vinyl chloride | 22. Propionitrile | 41. Dibromomethane | 60. Ethylbenzene | 79. sec-Butylbenzene |
| 4. Bromomethane | 23. Methyl acrylate | 42. Bromodichloromethane | 61. m-Xylene | 80. 1,3-Dichlorobenzene |
| 5. Chloroethane | 24. Methacrylonitrile | 43. 2-Nitropropane | 62. p-Xylene | 81. p-Isopropyltoluene |
| 6. Trichlorofluoromethane | 25. Bromochloromethane | 44. Chloroacetone | 63. o-Xylene | 82. 1,4-Dichlorobenzene |
| 7. Diethyl ether | 26. Tetrahydrofuran | 45. cis-1,3-Dichloropropene | 64. Styrene | 83. n-Butylbenzene |
| 8. 1,1-Dichloroethene | 27. Chloroform | 46. 4-Methyl-2-pentanone | 65. Bromoform | 84. 1,2-Dichlorobenzene |
| 9. Acetone | 28. Pentafluorobenzene (IS) | 47. 1,1-Dichloro-2-propanone | 66. Isopropylbenzene | 85. Hexachloroethane |
| 10. Iodomethane | 29. 1,1,1-Trichloroethane | 48. Toluene | 67. 4-Bromofluorobenzene (SS) | 86. 1,2-Dibromo-3-chloropropane |
| 11. Carbon disulfide | 30. 1-Chlorobutane | 49. trans-1,3-Dichloropropene | 68. 1,1,2,2-Tetrachloroethane | 87. Nitrobenzene |
| 12. Allyl chloride | 31. 1,1-Dichloropropene | 50. Ethyl methacrylate | 69. Bromobenzene | 88. 1,2,4-Trichlorobenzene |
| 13. Methylene chloride | 32. Carbon tetrachloride | 51. 1,1,2-Trichloroethane | 70. 1,2,3-Trichloropropane | 89. Hexachlorobutadiene |
| 14. Acrylonitrile | 33. Benzene | 52. Tetrachloroethene | 71. trans-1,4-Dichloro-2-butene | 90. Naphthalene |
| 15. Methyl-tert-butyl ether | 34. 1,2-Dichloroethane | 53. 1,3-Dichloropropane | 72. n-Propylbenzene | 91. 1,2,3-Trichlorobenzene |
| 16. trans-1,2-Dichloroethene | 35. 2,2-Dimethylhexane | 54. 2-Hexanone | 73. 2-Chlorotoluene | |
| 17. Hexane | 36. Fluorobenzene (IS) | 55. Dibromochloromethane | 74. 1,3,5-Trimethylbenzene | |
| 18. 1,1-Dichloroethane | 37. 1,4-Difluorobenzene (IS) | 56. 1,2-Dibromoethane | 75. 4-Chlorotoluene | |
| 19. 2-Butanone | 38. Trichloroethene | 57. 1-Chloro-3-fluorobenzene (IS) | 76. tert-Butylbenzene | |

Note: Some compounds not present in both chromatograms

Environmental Applications, Air Analysis

EPA Air Analysis Compendium Method TO-14 Standard

Column: DB-1
123-1063
60 m x 0.32 mm, 1.00 μ m

Carrier: Helium at 25 cm/s measured off of CO₂ at 35 °C
constant flow mode

Oven: 35 °C for 5 min
35-120 °C at 5 °C/min
120-220 °C at 30 °C/min
220 °C for 5 min

Injection: Entech 7100 cryogenic sample preconcentrator

Detector: MSD
Full scan of m/z 40-250

Sample: 400 mL of a 10 ppbV TO-14 standard
and 100 mL of a 20 ppbv IS/SS standard


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

- | | |
|---|---------------------------------------|
| 1. CO ₂ | 14. Bromochloromethane (IS) |
| 2. Freon 12 (dichlorodifluoromethane) | 15. Chloroform |
| 3. Chloromethane | 16. 1,2-Dichloroethane |
| 4. Freon 114 (1,2-dichloro-1,1,2,2-tetrafluoroethane) | 17. 1,1,1-Trichloroethane |
| 5. Vinyl chloride | 18. Benzene |
| 6. Bromomethane | 19. Carbon tetrachloride |
| 7. Chloroethane | 20. 1,4-Difluorobenzene (IS) |
| 8. Freon 11 (trichlorofluoromethane) | 21. 1,2-Dichloropropane |
| 9. 1,1-Dichloroethene | 22. Trichloroethene |
| 10. Methylene chloride | 23. cis-1,3-Dichloropropene |
| 11. Freon 113 (1,1,2-trichloro-1,2,2-trifluoroethane) | 24. trans-1,3-Dichloropropene |
| 12. 1,1-Dichloroethane | 25. 1,1,2-Trichloroethane |
| 13. cis-1,2-Dichloroethene | 26. Toluene-d ₈ (SS) |
| | 27. Toluene |
| | 28. 1,2-Dibromoethane |
| | 29. Tetrachloroethene |
| | 30. Chlorobenzene-d ₅ (SS) |
| | 31. Chlorobenzene |
| | 32. Ethylbenzene |
| | 33. m-Xylene |
| | 34. p-Xylene |
| | 35. Styrene |
| | 36. 1,1,2,2-Tetrachloroethane |
| | 37. o-Xylene |
| | 38. 4-Bromofluorobenzene (SS) |
| | 39. 1,3,5-Trimethylbenzene |
| | 40. 1,2,4-Trimethylbenzene |
| | 41. 1,3-Dichlorobenzene |
| | 42. 1,2-Dichlorobenzene |
| | 43. 1,4-Dichlorobenzene |
| | 44. 1,2,4-Trichlorobenzene |
| | 45. 1,2-Dibromobenzene (IS) |
| | 46. Hexachloro-1,3-butadiene |

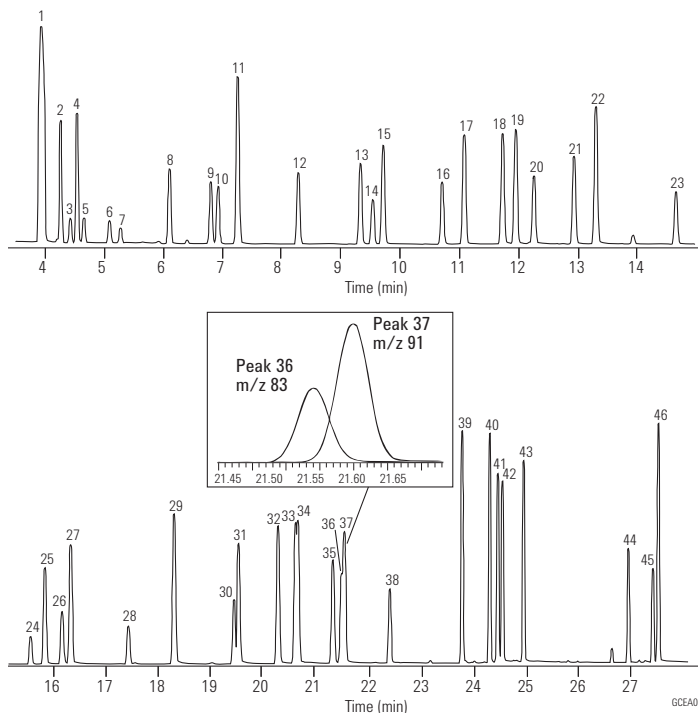

Agilent wishes to thank Entech Instruments for providing this chromatogram.

Formaldehyde, 50 ppb

Column: DB-5ms
123-5563
60 m x 0.32 mm, 1.00 µm

Carrier: Helium, 1.5 mL/min

Oven: 35 °C for 5 min
35-85 °C at 10 °C/min

Sampler: Entech 7100 cryogenic sample preconcentrator

Detector: GC/MS 6890/5973N
Scan 29-180 amu 0-6 min
33-280 amu 6-30 min
Electron impact 70 eV


Sample: 100 cc 50 ppb Formaldehyde/20 ppb others

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

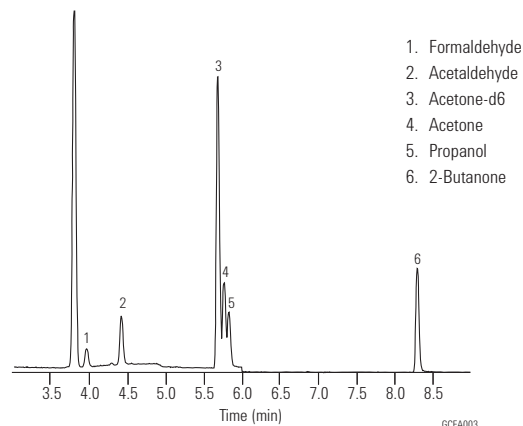

Agilent wishes to thank Entech Instruments for providing this chromatogram.

Sulfur in Air

Column: DB-5ms
123-5563
60 m x 0.32 mm, 1.00 µm


Carrier: Helium, 1.5 mL/min

Oven: 35 °C for 5 min
35-140 °C at 6 °C/min
140-220 °C at 15 °C/min
220 °C for 3 min

Sampler: Entech 7100 cryogenic sample preconcentrator

Detector: GC/MS 6890/5973N
Scan 29-180 amu 0-6 min
33-280 amu 6-30 min
Electron impact 70 eV

Sample: 400 cc 10 ppb sulfurs

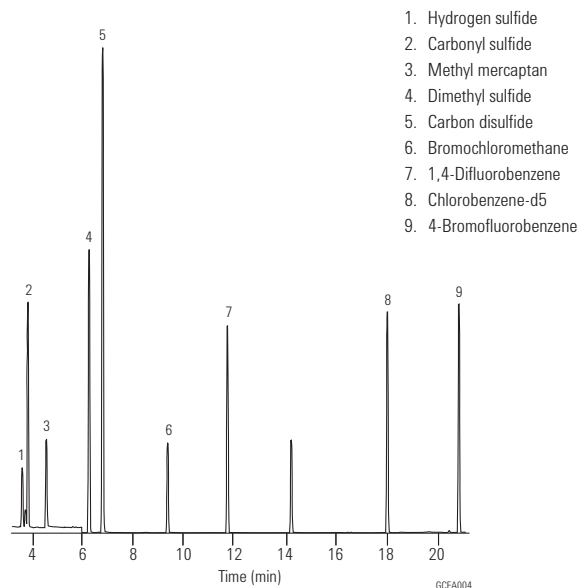

Agilent wishes to thank Entech Instruments for providing this chromatogram.

N₂O I

Column: HP-PLOT Q
19095P-Q04
30 m x 0.53 mm, 40.00 µm


Carrier: Helium, 5 psi (approximately 8 mL/min)

Oven: 35 °C isothermal

Injection: 250 µL, injected
Split ratio 1:3

Detector: TCD, 200 °C

Sample: Approximately 200 ppmv methane
200 ppmv CO₂
250 ppmv N₂O (nitrogen balance gas)

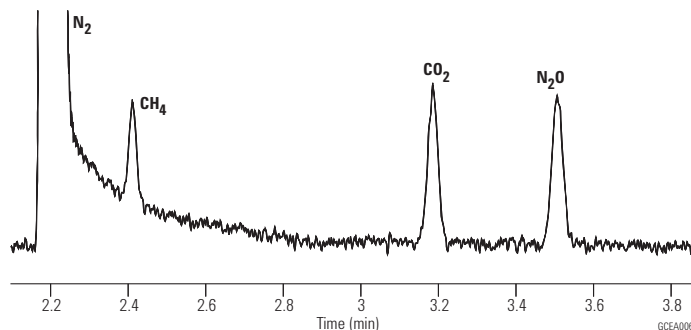**N₂O II**

Column: HP-PLOT Molesieve
19095P-MS6
30 m x 0.53 mm, 25.00 µm


Carrier: Helium, 6 psi (approximately 10 mL/min)

Oven: 50 °C (5 min), 25 °C/min to 200 °C and hold

Injection: 250 µL injected
Split ratio 1:4

Detector: TCD, 250 °C
Column compensation on

Sample: Approximately 200 ppmv methane
200 ppmv CO₂
250 ppmv N₂O (nitrogen balance gas)

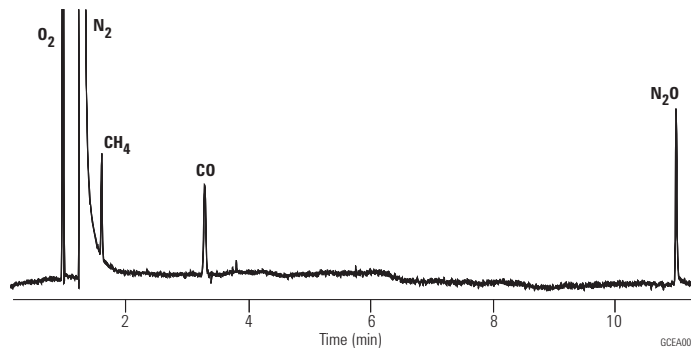**N₂O III**

Column: GS-CarbonPLOT
113-3133
30 m x 0.32 mm, 3.00 µm


Carrier: Helium, 12 psi (approximately 3 mL/min)

Oven: 35 °C isothermal

Injection: 250 µL injected
Split ratio 1:4

Detector: TCD, 200 °C


Sample: Approximately 200 ppmv methane
200 ppmv CO₂
250 ppmv N₂O (nitrogen balance gas)

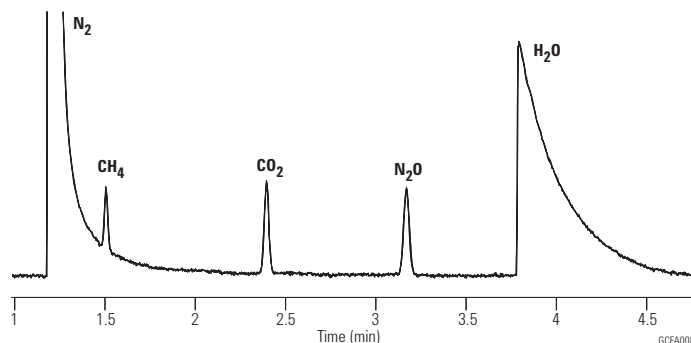

Permanent Gases on a Thick Film Molsieve Column


Column: CP-Molsieve 5Å
CP7538
25 m x 0.53 mm, 50.00 µm

Sample: 10 µL
Sample Conc: % range
Carrier: H₂
Oven: 30 °C
Injection: Split, 100 mL/min
Detector: TCD

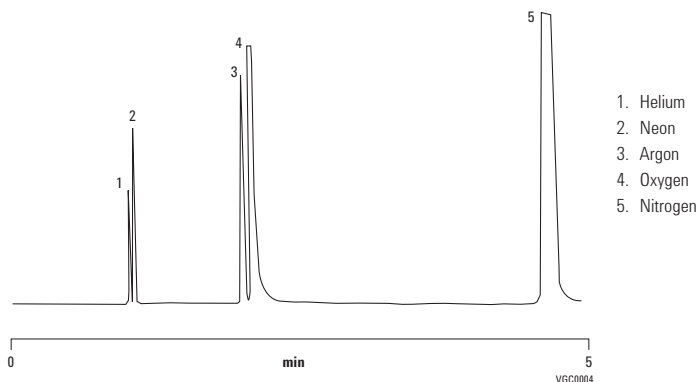**Fast Analysis of Permanent Gases and CO₂ using Tandem PLOT Columns**

Column: Select for Permanent Gases/CO₂
CP7429

Sample: 10 µL
Sample Conc: % level
Carrier: H₂, 60 kPa
Oven: 45 °C
Injection: Split, 50 mL/min
Detector: µ-TCD

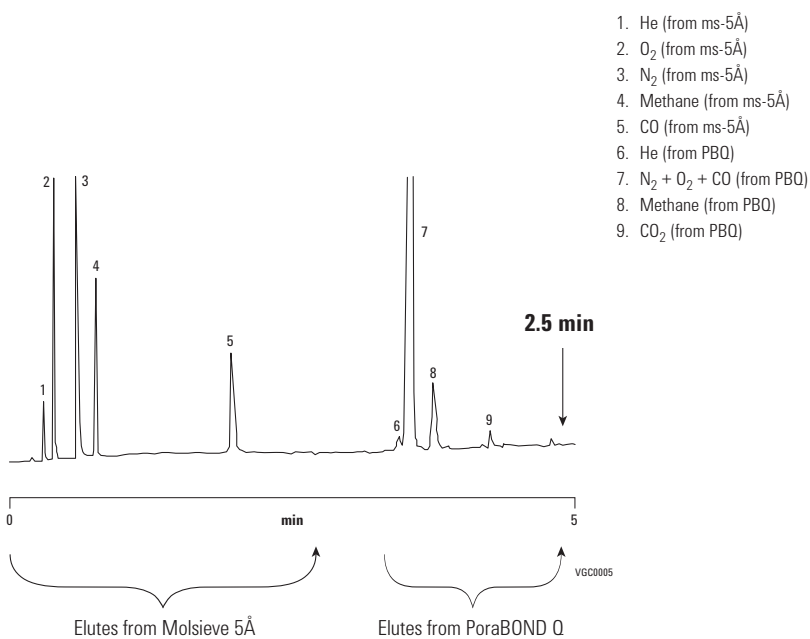

**EPA Air Analysis Method TO-15
(1 ppbv standard)**

Column: DB-5ms
123-5563
60 m x 0.32 mm, 1.00 µm

Carrier: Helium, 1.5 mL/min

Oven: 35 °C for 5 min
35-140 °C at 6 °C/min
140-220 °C at 15 °C/min
220 °C for 3 min

Sampler: Entech 7100 cryogenic sample preconcentrator


Detector: GC/MS 6890/5973N
Scan 29-180 amu 0-6 min
33-280 amu 6-30 min
Electron impact 70 eV

Sample: 400 mL sample load
All compounds at 10 ppbv except formaldehyde (50 ppbv),
acetaldehyde (20 ppbv), propanol (20 ppbv), acetone (30 ppbv),
2-butanone (30 ppbv)

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct, 1.5 mm id, 18740-80200
Seal: Gold plated seal, 18740-20885

Quantitation Ion		Quantitation Ion		Quantitation Ion	
1. Formaldehyde	30	26. n-Hexane	57	51. Tetrachloroethene	166
2. Propene	41	27. cis-1,2-Dichloroethene	96	52. 1,2-Dibromoethane	107
3. Dichlorodifluoromethane	85	28. Ethyl acetate	43	53. Chlorobenzene-d5 (IS)	117
4. Chloromethane	50	29. Bromochloromethane (IS)	128	54. Chlorobenzene	112
5. Dichlorotetrafluoroethane	85	30. Chloroform	83	55. Ethylbenzene	91
6. Acetaldehyde	29	31. Tetrahydrofuran	42	56. m-Xylene	91
7. Vinyl chloride	62	32. 1,1,1-Trichloroethane	97	57. p-Xylene	91
8. 1,3-Butadiene	39	33. 1,2-Dichloroethane	62	58. Styrene	104
9. Bromomethane	94	34. Benzene	78	59. o-Xylene	91
10. Chloroethane	64	35. Carbon tetrachloride	117	60. Bromoform	173
11. Bromoethene	106	36. Cyclohexane	56	61. 1,1,2,2-Tetrachloroethane	83
12. Trichlorofluoromethane	101	37. 1,4-Difluorobenzene (IS)	114	62. 4-Bromofluorobenzene	95
13. Acetone	58	38. 2,2,4-Trimethylpentane (isooctane)	57	63. 4-Ethyltoluene	105
14. Propanal	29	39. n-Heptane	41	64. 1,3,5-Trimethylbenzene	105
15. Isopropyl alcohol	45	40. Trichloroethene	130	65. 1,2,4-Trimethylbenzene	105
16. 1,1-Dichloroethene	61	41. 1,2-Dichloropropane	63	66. 1,3-Dichlorobenzene	146
17. 1,1,2-Trichloro-1,2,2-trifluoroethane	101	42. 1,4-Dioxane	88	67. Benzyl chloride	91
18. Methylene chloride	49	43. Bromodichloromethane	83	68. 1,4-Dichlorobenzene	146
19. 3-Chloro-1-propene (allyl chloride)	76	44. 4-Methyl-2-pentanone (MIBK)	43	69. 1,2-Dichlorobenzene	146
20. Carbon disulfide	76	45. cis-1,3-Dichloropropene	75	70. 1,2,4-Trichlorobenzene	180
21. trans-1,2-Dichloroethene	96	46. trans-1,3-Dichloropropene	75	71. Hexachlorobutadiene	225
22. tert-Butyl methyl ether (MTBE)	73	47. Toluene	91		
23. 1,1-Dichloroethane	63	48. 1,1,2-Trichloroethane	97		
24. Vinyl acetate	43	49. 2-Hexanone	43		
25. 2-Butanone (MEK)	72	50. Dibromochloromethane	129		

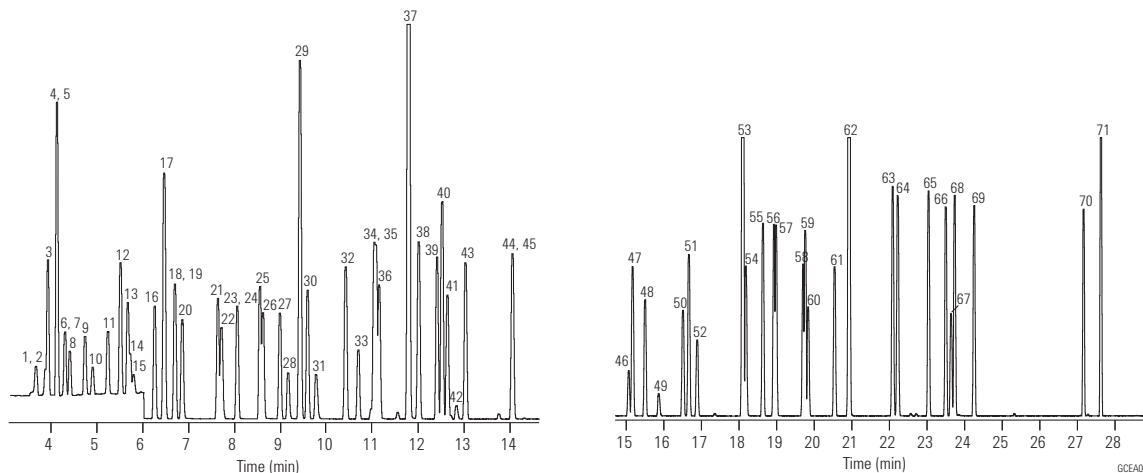

Agilent wishes to thank Entech Instruments for providing this chromatogram.

Food, Flavor, and Fragrance Applications

DB-624UI 1 µL/L Fermented Beverage Standard Mix

Column: DB-624 Ultra Inert

123-1334UI

30 m x 0.32 mm, 1.80 µm

Carrier: Helium, 2.3 mL/min, constant flow set at 35 °C

Oven: 35 °C for 5 min
10 °C/min to 100 °C for 1.5 min
15 °C/min to 220 °C for 3.0 min
25 °C/min to 250 °C for 2.8 min


Inlet: Split/splitless, 220 °C, 1 µL, split 20:1

MSD Restrictor: Scan mode 30-400 amu, source temp 230 °C,
quad temp 150 °C, transfer line temp 260 °C

Instrument: Agilent 7890/5975C equipped with MMI and FID

Sampler: Agilent 7697A headspace with 111 position tray,
1 mL sample loop

- | | |
|----------------------------|-----------------------------|
| 1. Acetyl aldehyde | 18. Isoamyl alcohol |
| 2. Methanol | 19. Active amyl alcohol |
| 3. Ethanol | 20. Isobutyl acetate |
| 4. Acetone | 21. 1-Pentanol |
| 5. Isopropanol | 22. Ethyl butanoate |
| 6. Isobutyl aldehyde | 23. Hexanal |
| 7. 1-Propanol | 24. Isoamyl acetate |
| 8. Butyl aldehyde | 25. Active amyl acetate |
| 9. 2,3 Butanedione (VDK) | 26. 1-Hexanol |
| 10. Ethyl acetate | 27. Heptanal |
| 11. 2-Butanol | 28. Octanal |
| 12. Isobutyl alcohol | 29. 1,3,5-Trioxane impurity |
| 13. 1-Butanol | 30. 1,3,5-Trioxane impurity |
| 14. 2,3 Pentanedione (VDK) | 31. Ethyl caprylate |
| 15. Ethyl propanoate | 32. 1-Phenyl ethyl acetate |
| 16. Propyl acetate | 33. Benzaldehyde, 3 methoxy |
| 17. 3-Pentanol | 34. Ethyl caprate |

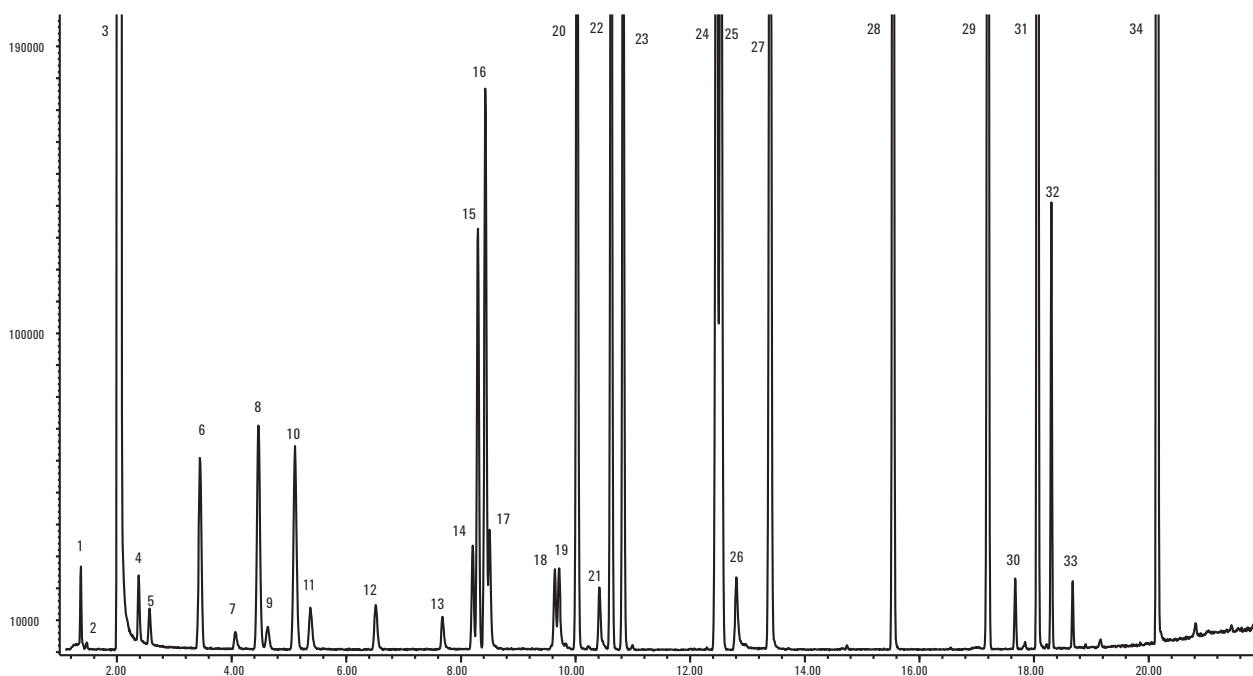

Spearmint Oil

Column A: DB-1
122-1032
30 m x 0.25 mm, 0.25 μ m


Column B: DB-1
121-1022
20 m x 0.18 mm, 0.18 μ m

Carrier: A: Helium 25 cm/s measured at 40 °C
B: Hydrogen 47 cm/s measured at 40 °C

Oven: A: 40 °C hold 1 min, 5 °C/min to 290 °C
B: 40 °C hold 0.38 min, 13 °C/min to 290 °C
hold 13.09 min


Injection: 250 °C, Split 40:1, 1 μ L injection

**Original method with a
DB-1, 30 m x 0.25 mm, 0.25 μ m column
and helium carrier**

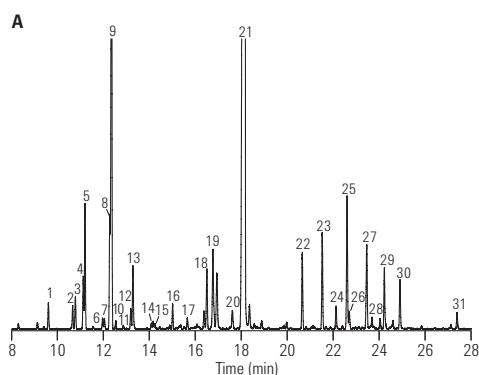

1. α -Pinene
2. Sabinene
3. β -Pinene
4. 3-Octanol
5. Myrcene
6. α -Terpinene
7. *p*-Cymene
8. 1,8-Cineol
9. Limonene
10. *cis*-Ocimene
11. *trans*-Ocimene
12. γ -Terpinene
13. *trans*-Sabinene hydrate
14. Terpinolene
15. Linalool
16. 3-Octyl acetate
17. Isomenthone
18. Terpinen-4-ol
19. Dihydro carvone
20. *trans*-Carveol
21. *l*-Carvone
22. *trans*-Dihydro carveol acetate
23. *cis*-Carvyl acetate
24. *cis*-Jasmone
25. β -Bourbonene
26. α -Bourbonene
27. β -Caryophyllene
28. α -Copaene
29. *trans*- β -Farnesene
30. Germacrene-d
31. Viridiflorol

**Faster method with a high efficiency
DB-1, 20 m x 0.18 mm, 0.18 μ m column
and hydrogen carrier**

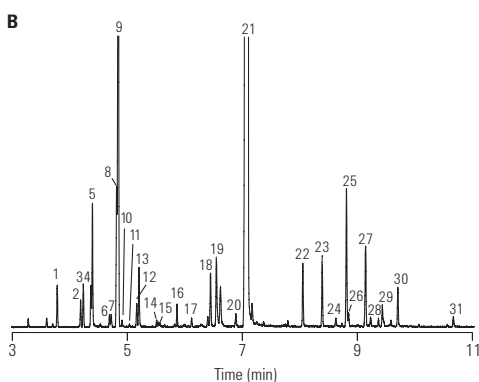

Using hydrogen as a carrier gas in conjunction with the high efficiency column resulted in an overall speed gain of 61% compared to the original method. In addition, the resolution was well maintained throughout the method translation process.

SPEARMINT

Lavender Oil Characterization

Column: DB-1ms Ultra Inert
122-0132UI
30 m x 0.25 mm, 0.25 µm

Instrument: Agilent 7890A/5975B MSD
and a 6890N FID equipped

Sampler: Agilent 7683B, 5.0 µL syringe (p/n 5188-5246),
1.0 µL injection

Carrier: Helium 40 cm/s, constant flow MSD system,
35 cm/s FID system

Inlet: 200:1 split

Oven: 62 °C 12.5 min hold, 3 °C/min to 92 °C,
then 5 °C/min to 165 °C,
then 100 °C/min to 310 °C, 2.5 min hold


Detector: MSD source at 300 °C, quadrupole at 180 °C,
transfer line at 280 °C, scan range 45-450 amu

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Single taper, MS certified liner with restriction to hold glass wool, 5188-6576

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

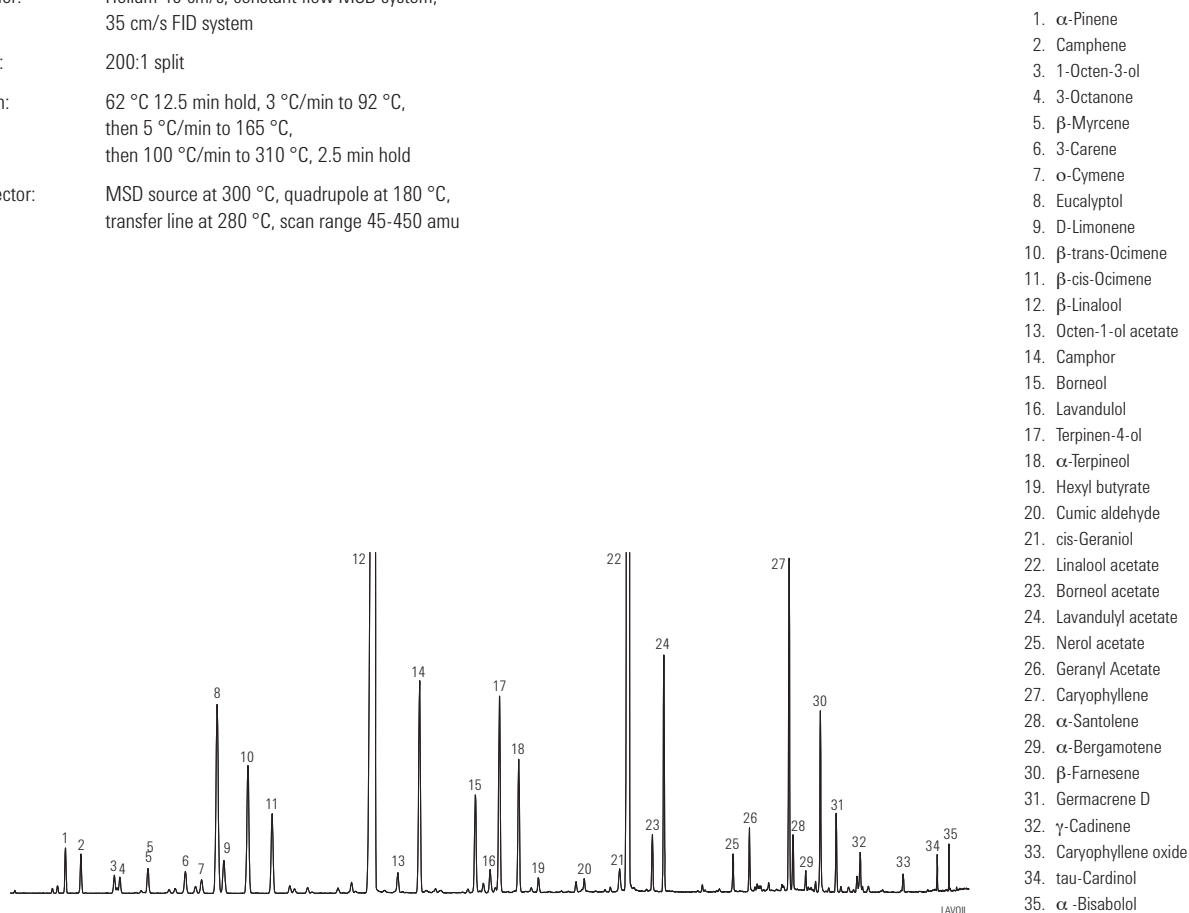

GC/MS total ion chromatogram of lavender oil sample on an Agilent J&W DB-1ms Ultra Inert 30 m x 0.25 mm, 0.25 µm capillary GC column (p/n 122-0132UI). The well-resolved, sharp peaks observed on the column ensure reliable analysis and fingerprinting of lavender oils.

Essential Oils

Column: DB-WAX
121-7022
20 m x 0.18 mm, 0.18 μ m


Carrier: Hydrogen at 44.3 cm/s
Measured at 45 °C

Oven: 45 °C hold 0.77 min
7.79 °C/min to 250 °C

Injection: Split 1:30, 250 °C
1 μ L of 1:35 oil in acetone


Detector: MSD full scan at m/z 40-500
250 °C transfer line

Wild chamomile

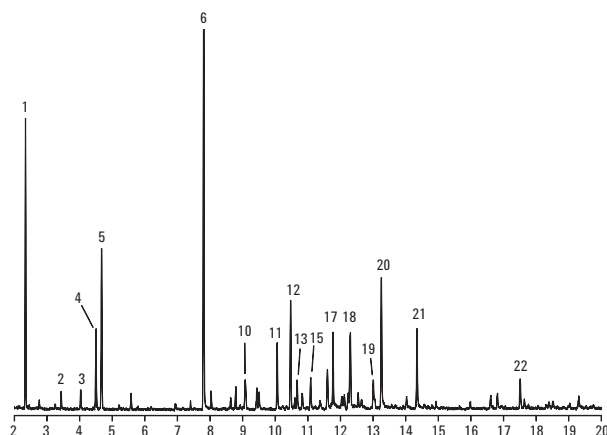

1. α -Pinene
2. β -Pinene
3. β -Myrcene
4. D-Limonene
5. Eucalyptol
6. 2,4-Hexadienal
7. Menthone
8. γ -Terpinene
9. Menthofuran
10. Iso-menthone
11. Δ -Carane
12. Bornyl acetate
13. β -Caryophyllene
14. Isomenthol
15. Citronellyl formate
16. Menthol
17. t- β -Farnesene
18. γ -Cadinene
19. δ -Cadinene
20. Citronellol
21. Nerol
22. β -Maaliene

Peppermint

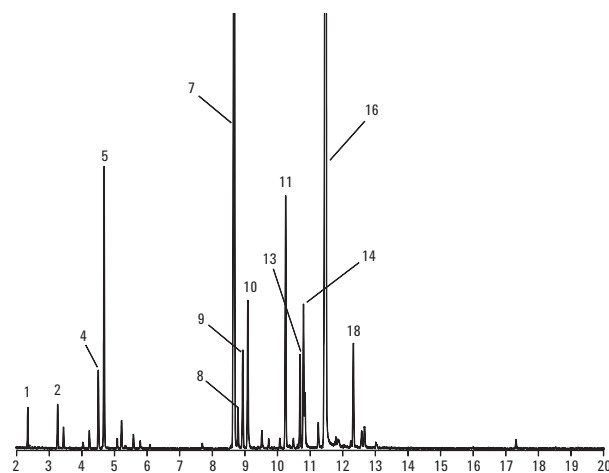

Fragrance Reference Standard

Column: DB-1
122-1032
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 25 cm/s, measured at 150 °C

Oven: 40 °C for 1 min
40-290 °C at 5 °C/min

Injection: Split, 250 °C
Split ratio 1:50

Detector: MSD, 300 °C transfer line


Sample: 1 µL of a 1:20 dilution of neat sample in acetone

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Split, single taper, low pressure drop, glass wool, 5183-4647
Seal: Gold plated seal, 18740-20885
Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

Many thanks to Carl Frey, Manager of Analytical Services, Dragoco, and Kevin Myung, Director of Flavor and Perfumery Research, Bush Boake Allen, Inc. for contributing to this work.

- | | | | |
|--|-----------------------------|----------------------------|--|
| 1. Acetone | 26. Hexyl acetate | 53. Ethyl octanoate | 90. Eugenyl acetate |
| 2. 2,3-Butanedione (diacetyl) | cis-Linalool oxide | 54. Octyl acetate | 91. Frambinone (raspberry ketone) |
| 3. Ethyl acetate | Methyl benzoate | 56. Fenchyl acetate | 93. Isoamyl salicylate |
| 4. 2,3-Pentanedione (acetyl propionyl) | trans-Linalool oxide | 57. Citronellol | 94. δ-Cadinene |
| 5. Ethyl propionate | 28. Methyl-cresol | 58. Neral | 95. cis-Nerolidol |
| 6. Methyl butyrate | 29. Benzyl alcohol | 59. Carvone | 96. Rosatol (rosetone) |
| 7. 3-Methylbutyl alcohol | 30. para-Cymene | Phenylethyl acetate | Geranyl butyrate |
| 8. 2-Methylbutyl alcohol | 31. 1,8-Cineol | 60. Geraniol | 97. trans-Nerolidol |
| 9. Isobutyl acetate | 32. Limonene | 61. Linalyl acetate | 98. n-Amyl salicylate |
| 10. Ethyl butyrate | 33. 2,6-Dimethylhept-5-enal | 62. Geranial | 99. Phenyl ethyl tiglate |
| 11. Furfural | 34. γ-Terpinene | 63. Hydroxycitronellal | 100. Ethyl dodecanoate |
| 12. Ethyl isovalerate | 35. Octanol | 64. Citronellyl formate | 101. Benzophenone |
| 13. Hexanol | 37. Ethyl heptanoate | 66. Bornyl acetate | 102. Dibenzyl ether |
| 14. Allyl butyrate | 38. Linalool | 67. Vertenex (isomer 1) | 103. γ-Dodecalactone |
| 15. Ethyl pentanoate | 39. Benzene ethanol | 68. Ethyl nonanoate | 104. Citronellyl tiglate |
| 16. Hexylene glycol | 41. Rose oxide, cis-rose | 69. Geranyl formate | 105. Evermyl |
| 17. α-Thujone | 42. Rose oxide, trans-rose | 70. Vertenex (isomer 2) | 106. Geranyl tiglate |
| 18. Benzaldehyde | 43. Camphor | 71. γ-Nonalactone | 107. Geranyl-2-methyl valerate |
| 19. α-Pinene | 44. Citronellal | 72. Citronellyl acetate | 108. Celestolide |
| 20. Camphene | 45. Benzyl acetate | 73. Neryl acetate | 109. Heptadec-1-ene |
| 21. 3,5,5-Trimethylhexanol | 46. Menthone | 74. Geranyl acetate | 110. Benzyl benzoate |
| 22. Sabinene | 47. Isoborneol | 76. Diphenyl oxide | 111. Ethyl tetradecanoate |
| 23. β-Pinene | 48. Isomenthone | 78. Ethyl decanoate | 112. Benzyl salicylate |
| 24. Ethyl hexanoate | 49. Borneol | 79. α-Copaene | 113. Tonalid |
| 25. Myrcene | 51. Terpinen-4-ol | 80. Florazone (isomer 1) | 114. Nonadec-1-ene |
| | 52. α-Terpineol | 81. Florazone (isomer 2) | 115. Isopropylmyristate |
| | | 82. β-Caryophyllene | 116. Ethyl pentadecanoate |
| | | 83. Citronellyl propionate | Nonadecane |
| | | 85. 3,7-Guaiadiene | 117. Ethyl hexadecanoate |
| | | 88. Dodecanol | 118. Musk T (ethylene brassylate) |
| | | 89. Ethyl undecanoate | 119. Eicosane |
| | | | 120. Cinnamyl phenyl acetate |
| | | | 121. Heneicosane |
| | | | 122. Phenyl ethyl cinnamate |
| | | | 123. Ethyl octadecanoate |
| | | | 124. Herculyn D (tetrahydro & dihydro methyl abietate) |
| | | | 125. Cinnamyl cinnamate |
| | | | 126. Cetearyl octanoate |
| | | | 127. Cetearyl decanoate |

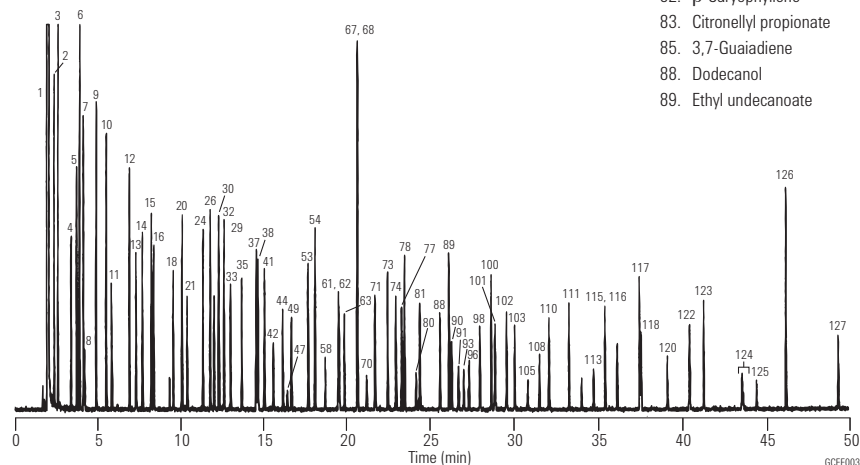

Fragrance Reference Standard

Column: DB-WAX
122-7032
30 m x 0.25 mm, 0.25 μ m

Carrier: Helium at 25 cm/s,
measured at 150 °C

Oven: 45 °C for 2 min
45-250 °C at 3 °C/min
250 °C for 34 min

Injection: Split, 250 °C
Split ratio 1:50

Detector: MSD, 250 °C transfer line

Sample: 1 μ L of a 1:20 dilution of neat sample in acetone

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759


Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 μ L tapered, FN 23-26s/42/HP, 5181-1273

Many thanks to Carl Frey, Manager of Analytical Services, Dragoco, and Kevin Myung, Director of Flavor and Perfumery Research, Bush Boake Allen, Inc. for contributing to this work.

- | | | | |
|--|----------------------------|---------------------------|------------------------------------|
| 1. Acetone | 28. Rose oxide, cis-rose | 55. Neral | 83. Ethyl tetradecanoate |
| 2. Ethyl acetate | 29. Hexanol | 56. α -Terpineol | 84. n-Amyl salicylate |
| 3. Ethyl propionate | 30. Rose oxide, trans-rose | 57. Geranyl formate | 85. Geranyl tiglate |
| 4. 2,3-Butanedione (diacetyl) | 31. Methyl-para-cresol | 58. Borneol | 86. Ethyl pentadecanoate |
| 5. Methyl butyrate | 32. Ethyl octanoate | 59. β -Bisabolene | 87. Isopropylmyristate |
| 6. Isobutyl acetate | 33. cis-Linalool oxide | 60. Benzyl acetate | 90. Phenyl ethyl tiglate |
| 7. α -Pinene | 34. Menthone | 61. Neryl acetate | 91. Rosatol (rosetone) |
| 8. Ethyl butyrate | 35. Furfural | 62. Geranial | 92. Eugenyl acetate |
| 9. 2,3-Pentanedione (acetyl propionyl) | 36. trans-Linalool oxide | 63. Ethyl undecanoate | 93. Ethyl hexadecanoate |
| 10. Camphene | 37. Octyl acetate | 64. δ -Cadinene | 94. γ -Dodecalactone |
| 11. Ethyl isovalerate | 38. Isomenthone | 65. Geranyl acetate | 95. Dibenzyl ether |
| 12. β -Pinene | 39. α -Copaene | 66. Citronellol | 96. Tonalid |
| 13. Ethyl pentanoate | 40. Camphor | 67. Ethyl dodecanoate | 97. Ethyl octadecanoate |
| 14. Myrcene | 41. Benzaldehyde | 68. Geraniol | 98. Benzophenone |
| 15. Allyl butyrate | 42. Ethyl nonanoate | 69. Benzyl alcohol | 99. Benzyl benzoate |
| 16. Limonene | 43. Linalool | 70. Geranyl butyrate | 100. Cetearyl octanoate |
| 17. 1,8-Cineol | 44. Linalyl acetate | 71. Nonadecane | 101. Musk T (ethylene brassylate) |
| 18. 3,5,5-Trimethylhexanol | 45. Vertenex (isomer 1) | 72. Benzene ethanol | 102. Cetearyl decanoate |
| 19. 3-Methylbutyl alcohol | 46. Octanol | 73. Nonadec-1-ene | 103. Frambinone (raspberry ketone) |
| 20. 2-Methylbutyl alcohol | 47. β -Caryophyllene | 74. Florazone (isomer 1) | 104. Cinnamyl phenyl acetate |
| 21. Ethyl hexanoate | 48. Vertenex (isomer 2) | 75. Florazone (isomer 2) | 105. Phenyl ethyl cinnamate |
| 22. γ -Terpinene | 49. Terpinen-4-ol | 76. Hydroxycitronellal | 106. Cinnamyl cinnamate |
| 23. p-Cymene | 50. Methyl benzoate | 77. Dodecanol | |
| 24. Hexyl acetate | 51. Hexylene glycol | 78. Diphenyl oxide | |
| 25. Terpinolene | 52. Ethyl decanoate | 79. Citronellyl tiglate | |
| 26. Ethyl heptanoate | 53. Citronellyl acetate | 80. Eugenyl methyl ether | |
| 27. 2,6-Dimethylhept-5-enal (melonal) | 54. Isoborneol | 81. γ -Nonalactone | |

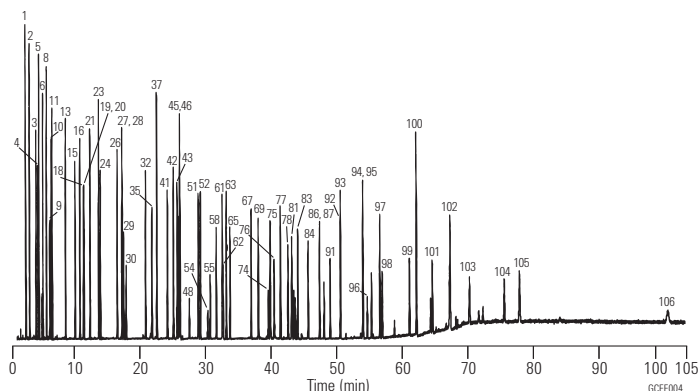

Perfume


Column: HP-INNOWax
19091N-133
30 m x 0.25 mm, 0.25 µm

Carrier: Helium, 30 cm/s
0.9 mL/min constant flow

Oven: 80 °C for 1 min
80-250 °C at 5 °C/min
250 °C for 2 min

Injection: Split, 250 °C
Split ratio 20:1

Detector: MSD, 280 °C

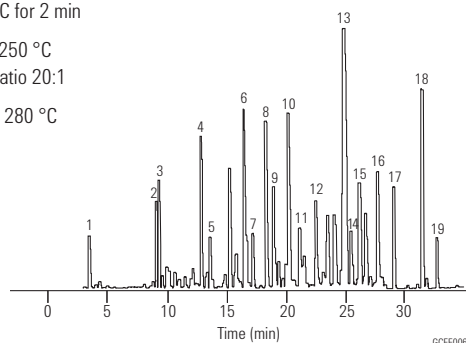

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Split, single taper, low pressure drop, glass wool, 5183-4647
Seal: Gold plated seal, 18740-20885
Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

- | | |
|---------------------------|-----------------------|
| 1. Limonene | 11. Commamyl acetate |
| 2. Linalool | 12. Acetyl cedrene |
| 3. Linalyl acetate | 13. Diethyl phthalate |
| 4. Benzyl acetate | 14. Tonalid |
| 5. Citronellol | 15. Coumarin |
| 6. Benzene ethanol | 16. Musk xylene |
| 7. α-Methyl ionone | 17. Benzyl benzoate |
| 8. Carvacrol and geraniol | 18. Benzyl salicylate |
| 9. Isoamyl salicylate | 19. Musk ketone |
| 10. n-Amyl salicylate | |

Chiral Compounds in Essential Oils and Fragrances

Column: HP-Chiral 20β
19091G-B233
30 m x 0.25 mm, 0.25 µm


Carrier: Hydrogen, 39 cm/s,
constant pressure

Oven: 65 °C for 1 min
65-170 °C at 5 °C/min

Injection: Split, 250 °C
Split ratio 30:1

Detector: FID, 300 °C

Sample: 1 µL
0.25 ng/µL each
analyte in Hexane

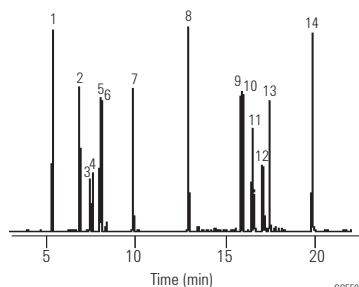

1. 1,2-Dimethylbenzene
2. Myrcene
3. (-)-Camphene
4. (+)-Camphene
5. (+)-β-Pinene
6. 1S-(-)-β-Pinene
7. Cineole
8. (R)-(-)-Citronellal
9. 1S,2R,5S-(+)-Menthhol
10. 1R,2S,5R-(-)-Menthhol
11. α-Terpineol
12. (+/-)-Isoborneol
13. (+)-Borneol
14. trans-Cinnamaldehyde

Menthhol

Column: Cyclodex-B
112-2532
30 m x 0.25 mm, 0.25 µm


Carrier: Hydrogen, 55 cm/s

Oven: 105 °C isothermal

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 1 µg/µL each chloroform

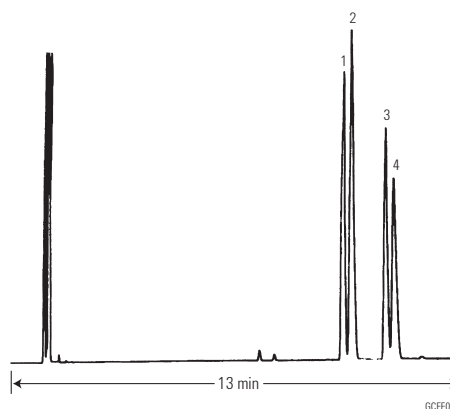

1. (+)-Neomenthol
2. (-)-Neomenthol
3. (+)-Menthhol
4. (-)-Menthhol

FAMEs

Column: DB-23
122-2362
60 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 43 cm/s,
constant pressure mode

Oven: 130 °C for 1.0 min
130-170 °C at 6.5 °C/min
170-215 °C at 2.75 °C/min
215 °C for 12 min
215-230 °C at 40 °C/min
230 °C for 3 min


Injection: Split, 270 °C
Split ratio 50:1

Detector: FID, 280 °C

Chromatogram provided courtesy of Steve Watkins and Jeremy Ching,
FAME Analytics, <http://www.fameanalytics.com>

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Split, single taper, low pressure drop,
glass wool, 5183-4647
Seal: Gold plated seal, 18740-20885
Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

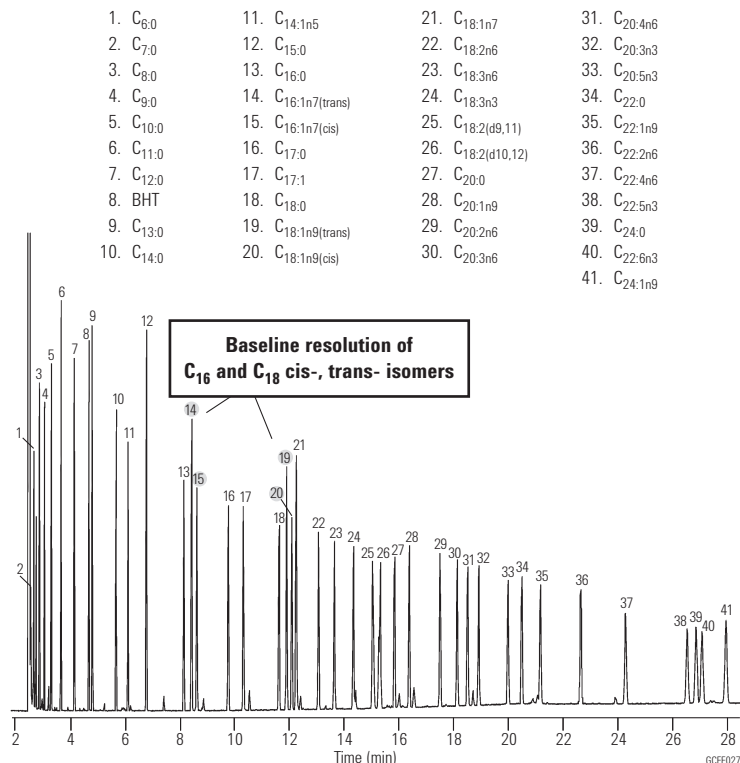**Analysis of Fragrance and Allergens**

Column: VF-WAXms
CP9205
30 m x 0.25 mm, 0.25 µm


Oven: 100 °C to 250 °C with 10 °C/min

Carrier: Helium, 1.0 mL/min

Injection: Split 1:30, T=250 °C

Detector: GC/MS Ion Trap
Trap: 200 °C
Manifold: 60 °C

Sample: 0.1 µL, Fragrances mixture (500 ppm)

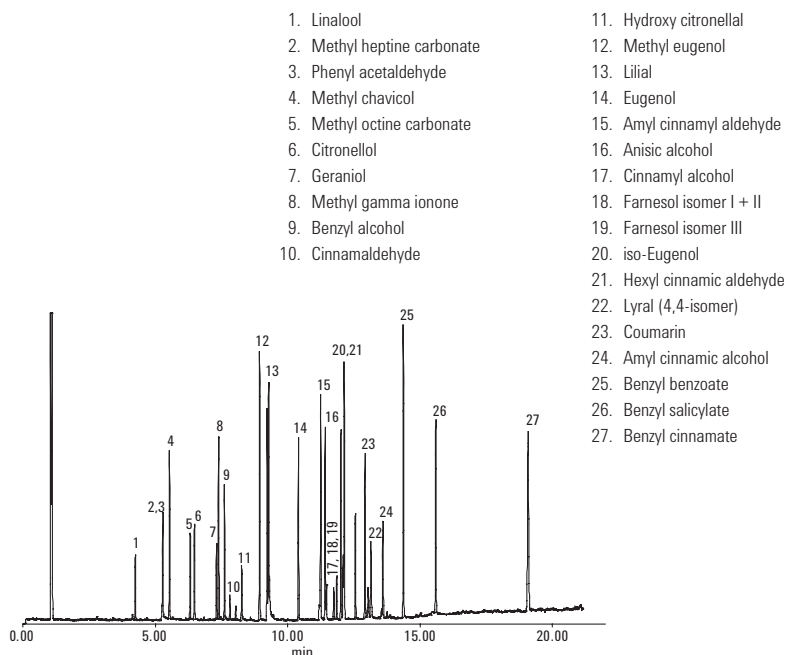

Organophosphorus Pesticide Residues in Olive Oil Extract

Column: DB-35ms Ultra Inert
122-3832UI
30 m x 0.25 mm, 0.25 µm

Instrument: Agilent 7890/5975C

Sampler: Agilent 7683B, 5.0 µL syringe (p/n 5181-1273)

CFT Device: Purged 2-way splitter (p/n G3180B)
Split ratio MSD:FPD = 1:1

MSD Restrictor: 1.43 m x 0.18 mm id deactivated fused silica tubing

FPD Restrictor: 0.53 m x 0.18 mm id deactivated fused silica tubing

Aux EPC: 3.8 psi constant pressure

Inlet: 2 µL splitless; 250 °C, purge flow 60 mL/min at 0.25 min,
gas saver on at 2 min 20 mL/min

Carrier: Helium, constant pressure 28.85 psi at 95 °C


Oven: 95 °C (0.5 min), 25 °C/min to 210 °C, 10 °C/min to 250 °C (0.5 min),
20 °C to 290 °C (4.5 min)

Postrun: 7.5 min at 290 °C, Aux EPC pressure 54 psi during backflush,

Backflush: 2 psi inlet pressure during backflush

Detector: MSD: 300 °C transfer line, 300 °C source, 150 °C quad
FPD: 230 °C, hydrogen 75 mL/min, air 100 mL/min,
carrier + makeup (N₂) 60 mL/min

- | | |
|----------------------|---|
| 1. Methamidophos | 10. Fenitrothion |
| 2. Acephate | 11. Parathion |
| 3. Omethoate | 12. Fenthion |
| 4. Diazinon | 13. Methidathion |
| 5. Dimethoate | 14. Carbophenothion |
| 6. Pirimiphos-methyl | 15. Triphenyl-phosphate (surrogate std) |
| 7. Parathion-methyl | 16. Azinphos-methyl |
| 8. Malathion | 17. Azinphos-ethyl |
| 9. Chlorpyrifos | |

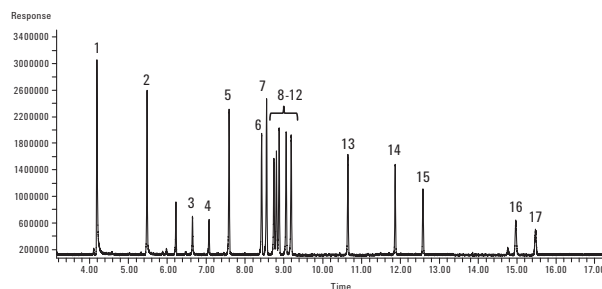

GC/FPD chromatogram of a 100 ng/mL matrix-matched organophosphorus pesticide standard with analyte protectant analyzed on an Agilent J&W DB-35ms UI GC column.

**TIPS & TOOLS**

View the latest GC column focused applications, products and educational resources at www.agilent.com/chem/myGCColumns

Fragrance Allergens

Column: HP-5ms
19091S-433
30 m x 0.25 mm, 0.25 μ m

Carrier: Helium, 1.2 mL/min,
constant pressure of 70 kPa

Oven: 50 °C in 1 min, 8 °C/min to 250 °C,
250-300 °C at 35 °C/min
300 °C hold, 5 min
5973N MSD in scan (40-350 amu)
Solvent delay, 3.0 min

Injection: Split, 250 °C
Split ratio 1:50

Sample: 1 μ L, 50 ppm standard


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 μ L tapered, FN 23-26s/42/HP, 5181-1273

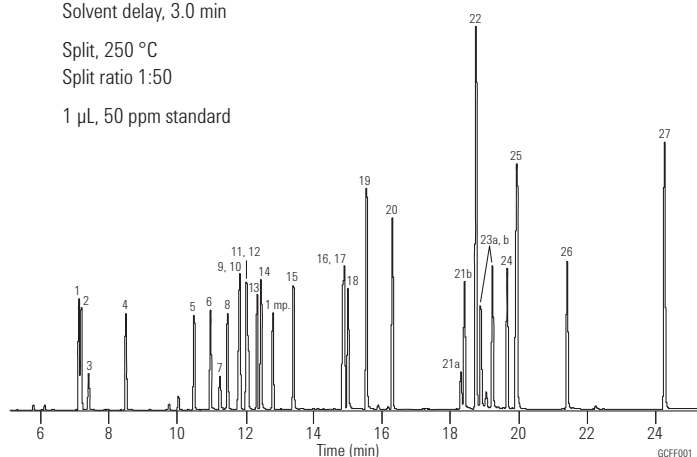

- | | |
|-----------------------------|----------------------------|
| 1. Limonene | 16. Coumarin |
| 2. Benzyl alcohol | 17. Cinnamyl acetate |
| 3. Phenyl acetaldehyde | 18. Isoeugenol |
| 4. Linalool | 19. Alpha isomethyl ionone |
| 5. Methyl heptene carbonate | 20. Lilial (BMHCA) |
| 6. Citronellol | 21a. Lyrar 1 |
| 7. Neral | 21b. Lyrar 2 |
| 8. Geraniol | 22. Amyl cinnamyl alcohol |
| 9. Citral (geranial) | 23a. Farnesol 1 |
| 10. Cinnamaldehyde | 23b. Farnesol 1 |
| 11. Anisyl alcohol | 24. Hexyl cinnamaldehyde |
| 12. Hydroxy citronellal | 25. Benzyl benzoate |
| 13. Methyl octene carbonate | 26. Benzyl salicylate |
| 14. Cinnamic alcohol | 27. Benzyl cinnamate |
| 15. Eugenol | |

Flavor Mixture

Column: Ultra 2
19091B-112
25 m x 0.32 mm, 0.52 μ m

Carrier: Helium, 90 kPa, 2.2 mL/min constant flow

Oven: 80 °C for 1 min
80-210 °C at 8 °C/min
210 °C for 2 min

Injection: Split, 250 °C
Split ratio 20:1

Detector: IRD, 280 °C
Wide Band MCT, 550 to 4000 cm^{-1}


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 μ L tapered, FN 23-26s/42/HP, 5181-1267

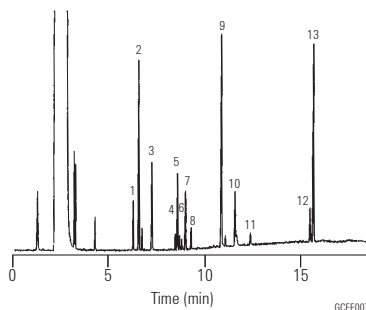

- | |
|-------------------------------|
| 1. Fenchone |
| 2. Thujone |
| 3. Benzaldehyde |
| 4. trans-Carveol |
| 5. Farnesol |
| 6. cis-Carveol |
| 7. trans-Geraniol |
| 8. Citral |
| 9. Eugenol |
| 10. Vanillin |
| 11. trans-Isoeugenol |
| 12. trans-Citronellyl tiglate |
| 13. cis-Citronellyl tiglate |

Lemon Oil

Column: DB-5
127-5022
20 m x 0.10 mm, 0.10 µm

Carrier: Hydrogen at 60 cm/s, measured at 40 °C

Oven: 40 °C for 3 min
40-185 °C at 30 °C/min
185 °C for 3 min

Injection: Split, 275 °C
Split ratio 1:275

Detector: Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

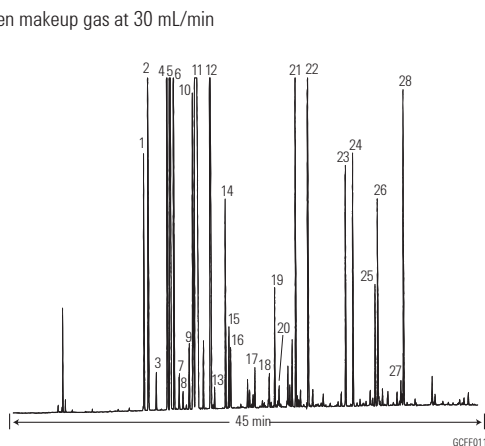

- | | |
|-------------------|-------------------------|
| 1. α-Thujone | 15. Linalool |
| 2. β-Thujone | 16. Nonanal |
| 3. Camphene | 17. Citronellal |
| 4. Sabinene | 18. Terpinen-4-ol |
| 5. β-Pinene | 19. α-Terpineol |
| 6. Myrcene | 20. Decanal |
| 7. Octanal | 21. Neral |
| 8. α-Phellandrene | 22. Geranial |
| 9. α-Terpinene | 23. Neryl acetate |
| 10. r-Cymene | 24. Geranyl acetate |
| 11. δ-Limonene | 25. β-Caryophyllene |
| 12. γ-Terpinene | 26. trans-α-Bergamotene |
| 13. Octanol | 27. α-Humulene |
| 14. Terpinolene | 28. β-Bisabolene |

Cold-pressed Orange Oil

Column: DB-5
127-5022
20 m x 0.10 mm, 0.10 µm

Carrier: Hydrogen at 60 cm/s, measured at 70 °C

Oven: 70 °C for 1 min
70-250 °C at 30 °C/min
250-310 °C at 20 °C/min
310 °C for 2 min

Injection: Split, 275 °C
Split ratio 1:275

Detector: FID, 350 °C
Nitrogen makeup gas at 30 mL/min

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

Chromatogram courtesy of Tastemaker

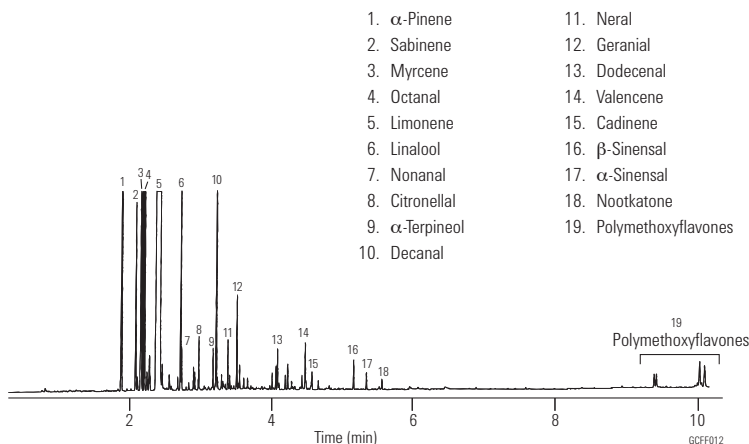

Peppermint Oil

Column: DB-WAX
122-7062
60 m x 0.25 mm, 0.25 μ m

Carrier: Helium at 25 cm/s (0.73 mL/min)

Oven: 75 °C for 8 min
75-200 °C at 4 °C/min
200 °C for 5 min

Injection: Split, 270 °C
Split ratio 1:150


Detector: FID, 270 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 μ L neat

Thanks to William Faas of A.M. Todd Company for providing the sample and assisting with peak identification.

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Split, single taper, low pressure drop, glass wool, 5183-4647
Seal: Gold plated seal, 18740-20885
Syringe: 5 μ L tapered, FN 23-26s/42/HP, 5181-1273

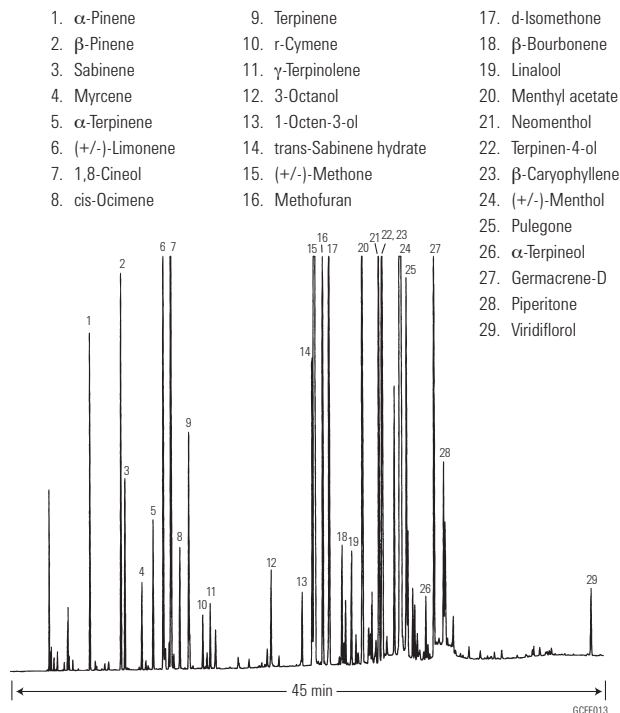**Spearmint Oil (Western)**

Column: DB-WAX
122-7062
60 m x 0.25 mm, 0.25 μ m

Carrier: Helium at 25 cm/s (0.73 mL/min)

Oven: 75 °C for 8 min
75-200 °C at 4 °C/min
200 °C for 5 min

Injection: Split, 270 °C
Split ratio 1:150


Detector: FID, 270 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 μ L neat

Thanks to William Faas of A.M. Todd Company for providing the sample and assisting with peak identification.

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Split, single taper, low pressure drop, glass wool, 5183-4647
Seal: Gold plated seal, 18740-20885
Syringe: 5 μ L tapered, FN 23-26s/42/HP, 5181-1273

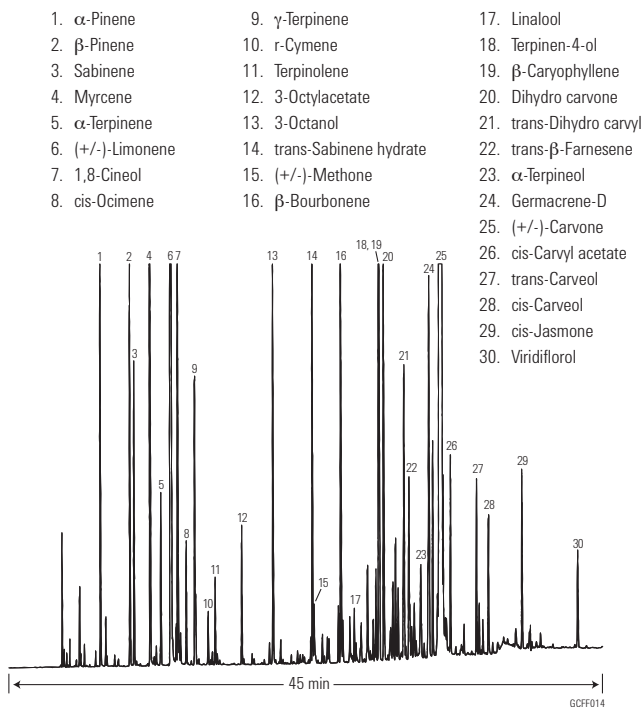

Ylang Ylang Oil

Column: DB-XLB
122-1232
30 m x 0.25 mm, 0.25 μ m

Carrier: Helium at 34 cm/s, measured at 50 °C


Oven: 50 °C for 1 min
50-250 °C at 3.5 °C/min

Injection: Split, 250 °C
Split ratio 1:125

Detector: MSD, 310 °C transfer line
full scan at m/z 35-550

Sample: 1 μ L of 10% oil in methylene chloride

- | | |
|--------------------|---------------------------|
| 1. r-Methylansiole | 7. β -Caryophyllene |
| 2. Linalool | 8. Cinnamyl acetate |
| 3. Methylbenzoate | 9. Germacrene-D |
| 4. Benzylacetate | 10. Benzyl benzoate |
| 5. Geraniol | 11. Farnesol acetate |
| 6. Geranyl acetate | 12. Benzyl salicylate |

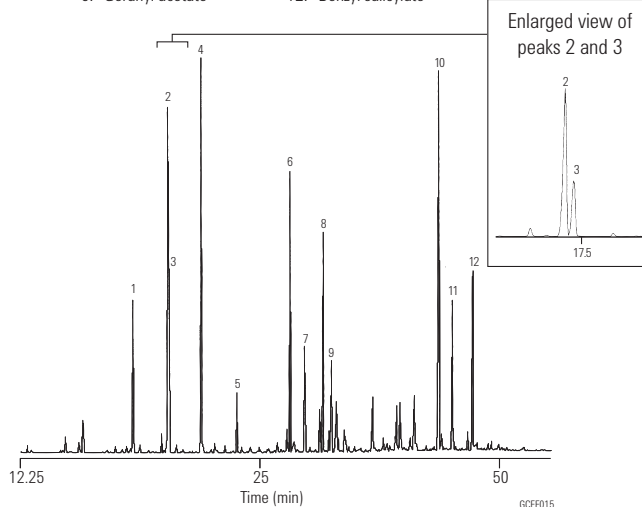

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Split, single taper, low pressure drop, glass wool, 5183-4647
Seal: Gold plated seal, 18740-20885
Syringe: 5 μ L tapered, FN 23-26s/42/HP, 5181-1273

Ylang Ylang Oil

Column: DB-WAX
121-7022
20 m x 0.18 mm, 0.18 μ m

Carrier: A: Helium 26.3 cm/s measured at 45 °C
B: Hydrogen 44.3 cm/s measured at 45 °C

Oven: A: 45 °C hold 1.28 min
4.68 °C/min to 250 °C hold 21.81 min
B: 45 °C hold 0.77 min
7.79 °C/min to 250 °C hold 13.09 min

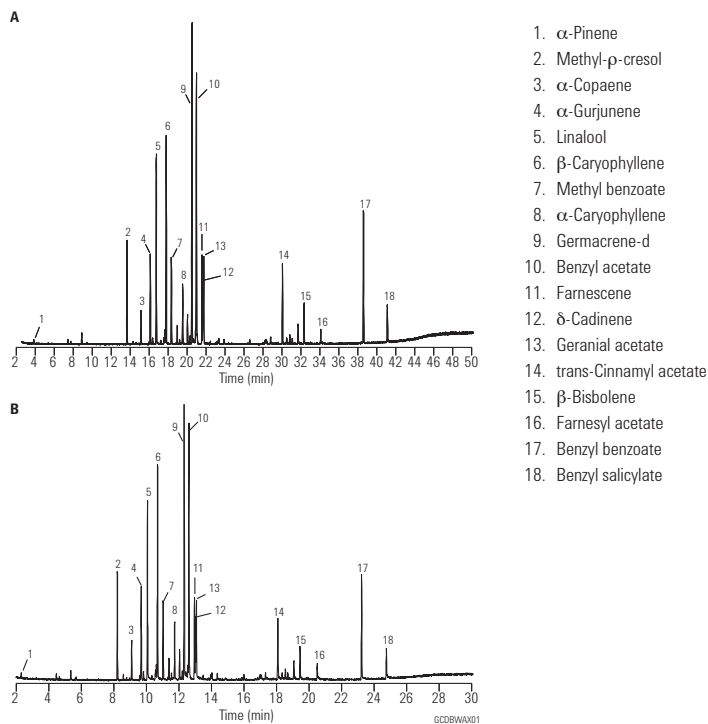

Rosemary Oil

Column: CycloSil-B
112-6632
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 40 cm/s, measured at 60 °C

Oven: 55 °C for 1 min
50-180 °C at 5 °C/min

Injection: Split, 250 °C
Split ratio 50:1

Detector: FID, 340 °C


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool,
5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

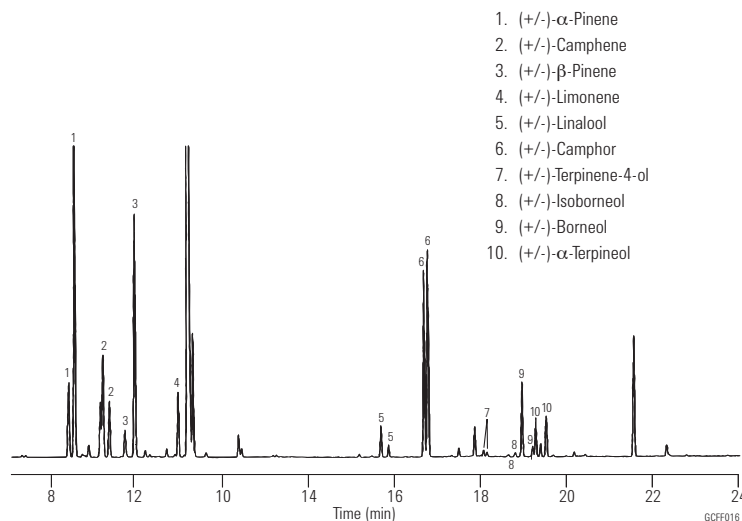

Citrus Flavored Carbonated Beverage (Soda)

Column: CycloSil-B
112-6632
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 37 cm/s,
measured at 40 °C

Oven: 40-190 °C at 2 °C/min

Sampler: Headspace
No stir, NaCl 1g/10 mL sample
Adsorption: 27 °C for 68 min
Desorption: 250 °C for 15 min

Injection: Split, 1:5
Polyacrylate fiber, 85 µm

Detector: MSD, 280 °C transfer line


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool,
5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

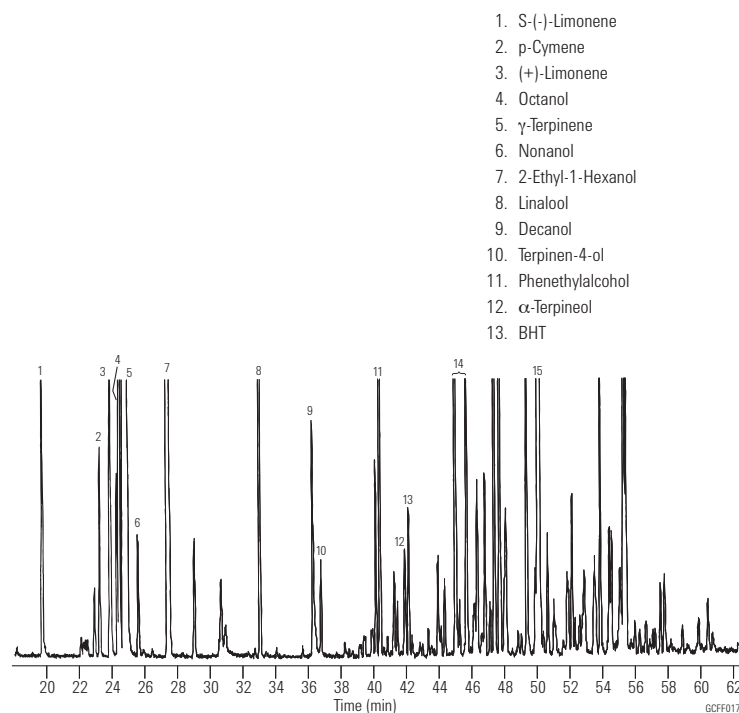


Alcohol Beverage Standard

Column: HP-FFAP
19091F-105
50 m x 0.20 mm, 0.33 μ m

Carrier: Hydrogen

Oven: 60 °C for 4 min
60-200 °C at 6 °C/min
200 °C for 2 min

Detector: FID

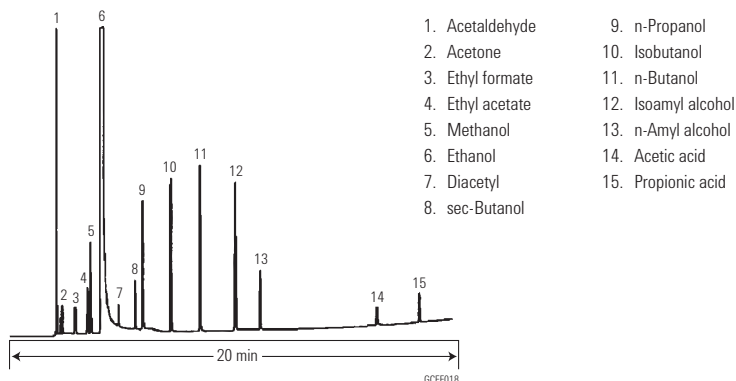

Bourbon

Column: HP-INNOWax
19091N-133
30 m x 0.25 mm, 0.25 μ m


Carrier: Helium, 33 cm/s, 15.5 psi (35 °C)
1.5 mL/min constant flow

Oven: 35 °C for 5 min
35-150 °C at 5 °C/min
150-250 °C at 20 °C/min
250 °C for 2 min

Injection: Split, 220 °C
Split ratio 25:1

Detector: FID, 280 °C

Sample: 1 μ L

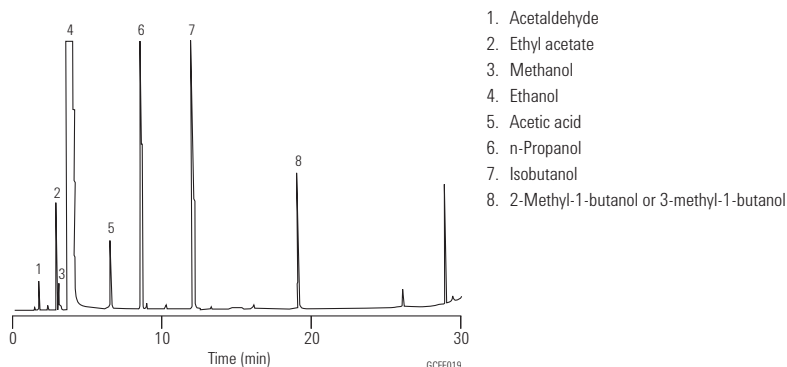

Alditol Acetates

Column: DB-225
122-2231
30 m x 0.25 mm, 0.15 μ m


Carrier: Hydrogen at 36.5 cm/s

Oven: 220 °C isothermal

Injection: Split, 225 °C
Split ratio 1:50

Detector: FID, 250 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 μ L

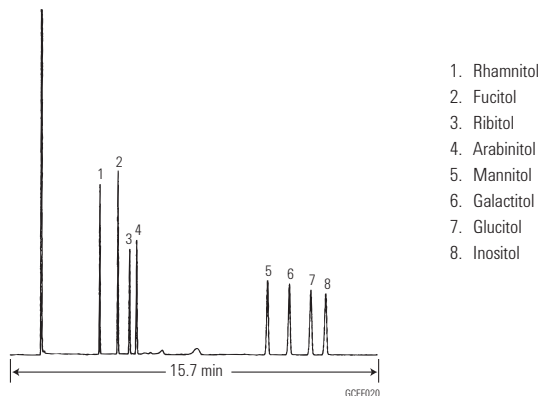

Strawberry Syrup


Column: HP-INNOWax
19091N-213
30 m x 0.32 mm, 0.50 µm

Carrier: Helium, 40 cm/s, 11.7 psi (60 °C)
2.5 mL/min constant flow

Oven: 60 °C for 1 min
60-250 °C at 10 °C/min
250 °C for 2 min

Injection: Split, 220 °C
Split ratio 60:1

Detector: FID, 275 °C

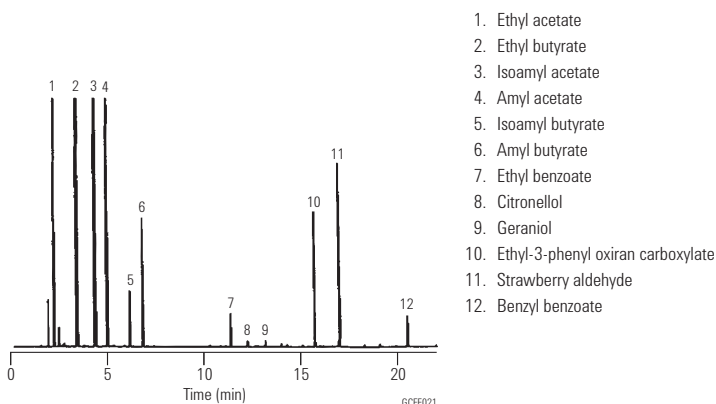**Separation of TMS-derivatized Sugars using VF-1ms**

Column: VF-1ms
CP8912
30 m x 0.25 mm, 0.25 µm

Sample: 5 µL, splitless 1 µL


Sample Conc: 40 ppb

Carrier: He, 1.0 mL/min

Oven: 105 °C to 240 °C,
4 °C/min to 300 °C,
20 °C/min

Injection: Split; 1:15

Detector: MS

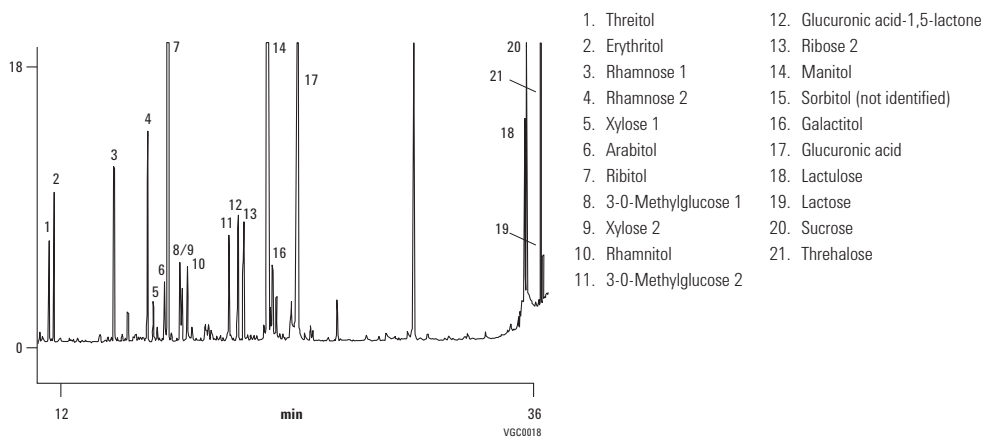

Organic Acids

Column: DB-FFAP
122-3232
30 m x 0.25 mm, 0.25 µm


Carrier: Helium at 40 cm/s, measured at 100 °C

Oven: 100 °C for 5 min
100-250 °C at 10 °C/min
250 °C for 12 min

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

- | | |
|----------------------------------|--------------------------------------|
| 1. Acetone | 10. Caproic acid (hexanoic acid) |
| 2. Formic acid | 11. Heptanoic acid |
| 3. Acetic acid | 12. Octanoic acid |
| 4. Propionic acid | 13. Decanoic acid |
| 5. Isobutyric acid | 14. Dodecanoic acid |
| 6. Butyric acid | 15. Tetradecanoic acid |
| 7. Isovaleric acid | 16. Hexadecanoic acid |
| 8. Valeric acid (pentanoic acid) | 17. Octadecanoic acid |
| 9. Isocaproic acid | 18. Arachidic acid (eicosanoic acid) |

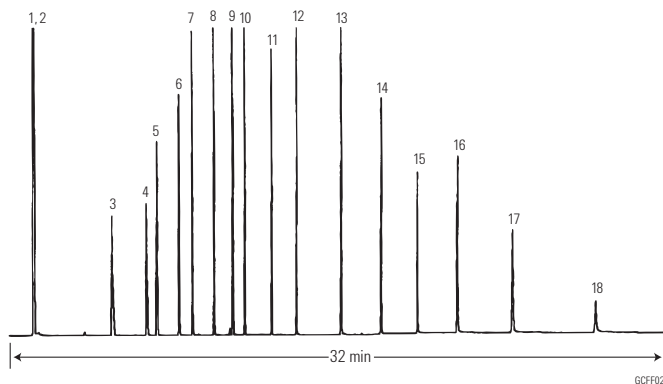

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

Acids

Column: VF-WAXms
CP9205
30 m x 0.25 mm, 0.25 µm


Sample: Acid sample, 0.1% (Cyclohexane), 1.0 µL

Carrier: Hydrogen, 75 kPa

Oven: 60 °C to 200 °C, 5 °C/min

Injection: 250 °C, split 100 mL/min

VF-WAXms

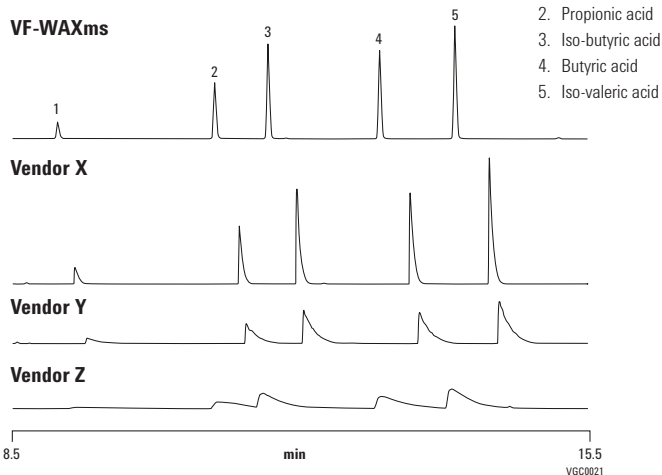

- | |
|---------------------|
| 1. Acetic acid |
| 2. Propionic acid |
| 3. Iso-butyric acid |
| 4. Butyric acid |
| 5. Iso-valeric acid |

Vendor X

Vendor Y

Vendor Z

Bacterial Fatty Acid Methyl Esters

Column: DB-5
122-5032
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 42 cm/s

Oven: 150 °C for 4 min
150-250 °C at 4 °C/min

Injection: Split ratio 1:100

Detector: FID
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

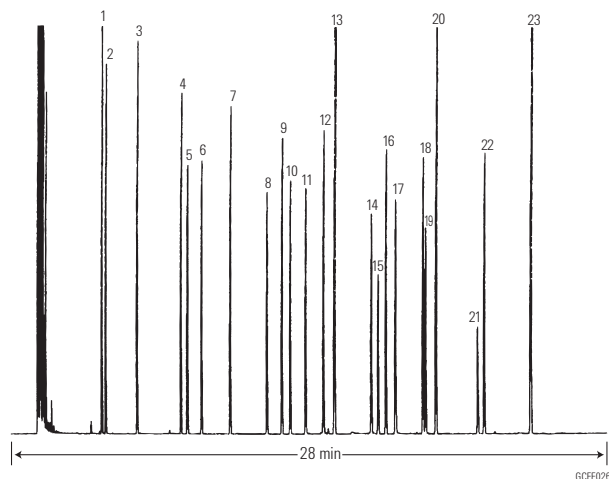

- | | |
|---------------------------------|---|
| 1. C _{11:0} | Methyl undecanoate |
| 2. 2-OH C _{10:0} | Methyl 2-hydroxydecanoate |
| 3. C _{12:0} | Methyl laurate |
| 4. C _{13:0} | Methyl tridecanoate |
| 5. 2-OH C _{12:0} | Methyl 2-hydroxydodecanoate |
| 6. 3-OH C _{12:0} | Methyl 3-hydroxydodecanoate |
| 7. C _{14:0} | Methyl myristate |
| 8. 12-Me C _{14:0} | Methyl 12-methyltetradecanoate |
| 9. C _{15:0} | Methyl pentadecanoate |
| 10. 2-OH C _{14:0} | Methyl 2-hydroxytetradecanoate |
| 11. 3-OH C _{14:0} | Methyl 3-hydroxytetradecanoate |
| 12. C _{16:1} | Methyl palmitoleate |
| 13. C _{16:0} | Methyl palmitate |
| 14. 14-Me C _{16:0} | Methyl 14-methylhexadecanoate |
| 15. 9,10-diMe C _{16:0} | Methyl cis-9,10-methyl hexadecanoate |
| 16. C _{17:0} | Methyl heptadecanoate |
| 17. 2-OH C _{16:0} | Methyl 2-hydroxyhexadecanoate |
| 18. C _{18:1} | Methyl oleate |
| 19. C _{18:1} | Methyl elaidate |
| 20. C _{18:0} | Methyl stearate |
| 21. 9,10-diMe C _{18:0} | Methyl cis-9,10-methylene octadecanoate |
| 22. C _{19:0} | Methyl nonadecanoate |
| 23. C _{20:0} | Methyl arachidate |

Separation of cis-trans FAME Isomers

Column: Select FAME
CP7421
200 m x 0.25 mm

Sample: 0.5 µL


Sample Conc: 5 ng approx. per component on the column

Carrier: Helium, 520 kPa

Oven: 185 °C

Injection: Split, 1:20

Detector: FID

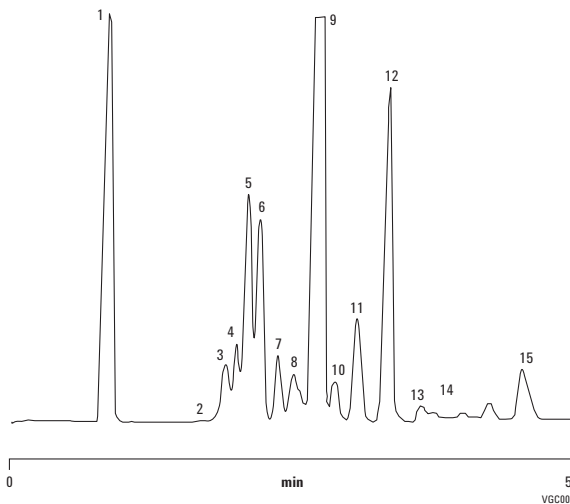

- | |
|-----------------------------------|
| 1. C _{18:0} |
| 2. C _{18:1} 7 trans |
| 3. C _{18:1} 8 trans |
| 4. C _{18:1} 9 trans |
| 5. C _{18:1} 10 trans |
| 6. C _{18:1} 11 trans |
| 7. C _{18:1} 12 trans |
| 8. C _{18:1} 13 trans + ? |
| 9. C _{18:1} 9 cis |
| 10. C _{18:1} 10 cis |
| 11. C _{18:1} 11 cis |
| 12. C _{18:1} 12 cis |
| 13. C _{18:1} 13 cis |
| 14. C _{18:1} 14 cis |
| 15. C _{18:1} 15 cis |

69 Component FAME Mix

Column: HP-88
112-8867
60 m x 0.25 mm, 0.20 μ m


Carrier: He at 1.4 mL/min constant flow

Oven: 125 °C
125 °C to 145 °C at 8 °C/min
145 °C for 26 min
145 °C to 220 °C at 2 °C/min
220 °C for 1 min

Injection: Split, 250 °C
Split ratio 50:1
1 μ L of 70 ppm each in CHCl₃

Detector: FID, 260 °C

- | | | | |
|-----------------|---------------------|------------------------|----------------------------|
| 1. nC6:0 | 16. C15:1 (14c) | 31. C19:1 (10t) | 50. C20:3 (8c,11c,14c) |
| 2. nC7:0 | 17. nC16:0 | 32. nC19:0 | 51. nC22:0 |
| 3. nC8:0 | 18. C16:1 (9t) | 33. C19:1 (7t) | 52. C22:1 (13t) |
| 4. nC9:0 | 19. C16:1 (9c) | 34. C18:2 (9c,12c) | 53. C20:4 (5c,8c,11c,14c) |
| 5. nC10:0 | 20. nC17:0 | 35. C19:1 (7c) | 54. C20:3 (11c,14c,17c) |
| 6. nC11:0 | 21. C17:1 (10t) | 36. C19:1 (10c) | 55. C21:2 (12c,15c) |
| 7. nC12:0 | 22. C17:1 (10c) | 37. C18:3 g(6c,9c,12c) | 56. C22:1 (13c) |
| 8. C12:1 (11c) | 23. nC18:0 | 38. nC20:0 | 57. nC23:0 |
| 9. nC13:0 | 24. C18:1 (6t) | 39. C18:3 (9c,12c,15c) | 58. C20:5 (EPA) |
| 10. nC14:0 | 25. C18:1 (9t) | 40. C20:1 (5c) | 59. C22:2 (13c,16c) |
| 11. C14:1 (9t) | 26. C18:1 (11t) | 41. C19:2 (10c,13c) | 60. C23:1 (14c) |
| 12. C14:1 (9c) | 27. nC18:1 (6c) | 42. C20:1 (11t) | 61. nC24:0 |
| 13. nC15:0 | 28. C18:1 (9c) | 43. C18:2 CONJ | 62. C22:3 (13c,16c,19c) |
| 14. C15:1 (10t) | 29. C18:1 (11c) | 44. C20:1 (8c) | 63. C22:4 (7c,10c,13c,16c) |
| 15. C15:1 (10c) | 30. nC18:2 (9t,12t) | 45. C20:1 (11c) | 64. C24:1 (15c) |
| | | 46. C18:2 (10t,12c) | 65. C22:5 (DPA) |
| | | 47. nC21:0 | 66. C22:6 (DHA) |
| | | 48. C20:2 (11c,14c) | 67. C18:1-12 Hydroxy (9t) |
| | | 49. C21:1 (12c) | 68. C18:0 12 Hydroxy |
| | | | 69. C18:1-12 Hydroxy (9c) |

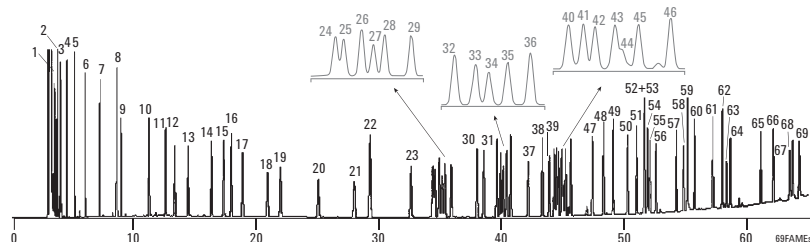

FAME Standard

Column: DB-WAX
127-7012
10 m x 0.10 mm, 0.10 μ m

Carrier: Hydrogen at 77 cm/s,
measured at 40 °C

Oven: 40 °C for 0.5 min
40-195 °C at 25 °C/min
195-205 °C at 3 °C/min
205-230 °C at 8 °C/min
230 °C for 1 min

Injection: Split, 250 °C
Split ratio 1:30

Detector: FID, 250 °C

- | | |
|---|--|
| 1. Butyric acid methyl ester (C _{4:0}) | 20. Linolelaidic acid methyl ester (C _{18:2n6t}) |
| 2. Caproic acid methyl ester (C _{6:0}) | 21. γ -Linolenic acid methyl ester (C _{18:3n6}) |
| 3. Caprylic acid methyl ester (C _{8:0}) | 22. Linolenic acid methyl ester (C _{18:3n3}) |
| 4. Capric acid methyl ester (C _{10:0}) | 23. Arachidic acid methyl ester (C _{20:0}) |
| 5. Undecanoic acid methyl ester (C _{11:0}) | 24. cis-11,14-Eicosenoic acid methyl ester (C _{20:1}) |
| 6. Lauric acid methyl ester (C _{12:0}) | 25. cis-11,14-Eicosadienoic acid methyl ester (C _{20:2}) |
| 7. Tridecanoic acid methyl ester (C _{13:0}) | 26. cis-8,11,14-Eicosatrienoic acid methyl ester (C _{20:3n6}) |
| 8. Myristic acid methyl ester (C _{14:0}) | 27. Heneicosanoic acid methyl ester (C _{21:0}) |
| 9. Myristoleic acid methyl ester (C _{14:1}) | 28. cis-11,14,17-Eicosatrienoic acid methyl ester (C _{20:3n3}) |
| 10. Pentadecanoic acid methyl ester (C _{15:0}) | 29. Arachidonic acid methyl ester (C _{20:4n6}) |
| 11. cis-10-Pentadecenoic acid methyl ester (C _{15:1}) | 30. cis-5,8,11,14,17-Eicosapentaenoic acid methyl ester (C _{20:5n3}) |
| 12. Palmitic acid methyl ester (C _{16:0}) | 31. Behenic acid methyl ester (C _{22:0}) |
| 13. Palmitoleic acid methyl ester (C _{16:1}) | 32. Erucic acid methyl ester (C _{22:1n9}) |
| 14. Heptadecanoic acid methyl ester (C _{17:0}) | 33. cis-13,16-Docosadienoic acid methyl ester (C _{22:2}) |
| 15. cis-10-Heptadecenoic acid methyl ester (C _{17:1}) | 34. Tricosanoic acid methyl ester (C _{23:0}) |
| 16. Stearic acid methyl ester (C _{18:0}) | 35. Lignoceric acid methyl ester (C _{24:0}) |
| 17. Oleic acid methyl ester (C _{18:1n9c}) | 36. cis-4,7,10,13,16,19-Docosahexaenoic acid methyl ester (C _{22:6n3}) |
| 18. Elaidic acid methyl ester (C _{18:1n9t}) | 37. Nervonic acid methyl ester (C _{24:1}) |
| 19. Linoleic acid methyl ester (C _{18:2n6c}) | |


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop,
glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 μ L tapered, FN 23-26s/42/HP,
5181-1273

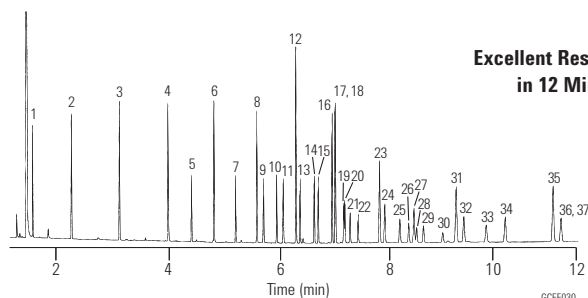

FAME Standard

Column: DB-225
127-2222
20 m x 0.10 mm, 0.10 µm

Carrier: Hydrogen at 59.3 cm/s,
measured at 35 °C

Oven: 35 °C for 0.5 min
35-195 °C at 25 °C/min
195-205 °C at 3 °C/min
205-230 °C at 8 °C/min
230 °C for 1 min

Injection: Split, 250 °C
Split ratio 1:30

Detector: FID, 250 °C

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

- | | |
|--|---|
| 1. Butyric acid methyl ester (C4:0) | 20. Linolelaidic acid methyl ester (C18:2n6t) |
| 2. Caproic acid methyl ester (C6:0) | 21. γ-Linolenic acid methyl ester (C18:3n6) |
| 3. Caprylic acid methyl ester (C8:0) | 22. Linolenic acid methyl ester (C18:3n3) |
| 4. Capric acid methyl ester (C10:0) | 23. Arachidic acid methyl ester (C20:0) |
| 5. Undecanoic acid methyl ester (C11:0) | 24. cis-11-Eicosenoic acid methyl ester (C20:1) |
| 6. Lauric acid methyl ester (C12:0) | 25. cis-11,14-Eicosadienoic acid methyl ester (C20:2) |
| 7. Tridecanoic acid methyl ester (C13:0) | 26. cis-8,11,14-Eicosatrienoic acid methyl ester (C20:3n6) |
| 8. Myristic acid methyl ester (C14:0) | 27. Heneicosanoic acid methyl ester (C21:0) |
| 9. Myristoleic acid methyl ester (C14:1) | 28. cis-11,14,17-Eicosatrienoic acid methyl ester (C20:3n3) |
| 10. Pentadecanoic acid methyl ester (C15:0) | 29. Arachidonic acid methyl ester (C20:4n6) |
| 11. cis-10-Pentadecenoic acid methyl ester (C15:1) | 30. cis-5,8,11,14,17-Eicosapentaenoic acid methyl ester (C20:5n3) |
| 12. Palmitic acid methyl ester (C16:0) | 31. Behenic acid methyl ester (C22:0) |
| 13. Palmitoleic acid methyl ester (C16:1) | 32. Erucic acid methyl ester (C22:1n9) |
| 14. Heptadecanoic acid methyl ester (C17:0) | 33. cis-13,16-Docosadienoic acid methyl ester (C22:2) |
| 15. cis-10-Heptadecenoic acid methyl ester (C17:1) | 34. Tricosanoic acid methyl ester (C23:0) |
| 16. Stearic acid methyl ester (C18:0) | 35. Lignoceric acid methyl ester (C24:0) |
| 17. Oleic acid methyl ester (C18:1n9c) | 36. cis-4,7,10,13,16,19-Docosahexaenoic acid methyl ester (C22:6n3) |
| 18. Elaidic acid methyl ester (C18:1n9t) | 37. Nervonic acid methyl ester (C24:1) |
| 19. Linoleic acid methyl ester (C18:2n6c) | |

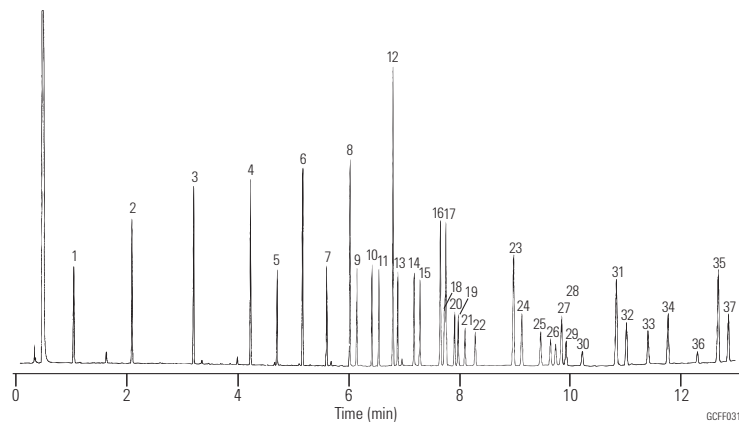

Canola Oil Margarine Partially Hydrogenated FAMES AOCs Method 1c-89

Column: DB-23
122-2362
60 m x 0.25 mm, 0.25 µm


Carrier: Helium at 15 cm/s (0.44 mL/min),
measured at 150 °C

Oven: 150-200 °C at 1.3 °C/min
200 °C for 10 min

Injection: Split, 210 °C
Split 1:100

Detector: FID, 210 °C

Sample: 1 µL

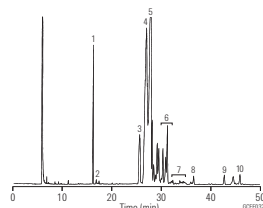

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop,
glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP,
5181-1273

1. C16:0 Methyl palmitate
2. C16:1 Methyl palmitoleate
3. C18:0 Methyl stearate
4. C18:1 trans-Methyl elaidate and multiple isomers
5. C18:1 cis-Methyl oleate and multiple isomers
6. C18:2 trans-Multiple isomers
7. C18:2 cis-Multiple isomers
8. C18:3 Methyl linolenate
9. C20:0 Methyl arachidate
10. C20:1 Methyl 11-eicosanoate

Butter Triglycerides I

Column: DB-5ht
123-5731
30 m x 0.32 mm, 0.10 µm


Carrier: Hydrogen at 55 cm/s, measured at 250 °C

Oven: 35-250 °C at 70 °C/min
250-400 °C at 5 °C/min
400 °C for 20 min

Injection: Cool on-column

Detector: FID, 400 °C
Nitrogen makeup gas at 30 mL/min
Baseline corrected

Sample: 1 µL of 9 µg/µL in toluene
(approximately 1% w/w solution)

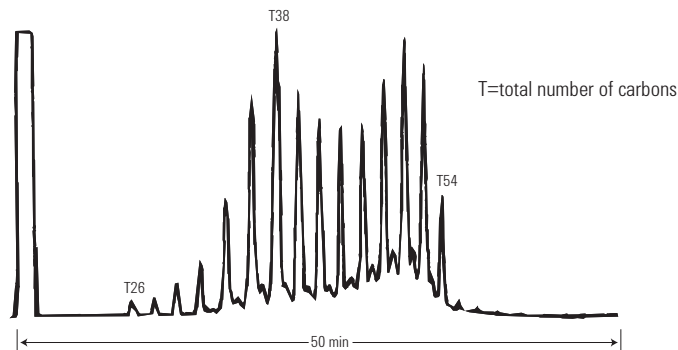

Butter Triglycerides II

Column: DB-17ht
123-1831
30 m x 0.32 mm, 0.15 µm


Carrier: Hydrogen at 40 cm/s

Oven: 250-365 °C at 5 °C/min
365 °C for 1 min

Injection: Cool on-column

Detector: FID, 400 °C
Nitrogen makeup gas at 30 mL/min
Baseline corrected

Sample: 1 µL of 9 µg/µL in toluene
(approximately 1% w/w solution)

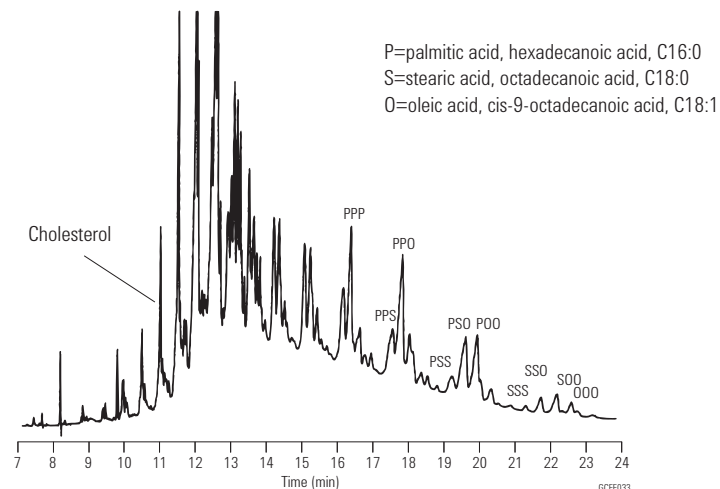

Fast Screening of FAME Isomers in Butter

Column: VF-23ms
CP8822
30 m x 0.25 mm, 0.25 µm


Sample: 0.5 µL ca. 5 ng per component on column

Carrier: Hydrogen, 70 kPa

Oven: 185 °C

Injection: Split, 1:100
T=275 °C

Detector: FID

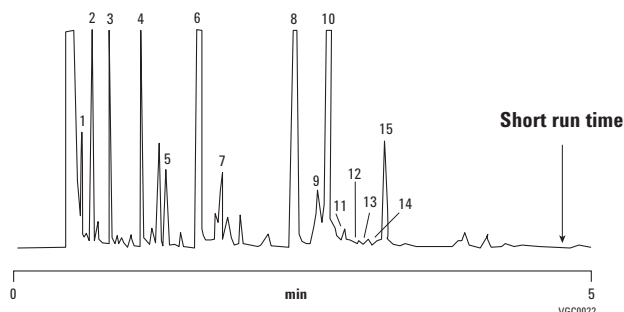

1. C8:0
2. C10:0
3. C12:0
4. C14:0
5. C14:1
6. C14:1
7. C16:1 9-cis
8. C16:1 9-cis
9. C18:1 trans
10. C18:1 9-cis
11. C18:1 13-cis
12. C18:2 9-trans, 12-trans
13. C18:2 9-cis, 12-trans
14. C18:2 9-trans, 12-cis
15. C18:2 9-cis, 12-cis

Pesticides in Sunflower Oil

Column: VF-5ms
CP8960
60 m x 0.25 mm, 0.25 µm

Sample: 5 µL, splitless

Sample Conc: 40 ppb

Carrier: He, 1.2 mL/min, constant flow



Oven: 70 °C (3.0 min), 25 °C to 190 °C/min (0.0 min) to
10 °C/min to 320 °C (10 min)

Injection: 1079 with carbofrit liner

Detector: A: Ion Trap in MS/MS, full scan
B: MS/MS

- | | |
|----------------------|------------------------|
| 1. β-HCH | 10. Bromofos |
| 2. γ-HCH | 11. o,p'-DDE |
| 3. δ-HCH | 12. α-Endosulfan |
| 4. + Vinclozolin | 13. p,p'-DDE |
| 5. Pyrimiphos methyl | 14. o,p'-DDD |
| 6. + Malathion | 15. Dieldrin |
| 7. Chloropyrifos | 16. p,p'-DDD |
| 8. Ethyl parathion | 17. b Endosulfan |
| 9. Pyrimiphos ethyl | 18. p,p'-DDT |
| | 19. Endosulfan sulfate |

- | | |
|----------------------|------------------------|
| 1. β-HCH | 10. Promofos |
| 2. γ-HCH | 11. o,p'-DDE |
| 3. δ-HCH | 12. α-Endosulfan |
| 4. + Vinclozolin | 13. p,p'-DDE |
| 5. Methyl parathion | 14. o,p'-DDD |
| 6. Pyrimiphos methyl | 15. Dieldrin |
| 7. + Fenitrothion | 16. p,p'-DDD |
| 8. Chloropyrifos | 17. b Endosulfan |
| 9. Pyrimiphos ethyl | 18. p,p'-DDT |
| | 19. Endosulfan sulfate |

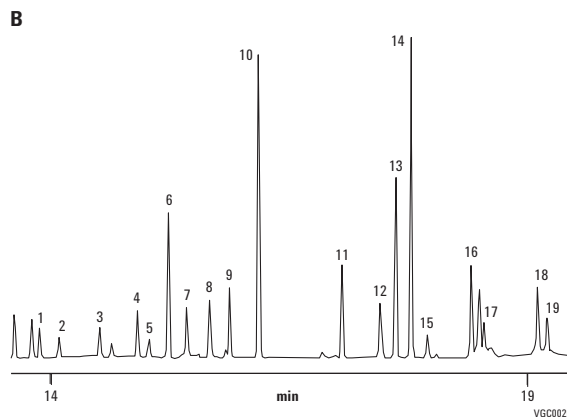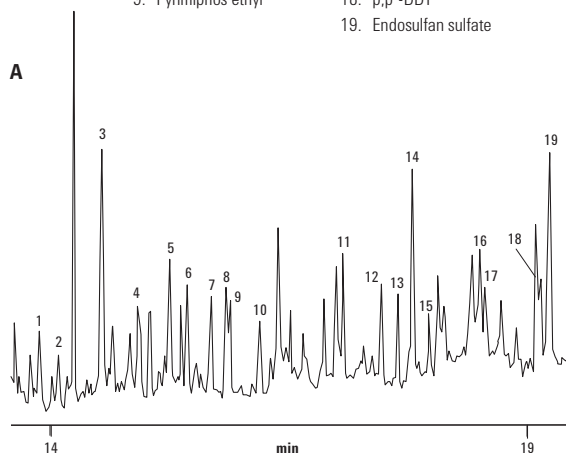

Energy and Fuels Applications

Fast Analysis of Aromatic Solvent


Column: HP-INNOWax
19091N-216
60 m x 0.32 mm, 0.50 μ m

Carrier: Helium at 20 psi constant pressure mode
Oven: 75 °C (10 min); 3 °C/min to 100 °C (0 min)
10 °C/min to 145 °C (0 min)
Injection: Split/splitless at 250 °C
100:1 split ratio
Detector: FID at 250 °C
Sample: 1.0 μ L


Column: HP-INNOWax
19091N-577
20 m x 0.18 mm, 0.18 μ m

Carrier: Helium at 33 psi constant pressure mode
Oven: 70 °C (3 min); 45 °C/min to 145 °C (1 min)
Injection: Split/splitless at 250 °C
100:1 to 600:1 split ratio
Detector: FID at 250 °C
Sample: 0.2 to 1.0 μ L

Unified aromatic solvent method

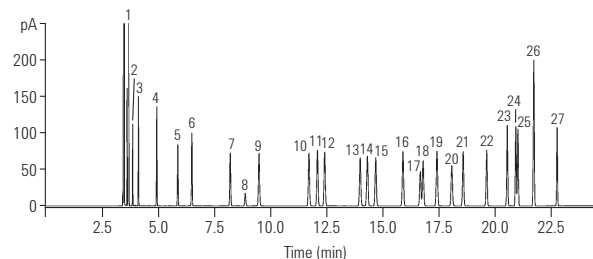

Optimized unified aromatic solvent method

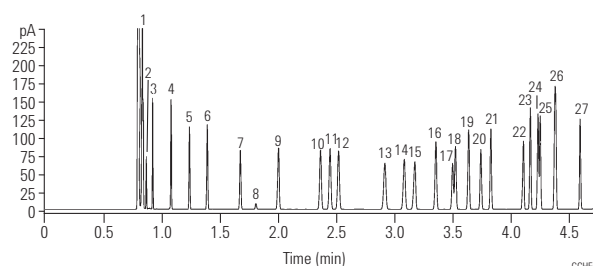

1. Heptane
2. Cyclohexane
3. Octane
4. Nonane
5. Benzene
6. Decane
7. Toluene
8. 1,4-Dioxane
9. Undecane
10. Ethylbenzene
11. p-Xylene
12. m-Xylene
13. Cumene
14. Dodecane
15. o-Xylene
16. Propylbenzene
17. p-Ethyltoluene
18. m-Ethyltoluene
19. t-Butylbenzene
20. s-Butylbenzene
21. Styrene
22. Tridecane
23. 1,3-Diethylbenzene
24. 1,2-Diethylbenzene
25. n-Butylbenzene
26. a-Methylstyrene
27. Phenylacetylene

This application showcases the practicality using high efficiency GC columns in daily aromatic solvent analysis. The result: a four-fold reduction in run time (compared to a 0.32 mm id column) with no compromise in resolution.

Refinery Gas I

Column: HP-PLLOT Q
19095P-Q04
30 m x 0.53 mm, 40.00 µm

Carrier: Helium p=9.0 psi at 60 °C

Oven: 60 °C for 5 min
60-200 °C at 20 °C/min
200 °C for 1 min

Injection: Split, 250 °C
Split flow 100 mL/min
0.25 cc valve

Detector: TCD, 250 °C


Sample: Refinery gas and others

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

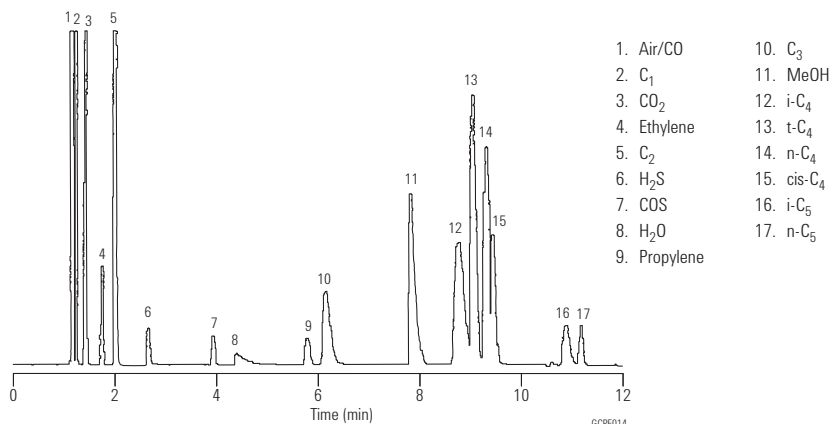**Unleaded Gasoline**

Column: DB-Petro
122-10A6
100 m x 0.25 mm, 0.50 µm

Carrier: Helium at 25.6 cm/s

Oven: 0 °C for 15 min
0-50 °C at 1 °C/min
50-130 °C at 2 °C/min
130-180 °C at 4 °C/min
180 °C for 20 min

Injection: Split, 200 °C
Split ratio 1:300


Detector: FID, 250 °C
Nitrogen makeup gas
at 30 mL/min

Sample: 1 µL of neat sample

- | | | |
|-----------------------|----------------------------|--------------------------------|
| 1. Methane | 12. 2,3,3-Trimethylpentane | 23. 1,2,4-Trimethylbenzene |
| 2. n-Butane | 13. 2-Methylheptane | 24. Isobutylbenzene |
| 3. Isopentane | 14. 4-Methylheptane | 25. sec-Butylbenzene |
| 4. n-Pentane | 15. n-Octane | 26. n-Decane |
| 5. n-Hexane | 16. Ethylbenzene | 27. 1,2,3-Trimethylbenzene |
| 6. Methylcyclopentane | 17. m-Xylene ** | 28. Butylbenzene |
| 7. Benzene | 18. p-Xylene | 29. n-Undecane |
| 8. Cyclohexane | 19. o-Xylene | 30. 1,2,4,5-Tetramethylbenzene |
| 9. Isooctane | 20. n-Nonane | 31. Naphthalene |
| 10. n-Heptane | 21. Isopropylbenzene | 32. Dodecane |
| 11. Toluene * | 22. Propylbenzene | 33. Tridecane |

*Valley point with 12 = 78%

**Valley point with 18 = 87%

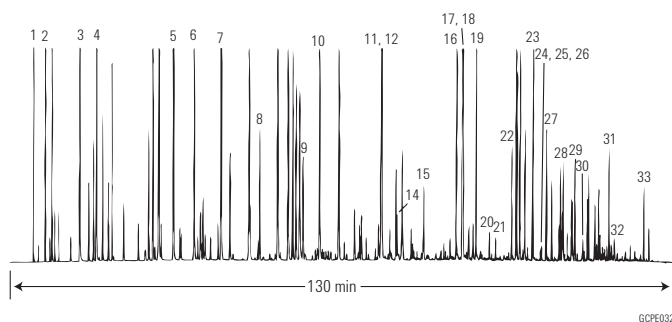

n-Paraffin Standard

Column: DB-HT Sim Dis
145-1001
5 m x 0.53 mm, 0.15 µm


Carrier: Helium at 18 mL/min, measured at 35 °C

Oven: -30-430 °C at 10 °C/min

Injection: OPTIC PTV
55-450 °C at 2 °C/s

Detector: FID, 450 °C
Nitrogen makeup gas at 15 mL/min

Sample: 0.5 µL of about 2% n-paraffins in CS₂

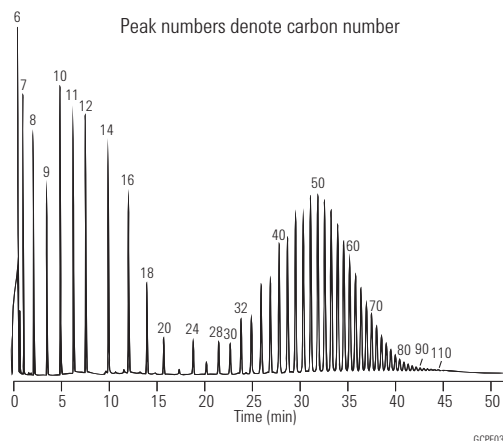**Sulfur Standards in Toluene**

Column: DB-Sulfur SCD
G3903-63001
60 m x 0.32 mm, 4.20 µm

Inlet: 275 °C, Split ratio 10:1
(Inert Flow Path split/splitless inlet)

Carrier: Helium, constant flow mode, 2.8 mL/min

Oven: 35 °C for 3 min,
35 °C to 250 °C at 10 °C/min,
250 °C for 10 min

Injection: 1 µL

Burner temperature: 800 °C

Vacuum of burner: 364 torr

Vacuum of reaction cell: 5 torr

Hydrogen: 40 mL/min

Air: 60 mL/min


Suggested Supplies

Septum: Non-stick bleed and temperature optimized (BTO) septa, 11 mm, 50/pk, 5183-4757

Liner: Low pressure drop, Ultra Inert Liner with glass wool, 5190-2295

Seal: Ultra Inert gold plated seal and washer, 5190-6144

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

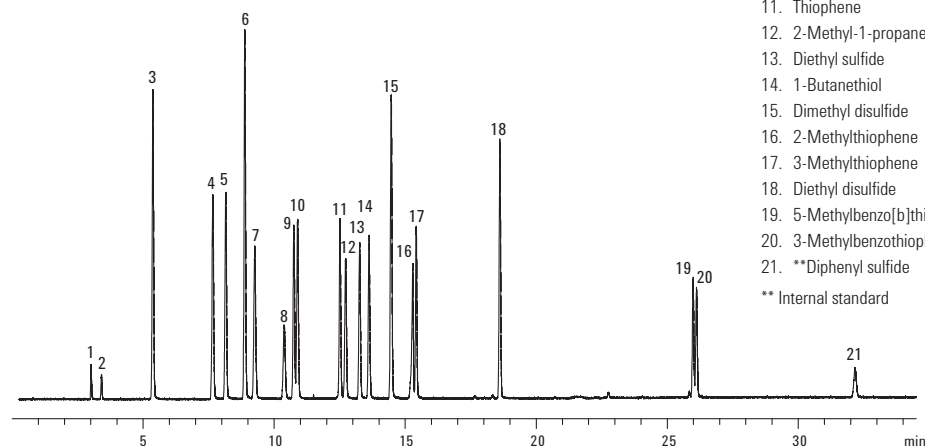

CAS No.	Formula	Concentration (mg/kg)
1. Hydrogen sulfide	H ₂ S	2000
2. Carbonyl sulfide	COS	2000
3. Methanethiol	CH ₃ SH	2000
4. Ethanethiol	C ₂ H ₅ SH	2000
5. Dimethyl sulfide	(CH ₃) ₂ S	2000
6. Carbon disulfide	CS ₂	2000
7. 2-Propanethiol	C ₃ H ₇ S	2000
8. 2-Methyl-2-propanethiol	C ₄ H ₁₀ S	2000
9. 1-Propanethiol	C ₃ H ₇ S	2000
10. Ethyl methyl sulfide	C ₂ H ₅ SCH ₃	2000
11. Thiophene	C ₄ H ₄ S	2000
12. 2-Methyl-1-propanethiol	C ₄ H ₁₀ S	2000
13. Diethyl sulfide	(C ₂ H ₅) ₂ S	2000
14. 1-Butanethiol	C ₄ H ₁₀ S	2000
15. Dimethyl disulfide	(CH ₃) ₂ S ₂	2000
16. 2-Methylthiophene	C ₅ H ₆ S	2000
17. 3-Methylthiophene	C ₅ H ₆ S	2000
18. Diethyl disulfide	(C ₂ H ₅) ₂ S ₂	2000
19. 5-Methylbenzo[b]thiophene	C ₉ H ₈ S	2000
20. 3-Methylbenzothiophene	C ₉ H ₈ S	2000
21. **Diphenyl sulfide	C ₁₂ H ₁₀ S	2000

** Internal standard

Sulfur Compounds in Propylene (1 ppm)

Column: GS-GasPro
113-4332
30 m x 0.32 mm


Oven: 60 °C for 2.5 min
60-250 °C at 10 °C/min

Injection: OI Analytical Volatiles Inlet
Split ratio 5:1
200 µL gas sampling valve

Detector: OI Analytical Model 5380 PFPD

Sample: 1 ppm sulfur compounds in propylene

Chromatogram courtesy of OI Analytical

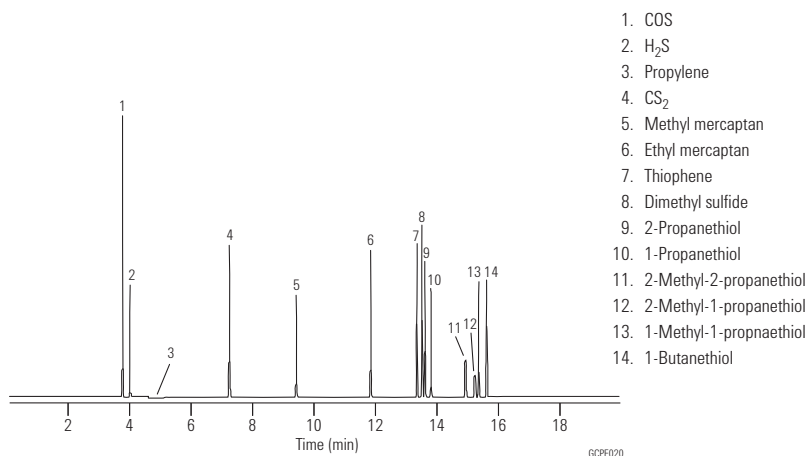**Sulfur Impurities in Propylene**

Column: Select Low Sulfur
CP8575
60 m x 0.32 mm

Oven: 65 °C for 4 min, 30 °C/min to 120 °C for 5 min


Carrier: Helium, constant flow, 2.0 mL/min

Injection: Gas sampling valve
220 °C, split 1:10

Detector: SCD, 200 °C

Sample: Polypropylene matrix containing
~300 ppb H₂S and CH₃SH, ~500 ppb COS

Injection Volume: 1 mL

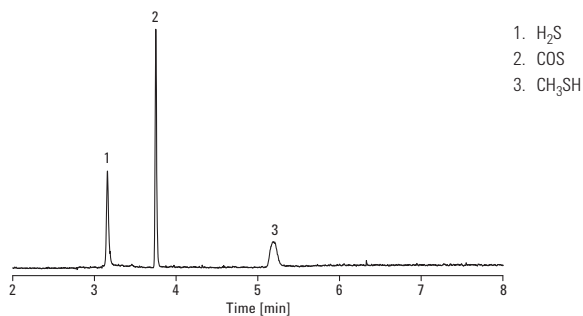

C₁ to C₄ Hydrocarbon Mix

Column: PoraPLOT Q PT
CP7550PT
10 m x 0.32 mm, 10.00 µm

Carrier: Helium, 1 mL/min in constant flow mode


Oven: 50 °C (5 min) then to 120 °C at 50 °C/min, hold 4.6 min

Sampler: Headspace unit
Oven 40 °C, valve 50 °C, transfer line 60 °C

Detector: FID or TCD at 250 °C

Injection Volume: 0.1 mL loop fitted to inlet valve of headspace unit

Inlet: Split mode at 5:1, typically at 70 °C or higher depending on column oven initial conditions

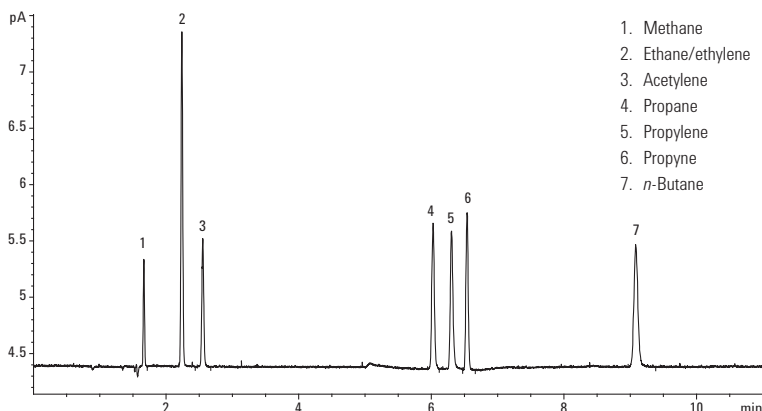

PoraPLOT Q PT, 10 m x 0.32 mm, with attached manufacturer-prepared integrated dual-ended particle trap, showing the absence of particles or spikes on FID.

Column: PoraPLOT U PT
CP7584PT
25 m x 0.53 mm, 20.00 µm

Carrier: Helium, 2 mL/min in constant flow mode


Oven: 85 °C isothermal

Sampler: Headspace unit
Oven 40 °C, valve 50 °C, transfer line 60 °C

Detector: FID or TCD at 250 °C

Injection Volume: 0.1 mL loop fitted to inlet valve of headspace unit

Inlet: Split mode at 5:1, typically at 70 °C or higher depending on column oven initial conditions

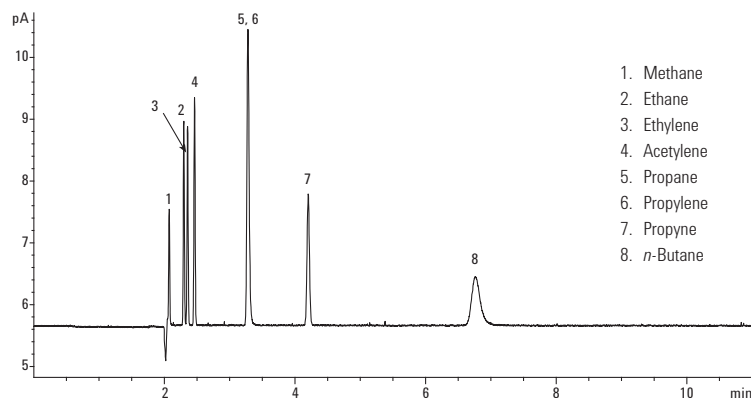

PoraPLOT U PT, 25 m x 0.53 mm, 20 µm film, with attached manufacturer-prepared integrated dual-ended particle trap, showing the lack of particles or spikes on FID.

Column: HP-PLOT Al₂O₃ KCI PT
19095P-K25PT
50 m x 0.53 mm, 15.00 µm

Carrier: Helium, 3 mL/min in constant flow mode


Oven: 100 °C (10 min) then to 120 °C at 30 °C/min, hold 3 min

Sampler: Headspace unit
Oven 40 °C, valve 50 °C, transfer line 60 °C

Detector: FID or TCD at 250 °C

Injection Volume: 0.1 mL loop fitted to inlet valve of headspace unit

Inlet: Split mode at 5:1, typically at 70 °C or higher depending on column oven initial conditions

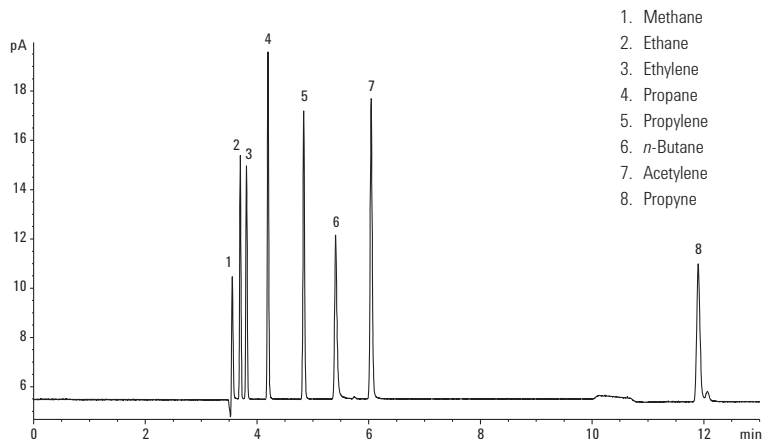


HP-PLOT Al₂O₃ KCI PT, 50 m x 0.53 mm, 15 µm film, with integrated dual-ended particle trap, showing lack of particles or spikes on FID.

Trace Sulfur Compounds in Methane (50 ppbv)

Column: Select Low Sulfur
CP8575
60 m x 0.32 mm

Oven: 40 °C (6 min), to 120 °C at 6 °C/min,
to 180 °C (5 min) at 10 °C/min

Sample: 1 mL, split ratio: 3:1

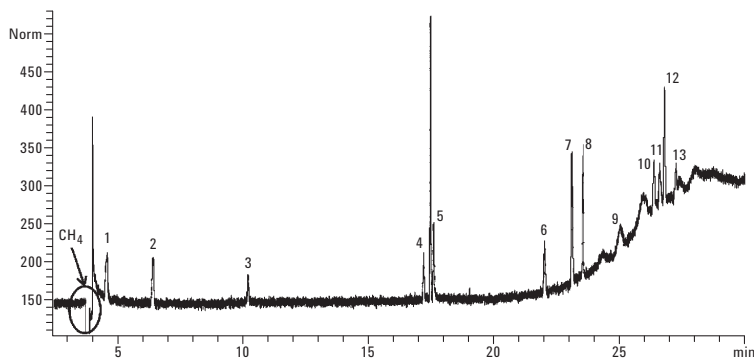

Compound	Signal/noise
1. Hydrogen sulfide	3.8
2. Carbonyl sulfide	4.0
3. Methylmercaptan	2.2
4. Ethylmercaptan	3.8
5. Dimethyl sulfide	6.3
6. 2-Propanethiol	4.3
7. Methyl ethyl sulfide	11
8. Thiophene	11
9. tert-Butyl mercaptan	2.1
10. 2-Butanethiol	4.5
11. 2-Methyl-1 propanethiol	3.7
12. Diethyl sulfide	9.8
13. 1-Butanethiol	2.4

Trace Oxygenates
in Light Hydrocarbon Matrices

Column: DB-1
125-102J
25 m x 0.53 mm, 1.00 µm


Column: GS-OxyPLOT
115-4912
10 m x 0.53 mm

Carrier: Helium (tm = 0.96 min at 50 °C)

Oven: 50 °C for 5 min
50 °C to 240 °C

Injection: Split

Detector: FID

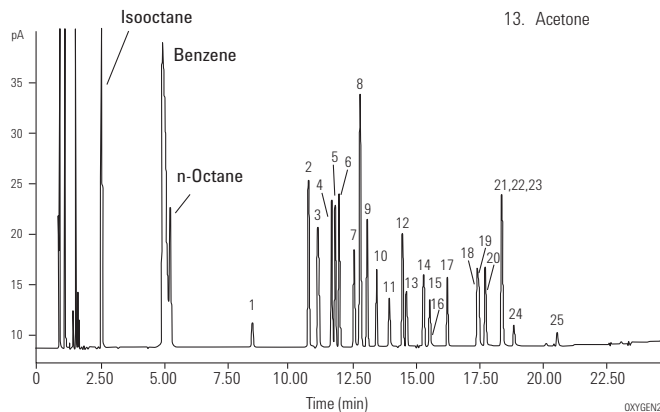

1. Dimethyl ether	14. Isovaleraldehyde
2. Diethyl ether	15. Valeraldehyde
3. Acetaldehyde	16. Methyl ethyl ketone
4. Ethyl t-butyl ether	17. Ethanol
5. Methyl t-butyl ether	18. n-Propanol
6. Diisopropyl ether	19. Isopropanol
7. Propionaldehyde	20. Allyl alcohol
8. Tert-amyl methyl ether	21. Isobutanol
9. Propyl ether	22. t-Butyl alcohol
10. Isobutylaldehyde	23. s-Butyl alcohol
11. Butylaldehyde	24. n-Butyl alcohol
12. Methanol	25. 2-Methyl-2 pentanol
13. Acetone	

Selected Oxygenates


Column: GS-OxyPLOT
115-4912
10 m x 0.53 mm

Carrier: Helium at 41 cm/s

Oven: 150 °C isothermal

Injection: Split, 1:40, 250 °C

Detector: FID, 290 °C

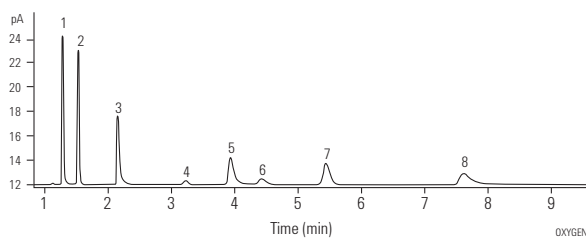

1. n-Dodecane
2. Methyl t-butyl ether
3. n-Tridecane
4. Iso-Butyraldehyde
5. n-Tetradecane
6. Methanol
7. Acetone
8. n-Pentadecane

Noble Gases

Column: HP-PLOT Molesieve
19095P-MS0
30 m x 0.53 mm, 50.00 µm

Carrier: Helium, 4 mL/min

Oven: 35 °C for 3 min
35-120 °C at 25 °C/min
120 °C for 5 min

Injection: Split ratio 50:1

Detector: TCD


Sample: 250 µL

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

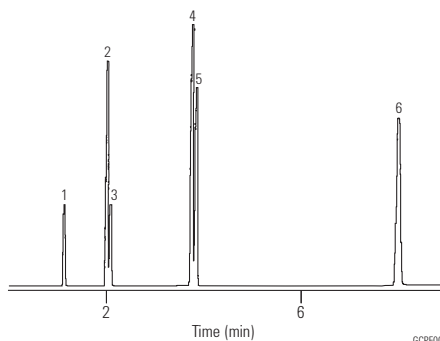

1. Neon
2. Argon
3. Oxygen
4. Nitrogen
5. Krypton
6. Xenon

Permanent Gases

Column: HP-PLOT Molesieve
19091P-MS4
30 m x 0.32 mm, 12.00 µm

Carrier: Helium, 2 mL/min

Oven: 40 °C isothermal

Injection: Split ratio 75:1

Detector: TCD


Sample: 250 µL

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

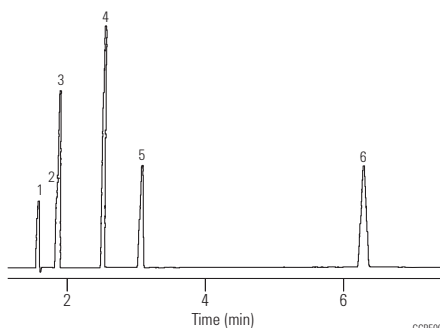

1. Neon
2. Argon
3. Oxygen
4. Nitrogen
5. Methane
6. Carbon monoxide

Baseline Resolution of Air/CO, CO₂, and Methane in a Natural Gas Sample

Column: HP-PLOT Q
19095P-Q04
30 m x 0.53 mm, 40.00 µm

Carrier: Helium (8.6 mL/min at 60 °C)

Oven: 60 °C for 2 min
60-240 °C at 30 °C/min
240 °C for 1 min

Injection: Split ratio 12:1

Detector: TCD, 250 °C


Sample: 0.25 cc natural gas sample, methane, 80%+

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

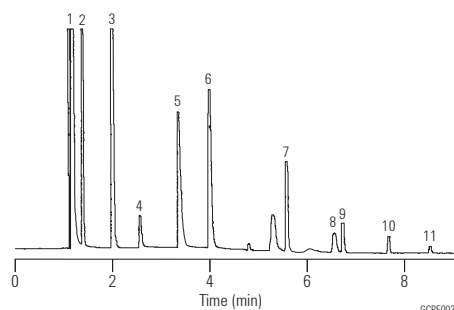

1. Air/CO
2. CO₂
3. Ethane
4. H₂S
5. Water
6. C₃
7. i-C₄/n-C₄
8. neo-C₅
9. i-C₅/n-C₅
10. C₆
11. C₇

Natural Gas

Column: HP-PLOT Al₂O₃ S
19095P-S21
15 m x 0.53 mm, 15.00 µm

Carrier: Helium, 50 cm/s (100 °C), 6 mL/min

Oven: 100 °C for 1.5 min
100-180 °C at 30 °C/min

Injection: Split, 250 °C
Split ratio 50:1

Detector: FID, 250 °C


Sample: 5 µL natural gas, p/n 5080-8756

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

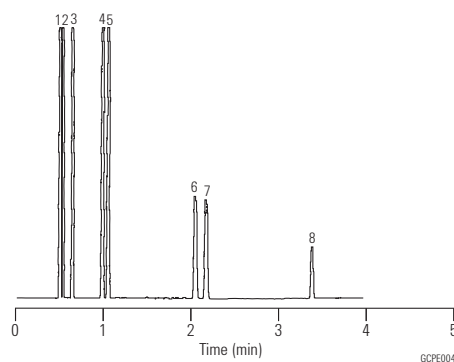

1. Methane
2. Ethane
3. Propane
4. iso-Butane
5. n-Butane
6. iso-Pentane
7. n-Pentane
8. n-Hexane

Ethylene

Column: HP-PLLOT Al₂O₃ S
19095P-S25
50 m x 0.53 mm, 15.00 µm


Carrier: Helium, 50 cm/s (35 °C),
7 mL/min constant flow

Oven: 35 °C for 2 min
35-100 °C at 5 °C/min

Injection: Split, 250 °C
Split ratio 65:1

Detector: FID, 250 °C

Sample: 5 µL
ethylene 98.4%

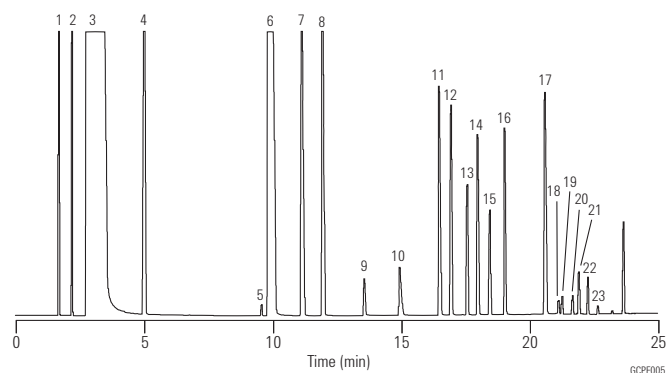

1. Methane
2. Ethane
3. Ethylene
4. Propane
5. Cyclopropane
6. Propylene
7. Isobutane
8. n-Butane
9. Propadiene
10. Acetylene
11. trans-2-Butene
12. Butene-1
13. Isobutylene
14. cis-2-Butene
15. Isopentane
16. n-Pentane
17. 1,3-Butadiene
18. Propyne
19. trans-2-Pentene
20. 2-Methyl-2-butene
21. Pentene-1
22. cis-2-Pentene
23. n-Hexane

Impurities in Ethylene

Column: GS-Alumina KCl
115-3352
50 m x 0.53 mm


Carrier: Helium at 8 mL/min, measured at 35 °C

Oven: 35 °C for 2 min
35-190 °C at 10 °C/min
190 °C for 3 min

Injection: Split, 200 °C
Split ratio 1:40

Detector: FID, 200 °C
Nitrogen makeup gas at 20 mL/min

Sample: 0.2 mL of trace hydrocarbons
in ethylene

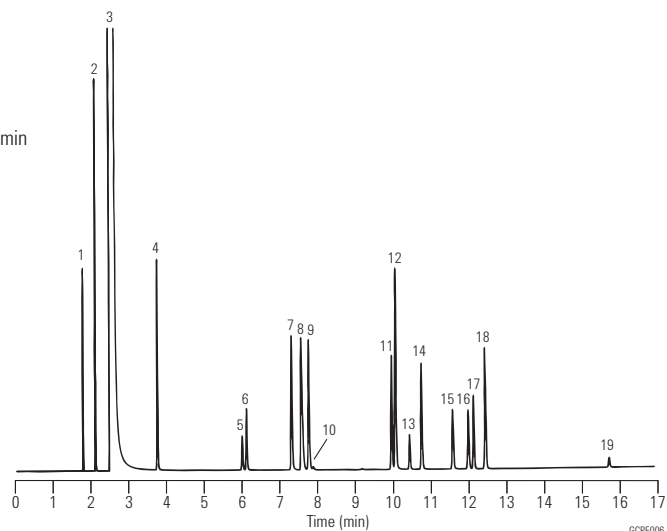

1. Methane
2. Ethane
3. Ethylene
4. Propane
5. Cyclopropane
6. Propylene
7. Isobutane
8. Acetylene
9. n-Butane
10. Propadiene
11. trans-2-Butene
12. 1-Butene
13. Isobutylene
14. cis-2-Butene
15. Isopentane
16. n-Pentane
17. Propyne
18. 1,3-Butadiene
19. 1-Pentene

Impurities in Propylene

Column: GS-Alumina KCl
115-3352
50 m x 0.53 mm

Carrier: Helium at 10 mL/min,
measured at 35 °C

Oven: 35 °C for 2 min
35-190 °C at 10 °C/min
190 °C for 3 min

Injection: Split, 200 °C
Split ratio 1:30

Detector: FID, 200 °C
Nitrogen makeup gas
at 20 mL/min


Sample: 0.2 mL of trace
hydrocarbons in propylene

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

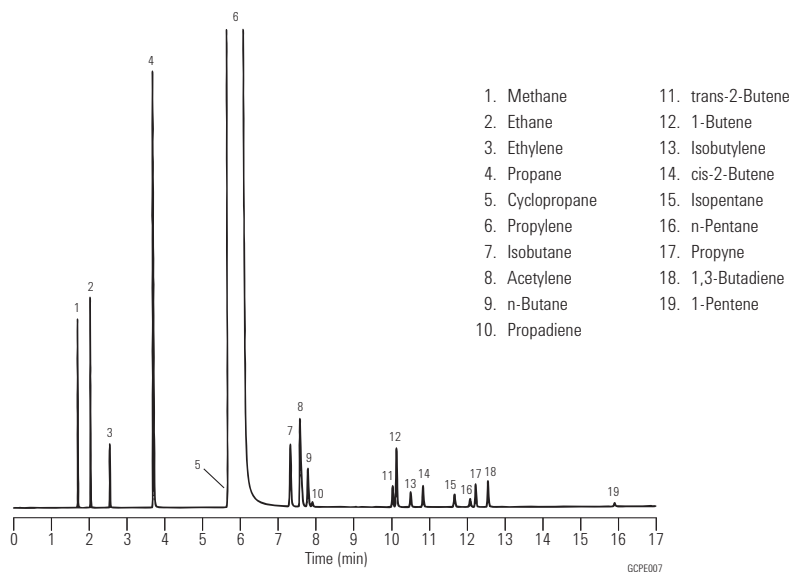**Propylene**

Column: GS-Alumina
115-3552
50 m x 0.53 mm

Carrier: Helium at 10 mL/min,
measured at 35 °C

Oven: 35 °C for 2 min
35-190 °C at 10 °C/min
190 °C for 3 min

Injection: Split, 200 °C
Split ratio 1:30

Detector: FID, 200 °C
Nitrogen makeup gas
at 20 mL/min


Sample: 0.2 mL of trace
hydrocarbons in propylene

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

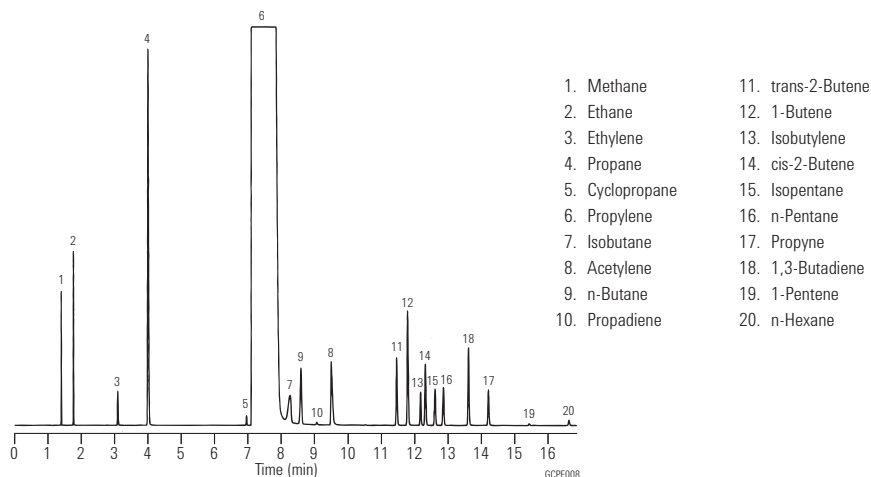

1,3-Butadiene

Column: DB-624
128-1324
25 m x 0.20 mm, 1.12 µm

Carrier: Helium at 1.0 mL/min

Oven: -20 °C for 3 min
-20 °C to 20 °C at 4 °C/min
20 °C to 200 °C at 8 °C/min
200 °C for 10 min

Injection: Split, 250 °C
Split ratio 1:150

Detector: FID, 250 °C

Sample: 0.5 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

*Agilent Technologies wishes to thank DCG Industries
(Pearland, TX) for providing this chromatogram.*

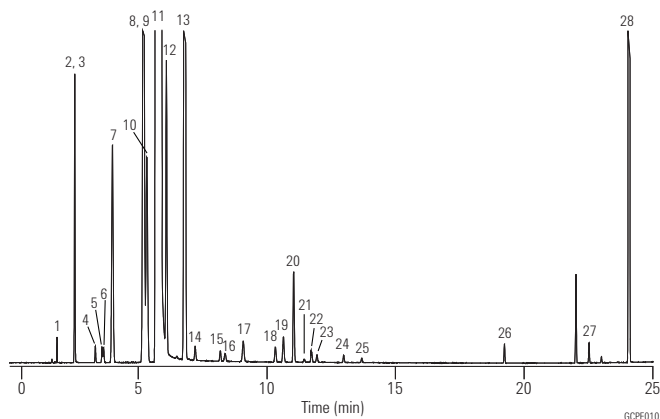**Refined Butadiene Standard
Component****Gravimetric concentration
(PPM)**

1. Acetylene	20.7
2. Propane	19.8
3. Propylene	296
4. Propadiene (allene)	21.1
5. Propyne (methylacetylene)	21
6. Cyclopropane	20
7. Isobutane	506
8. Butene-1	999
9. Isobutylene	495
10. n-Butane	494
11. 1,3-Butadiene	balance
12. trans-2-Butene	442
13. cis-2-Butene	1946
14. 1-Butyne (ethylacetylene)	20.2
15. 1,2-Butadiene	28.9
16. 3-Methyl-1-butene	19.8
17. Isopentane	50.1
18. Pentene-1	29.8
19. n-Pentane	50.1
20. 2-Butyne (dimethylacetylene)	150
21. trans-2-Pentene	5.57
22. Isoprene	20
23. cis-2-Pentene	13.9
24. trans-1,3-Pentadiene	13.8
25. cis-1,3-Pentadiene	7.73
26. Benzene	20.3
27. Toluene	20.2
28. Dimer (4-vinylcyclohexene-1)	

1,3-Butadiene Purity

Column: GS-Alumina
115-3552
50 m x 0.53 mm

Carrier: Helium, 6.0 mL/min
(constant flow mode)

Oven: 45 °C for 3 min
6 °C/min to 195 °C
195 °C for 15 min

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 250 °C

Sample: 0.5 µL

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

Refined Butadiene Standard Component**Gravimetric concentration (PPM)**

1. Acetylene	20.7
2. Propane	19.8
3. Propylene	296
4. Propadiene (allene)	21.1
5. Propyne (methylacetylene)	21
6. Cyclopropane	20
7. Isobutane	506
8. Butene-1	999
9. Isobutylene	495
10. n-Butane	494
11. 1,3-Butadiene	Balance
12. trans-2-Butene	442
13. cis-2-Butene	1946
14. 1-Butyne (ethylacetylene)	20.2
15. 1,2-Butadiene	28.9
16. 3-Methyl-1-butene	19.8
17. Isopentane	50.1
18. Pentene-1	29.8
19. n-Pentane	50.1
20. 2-Butyne (dimethylacetylene)	150
21. trans-2-Pentene	5.57
22. Isoprene	20
23. cis-2-Pentene	13.9
24. trans-1,3-Pentadiene	13.8
25. cis-1,3-Pentadiene	7.73
26. Benzene	20.3
27. Toluene	20.2
28. Dimer (4-vinylcyclohexene-1)	

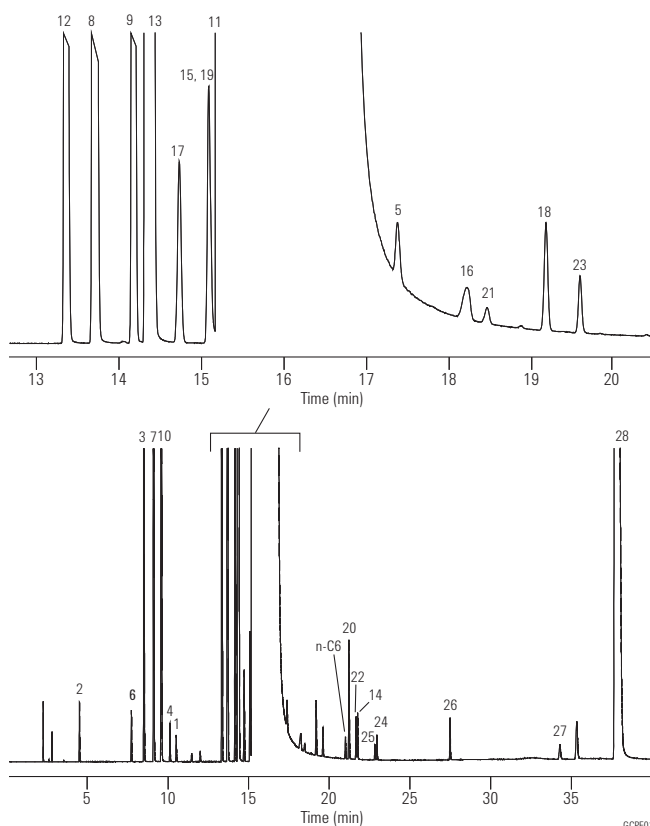

GCPED011

Extended Hydrocarbon Analysis I

Column: GS-Alumina
115-3532
30 m x 0.53 mm

Carrier: Helium at 52 cm/s (6.7 mL/min),
measured at 100 °C

Oven: 100 °C for 1 min
100-140 °C at 8 °C/min
140 °C for 0.5 min
140-200 °C at 30 °C/min

Injection: Split, 250 °C
Split ratio 1:8

Detector: FID, 275 °C
Nitrogen makeup gas at 29 mL/min


Sample: 300 µL injection of 100 ppmv
SUMMA canister mixture

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

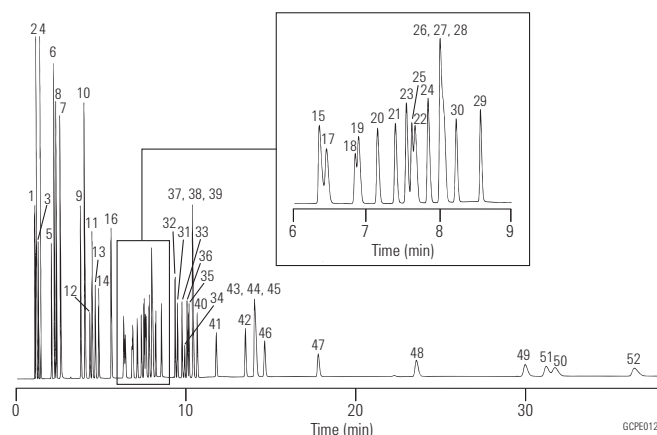

- | | |
|------------------------|--|
| 1. Methane | 27. 2-Methylpentane |
| 2. Ethane | 28. 3-Methylpentane |
| 3. Ethylene | 29. Isoprene |
| 4. Propane | 30. n-Hexane |
| 5. Propylene | 31. 4-Methyl-1-pentene |
| 6. Isobutane | 32. trans-2-Hexene |
| 7. Acetylene | 33. 2-Methyl-1-pentene |
| 8. n-Butane | 34. cis-2-Hexene |
| 9. trans-2-Butene | 35. 2,4-Dimethylpentane |
| 10. 1-Butene | 36. Methylcyclohexane |
| 11. cis-2-Butene | 37. 2,3-Dimethylpentane |
| 12. Cyclopentane | 38. 2-Methylhexane |
| 13. Isopentane | 39. 3-Methylhexane |
| 14. n-Pentane | 40. n-Heptane |
| 15. Propyne | 41. Benzene |
| 16. 1,3-Butadiene | 42. Isooctane (2,2,4-trimethylpentane) |
| 17. Cyclopentene | 43. 2,3,4-Trimethylpentane |
| 18. 3-Methyl-1-butene | 44. 3-Methylheptane |
| 19. trans-2-Pentene | 45. 2-Methylheptane |
| 20. 2-Methyl-2-butene | 46. n-Octane |
| 21. 1-Pentene | 47. Toluene |
| 22. cis-2-Pentene | 48. n-Nonane |
| 23. Methylcyclopentane | 49. Ethylbenzene |
| 24. 2,2-Dimethylbutane | 50. m-Xylene |
| 25. Cyclohexane | 51. p-Xylene |
| 26. 2,3-Dimethylbutane | 52. o-Xylene |

Extended Hydrocarbon Analysis II

Column: GS-GasPro
113-4362
60 m x 0.32 mm

Carrier: Helium at 40 cm/s (3.3 mL/min),
measured at 80 °C

Oven: 80 °C for 0.5 min
80-175 °C at 25 °C/min
175 °C for 2 min
175-250 °C at 25 °C/min

Injection: Split, 250 °C
Split ratio 1:17

Detector: FID, 275 °C
Nitrogen makeup gas at 32 mL/min


Sample: 500 µL injection of 100 ppmv
SUMMA canister mixture

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

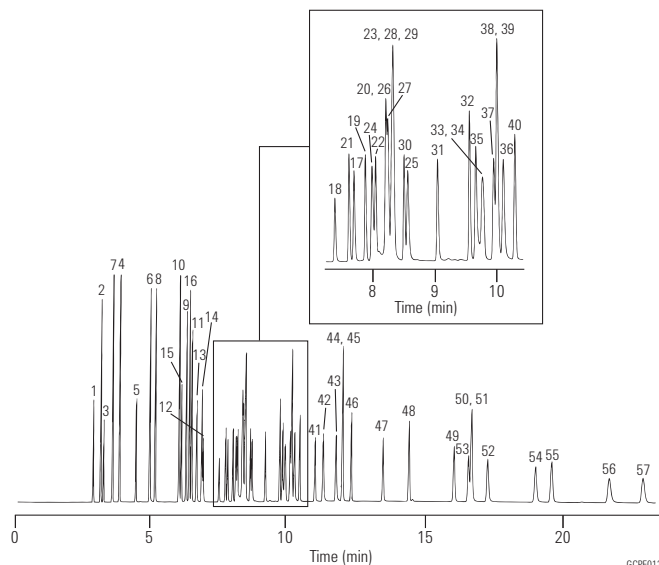

- | | |
|------------------------|--|
| 1. Methane | 30. n-Hexane |
| 2. Ethane | 31. 4-Methyl-1-pentene |
| 3. Ethylene | 32. trans-2-Hexene |
| 4. Propane | 33. 2-Methyl-1-pentene |
| 5. Propylene | 34. cis-2-Hexene |
| 6. Isobutane | 35. 2,4-Dimethylpentane |
| 7. Acetylene | 36. Methylcyclohexane |
| 8. n-Butane | 37. 2,3-Dimethylpentane |
| 9. trans-2-Butene | 38. 2-Methylhexane |
| 10. 1-Butene | 39. 3-Methylhexane |
| 11. cis-2-Butene | 40. n-Heptane |
| 12. Cyclopentane | 41. Benzene |
| 13. Isopentane | 42. Isooctane (2,2,4-trimethylpentane) |
| 14. n-Pentane | 43. 2,3,4-Trimethylpentane |
| 15. Propyne | 44. 3-Methylheptane |
| 16. 1,3-Butadiene | 45. 2-Methylheptane |
| 17. Cyclopentene | 46. n-Octane |
| 18. 3-Methyl-1-butene | 47. Toluene |
| 19. trans-2-Pentene | 48. n-Nonane |
| 20. 2-Methyl-2-butene | 49. Ethylbenzene |
| 21. 1-Pentene | 50. m-Xylene |
| 22. cis-2-Pentene | 51. p-Xylene |
| 23. Methylcyclopentane | 52. o-Xylene |
| 24. 2,2-Dimethylbutane | 53. Styrene |
| 25. Cyclohexane | 54. Isopropylbenzene (cumene) |
| 26. 2,3-Dimethylbutane | 55. n-Propylbenzene |
| 27. 2-Methylpentane | 56. 1,3,5-Trimethylbenzene |
| 28. 3-Methylpentane | 57. 1,2,4-Trimethylbenzene |
| 29. Isoprene | |

Refinery Gas

Column: HP-PLLOT Al₂O₃ S
19095P-S25
50 m x 0.53 mm, 15.00 µm

Carrier: Helium 7 mL/min

Oven: 100 °C isothermal

Injection: Split, 250 °C
Split ratio 100:1

Detector: FID, 250 °C


Sample: 5 µL

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

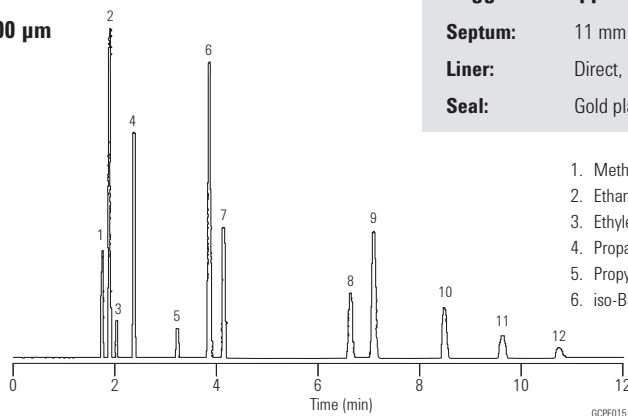

- | | |
|---------------|-------------------|
| 1. Methane | 7. n-Butane |
| 2. Ethane | 8. trans-2-Butene |
| 3. Ethylene | 9. 1-Butene |
| 4. Propane | 10. cis-2 Butene |
| 5. Propylene | 11. iso-Pentane |
| 6. iso-Butane | 12. n-Pentane |

**Sulfur Gas Analysis
in Light Hydrocarbon Streams I**

Column: GS-GasPro
113-4332
30 m x 0.32 mm


Carrier: Helium, 10 psig, 2.0 mL/min at 60 °C

Oven: 60 °C for 2 min, 20 °C/min to 260 °C
and hold

Injection: Split, 200 °C
Split ratio 1:20

Detector: Two separate analyses under identical
conditions on FID and PFPD


PFPD

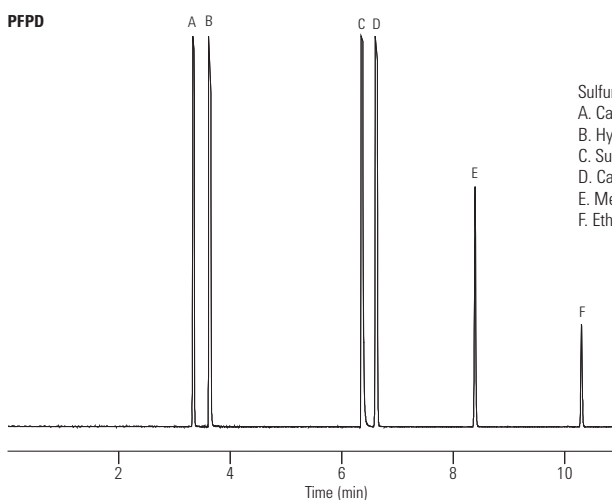

Sulfur compounds (PFPD)

- A. Carbonyl sulfide
B. Hydrogen sulfide
C. Sulfur dioxide
D. Carbon disulfide
E. Methyl mercaptan
F. Ethyl mercaptan

FID

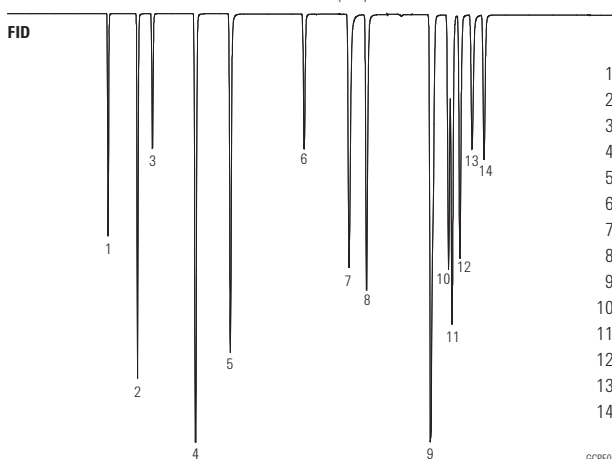**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

- | |
|------------------------------|
| 1. Methane |
| 2. Ethane |
| 3. Ethylene |
| 4. Acetylene |
| 5. Propane |
| 6. Propylene |
| 7. iso-Butane |
| 8. n-Butane |
| 9. 1-Butene/methyl acetylene |
| 10. trans-2-Butene |
| 11. 1,3-Butadiene |
| 12. cis-2-Butene |
| 13. iso-Pentane |
| 14. n-Pentane |

Sulfur Gas Analysis in Light Hydrocarbon Streams II

Column: GS-Q
113-3432
30 m x 0.32 mm, 0.20 µm


Carrier: Helium, 10 psig, 1.7 mL/min at 100 °C

Oven: 100 °C for 2 min, 20 °C/min to 250 °C and hold

Injection: Split, 200 °C
Split ratio 1:20

Detector: Two separate analyses under identical conditions on FID and PFPD

PFPD

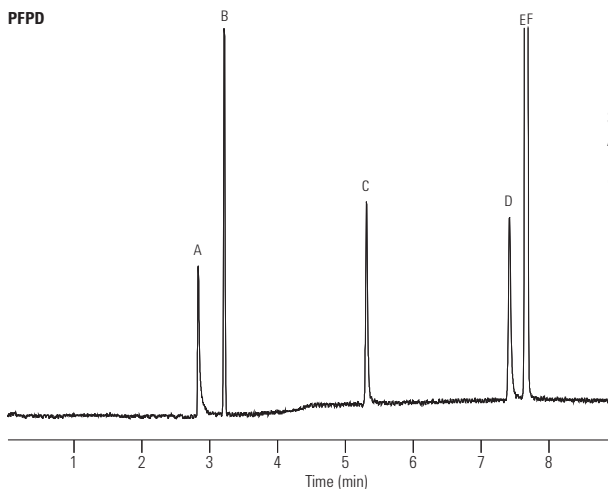

Sulfur compounds (PFPD)
A. Hydrogen sulfide
B. Carbonyl sulfide
C. Methyl mercaptan
D. Ethyl mercaptan
E. Carbon disulfide


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

FID

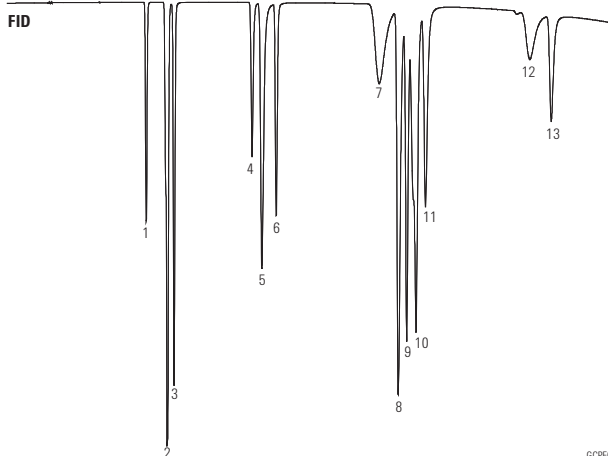

1. Methane
2. Ethylene/acetylene
3. Ethane
4. Propylene
5. Propane
6. Methyl acetylene
7. iso-Butane
8. 1-Butene
9. 1,3-Butadiene
10. n-Butane/cis-2-butene
11. trans-2-Butene
12. iso-Pentane
13. n-Pentane

GCPE019

Sulfur Compounds in Propylene (1 ppm)


Column: GS-GasPro
113-4332
30 m x 0.32 mm

Oven: 60 °C for 2.5 min
60-250 °C at 10 °C/min

Injection: OI Analytical Volatiles Inlet
Split ratio 5:1
200 µL gas sampling valve

Detector: OI Analytical Model 5380 PFPD

Sample: 1 ppm sulfur compounds in propylene

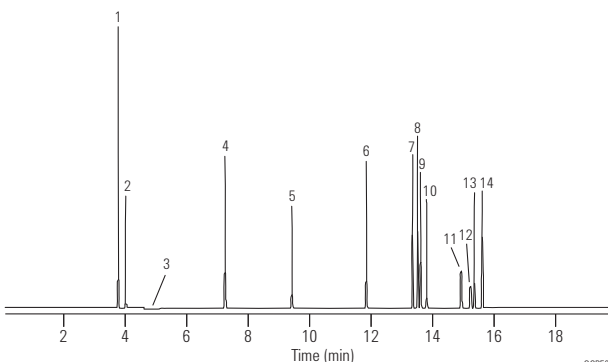

1. COS
2. H₂S
3. Propylene
4. CS₂
5. Methyl mercaptan
6. Ethyl mercaptan
7. Thiophene
8. Dimethyl sulfide
9. 2-Propanethiol
10. 1-Propanethiol
11. 2-Methyl-2-propanethiol
12. 2-Methyl-1-propanethiol
13. 1-Methyl-1-propanethiol
14. 1-Butanethiol

Chromatogram courtesy of OI Analytical

GCPE020

Mercaptans

Column: GS-GasPro
113-4332
30 m x 0.32 mm

Carrier: Helium at 25 cm/s

Oven: 175 °C for 2 min
175-260 °C at 10 °C/min

Injection: Split
Split flow 80 mL/min

Detector: FID


Sample: 0.2 mL

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

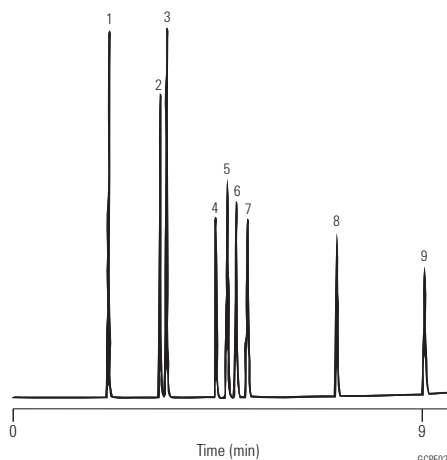

1. Ethyl mercaptan
2. 2-Propyl mercaptan
3. 1-Propyl mercaptan
4. 2-Methyl-2-propyl mercaptan
5. 2-Methyl-1-propyl mercaptan
6. 1-Methyl-1-propyl mercaptan
7. 1-Butyl mercaptan
8. 1-Pentyl mercaptan
9. 1-Hexyl mercaptan

**Sulfur Compounds in Natural Gas –
Synthetic Mixture**

Column: HP-1
19091Z-205
50 m x 0.20 mm, 0.50 µm

Carrier: Helium

Oven: 35 °C for 10 min
35-300 °C at 7 °C/min

Injection: Split 100:1

Detector: FPD


Sample: 0.5 mL

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

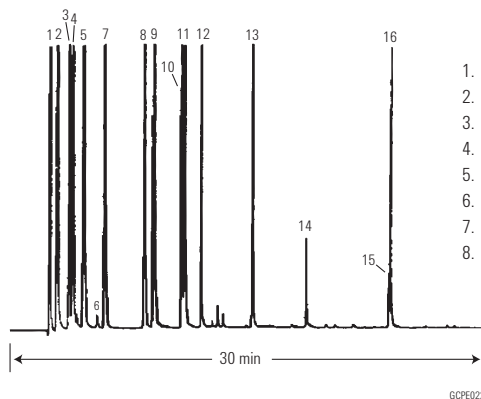

- | | |
|--------------------------------------|----------------------------|
| 1. Hydrogen sulfide | 9. Isobutyl mercaptan |
| 2. Methyl mercaptan | 10. n-Butyl mercaptan |
| 3. Ethyl mercaptan | 11. tert-Amyl mercaptan |
| 4. Dimethyl sulfide | 12. Isoamyl mercaptan |
| 5. Isopropyl mercaptan | 13. n-Amyl mercaptan |
| 6. tert-Butyl mercaptan | 14. n-Hexyl mercaptan |
| 7. n-Propyl mercaptan | 15. tert-Dibutyl disulfide |
| 8. Thiophene and sec-butyl mercaptan | 16. n-Octyl mercaptan |

Sulfur Compounds in Naphtha

Column: HP-PONA
19091S-001
50 m x 0.20 mm, 0.50 μ m


Carrier: Helium, 26 cm/s

Oven: 35 °C for 15 min
35-70 °C at 8 °C/min
70-130 °C at 15 °C/min

Injection: Split ratio 400:1

Detector: FPD

Sample: 3 μ L

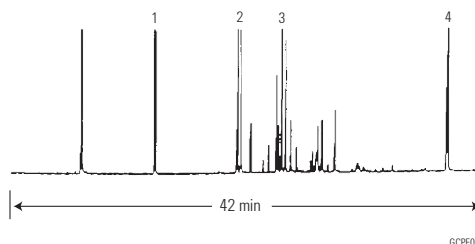

1. Thiophene
2. Methyl thiophenes
3. Ethyl and dimethyl thiophenes
4. Benzothiophene

GCPE023

Aromatics Analysis – ASTM D16 Analytes


Column: HP-INNOWax
19091N-216
60 m x 0.32 mm, 0.50 μ m

Carrier: Helium at 20 psi, constant pressure mode

Oven: 75 °C for 10 min
3 °C/min to 100 °C
10 °C/min to 145 °C

Injection: Split, 250 °C
Split ratio 100:1 to 400:1

Detector: FID, 250 °C
Data acquisition rate at 20 Hz

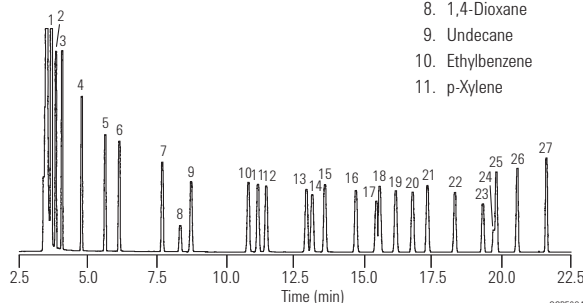

- | | |
|------------------|-----------------------------|
| 1. Heptane | 12. m-Xylene |
| 2. Cyclohexane | 13. Cumene |
| 3. Octane | 14. Dodecane |
| 4. Nonane | 15. o-Xylene |
| 5. Benzene | 16. Propylbenzene |
| 6. Decane | 17. p-Ethyltoluene |
| 7. Toluene | 18. m-Ethyltoluene |
| 8. 1,4-Dioxane | 19. tert-Butylbenzene |
| 9. Undecane | 20. sec-Butylbenzene |
| 10. Ethylbenzene | 21. Styrene |
| 11. p-Xylene | 22. Tridecane |
| | 23. Diethylbenzene isomer |
| | 24. Diethylbenzene isomer |
| | 25. n-Butylbenzene |
| | 26. α -Methylstyrene |
| | 27. Phenylacetylene |

GCPE024

Aromatics Analysis – Ethylbenzene Impurities


Column: HP-INNOWax
19091N-216
60 m x 0.32 mm, 0.50 μ m

Carrier: Helium at 20 psi, constant pressure mode

Oven: 75 °C for 10 min
3 °C/min to 100 °C
10 °C/min to 145 °C

Injection: Split, 250 °C
Split ratio 100:1 to 400:1

Detector: FID, 250 °C
Data acquisition rate at 20 Hz

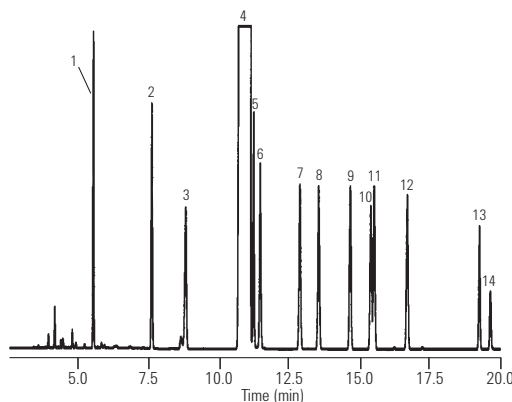

1. Benzene
2. Toluene
3. Undecane
4. Ethylbenzene
5. p-Xylene
6. m-Xylene
7. Isopropylbenzene
8. o-Xylene
9. n-Propylbenzene
10. p-Ethyltoluene
11. m-Ethyltoluene
12. s-Butylbenzene
13. Diethylbenzene
14. Diethylbenzene

GCPE025

Impurities in p-Xylene – ASTM D3798

Column: HP-INNOWax
19091N-216
60 m x 0.32 mm, 0.50 µm


Carrier: Helium, 32 cm/s, 19.9 psi (60 °C),
2.5 mL/min constant flow

Oven: 60 °C for 1 min
60-92 °C at 4 °C/min
92 °C for 4.5 min
92-220 °C at 20 °C/min
220 °C for 5 min

Injection: Split, 220 °C
Split ratio 100:1

Detector: FID, 270 °C

Sample: 0.5 µL
Neat, 99%+

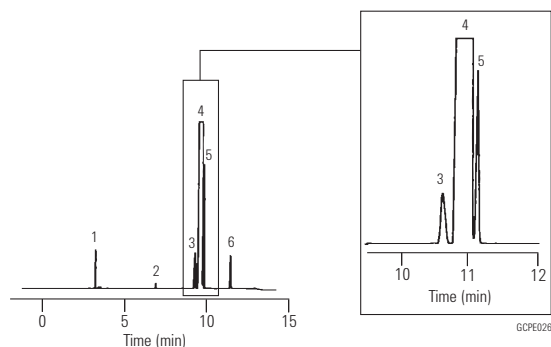

1. Non-aromatic hydrocarbon
2. Toluene
3. Ethylbenzene
4. p-Xylene
5. m-Xylene
6. o-Xylene

Ethylene Oxide Synthetic Standard

Column: HP-PLLOT Q
19095P-Q04
30 m x 0.53 mm, 40.00 µm

Carrier: Helium, 25 psi

Oven: 50 °C for 2 min
50-250 °C at 15 °C/min

Injection: Split ratio 40:1

Detector: FID

Sample: 1 µL liquid injection
sample 2000 ppm v/v


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

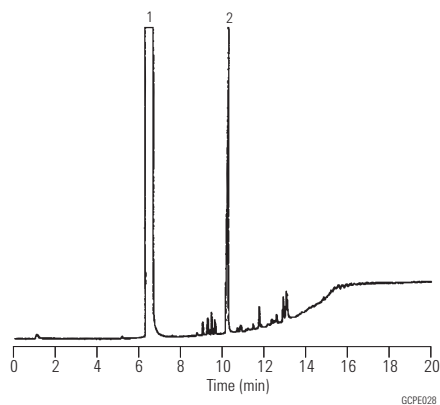

1. Ethylene oxide
2. 2-Chloropropene

Analysis of Oxygenates in Mixed C4 Streams

Column: PoraBOND Q PT
CP7351PT
25 m x 0.32 mm, 5.00 µm

Instrument: Agilent 7890A Series


Carrier: Helium, constant flow mode, 35 cm/s, 45 °C

Oven: 45-90 °C at 6 °C/min, 90-240 °C at 15 °C/min,
240 °C for 10 min

Injection: 200 °C, split ratio 30:1, 200 µL gas sampling valve

Detector: FID at 250 °C

Sample: 50-100 mg/L oxygenates in mixed C4

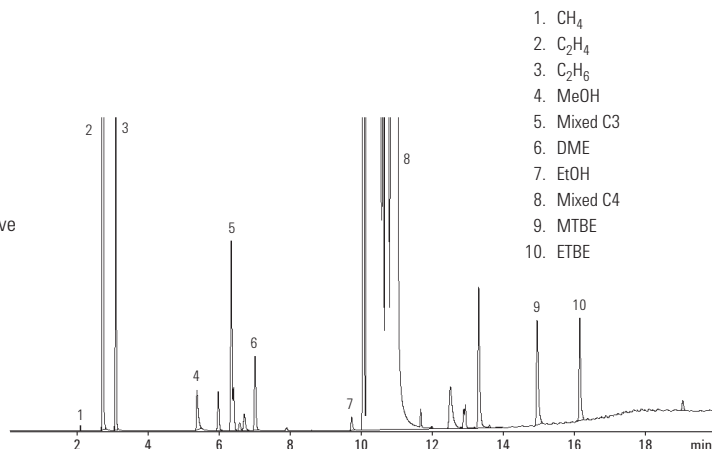

Oxygenates in Gasoline ASTM D5599 (GC-OFID)

Column: HP-1
19091Z-236
60 m x 0.25 mm, 1.00 µm


Carrier: Helium, 30 cm/s constant flow

Oven: 40 °C for 6 min
40-50 °C at 5 °C/min
50 °C for 4 min
50-175 °C at 25 °C/min
175 °C for 5 min

Injection: Split ratio 150:1

Detector: Wasson ECE OFID

Sample: 0.5 µL

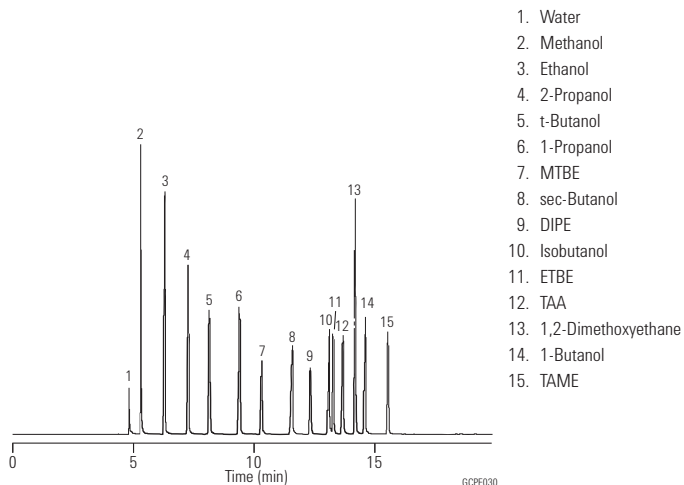

Denatured Fuel Ethanol – ASTM D5501

Column: HP-1
19091Z-530
100 m x 0.25 mm, 0.50 µm


Carrier: Helium 24 cm/s

Oven: 15 °C for 12 min
15-250 °C at 19 °C/min
250 °C for 20 min

Injection: Split ratio 200:1

Detector: FID, 250 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.5 µL

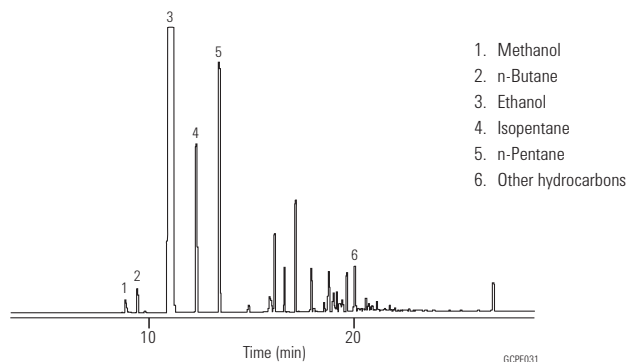**PONA Mix as Specified by AFNOR Method #2**

Column: DB-Petro
128-1056
50 m x 0.20 mm, 0.50 µm


Carrier: Helium at 16.7 cm/s, measured at 35 °C

Oven: 10 °C for 15 min
10-70 °C at 1.3 °C/min
70-250 °C at 1.7 °C/min

Injection: Split, 250 °C
Split ratio 1:200

Detector: FID, 250 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.3 µL petroleum reformat

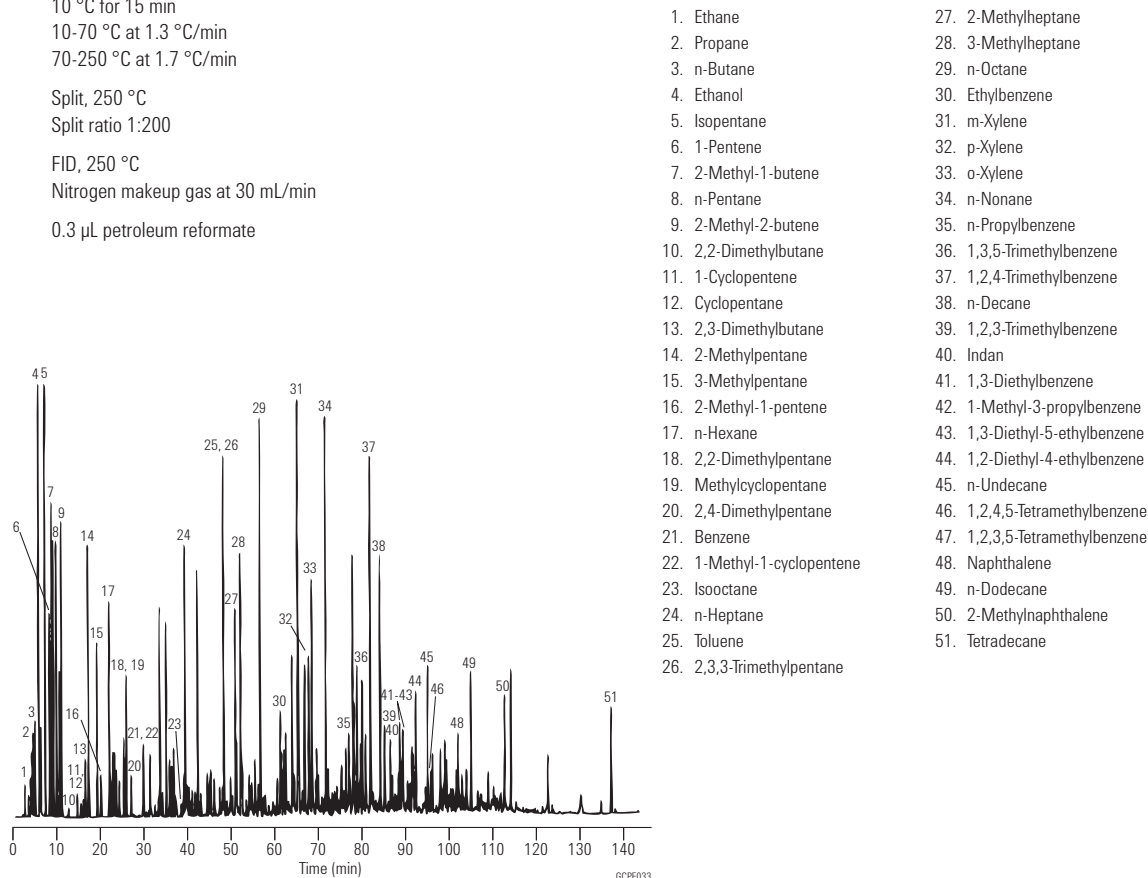

Aromatics in Finished Gasoline – ASTM Method D5769

Column: DB-1
122-1063
60 m x 0.25 mm, 1.00 µm

Carrier: Helium at 35 cm/s,
measured at 50 °C


Oven: 50 °C for 1 min
50-190 °C at 2 °C/min
190 °C for 1 min

Injection: Split, 250 °C
Split ratio 1:100

Detector: MSD

Sample: 0.3 µL unleaded gasoline
Calibration standard: ASTM/EPA gasoline
refinery aromatics
(AccuStandard M-GRA-CAL/IS-SET)

- | | | |
|-----------------------------------|-----------------------------|--------------------------------|
| 1. Methyl-tert-butyl-ether (MTBE) | 12. p-Xylene | 23. 1,4-Diethylbenzene |
| 2. n-Hexane | 13. o-Xylene | 24. n-Butylbenzene (valley) |
| 3. Benzene-d6 (IS) | 14. n-Propylbenzene | 25. 1,2-Diethylbenzene |
| 4. Benzene | 15. 1-Methyl-3-ethylbenzene | 26. 1,2,4,5-Tetramethylbenzene |
| 5. Isooctane | 16. 1-Methyl-4-ethylbenzene | 27. 1,2,3,5-Tetramethylbenzene |
| 6. n-Heptane | 17. 1,3,5-Trimethylbenzene | 28. Naphthalene-d8 (IS) |
| 7. Toluene | 18. 1-Methyl-2-ethylbenzene | 29. Naphthalene |
| 8. n-Octane | 19. 1,2,4-Trimethylbenzene | 30. n-Dodecane |
| 9. Ethylbenzene-d10 (IS) | 20. n-Decane | 31. Pentamethylbenzene |
| 10. Ethylbenzene | 21. 1,2,3-Trimethylbenzene | 32. 2-Methylnaphthalene |
| 11. m-Xylene | 22. Indan | 33. 1-Methylnaphthalene |

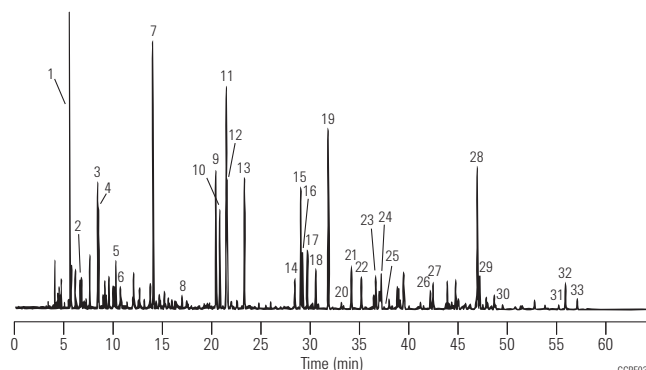

GCPE034

Simulated Distillation

Column: DB-2887
125-2814
10 m x 0.53 mm, 3.00 µm

Carrier: Helium at 7 mL/min

Oven: 35-350 °C at 15 °C/min

Injection: Direct

Detector: FID
Nitrogen makeup gas
at 30 mL/min


Suggested Supplies

Septum: Non-stick bleed and temperature optimized (BTO) septa, 11 mm, 50/pk, 5183-4757

Liner: Direct connect, dual taper, deactivated, 4 mm id, G1544-80700

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

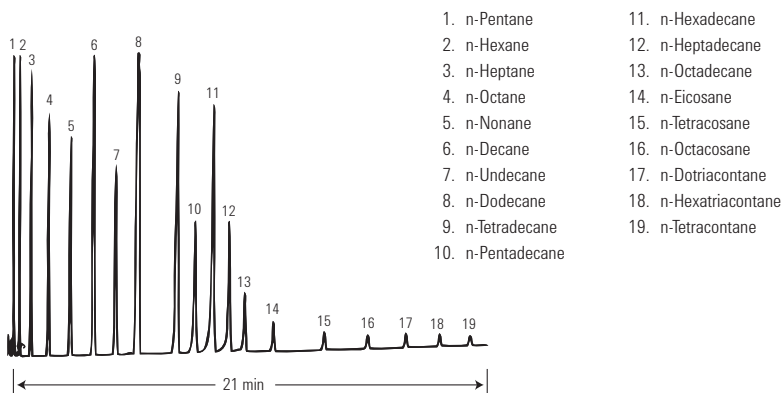


GCPE035

Reference Gas Oil

Column: DB-2887
125-2814
10 m x 0.53 mm, 3.00 μ m

Carrier: Helium at 7 mL/min
Oven: 35-350 °C at 15 °C/min
Injection: Direct
Detector: FID
Nitrogen makeup gas
at 30 mL/min

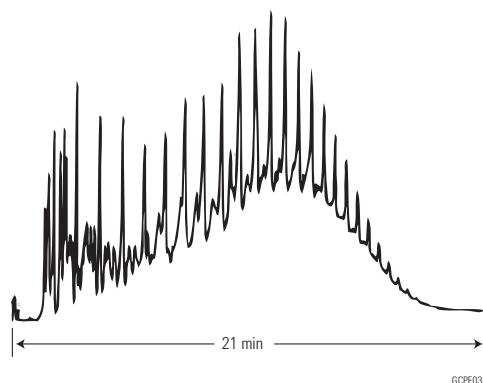
**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, dual taper, deactivated, 4 mm id, G1544-80700
Seal: Gold plated seal, 18740-20885
Syringe: 10 μ L tapered, FN 23-26s/42/HP, 5181-1267

Regular Unleaded Gasoline (California Phase 1) – "Normal" GC Run I

Column: DB-Petro
122-10A6
100 m x 0.25 mm, 0.50 μ m

Carrier: Hydrogen at 31 cm/s
Oven: 35 °C for 9.5 min
35-45 °C at 13.3 °C/min
45 °C for 11 min
45-60 °C at 1.4 °C/min
60 °C for 11 min
60-220 °C at 2.7 °C/min
220 °C for 3.6 min
Injection: Split ratio 1:200
Detector: FID, 300 °C
Sample: 0.2 μ L

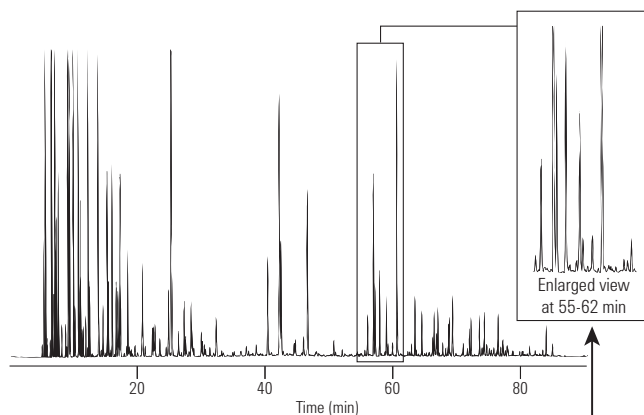


Compare Resolution

Regular Unleaded Gasoline (California Phase 1) – "Normal" GC Run II

Column: DB-1
127-1046
40 m x 0.10 mm, 0.20 μ m

Carrier: Hydrogen at 34.8 cm/s
Oven: 35 °C for 3.6 min
35-45 °C at 36.1 °C/min
45 °C for 4.2 min
45-60 °C at 3.9 °C/min
60 °C for 4.2 min
60-220 °C at 6.9 °C/min
220 °C for 1.4 min
Injection: Split ratio 1:400
Detector: FID, 300 °C
Sample: 0.2 μ L

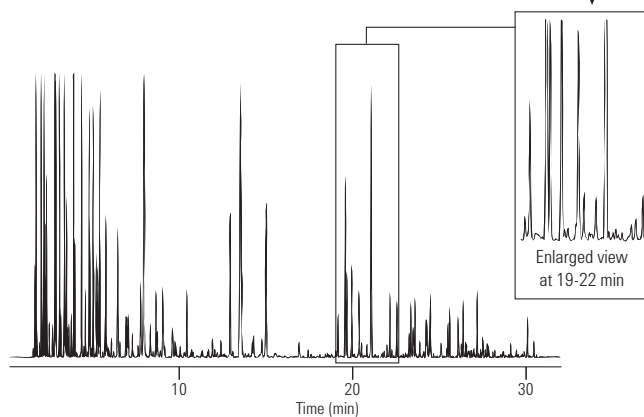

Gasoline Unleaded ASTM D5769

Column: CP-Sil PONA CB
CP7530
100 m x 0.25 mm, 0.50 μ m


Sample: 0.1 μ L

Carrier: Helium, 240 kPa (2.4 bar, 34 psi)

Oven: 35 $^{\circ}$ C (7 min) to 250 $^{\circ}$ C, 3 $^{\circ}$ C/min

Injection: Split, 80 mL/min

Detector: FID

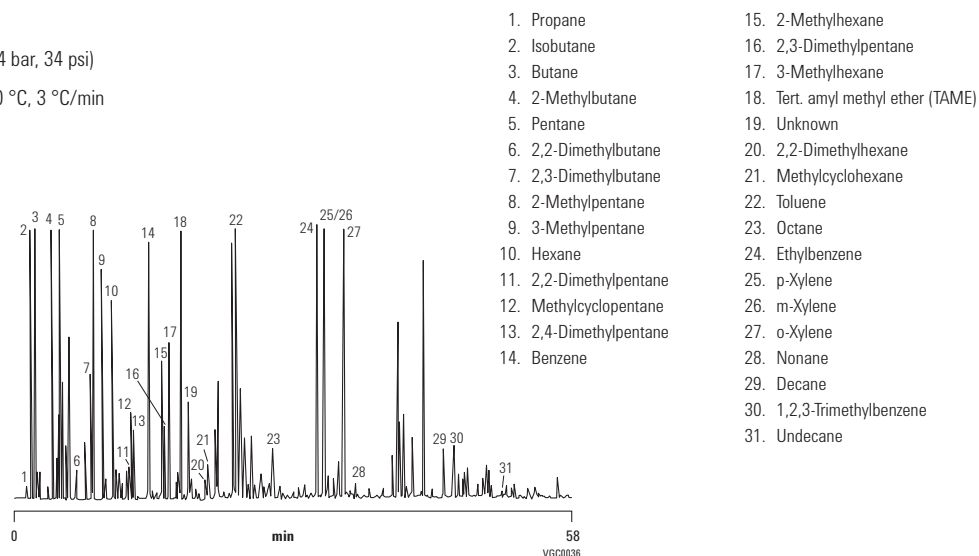**Polyethylene**

Column: DB-1
125-1011
15 m x 0.53 mm, 0.15 μ m

Carrier: Helium at 8 mL/min

Oven: 120-360 $^{\circ}$ C at 10 $^{\circ}$ C/min


Injection: Split ratio 1:500

Detector: FID, 300 $^{\circ}$ C

Nitrogen makeup gas at 30 mL/min

Sample: 0.5 μ L

3% solution in CS₂

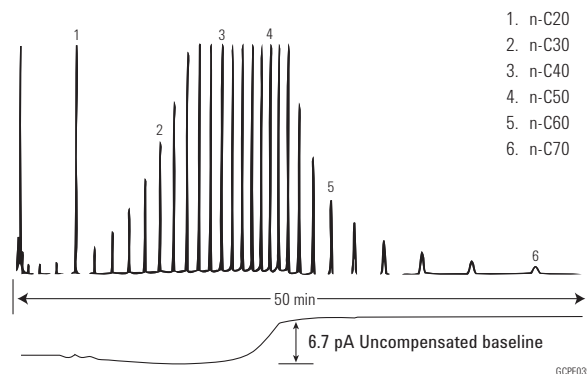


Diesel Analysis

Column: VF-5ht Fused Silica
CP9047
15 m x 0.32 mm, 0.10 µm

Carrier: H₂, 60 kPa, 0.6 bar, 8.6 psi

Oven: 50 °C (1 min), 15 °C to 180 °C,
7 °C to 230 °C, 30 °C to 380 °C

Detector: FID

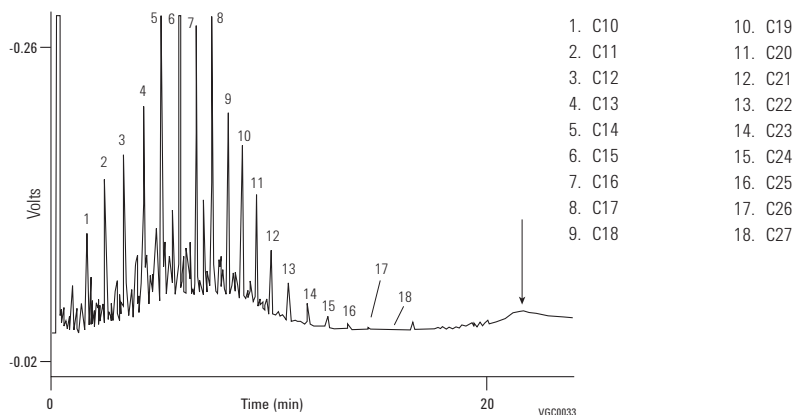**Analysis of Oxygenates
in a C1 to C5 Hydrocarbon Mix**

Column: Lowox
CP8587
10 m x 0.53 mm, 10.00 µm

Sample: 1 µL

Sample Conc: 0.01% per compound

Solvent: Cyclohexane


Carrier: He, 28.8 kPa (0.288 bar, 4.1 psi)

Oven: 50 °C (5 min) to 240 °C, 10 °C/min

Injection: Split, T=250 °C

Detector: FID, T=250 °C

- | | |
|-------------------------------|--|
| 1. Acetaldehyde | 11. Methanol |
| 2. Diethyl ether | 12. Acetone |
| 3. Ethyl tert-butyl ether | 13. Isovaleraldehyde |
| 4. Methyl tert-butyl ether | 14. Valeraldehyde |
| 5. Diisopropyl ether | 15. 2-Butanone |
| 6. Propionaldehyde (propanol) | 16. Ethanol |
| 7. Tert-amyl methyl ether | 17. 1-Propanol |
| 8. Dipropyl ether | 18. 2-Methyl-1-propanol (isobutanol) |
| 9. Isobutyraldehyde | 19. 2-Methyl-2-propanol (tert-butanol) |
| 10. Butyraldehyde | 20. 1-Butanol |

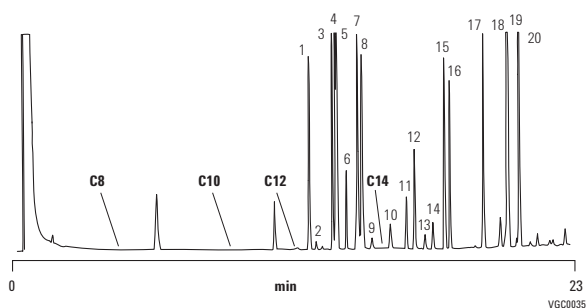

Analysis of Process Gas

Column: HP-PLLOT Q PT
19095P-Q04PT
30 m x 0.53 mm, 40.00 µm


Instrument: Agilent 7890A

Carrier: Hydrogen, constant flow mode, 40 cm/s, 32 °C

Oven: 32 °C for 5 min, 32 °C to 70 °C at 30 °C/min,
70 °C for 5 min, 70 to 160 °C at 10 °C/min

Injection: 170 °C, split ratio 5:1, 250 µL gas sampling loop

Detector: TCD at 250 °C

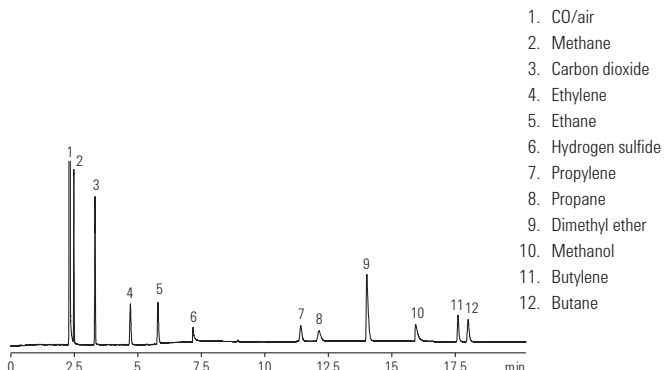**Detailed Hydrocarbon Analysis of Petroleum Naphthas Through N-nonane Using ASTM D5134**

Column: CP-Sil PONA for ASTM D5134
CP7531
50 m x 0.21 mm, 0.50 µm


Sample: 0.2 µL

Carrier: Helium

Oven: 35 °C (30 min) at 2 °C/min to 200 °C (10 min)

Injection: Split/splitless 1177, full EFC control,
250 °C, split 200 mL/min

Detector: FID, 250 °C

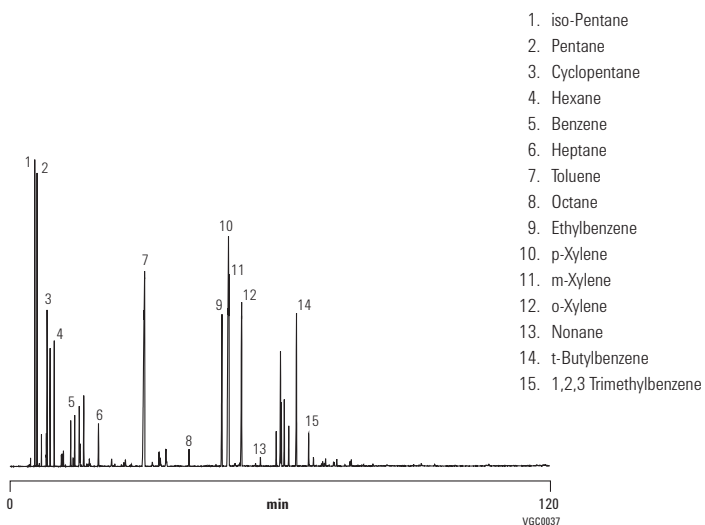

Industrial Chemical Applications

Alcohols I

Column: DB-624
125-1334
30 m x 0.53 mm, 3.00 μ m

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min
260 °C for 3 min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 μ L of 0.01-0.05% each solvent in CS₂

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal kit, 5188-5367

Syringe: 5 μ L tapered, FN 23-26s/42/HP, 5181-1273

- | | |
|--|--|
| 1. Methanol | 27. 2-Penten-1-ol |
| 2. Ethanol | 28. 3-Methyl-2-buten-1-ol |
| 3. Isopropanol | 29. Cyclopentanol |
| 4. tert-Butanol | 30. 3-Hexanol |
| 5. 2-Propen-1-ol (allyl alcohol) | 31. 2-Hexanol |
| 6. 1-Propanol | 32. 4-Hydroxy-4-methyl-2-pentanone |
| 7. 2-Propyn-1-ol (propargyl alcohol) | 33. Furfuryl alcohol |
| 8. sec-Butanol | 34. cis-3-Hexen-1-ol |
| 9. 2-Methyl-3-buten-2-ol | 35. 1-Hexanol |
| 10. Isobutanol | 36. cis-2-Hexen-1-ol |
| 11. 2-Methoxyethanol (methyl cellosolve) | 37. Cyclohexanol |
| 12. 3-Buten-1-ol | 38. 3-Heptanol |
| 13. 2-Methyl-2-butanol (tert-amyl alcohol) | 39. 2-Heptanol |
| 14. 1-Butanol | 40. 2-Butoxyethanol (butyl cellosolve) |
| 15. 2-Buten-1-ol (crotyl alcohol) | 41. cis-4-Hepten-1-ol |
| 16. Ethylene glycol | 42. trans-2-Hepten-1-ol |
| 17. 1-Penten-3-ol | 43. 1-Heptanol |
| 18. 2-Pentanol | 44. Benzyl alcohol |
| 19. Glycidol | 45. 2-Ethyl-1-hexanol |
| 20. 3-Pentanol | 46. α -Methylphenyl alcohol |
| 21. 2-Ethoxyethanol (cellosolve) | 47. 1-Octanol |
| 22. Propylene glycol | 48. 1-Nonanol |
| 23. 3-Methyl-1-butanol (isoamyl alcohol) | 49. 2-Phenoxyethanol |
| 24. 2-Methyl-1-butanol (active amyl alcohol) | 50. α -Ethylphenethyl alcohol |
| 25. 4-Methyl-2-pentanol | 51. β -Ethylphenethyl alcohol |
| 26. 1-Pentanol | 52. 1-Decanol |

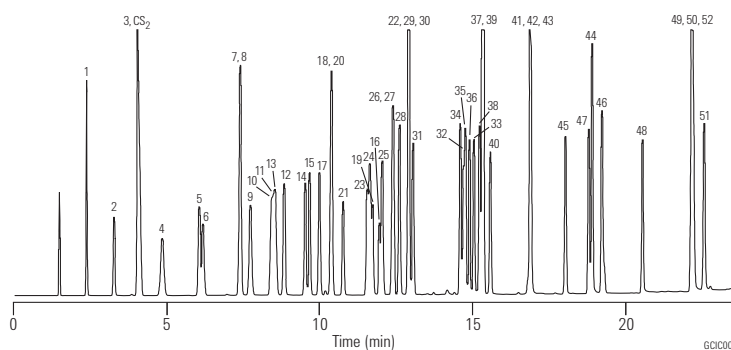

Halogenated Hydrocarbons I

Column: DB-624
123-1334
30 m x 0.32 mm, 1.80 µm

Carrier: Helium at 35 cm/s

Oven: 35 °C for 5 min
35-245 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal kit, 5188-5367

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

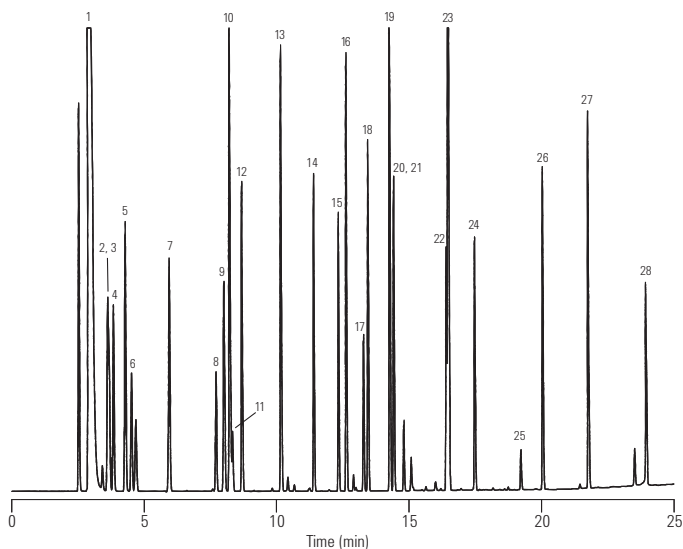

1. Pentane
2. Iodomethane
3. 1,1-Dichloroethene
4. 1,1,2-Trichlorotrifluoroethane (freon 113)
5. 3-Chloropropene (allyl chloride)
6. Methylene chloride
7. 1,1-Dichloroethane
8. Chloroform
9. 1,1,1-Trichloroethane
10. 1-Chlorobutane
11. Carbon tetrachloride
12. 1,2-Dichloroethane
13. 1,2-Dichloropropane
14. cis-1,2-Dichloropropene
15. trans-1,2-Dichloropropene
16. 1,1,2-Trichloroethane
17. 1,1,1,2-Tetrachloroethane
18. 1,2-Dibromoethane (EDB)
19. 1-Chlorohexane
20. trans-1,4-Dichloro-2-butene
21. Iodoform
22. Hexachlorobutadiene
23. 1,2,3-Trichloropropane
24. 1,1,2,2-Tetrachloroethane
25. Pentachloroethane
26. 1,2-Dibromo-3-chloropropane (DBCP)
27. Hexachloroethane
28. Hexachlorocyclopentadiene

GCIC034

Aromatic Solvents

Column: DB-200
122-2032
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 31 cm/s

Oven: 50 °C for 5 min
50-160 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.5 µL of 0.5 µg/µL
standard in hexane


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal kit, 5188-5367

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

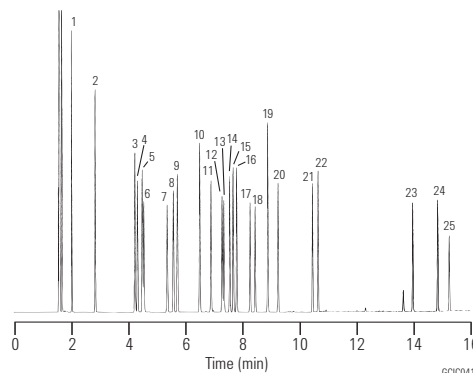

- | | |
|---------------------|----------------------------|
| 1. Benzene | 14. tert-Butylbenzene |
| 2. Toluene | 15. sec-Butylbenzene |
| 3. Ethylbenzene | 16. Isobutylbenzene |
| 4. Chlorobenzene | 17. 1,3-Dichlorobenzene |
| 5. p-Xylene | 18. 1,4-Dichlorobenzene |
| 6. m-Xylene | 19. n-Butylbenzene |
| 7. o-Xylene | 20. 1,2-Dichlorobenzene |
| 8. Styrene | 21. 1,3-Diisopropylbenzene |
| 9. Isopropylbenzene | 22. 1,4-Diisopropylbenzene |
| 10. n-Propylbenzene | 23. 2-Nitrotoluene |
| 11. 2-Chlorotoluene | 24. 3-Nitrotoluene |
| 12. 3-Chlorotoluene | 25. 4-Nitrotoluene |
| 13. 4-Chlorotoluene | |

Phenols I

Column: HP-5ms
19091S-433
30 m x 0.25 mm, 0.25 µm

Carrier: Helium, 33 cm/s, constant flow

Oven: 35 °C for 5 min
35-220 °C at 8 °C/min

Injection: Splitless, 250 °C

Detector: FID, 300 °C

Sample: 1 µL
20 µg/mL phenols in methylene chloride


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Seal: Gold plated seal kit, 5188-5367

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

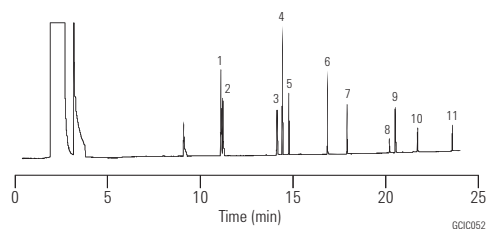

- | |
|--------------------------------|
| 1. Phenol |
| 2. 2-Chlorophenol |
| 3. 2-Nitrophenol |
| 4. 2,4-Dimethylphenol |
| 5. 2,4-Dichlorophenol |
| 6. 4-Chloro-3-methylphenol |
| 7. 2,4,6-Trinitrophenol |
| 8. 2,4-Dinitrophenol |
| 9. 4-Nitrophenol |
| 10. 2-Methyl-4,6-dinitrophenol |
| 11. Pentachlorophenol |

Inorganic Gases

Column: GS-GasPro
113-4332
30 m x 0.32 mm


Carrier: Helium at 53 cm/s

Oven: 25 °C for 3 min
25-200 °C at 10 °C/min
200 °C hold

Injection: Split, 200 °C
Split ratio 1:50

Detector: TCD, 250 °C

Sample: 50 µL

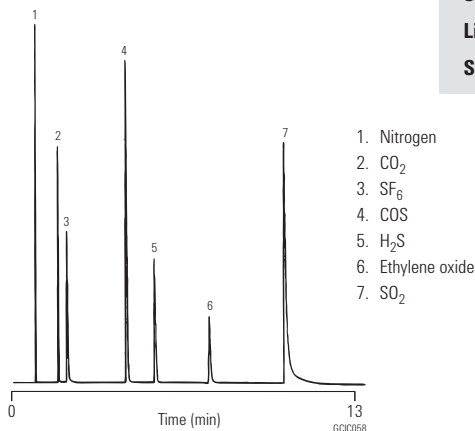**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal kit, 5188-5367

Alcohols II

Column: DB-WAXetr
123-7354
50 m x 0.32 mm, 1.00 µm

Carrier: Helium at 50 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-230 °C at 10 °C/min
230 °C for 5 min

Injection: Split, 250 °C
Split ratio 1:5


Detector: FID, 250 °C
Nitrogen makeup gas
at 35 mL/min

Sample: 1 µL of 0.15%
each solvent in CS₂

1. Methanol
2. Ethanol
3. Isopropanol
4. tert-Butanol
5. 2-Propen-1-ol (allyl alcohol)
6. 1-Propanol
7. 2-Propyn-1-ol (propargyl alcohol)
8. sec-Butanol
9. 2-Methyl-3-buten-2-ol
10. Isobutanol
11. 2-Methoxyethanol (methyl cellosolve)
12. 3-Buten-1-ol

13. 2-Methyl-2-butanol (tert-amyl alcohol)
14. 1-Butanol
15. 2-Buten-1-ol (crotyl alcohol)
16. Ethylene glycol
17. 1-Penten-3-ol
18. 2-Pentanol
19. Glycidol
20. 3-Pentanol
21. 2-Ethoxyethanol (cellosolve)
22. Propylene glycol
23. 3-Methyl-1-butanol (isoamyl alcohol)
24. 2-Methyl-1-butanol (active amyl alcohol)
25. 4-Methyl-2-pentanol
26. 1-Pentanol
27. 2-Penten-1-ol
28. 3-Methyl-2-buten-1-ol
29. Cyclopentanol
30. 3-Hexanol
31. 2-Hexanol
32. 4-Hydroxy-4-methyl-2-pentanone

33. Furfuryl alcohol
34. cis-3-Hexen-1-ol
35. cis-2-Hexen-1-ol
36. Cyclohexanol
37. 3-Heptanol
38. 2-Heptanol
39. 2-Butoxyethanol (butyl cellosolve)
40. cis-4-Hepten-1-ol
41. trans-2-Hepten-1-ol
42. 1-Heptanol
43. Benzyl alcohol
44. 2-Ethyl-1-hexanol
45. 1-Octanol
46. 1-Nonanol
47. 2-Phenoxyethanol
48. a-Ethylphenethyl alcohol
49. b-Ethylphenethyl alcohol
50. 1-Decanol

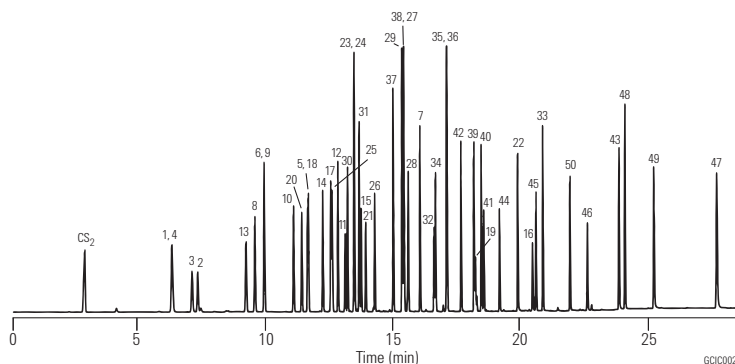

Alcohols III

Column: HP-INNOWax
19095N-123
30 m x 0.53 mm, 1.00 µm


Carrier: Helium, 29 cm/s, 3.0 psi (45 °C)

Oven: 45 °C for 1 min
45-150 °C at 10 °C/min
4 mL/min constant flow

Injection: Split, 250 °C
Split ratio 25:1

Detector: FID, 250 °C

Sample: 1 µL

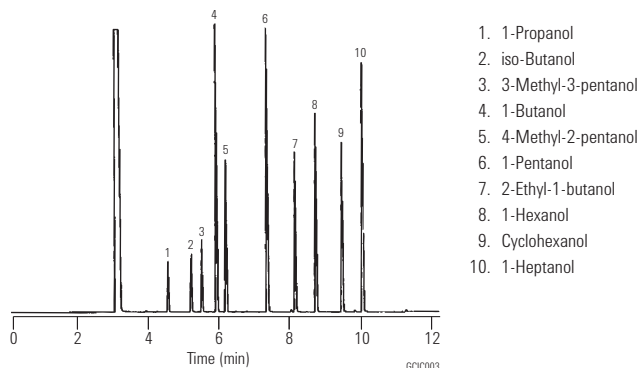

Analysis of Amino Alcohols in Water

Column: CP-Sil 5 CB
CP7640
50 m x 0.53 mm, 2.00 µm

Sample: 0.2 µL

Sample Conc: 1 ppm


Solvent: Water

Carrier: He, 0.7 mL/min, 70 kPa (0.7 bar, 9 psi)

Oven: 65 °C to 100 °C, 10 °C/min

Injection: Splitless

Detector: MS

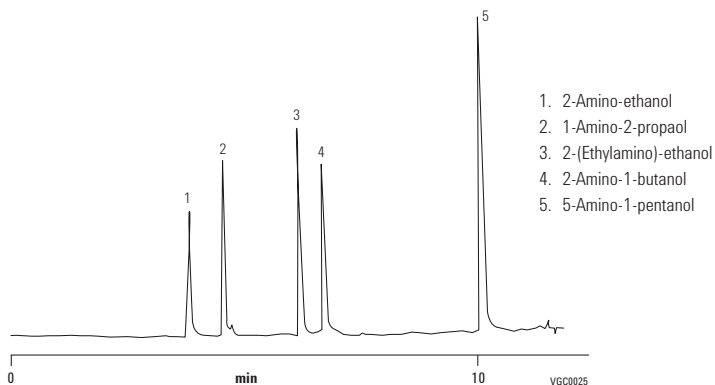

*Courtesy of Victor Berezkin and Aleksey B. Lapin,
Institute of Petrochemical Synthesis, Russian Academy of Science, Moscow, Russia*

Amines and Alcohols

Column: CP-Volamine
CP7446
15 m x 0.32 mm

Sample: 0.5 µL

Sample Conc: 1000 ppm, approx. 5 ng per component
on the column


Solvent: Methanol

Carrier: Helium, 50 kPa, 55 cm/s

Oven: 35 °C (0.5 min) to 240 °C, 30 °C/min

Injection: Split

Detector: MS

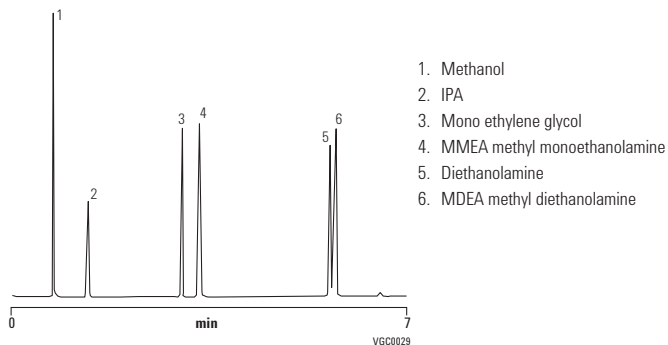

Courtesy of J. Luong, Dow Chemical Canada

Analysis of Ethanolamines

Column: CP-Sil 8 CB for Amines
CP7596
30 m x 0.32 mm, 1.00 µm

Sample Conc: 5-10 ng per component on the column


Solvent: Methanol

Carrier: Helium, 50 kPa (0.5 bar, 7 psi)

Oven: 60 °C (5 min) to 220 °C, 6 °C/min

Injection: Split

Detector: FID

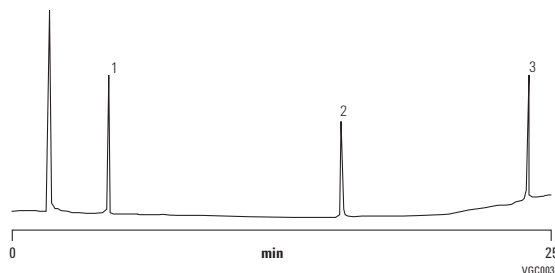

1. MEA (mono-ethanolamine)
2. DEA (di-ethanolamine)
3. TEA (tri-ethanolamine)

Ethoxyethanol

Column: HP-FFAP
19095F-123
30 m x 0.53 mm, 1.00 µm

Carrier: Helium, 10 mL/min

Oven: 60 °C for 1 min
60-100 °C at 5 °C/min
100-210 °C at 10 °C/min

Injection: Split ratio 10:1

Detector: TCD


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

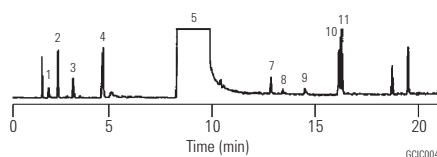

- | | |
|--------------------------|---------------------------------|
| 1. Ethylene oxide | 7. Hydroxy acetate |
| 2. Ethyl formate | 8. Acetic acid |
| 3. Ethyl alcohol | 9. Formic acid |
| 4. Water | 10. Ethylene glycol/monoformate |
| 5. 2-Ethoxyethanol | 11. Ethylene glycol/monoacetate |
| 6. 2-Ethoxyethyl acetate | |

Organic Acids

Column: DB-WAXetr
125-7332
30 m x 0.53 mm, 1.00 µm

Carrier: Helium at 37 cm/s,
measured at 40 °C

Oven: 125 °C for 5 min
125-180 °C at 15 °C/min
180 °C for 12 min

Injection: Split, 250 °C

Detector: FID, 250 °C


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

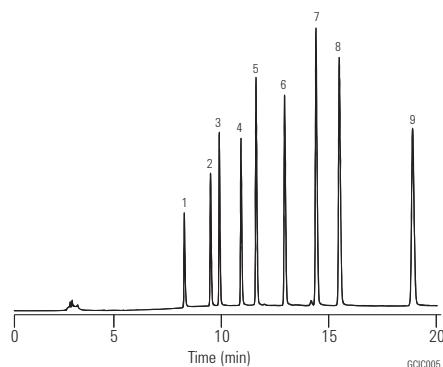

- | | |
|--------------------|----------------------------------|
| 1. Acetic acid | 6. Valeric acid (pentanoic acid) |
| 2. Propionic acid | 7. Isocaproic acid |
| 3. Isobutyric acid | 8. Caproic acid (hexanoic acid) |
| 4. Butyric acid | 9. Heptanoic acid |
| 5. Isovaleric acid | |

Free Organic Acids/C₄-C₅ Isomers

Column: HP-INNOWax
19091N-133
30 m x 0.25 mm, 0.25 µm


Carrier: Helium 42 cm/s, 24 psi (120 °C)
1.8 mL/min constant flow

Oven: 110 °C for 1 min
110-133 at 2 °C/min
133-160 °C at 3 °C/min

Injection: Split, 250 °C
Split ratio 40:1

Detector: FID, 300 °C

Sample: 1 µL

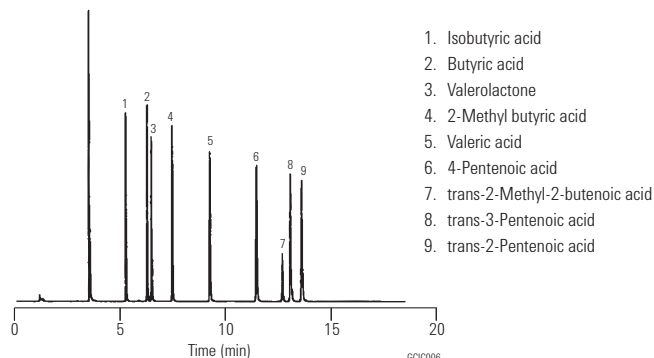

Volatile Amines


Column: DB-1
125-1035
30 m x 0.53 mm, 5.00 µm

Oven: 30 °C isothermal

Sampler: Headspace

Injection: Split ratio 1:10

Detector: FID
Nitrogen makeup gas at 30 mL/min

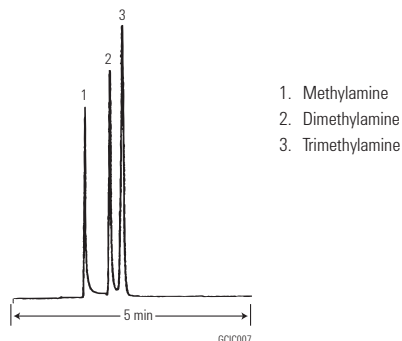

Trace Active Amines, 10 ng on-column

Column: HP-5ms
19091S-213
30 m x 0.32 mm, 1.00 µm


Carrier: Helium, constant pressure 9.79 psi

Oven: 75 °C for 0.5 min
75-250 °C at 10 °C/min
250-320 °C at 25 °C/min
320 °C for 5 min


Injection: On-column
Oven tracking mode

Detector: FID, 300 °C

Sample: 0.5 µL of each standard in methanol

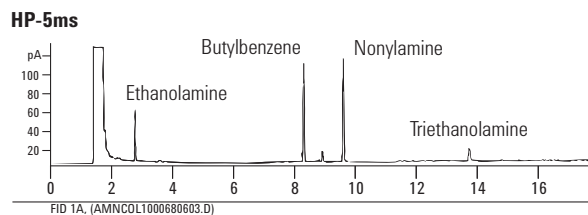

non-Agilent 5% phenyl amines column

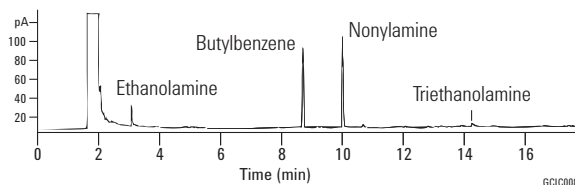

Primary Amines


Column: CAM
112-2132
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 40 cm/s

Oven: 110 °C isothermal

Injection: Split

Detector: FID
 Nitrogen makeup gas at 30 mL/min

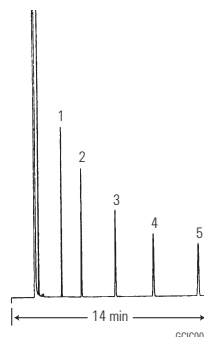

1. n-Octylamine
2. n-Nonylamine
3. n-Decylamine
4. Benzylamine
5. Dicyclohexylamine

Polyethyleneamines

Column: DB-5ms
122-5536
30 m x 0.25 mm, 0.50 µm


Carrier: Helium at 30 cm/s, measured at 100 °C

Oven: 100 °C for 1 min
 100-320 °C at 20 °C/min

Injection: Split, 250 °C
 Split ratio 1:50

Detector: FID, 300 °C
 Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 100 ng/µL standard in methanol

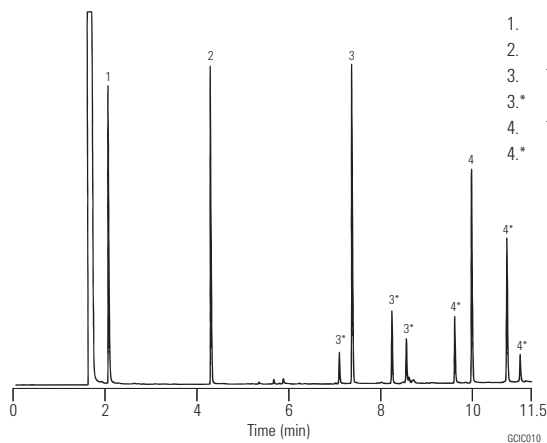

1. Ethylenediamine
2. Diethylenetriamine
3. Triethylenetetramine
- 3.* Branched and piperazine analogs of peak 3
4. Tetraethylenepentamine
- 4.* Branched and piperazine analogs of peak 4

Amines and Nitriles

Column: DB-5ms
122-5536
30 m x 0.25 mm, 0.50 μ m

Carrier: Helium at 22 cm/s, measured at 40 °C

Oven: 40 °C for 1 min
40-260 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 μ L of 100 ng/ μ L standard in methanol


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 μ L tapered, FN 23-26s/42/HP, 5181-1273

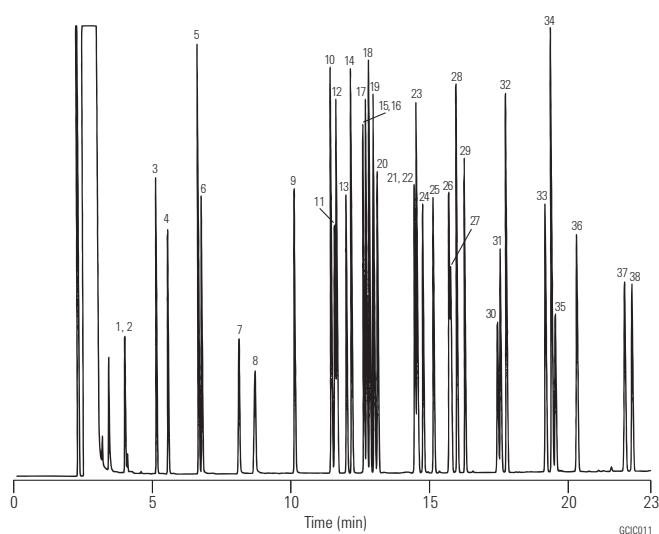

- | | |
|-----------------------------|--------------------------|
| 1. Diethylamine | 20. 2-Cyanopyridine |
| 2. Propionitrile | 21. 2-Chloroaniline |
| 3. Diisopropylamine | 22. n-Nonylamine |
| 4. Triethylamine | 23. 2,4-Dimethylaniline |
| 5. Pyridine | 24. 4-Chlorobenzonitrile |
| 6. Pyrimidine | 25. 2,6-Dimethylaniline |
| 7. Pyrazole | 26. 3-Chloroaniline |
| 8. Acrylamide | 27. 4-Chloroaniline |
| 9. Pyridazine | 28. N,N-Diethylaniline |
| 10. Aniline | 29. n-Decylamine |
| 11. 3-Bromopyridine | 30. 4-Bromoaniline |
| 12. Benzonitrile | 31. 3,4-Diaminotoluene |
| 13. 3-Cyanopyridine | 32. 2,6-Diethylaniline |
| 14. Benzylamine | 33. 2-Nitroaniline |
| 15. n-Octylamine | 34. Dicyclohexylamine |
| 16. 1-Methyl-2-pyrrolidine | 35. 3,4-Dichloroaniline |
| 17. N,N-Dimethylbenzylamine | 36. 3-Nitroaniline |
| 18. Phenylethylamine | 37. 4-Nitroaniline |
| 19. N-Benzylmethylamine | 38. Diphenylaniline |

Amines in Water


Column: CAM
112-2132
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 38 cm/s

Oven: 120-220 °C at 10 °C/min

Injection: Split

Detector: FID
Nitrogen makeup gas at 30 mL/min

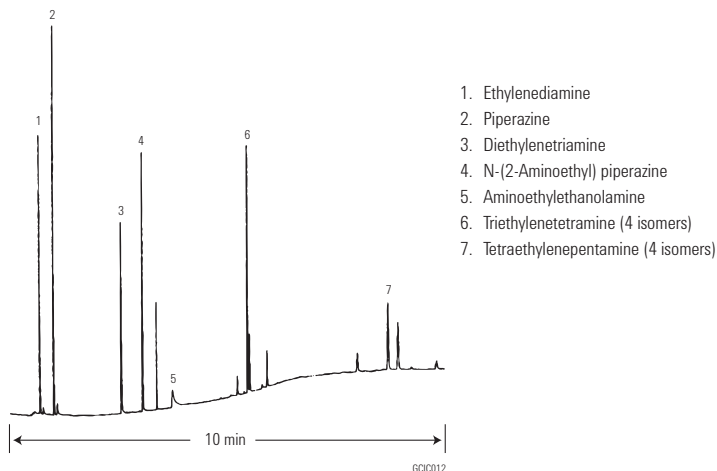**Aldehydes and Acids**

Column: HP-INNOWax
19091N-213
30 m x 0.32 mm, 0.50 µm

Carrier: Helium, 40 cm/s, 11.7 psi (60 °C)

Oven: 60 °C for 1 min
60-250 °C at 10 °C/min
2.5 mL/min constant flow

Injection: Split, 250 °C
Split ratio 40:1

Detector: FID, 275 °C

Sample: 0.5 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

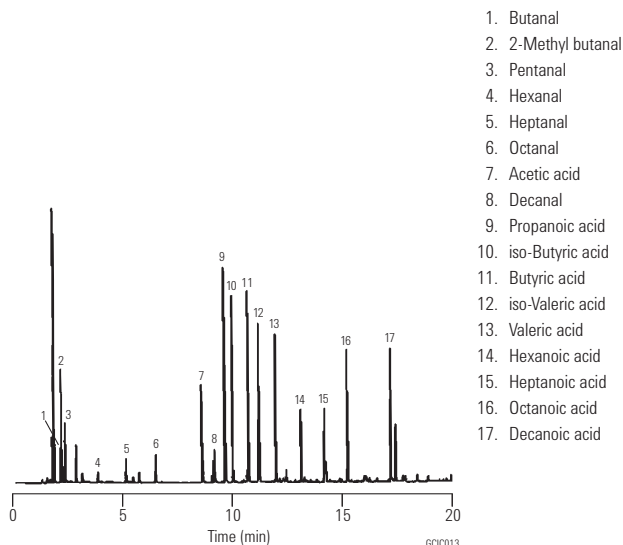

Aldehydes and Ketones

Column: DB-1
123-1034
30 m x 0.32 mm, 3.00 µm

Column: DB-WAX
123-7033
30 m x 0.32 mm, 0.50 µm

Carrier: Helium at 32 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-210 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273


DB-1

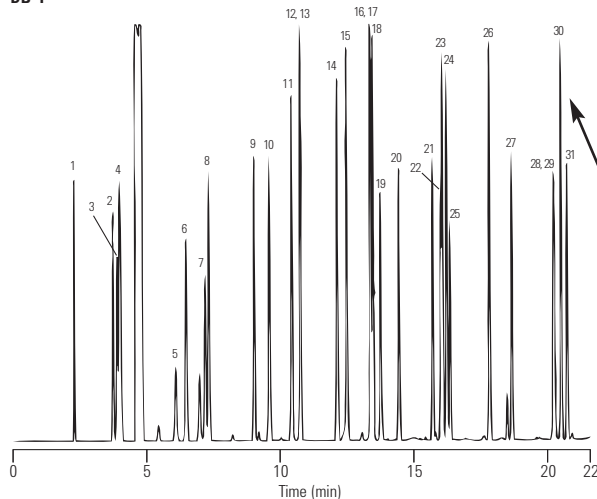

- | | |
|---------------------------------|--------------------|
| 1. Acetaldehyde | 17. Cyclopentanone |
| 2. Acrolein | 18. 2-Hexanone |
| 3. Acetone | 19. Hexanal |
| 4. Propionaldehyde | 20. Furfural |
| 5. Isobutyraldehyde | 21. 4-Heptanone |
| 6. Methacrolein | 22. 3-Heptanone |
| 7. Butyraldehyde | 23. 2-Heptanone |
| 8. 2-Butanone (MEK) | 24. Cyclohexanone |
| 9. Crotonaldehyde | 25. Heptanal |
| 10. 3-Methyl-2-butanone | 26. Benzaldehyde |
| 11. 2-Pentanone | 27. Octyl aldehyde |
| 12. 3-Pentanone | 28. o-Tolualdehyde |
| 13. Valeraldehyde (pentanal) | 29. m-Tolualdehyde |
| 14. 4-Methyl-2-pentanone (MIBK) | 30. p-Tolualdehyde |
| 15. 2-Methyl-3-pentanone | 31. Nonyl aldehyde |
| 16. 3-Hexanone | |

DB-1 provides the best overall resolution;
however, DB-WAX provides better resolution
of o- and m-tolualdehyde.

DB-WAX

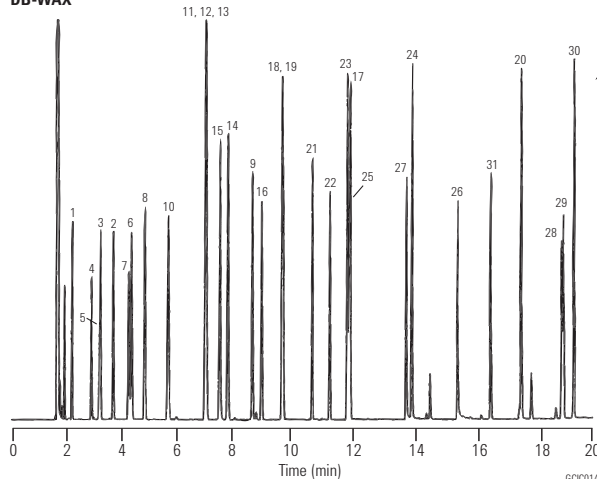

Formaldehyde Underivatized


Column: DB-WAX
123-7033
30 m x 0.32 mm, 0.50 µm

Carrier: Helium at 36 cm/s,
measured at 35 °C

Oven: 35 °C isothermal

Injection: Split, 200 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

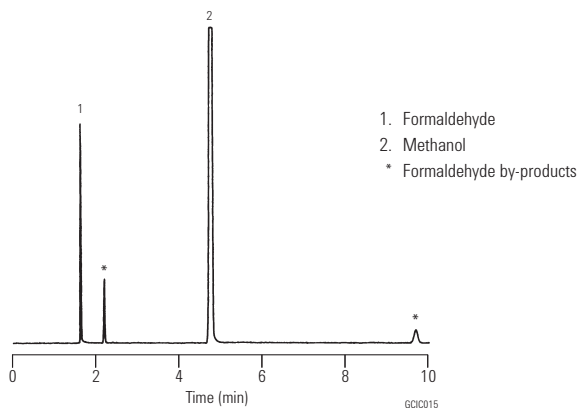**Formaldehyde-DNPH Derivative**

Column: DB-1
123-1012
15 m x 0.32 mm, 0.25 µm

Carrier: Helium at 35 cm/s,
measured at 150 °C

Oven: 150-250 °C at 20 °C/min

Injection: Split, 300 °C
Split ratio 1:100

Detector: ECD, 375 °C
Nitrogen makeup gas at 35 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

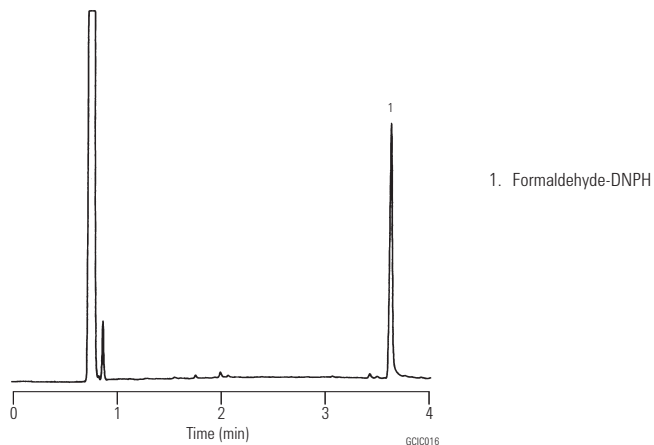

PFBHA Derivative


Column: DB-1
123-1012
15 m x 0.32 mm, 0.25 µm

Carrier: Helium at 40 cm/s,
measured at 60 °C

Oven: 60-100 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 375 °C
Nitrogen makeup gas at 35 mL/min

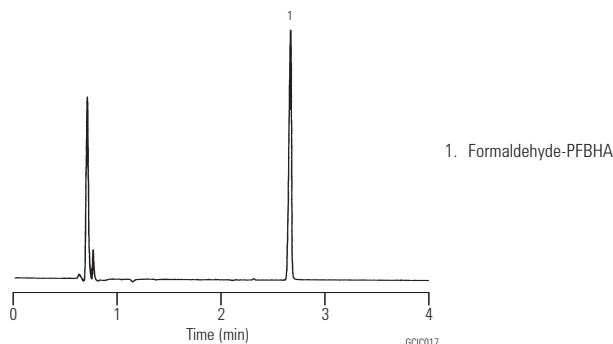

Aromatics I

Column: DB-1
125-1034
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

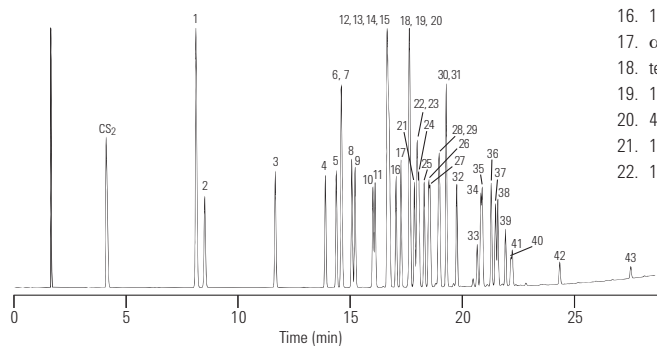

- | | |
|---|--|
| 1. Benzene | 23. Isobutylbenzene |
| 2. Fluorobenzene | 24. sec-Butylbenzene |
| 3. Toluene | 25. 1,2,3-Trimethylbenzene (hemimellitene) |
| 4. Chlorobenzene | 26. 1,2-Dichlorobenzene |
| 5. Ethylbenzene | 27. Iodobenzene |
| 6. m-Xylene | 28. Styrene oxide |
| 7. p-Xylene | 29. Butylbenzene |
| 8. Styrene | 30. 4-Chlorostyrene |
| 9. o-Xylene | 31. Nitrobenzene |
| 10. Isopropylbenzene (cumene) | 32. 4-tert-Butyltoluene |
| 11. Bromobenzene | 33. 1,3,5-Trichlorobenzene |
| 12. Propylbenzene | 34. 2-Nitrotoluene |
| 13. 2-Chlorotoluene | 35. 1,3-Diisopropylbenzene |
| 14. 3-Chlorotoluene | 36. 1,4-Diisopropylbenzene |
| 15. 4-Chlorotoluene | 37. 1,2,4-Trichlorobenzene |
| 16. 1,3,5-Trimethylbenzene (mesitylene) | 38. 3-Nitrotoluene |
| 17. α-Methylstyrene | 39. 4-Nitrotoluene |
| 18. tert-Butylbenzene | 40. 1,2,3-Trichlorobenzene |
| 19. 1,2,4-Trimethylbenzene (pseudocumene) | 41. 1-Chloro-4-nitrobenzene |
| 20. 4-Methylstyrene | 42. 1,2,4,5-Tetrachlorobenzene |
| 21. 1,3-Dichlorobenzene | 43. Pentachlorobenzene |
| 22. 1,4-Dichlorobenzene | |

Aromatics II

Column: DB-WAX
125-7032
30 m x 0.53 mm, 1.00 µm

Carrier: Helium at 30 cm/s, measured at 40 °C

Oven: 40 °C for 5 min
40-230 °C at 10 °C/min
230 °C for 7 min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

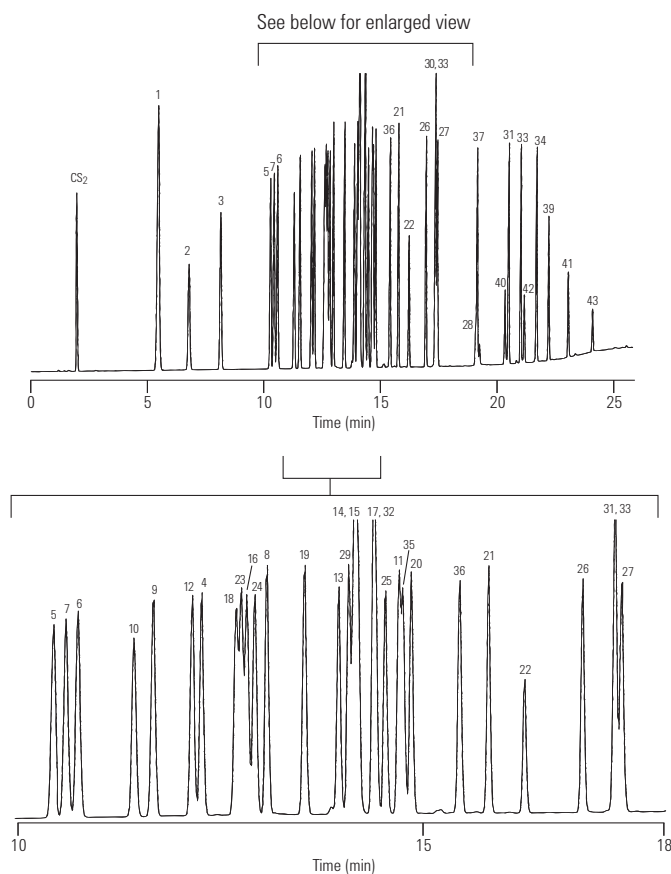

1. Benzene
2. Fluorobenzene
3. Toluene
4. Chlorobenzene
5. Ethylbenzene
6. m-Xylene
7. p-Xylene
8. Styrene
9. o-Xylene
10. Isopropylbenzene (cumene)
11. Bromobenzene
12. Propylbenzene
13. 2-Chlorotoluene
14. 3-Chlorotoluene
15. 4-Chlorotoluene
16. 1,3,5-Trimethylbenzene (mesitylene)
17. α-Methylstyrene
18. tert-Butylbenzene
19. 1,2,4-Trimethylbenzene (pseudocumene)
20. 4-Methylstyrene
21. 1,3-Dichlorobenzene
22. 1,4-Dichlorobenzene
23. Isobutylbenzene
24. sec-Butylbenzene
25. 1,2,3-Trimethylbenzene (hemimellitene)
26. 1,2-Dichlorobenzene
27. Iodobenzene
28. Styrene oxide (peak not shown)
29. Butylbenzene
30. 4-Chlorostyrene
31. Nitrobenzene
32. 4-tert-Butyltoluene
33. 1,3,5-Trichlorobenzene
34. 2-Nitrotoluene
35. 1,3-Diisopropylbenzene
36. 1,4-Diisopropylbenzene
37. 1,2,4-Trichlorobenzene
38. 3-Nitrotoluene
39. 4-Nitrotoluene
40. 1,2,3-Trichlorobenzene
41. 1-Chloro-4-nitrobenzene
42. 1,2,4,5-Tetrachlorobenzene
43. Pentachlorobenzene

GIC019

Impurities in Styrene

Column: DB-WAXetr
123-7363
60 m x 0.32 mm, 0.50 µm

Carrier: Helium at 29.4 cm/s, measured at 70 °C

Oven: 80 °C isothermal

Injection: Split, 230 °C
Split ratio 1:150

Detector: FID, 240 °C


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

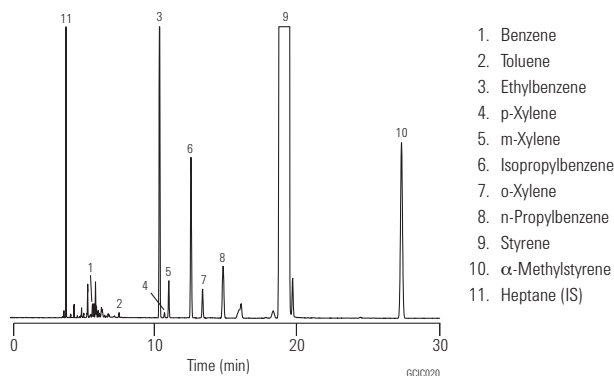**Impurities in Ethylbenzene**

Column: HP-INNOWax
19091N-216
60 m x 0.32 mm, 0.50 µm


Carrier: Helium, 32 cm/s, 19.9 psi (60 °C)
2.5 mL/min constant flow

Oven: 60 °C for 1 min
60-92 °C at 4 °C/min
92 °C for 4.5 min
92-220 °C at 20 °C/min

Injection: Split, 220 °C
Split ratio 100:1
ASTM Method D5060

Detector: FID, 270 °C

Sample: 0.5 µL
Neat, 99%+

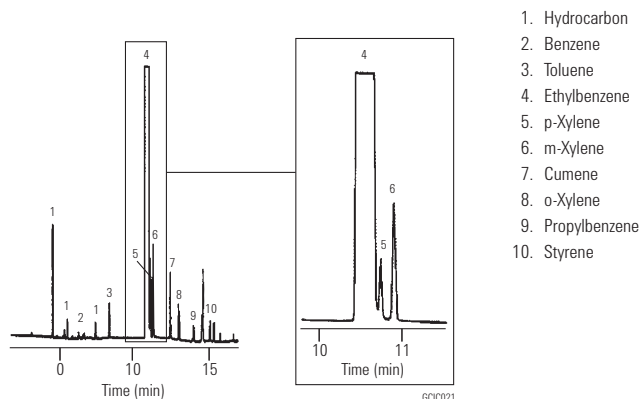

Pyrolysates of Polystyrene

Column: Ultra 1
19091A-105
50 m x 0.20 mm, 0.33 µm


Carrier: Helium, 30 psi, 12 mL/min

Oven: 0-280 at 5 °C/min

Injection: Split, 280 °C
Split ratio 30:1
Pyrolyzer 600 °C

Detector: FID, 300 °C

Sample: 100 mg pyrolyzed

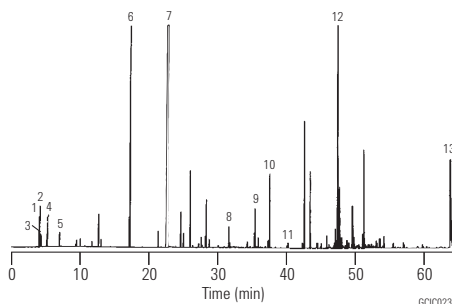

1. Propylene
2. Propane
3. 1-Butene
4. Butene
5. Pentane
6. Toluene
7. Styrene
8. $C_2H_5-C(Ph)=CH_2$
9. $C_4H_9-CH_2CH_2-Ph$
10. $C_4H_9-C(Ph)=CH_2$
11. $C_4H_9-CH=C(Ph)CH_3$
12. Styrene dimer
13. Styrene trimer

Esters I


Column: DB-1
125-1034
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

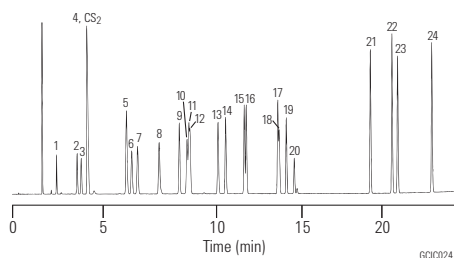**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711
Seal: Gold plated seal, 18740-20885
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

- | | |
|------------------------|---------------------------|
| 1. Methyl formate | 13. sec-Butyl acetate |
| 2. Ethyl formate | 14. Isobutyl acetate |
| 3. Methyl acetate | 15. Propyl propionate |
| 4. Vinyl acetate | 16. Butyl acetate |
| 5. Ethyl acetate | 17. Isoamyl acetate |
| 6. Propyl formate | 18. Amyl acetate |
| 7. Methyl propionate | 19. 2-Ethoxyethyl acetate |
| 8. Isopropyl acetate | 20. 2-Methylbutyl acetate |
| 9. Ethyl acrylate | 21. Methyl benzoate |
| 10. tert-Butyl acetate | 22. Benzyl acetate |
| 11. Ethyl propionate | 23. Ethyl benzoate |
| 12. Propyl acetate | 24. Propyl benzoate |

Esters II

Column: DB-624
125-1334
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min
260 °C for 3 min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

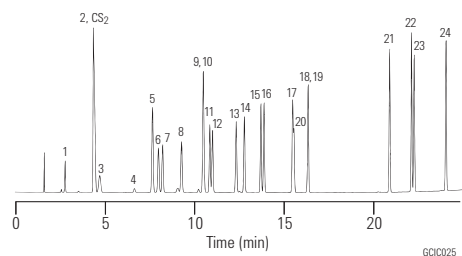

- | | |
|------------------------|---------------------------|
| 1. Methyl formate | 13. sec-Butyl acetate |
| 2. Ethyl formate | 14. Isobutyl acetate |
| 3. Methyl acetate | 15. Propyl propionate |
| 4. Vinyl acetate | 16. Butyl acetate |
| 5. Ethyl acetate | 17. Isoamyl acetate |
| 6. Propyl formate | 18. Amyl acetate |
| 7. Methyl propionate | 19. 2-Ethoxyethyl acetate |
| 8. Isopropyl acetate | 20. 2-Methylbutyl acetate |
| 9. Ethyl acrylate | 21. Methyl benzoate |
| 10. tert-Butyl acetate | 22. Benzyl acetate |
| 11. Ethyl propionate | 23. Ethyl benzoate |
| 12. Propyl acetate | 24. Propyl benzoate |

Esters III

Column: HP-INNOWax
19095N-123
30 m x 0.53 mm, 1.00 µm

Carrier: Helium 29 cm/s, 3.0 psi (45 °C)
4 mL/min constant flow

Oven: 45 °C for 1 min
45-200 °C at 5 °C/min

Injection: Split, 250 °C
Split ratio 25:1

Detector: FID, 250 °C

Sample: 1 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

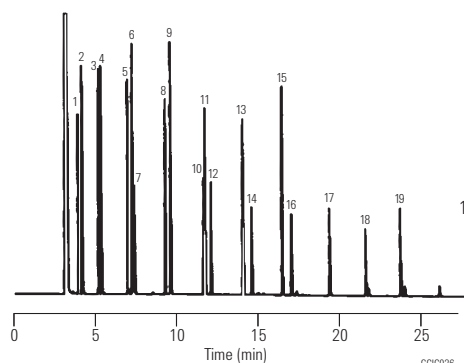

- | | |
|----------------------|---------------------------|
| 1. Ethyl propionate | 11. Propyl caproate |
| 2. Propyl acetate | 12. Methyl decanoate |
| 3. Ethyl butyrate | 13. Butyl caproate |
| 4. Propyl propionate | 14. Methyl dodecanoate |
| 5. Propyl butyrate | 15. Butyl heptanoate |
| 6. Ethyl valerate | 16. Methyl tetradecanoate |
| 7. Butyl propionate | 17. Methyl hexadecanoate |
| 8. Propyl valerate | 18. Methyl octadecanoate |
| 9. Ethyl caproate | 19. Methyl eicosenoate |
| 10. Butyl valerate | |

Ethers

Column: DB-624
125-1334
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min
260 °C for 3 min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

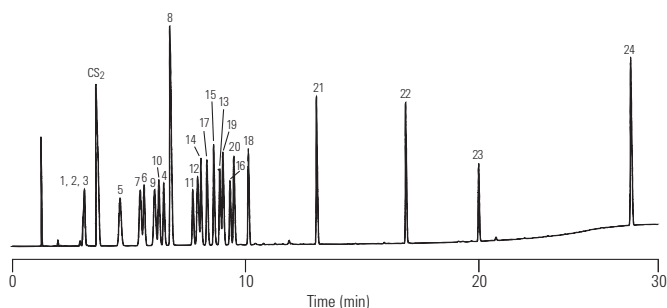

GCIC027

- | | |
|---|--|
| 1. Furan | 13. Diglyme (diethylene glycol dimethyl ether) |
| 2. Ethyl vinyl ether | 14. Propyl ether |
| 3. Ethyl ether | 15. Allyl ether |
| 4. 1,3-Dioxalane | 16. 1,4-Dioxane |
| 5. Methyl-tert-butyl ether (MTBE) | 17. Butyl ethyl ether |
| 6. Allyl ethyl ether | 18. Epichlorohydrin |
| 7. Isopropyl ether | 19. Tetrahydropyran |
| 8. Tetrahydrofuran (THF) | 20. Acetal (acetaldehyde diethyl acetal) |
| 9. tert-Amyl methyl ether | 21. Butyl ether |
| 10. Butyl methyl ether | 22. Pentyl ether |
| 11. Glyme (propylene glycol dimethyl ether) | 23. Triglyme (triethylene glycol dimethyl ether) |
| 12. tert-Amyl methyl ether | 24. Benzyl ether |

Glycols I

Column: DB-WAX
124-7032
30 m x 0.45 mm, 0.85 µm

Carrier: Helium at 35 cm/s,
measured at 50 °C

Oven: 50 °C for 2 min
50-220 °C at 10 °C/min

Injection: Megabore direct, 250 °C

Detector: FID, 280 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

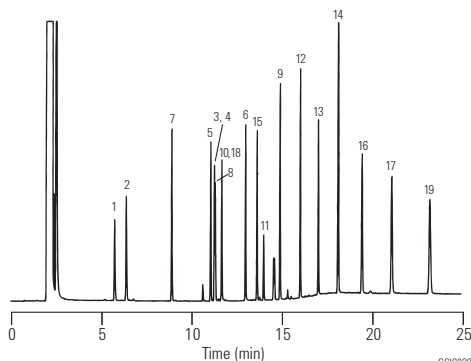

GCIC028

- | | |
|---------------------------------------|---------------------------------------|
| 1. Ethylene glycol monomethyl ether | 11. Dipropylene glycol |
| 2. Ethylene glycol monoethyl ether | 12. 1,5-Pentandiol |
| 3. 1,3-Propandiol | 13. 1,6-Hexandiol |
| 4. 1,2-Propandiol (propylene glycol) | 14. 1,7-Heptandiol |
| 5. 2,3-Butandiol | 15. Diethylene glycol monobutyl ether |
| 6. 1,3-Butandiol | 16. 1,8-Octandiol |
| 7. Ethylene glycol monobutyl ether | 17. 1,9-Nonandiol |
| 8. Diethylene glycol monomethyl ether | 18. Ethylene glycol |
| 9. 1,4-Butandiol | 19. 1,10-Decandiol |
| 10. Diethylene glycol monoethyl ether | |

Glycols II

Column: DB-624
125-1334
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min
260 °C for 3 min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, dual taper, deactivated, 4 mm id, G1544-80700

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

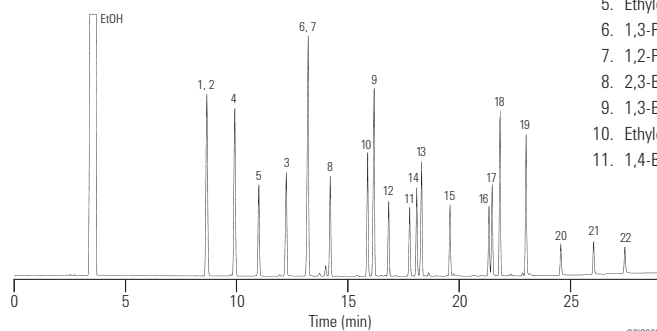

- | | |
|---------------------------------------|--|
| 1. Ethylene glycol monomethyl ether | 12. Diethylene glycol monomethyl ether |
| 2. Glyme | 13. Diethylene glycol |
| 3. Ethylene glycol | 14. Diethylene glycol monoethyl ether |
| 4. Diglyme | 15. 1,5-Pentanediol |
| 5. Ethylene glycol monoethyl ether | 16. 1,6-Hexanediol |
| 6. 1,3-Propanediol | 17. Diethylene glycol monobutyl ether |
| 7. 1,2-Propanediol (propylene glycol) | 18. Triglyme |
| 8. 2,3-Butanediol | 19. 1,7-Heptanediol |
| 9. 1,3-Butanediol | 20. 1,8-Octanediol |
| 10. Ethylene glycol monobutyl ether | 21. 1,9-Nonanediol |
| 11. 1,4-Butanediol | 22. 1,10-Decanediol |

Glycols III

Column: DB-1
124-1032
30 m x 0.45 mm, 1.27 µm

Carrier: Helium at 35 cm/s,
measured at 50 °C

Oven: 50 °C for 2 min
50-260 °C at 10 °C/min

Injection: Split, 250 °C

Detector: FID, 280 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, dual taper, deactivated, 4 mm id, G1544-80700

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

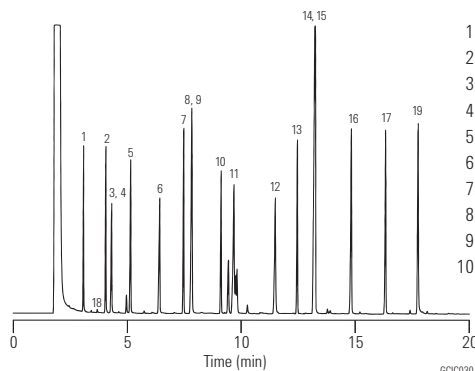

- | | |
|---------------------------------------|---------------------------------------|
| 1. Ethylene glycol monomethyl ether | 11. Dipropylene glycol |
| 2. Ethylene glycol monoethyl ether | 12. 1,5-Pentanediol |
| 3. 1,3-Propanediol | 13. 1,6-Hexanediol |
| 4. 1,2-Propanediol | 14. 1,7-Heptanediol |
| 5. 2,3-Butanediol | 15. Diethylene glycol monobutyl ether |
| 6. 1,3-Butanediol | 16. 1,8-Octanediol |
| 7. Ethylene glycol monobutyl ether | 17. 1,9-Nonanediol |
| 8. Diethylene glycol monomethyl ether | 18. Ethylene glycol |
| 9. 1,4-Butanediol | 19. 1,10-Decanediol |
| 10. Diethylene glycol monoethyl ether | |

Triethylene Glycol and Impurities

Column: DB-1
124-1032
30 m x 0.45 mm, 1.27 µm

Carrier: Helium at 35 cm/s,
measured at 50 °C

Oven: 170 °C isothermal

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 280 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.5 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

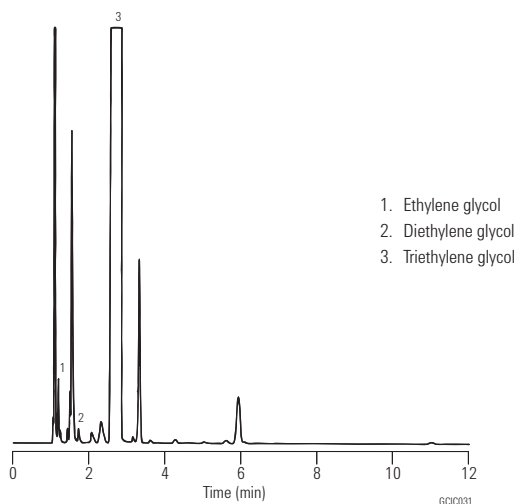**Ethylene Glycol Mixture**

Column: Ultra 1
19091A-101
12 m x 0.20 mm, 0.33 µm

Carrier: Helium, 25 cm/s

Oven: 100 °C for 0.5 min
100-200 °C at 20 °C/min

Injection: Split, 250 °C
Split ratio 100:1

Detector: FID

Sample: 1 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Liner, splitless, single-taper, glass wool, deactivated, 5062-3587

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

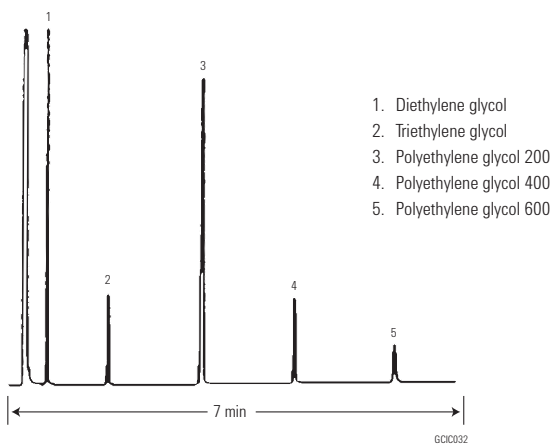

Glycols/Diols

Column: HP-1
19095Z-023
30 m x 0.53 mm, 0.88 µm


Carrier: Helium

Oven: 50 °C for 3 min
50-180 °C at 8 °C/min

Injection: On-column

Detector: FID, 250 °C

Sample: 1 µL

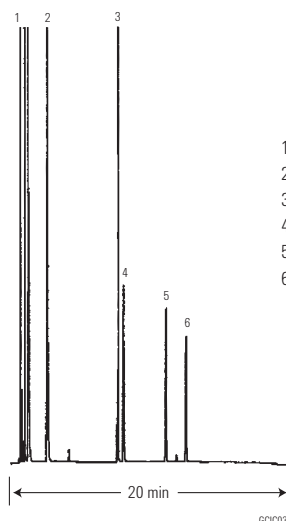

1. Ethylene glycol
2. 1,3-Butandiol
3. Ethylene glycol phenyl ether
4. 1,7-Hepatanediol
5. 1,9-Nonanediol
6. 1,10-Decanediol

Halogenated Hydrocarbons II

Column: DB-1
123-1034
30 m x 0.32 mm, 3.00 µm

Carrier: Helium at 35 cm/s, measured at 35 °C

Oven: 35 °C for 5 min
35-245 °C at 10 °C/min
245 °C for 2 min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: In pentane


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

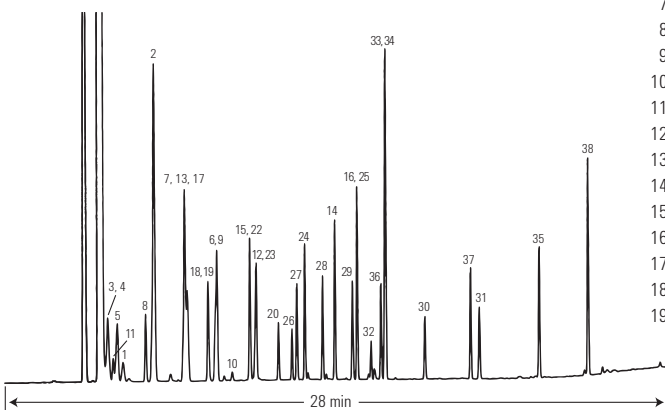

1. 1,1,2-Trichlorotrifluoroethane (freon 113)
2. 1,1-Dichloroethene
3. Bromoethane (ethyl bromide)
4. Iodomethane
5. 3-Chloropropene (allyl chloride)
6. 1-Chlorobutane
7. 2,2-Dichloropropane
8. trans-1,2-Dichloroethene
9. 1,1,1-Trichloroethane
10. Carbon tetrachloride
11. Methylene chloride
12. Trichloroethene
13. Chloroform
14. Tetrachloroethene
15. 1,2-Dichloropropane
16. 1-Chlorohexane
17. Bromochloromethane
18. 1,1-Dichloroethane
19. 1,2-Dichloroethane
20. Iodoform
21. cis-1,3-Dichloropropene
22. Dibromomethane
23. Bromodichloromethane
24. 1,3-Dichloropropane
25. 1,1-Dichloropropane
26. trans-1,3-Dichloropropene
27. 1,1,2-Trichloroethane
28. 1,2-Dibromoethane (EDB)
29. 1,1,1,2-Tetrachloroethane
30. Pentachloroethane
31. Hexachloroethane
32. Bromoform
33. trans-1,4-Dichloro-2-butene
34. 1,2,3-Trichloropropane
35. Hexachlorobutadiene
36. 1,1,2,2-Tetrachloroethane
37. 1,2-Dibromo-3-chloropropane (DBCP)
38. Hexachlorocyclopentadiene

Chlorinated Isooctane

Column: HP-INNOWax
19091N-136
60 m x 0.25 mm, 0.25 µm

Carrier: Helium, 33 cm/s, 35.7 psi (80 °C) 2 mL/min

Oven: 80 °C isothermal

Injection: Split, 250 °C
Split ratio 150:1

Detector: FID, 300 °C

Sample: Monochloro isomers, 0.5 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

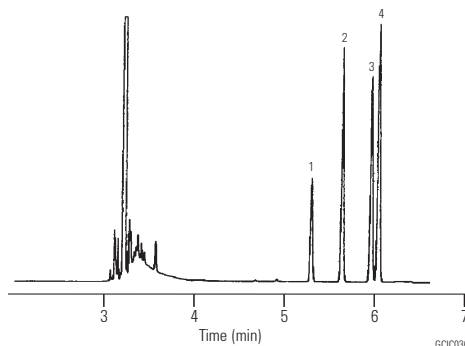

1. 1-Chloro isooctane
2. 4-Chloromethyl 2,2'-dimethyl pentane
3. 3-Chloro isooctane
4. 4-Chloro isooctane

Solvents I

Column: DB-WAXetr
125-7332
30 m x 0.53 mm, 1.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-140 °C at 5 °C/min

Injection: Split, 250 °C

Detector: FID, 250 °C


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

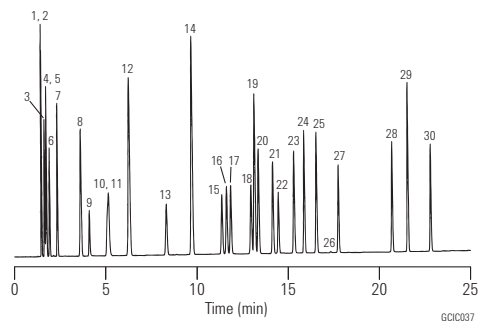

- | | |
|-----------------------------------|----------------------------|
| 1. 3-Methylpentane | 16. p-Xylene |
| 2. Hexane | 17. m-Xylene |
| 3. Isooctane | 18. Cumene |
| 4. Methyl-tert-butyl ether (MTBE) | 19. Dodecane |
| 5. Heptane | 20. o-Xylene |
| 6. Cyclohexane | 21. Propylbenzene |
| 7. Octane | 22. Chlorobenzene |
| 8. Nonane | 23. Mesitylene |
| 9. Methanol | 24. Styrene |
| 10. Ethanol | 25. 1,2,4-Trimethylbenzene |
| 11. Benzene | 26. Naphthalene |
| 12. Decane | 27. 4-Chlorotoluene |
| 13. Toluene | 28. 1,3-Dichlorobenzene |
| 14. Undecane | 29. 1,4-Dichlorobenzene |
| 15. Ethylbenzene | 30. 1,2-Dichlorobenzene |

Solvents II

Column: DB-WAXetr
123-7354
50 m x 0.32 mm, 1.00 µm

Carrier: Helium at 41 cm/s, measured at 50 °C

Oven: 50 °C for 5 min
50-170 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 280 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

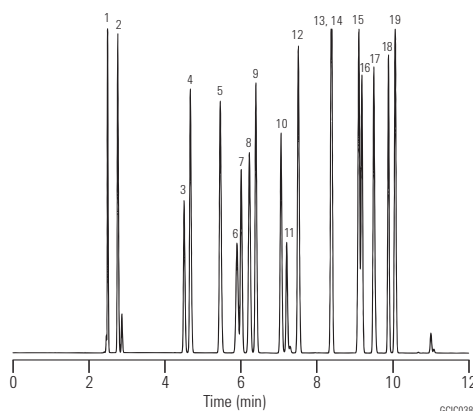

1. Hexane
2. Isooctane
3. Acetone
4. Ethyl formate
5. Tetrahydrofuran
6. Trichloroethane
7. Ethyl acetate
8. Isopropyl acetate
9. Methyl ethyl ketone
10. Isopropyl alcohol
11. Methylene chloride
12. Benzene
13. 2-Pentanone
14. Methyl isobutyl ketone
15. Isobutyl acetate
16. Chloroform
17. sec-Butyl alcohol
18. Toluene
19. n-Propanol

Solvents III

Column: DB-200
122-2033
30 m x 0.25 mm, 0.50 µm

Carrier: Helium at 31 cm/s

Oven: 45 °C for 7 min
45-145 °C at 20 °C/min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.5 µL of 0.5-1.0 µg/µL
standard in water


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

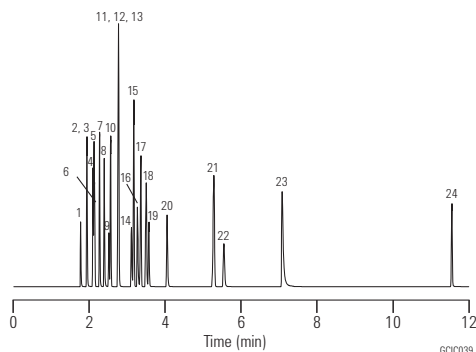

- | | |
|-----------------------|-------------------------------|
| 1. Methanol | 13. Acetone |
| 2. Ethanol | 14. Acetonitrile |
| 3. Ethyl ether | 15. Benzene |
| 4. Isopropanol | 16. Tetrahydrofuran (THF) |
| 5. n-Hexane | 17. Trichloroethylene |
| 6. Methylene chloride | 18. n-Butanol |
| 7. tert-Butanol | 19. Ethyl acetate |
| 8. n-Propanol | 20. Methyl ethyl ketone (MEK) |
| 9. Chloroform | 21. Toluene |
| 10. Cyclohexane | 22. 1,4-Dioxane |
| 11. sec-Butanol | 23. Pyridine |
| 12. n-Heptane | 24. Dimethylformamide (DMF) |

Solvents IV

Column: HP-1
19091Z-205
50 m x 0.20 mm, 0.50 µm

Carrier: Helium, 30 psi

Oven: 70-200 °C at 5 °C/min
200 °C for 2 min

Injection: Split

Detector: TCD

Sample: 1 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

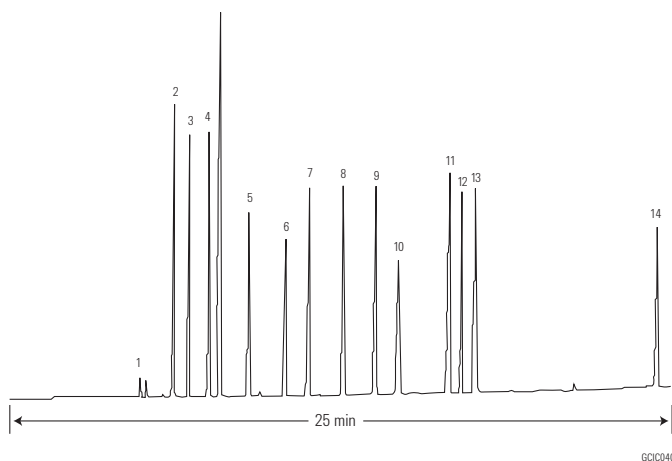

1. Isopropanol
2. Methyl ethyl ketone
3. Ethyl acetate
4. n-Butyl alcohol
5. Ethyl cellosolve
6. Methyl isobutyl ketone
7. Toluene
8. n-Butyl acetate
9. Diacetone alcohol
10. p-Xylene
11. Cellosolve acetate
12. o-Xylene
13. Butyl cellosolve
14. Butyl cellosolve acetate

Solvents

Column: PoraBOND Q PT
CP7348PT
25 m x 0.25 mm, 3.00 µm


Carrier: Helium, 1.5 mL/min

Oven: 90 °C to 140 °C at 10 °C/min
140 °C for 5 min
140 °C to 210 °C at 4 °C/min
210 °C for 6 min

Injection: Split, 250 °C, split ratio 1:150

Detector: MSD, 280 °C transfer line
Full scan at m/z 30-350

Sample: 1 µL

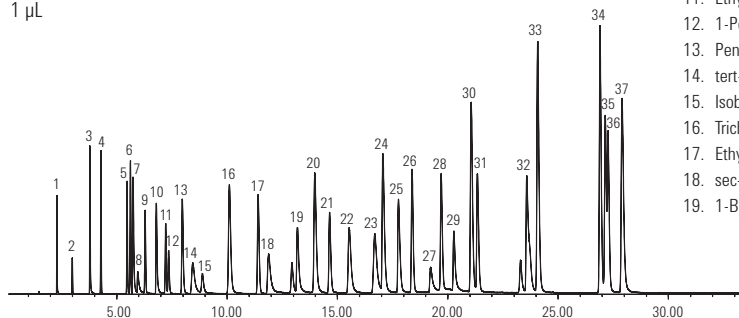

- | | |
|------------------------|----------------------------|
| 1. Methyl alcohol | 20. Benzene |
| 2. Acetaldehyde | 21. Hexane |
| 3. Ethanol | 22. 1,4-Dioxane |
| 4. Acetonitrile | 23. Ethyl tert-butyl ether |
| 5. Acetone | 24. Pyridine |
| 6. Methylene chloride | 25. N,N-dimethylformamide |
| 7. Isopropyl alcohol | 26. N-Propyl acetate |
| 8. 2-Propanamine | 27. 3-Methyl-1-butanol |
| 9. Ethyl formate | 28. n-Propyl ether |
| 10. 1-Propanol | 29. 1-Pentanol |
| 11. Ethyl ether | 30. Toluene |
| 12. 1-Pentene | 31. Heptane |
| 13. Pentane | 32. N,N-dimethylacetamide |
| 14. tert-Butyl alcohol | 33. Chlorobenzene |
| 15. Isobutyraldehyde | 34. Ethylbenzene |
| 16. Trichloromethane | 35. m-Xylene |
| 17. Ethyl acetate | 36. p-Xylene |
| 18. sec-Butyl alcohol | 37. o-Xylene |
| 19. 1-Butanol | |

Analysis of Solvents

Column: PoraBOND Q
CP7354
25 m x 0.53 mm, 10.00 µm

Sample: 5 µL

Sample Conc: 0.1% per compound


Solvent: DMSO

Carrier: He, 25 kPa (0.25 bar, 3.5 psi)

Oven: 100 °C (2 min) to 300 °C, 5 °C/min

Injection: Split, T=250 °C

Detector: FID, T=250 °C

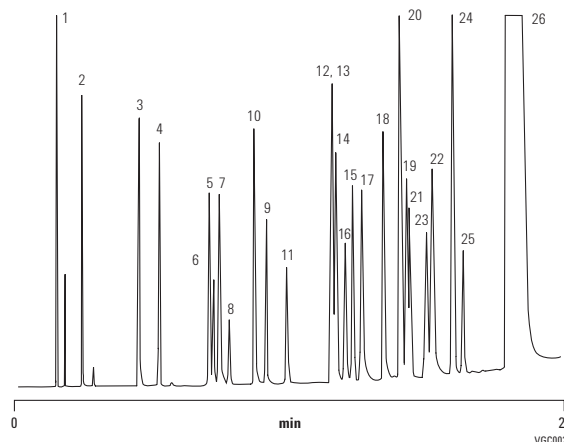

1. Methane
2. Methanol
3. Ethanol
4. Acetonitrile
5. Acetone
6. Dichloromethane
7. 2-Propanol
8. Dimethyl sulfide
9. Diethyl ether
10. 1-Propanol
11. Pentane
12. 2-Butanone
13. Trichloromethane
14. Tetrahydrofuran
15. Ethyl acetate
16. 2-Methoxyethanol
17. Isobutanol
18. Butanol
19. Hexane
20. Benzene
21. Trichloroethylene
22. Cyclohexane
23. 1,4-Dioxane
24. Pyridine
25. N,N-dimethylformamide
26. Dimethyl sulfoxide

Nitrogen-based Solvents I

Column: DB-1
125-1034
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

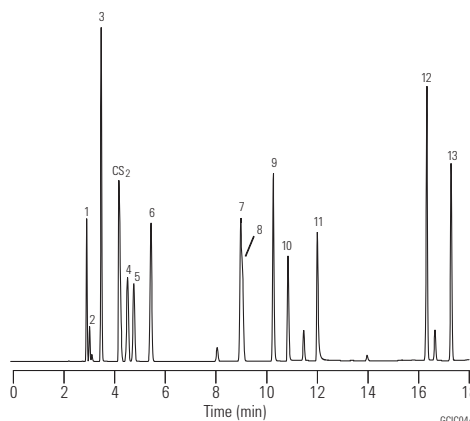

1. Acetonitrile
2. Acrolein
3. Acrylonitrile
4. Propionitrile
5. Methacrolein
6. Methacrylonitrile
7. Triethylamine
8. Ethyl acrylate
9. Pyridine
10. DMF (dimethylformamide)
11. DMSO (dimethyl sulfoxide)
12. Benzonitrile
13. 1-Methyl-2-pyrrolidinone

Nitrogen-based Solvents II

Column: DB-624
125-1334
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 30 cm/s,
measured at 40 °C

Oven: 40 °C for 5 min
40-260 °C at 10 °C/min
260 °C for 3 min

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

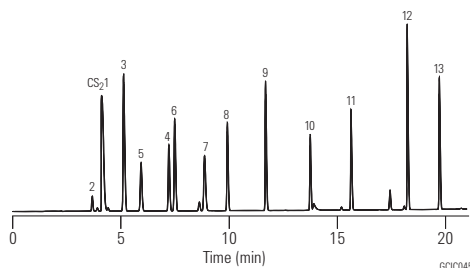

1. Acetonitrile
2. Acrolein
3. Acrylonitrile
4. Propionitrile
5. Methacrolein
6. Methacrylonitrile
7. Triethylamine
8. Ethyl acrylate
9. Pyridine
10. DMF (dimethylformamide)
11. DMSO (dimethyl sulfoxide)
12. Benzonitrile
13. 1-Methyl-2-pyrrolidinone

Acrylate Impurities I

Column: DB-200
125-2032
30 m x 0.53 mm, 1.00 µm

Carrier: Helium at 34.5 cm/s,
measured at 35 °C

Oven: 35 °C for 5 min,
35-200 °C at 10 °C/min

Injection: Split, 230 °C
Split ratio 1:10

Detector: FID, 250 °C


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

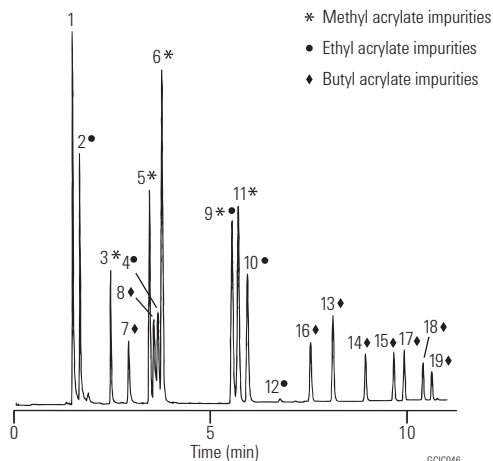

- * Methyl acrylate impurities
- Ethyl acrylate impurities
- ♦ Butyl acrylate impurities

1. Methanol
2. Ethanol
3. Methyl acetate
4. Ethyl acetate
5. Methyl acrylate
6. Methyl propionate
7. Isobutanol
8. Butanol
9. Ethyl acrylate
10. Ethyl propionate
11. Methyl methacrylate
12. Isopropyl acrylate
13. Isobutyl acetate
14. Butyl acetate
15. Isobutyl acrylate
16. Dibutyl ether
17. Isobutyl propionate
18. Butyl acrylate
19. Butyl propionate

Acrylate Impurities II

Column: DB-1701
125-0732
30 m x 0.53 mm, 1.00 µm

Carrier: Helium at 36.8 cm/s,
measured at 35 °C

Oven: 35 °C for 5 min,
35-200 °C at 10 °C/min

Injection: Split, 230 °C
Split ratio 1:10

Detector: FID, 250 °C


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

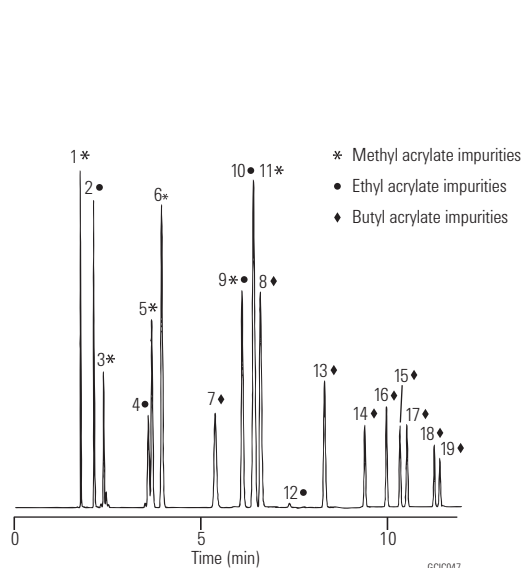

1. Methanol
2. Ethanol
3. Methyl acetate
4. Ethyl acetate
5. Methyl acrylate
6. Methyl propionate
7. Isobutanol
8. Butanol
9. Ethyl acrylate
10. Ethyl propionate
11. Methyl methacrylate
12. Isopropyl acrylate
13. Isobutyl acetate
14. Butyl acetate
15. Isobutyl acrylate
16. Dibutyl ether
17. Isobutyl propionate
18. Butyl acrylate
19. Butyl propionate

Acrylates

Column: HP-FFAP
19095F-121
10 m x 0.53 mm, 1.00 µm


Carrier: Hydrogen

Oven: 35 °C for 1 min
35-60 °C at 10 °C/min
60-160 °C at 15 °C/min

Injection: On-column

Detector: FID

Sample: 1 µL

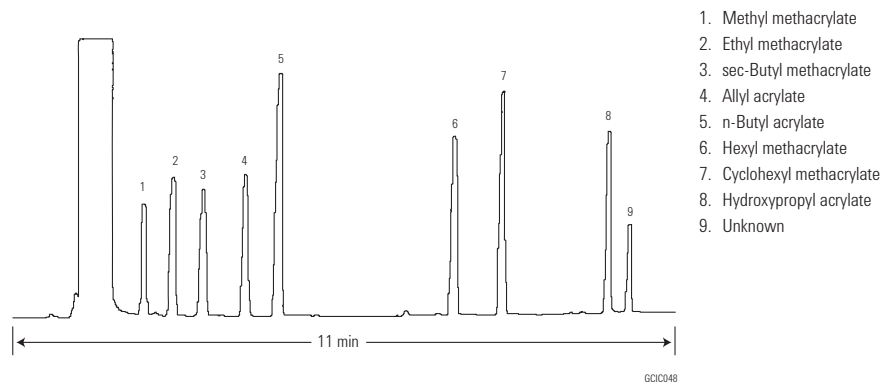

1. Methyl methacrylate
2. Ethyl methacrylate
3. sec-Butyl methacrylate
4. Allyl acrylate
5. n-Butyl acrylate
6. Hexyl methacrylate
7. Cyclohexyl methacrylate
8. Hydroxypropyl acrylate
9. Unknown

Anilines

Column: DB-35ms
128-3822
25 m x 0.20 mm, 0.33 µm

Carrier: Helium at 35 cm/s,
measured at 50 °C

Oven: 50 °C for 2 min
50-340 °C at 20 °C/min
340 °C for 10 min

Injection: Splitless, 280 °C
0.50 min purge activation time

Detector: FID, 320 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 5 ng
on-column per component


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

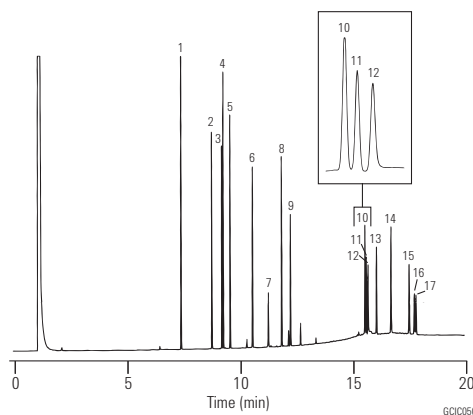

1. o-Toluidine
2. 4-Chloroaniline
3. 2-Methoxy-5-methylaniline
4. 2,4,5-Trimethylaniline
5. 4-Chloro-2-methylaniline
6. 2,4-Diaminotoluene
7. 2,4-Diaminoanisole
8. 2-Aminonaphthalene
9. 2-Methyl-5-nitroaniline
10. 4,4'-Oxydianiline
11. 4,4'-Methylenedianiline
12. Benzidine
13. 2-Aminoazotoluene
14. o-Tolidine
15. 4,4'-Thiodianiline
16. 3,3'-Dimethoxybenzidine
17. 3,3'-Dichlorobenzidine

Substituted Anilines

Column: DB-5ms
122-5536
30 m x 0.25 mm, 0.50 µm

Carrier: Helium at 33.3 cm/s,
measured at 150 °C

Oven: 40 °C for 5 min
40-290 °C at 12 °C/min
290 °C for 10 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: MSD, 325 °C transfer line

Sample: 1 µL of 25 ng/µL standard


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

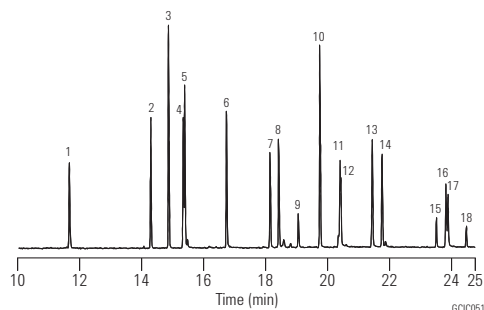

- | | m/z |
|---------------------------------|-----|
| 1. Aniline | 93 |
| 2. 2-Chloroaniline | 127 |
| 3. 2,6-Dimethylaniline | 121 |
| 4. 3-Chloroaniline | 127 |
| 5. 4-Chloroaniline | 127 |
| 6. 4-Bromoaniline | 171 |
| 7. 2-Nitroaniline | 138 |
| 8. 3,4-Dichloroaniline | 161 |
| 9. 3-Nitroaniline | 65 |
| 10. 2,4,5-Trichloroaniline | 195 |
| 11. 4-Chloro-2-nitroaniline | 172 |
| 12. 4-Nitroaniline | 138 |
| 13. 2-Chloro-4-nitroaniline | 172 |
| 14. 2,6-Dichloro-4-nitroaniline | 176 |
| 15. 2-Chloro-4,6-dinitroaniline | 217 |
| 16. 2,6-Dibromo-4-nitroaniline | 266 |
| 17. 2,4-Dinitroaniline | 183 |
| 18. 2-Bromo-4,6-dinitroaniline | 261 |

Phenols II

Column: DB-5ms
122-5536
30 m x 0.25 mm, 0.50 µm

Carrier: Helium at 22 cm/s,
measured at 100 °C

Oven: 100 °C for 1 min
100-270 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 50 ng/µL standard
in toluene/p-xylene


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

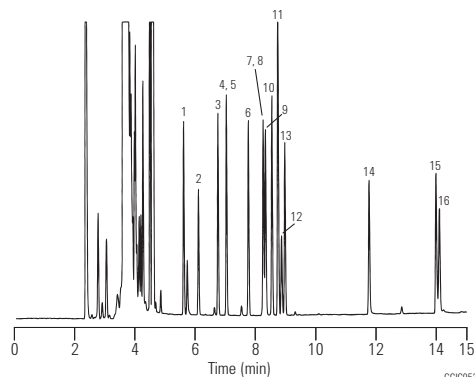

1. Phenol
2. 2-Chlorophenol
3. o-Cresol
4. m-Cresol
5. p-Cresol
6. 2,6-Xylenol
7. 2,4-Xylenol
8. 2,5-Xylenol
9. 2-Nitrophenol
10. 3,5-Xylenol
11. 2,3-Xylenol
12. 2,4-Dichlorophenol
13. 3,4-Xylenol
14. 2,4,6-Trichlorophenol
15. 2,4-Dinitrophenol
16. 1-Naphthol

Phenols III

Column: DB-WAX
122-7032
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 43 cm/s

Oven: 165 °C isothermal

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

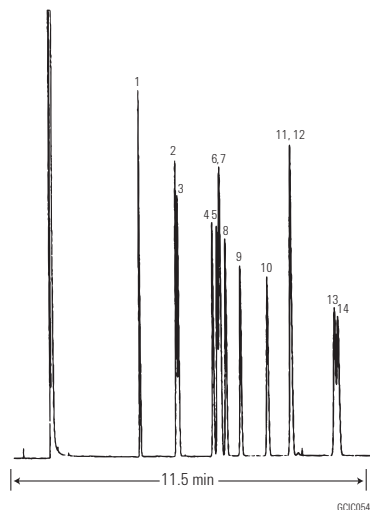

1. 2,6-Xylenol
2. 2-Cresol
3. Phenol
4. 2-Ethylphenol
5. 2,5-Xylenol
6. 4-Cresol
7. 2,4-Xylenol
8. 3-Cresol
9. 2-Isopropylphenol
10. 2,3-Xylenol
11. 3,5-Xylenol
12. 4-Ethylphenol
13. 3,4-Xylenol
14. 2,3,5-Trimethylphenol

Halocarbons

Column: GS-GasPro
113-4332
30 m x 0.32 mm

Carrier: Helium at 30 cm/s

Oven: 130 °C for 4 min
130-225 °C at 10 °C/min
225 °C hold

Injection: Split, 250 °C
Split ratio 1:67

Detector: FID, 250 °C


Sample: 1 µL

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

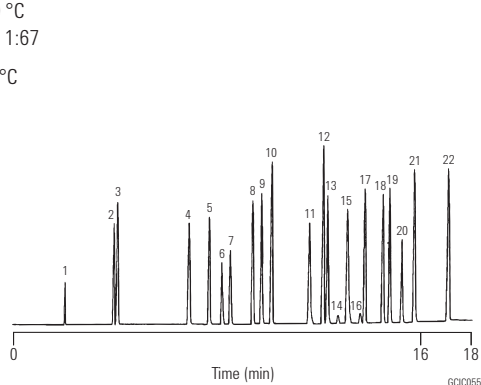

- | | |
|---|---|
| 1. CH ₄ | 12. cis-ClCH=CHCl |
| 2. CHClF ₂ (Freon 22) | 13. CHCl ₃ |
| 3. CCl ₂ F ₂ (Freon 12) | 14. CCl ₄ |
| 4. ClCF ₂ CF ₂ Cl (Freon 114) | 15. CCl ₄ |
| 5. CHCl ₂ F (Freon 21) | 16. CCl ₄ |
| 6. CCl ₃ F (Freon 11) | 17. CH ₃ CH ₂ I |
| 7. CF ₂ Br ₂ (Freon 12B2) | 18. CH ₂ Br ₂ |
| 8. CH ₃ I | 19. CHCl ₂ Br |
| 9. CH ₂ Cl ₂ | 20. C ₄ F ₉ I |
| 10. trans-ClCH=CHCl | 21. CHClBr ₂ |
| 11. CF ₃ CCl ₃ (Freon 113) | 22. CH ₃ CH ₂ CH ₂ I |

Ethylene Oxide

Column: DB-WAX
122-7032
30 m x 0.25 mm, 0.25 µm


Carrier: Helium at 1 mL/min

Oven: 60 °C for 2 min
60-180 °C at 16 °C/min

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Courtesy of J. Chromatogr. Sci., 28:97 [1990]

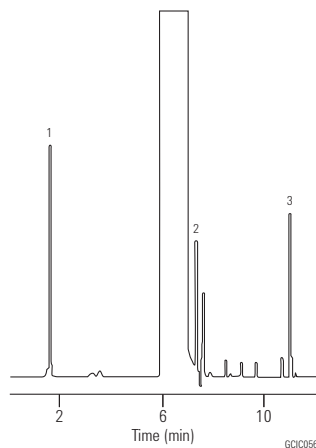

1. Ethylene oxide
2. 2-Chloroethanol
3. Ethylene glycol (solvent: dimethylformamide)

Impurities in Mixed Xylenes

Column: DB-WAXetr
123-7362
60 m x 0.32 mm, 0.25 µm


Carrier: Helium at 20 cm/s,
measured at 145 °C

Oven: 60 °C for 10 min
60-150 °C at 5 °C/min
150 °C for 10 min

Injection: Split, 230 °C
Split ratio 1:150

Detector: FID, 240 °C

ASTM D2360

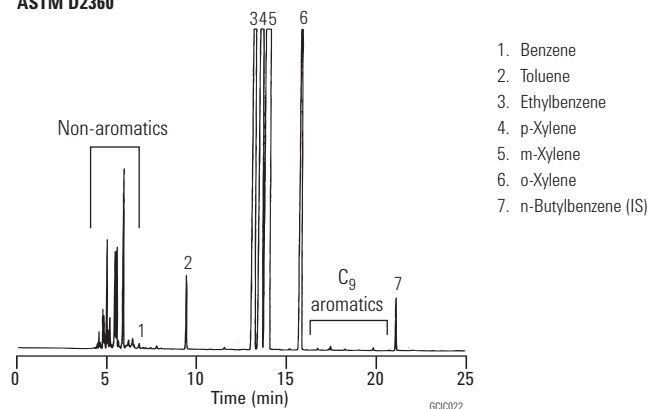**High Resolution Separation of Xylene Isomers**

Column: CP-Chirasil-Dex CB
CP7502
25 m x 0.25 mm, 0.25 µm

Sample: 0.5 µL


Sample Conc: 10-20%

Carrier: Helium, 40 kPa, 6 psi

Oven: 80 °C, (6 min) to 130 °C, 25 °C/min

Injection: Split, T=210 °C, 1:20

Detector: FID, T=230 °C

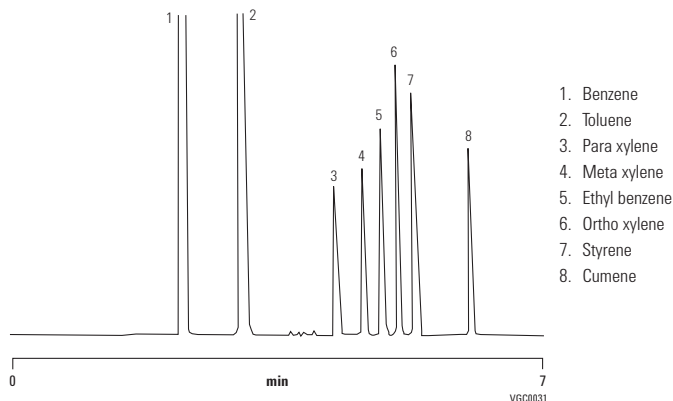**Halothane**

Column: GS-GasPro
113-4312
15 m x 0.32 mm


Carrier: Helium at 45 cm/s

Oven: 240 °C isothermal

Injection: Split, 200 °C
Split ratio 1:100

Detector: FID, 200 °C

Sample: 0.2 µL

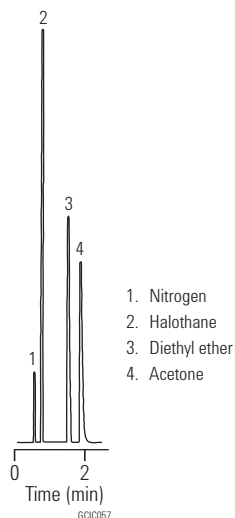**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

Inorganic Hydride Gases

Column: HP-1
19091Z-205
50 m x 0.20 mm, 0.50 μ m


Carrier: Helium, 35 cm/s

Oven: A: -70 °C isothermal
B: 30 °C isothermal

Injection: Split ratio 25:1

Detector: FPD, 535 μ m filter

Sample: 1 μ L

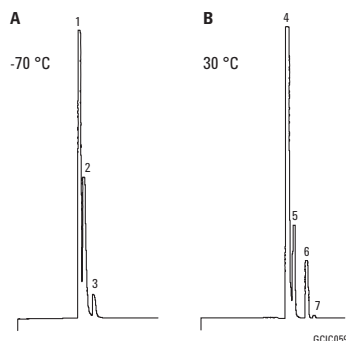

1. Arsine 0.1%
2. Phosphine 0.1%
3. Selenide 0.1%
4. Diborane 0.10 ppm
5. Tetraborane 0.10 ppm
6. Pentaborane 0.10 ppm
7. Dihydropentaborane 0.60 ppm


Polymer Additives

Column: HP-35 (use only 10 m)
19091G-013
30 m x 0.32 mm, 0.15 μ m

Carrier: Helium, 6 psi (4 mL/min at 50 °C) hold for 5 min,
ramp to 50 psi (21 mL/min at 350 °C) at 5 psi/min

Injection: EPC on-column, oven track 0.5 μ L injection

Detector: FID

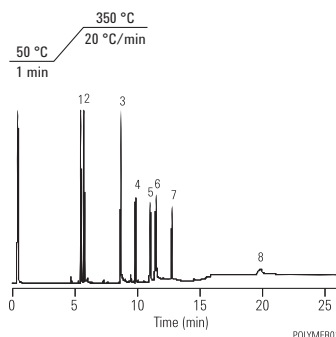

1. BHT
2. BHEB
3. Tinuvin P
4. Isonox 129
5. Irgafos 168
6. Irganox 1076
7. MD 1024
8. Irganox 1010

Fast Separation of Silanes


Column: VF-200ms
CP8860
30 m x 0.25 mm, 1.00 μ m

Carrier: Hydrogen, ca 1.0 mL/min, 60 kPa

Oven: 50 °C

Injection: Split/splitless, in split mode, 1:100

Detector: FID

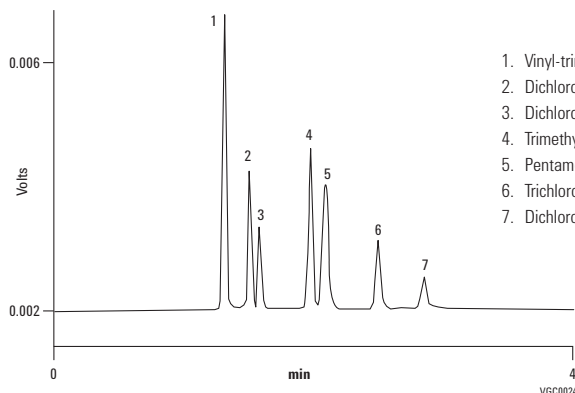

1. Vinyl-trimethyl silane
2. Dichloromethyl silane
3. Dichloromethane
4. Trimethylchloro silane
5. Pentamethyl disiloxane
6. Trichloromethyl silane
7. Dichlorodimethyl silane

Sulfur Gases

Column: PoraPLOT U
CP7584
25 m x 0.53 mm, 20.00 μ m


Sample: ± 100 ppm

Carrier: H_2

Oven: 50 $^{\circ}C$

Injection: 100 mL/min

Detector: FPD

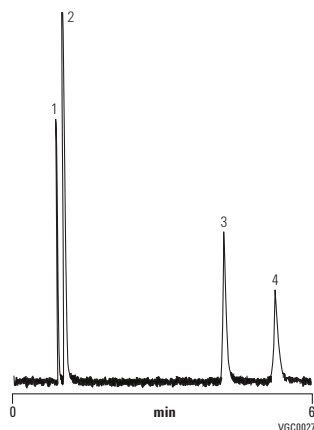

1. Hydrogen sulfide
2. Carbonyl sulfide
3. Sulfur dioxide
4. Methyl sulfide

Analysis of Acetylenes' Mixture

Column: Select Al_2O_3
CP7432
50 m x 0.53 mm

Sample Conc: Approx 100 ppm in nitrogen, synthetic standard


Carrier: Helium, 4 psig, 4 min to 11 psig, 0.5 psig/min, 2 min

Oven: 40 $^{\circ}C$, 5 min to 160 $^{\circ}C$, 10 $^{\circ}C$ /min to 200 $^{\circ}C$,
20 $^{\circ}C$ /min, hold 1 min

Injection: Split, 60 mL/min

Detector: FID

Courtesy of J. Luong, Dow Chemical Canada

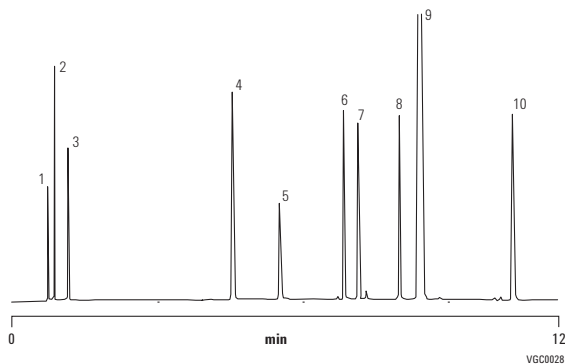

1. Methane
2. Ethane
3. Ethylene
4. n-Butane
5. Propadiene
6. 1-Butene
7. Iso-butene
8. 1,2-Butadiene
9. 1,3-Butadiene
10. Ethyl acetylene

Forensic Toxicology and Pharma Applications

DB-Select 624 UI for <467>

Megabore

Early Eluting Peaks

Column: DB-Select 624 Ultra Inert
125-0334UI
30 m x 0.53 mm, 3.00 µm

Carrier: Helium 44 cm/s (approx. 6 mL/min) set at 40 °C,
EPC – Constant Flow


Oven: 40 °C 20 min hold, then 10°/min to 170 °C

Injection: 20 Hz

Detector: FID at 240 °C, H₂ at 30 mL/min
Air at 400 mL/min
N₂ makeup at 35 mL/min
(constant column + makeup)

Sample: FID signal

- | | |
|--------------------------------|--|
| 1. Methane | 18. 1,1,1-Trichloroethane |
| 2. Methanol | 19. Cyclohexane |
| 3. Ethanol | 20. Carbon tetrachloride |
| 4. Diethyl ether | 21. Benzene |
| 5. 1,1-Dichloroethylene | 22. 1,2-Dichloroethane |
| 6. 2-Propanol | 23. Isooctane (2,2,4-trimethylpentane) |
| 7. Acetonitrile | 24. 3-Methyl-2-butanone |
| 8. Methyl acetate | 25. n-Heptane |
| 9. Dichloromethane | 26. Trichloroethylene |
| 10. trans-1,2-Dichloroethylene | 27. Methylcyclohexane |
| 11. n-Hexane | 28. 1,4-Dioxane |
| 12. 1-Propanol | 29. Propyl acetate |
| 13. Nitromethane | 30. 2-Ethoxyethanol |
| 14. cis-1,2-Dichloroethylene | 31. Pyridine excellent peak shape |
| 15. Ethyl acetate | 32. Toluene |
| 16. 2-Butanol | 33. 3-Methyl-1-butanol |
| 17. Chloroform | |

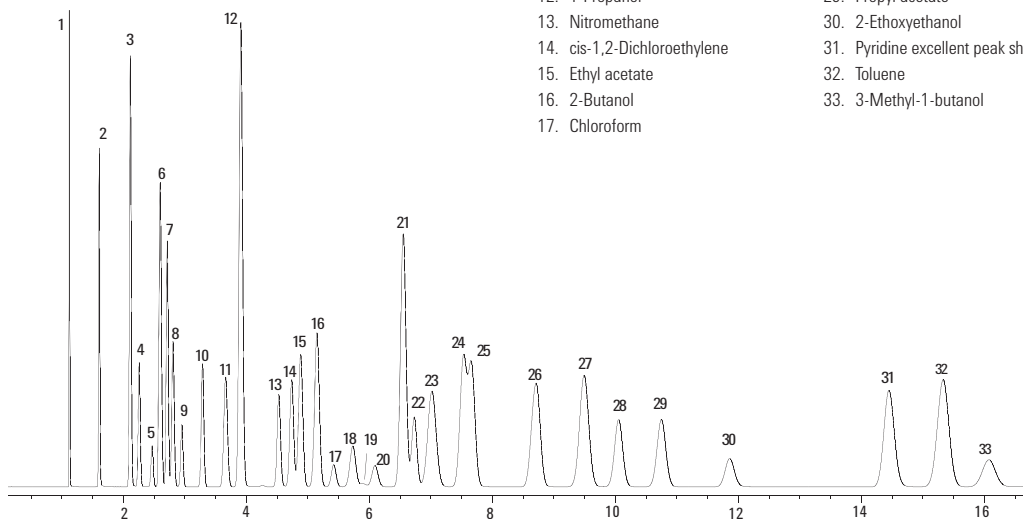

Benzodiazepines I

Column: DB-5ms Ultra Inert
122-5532UI
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen, 53 cm/s, constant flow
1.6 for 11 min
1.6-2.4 at 60 mL/min, hold 2 min
2.4-5.0 at 50 mL/min, hold 9 min


Oven: 170 °C for 3.2 min
170-250 °C at 24.7 °C/min, hold 5.3 min
250-280 °C at 18.6 °C/min, hold 4.0 min
280-325 °C at 50.0 °C/min, hold 4.0 min

Injection: Pulsed splitless, 280 °C
20 psi pulse pressure for 0.38 min
50 mL/min purge at 0.40 min
Direct connect liner (p/n G1544-80730)

Detector: FID, 350 °C

Sample: 1 µL of 5-10 ppm

Analysis of benzodiazepines and other drugs is particularly challenging because of their high level of activity. For this reason, all aspects of the sample path – particularly the GC Column – must be as inert as possible.

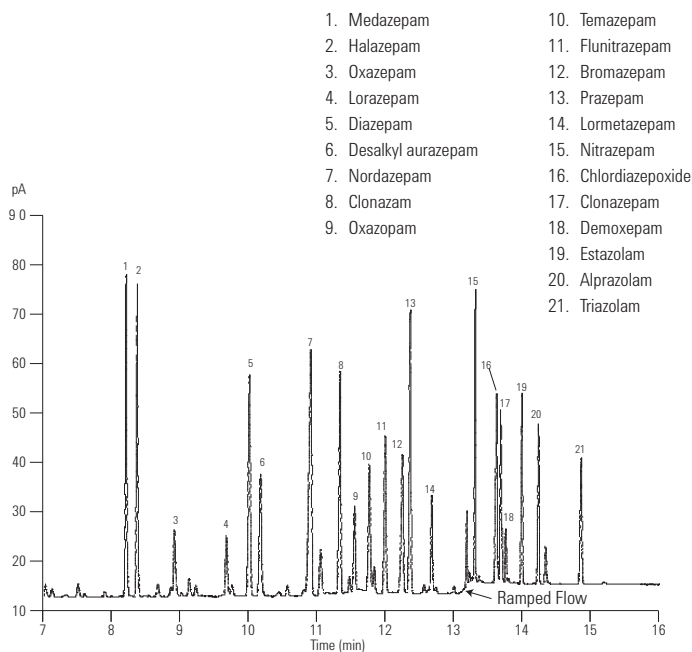**Amphetamines and Precursors – TMS Derivatives**

Column: DB-5
121-5023
20 m x 0.18 mm, 0.40 µm

Carrier: Helium at 39 cm/s, measured at 100 °C

Oven: 100-240 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 2 µg/µL each in pyridine


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

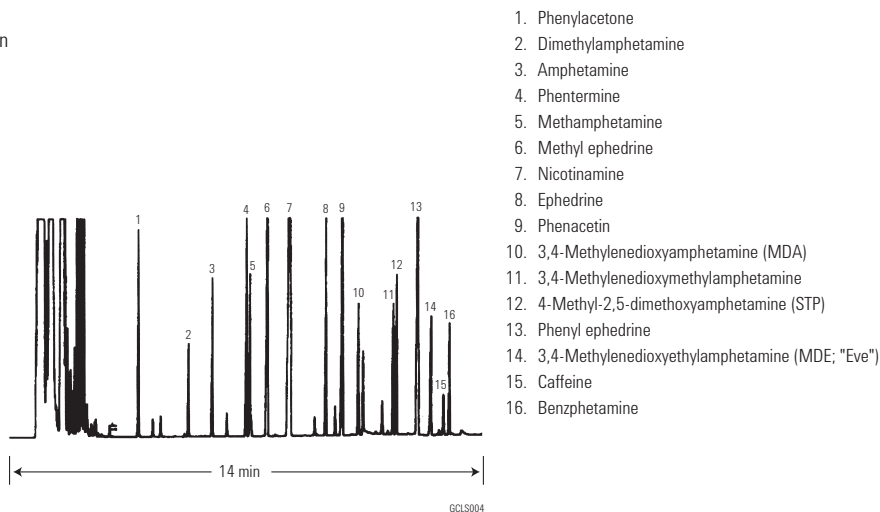

Barbiturates

Column: DB-35ms
122-3832
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 31 cm/s, measured at 50 °C

Oven: 50 °C for 0.5 min
50-150 °C at 25 °C/min
150-300 °C at 10 °C/min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: MSD, 280 °C transfer line
full scan at m/z 40-270


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

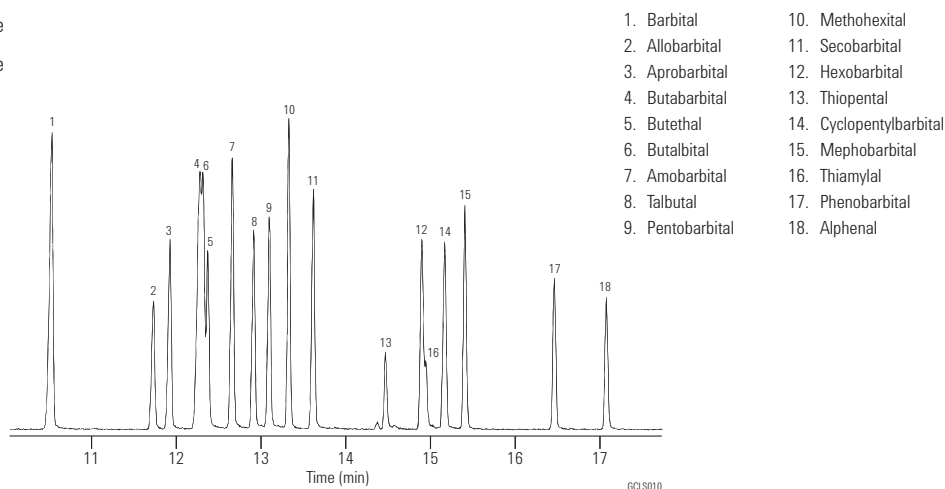

Narcotics

Column: DB-5ms
122-5532
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 31 cm/s, measured at 50 °C

Oven: 50 °C for 0.5 min
50-150 °C at 25 °C/min
150-325 °C at 10 °C/min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: MSD, 300 °C transfer line
full scan at m/z 40-380


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

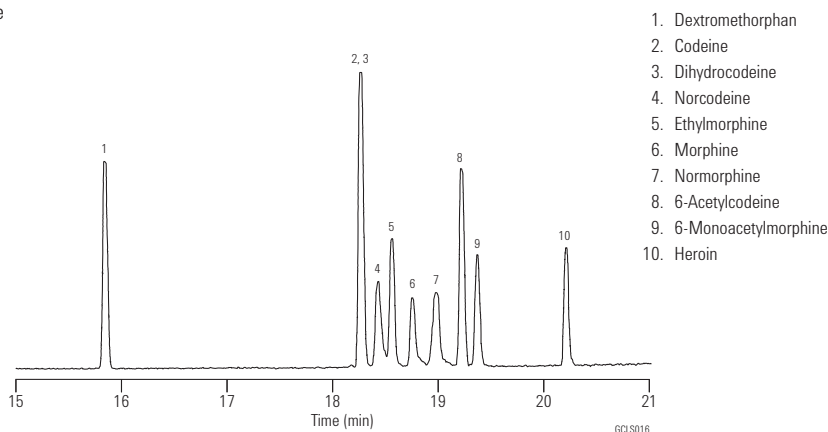

Blood Alcohols I (Static Headspace/Split)

Column: DB-ALC1
125-9134
30 m x 0.53 mm, 3.00 µm

Carrier: Helium at 80 cm/s,
measured at 40 °C

Oven: 40 °C isothermal

Sampler: Headspace

Injection: Split, 250 °C
Split ratio 1:10


Detector: FID, 300 °C
Nitrogen makeup gas
at 23 mL/min

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

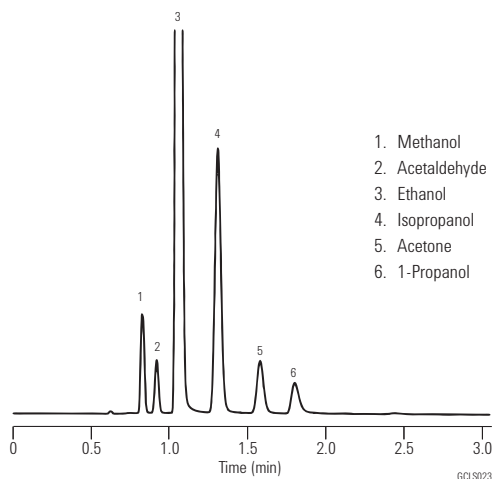

Blood Alcohols II (Static Headspace/Split)

Column: DB-ALC2
125-9234
30 m x 0.53 mm, 2.00 µm

Carrier: Helium at 80 cm/s,
measured at 40 °C

Oven: 40 °C isothermal

Sampler: Headspace

Oven: 70 °C
Loop: 80 °C
Transfer line: 90 °C
Vial equil. time: 10 min
Pressurization time: 0.20 min
Loop fill time: 0.20 min
Loop equil. time: 0.05 min
Inject time: 0.1-0.2 min
Sample loop size: 1.0 mL

Injection: Split, 250 °C
Split ratio 1:10

Detector: FID, 300 °C
Nitrogen makeup gas
at 23 mL/min


Sample: 0.1% Ethanol,
0.001% Others

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

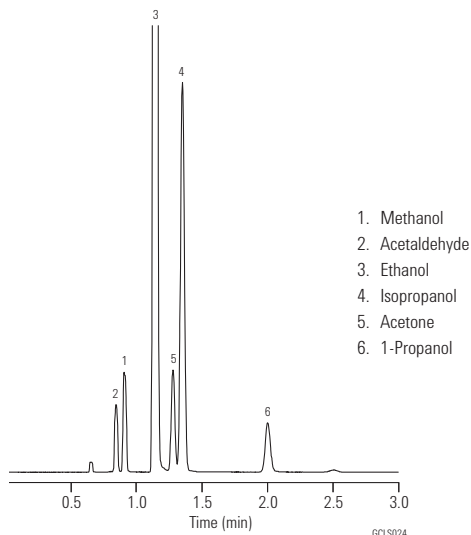

Residual Solvents, DMI Diluent

Column: DB-624
123-1364
60 m x 0.32 mm, 1.80 µm

Oven: 50-60 °C, 1 °C/min
60-115 °C, 9.2 °C/min
115-220 °C, 35 °C/min
220 °C – hold 6 min


Sampler: Headspace
Plate 140 °C
Transfer line, valve 250 °C
Sample loop 2 mL

Injection: Split, 250 °C
Split ratio 1:18

Detector: FID, 270 °C
Nitrogen makeup

Sample: 5,000 ppm standard

- | | | |
|---------------------------------------|--------------------------|--|
| 1. Methanol | 12. 2-Butanone (MEK) | 23. MIBK (2-Pentanone) |
| 2. Ethanol | 13. Ethyl acetate | 24. Toluene |
| 3. Acetone | 14. 2-Butanol | 25. 1-Pentanol |
| 4. 2-Propanol | 15. Tetrahydrofuran | 26. n,n-Dimethylformamide (DMF) |
| 5. Acetonitrile | 16. Cyclohexane | 27. Ethyl benzene |
| 6. Methylene chloride | 17. Isopropyl acetate | 28. m,p-Xylene |
| 7. 2-Methyl-2-propanol (tert-butanol) | 18. 1,2-Dimethoxyethane | 29. o-Xylene |
| 8. MTBE | 19. Heptane | 30. Dimethyl sulfoxide (DMSO) |
| 9. Hexane | 20. 1-Methoxy-2-propanol | 31. n,n-Dimethylacetamide |
| 10. 1-Propanol | 21. Methylcyclohexane | 32. n-Methylpyrrolidone |
| 11. DMI impurity | 22. 2-Ethoxyethanol | 33. 1,3-Dimethyl-2-imidazolidinone (DMI) |

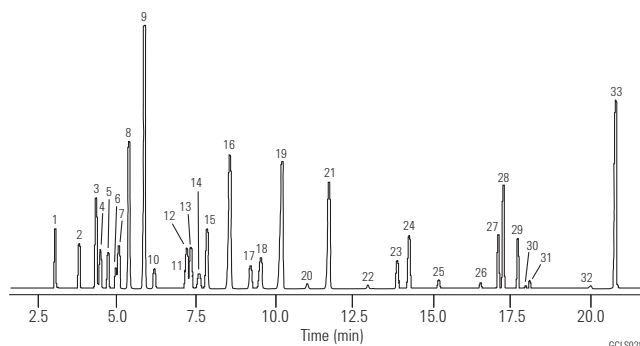

Special thanks to Julie Kancler, Brian Wallace, Teledyne.

Suggested Supplies

Septum: 11 mm Advanced Green septa,
5183-4759

Liner: Direct, 1.5 mm id,
18740-80200

Seal: Gold plated seal, 18740-20885

Underivatized Drugs of Abuse – Agilent Fast Toxicology Analyzer

Column: DB-35ms Ultra Inert
122-3812UI
15 m x 0.25 mm, 0.25 µm

Carrier: Helium, fixed pressure 35.0 psi

Injection: Splitless 1 µL 280 °C, total flow 56.4 mL/min,
3 mL/min switched septum purge, gas saver off,
50 mL/min after 0.4 min

Liner: Splitless, dual taper, deactivated, 4 mm id,
(p/n 5181-3315)

Sample: Agilent GC/MS toxicology checkout mixture
(p/n 5190-0471)


Backflush: Post run: 1 min 1 psi inlet, 75 psi aux EPC

Oven: 100 °C (0.25 min) to 345 °C
(40 °C/min, 2.25 min hold)

Detector: MSD: Transfer line 300 °C, source 300 °C
Quadrupole: 180 °C scan mode
NPD: Bloss bead 300 °C H₂ 3 mL/min, 60 mL/min air,
11 mL/min makeup and column flow

CFT Device: 2-Way splitter with solvent venting between
MSD and NPD

- | | | |
|---|-----------------------------|----------------------|
| 1. Amphetamine | 9. Phencyclidine | 19. Oxycodone |
| 2. Phentermine | 10. Methadone | 20. Temazepam |
| 3. Methamphetamine | 11. Cocaine | 21. Diacetylmorphine |
| 4. Nicotine | 12. SKF-525a (RTL compound) | 22. Flunitrazepam |
| 5. Methylenedioxyamphetamine (MDA) | 13. Oxazepam | 23. Nitrazepam |
| 6. Methylenedioxymethamphetamine (MDMA) | 14. Tetrahydrocannabinol | 24. Clonazepam |
| 7. Methylenedioxyethylamphetamine | 15. Codeine | 25. Alprazolam |
| 8. Meperidine | 16. Lorazepam | 26. Verapamil |
| | 17. Diazepam | 27. Strychnine |
| | 18. Hydrocodone | 28. Trazodone |

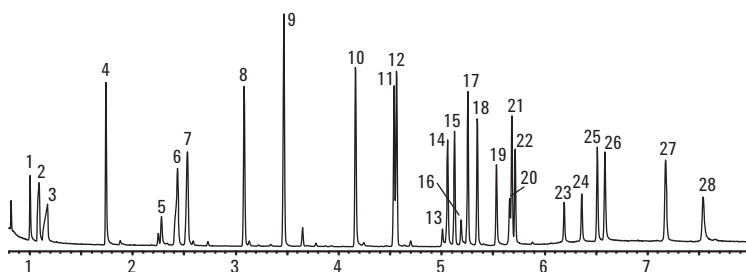

Example NPD chromatogram of underivatized drugs of abuse 5 ng/component on an Agilent J&W DB-35ms UI column. Component number 12 is used for retention time locking in the deconvolution reporting software database.

Benzodiazepines II

Column: DB-35ms
122-3832
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 31 cm/s, measured at 50 °C

Oven: 50 °C for 0.5 min
50-150 °C at 25 °C/min
150-340 °C at 10 °C/min
340 °C for 6 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: MSD, 280 °C transfer line
full scan at m/z 40-400

Suggested Supplies


Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

- | | |
|----------------------|-------------------|
| 1. Medazepam | 10. Flunitrazepam |
| 2. Halazepam | 11. Delorazepam |
| 3. Oxazepam | 12. Bromazepam |
| 4. Lorazepam | 13. Prazepam |
| 5. Diazepam | 14. Flurazepam |
| 6. Demoxepam | 15. Clonazepam |
| 7. Desmethyldiazepam | 16. Alprazolam |
| 8. Clobazam | 17. Triazolam |
| 9. Temazepam | |

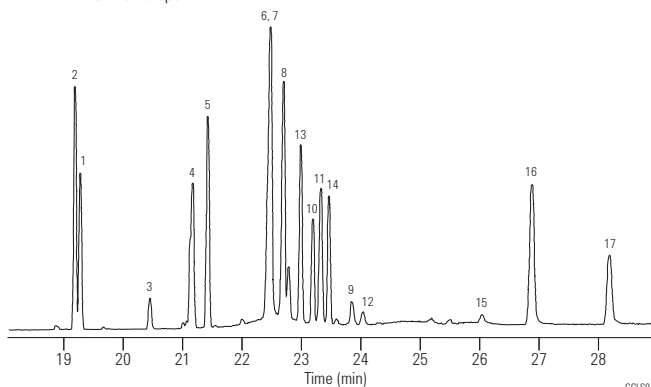

GCL5011

Drug Screen

Column: DB-1ms
122-0132
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 40 cm/s,
measured at 50 °C


Oven: 50 °C for 1.0 min
50-125 °C at 25 °C/min
125-325 °C at 10 °C/min
325 °C for 5 min

Injection: Cold splitless
Optic II injector, 50-250 °C at 10 °C/s
45 s purge activation time

Detector: FID, 300 °C

Sample: 1 µL injection of 50-150 ppm standard

- | | |
|---------------------------------|-------------------|
| 1. Nicotine | 10. Cocaine |
| 2. Caffeine | 11. Desipramine |
| 3. Glutethimide | 12. Carbamazepine |
| 4. Lidocaine | 13. Trimipramine |
| 5. PCP | 14. Heroin |
| 6. Phenobarbital | 15. Fentanyl |
| 7. Methadone primary metabolite | 16. Ibogaine |
| 8. Methaqualone | 17. Triazolam |
| 9. Methadone | 18. LSD |

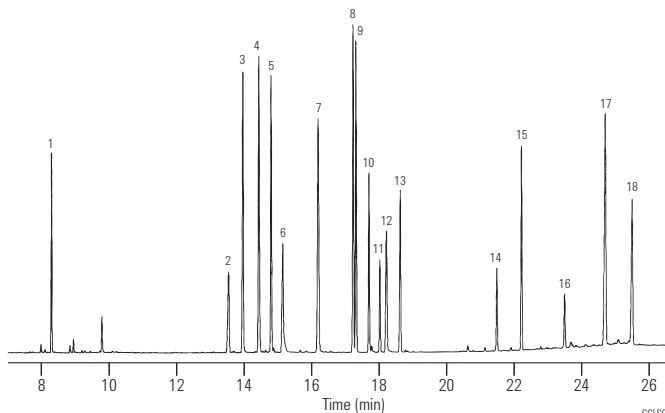

GCL5002

Common Drug Screen

Column: DB-5
122-5032
30 m x 0.25 mm, 0.25 µm

Column: DB-17
122-1732
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 41 cm/s,
measured at 80 °C

Oven: 80 °C for 1 min
80-280 °C at 10 °C/min
280 °C for 9 min

Injection: Split, 250 °C
Split ratio 1:40

Detector: FID, 300 °C

Suggested Supplies



Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper,
glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP,
5181-1267

	DB-17 Time	DB-5 Time		DB-17 Time	DB-5 Time
1. Nicotine	9.87	8.57	18. Hexobarbital	17.52	15.22
2. Phenmetrazine	11.8	9.95	19. Doxylamine	17.69	15.87
3. Ibuprofen	12.06	10.64	20. Caffeine	18.05	13.11
4. Procaine	13.48	14.82	21. Chlorpheniramine	18.47	16.35
5. Allobarbitol	13.91	12.02	22. Methapyrilene	18.72	16.68
6. Aprobital	14.14	12.27	23. Theryldiamine	18.87	16.85
7. Butabarbital	14.56	12.76	24. Phenobarbital	19.11	16.29
8. Secobarbital	14.87	14.31	25. Bromopheniramine	19.71	17.39
9. Pentobarbital	15.41	13.73	26. Chlorcyclizine	20.75	19.13
10. Phenacetin	15.72	12.94	27. Cocaine	21.32	18.88
11. Amobarbital	15.87	13.43	28. Pyrrobutamine	22.79	20.89
12. Benzphetamine	16.14	14.96	29. Codeine	24.27	20.66
13. Acetaminophen	16.34	11.12	30. Diazepam	25.27	21.13
14. Hydroxyphenamate	16.47	15.31	31. Morphine	25.36	21.12
15. Dimenhydrinate	16.93	13.79	32. Hydrocodone	25.98	21.26
16. Meprobamate	17.12	14.44	33. Oxymorphone	28.27	22.21
17. Benactyzine	17.26	14.71	34. Heroin	29.32	23.14

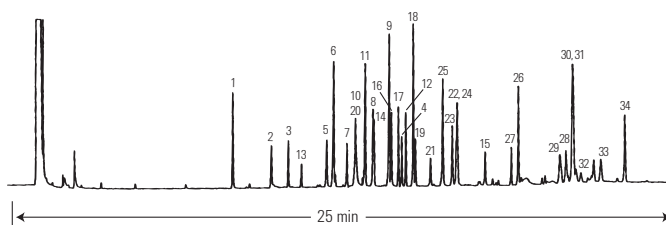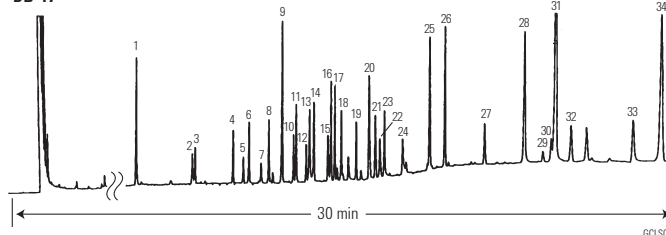DB-5**DB-17**

GCL5001

Urine Drug Screen


Column: Ultra 2
19091B-115
50 m x 0.32 mm, 0.52 µm

Carrier: Hydrogen, 80 cm/s

Oven: 45 °C for 1.5 min
45-300 °C at 6 °C/min

Injection: Splitless

Detector: FID

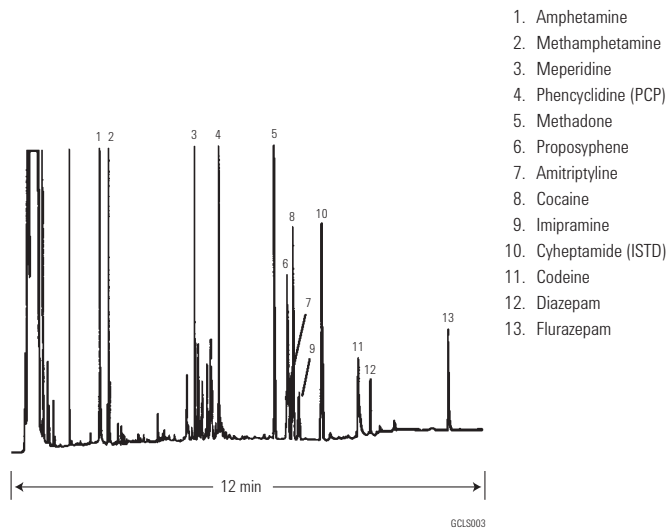

Analysis of Drugs of Abuse in Urine via GC/MS

Column: VF-DA
CP8964
12 m x 0.20 mm, Optimized µm

Sample: 1 µL

Solvent: Methanol

Carrier: He, ca 1.0 mL/min


Oven: 70 °C, 1.2 min to 200 °C,
20 °C/min to 270 °C,
7 °C/min to 320 °C, 20 °C/min

Pressure: 58.7 kPa, 2.2 min to 97 kPa, 58 kPa/min to 132 kPa,
3 kPa/min to 180 kPa, 12 kPa/min

Injection: Splitless

Detector: MS

Derivatization: Acetic acid anhydride to form acetates

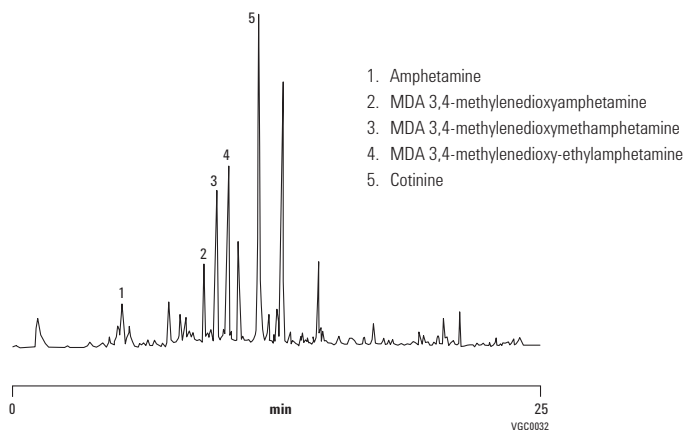

Anesthetics

Column: DB-5ms EVDX
128-8522
25 m x 0.20 mm, 0.33 µm

Carrier: Helium at 35 cm/s, measured at 55 °C


Oven: 55 °C for 1 min
55-130 °C at 25 °C/min
130-325 °C at 15 °C/min

Injection: Splitless, 250 °C
45 s purge activation time

Detector: MSD, 280 °C transfer line
full scan at m/z 35-400

Sample: 1 µL of 50-100 ng/µL standard in methanol

1. Salicylamide
2. Benzocaine
3. Lidocaine
4. Procaine
5. Nefopam
6. Mepivacaine
7. Tetracaine
8. Butacaine
9. Dibucaine

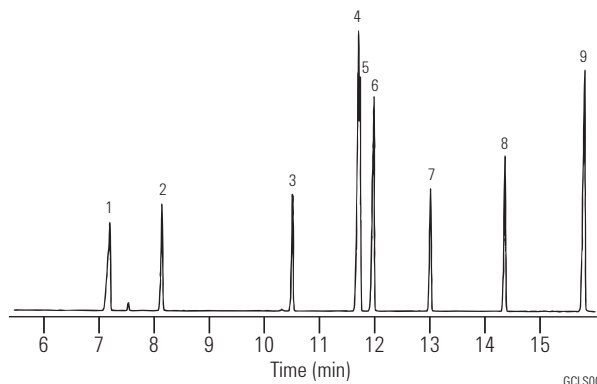**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316
Seal: Gold plated seal, 18740-20885
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Anticonvulsants

Column: DB-1
125-1032
30 m x 0.53 mm, 1.50 µm

Carrier: Helium at 8 mL/min


Oven: 160 °C for 2 min
160-275 °C at 15 °C/min

Injection: Megabore direct, 250 °C

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 100 ng/µL in methanol

1. Ethosuximide
2. Methsuximide
3. Phensuximide
4. N-Desmethyl methsuximide
5. Phenylethylmalonamide
6. Phenobarbital
7. Primidone
8. Carbamazepine
9. Phenytoin
10. 5-Methyl-5-phenylhydantoin

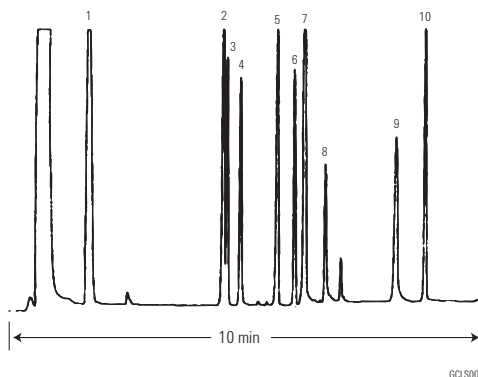**Suggested Supplies**

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730
Seal: Gold plated seal, 18740-20885
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Antihistamines

Column: DB-5
123-5032
30 m x 0.32 mm, 0.25 µm

Carrier: Helium at 40 cm/s, measured at 55 °C

Oven: 55 °C for 1 min
55-175 °C at 30 °C/min
175-320 °C at 10 °C/min
320 °C for 1 min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 50 ng/µL each in methanol

- | | |
|----------------------|--------------------|
| 1. Pheniramine | 13. Thonzylamine |
| 2. Dimenhydrinate | 14. Chlorcyclizine |
| 3. Diphenhydramine | 15. Pyrilamine |
| 4. Doxylamine | 16. Triprolidine |
| 5. Phenyltoloxamine | 17. Promethazine |
| 6. Triphenennamine | 18. Antazoline |
| 7. Methapyrilene | 19. Clemizole |
| 8. Chlorpheniramine | 20. Hydroxyzine |
| 9. Cyclizine | 21. Meclizine |
| 10. Carbinoxamine | 22. Cinnanzine |
| 11. Diphenylpyraline | 23. Bucizine |
| 12. Bromopheniramine | |


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Splitless, single taper, deactivated, 4 mm id, 5181-3316

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

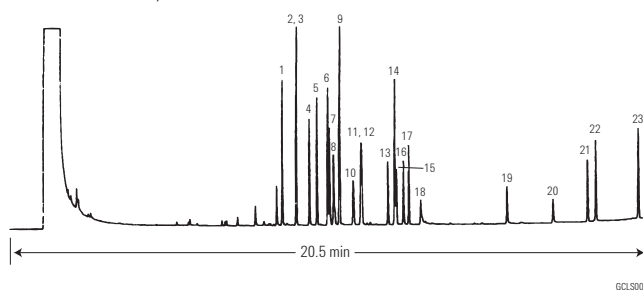

GCL5007

Antiepileptic Drugs


Column: Ultra 2
19091B-012
25 m x 0.32 mm, 0.17 µm

Carrier: Helium, 14 psi

Oven: 100-230 °C at 15 °C/min

Injection: Split ratio 35:1

Detector: NPD

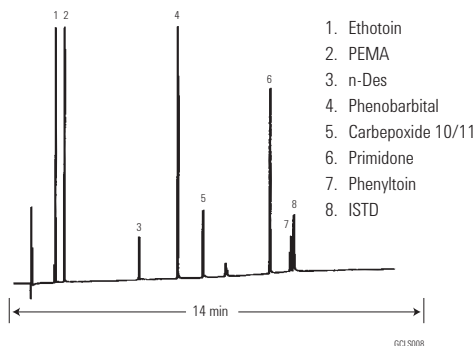

GCL5008

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Tricyclic Antipsychotics


Column: Ultra 2
19091B-101
12 m x 0.20 mm, 0.33 µm

Carrier: Hydrogen, 106 cm/s

Oven: 250 °C for 3 min
250-290 °C at 10 °C/min
290 °C for 10 min

Injection: Split ratio 75:1

Detector: FPD

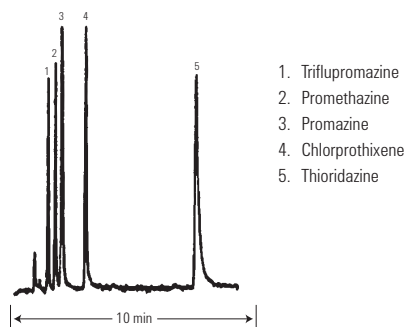

GCL5009

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: General purpose split/splitless liner, taper, glass wool, 5183-4711

Seal: Gold plated seal, 18740-20885

Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

Fentanyl

Column: DB-1701
125-0732
30 m x 0.53 mm, 1.00 µm

Carrier: Hydrogen at 15 mL/min

Oven: 270 °C isothermal

Injection: Split, 250 °C
Split ratio 1:5

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 0.8 µL


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

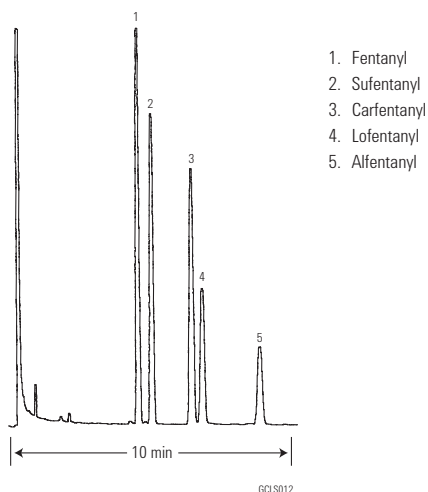**Tocopherols**

Column: DB-17ms
122-4732
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 40 cm/s,
measured at 150 °C

Oven: 300 °C for 1 min
300-320 °C at 25 °C/min
320 °C for 4 min

Injection: Split, 310 °C
Split ratio 1:25

Detector: MSD, 310 °C transfer line
full scan at m/z 45-550

Sample: 1 µL of 1-10 ng/µL in isooctane


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Split, single taper, low pressure drop, glass wool, 5183-4647

Seal: Gold plated seal, 18740-20885

Syringe: 5 µL tapered, FN 23-26s/42/HP, 5181-1273

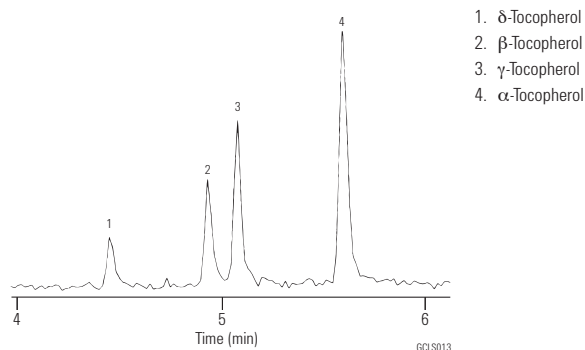

Hallucinogens

Column: DB-17ms
122-4732
30 m x 0.25 mm, 0.25 µm

Carrier: Helium at 30 cm/s, measured at 50 °C

Oven: 50 °C for 0.5 min
50-125 °C at 25 °C/min
125-255 °C at 10 °C/min
255-320 °C at 25 °C/min
320 °C for 16 min


Injection: Splitless, 250 °C
30 s purge activation time

Detector: MSD, 300 °C transfer line
full scan at m/z 40-350

Sample: 1 µL of 10-50 ng/µL standard in methanol

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730
Seal: Gold plated seal, 18740-20885
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

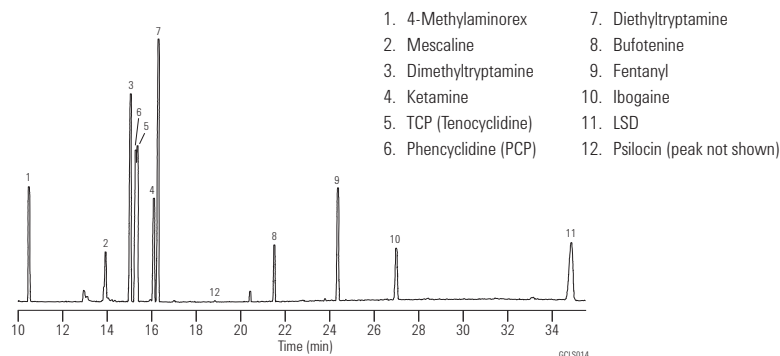

Sedative Hypnotics

Column: DB-5ms EVDX
128-8522
25 m x 0.20 mm, 0.33 µm

Carrier: Helium at 35 cm/s, measured at 55 °C

Oven: 55 °C for 1 min
55-130 °C at 25 °C/min
130-325 °C at 15 °C/min
325 °C for 4 min


Injection: Splitless, 250 °C
45 s purge activation time

Detector: MSD, 280 °C transfer line
full scan at m/z 35-400

Sample: 1 µL of 50-100 ng/µL standard in methanol

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759
Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730
Seal: Gold plated seal, 18740-20885
Syringe: 10 µL tapered, FN 23-26s/42/HP, 5181-1267

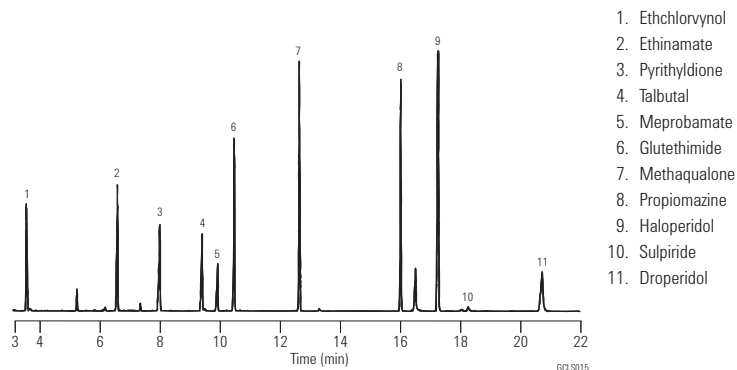

Narcotics and Adulterants

Column: DB-5
123-5032
30 m x 0.32 mm, 0.25 µm


Carrier: Helium at 40 cm/s, measured at 140 °C

Oven: 140-320 °C at 12 °C/min
320 °C for 4 min

Injection: Split, 250 °C
Split ratio 1:75

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 0.5 µg/µL each in methanol

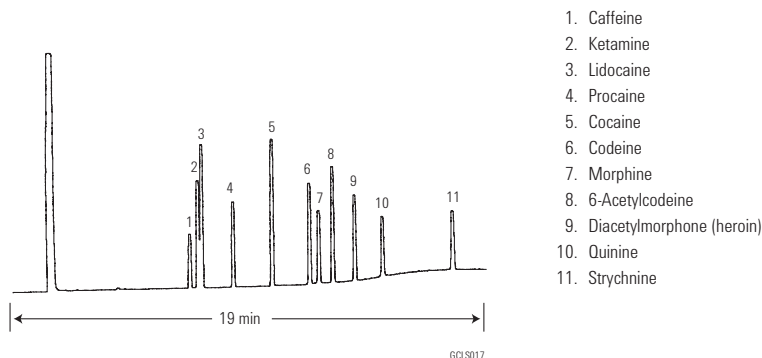**Over-the-Counter Pain Killers – TMS Derivatives**

Column: DB-5
121-5023
20 m x 0.18 mm, 0.40 µm


Carrier: Helium at 39 cm/s, measured at 100 °C

Oven: 100-240 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

Sample: 1 µL of 2 µg/µL each in pyridine

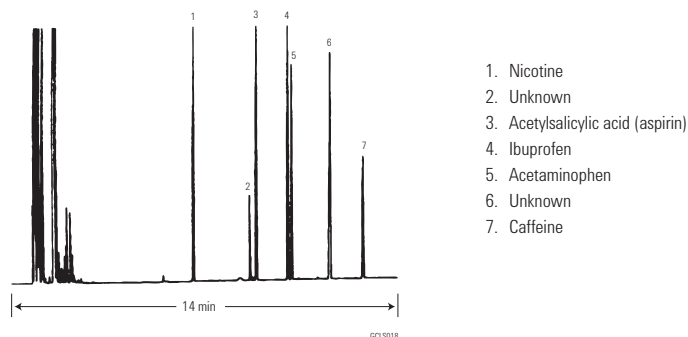**Aspirin and Ibuprofen in Methanol**


Column: DB-FFAP
122-3232
30 m x 0.25 mm, 0.25 µm

Carrier: Hydrogen at 24 cm/s, measured at 180 °C

Oven: 180 °C isothermal

Injection: Split, 250 °C
Split ratio 1:50

Detector: FID, 300 °C
Nitrogen makeup gas at 30 mL/min

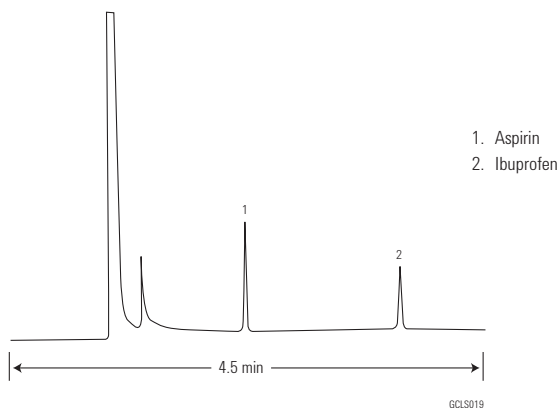

Free Steroids

Column: DB-17
122-1731
30 m x 0.25 mm, 0.15 µm


Carrier: Hydrogen at 44 cm/s

Oven: 260 °C isothermal

Injection: Split, 250 °C
Split ratio 1:100

Detector: FID, 300 °C
Nitrogen makeup gas at
30 mL/min

Sample: 1 µL

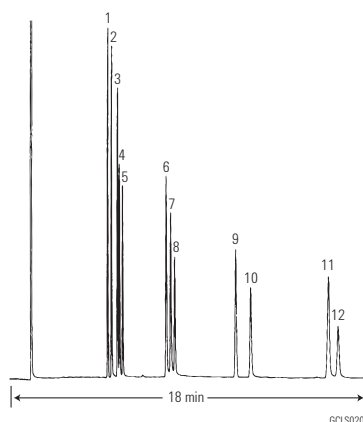

1. Coprostan (5- β -cholestane)
2. 5- β -Androsterone
3. 5- α -Cholestane
4. Androsterone
5. Epiandrosterone (trans-androsterone)
6. 17- α -Estradiol
7. β -Estradiol
8. Estrone
9. Progesterone
10. Cholesterol
11. Estriol
12. Stigmastrol

Anabolic Steroids

Column: DB-1
122-1031
30 m x 0.25 mm, 0.10 µm


Carrier: Helium at 40 cm/s, measured at 180 °C

Oven: 180-320 °C at 10 °C/min
320 °C for 4 min

Injection: Split ratio 1:40

Detector: FID, Nitrogen makeup gas at 30 mL/min

Sample: 2 µL of 0.125 µg/µL each in methanol

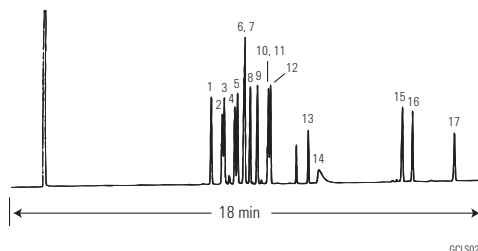

1. Dehydroisoandrosterone (prasterone)
2. 5 α -Androstan-17 α -ol-3-one (stanolone)
3. 19-Nortestosterone (nandrolone)
4. Mesterolone
5. Testosterone
6. 1-Dehydrotestosterone (boldenone)
7. 17 α -Methyltestosterone
8. 1-Dehydro-17- α -methyltestosterone (methandrostenolone)
9. Norethandrolone
10. 1-Dehydrotestosterone acetate
11. Oxymetholone
12. 19-Nortestosterone-17-propionate
13. 4-Chlortestosterone-17-acetate (clostebol)
14. Stanozolol
15. 1-Dehydrotestosterone benzoate
16. 19-Nortestosterone-17-decanoate
17. 1-Dehydrotestosterone undecylenate

Marijuana (Δ^9 -THC) and Major Metabolites – TMS Derivatives

Column: DB-5
123-5032
30 m x 0.32 mm, 0.25 μ m

Carrier: Helium at 40 cm/s, measured at 100 °C

Oven: 100 °C for 1 min
100-175 °C at 30 °C/min
175-295 °C at 12 °C/min

Injection: Splitless, 250 °C
30 s purge activation time

Detector: FID, 300 °C
Nitrogen makeup
gas at 30 mL/min

Sample: 1 μ L of 0.1 μ g/ μ L each in pyridine


Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct connect, single taper, deactivated, 4 mm id, G1544-80730

Seal: Gold plated seal, 18740-20885

Syringe: 10 μ L tapered, FN 23-26s/42/HP, 5181-1267

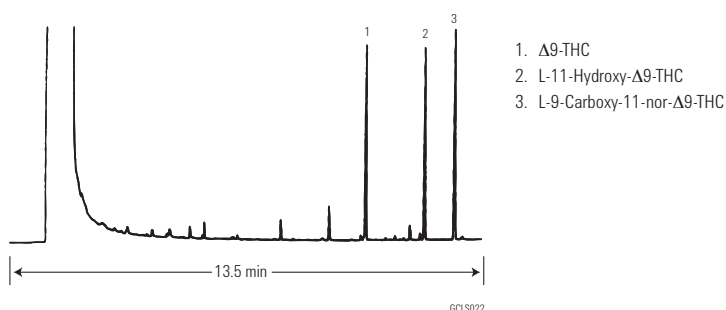

Blood Pollutants I

Column: DB-ALC1
125-9134
30 m x 0.53 mm, 3.00 μ m

Carrier: Helium, 36 cm/s, measured at 40 °C

Oven: 40 °C for 5 min
40-210 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:10


Detector: FID, 300 °C

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

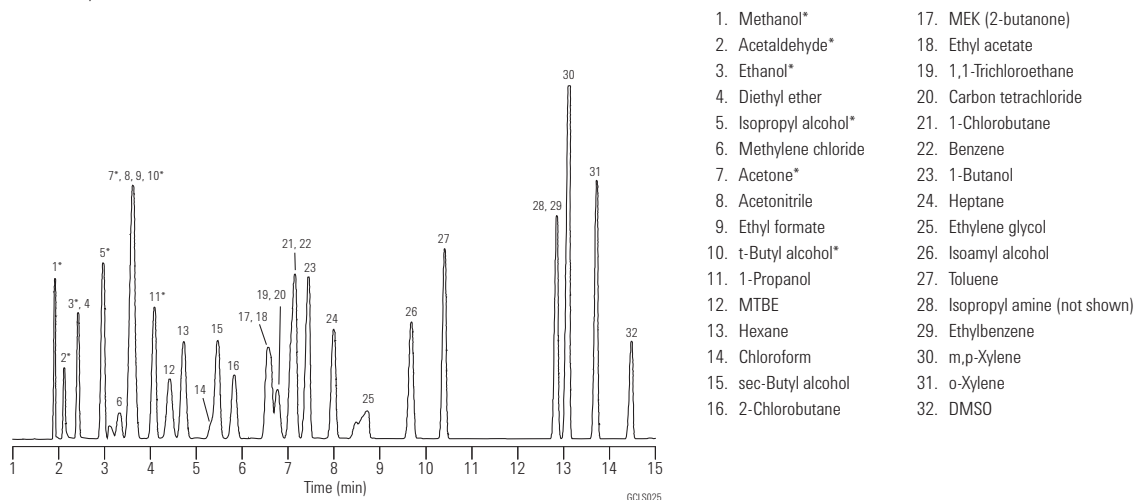

Blood Pollutants II

Column: DB-ALC2
125-9234
30 m x 0.53 mm, 2.00 µm

Carrier: Helium, 36 cm/s, measured at 40 °C

Oven: 40 °C for 5 min
40-210 °C at 10 °C/min

Injection: Split, 250 °C
Split ratio 1:10


Detector: FID, 300 °C

Suggested Supplies

Septum: 11 mm Advanced Green septa, 5183-4759

Liner: Direct, 1.5 mm id, 18740-80200

Seal: Gold plated seal, 18740-20885

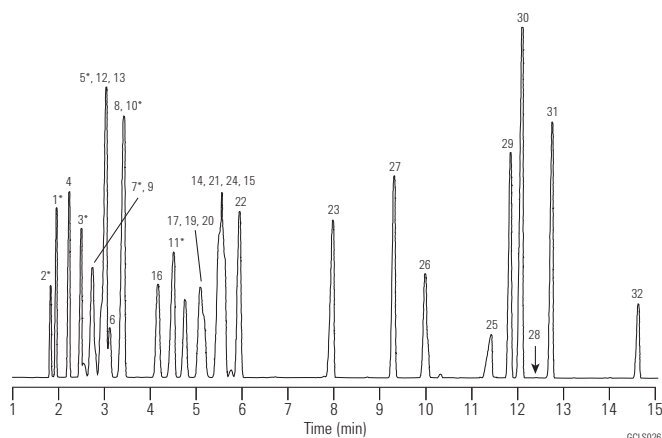

- | | |
|-----------------------|---------------------------------|
| 1. Methanol* | 17. MEK (2-butanone) |
| 2. Acetaldehyde* | 18. Ethyl acetate |
| 3. Ethanol* | 19. 1,1-Trichloroethane |
| 4. Diethyl ether | 20. Carbon tetrachloride |
| 5. Isopropyl alcohol* | 21. 1-Chlorobutane |
| 6. Methylene chloride | 22. Benzene |
| 7. Acetone* | 23. 1-Butanol |
| 8. Acetonitrile | 24. Heptane |
| 9. Ethyl formate | 25. Ethylene glycol |
| 10. t-Butyl alcohol* | 26. Isoamyl alcohol |
| 11. 1-Propanol | 27. Toluene |
| 12. MTBE | 28. Isopropyl amine (not shown) |
| 13. Hexane | 29. Ethylbenzene |
| 14. Chloroform | 30. m,p-Xylene |
| 15. sec-Butyl alcohol | 31. o-Xylene |
| 16. 2-Chlorobutane | 32. DMSO |

Residual Solvents, USP 467

Column: DB-624
125-1334
30 m x 0.53 mm, 3.00 µm


Carrier: Helium at 35 cm/s, measured at 40 °C

Oven: 40 °C for 20 min
40-90 °C at 5 °C/min

Injection: Megabore direct, 250 °C
5 m phenylmethylsilane deactivated
retention gap

Detector: FID, 300 °C
Nitrogen makeup gas at
30 mL/min

- | | |
|-------------------------------|-----------------------------|
| 1. Methanol | 13. Tetrahydrofuran (THF) |
| 2. Ethanol | 14. Chloroform |
| 3. Ethyl ether | 15. sec-Butanol |
| 4. Acetone | 16. Cyclohexane |
| 5. Isopropanol | 17. Benzene |
| 6. Acetonitrile | 18. n-Heptane |
| 7. Methylene chloride | 19. Trichloroethylene |
| 8. tert-Butanol | 20. n-Butanol |
| 9. n-Hexane | 21. 1,4-Dioxane |
| 10. n-Propanol | 22. Pyridine |
| 11. Methyl ethyl ketone (MEK) | 23. Toluene |
| 12. Ethyl acetate | 24. Dimethylformamide (DMF) |

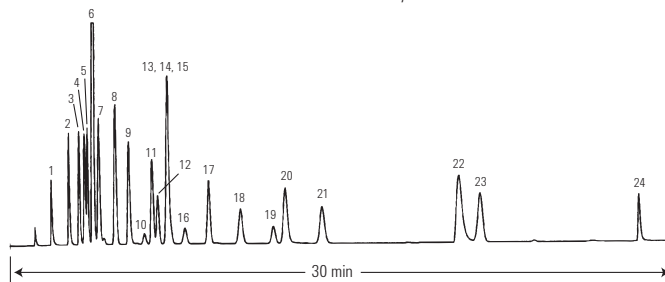

Column Performance for USP <467> Standards

Column: DB-Select 624 Ultra Inert
123-0334UI
30 m x 0.32 mm, 1.80 µm

Carrier: Helium, 2.2 mL/min
constant flow at 40 °C

Oven: 40 °C for 20 min, then
10 °C/min to 240 °C 5 min

Inlet: MMI, 140 °C, 1 µL split 5:1

Inlet liner: 1 mm straight single taper Ultra Inert liner

Sample Conc: 1.0 mL loop


Detector: FID: 250 °C, H₂ 30 mL/min, air 400 mL/min,
N₂ constant col + makeup = 30 mL/min

Suggested Supplies

Septum: Non-stick bleed and temperature optimized (BTO) septa, 11 mm, 50/pk, 5183-4757


Liner: Liner, GC, Ultra Inert, straight, 1 mm id, 5190-4047

Seal: Certified gold plated seal kit, includes washer, 10/pk, 5190-2209

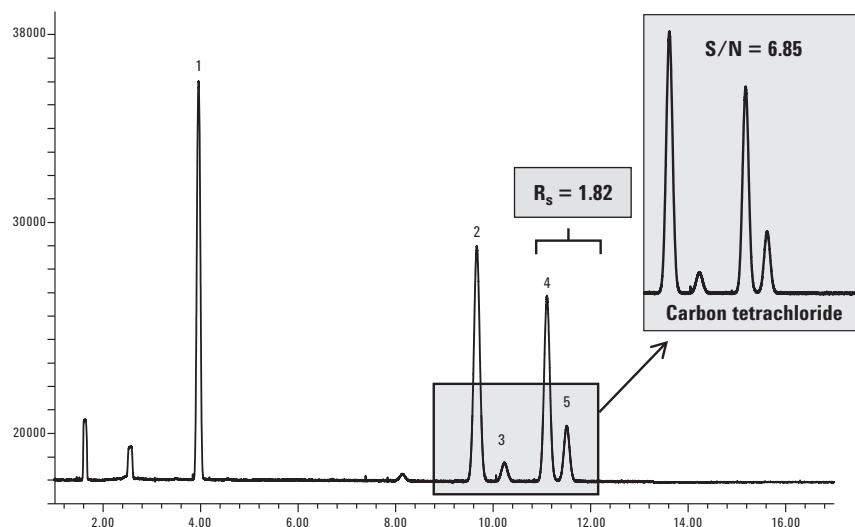

FID trace of Class 1 solvent standard at USP <467> specified limits on an Agilent J&W DB-Select 624UI for USP <467>, 30 m x 0.32 mm, 1.80 µm column

1. 1,1-Dichloroethene
2. 1,1,1-Trichloroethane
3. Carbon tetrachloride
4. Benzene
5. 1,2-Dichloroethane

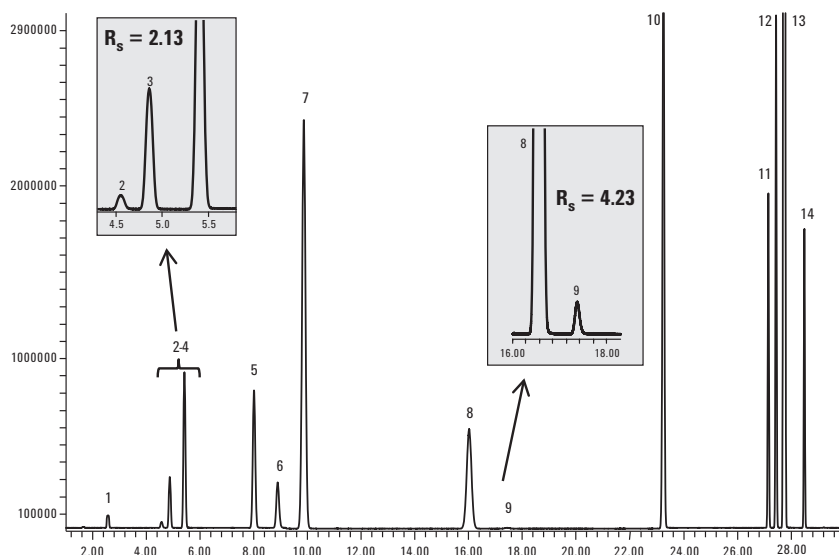

FID trace of Class 2A solvent standard at USP <467> specified limits on an Agilent J&W DB-Select 624UI for USP <467>, 30 m x 0.32 mm, 1.80 µm column

1. Methanol
2. Acetonitrile
3. Dichloromethane
4. *trans*-1,2-Dichloroethane
5. *cis*-1,2-Dichloroethane
6. Tetrahydrofuran
7. Cyclohexane
8. Methylcyclohexane
9. 1,4-Dioxane
10. Toluene
11. Chlorobenzene
12. Ethylbenzene
13. *m/p*-Xylene
14. *O*-Xylene

Part Number Index

0100-0057.....	113, 131	05971-20143.....	154-155, 157	115-3352E.....	429	121-5522.....	126, 300
0100-0161.....	127	05971-60571.....	163, 176	115-3352PT.....	429	121-5522E.....	300
0100-0549.....	169	05971-80103.....	165	115-3412.....	426	121-5522LTM.....	300, 456
0100-1324.....	41	05980-20018.....	146, 176-177	115-3422.....	426	121-5522UI.....	266, 290
0100-1325.....	41	05980-60051.....	146, 174	115-3432.....	426	121-5522UILTM.....	290, 457
0100-1326.....	41	05988-20066.....	38, 40, 148, 172	115-3432E.....	426	121-5523.....	300
0100-1331.....	41	07673-20570.....	51	115-3432PT.....	421, 426	121-5523LTM.....	300, 456
0100-1332.....	41, 91, 99, 105, 113, 129, 131, 136-137	07673-40180.....	50	115-34H2.....	426	121-5523UI.....	290
0100-1342.....	41	07673-60840.....	51	115-3532.....	433	121-5523UILTM.....	290, 457
0100-1344.....	41	0960-0897.....	162	115-3532PT.....	278-279, 421, 433	121-5542.....	300
0100-1365.....	41	1000-1437.....	113-115, 131	115-3552.....	433	121-5621.....	394
0100-1375.....	41, 85, 87	1000-1438.....	114	115-3552PT.....	279, 285, 421, 433	121-5622.....	394
0100-1378.....	41	112-2032.....	361	115-4912.....	366	121-7012.....	352
0100-1379.....	41	112-2032LTM.....	361, 453	115-4912E.....	366	121-7012LTM.....	352, 458
0100-1381.....	41	112-2112.....	442	121-0122.....	295	121-7013LTM.....	458
0100-1389.....	80	112-2132.....	442	121-0122LTM.....	295, 455	121-7022.....	352
0100-1597.....	169	112-2133.....	442	121-0122UI.....	289	121-7022LTM.....	352, 458
0100-2138.....	122-123	112-2133LTM.....	442, 453	121-0713.....	349	121-7023.....	352
0100-2430.....	122-123	112-2162.....	442	121-0722.....	266, 349	121-7023LTM.....	352, 458
0100-2594.....	53	112-2532.....	411	121-0722LTM.....	349, 454	121-7042.....	352
0100-2595.....	53	112-2532E.....	411	121-1012.....	319	121-7042E.....	352
0101-0282.....	48	112-2532LTM.....	411, 453	121-1012E.....	319	121-7043.....	352
0101-0299.....	48	112-2562.....	411	121-1012LTM.....	319, 453	121-9627.....	390
0101-0300.....	48	112-5432.....	444	121-1013.....	319	121-9723.....	264, 276, 292, 392
0101-0301.....	48	112-5462.....	444	121-1013LTM.....	319, 453	122-0112.....	295
0101-0302.....	48	112-6632.....	410	121-101A.....	319	122-0112E.....	295
0101-0303.....	48	112-6632LTM.....	410, 453	121-101ALTM.....	319, 453	122-0112LTM.....	295, 455
0101-0304.....	48	112-8837.....	126, 405	121-1022.....	319	122-0112UI.....	289
0101-0355.....	48	112-8837E.....	405	121-1022E.....	319	122-0131.....	295
0101-0403.....	48	112-8837LTM.....	405, 460	121-1022LTM.....	319, 453	122-0132.....	126, 273, 295
0101-0532.....	48	112-8867.....	405	121-1023.....	319	122-0132E.....	295
0101-0584.....	47	112-8867E.....	405	121-1023LTM.....	319, 453	122-0132LTM.....	455
0101-0585.....	47	112-88A7.....	405	121-1043.....	319	122-0132UI.....	289
0101-0633.....	48	112-88A7E.....	405	121-1222.....	304	122-0132UIE.....	289
0101-0636.....	47	113-2032.....	361	121-1232.....	304	122-0162.....	295
0101-0637.....	47	113-2132.....	442	121-1324.....	267, 271, 275, 284, 398	122-0162UI.....	289
0101-0638.....	47	113-2132LTM.....	442, 453	121-1324E.....	398	122-0212.....	344
0101-0639.....	47	113-2133.....	442	121-1324LTM.....	398, 457	122-0232.....	344
0101-0666.....	48	113-2532.....	411	121-1324UI.....	275, 291, 397	122-0232E.....	344
0101-0667.....	48	113-2532E.....	411	121-1524.....	267, 271, 275, 284, 400	122-0233.....	267, 344
0101-0946.....	47	113-3032.....	444	121-1524LTM.....	400, 457	122-0334UI.....	292, 416
0101-0947.....	47	113-3112.....	437	121-1544.....	400	122-0364UI.....	292, 416
0101-0948.....	48	113-3132.....	437	121-1544E.....	400	122-0712.....	349
0101-0954.....	48	113-3133.....	437	121-1722.....	339	122-0713.....	349
0101-0955.....	48	113-3133LTM.....	437, 458	121-1722LTM.....	339, 454	122-0713LTM.....	349, 454
0101-0956.....	48	113-3162.....	437	121-1723.....	339	122-0731.....	349
0101-0957.....	48	113-3432.....	426	121-2223.....	345	122-0732.....	126, 266, 268, 349
0101-1001.....	48	113-3432E.....	426	121-2323.....	342	122-0732E.....	349
0101-1472.....	47	113-3432LTM.....	426, 458	121-3822.....	306	122-0732LTM.....	349, 454
0101-1473.....	47	113-4302.....	434	121-3822UI.....	291	122-0733.....	274, 349
03396-61010.....	127	113-4312.....	434	121-4722.....	308	122-0733E.....	349
0515-0680.....	113, 115, 131	113-4332.....	283, 434	121-4722LTM.....	308, 454	122-0733LTM.....	349, 454
0515-0683.....	114	113-4362.....	278, 434	121-5012.....	329	122-0761.....	349
0515-2495.....	120, 133	113-5432.....	444	121-5012E.....	329	122-0762.....	349
0515-2712.....	128	113-6632.....	410	121-5012LTM.....	329, 455	122-0763.....	349
0515-2726.....	120, 133	113-6632LTM.....	410, 453	121-5013.....	329	122-0763E.....	349
0535-0071.....	157	115-2132.....	282, 442	121-5013LTM.....	329, 455	122-0766.....	349
05890-61525.....	50	115-2132LTM.....	442, 453	121-5022.....	329	122-1011.....	319
05890-80660.....	49	115-3113.....	437	121-5022E.....	329	122-1012.....	319
05921-21170.....	39-40	115-3133.....	437	121-5022LTM.....	329, 455	122-1012LTM.....	319, 453
05970-60045.....	182	115-3133LTM.....	437, 458	121-5023.....	329	122-1013.....	319
05971-20134.....	151, 154	115-3332.....	429	121-5023LTM.....	329, 455	122-1022.....	319
		115-3352.....	429	121-5042.....	329	122-1022LTM.....	319, 453

122-1022LTM	319, 453	122-1811	364	122-5032E	330	122-6832	269-270, 388
122-1031	319	122-1831	364	122-5032LTM	330, 456	122-7012	353
122-1032	319	122-1831LTM	364, 454	122-5033	272, 330	122-7012E	353
122-1032E	319	122-1932	337	122-5033E	330	122-7012LTM	353, 458
122-1032G	445	122-1962	337	122-5033LTM	330, 456	122-7013	272, 353
122-1032LTM	319, 453	122-2032	343	122-503E	330	122-7013LTM	353, 458
122-1033	267, 319	122-2032LTM	343, 455	122-503ELTM	330, 456	122-7031	353
122-1033E	319	122-2033	343	122-5052	330	122-7032	278, 353
122-1033LTM	319, 453	122-2033LTM	343, 455	122-5061	330	122-7032E	353
122-103E	319	122-2212	345	122-5062	330	122-7032LTM	353, 458
122-103ELTM	319, 453	122-2212LTM	345, 455	122-5063	330	122-7033	272, 276, 353
122-1052	319	122-2231	345	122-506E	330	122-7033E	353
122-1061	319	122-2232	345	122-5511	300	122-7033LTM	353, 458
122-1062	319	122-2232LTM	345, 455	122-5511LTM	300, 456	122-7061	353
122-1063	319	122-2312	342	122-5512	300	122-7062	353
122-106E	278, 319	122-2331	342	122-5512LTM	300, 456	122-7062E	353
122-10A6	369	122-2332	342	122-5512UI	267, 290	122-7063	353
122-10A6E	283, 369	122-2332E	342	122-5512UILTM	290, 457	122-7063E	353
122-10AE	319	122-2332LTM	342, 455	122-5513	300	122-7332	354
122-10G3	319	122-2361	342	122-5513UI	290	122-7332E	354
122-1111	362	122-2361E	342	122-5516	300	122-7332LTM	354, 458
122-1111E	362	122-2362	285, 342	122-5522	300	122-7333	354
122-1131	362	122-2362E	342	122-5522LTM	300, 456	122-7362	354
122-1131LTM	455	122-2461	396	122-5522UI	290	122-7363	354
122-1211	304	122-2462	396	122-5522UILTM	290, 457	122-7732	283-284, 386
122-1211LTM	304, 458	122-2912	312	122-5531	300	122-7732LTM	386, 454
122-1212	304	122-2912LTM	312, 455	122-5532	265, 267-268, 272-274, 300	122-96L2	390
122-1231	304	122-2932	283, 312	122-5532E	300	122-9732	276, 292, 392
122-1232	267, 269, 272, 282-285, 304	122-2932E	312	122-5532G	445	122-9736	292, 392
122-1232G	445	122-2932LTM	312, 455	122-5532LTM	300, 456	123-0112	295
122-1232LTM	304, 458	122-2962	312	122-5532UI	266, 268-271, 290	123-0112UI	289
122-1233	304	122-3212	358	122-5532UIE	290	123-0131	295
122-1236	265, 304	122-3232	358	122-5532UILTM	290, 457	123-0132	295
122-1262	304	122-3232E	358	122-5533	267, 273-274, 276, 300	123-0132UI	289
122-1332	347	122-3232LTM	358, 457	122-5533E	300	123-0162	295
122-1332E	347	122-3233	358	122-5533G	445	123-0213	344
122-1333	267, 347	122-3262	358	122-5533LTM	300, 456	123-0232	344
122-1334	265, 269, 272, 398	122-3262E	358	122-5533UI	266, 274, 290	123-0233	344
122-1334E	398	122-3263	358	122-5533UILTM	290, 457	123-0334UI	292, 416
122-1334LTM	398, 457	122-3812	306	122-5536	270, 300	123-0364UI	292, 416
122-1334UI	265, 272, 291, 397	122-3812UI	291	122-5536E	300	123-0712	349
122-1362	347	122-3831	306	122-5536G	445	123-0712LTM	349, 454
122-1363	347	122-3832	265, 270, 274, 283-284, 306	122-5536UI	290	123-0713	349
122-1363E	347	122-3832E	306	122-5536UILTM	290, 457	123-0731	349
122-1364	265, 267, 269, 271, 284, 398	122-3832LTM	306, 455	122-5552	300	123-0732	349
122-1364E	398	122-3832UI	271, 291	122-5552UI	290	123-0732E	349
122-1364UI	265, 267, 275, 291, 397	122-3862	306	122-5561	300	123-0733	349
122-1534	265, 269, 272, 400	122-4711	308	122-5562	300	123-0733E	349
122-1534LTM	400, 457	122-4711LTM	308, 454	122-5562E	300	123-0753	349
122-1564	265, 267, 269, 271-272, 275, 284, 400	122-4712	308	122-5562G	445	123-0762	349
122-1564E	400	122-4712LTM	308, 454	122-5562UI	270-271, 276, 290	123-0763	349
122-1712	339	122-4731	308	122-5563	300	123-0763E	349
122-1713	339	122-4732	270, 273, 308	122-5563UI	290	123-100ALTM	453
122-1713E	339	122-4732E	308	122-5631	394	123-1011	320
122-1731	339	122-4732LTM	308, 454	122-5631G5	445	123-1011LTM	320, 453
122-1731E	339	122-4762	308	122-5632	394	123-1012	320
122-1732	339	122-5002LTM	456	122-5633	394	123-1012LTM	320, 453
122-1732E	339	122-5011	330	122-5661	394	123-1013	320
122-1732LTM	339, 454	122-5012	283, 330	122-5711	363	123-1014	320
122-1733	339	122-5012LTM	330, 456	122-5711E	363	123-1015	320
122-1762	339	122-5013	330	122-5711LTM	363, 456	123-1015LTM	320, 453
122-1801	364	122-501E	330	122-5731	363	123-1022	283, 320
122-1801LTM	364, 454	122-5022	330	122-5731LTM	363, 456	123-1026	320
		122-5031	330	122-6432	443	123-1027	320
		122-5032	284, 330	122-6462	443	123-102F	320
						123-1031	320

PART NUMBER INDEX

123-1032	320	123-1861	364	123-5512LTM	300, 456	123-BD01	378
123-1032LTM	320, 453	123-1932	337	123-5513	300	123-BD11	378
123-1033	282, 320	123-1933	337	123-5513LTM	300, 456	123-BD34	378
123-1033E	320	123-1933E	337	123-5526	300	124-0034	402
123-1033LTM	320, 453	123-1933LTM	337, 455	123-5531	300	124-1032	320
123-1034	320	123-2032	343	123-5532	270, 273, 300	124-1034	320
123-1035	320	123-2033	343	123-5532E	300	124-1334	280, 283, 398
123-1035LTM	320, 453	123-2232	345	123-5532UI	275, 290	124-1334LTM	398, 457
123-103B	320	123-2232E	345	123-5532UIE	290	124-1374	269, 398
123-103BLTM	320, 453	123-2332	342	123-5533	276, 300	124-1534	281, 400
123-103E	270, 274, 320	123-2332E	342	123-5533LTM	300, 456	124-1574	281, 400
123-103ELTM	320, 453	123-2362	342	123-5533UI	290	124-3232	282, 358
123-1052	320	123-2932	312	123-5536	270, 274, 300	125-032	331
123-1055	320	123-3212	358	123-5536LTM	300, 456	124-5037	331
123-1056	320	123-3223	358	123-5536UI	290	124-7032	280, 282, 353
123-105C	320	123-3232	281, 358	123-5561	300	125-0212	279, 344
123-105F	320	123-3232E	358	123-5562	300	125-0232	344
123-1061	320	123-3232LTM	358, 457	123-5563	300	125-0232LTM	344, 455
123-1062	320	123-3233	358	123-5563UI	290	125-0334UI	292, 416
123-1062E	320	123-3233LTM	358, 457	123-5566	300	125-0712	349
123-1063	320	123-3234	358	123-5631	394	125-0712E	349
123-1063E	320	123-3234LTM	358, 457	123-5632	394	125-0712LTM	349, 454
123-1064	320	123-3253	358	123-5701	363	125-0731	349
123-1064E	320	123-3262	358	123-5701LTM	363, 456	125-0732	349
123-1065	320	123-3263	358	123-5711	363	125-0732E	349
123-1065E	320	123-3264	358	123-5711E	363	125-0732G	445
123-106B	320	123-3812	306	123-5731	363	125-0733	349
123-106BE	320	123-3832	273-274, 306	123-5731E	363	125-0737	349
123-106E	279, 320	123-3832E	306	123-6133	443	125-0762	349
123-106G	320	123-3832UI	266, 268, 270, 275, 291	123-6412	443	125-0762E	349
123-1102LTM	455	123-4712	308	123-6432	443	125-1002	321
123-1111	362	123-4732	270, 308	123-7012	353	125-1005	321
123-1111LTM	362, 455	123-4732LTM	308, 454	123-7012LTM	353, 458	125-1005LTM	321, 454
123-1131	362	123-5003LTM	456	123-7013	353	125-100B	321
123-1131E	362	123-500E	330	123-7013LTM	353, 458	125-1011	321
123-1232	304	123-500ELTM	330, 456	123-7031	353	125-1011E	321
123-1236	266, 268, 270, 273, 304	123-5011	330	123-7032	353	125-1011LTM	321, 454
123-1262	304	123-5011LTM	330, 456	123-7032E	353	125-1012	321
123-1332	347	123-5012	330	123-7032LTM	353, 458	125-1012E	321
123-1333	347	123-5012E	330	123-7033	353	125-1012LTM	321, 454
123-1334	398	123-5012LTM	330, 456	123-7033E	353	125-1014	321
123-1334E	398	123-5013	330	123-7033LTM	353, 458	125-1015	321
123-1334LTM	398, 457	123-5013E	330	123-7062	281, 353	125-1015LTM	321, 454
123-1334UI	291, 397	123-5013LTM	330, 456	123-7063	353	125-1017	321
123-1363	347	123-5022	330	123-7063E	353	125-101J	321
123-1364	398	123-5022LTM	330, 456	123-7312	354	125-101K	321
123-1364E	398	123-5026	330	123-7314	354	125-1025	321
123-1364UI	275, 291, 397	123-502D	330	123-7332	354	125-1025LTM	321, 454
123-1464	402	123-502F	330	123-7333	354	125-102J	321
123-1534	283, 400	123-5031	330	123-7334	354	125-1032	321
123-1564	400	123-5032	330	123-7334LTM	354, 458	125-1032LTM	321, 454
123-1632	284, 403	123-5032E	330	123-7354	354	125-1034	283, 321
123-1711	339	123-5032LTM	330, 456	123-7354E	354	125-1034E	321
123-1712	339	123-5033	282, 330	123-7362	354	125-1034LTM	321, 454
123-1713	339	123-5033E	330	123-7363	354	125-1035	278, 284-285, 321
123-1730	266, 282, 388	123-503B	330	123-7364	285, 354	125-1035E	321
123-1730LTM	388, 457	123-503BLTM	330, 456	123-7722	386	125-1035LTM	321, 454
123-1731	339	123-503E	330	123-7732	386	125-1037	321
123-1732	339	123-503ELTM	330, 456	123-8232	265-266, 268, 273, 275, 384	125-1039	321
123-1732E	339	123-5052	330	123-8336	265-266, 268, 273, 275, 384	125-103B	321
123-1732LTM	339, 454	123-5053	330	123-9134	414	125-103J	321
123-1733	339	123-5056	330	123-9134LTM	414, 457	125-103JLTM	321, 454
123-1733E	339	123-5062	330	123-9234	414	125-103K	321
123-1762	339	123-5063	330	123-9234E	414	125-103KE	321
123-1811	364	123-5511	300			125-103KLTM	321, 454
123-1831	364	123-5512	300			125-1055	321

125-1062	321	125-3237	358	127-0113	295	128-3812	306
125-1062E	321	125-3262	358	127-0122	295	128-3822	306
125-1064	321	125-32H2	358	127-0123	295	128-5012	330
125-1065	321	125-3832	306	127-0123LTM	295, 455	128-5022	330
125-1065E	321	125-3837	306	127-0722	349	128-5022LTM	330, 455
125-106J	321	125-5012	273, 331	127-0723	349	128-5052	330
125-106JE	321	125-5012E	331	127-100A	318	128-50H7	330
125-10B5	321	125-5012LTM	331, 456	127-100ALTM	318, 453	128-5512	300
125-10H5	321	125-5017	331	127-1012	318	128-5522	300
125-10HB	321	125-501J	331	127-1012E	318	128-5522LTM	300, 456
125-10HBE	321	125-501K	331	127-1013	318	128-5552	300
125-10HBLTM	321, 454	125-5025	331	127-1013E	318	128-7022	352
125-1131	362	125-5032	282, 331	127-1013LTM	318, 453	128-7032	352
125-1212	304	125-5032E	331	127-1022	318	128-7032LTM	352, 458
125-1232	304	125-5032LTM	331, 456	127-1022E	318	128-7052	352
125-1312	347	125-5034	331	127-1023	318	128-7323	354
125-1314	398	125-5035	331	127-1023LTM	318, 453	128-8522	415
125-1332	347	125-5035E	331	127-1043	318	12A-1015	319
125-1333	347	125-5035LTM	331, 456	127-1046	318	12A-1015LTM	319, 453
125-1333LTM	347, 454	125-5037	331	127-1046E	318	12A-5015	329
125-1334	279-280, 285, 398	125-503B	331	127-1712	339	12A-5015LTM	329, 455
125-1334E	398	125-503D	331	127-1712LTM	339, 454	1300502506	54
125-1334G5	445	125-503J	331	127-1713	339	14-3823-000	62
125-1334LTM	398, 457	125-503K	331	127-1722	339	14-6855-000	62
125-1334UI	291, 397	125-5062	331	127-2222	345	14-8911-003	63
125-1364	398	125-5065	331	127-3212	358	1400-0015	120, 133
125-1364E	398	125-5065E	331	127-3212LTM	358, 457	145-1001	285, 370, 417-418
125-1374	269, 398	125-50HB	331	127-32H2	358	145-1009	370, 417-418
125-1374E	398	125-5512	300	127-32H2LTM	358, 457	145-1011	283, 418
125-1374UI	291, 397	125-5532	300	127-5012	329	145-1032	418
125-14A4	402	125-5532LTM	300, 456	127-5012E	329	145-1334	419
125-1704	339	125-5537	300	127-5012LTM	329, 455	145-2814	417-418
125-1711	339	125-5537G	445	127-5013	329	145-7032	419
125-1712	339	125-553J	300	127-5013LTM	329, 455	1460-1160	113, 115, 131
125-1712LTM	339, 454	125-553JLTM	300, 456	127-501E	329	1460-1914	126
125-1713	339	125-6837	272-273, 275, 388	127-501ELTM	329, 455	1500334701	218, 220
125-1713LTM	339, 454	125-7012	353	127-501N	329	1535-4952	47, 52
125-1717	339	125-7012E	329	127-5022	329	1535-4954	47
125-1730	282, 388	125-7017	353	127-5023	329	160-1010-1	464
125-1731	339	125-7031	353	127-7012	352	160-1010-10	464
125-1732	339	125-7031LTM	353, 458	127-7012E	352	160-1010-5	464
125-1732E	339	125-7032	279-280, 353	127-7012LTM	352, 458	160-2200-10	468
125-1732LTM	339, 454	125-7032E	353	127-7013	352	160-2200-5	468
125-1733	339	125-7032LTM	353, 458	127-7013LTM	352, 458	160-2205-1	464
125-1737	339	125-7037	353	127-7022	352	160-2205-10	464
125-1762	339	125-7062	353	127-7022LTM	352, 458	160-2205-5	464
125-1912	337	125-7062E	353	127-7023	352	160-2250-10	468
125-1932	337	125-7312	354	127-7023E	352	160-2250-5	468
125-1932LTM	337, 455	125-7314	354	127-7023FF	353	160-2255-1	465
125-1937	337	125-7332	354	127-7023LTM	352, 458	160-2255-10	465
125-2032	343	125-7332E	354	128-0112	295	160-2255-30	465
125-2212	273, 345	125-7333	354	128-0122	295	160-2255-5	43, 53, 465
125-2232	345	125-7333LTM	354, 458	128-0122E	295	160-2320-10	468
125-2237	345	125-7334	354	128-0122LTM	295, 455	160-2320-5	468
125-2312	342	125-7334E	354	128-1012	319	160-2325-1	465
125-2332	342	125-7362	354	128-1012LTM	319, 453	160-2325-10	465
125-2814	279-280, 283, 370	125-7732	386	128-1022	319	160-2325-30	465
125-2814E	370	125-9134	414	128-1022LTM	319, 453	160-2325-5	43, 53, 465
125-2814LTM	370, 455	125-9134E	414	128-1034	319	160-2455-1	465
125-3212	358	125-9234	414	128-1052	319	160-2455-10	465
125-3217	358	126-1012	318	128-1056	284, 369	160-2455-5	53, 465
125-3217LTM	358, 457	126-1013	318	128-1222	304	160-2530-10	468
125-3231	358	126-7012	352	128-1314LTM	457	160-2530-5	468
125-3232	358	126-7013	352	128-1324	398	160-2535-1	465
125-3232E	358	127-0112	295	128-1324E	398	160-2535-10	465
125-3233	358	127-0112LTM	295, 455	128-1324LTM	398, 457	160-2535-30	465

PART NUMBER INDEX

160-2535-5.....	43, 53, 465	19091A-115.....	327	19091J-216E.....	333	19091N-216E.....	355
160-2610-10.....	468	19091B-002.....	328	19091J-231.....	333	19091N-231.....	272, 355
160-2610-5.....	468	19091B-005.....	328	19091J-233.....	272, 333	19091N-233.....	276, 355
160-2615-1.....	464	19091B-012.....	328	19091J-233LTM.....	333, 459	19091N-233E.....	355
160-2615-10.....	464	19091B-012E.....	328	19091J-236.....	333	19091N-236.....	355
160-2615-5.....	464	19091B-015.....	328	19091J-313.....	333	19091N-331.....	355
160-2625-1.....	43, 464	19091B-101.....	328	19091J-330.....	333	19091N-577.....	355
160-2625-10.....	43, 464	19091B-101LTM.....	328, 461	19091J-330LTM.....	333, 459	19091N-577E.....	355
160-2625-5.....	43, 464	19091B-102.....	328	19091J-333.....	333	19091N-577LTM.....	355, 460
160-2634-10.....	468	19091B-102E.....	328	19091J-411.....	333	19091P-K15.....	428
160-2634-5.....	468	19091B-102LTM.....	328, 461	19091J-411LTM.....	333, 459	19091P-K15E.....	428
160-2635-1.....	464	19091B-105.....	328	19091J-413.....	333	19091P-K15PT.....	421, 428
160-2635-10.....	464	19091B-105E.....	328	19091J-413E.....	333	19091P-K33.....	428
160-2635-5.....	464	19091B-112.....	328	19091J-413LTM.....	333, 459	19091P-K33LTM.....	428, 461
160-2644-10.....	468	19091B-112LTM.....	328, 461	19091J-416.....	333	19091P-M15.....	434
160-2644-5.....	468	19091B-115.....	328	19091J-431.....	333	19091P-M15E.....	434
160-2650-10.....	468	19091B-115E.....	328	19091J-431E.....	333	19091P-MS4.....	438
160-2650-5.....	468	19091F-102.....	359	19091J-433.....	333	19091P-MS4E.....	438
160-2655-1.....	464	19091F-102E.....	359	19091J-433E.....	333	19091P-MS7.....	438
160-2655-10.....	464	19091F-102LTM.....	359, 460	19091J-433LTM.....	333, 459	19091P-MS7LTM.....	438, 461
160-2655-5.....	464	19091F-105.....	359	19091J-436.....	333	19091P-MS8.....	438
160-2660-5.....	468	19091F-105E.....	359	19091J-436E.....	333	19091P-MS8LTM.....	438, 461
160-2815-5.....	465	19091F-112.....	359	19091J-577.....	332	19091P-Q03.....	425
160-2825-5.....	465	19091F-112E.....	359	19091J-577E.....	332	19091P-Q03LTM.....	425, 461
160-2845-10.....	465	19091F-112LTM.....	359, 460	19091J-577LTM.....	332, 459	19091P-Q03PT.....	421, 425
160-2845-5.....	465	19091F-115.....	359	19091L-101.....	340	19091P-Q04.....	425
160-2855-10.....	465	19091F-115E.....	359	19091L-113.....	340	19091P-Q04E.....	425
160-2855-5.....	465	19091F-413.....	359	19091L-113E.....	340	19091P-Q04LTM.....	425, 461
160-2865-10.....	465	19091F-433.....	359	19091L-133.....	340	19091P-Q04PT.....	421, 425
160-2865-5.....	465	19091F-433E.....	359	19091L-330.....	340	19091P-S12.....	432
160-4535-5.....	53, 469	19091F-433LTM.....	359, 460	19091L-330LTM.....	340, 459	19091P-S12LTM.....	432, 461
169-0013-HSP.....	55	19091G-113.....	338	19091L-333.....	340	19091P-S12PT.....	421, 432
18596-40015.....	51	19091G-131.....	338	19091L-413.....	340	19091P-S15.....	432
18710-20119.....	96, 119, 130, 135	19091G-131E.....	338	19091L-413E.....	340	19091P-S15E.....	432
18710-60170.....	105, 136	19091G-131LTM.....	338, 459	19091L-416.....	340	19091P-S15PT.....	421, 432
18711-60060.....	105, 136	19091G-133.....	338	19091L-431.....	340	19091P-S33.....	432
18713-60040.....	102, 132	19091G-213.....	338	19091L-431LTM.....	340, 459	19091P-U04.....	428
18713-60050.....	102, 132	19091G-B133.....	411	19091L-433.....	340	19091P-U04E.....	428
18740-20800.....	75, 139	19091G-B213.....	411	19091L-433LTM.....	340, 459	19091P-U04LTM.....	428, 461
18740-20880.....	67, 139	19091G-B233.....	411	19091N-013.....	355	19091R-303.....	401
18740-60830.....	74, 76, 138	19091G-B233E.....	411	19091N-013LTM.....	355, 460	19091R-306.....	267, 271, 275, 401
18740-60835.....	74, 76, 83, 85, 90-91, 137-138	19091J-002.....	332	19091N-030.....	355	19091R-316.....	401
18740-60840.....	31	19091J-005.....	332	19091N-030LTM.....	355, 460	19091R-319.....	401
18740-80190.....	31	19091J-012.....	333	19091N-033.....	355	19091S-001.....	282-283, 367
18740-80200.....	30	19091J-015.....	333	19091N-036.....	355	19091S-001E.....	367
18740-80220.....	30	19091J-101.....	332	19091N-102.....	355	19091S-010.....	388
18789-60060.....	120, 133	19091J-102.....	332	19091N-102LTM.....	355, 460	19091S-101.....	301
18789-80070.....	96, 130	19091J-102E.....	332	19091N-105.....	355	19091S-101LTM.....	301, 460
18900-21000.....	48	19091J-105.....	332	19091N-105E.....	355	19091S-102.....	301
18900-60640.....	124	19091J-105E.....	332	19091N-111.....	355	19091S-102E.....	301
19091-21050.....	468	19091J-108.....	332	19091N-113.....	282, 355	19091S-102LTM.....	301, 460
19091-60312.....	322	19091J-112.....	333	19091N-113E.....	355	19091S-105.....	301
19091-60620E.....	464	19091J-112E.....	333	19091N-116.....	278, 285, 355	19091S-111LTM.....	460
19091-63000.....	80	19091J-113.....	333	19091N-131.....	355	19091S-112.....	301
19091A-002.....	327	19091J-113E.....	333	19091N-131E.....	355	19091S-112E.....	301
19091A-005.....	327	19091J-113LTM.....	333, 459	19091N-133.....	355	19091S-113.....	270, 274, 301
19091A-008.....	327	19091J-115.....	333	19091N-133E.....	355	19091S-133.....	267, 270, 273-274, 276, 301
19091A-012.....	327	19091J-115E.....	333	19091N-133LTM.....	355, 460	19091S-133UI.....	271, 274, 290
19091A-015.....	327	19091J-202.....	332	19091N-136.....	278, 355	19091S-133UILTM.....	290, 460
19091A-101.....	327	19091J-205.....	332	19091N-136E.....	355	19091S-139.....	395
19091A-102.....	327	19091J-212.....	333	19091N-202.....	355	19091S-213.....	276, 282, 301
19091A-102E.....	327	19091J-213.....	333	19091N-205.....	355	19091S-213UI.....	290
19091A-105.....	327	19091J-213E.....	333	19091N-205E.....	355	19091S-213UILTM.....	290, 460
19091A-108.....	327	19091J-215.....	333	19091N-213.....	355	19091S-231.....	301
19091A-112.....	327	19091J-215E.....	333	19091N-213E.....	355	19091S-233.....	301
		19091J-216.....	333	19091N-216.....	281-282, 355		

19091S-233E.....	301	19091W-012LTM.....	361, 459	19095F-121LTM.....	359, 460	19095S-100.....	323
19091S-233UI.....	266, 290	19091W-015.....	361	19095F-123.....	359	19095S-100E.....	323
19091S-233UILTM.....	290, 460	19091W-015E.....	361	19095F-123E.....	359	19095S-200.....	369
19091S-313.....	301	19091W-102.....	361	19095F-123LTM.....	359, 460	19095S-205.....	369
19091S-331.....	301	19091W-105.....	361	19095J-023.....	333	19095V-420.....	416
19091S-331LTM.....	301, 460	19091Y-012.....	441	19095J-023E.....	333	19095V-420E.....	416
19091S-333.....	301	19091Y-012E.....	441	19095J-121.....	333	19095V-420LTM.....	416, 460
19091S-336.....	301	19091Y-015.....	441	19095J-121E.....	333	19095W-121.....	361
19091S-413.....	285, 301	19091Y-102.....	441	19095J-121LTM.....	333, 459	19095W-123.....	361
19091S-413E.....	301	19091Z-002.....	322	19095J-123.....	333	19095Z-020.....	279, 323
19091S-413LTM.....	301, 460	19091Z-002LTM.....	322, 459	19095J-123E.....	333	19095Z-021.....	279, 323
19091S-413UI.....	275, 290	19091Z-005.....	322	19095J-321.....	273, 333	19095Z-021E.....	323
19091S-413UILTM.....	290, 460	19091Z-008.....	322	19095J-323.....	333	19095Z-021LTM.....	323, 459
19091S-416.....	301	19091Z-012.....	323	19095J-323E.....	333	19095Z-023.....	323
19091S-431.....	301	19091Z-012LTM.....	323, 459	19095J-621.....	333	19095Z-023E.....	323
19091S-431LTM.....	301, 460	19091Z-015.....	323	19095J-623.....	285, 333	19095Z-023LTM.....	323, 459
19091S-431UI.....	267, 290	19091Z-102.....	322	19095J-623E.....	333	19095Z-121.....	323
19091S-431UILTM.....	290, 460	19091Z-102E.....	322	19095L-021.....	340	19095Z-121E.....	323
19091S-433.....	267, 283-284, 301	19091Z-105.....	322	19095L-021LTM.....	340, 459	19095Z-121LTM.....	323, 459
19091S-433E.....	301	19091Z-112.....	323	19095L-023.....	340	19095Z-123.....	323
19091S-433LTM.....	301, 460	19091Z-112E.....	323	19095L-023E.....	340	19095Z-123E.....	323
19091S-433UI.....	266, 268, 271, 290	19091Z-115.....	281, 285, 323	19095L-121.....	441	19095Z-123LTM.....	323, 459
19091S-433UIE.....	290	19091Z-115E.....	323	19095L-523.....	340	19095Z-220.....	323
19091S-433UILTM.....	290, 460	19091Z-202.....	322	19095N-121.....	355	19095Z-221.....	323
19091S-436.....	301	19091Z-202LTM.....	322, 459	19095N-123.....	355	19095Z-221E.....	323
19091S-436E.....	301	19091Z-205.....	322, 367	19095N-123E.....	355	19095Z-321.....	323
19091S-436UI.....	290	19091Z-205E.....	367	19095N-123LTM.....	355, 460	19095Z-323.....	323
19091S-577.....	301	19091Z-211.....	323	19095N-126.....	355	19095Z-323E.....	323
19091S-577LTM.....	301, 460	19091Z-212.....	323	19095P-K23.....	428	19095Z-421.....	323
19091S-577UI.....	266, 290	19091Z-213.....	323	19095P-K23LTM.....	428, 461	19095Z-423.....	323
19091S-577UILTM.....	290, 460	19091Z-213E.....	323	19095P-K23PT.....	278, 421, 428	19095Z-423E.....	323
19091S-602.....	296	19091Z-215.....	323	19095P-K25.....	428	19095Z-621.....	282, 323
19091S-602E.....	296	19091Z-216.....	323	19095P-K25E.....	428	19095Z-623.....	323
19091S-612.....	296	19091Z-216E.....	323	19095P-K25PT.....	285, 421, 428	19095Z-623E.....	323
19091S-612UI.....	289	19091Z-231.....	322	19095P-M23.....	434	19095Z-623LTM.....	323, 459
19091S-633.....	296	19091Z-233.....	322	19095P-M25.....	434	19095Z-626.....	285, 323
19091S-633LTM.....	296, 459	19091Z-233E.....	322	19095P-M25PT.....	421, 434	19095Z-627.....	283, 323
19091S-633UI.....	289	19091Z-236.....	284-285, 322	19095P-MS0.....	278, 438	1909BD-113.....	378
19091S-677.....	296	19091Z-236E.....	322	19095P-MS0E.....	438	19231-20910.....	98
19091S-677LTM.....	296, 459	19091Z-313.....	323	19095P-MS0LTM.....	438, 461	19231-20940.....	98
19091S-677UI.....	289	19091Z-313LTM.....	323, 459	19095P-MS5.....	438	19231-20980.....	98
19091S-713.....	280, 283, 296	19091Z-331.....	322	19095P-MS6.....	438	19231-21050.....	128
19091S-713LTM.....	296, 459	19091Z-333.....	322	19095P-MS6E.....	438	19231-21060.....	98
19091S-713UI.....	289	19091Z-411.....	323	19095P-MS9.....	278, 438	19231-60680.....	98
19091S-733.....	296	19091Z-413.....	323	19095P-Q03.....	425	19231-80520.....	99, 121, 129, 134
19091S-733E.....	296	19091Z-413E.....	323	19095P-Q03E.....	425	19231-80530.....	99, 121, 129, 134
19091S-733LTM.....	296, 459	19091Z-416.....	323	19095P-Q03LTM.....	425, 461	19233-20755.....	102, 132
19091S-733UI.....	289	19091Z-431.....	322	19095P-Q03PT.....	278, 282, 284, 421, 425	19234-60700.....	90, 99, 102, 129, 132, 137
19091S-833.....	296	19091Z-433.....	322	19095P-Q04.....	425	19234-60715.....	90, 99, 102, 129, 132, 137
19091S-833LTM.....	296, 459	19091Z-433E.....	322	19095P-Q04E.....	425	19234-60720.....	91, 121, 134
19091S-911.....	296	19091Z-436.....	322	19095P-Q04LTM.....	425, 461	19243-00070.....	75, 139
19091S-911UI.....	289	19091Z-513.....	323	19095P-Q04PT.....	283-284, 421, 425	19243-80530.....	90-91, 137
19091S-913.....	296	19091Z-513E.....	323	19095P-S21.....	432	19243-80540.....	90-91, 137
19091S-913E.....	296	19091Z-530.....	283, 322, 367	19095P-S23.....	432	19243-80570.....	90-91, 137
19091S-913UI.....	289	19091Z-530E.....	322, 367	19095P-S23PT.....	278, 421, 432	19244-80540.....	90-91, 137
19091S-916.....	296	19091Z-577.....	322	19095P-S25.....	432	19244-80560.....	96, 119, 130, 135
19091S-931.....	296	19091Z-577E.....	322	19095P-S25E.....	432	19244-80610.....	99, 121, 129, 134
19091S-931UI.....	289	19091Z-613.....	284, 323	19095P-S25PT.....	281, 421, 432	19244-80620.....	96, 119, 130, 135
19091S-933.....	296	19091Z-613LTM.....	323, 459	19095P-U03.....	428	19245-20510.....	80
19091S-933E.....	296	19091Z-713.....	323	19095P-U04.....	428	19245-20515.....	80
19091S-933LTM.....	296, 459	19091Z-713E.....	323	19095P-U04E.....	428	19245-20525.....	80
19091S-933UI.....	289	19091Z-713LTM.....	323, 459	19095P-U04LTM.....	428, 461	19245-20580.....	80
19091S-936.....	296	19091Z-716.....	323	19095P-U04PT.....	284, 421, 428	19245-20780.....	80
19091S-936E.....	296	19095F-120.....	359	19095R-420.....	401	19245-40050.....	80
19091W-012.....	361	19095F-120E.....	359	19095R-429.....	401		
19091W-012E.....	361	19095F-121.....	359				

PART NUMBER INDEX

19245-60760.....	80	222-4712LTM.....	463	390607400.....	220	393050492.....	180
19245-80521.....	80	222-4732LTM.....	463	390607401.....	220	393050493.....	180
19251-60540.....	29	222-5512LTM.....	462	390607900.....	220	393053501.....	180
19251-80680.....	126	222-5512UILTM.....	462	390812700.....	216-217	393053502.....	179-180
19256-00090.....	114	222-5532LTM.....	462	390820601.....	213-217	393060191.....	180
19256-00200.....	114	222-5532UILTM.....	462	390842300.....	213-217	393065201.....	182
19256-00320.....	114	222-5711LTM.....	462	391821100.....	215	393082491.....	179
19256-20690.....	113, 115, 131	222-5731LTM.....	462	391866306.....	215	393101291.....	179
19256-20705.....	113, 131	222-7013LTM.....	463	391866308.....	213	393112601.....	179
19256-20900.....	114	222-7033LTM.....	463	391867600.....	214	393112702.....	179
19256-20910.....	113-115, 131	2302533140.....	54	392017401.....	179	393113001.....	179
19256-21140.....	113, 131	2307230001.....	54	392027300.....	179-180	393161001.....	179
19256-21150.....	113, 131	2307232901.....	54	392030500.....	179-180	393164493.....	179
19256-60510.....	114	232-2790010-EHS.....	55	392035300.....	179-180	393167593.....	179
19256-60700.....	113, 131	232-2790012-EHS.....	55	392043700.....	180	393171201.....	179
19256-60750.....	115	2321700003.....	54-55	392047100.....	180	393175101.....	179
19256-60800.....	113, 131	2321700004.....	54-55	392511901.....	219-220	393706201.....	146, 174
19256-80000.....	114	2322590004.....	54-55	392512800.....	220	393847701.....	179
19256-80010.....	113-115, 131	2322590005.....	54-55	392513800.....	219	394958700.....	218
19256-80640.....	113, 115, 131	2322700011.....	54	392514300.....	219-220	394966601.....	214-215
19258-20830.....	55	2710100200.....	179	392514500.....	219-220	410105017.....	54
19258-20870.....	55	2710100400.....	179	392515101.....	219	4177-0607.....	53
19298-60500.....	124	2710100500.....	179	392515102.....	219	430-1020.....	126
19301-60660.....	128	2722990700.....	180	392515103.....	219	450-1000.....	86, 90, 137
19320-80625.....	80	2735000500.....	179	392515104.....	219	480-0003.....	86, 90, 137
19325-60660.....	48	2740236100.....	220	392515105.....	219	500-2114.....	37
19354-60510.....	124	2740292400.....	219	392515500.....	219-220	500-2118.....	37
200-0010.....	482	2740928202.....	220	392517100.....	219	5020-8292.....	40
200-0032.....	482	29091N-133LTM.....	463	392517600.....	219	5020-8293.....	40
200-0070.....	482	29091N-577LTM.....	463	392517700.....	219	5020-8294.....	40
200-0110.....	482	29091R-303LTM.....	463	392517800.....	219	5021-7107.....	127
200-0113.....	482	29091S-431UILTM.....	462	392517901.....	219	5021-7133.....	45
200-0185.....	482	29091S-433LTM.....	463	392519200.....	219	5021-7134.....	45
200-0310.....	482	29091S-433UILTM.....	462	392543101.....	217	5021-7136.....	45
200-0370.....	482	29091S-577UILTM.....	462	392544001.....	214	5021-7137.....	45
200187500.....	218	29091S-677LTM.....	462	392544011.....	214	5021-7146.....	45
200187600.....	220	29091S-833LTM.....	462	392544391.....	219	5021-7148.....	45
200193800.....	218, 220	29091S-931LTM.....	462	392548201.....	217	5021-7164.....	45
210-3003.....	30	301-011-HSP.....	55	392548301.....	216	5021-7166.....	45
210-3003-5.....	30	301-015-HSP.....	55	392560591.....	218	5021-7168.....	45
210-4004-5.....	29	301-016-HSP.....	55	392561290.....	218	5021-7169.....	45
210-4022-5.....	29	301-017-HSP.....	55	392567111.....	179	5021-7170.....	45
2100003100.....	220	301-152-HSP.....	55	392585291.....	218	5040-4667.....	61
2100003200.....	218, 220	301-169-HSP.....	55	392585292.....	218	5040-4668.....	61
221-0122LTM.....	462	301-170-HSP.....	55	392595501.....	216-217	5040-4669.....	61
221-1324LTM.....	463	301-205-HSP.....	55	392597101.....	213	5040-4671.....	61
221-1524LTM.....	463	301-211-HSP.....	55	392597301.....	213	5060-9086.....	124
221-3822LTM.....	463	301-212-HSP.....	55	392597302.....	213	5061-5869.....	75, 139
221-4722LTM.....	463	3050-0891.....	154-155, 157	392597303.....	213	5061-5886.....	114
221-5522LTM.....	462	3050-1246.....	98	392597501.....	213	5061-5890.....	114
221-5522UILTM.....	462	3050-1301.....	154-155, 157	392599401.....	213	5061-5896.....	146, 174
221-6822LTM.....	463	3050-1374.....	157	392599411.....	213	5062-3506.....	37-38, 148, 172
222-0112LTM.....	462	3050-1375.....	154-155, 159	392599501.....	213	5062-3507.....	37, 172
222-0132LTM.....	462	3150-0602.....	84	392609901.....	23	5062-3508.....	37-38, 148
222-0732LTM.....	463	321-002-HSP.....	55	392609902.....	22	5062-3511.....	37
222-1032LTM.....	462	321-055-HSP.....	55	392609903.....	23	5062-3512.....	37
222-1111LTM.....	462	321-056-HSP.....	55	393001991.....	180	5062-3513.....	37
222-1131LTM.....	462	321-057-HSP.....	55	393010918.....	179	5062-3514.....	37-39, 148
222-1334LTM.....	463	325-062-HSP.....	55	393010920.....	179	5062-3515.....	37
222-1534LTM.....	463	325-132-HSP.....	55	393010924.....	179	5062-3516.....	37
222-2912LTM.....	463	325-185-HSP.....	55	393011391.....	180	5062-3519.....	47
222-2932LTM.....	463	35900-60670.....	127	393031501.....	180	5062-3525.....	45, 84, 86
222-3212LTM.....	463	35900-60800.....	127	393050292.....	180	5062-3538.....	37, 115
222-3232LTM.....	463	35900-60920.....	127	393050293.....	180	5062-3580.....	40
222-3812LTM.....	463	3600500001.....	55	393050392.....	180	5062-3581.....	40
222-3832LTM.....	463	3600500002.....	55	393050393.....	180	5062-3587.....	30

5062-9508.....	47	5181-7460.....	47	5182-9747.....	84-86	5183-4757.....	18, 74, 76,
5062-9509.....	47	5181-8806.....	190	5182-9748.....	84-85, 87		90-91, 137-138
5062-9510.....	47	5181-8808.....	190	5182-9749.....	84-85	5183-4757-100.....	18, 74, 76, 138
5062-9511.....	47	5181-8809.....	190	5182-9754.....	83, 86	5183-4758.....	18, 80
5067-0226.....	56	5181-8810.....	189	5182-9756.....	83, 86	5183-4759.....	19, 85, 90, 137
5067-0227.....	56	5181-8811.....	190	5182-9757.....	84, 87	5183-4759-100.....	19, 85
5067-0234.....	58	5181-8813.....	190	5182-9758.....	84, 87	5183-4760.....	19, 80
5080-5400.....	86, 90, 137, 146, 174	5181-8815.....	22, 74, 76, 90	5182-9759.....	84, 87	5183-4761.....	19, 74, 76, 85,
5080-8716.....	124	5181-8816.....	22, 74, 76	5182-9760.....	84-85, 87		90, 137-138
5080-8728-100.....	21	5181-8818.....	30	5182-9761.....	83, 86	5183-4761-100.....	19, 74, 76, 85, 138
5080-8728-50.....	21	5181-8830.....	39-40, 91, 99, 105,	5182-9762.....	83, 86	5183-4762.....	19, 80
5080-8732.....	90-91, 137		121, 129, 132, 134, 137	5182-9763.....	83, 86	5188-1181.....	169
5080-8750.....	91, 99, 127, 129, 137	5181-8833.....	22	5182-9768.....	83, 86	5188-1447.....	63
5080-8751.....	91, 99, 129, 134	5181-8836.....	53, 126	5182-9769.....	83, 86	5188-1448.....	63
5080-8752.....	121, 134	5181-8839.....	22	5182-9770.....	83, 86	5188-2717.....	68, 138
5080-8753.....	121, 134	5181-8863.....	146	5182-9775.....	84, 87	5188-2753.....	59
5080-8755.....	124	5181-8866.....	191	5182-9799.....	191	5188-2759.....	59
5080-8756.....	124	5182-0551.....	50	5183-0314.....	189	5188-5241.....	75, 139
5080-8759.....	124	5182-0773.....	63	5183-0316.....	190	5188-5245.....	113, 115, 131
5080-8761.....	124	5182-0774.....	63	5183-0318.....	87, 190	5188-5246.....	190
5080-8768.....	124	5182-0775.....	63	5183-2007.....	116, 120, 133	5188-5247.....	190
5080-8769.....	124	5182-0781.....	63	5183-2036.....	32, 85, 87	5188-5311.....	33
5080-8773.....	105	5182-0783.....	63	5183-2037.....	32, 87	5188-5312.....	40
5080-8774.....	90-91, 99, 102, 105, 121,	5182-0794.....	62	5183-2038.....	32, 87	5188-5313.....	32, 85, 87
	129, 132, 134, 136-137	5182-0795.....	62	5183-2042.....	189	5188-5314.....	40
5080-8853.....	37, 105	5182-0796.....	62	5183-2058.....	87, 190	5188-5315.....	40-41
5080-8858.....	482	5182-0830.....	190	5183-4474.....	56	5188-5316.....	124
5080-8894-100.....	21	5182-0831.....	80, 191	5183-4475.....	56	5188-5317.....	124
5080-8896-50.....	21, 90	5182-0832.....	189, 191	5183-4477.....	58	5188-5347.....	172, 182
5080-8898.....	90-91, 137	5182-0833.....	80, 191	5183-4478.....	58	5188-5348.....	172, 182
5180-4103.....	99, 105-106,	5182-0834.....	190	5183-4479.....	58	5188-5356.....	32, 85, 87
	129, 136-137	5182-0835.....	189	5183-4480.....	58	5188-5357.....	172, 182
5180-4105.....	90-91, 99, 102, 105, 121,	5182-0836.....	80, 191	5183-4641.....	127	5188-5365.....	33, 74, 76, 138, 163, 176
	129, 132, 134, 136-137	5182-0837.....	56	5183-4642.....	127	5188-5366.....	33, 68, 138
5180-4124.....	136	5182-0838.....	56	5183-4644.....	127	5188-5367.....	67, 75, 139
5180-4150.....	128	5182-0839.....	58	5183-4645.....	127	5188-5370.....	190
5180-4152.....	128	5182-0840.....	58	5183-4647.....	28-30	5188-5371.....	190
5180-4165.....	98, 128	5182-0844.....	62	5183-4691.....	29	5188-5372.....	128
5180-4168.....	33, 74, 76, 138	5182-0845.....	62	5183-4692.....	29	5188-5379.....	482
5180-4173.....	33, 74, 76, 138	5182-0846.....	62	5183-4693.....	30	5188-5392.....	59
5180-4196.....	127	5182-0847.....	62	5183-4694.....	30	5188-5953.....	113, 115, 131
5181-1260.....	21, 80	5182-0848.....	62	5183-4695.....	30	5188-6471.....	31
5181-1261.....	21	5182-0849.....	62	5183-4696.....	30	5188-6493.....	74, 138
5181-1267.....	189	5182-0850.....	62	5183-4697.....	31	5188-6495.....	140
5181-1273.....	189	5182-0851.....	62	5183-4698.....	31	5188-6496.....	74, 138
5181-1291.....	48-49	5182-0852.....	62	5183-4699.....	31	5188-6497.....	74, 138
5181-1292.....	48-49	5182-0853.....	22	5183-4700.....	31	5188-6498.....	90-91, 137
5181-3308.....	37, 40, 172	5182-0875.....	190	5183-4701.....	28-31	5188-6537.....	59
5181-3315.....	30	5182-3442.....	22, 74, 76	5183-4702.....	28-31	5188-6538.....	59
5181-3316.....	30	5182-3444.....	22, 74, 76, 83, 85, 90-91	5183-4703.....	30	5188-8813.....	63
5181-3316i.....	30	5182-3445.....	22, 74, 76, 83, 85, 90-91	5183-4704.....	30	5188-8814.....	63
5181-3319.....	189	5182-3466.....	147	5183-4705.....	30	5188-8815.....	63
5181-3321.....	189	5182-3477.....	106, 136	5183-4706.....	30	5188-8816.....	63
5181-3322.....	37	5182-9622.....	191	5183-4707.....	30	5188-8817.....	63
5181-3323.....	37-39, 148, 172	5182-9626.....	191	5183-4708.....	30	5188-8818.....	63
5181-3351.....	440	5182-9633.....	80	5183-4709.....	30	5188-8819.....	63
5181-3352.....	440	5182-9645.....	80	5183-4710.....	30	5188-8820.....	63
5181-3354.....	189	5182-9651.....	87	5183-4711.....	29	5188-8821.....	63
5181-3356.....	189	5182-9652.....	67, 75, 139	5183-4712.....	29	5188-8822.....	63
5181-3358.....	189-190	5182-9673.....	106, 136	5183-4713.....	29	5190-0468.....	181
5181-3360.....	189	5182-9676.....	106, 136	5183-4728.....	190	5190-0471.....	181
5181-3361.....	189	5182-9677.....	106, 136	5183-4729.....	190	5190-0472.....	181
5181-3365.....	189	5182-9679.....	106, 136	5183-4730.....	190	5190-0473.....	181
5181-3382.....	90-91, 137	5182-9722.....	120, 133	5183-4732.....	40, 105	5190-0490.....	181
5181-3388.....	40	5182-9733.....	52, 54-55	5183-4741.....	64	5190-0491.....	181
5181-7459.....	47	5182-9734.....	191			5190-0492.....	181

PART NUMBER INDEX

5190-0493.....	181	5190-3193.....	190	8001-0159.....	227-229	8003-0106.....	221
5190-0494.....	181	5190-3976.....	64	8001-0160.....	227-229	8003-0107.....	221
5190-0513.....	181	5190-3978.....	64	8001-0162.....	228	8003-0108.....	222
5190-0531.....	177	5190-3983.....	28	8001-0163.....	228	8003-0109.....	222
5190-0585.....	172	5190-4006.....	28, 31	8001-0201.....	230	8003-0110.....	221
5190-1407.....	378	5190-4007.....	28	8001-0202.....	230	8003-0111.....	221
5190-1408.....	378	5190-4047.....	28, 31	8001-0203.....	230	8003-0151.....	221
5190-1409.....	378	5190-4048.....	28, 31	8001-0211.....	231	8003-0153.....	221
5190-1410.....	378	5190-4054.....	37, 40	8001-0212.....	231	8003-0154.....	222
5190-1426.....	32, 85, 87	5190-4061.....	60	8001-0213.....	231	8003-0155.....	222
5190-1437.....	53	5190-4062.....	60	8001-0214.....	231	8003-0157.....	221
5190-1438.....	53	5190-4063.....	60	8001-0221.....	231	8003-0158.....	221
5190-1441.....	181	5190-4064.....	60	8001-0222.....	231	8003-0159.....	221
5190-1445.....	63	5190-4065.....	60	8001-0223.....	231	8003-0160.....	221
5190-1446.....	63	5190-4066.....	60	8001-0224.....	231	8003-0162.....	221
5190-2209.....	16, 67, 75, 139	5190-4067.....	60	8001-0311.....	232	8003-0163.....	222
5190-2231.....	63	5190-5233.....	38, 40, 148, 172	8001-0312.....	232	8003-0165.....	221
5190-2232.....	63	5190-6144.....	67, 75	8002-0001.....	236	8003-0166.....	221
5190-2233.....	63	5190-6145.....	67, 75	8002-0002.....	236	8003-0202.....	222
5190-2234.....	63	5190-6149.....	16, 67	8002-0003.....	236	8003-0203.....	222
5190-2235.....	63	5190-6168.....	28, 53	8002-0004.....	236, 239	8003-0204.....	222
5190-2238.....	56	5190-6194.....	39-40, 91	8002-0005.....	236, 239	8003-0205.....	222
5190-2239.....	56	5190-6979.....	44	8002-0101.....	234	8003-0211.....	224
5190-2257.....	58	5190-6980.....	44	8002-0102.....	234	8003-0212.....	224
5190-2258.....	58	5190-6981.....	47	8002-0103.....	234	8003-0216.....	223
5190-2265.....	54	5190-6982.....	47	8002-0104.....	234	8003-0219.....	224
5190-2266.....	54	5958-9441.....	51	8002-0105.....	234	8003-0221.....	224
5190-2268.....	16, 33, 68, 138	5958-9442.....	51	8002-0106.....	234	8003-0222.....	224
5190-2269.....	16, 33	5958-9443.....	51	8002-0107.....	234	8003-0223.....	224
5190-2270.....	16, 30	5958-9444.....	51	8002-0151.....	234	8003-0311.....	224
5190-2271.....	16, 30	5958-9445.....	51	8002-0152.....	234	8004-0001.....	210
5190-2272.....	16, 30	5958-9450.....	51	8002-0153.....	234	8004-0002.....	210
5190-2275.....	16, 29	5982-0024.....	181	8002-0154.....	234	8004-0003.....	210
5190-2285.....	56	5982-0025.....	181	8002-0155.....	234	8004-0004.....	210
5190-2286.....	56	6040-0289.....	173-174	8002-0156.....	234	8004-0005.....	210
5190-2287.....	56	6040-0798.....	170	8002-0157.....	234	8004-0006.....	210
5190-2288.....	56	6040-0809.....	170	8002-0159.....	234	8004-0007.....	210
5190-2292.....	28	6040-0834.....	170	8002-0160.....	234	8004-0101.....	206
5190-2293.....	28	6040-1361.....	170	8002-0161.....	234	8004-0102.....	206
5190-2294.....	28	621-9723.....	292, 392	8002-0201.....	235	8004-0103.....	206
5190-2295.....	28	622-9732.....	292, 392	8002-0203.....	235	8004-0104.....	206
5190-2296.....	31	6410090050.....	54	8002-0204.....	235	8004-0105.....	206
5190-2297.....	28, 31	7200008400.....	213-217	8002-0211.....	235	8004-0106.....	206
5190-3151.....	62	8001-0004.....	232	8002-0212.....	235	8004-0107.....	206
5190-3152.....	62	8001-0005.....	232	8002-0213.....	235	8004-0108.....	207
5190-3153.....	62	8001-0006.....	232	8002-0214.....	235	8004-0109.....	207
5190-3154.....	62	8001-0007.....	232	8002-0215.....	235	8004-0110.....	207
5190-3157.....	16, 18	8001-0010.....	232	8002-0216.....	235	8004-0111.....	208
5190-3158.....	16, 19	8001-0011.....	232	8002-0217.....	235	8004-0112.....	207
5190-3162.....	28	8001-0012.....	232	8002-0220.....	235	8004-0113.....	206
5190-3163.....	28	8001-0013.....	232	8002-0221.....	235	8004-0114.....	206
5190-3164.....	28	8001-0014.....	232	8002-0222.....	235	8004-0116.....	206
5190-3165.....	28	8001-0101.....	227-229	8002-0311.....	236	8004-0118.....	207
5190-3166.....	28	8001-0102.....	227, 229	8002-0312.....	236	8004-0119.....	206
5190-3167.....	28	8001-0103.....	227-229	8003-0001.....	225	8004-0151.....	206
5190-3168.....	28	8001-0104.....	228	8003-0002.....	225	8004-0152.....	206
5190-3169.....	28	8001-0105.....	230	8003-0003.....	225	8004-0153.....	206
5190-3170.....	16, 28	8001-0106.....	227-228	8003-0004.....	225	8004-0154.....	206
5190-3171.....	16, 28	8001-0151.....	227-229	8003-0005.....	225	8004-0155.....	206
5190-3172.....	16, 28	8001-0152.....	227, 229	8003-0006.....	225	8004-0156.....	206
5190-3173.....	16, 28	8001-0153.....	227-229	8003-0007.....	225	8004-0157.....	206
5190-3188.....	60	8001-0154.....	228	8003-0008.....	225	8004-0158.....	206
5190-3189.....	60	8001-0155.....	230	8003-0101.....	221	8004-0159.....	207
5190-3190.....	60	8001-0156.....	227-228	8003-0103.....	221	8004-0160.....	207
5190-3191.....	60	8001-0157.....	227-229	8003-0104.....	222	8004-0161.....	207
5190-3192.....	60	8001-0158.....	227-229	8003-0105.....	221	8004-0162.....	207, 215

8004-0163.....	208	8010-0244.....	204, 238	8121-0723.....	127	CP1308.....	474
8004-0164.....	207, 214	8010-0245.....	204, 238	8121-0940.....	127	CP1309.....	476
8004-0165.....	206, 213	8010-0246.....	204, 238	8121-1222.....	127	CP1483.....	475
8004-0167.....	207, 215	8010-0249.....	205, 212, 238	8121-1301.....	127	CP17973.....	147, 164
8004-0168.....	208, 216	8010-0250.....	205, 212, 238	8121-1787.....	127	CP17977.....	147
8004-0170.....	206	8010-0251.....	205, 212, 238	8500-0656.....	146, 177, 182	CP17988.....	147, 164
8004-0171.....	207	8010-0252.....	205, 212, 238	8500-1233.....	146	CP2045.....	475
8004-0173.....	206	8010-0253.....	205, 212, 238	8500-5440.....	172, 182	CP2046.....	475
8004-0176.....	207	8010-0254.....	205, 212, 238	8500-5808.....	182	CP2050.....	479
8004-0178.....	206	8010-0255.....	205, 212, 226, 238	8500-5851.....	182	CP2055.....	474
8004-0201.....	208, 211, 213, 226, 233, 237	8010-0256.....	205, 212, 226, 238	8500-5995.....	182	CP2058.....	472
8004-0202.....	208, 213	8010-0257.....	205, 212, 238	8500-6812.....	482	CP2059.....	476
8004-0203.....	208	8010-0258.....	205, 212, 238	8500-6813.....	482	CP2062.....	473
8004-0204.....	208, 214	8010-0259.....	205, 238	8500-8510.....	163, 176	CP2065.....	475
8004-0211.....	209	8010-0260.....	205, 238	8500-8510.....	182	CP2068.....	473
8004-0212.....	209	8010-0261.....	205, 238	8650-0029.....	146, 174	CP2073.....	478
8004-0213.....	209	8010-0262.....	205, 238	8650-0030.....	146, 172, 174	CP4016.....	440
8004-0214.....	209	8010-0263.....	205, 233	8710-0510.....	40, 147, 172	CP4017.....	440
8004-0215.....	209	8010-0264.....	205, 233	8710-0899.....	147	CP4018.....	440
8004-0216.....	209	8010-0301.....	209, 223	8710-0900.....	147	CP4788.....	440
8004-0217.....	210, 215-217	8010-0302.....	209, 223	8710-1220.....	147	CP4789.....	440
8004-0218.....	209	8010-0303.....	209, 223	8710-1346.....	74, 138	CP4795.....	440
8004-0219.....	209	8010-0304.....	209, 223	8710-1561.....	120, 133	CP4796.....	440
8004-0311.....	210, 213-217	8010-0305.....	210, 215-217, 224	8710-1615.....	147	CP5881.....	297
8004-0312.....	210, 215	8010-0306.....	209, 223	8710-1622.....	147	CP5882.....	309
8010-0201.....	203, 211, 213, 237	8010-0307.....	209, 223	8829951700.....	170, 179	CP5883.....	309
8010-0202.....	203, 211, 237	8010-0308.....	209, 223	8829953800.....	179	CP5884.....	309
8010-0203.....	203, 211, 216-217, 237	8010-0309.....	209, 223	9300003590.....	179	CP5887.....	307
8010-0204.....	203, 211, 237	8010-0310.....	223, 231, 235	9301-0658.....	80	CP5889.....	307
8010-0205.....	203, 211, 237	8010-0311.....	223, 231, 235	9301-0713.....	23, 87, 190	CP5891.....	311
8010-0206.....	203, 211, 237	8010-0312.....	223, 231	9301-0714.....	190	CP5892.....	311
8010-0207.....	203, 211, 226, 237	8010-0313.....	223, 231	9301-0718.....	58	CP6172.....	406
8010-0208.....	203, 211, 226, 237	8010-0314.....	224, 231	9301-0719.....	58	CP6173.....	406
8010-0209.....	203, 211, 214-215, 237	8010-0315.....	224, 231	9301-0721.....	58	CP6174.....	406
8010-0210.....	203, 211, 237	8010-0351.....	239	9301-0723.....	50	CP6175.....	406
8010-0211.....	203, 237-238	8010-0352.....	239	9301-0725.....	190	CP6530.....	379
8010-0212.....	203, 237-238	8010-0353.....	236	9301-0891.....	190	CP6540.....	372, 417-418
8010-0213.....	203, 237	8010-0354.....	236	9301-0892.....	23, 87, 190	CP6550.....	372, 417-418
8010-0214.....	203, 237	8010-0355.....	236, 239	9301-0976.....	58	CP6571.....	469
8010-0215.....	203, 233	8010-0356.....	239	9301-1031.....	50	CP6572.....	469
8010-0216.....	203, 233	8010-0357.....	239	9310-4828.....	146, 174	CP6573.....	469
8010-0217.....	202, 211-213, 237	8010-0358.....	239	998-0000053-EHS.....	55	CP6574.....	469
8010-0218.....	202-205, 211, 226, 233, 237-238	8010-0359.....	239	C-102SSC.....	66	CP6575.....	469
8010-0219.....	202, 211, 216-217, 237	8010-0360.....	239	C-AT010C.....	66	CP6576.....	469
8010-0220.....	202, 211, 237	8010-0361.....	239	C-B010M.....	66	CP6577.....	469
8010-0221.....	202, 211, 237	8010-0362.....	239	C-BIO10.....	66	CP6578.....	469
8010-0222.....	202, 211, 237	8010-0363.....	239	C-BTX1UG.....	66	CP6579.....	469
8010-0223.....	202, 211, 226, 237	8010-0364.....	239	C-CF020.....	66	CP6580.....	469
8010-0224.....	202, 211, 226, 237	8010-0365.....	239	C-CPL0K.....	66	CP6581.....	469
8010-0225.....	202, 211, 214-215, 237	8010-0366.....	239	C-DF010.....	66	CP6666.....	326, 418
8010-0226.....	202, 211, 237	8010-0367.....	236, 239	C-G1CM10.....	66	CP6670.....	326, 418
8010-0227.....	202, 237	8010-0368.....	236	C-GAT010C.....	66	CP6671.....	326, 418
8010-0228.....	202, 237	8010-0371.....	239	C-GT010.....	66	CP6680.....	335, 419
8010-0229.....	202, 237	8010-0401.....	222, 235	C-HY010C.....	66	CP6681.....	419
8010-0230.....	202, 237	8120-5342.....	127	C-PL010.....	66	CP67569.....	418
8010-0231.....	202, 233	8120-6360.....	127	C-TBE10.....	66	CP68186.....	467
8010-0232.....	202, 233	8120-6894.....	127	C-TBP1C1C.....	66	CP68187.....	467
8010-0233.....	204, 212-213, 238	8120-6903.....	127	C-TBP1CXC.....	66	CP6918.....	419, 431
8010-0234.....	204, 212, 238	8120-8619.....	127	C-TBP1TC.....	66	CP6937.....	419, 439
8010-0239.....	204, 212, 226, 238	8120-8620.....	127	C-TNXTA.....	66	CP6938.....	419, 439
8010-0240.....	204, 212, 226, 238	8120-8621.....	127	C-UN010C.....	66	CP6938I5.....	419, 439
8010-0241.....	204, 212, 214-215, 238	8120-8622.....	127	CP0031.....	483	CP6953.....	419, 424
8010-0242.....	204, 212, 238	8121-0070.....	127	CP108194.....	467	CP6954.....	419, 424
8010-0243.....	204, 238	8121-0161.....	127	CP1305.....	474	CP6968.....	419, 431
		8121-0675.....	127	CP1306.....	475	CP7128.....	357
		8121-0710.....	127	CP1307.....	473	CP7130.....	418

PART NUMBER INDEX

CP7135.....	326, 418	CP7440.....	391, 419	CP7529.....	283, 375	CP7584PT.....	421, 427
CP7140.....	326, 418	CP7441.....	324	CP7530.....	368	CP7586.....	431
CP7141.....	419	CP7442.....	325	CP7531.....	282, 368	CP7587.....	431
CP7148.....	357, 419	CP7443.....	324	CP7531I5.....	368	CP7588.....	270-271, 276, 396
CP7150.....	326, 418	CP7443I5.....	324	CP7532.....	279, 372, 417-418	CP7591.....	383
CP7160.....	326, 418	CP7444.....	325	CP7532I5.....	372, 417	CP7592.....	280, 372, 417-418
CP7168.....	357, 419	CP7446.....	381	CP7533.....	278, 439	CP7593.....	372, 417-418
CP7170.....	326, 418	CP7447.....	381	CP7534.....	439	CP7594.....	383
CP7178.....	357, 419	CP7447I5.....	381	CP7534I5.....	439	CP7595.....	382
CP7179.....	357, 419	CP7448.....	381	CP7534PT.....	421, 439	CP7595I5.....	382
CP7196.....	335, 419	CP7448I5.....	381	CP7535.....	439	CP7596.....	269, 382
CP7310.....	324	CP7451.....	334	CP7535I5.....	439	CP7596I5.....	382
CP7311.....	324	CP7452.....	335	CP7536.....	439	CP7597.....	273, 382
CP7334.....	356	CP7453.....	334	CP7536PT.....	421, 439	CP7598.....	274, 382
CP7335.....	356	CP7453I5.....	334	CP7537.....	278, 439	CP7598I5.....	382
CP7340.....	350	CP7454.....	335	CP7538.....	439	CP7599.....	382
CP7347.....	422	CP7461.....	390	CP7538I5.....	439	CP7608.....	326
CP7348.....	422	CP7462.....	390	CP7538PT.....	421, 439	CP7614.....	282, 360
CP7348PT.....	421-422	CP7463.....	409, 419	CP7539.....	439	CP7615.....	285, 409
CP7350.....	280, 422	CP7475.....	367	CP7539PT.....	421, 439	CP7617.....	409
CP7350I5.....	422	CP7476.....	395	CP7540.....	439	CP7619.....	336
CP7351.....	422	CP7477.....	395	CP7540I5.....	439	CP7620.....	326
CP7351I5.....	422	CP7481.....	387	CP7541.....	371	CP7621.....	335
CP7351PT.....	281-282, 421-422	CP7482.....	281, 393	CP7542.....	372, 417-418	CP7622.....	324
CP7352.....	422	CP7483.....	409	CP7543.....	439	CP7624.....	360
CP7352I5.....	422	CP7485.....	410	CP7544.....	439	CP7625.....	326
CP7352PT.....	421-422	CP7485I5.....	410	CP7548.....	424	CP7628.....	357
CP7353.....	422	CP7486.....	279, 410	CP7549.....	424	CP7629.....	336
CP7353I5.....	422	CP7487.....	407	CP7550.....	424	CP7630.....	326
CP7353PT.....	280, 421-422	CP7488.....	278-279, 407	CP7550I5.....	424	CP7631.....	335
CP7354.....	272, 422	CP7489.....	407	CP7550PT.....	421, 424	CP7635.....	326
CP7354I5.....	422	CP7491.....	404	CP7551.....	424	CP7636.....	335
CP7354PT.....	281-282, 421-422	CP749103.....	404	CP7551I5.....	424	CP7637.....	350
CP7355.....	422	CP749106.....	404	CP7551PT.....	421, 424	CP7638.....	357
CP7371.....	436	CP7491I5.....	404	CP7552.....	424	CP7640.....	326
CP7372.....	436	CP7494.....	412	CP7553.....	424	CP7641.....	335
CP7374.....	278, 436	CP7495.....	412	CP7554.....	279, 424	CP7645.....	326
CP7375.....	436	CP7495I5.....	412	CP7554I5.....	424	CP7646.....	335
CP7381.....	423	CP7498.....	396	CP7554PT.....	421, 424	CP7647.....	350
CP7406.....	387	CP7500.....	413	CP7555.....	424	CP7648.....	357
CP7407.....	387	CP7500I5.....	413	CP7556.....	424	CP7649.....	336
CP7409.....	387	CP7501.....	413	CP7557.....	424	CP7654.....	272, 360
CP7411.....	399	CP7502.....	412	CP7557PT.....	421, 424	CP7656.....	335
CP7412.....	399	CP7502I5.....	412	CP7558.....	278, 424	CP7657.....	350
CP7413.....	399	CP7503.....	412	CP7559.....	424	CP7658.....	357
CP7414.....	279, 399	CP7504.....	387	CP7562.....	372, 417-418	CP7658I5.....	357
CP7415.....	399	CP7511.....	430	CP7565.....	431	CP7659.....	336
CP7416.....	399	CP7512.....	372, 417-418	CP7565I5.....	431	CP7664.....	360
CP7416I5.....	399	CP7512I5.....	372	CP7565PT.....	421, 431	CP7665.....	360
CP7417.....	279, 399	CP7513.....	436	CP7567.....	281, 431	CP7666.....	335
CP7418.....	399	CP7514.....	436	CP7568.....	431	CP7667.....	350
CP7419.....	407	CP7515.....	430	CP7568PT.....	421, 431	CP7668.....	357
CP7419I5.....	407	CP7515I5.....	430	CP7569.....	372, 417-418	CP7670.....	324
CP7420.....	407	CP7515PT.....	279, 421, 430	CP7569I5.....	372, 417	CP7670I5.....	324
CP7421.....	407	CP7517.....	430	CP7570.....	279, 372, 417-418	CP7671.....	334
CP7422.....	382	CP7517PT.....	421, 430	CP7571.....	372, 417-418	CP7672.....	350
CP7424.....	382	CP7518.....	430	CP7572.....	372, 417-418	CP7673.....	356
CP7428.....	281	CP7518PT.....	278-279, 421, 430	CP7574.....	427	CP7673I5.....	356
CP7429.....	376	CP7519.....	430	CP7576.....	430	CP7675.....	326
CP7430.....	376	CP7521.....	371	CP7577.....	430	CP7676.....	335
CP7431.....	376	CP7522.....	371	CP7579.....	427	CP7678.....	357
CP7432.....	376	CP7522I5.....	371	CP7580.....	427	CP7680.....	325
CP7433.....	376	CP7525.....	280-282, 284, 373	CP7581.....	427	CP7680I5.....	325
CP7434.....	380	CP7525I5.....	373	CP7582.....	279, 372, 417-418	CP7681.....	335
CP7435.....	380	CP7527.....	408	CP7583.....	427	CP7682.....	324
CP7437.....	380	CP7528.....	403	CP7584.....	427	CP7684.....	324

CP7685.....	285, 326	CP7770.....	325	CP81073.....	474	CP8718.....	357
CP7685I5.....	326	CP7770I5.....	325	CP8120.....	408	CP8722.....	350
CP7686.....	410	CP7771.....	335	CP8121.....	389, 467	CP8723.....	356
CP7688.....	326	CP7772.....	350	CP8128.....	466	CP8730.....	326
CP7690.....	325	CP7773.....	357	CP8129.....	466	CP8735.....	278-281, 326
CP7690I5.....	325	CP7773I5.....	357	CP8131.....	389	CP8735I5.....	326
CP7691.....	335	CP7775.....	356	CP8132.....	389	CP8736.....	272, 335
CP7691I5.....	335	CP7777.....	360	CP8133.....	389	CP8736I5.....	335
CP7692.....	324	CP7778.....	360	CP8134.....	389, 466	CP8737.....	350
CP7693.....	324	CP7779.....	335	CP8135.....	466	CP8738.....	279-280, 285, 357
CP7694.....	324	CP7785.....	356	CP8186.....	467	CP8738I5.....	357
CP7695.....	326	CP7787.....	360	CP8187.....	467	CP8741.....	324
CP7696.....	335	CP7789.....	335	CP8430.....	280, 325	CP8741I5.....	324
CP7697.....	350	CP7791.....	356	CP8431.....	335	CP8742.....	325
CP7698.....	357	CP7792.....	356	CP8510.....	324	CP8743.....	324
CP7698I5.....	357	CP7797.....	360	CP8511.....	334	CP8744.....	325
CP7700.....	324	CP7809.....	350	CP8513.....	356	CP8745.....	356
CP7702.....	350	CP7813.....	336	CP8521.....	334	CP8746.....	356
CP7703.....	356	CP7820.....	341	CP8529.....	325	CP8748.....	356
CP7709.....	324	CP7821.....	341	CP8530.....	325	CP8750.....	334
CP7710.....	324	CP7822.....	341	CP8531.....	335	CP8751.....	276, 334
CP7711.....	334	CP7822I5.....	341	CP8533.....	357	CP8752.....	335
CP7712.....	350	CP7824.....	341	CP8540.....	325	CP8753.....	334
CP7713.....	356	CP7829.....	350	CP8541.....	335	CP8754.....	335
CP7713I5.....	356	CP7830.....	341	CP8542.....	350	CP8756.....	335
CP7714.....	373	CP7831.....	341	CP8543.....	357	CP8760.....	279, 281, 325
CP7715.....	346	CP7832.....	341	CP8550.....	325	CP8761.....	335
CP7717.....	360	CP7834.....	341	CP8553.....	357	CP8762.....	350
CP7717I5.....	360	CP7839.....	350	CP8560.....	325	CP8763.....	279-280, 285, 357
CP7719.....	324	CP7871.....	341	CP8562.....	350	CP8770.....	324
CP7720.....	324	CP7871I5.....	341	CP8564.....	435	CP8771.....	334
CP7721.....	334	CP7879.....	357	CP8566.....	435	CP8772.....	350
CP7722.....	350	CP7884.....	334	CP8567.....	435	CP8773.....	282, 357
CP7723.....	356	CP7889.....	357	CP8568.....	435	CP8774.....	326
CP7723I5.....	356	CP7906.....	336	CP8570.....	435	CP8775.....	284, 326
CP7725.....	346	CP7907.....	336	CP8570I5.....	435	CP8780.....	324
CP7727.....	360	CP7917.....	336	CP8571.....	435	CP8781.....	334
CP7730.....	325	CP7926.....	336	CP8575.....	375	CP8791.....	335
CP7731.....	335	CP7926I5.....	336	CP8587.....	366	CP8796.....	335
CP7732.....	350	CP7936.....	336	CP8587I5.....	366	CP8798.....	280, 357
CP7739.....	325	CP7937.....	336	CP8605.....	348	CP8799.....	326
CP7740.....	325	CP7945.....	368	CP8607.....	348	CP8801.....	305
CP7741.....	335	CP7946.....	336	CP8610.....	348	CP8805.....	305
CP7741I5.....	335	CP7947.....	336	CP8613.....	348	CP8806.....	305
CP7742.....	350	CP7977.....	336	CP8662.....	350	CP8807.....	305
CP7743.....	357	CP8003.....	465	CP8667.....	474	CP8809.....	305
CP7745.....	346	CP8007.....	466	CP86678.....	474	CP8813.....	305
CP7747.....	360	CP8008.....	466	CP8673.....	326	CP8816.....	305
CP7749.....	325	CP8009.....	466	CP8674.....	326	CP8821.....	310
CP7749I5.....	325	CP8014.....	335	CP8675.....	326	CP8822.....	310
CP7750.....	325	CP8015.....	404, 466	CP8676.....	326	CP8822I5.....	310
CP7750I5.....	325	CP8016.....	466	CP8677.....	280-281, 326	CP8824.....	310
CP7751.....	335	CP8017.....	466	CP8678.....	335	CP8824I5.....	310
CP7751I5.....	335	CP8018.....	466	CP8685.....	326	CP8827.....	310
CP7752.....	350	CP8019.....	466	CP8687.....	280-281, 325	CP8828.....	310
CP7753.....	357	CP8070.....	467	CP8687I5.....	325	CP8829.....	310
CP7754.....	373	CP8073.....	357	CP8688.....	325	CP8831.....	310
CP7757.....	360	CP8073I5.....	357	CP8688I5.....	325	CP8842.....	350
CP7758.....	325	CP8080.....	467	CP8689.....	325	CP8843.....	357
CP7760.....	325	CP8087.....	466	CP8690.....	325	CP8853.....	357
CP7761.....	335	CP8088.....	466	CP8690I5.....	325	CP8855.....	311
CP7762.....	350	CP8089.....	466	CP8710.....	324	CP8857.....	311
CP7763.....	357	CP8090.....	467	CP8712.....	350	CP8858.....	271, 311
CP7764.....	373	CP81025.....	475	CP8713.....	356	CP8859.....	311
CP7767.....	360	CP81069.....	474	CP8713I5.....	356	CP8859I5.....	311
CP7769.....	334	CP81071.....	476	CP8716.....	335	CP8860.....	270, 311

PART NUMBER INDEX

CP8861.....	311	CP8951.....	303	CP9053.....	316	CP9205.....	313
CP8864.....	311	CP8953.....	303	CP9054.....	267, 316	CP9205I5.....	313
CP8865.....	311	CP8955.....	265, 276, 303	CP9055.....	316	CP9206.....	269, 313
CP8867.....	311	CP8956.....	303	CP9056.....	316	CP9207.....	313
CP8868.....	274, 311	CP8957.....	267, 275, 303	CP9057.....	316	CP9210.....	314
CP8870.....	325	CP8958.....	303	CP9058.....	316	CP9211.....	314
CP8871.....	335	CP8959.....	303	CP9062.....	316	CP9212.....	314
CP8872.....	307	CP8960.....	271, 303	CP9063.....	316	CP9213.....	314
CP8873.....	307	CP8961.....	303	CP9064.....	316	CP9214.....	314
CP8874.....	307	CP8962.....	303	CP9066.....	316	CP9215.....	272, 274, 276, 314
CP8875.....	307	CP8963.....	303	CP9070.....	265-266, 270-271, 275, 386	CP9216.....	314
CP8877.....	270, 275, 307	CP8964.....	415	CP9071.....	386	CP9217.....	314
CP8878.....	307	CP8965.....	298	CP9072.....	386	CP9218.....	313
CP8879.....	273, 307	CP8967.....	273, 298	CP9073.....	385	CP9219.....	313
CP8880.....	307	CP8968.....	298	CP9074.....	265-267, 270-271, 275, 385	CP9220.....	313
CP8882.....	307	CP8969.....	274, 298	CP9075.....	385	CP9221.....	313
CP8883.....	307	CP8970.....	298	CP9076.....	379, 419	CP9222.....	313
CP8884.....	307	CP8971.....	303	CP9077.....	379, 419	CP9223.....	313
CP8888.....	274, 307	CP8974.....	303	CP9078.....	285, 379, 419	CP9225.....	314
CP8900.....	297	CP8975.....	274, 303	CP9079.....	379, 419	CP9226.....	272, 314
CP8901.....	297	CP8976.....	270, 272-274, 303	CP9080.....	379	CP9228.....	314
CP8902.....	297	CP8977.....	309	CP9083.....	379	CP9229.....	313
CP8903.....	297	CP8979.....	309	CP9090.....	365, 417	CP97638.....	408
CP8904.....	297	CP8980.....	309	CP9091.....	365, 418	CP97658.....	408
CP8905.....	297	CP8981.....	309	CP9092.....	365, 417	CP97713.....	408
CP8906.....	280-281, 284, 298	CP8982.....	270, 309	CP9093.....	365, 418	CP97721.....	408
CP8907.....	298	CP8982I5.....	309	CP9094.....	365, 417	CP97723.....	408
CP8908.....	298	CP8983.....	271, 309	CP9094I5.....	365, 417	CP97723I5.....	408
CP8908I5.....	298	CP8984.....	309	CP9095.....	365, 418	CP97743.....	408
CP8909.....	298	CP8986.....	309	CP9096.....	365, 417	CP97753.....	408
CP8910.....	298	CP8990.....	309	CP9097.....	365, 418	CP97753I5.....	408
CP8911.....	298	CP8991.....	265, 309	CP9100.....	315	CP97763.....	408
CP8912.....	298	CP8996.....	309	CP9101.....	315	CP97773.....	408
CP8912I5.....	298	CP8998.....	309	CP9102.....	265, 267, 269, 315	CR213105.....	209
CP8913.....	267, 298	CP9001.....	272, 309	CP9102I5.....	315	DY50295500.....	64
CP8914.....	298	CP9002.....	273, 309	CP9103.....	265, 267, 272, 275, 315	DY50296800.....	64
CP8915.....	298	CP9010.....	446	CP9103I5.....	315	DY50390600.....	64
CP8916.....	298	CP9011.....	446	CP9104.....	315	DY50402400.....	64
CP8917.....	298	CP9012.....	446	CP9104I5.....	315	DY50504700.....	64
CP8919.....	298	CP9013.....	446	CP9105.....	269, 271, 275, 315	DY50546100.....	64
CP8921.....	298	CP9014.....	446	CP9106.....	315	DY50546390.....	64
CP8922.....	298	CP9015.....	446	CP9106I5.....	315	DY50548400.....	64
CP8923.....	298	CP9016.....	446	CP9107.....	272, 315	DY50549100.....	64
CP8924.....	265, 272, 298	CP9019.....	446	CP9108.....	269, 271, 315	DY50549290.....	64
CP8925.....	298	CP9021.....	446	CP9109.....	315	DY50549500.....	64
CP8926.....	265, 298	CP9022.....	446	CP9110.....	315	DY50551400.....	64
CP8928.....	298	CP9023.....	446	CP9145.....	317	DY50559800.....	64
CP8929.....	298	CP9024.....	317	CP9150.....	317	DY50559900.....	64
CP8930.....	298	CP9025.....	446	CP9151.....	266, 268, 317	DY50572600.....	64
CP8934.....	302	CP9026.....	446	CP9151I5.....	317	DY50573990.....	64
CP8935.....	302	CP9027.....	446	CP9152.....	317	DY50574190.....	64
CP8936.....	302	CP9030.....	297	CP9152I5.....	317	DY50574390.....	64
CP8937.....	302	CP9031.....	297	CP9154.....	317	DY50574500.....	64
CP8938.....	303	CP9032.....	297	CP9156.....	317	DY70001990.....	64
CP8939.....	267, 303	CP9034.....	302	CP9162.....	317	DY70007691.....	64
CP8940.....	303	CP9035.....	302	CP9163.....	265, 317	DY70007701.....	64
CP8941.....	303	CP9036.....	302	CP9165.....	317	DY70007791.....	64
CP8943.....	270, 303	CP9037.....	302	CP9166.....	317	DY70008101.....	64
CP8944.....	266-271, 273-274, 276, 303	CP9038.....	302	CP9170.....	317	DY70008590.....	64
CP8944I5.....	303	CP9039.....	302	CP9171.....	272-274, 317	G1072-20008.....	154-155
CP8945.....	270, 276, 303	CP9041.....	305	CP9176.....	446	G1099-20030.....	38, 126, 147-148, 172
CP8946.....	266, 273, 275-276, 303	CP9044.....	365	CP9177.....	446	G1099-20130.....	154
CP8947.....	303	CP9045.....	365	CP9201.....	313	G1099-20132.....	154
CP8948.....	270, 303	CP9046.....	365	CP9202.....	313	G1099-20133.....	154-155, 157
CP8949.....	303	CP9047.....	365	CP9203.....	313	G1099-20136.....	154-155
CP8950.....	303	CP9048.....	365	CP9204.....	313	G1099-60566.....	146

G1099-80039.....	170	G2589-20044.....	155, 157	G3440-20003.....	48	G3591-70013.....	475
G1530-60560.....	127	G2589-20045.....	151, 155	G3440-20004.....	48	G3591-70014.....	478
G1530-60570.....	127	G2589-20100.....	151, 155	G3440-20005.....	48	G3591-70015.....	476
G1530-60930.....	127	G2591B.....	153	G3440-20007.....	48	G3591-70016.....	471
G1530-61230.....	126	G2591C.....	153	G3440-20008.....	48	G3591-70017.....	476
G1530-61580.....	126	G2591D.....	153	G3440-20033.....	49	G3591-70018.....	478
G1530-61610.....	126	G2617-20510.....	84	G3440-20035.....	49	G3591-74001.....	474
G1530-61640.....	126	G2617-60506.....	83	G3440-20037.....	49	G3591-80000.....	480
G1530-80650.....	123, 127	G2617-60507.....	83, 85, 87	G3440-60033.....	49	G3591-80001.....	472
G1531-20690.....	98	G2617-60508.....	83	G3440-60035.....	49	G3591-80002.....	474
G1531-20700.....	98	G2617-60510.....	84	G3440-60136.....	49	G3591-80003.....	476
G1531-20740.....	98	G2617-80540.....	84	G3440-60233.....	49	G3591-80004.....	474
G1531-21090.....	98	G2617-80550.....	85, 87	G3440-60234.....	49	G3591-80006.....	480
G1531-60690.....	98	G2618-80500.....	83, 85	G3440-60236.....	49	G3591-80007.....	479
G1531-80560.....	96, 119	G2619-60501.....	83	G3440-60300.....	49	G3591-80008.....	479
G1531-80620.....	96, 119	G2630-60710.....	123	G3440-60310.....	49	G3591-80011.....	473
G1532-20710.....	105, 113, 121, 131, 134, 136	G2630-61230.....	107	G3440-60333.....	49	G3591-80012.....	476
G1532-60675.....	107	G2646-60500.....	123	G3440-60334.....	49	G3591-80013.....	478
G1532-60685.....	107	G2647-60501.....	113, 131	G3440-60336.....	49	G3591-80014.....	479
G1532-60690.....	107	G2648-60501.....	113, 131	G3440-60600.....	49	G3591-80015.....	479
G1532-60695.....	107	G2855-20530.....	43, 49	G3440-60610.....	49	G3591-80016.....	478
G1532-80540.....	105, 136	G2855-20532.....	49	G3440-60620.....	49	G3591-80017.....	475
G1534-20530.....	120, 133	G2855-20555.....	38, 40, 43, 49, 148, 172	G3440-80217.....	39-40, 74, 76	G3591-80019.....	474
G1534-20590.....	120, 133	G2855-20590.....	43, 53	G3440-80218.....	39-40, 74, 76	G3591-80020.....	474
G1534-60570.....	116, 120, 133	G2855-40001.....	53	G3440-81664.....	49	G3591-80021.....	471
G1534-60640.....	120, 133	G2855-60200.....	38, 43, 53, 148, 172	G3440-81665.....	49	G3591-80022.....	475
G1534-80510.....	120, 133	G2855-60570.....	53	G3440-85007.....	181	G3591-80023.....	474
G1534-80580.....	119	G2933-85001.....	123	G3440-85009.....	181	G3591-80025.....	476
G1534-80590.....	119, 135	G2933-85003.....	123	G3440-85012.....	181	G3591-80026.....	479
G1535-00010.....	114	G3163-20530.....	154	G3440-85013.....	181	G3591-80027.....	472
G1535-00030.....	114	G3170-20126.....	154-155, 157, 159	G3440-85017.....	181	G3591-80028.....	476
G1535-60600.....	113, 115, 131	G3170-20530.....	154-155	G3440-85018.....	181	G3591-80029.....	479
G1535-60610.....	114	G3170-20540.....	159	G3440-85026.....	181	G3591-80030.....	472
G1535-80520.....	114	G3170-60053.....	157	G3440-85027.....	181	G3591-80031.....	476
G1540-30025.....	48	G3170-60204.....	163	G3440-85028.....	181	G3591-80032.....	472
G1540-30026.....	48	G3170-60416.....	159	G3440-85029.....	181	G3591-80033.....	474
G1540-30027.....	48	G3170-80001.....	162	G3440-85035.....	181	G3591-80034.....	480
G1540-30028.....	48	G3170-80002.....	163, 176-177	G3440-85036.....	181	G3591-80035.....	476
G1540-80013.....	90	G3170-80002.....	146, 174	G3451-80501.....	90, 137	G3591-80036.....	477
G1544-20590.....	75, 139	G3170-80008.....	165	G3452-20512.....	76	G3591-80037.....	473
G1544-60585.....	74	G3170-80100.....	165	G3452-60570.....	74	G3591-80039.....	472
G1544-80700.....	32	G3170-80103.....	165, 173-174	G3452-60586.....	74	G3591-80040.....	480
G1544-80730.....	32	G3182-61580.....	42	G3452-60730.....	74	G3591-80043.....	476
G1544-80731.....	32	G3182-61581.....	42	G3452-60835.....	53	G3591-80044.....	477
G1545-80520.....	80	G3184-60065.....	43	G3452-80570.....	74	G3591-80045.....	473
G1888-60701.....	55	G3188-20509.....	49	G3480-20002.....	48	G3591-80046.....	475
G1888-60702.....	54	G3188-27501.....	37-39, 41, 43, 49, 148, 172	G3480-60663.....	48	G3591-80047.....	474
G1888-60703.....	54	G3188-27502.....	37-39, 41, 43, 49, 148, 172	G3480-67585.....	74	G3591-80048.....	476
G1888-60704.....	54	G3188-27503.....	37, 41, 43, 49, 172	G3500-80000.....	87	G3591-80049.....	480
G1890-60000.....	54	G3188-27504.....	49, 172	G3504-20504.....	40	G3591-80051.....	478
G1960-80303.....	174	G3188-27505.....	49	G3504-60620.....	53	G3591-80052.....	478
G1999-20021.....	154-155, 159	G3188-27506.....	49	G3507-60660.....	48	G3591-80053.....	478
G1999-20022.....	154-155, 159	G3397A.....	162	G3510-20018.....	76	G3591-80054.....	476
G1999-20430.....	159	G3397B.....	162	G3520-20210.....	53	G3591-80055.....	471
G1999-20432.....	159	G3430-60011.....	74	G3591-70001.....	478	G3591-80056.....	475
G1999-20433.....	159	G3431-60680.....	98	G3591-70002.....	475	G3591-80058.....	476
G1999-20443.....	159	G3432-60220.....	107	G3591-70003.....	475	G3591-80059.....	477
G1999-20444.....	159	G3432-60221.....	107	G3591-70004.....	475	G3591-80060.....	473
G1999-20446.....	159	G3433-63000.....	102, 132	G3591-70005.....	474	G3591-80061.....	475
G1999-60412.....	159	G3434-60806.....	116, 120	G3591-70006.....	474	G3591-80062.....	475
G1999-60452.....	163, 176	G3435-60350.....	115	G3591-70007.....	474	G3591-80063.....	479
G1999-80410.....	147, 164	G3435-81330.....	115	G3591-70008.....	475	G3591-80064.....	475
G2397-20540.....	102, 132	G3435-81360.....	115	G3591-70009.....	478	G3591-80065.....	475
G2397-80520.....	102	G3440-05012.....	181	G3591-70010.....	478	G3591-80066.....	478
G2589-20043.....	155			G3591-70011.....	474	G3591-80067.....	472
				G3591-70012.....	471	G3591-80068.....	479

PART NUMBER INDEX

G3591-80071.....	479	G3591-81016.....	478	G3591-81098.....	476	G3591-82048.....	476
G3591-80072.....	477	G3591-81017.....	475	G3591-81099.....	471	G3591-82049.....	480
G3591-80074.....	475	G3591-81019.....	474	G3591-81100.....	474	G3591-82050.....	471
G3591-80075.....	475	G3591-81020.....	474	G3591-81101.....	476	G3591-82051.....	478
G3591-80077.....	475	G3591-81021.....	471	G3591-81102.....	474	G3591-82053.....	478
G3591-80082.....	472	G3591-81022.....	475	G3591-81103.....	475	G3591-82063.....	479
G3591-80083.....	478	G3591-81023.....	474	G3591-81104.....	475	G3591-82064.....	475
G3591-80084.....	470	G3591-81025.....	476	G3591-81105.....	470	G3591-82067.....	472
G3591-80085.....	476	G3591-81026.....	479	G3591-81106.....	478	G3591-82072.....	477
G3591-80086.....	477	G3591-81027.....	472	G3591-81119.....	477	G3591-82082.....	472
G3591-80087.....	477	G3591-81028.....	476	G3591-81120.....	478	G3591-82084.....	470
G3591-80088.....	473	G3591-81029.....	479	G3591-81121.....	474	G3591-82087.....	477
G3591-80090.....	475	G3591-81030.....	472	G3591-81122.....	477	G3591-82088.....	473
G3591-80091.....	473	G3591-81031.....	476	G3591-81123.....	477	G3591-82093.....	476
G3591-80093.....	476	G3591-81032.....	472	G3591-81135.....	478	G3591-82095.....	470
G3591-80094.....	480	G3591-81033.....	474	G3591-81136.....	478	G3591-82099.....	471
G3591-80095.....	470	G3591-81034.....	480	G3591-81137.....	478	G3591-82101.....	476
G3591-80096.....	478	G3591-81035.....	476	G3591-81138.....	478	G3591-82102.....	474
G3591-80097.....	476	G3591-81036.....	477	G3591-81139.....	471	G3591-82104.....	475
G3591-80098.....	476	G3591-81037.....	473	G3591-81140.....	472	G3591-82106.....	478
G3591-80099.....	471	G3591-81039.....	472	G3591-81141.....	479	G3591-82119.....	477
G3591-80100.....	474	G3591-81040.....	480	G3591-81142.....	479	G3591-82121.....	474
G3591-80101.....	476	G3591-81043.....	476	G3591-81143.....	474	G3591-82122.....	477
G3591-80103.....	475	G3591-81044.....	477	G3591-81146.....	474	G3591-82123.....	477
G3591-80104.....	475	G3591-81045.....	473	G3591-81147.....	475	G3591-82135.....	478
G3591-80105.....	470	G3591-81046.....	475	G3591-81149.....	475	G3591-82136.....	478
G3591-80106.....	478	G3591-81047.....	474	G3591-81156.....	473	G3591-82137.....	478
G3591-80107.....	475	G3591-81048.....	476	G3591-81157.....	478	G3591-82139.....	471
G3591-80108.....	479	G3591-81049.....	480	G3591-81160.....	472	G3591-82140.....	472
G3591-80110.....	477	G3591-81050.....	471	G3591-81161.....	472	G3591-82159.....	474
G3591-80119.....	477	G3591-81051.....	478	G3591-81162.....	476	G3591-82176.....	479
G3591-80120.....	478	G3591-81052.....	478	G3591-81176.....	479	G3850-60014.....	147
G3591-80121.....	474	G3591-81053.....	478	G3591-81209.....	475	G3870-20021.....	157
G3591-80122.....	477	G3591-81054.....	476	G3591-81210.....	473	G3870-20135.....	154, 157
G3591-80123.....	477	G3591-81055.....	471	G3591-81211.....	473	G3870-20440.....	157
G3591-80124.....	474	G3591-81056.....	475	G3591-81212.....	473	G3870-20444.....	151, 157
G3591-80135.....	478	G3591-81058.....	476	G3591-81213.....	473	G3870-20445.....	157
G3591-80136.....	478	G3591-81059.....	477	G3591-81214.....	476	G3870-20446.....	157
G3591-80137.....	478	G3591-81060.....	473	G3591-81215.....	480	G3870-20448.....	151
G3591-80138.....	478	G3591-81061.....	475	G3591-82000.....	480	G3870-20449.....	151
G3591-80139.....	471	G3591-81062.....	475	G3591-82001.....	472	G3870-20530.....	157
G3591-80140.....	472	G3591-81063.....	479	G3591-82002.....	474	G3870-60171.....	157
G3591-80141.....	479	G3591-81064.....	475	G3591-82003.....	476	G3870-60172.....	154
G3591-80142.....	479	G3591-81065.....	475	G3591-82004.....	474	G3870-60179.....	155
G3591-80156.....	473	G3591-81066.....	478	G3591-82011.....	473	G3880-20030.....	38, 147-148
G3591-80157.....	478	G3591-81067.....	472	G3591-82012.....	476	G3880-80010.....	162
G3591-80158.....	473	G3591-81068.....	479	G3591-82013.....	478	G3880-80011.....	162
G3591-80160.....	472	G3591-81071.....	479	G3591-82014.....	479	G3900-63001.....	462
G3591-80161.....	472	G3591-81072.....	477	G3591-82015.....	479	G3900-63002.....	462
G3591-80162.....	476	G3591-81074.....	475	G3591-82016.....	478	G3900-63003.....	463
G3591-80170.....	479	G3591-81075.....	475	G3591-82017.....	475	G3900-63004.....	462
G3591-80171.....	479	G3591-81077.....	475	G3591-82020.....	474	G3900-63005.....	462
G3591-80176.....	479	G3591-81082.....	472	G3591-82022.....	475	G3900-63006.....	463
G3591-81000.....	480	G3591-81083.....	478	G3591-82023.....	474	G3900-63007.....	463
G3591-81001.....	472	G3591-81084.....	470	G3591-82025.....	476	G3900-63008.....	463
G3591-81002.....	474	G3591-81085.....	476	G3591-82026.....	479	G3900-63009.....	462
G3591-81003.....	476	G3591-81086.....	477	G3591-82029.....	479	G3900-63010.....	463
G3591-81004.....	474	G3591-81087.....	472	G3591-82032.....	472	G3900-63011.....	463
G3591-81006.....	480	G3591-81088.....	473	G3591-82033.....	474	G3900-63012.....	463
G3591-81007.....	479	G3591-81090.....	475	G3591-82035.....	476	G3900-63013.....	462
G3591-81008.....	479	G3591-81091.....	473	G3591-82036.....	477	G3900-63014.....	462
G3591-81011.....	473	G3591-81093.....	476	G3591-82037.....	473	G3900-63015.....	463
G3591-81012.....	476	G3591-81094.....	480	G3591-82039.....	472	G3900-63016.....	462
G3591-81013.....	478	G3591-81095.....	470	G3591-82040.....	480	G3900-63017.....	462
G3591-81014.....	479	G3591-81096.....	478	G3591-82043.....	476	G3900-63018.....	462
G3591-81015.....	479	G3591-81097.....	476	G3591-82044.....	477	G3900-63019.....	462

G3900-63020.....	463	G4513-80233.....	188	MKI-UTD-5064.....	66
G3900-63021.....	463	G4513-80234.....	187	MKI-Z-0285.....	65
G3900-63022.....	463	G4513-80235.....	187	MKI-Z-0351.....	65
G3900-63023.....	463	G4513-80236.....	187-188	RDT-1020.....	52, 140
G3900-63024.....	463	G4513-80239.....	187	RDT-1023.....	140
G3900-63025.....	463	G4513-80240.....	187	RMSN-2.....	164
G3900-63026.....	463	G4513-80241.....	188		
G3900-63027.....	463	G4513-80242.....	187		
G3900-63028.....	463	G4513-80243.....	187		
G3900-63029.....	463	G4513-80244.....	187		
G3900-63030.....	462	G4514-60610.....	127		
G3900-63031.....	462	G4514-60710.....	50		
G3900-63032.....	462	G4514-67505.....	50		
G3900-63033.....	462	G4525-60701.....	50		
G3900-63034.....	463	G4525-60702.....	50		
G3900-63035.....	463	G4525-60703.....	50		
G3900-63036.....	463	G4525-60704.....	50		
G3900-63037.....	463	G4556-60019.....	52		
G3900-63038.....	462	G4556-60125.....	52		
G3900-63039.....	462	G4556-67010.....	52		
G3900-63040.....	462	G4556-90500.....	52		
G3900-63041.....	462	G6012A.....	169		
G3900-63042.....	462	G6014A.....	169, 179		
G3900-64016.....	450	G6600-60037.....	123		
G3900-65001.....	448	G6600-60038.....	122		
G3900-65003.....	449	G6600-67007.....	122-123		
G3900-65004.....	449	G6600-67008.....	122-123		
G3903-61004.....	448	G6600-80018.....	41, 122-123		
G3903-63001.....	284, 374	G6600-80042.....	122		
G3903-63002.....	374	G6600-80043.....	122-123		
G3903-63003.....	283, 374	G6600-80044.....	122-123		
G3903-63004.....	374	G6600-80045.....	122		
G3903-65002.....	448	G6600-80050.....	122		
G4333-63000.....	132	G6600-80051.....	122-123		
G4513-20561.....	50	G6600-80063.....	122		
G4513-40525.....	50	G6600-80072.....	122-123		
G4513-40529.....	50	G6600-85000.....	123		
G4513-60560.....	188	G6600-85001.....	122-123		
G4513-60561.....	188	G6600-85002.....	122-123, 170		
G4513-80200.....	187	G7005-60061.....	154-155, 162, 174		
G4513-80201.....	187	G7005-60072.....	159, 162, 174		
G4513-80202.....	187	MKI-C-QSC10.....	65		
G4513-80203.....	188	MKI-MTD-1169.....	65		
G4513-80204.....	187	MKI-MTD-1204.....	65		
G4513-80205.....	187	MKI-SERUTD-5065.....	66		
G4513-80206.....	187	MKI-U-COV06.....	65		
G4513-80208.....	188	MKI-U-COV07.....	65		
G4513-80209.....	23, 187	MKI-U-COV10.....	65		
G4513-80210.....	188	MKI-U-DISK1.....	65		
G4513-80211.....	187	MKI-U-DISK3.....	65		
G4513-80212.....	187	MKI-U-T10CW-2S.....	65		
G4513-80213.....	23, 187	MKI-U-T11GPC.....	65		
G4513-80215.....	187	MKI-U-T11GPC-2S.....	65		
G4513-80216.....	187	MKI-U-T12ME.....	65		
G4513-80218.....	187	MKI-U-T12ME-2S.....	65		
G4513-80219.....	188	MKI-U-T13DHS.....	65		
G4513-80220.....	188	MKI-U-T13DHS-2S.....	65		
G4513-80221.....	188	MKI-U-T15ATA-2S.....	65		
G4513-80222.....	188	MKI-U-T16GHG-2S.....	65		
G4513-80223.....	188	MKI-U-T1703P-2S.....	65		
G4513-80224.....	187	MKI-U-T1HBL-2S.....	65		
G4513-80225.....	187-188	MKI-U-T3ATX.....	65		
G4513-80226.....	187	MKI-U-T3ATX-2S.....	65		
G4513-80227.....	188	MKI-U-T6SUL-2S.....	65		
G4513-80228.....	188	MKI-U-T9TNX.....	65		
G4513-80229.....	187	MKI-U-T9TNX-2S.....	65		

Product Index

Agilent J&W GC Columns

Application and Method Guides

ASTM Methods	278-285
Application Cross Reference Guide	262-264
EPA Methods	
Drinking Water	265-268
Solid Waste	272-276
Waste Water	269-271
United States Pharmacopoeia (USP)	
GC Phases	277

Capillary Columns

By Group

Foods, Flavors and Fragrances	405
High Temperature	362
Life Sciences	414
Low-bleed GC/MS	294
Non-Bonded	441
PAH	390
PLOT	420
Pesticides	384
Petroleum	366
Polyethylene Glycol (PEG)	351-413
Premium Polysiloxane	318
Semivolatiles	392
Ultra Inert	286
Volatiles	397

By Stationary Phase

Biodiesel	377-378
CAM	442
CarboBOND	436
CarboPLOT P7	436
Carbowax 20M and HP-20M	361
CP-1301	348
CP-Al2O3/KCl	430-431
CP-Al2O3/Na2SO4	430-431
CP-Carbowax 400 for Volatiles in Alcohol	408
CP-Chirasil Val	412
CP-Chirasil-Dex CB	412
CP-Cyclodextrin- β -2,3,6-M-19	413
CP-FFAP CB for Free Fatty Acids	
in Dairy Products	410
CP-Molsieve 5A	439
CP-Select 624 CB	399
CP-Select CB for MTBE	403
CP-Sil 13 CB	336
CP-Sil 19 CB	350
CP-Sil 19 CB for Pesticides	387
CP-Sil 2 CB	373
CP-Sil 24 CB	341
CP-Sil 43 CB	346
CP-Sil 5 CB	324-326
CP-Sil 5 CB for Formaldehyde	367
CP-Sil 5 CB for Sulfur	375
CP-Sil 5/C18 CB for PCB	395
CP-Sil 8 CB	334-335
CP-Sil 8 CB for Amines	382
CP-Sil 8 CB for PCB	393
CP-Sil 8 CB for Pesticides	387
CP-Sil 88	406
CP-Sil 88 for Dioxins	396
CP-Sil 88 for FAME	407
CP-Sil PAH CB UltiMetal	391

CP-Sil PONA CB	368
CP-Sil PONA for ASTM D5134	368
CP-SilicaPLOT	435
CP-SimDist	371-372
CP-TAP CB for Triglycerides	409
CP-TCEP for Alcohols in Gasoline	373
CP-Volamine	381
CP-Wax 52 CB	356-357
CP-Wax 57 CB	408
CP-Wax 57 CB for Glycols and Alcohols	409
CP-Wax 58 FFAP CB	360
CP-Wax for Volatile Amines and Diamines	382
CycloSil-B	410
Cyclodex-B	411
DB-1	318-321
DB-1301	347
DB-17	339
DB-1701	349
DB-1701P	386
DB-17ht	364
DB-17ms	308
DB-1ht	362
DB-1ms	295
DB-1ms Ultra Inert	289
DB-200	343
DB-210	344
DB-225	345
DB-225ms	312
DB-23	342
DB-2887	370
DB-35	337
DB-35ms	306
DB-35ms Ultra Inert	291
DB-5	329-331
DB-5.625	394
DB-502.2	402
DB-5ht	363
DB-5ms	299-300
DB-5ms EVDX	415
DB-5ms Ultra Inert	290
DB-608	388
DB-624	398
DB-624 Ultra Inert	291, 397
DB-ALC1	414
DB-ALC2	414
DB-CLP1	384
DB-CLP2	384
DB-Dioxin	396
DB-EUPAH	390
DB-FFAP	358
DB-HT SimDis	370
DB-MTBE	402
DB-Petro	369
DB-Select 624 UI for <467>	292, 416
DB-Sulfur SCD	374
DB-TPH	403
DB-UI 8270D Ultra Inert	292
DB-UI 8270D for Semivolatiles	392
DB-VRX	400
DB-WAX	351-353
DB-WAX FF	351-353
DB-WAXetr	354
DB-XLB	304

DX-1	443
DX-4	443
GS-Alumina	433
GS-Alumina KCl	429
GS-CarbonPLOT	437
GS-GasPro	434
GS-OxyPLOT	366
GS-Q	426
HP-1	322-323
HP-1 Aluminum Clad	369
HP-101	441
HP-17	441
HP-1ms	296
HP-1ms Ultra Inert	289
HP-20M	361
HP-35	338
HP-5	332-333
HP-50+	340
HP-5ms	301
HP-5ms Semivolatile	395
HP-5ms Ultra Inert	290
HP-88	405
HP-Chiral β	411
HP-FFAP	359
HP-Fast Residual Solvent	416
HP-INNOWax	355
HP-PAS5	388
HP-PLOT Al2O3 KCl	428
HP-PLOT Al2O3 M	434
HP-PLOT Al2O3 S	432
HP-PLOT Molesieve	438
HP-PLOT Q	425
HP-PLOT U	428
HP-PONA	367
HP-VOC	401
Lowox	366
PoraBOND Q	422
PoraBOND U	423
PoraPLOT Amines	383
PoraPLOT Q	424
PoraPLOT Q-HT	424
PoraPLOT S	427
PoraPLOT U	427
Rapid-MS	389
SE-30	444
SE-54	444
Select Al2O3 MAPD	376
Select Biodiesel	379
Select FAME	407
Select Low Sulfur	375
Select Mineral Oil	404
Select PAH	390
Select Silanes	380
Select for Permanent Gases	376
Ultra 1	327
Ultra 2	328
VF-1301ms	315-316
VF-1701 Pesticides	386
VF-1701ms	317
VF-17ms	309
VF-1ms	297-298
VF-200ms	311
VF-23ms	310

VF-35ms.....	307	Carboxen-1000.....	471	MSD Contamination.....	142-145
VF-5 Pesticides.....	385	Chromosorb 101.....	471	MSD Electron Multipliers and Replacement	
VF-5ht.....	365	Chromosorb 102.....	471	Horn.....	165
VF-5ms.....	302-303	HayeSep A.....	473	MSD Filaments.....	162
VF-624ms.....	315	HayeSep D.....	473	Maintenance Schedule.....	141
VF-DA.....	415	HayeSep DB.....	473	Quadrupole Mass Filter.....	164
VF-WAXms.....	313-314	HayeSep N.....	473	Recommended MS Interface Connections.....	148
VF-Xms.....	305	HayeSep N + HayeSep R 1:1.....	473	Standards	
Metal Columns.....	417	HayeSep P.....	473	Analyzer Kit Standards.....	181
Column Selection		HayeSep Q.....	474	Test and Performance Samples.....	182
Column Diameter.....	256-257	HayeSep R.....	474	Vacuum Systems and Pumps.....	166
Column Film Thickness.....	260-261	HayeSep T.....	474	Diffusion Pump.....	168
Column Length.....	258-259	MolSieve 13X.....	476	Foreline Pump.....	170
Gas-Solid or PLOT Columns.....	254	MolSieve 5A.....	475	Pressure Symptoms.....	166-167
Polarity.....	252-253	Porapak N.....	477	Quiet Cover.....	169
Selecting Stationary Phases.....	248-251	Porapak N + Porapak R 1:1.....	477	Vent Valve Supplies.....	163
Stationary Phase Selection Summary.....	255	Porapak Q.....	478	Inlet Liners.....	24-25
Custom Ordering.....	481	Porapak QS.....	478	Specialty Injection	
Fused Silica Tubing		Porapak R.....	478	Direct Connect.....	32
Deactivated.....	464-467	Porapak T.....	478	MultiMode Inlet Heavy Matrix.....	31
Undeactivated.....	468	Silica Gel.....	479	PTV.....	32
Guard Columns		Test Standards.....	482-483	Split, Original Deactivation.....	29
DuraGuard.....	445			Splitless, Original Deactivation.....	30
EZ-Guard.....	446			Ultra Inert.....	28
Installation and Troubleshooting				Inlet Septa.....	17
Causes of Column Performance				Advanced Green.....	19
Degradation.....	488-493			Bleed and Temperature Optimized (BTO).....	18
Evaluating the Problem.....	494-495			General Purpose.....	21
Quick Reference Guide.....	485-487			Long-Life.....	19
Troubleshooting Guides				Inlet Systems	
Baseline Instability or Disturbances.....	496			Cool On-Column.....	78-80
Change in Peak Size.....	498			Flip Top Inlet Sealing System.....	68
Excessive Baseline Noise.....	496			Multimode.....	76-77
Loss of Resolution.....	499			Programmable Temperature Vaporizer (PTV).....	81-87
Retention Time Shift.....	498			Purged Packed.....	88-91
Split Peaks.....	497			Split/Splitless Inlet Seals.....	67, 69-75
Tailing Peaks.....	497			Large Valve Oven.....	46
LTM Column Modules.....	447-463			Liner O-rings.....	33
Packed Columns				Maintenance Schedule.....	14-15
0.1% SP-1000.....	479			Markes Thermal Desorption.....	65-66
1.5% OV-101.....	476			Merlin Microseal.....	22-23
10% Carbowax 20M (G16, GS1).....	471			Oven Exhaust Deflector.....	123
10% Carbowax 20M (G16, GS1) + 2% KOH.....	471			Purge and Trap Systems	
10% OV-101.....	476			Archon.....	64
10% PEG-20M.....	477			Teledyne Tekmar.....	62-63
10% SE-30.....	478			Sample Introduction Systems	
10% UC W982.....	480			7693A Automatic Liquid Sampler.....	50
12% UC W982.....	480			7697A Headspace Sampler.....	52-53
15% Carbowax 1540.....	470			Automatic Liquid Sampler.....	51
15% Hallcomid M-18.....	472			G1883A Network Headspace Sampler.....	55
15% SP-2100.....	479			G1888A Network Headspace Sampler.....	54
20% Carbowax 20M (G16, GS1).....	471			G3520A XLSI Accessory.....	53
20% OV-101.....	476			Syringes.....	183
20% PEG-20M.....	477			7673/7683 On-Column Autosampler Syringes.....	191
20% Sebaconitrile.....	479			Blue Line Autosampler Syringes for 7693A ALS	
20% Sebaconitrile/2% H3PO4.....	479			Advanced Sample Enhancement.....	188
20% TCEP.....	480			Fitted Plungers.....	187
25% DC-200 (500 cSt).....	472			PTFE-Tipped Plungers.....	188
25% SP-2100.....	479			Gold Standard Autosampler Syringes	
30% DC 200/500.....	472			Straight Needle, 23 and 26s Gauge.....	190
30% DC-200 (500 cSt).....	472			Tapered Needle, 23-26s Gauge.....	189
35% DC-200 (500 cSt).....	472			HP 7670/71/72 Autosampler Syringes.....	191
5% Carbowax 20M (G16, GS1).....	470			Needle Gauge.....	184
7% Carbowax M + 3% Polyphenolether				Needle Termination.....	185
6 ring + 2% KOH.....	471			Needle Tip Design.....	185
Carbosieve S-II.....	470			Valves and Loops.....	47-48

Vials and Closures

Headspace	
20 mm Caps and Septa.....	58
CombiPAL 18 mm Screw Top Caps with Septa.....	59
CombiPAL Screw Top Vials.....	59
Convenience Packs.....	58
Crimp Top Vials.....	56
High Performance Septa.....	57-58
Electronic Crimpers and Decappers.....	60
Manual Crimpers and Decappers.....	61

Applications

Energy and Fuels.....	576
Environmental	
Air Analysis.....	549
Hydrocarbons.....	501
Pesticides and Herbicides.....	506
Semivolatiles.....	532
Volatiles.....	544
Food, Flavor and Fragrance.....	554
Industrial Chemical.....	602
Forensic Toxicology and Pharma.....	635

CrossLab Supplies

Bruker, Varian GC Systems	
Autosampler Syringes.....	210
Capillary Column Ferrules.....	209
Column Nuts.....	210
Detector Replacement	
Flame Ionization Detector (FID).....	218
Pulsed Flame Photometric Detector (PFPD).....	219-220
Thermal Conductivity Detector (TCD).....	218
Thermionic Specific Detector (TSD).....	220
Injector Replacement	
1041 Packed/Wide Bore	
On-Column (PWOC) injector.....	217
1061 Packed/530 µm	
Capillary Column Injector.....	216
1079 Large Volume Injector (LVI).....	214
1093 Cool On-Column (COC) Injector.....	215
1177 Split/Splitless Injector.....	213
Inlet Liners	
1060/1061 Injector.....	208
1075/1077 Injector.....	208
1078/1079 Injector.....	207-208
1093/1094 Injector.....	207-208
1177 Injector.....	206-220
Inlet Septa.....	211-212
Liner O-rings.....	208
Packed Column Ferrules.....	210
CTC GC Autosamplers	
Autosampler Syringes.....	239
Column Ferrules.....	195-196
Inlet Liners.....	193-194
Inlet Septa.....	199-201
Advanced Green.....	203
Bleed and Temperature Optimized (BTO).....	202
General Purpose.....	205
Long-Life.....	204
Liner O-rings.....	195

Perkin Elmer GC Systems

Autosampler Syringes.....	225
Capillary Column Ferrules.....	223
Column Nuts.....	224
Inlet Liners	
AutoSystem.....	221-222
AutoSystem XL.....	221-222
Clarus.....	221-222
Inlet Septa.....	226
Liner O-rings.....	222
Packed Column Ferrules.....	224

Shimadzu GC Systems

Autosampler Syringes.....	232
Capillary Column Ferrules.....	231
Column Nuts.....	232
Inlet Liners	
14 Systems.....	230
17A Systems.....	229
2010 and 2010 Plus Systems.....	228
2014 Systems.....	227
Inlet Septa.....	233
Liner O-rings.....	230
Packed Column Ferrules.....	231
Syringes.....	197-198

Thermo Scientific GC Systems

Autosampler Syringes.....	236
Capillary Column Ferrules.....	235
Column Nuts.....	236
Inlet Liners	
Focus Systems.....	234
Trace Systems.....	234
Inlet Septa.....	237-238
Liner O-rings.....	235

Application Title Index

1,3-Butadiene	586
1,3-Butadiene Purity	587
10 ng/ μ L Semivolatile Checkout Standard on a 20 m x 0.18 mm, 0.36 μ m Agilent J&W DB-UI 8270D Capillary GC Column using an Ultra Inert Liner with Wool	542
15+1 EU Priority PAHs	505
69 Component FAME Mix	572
78 Semi-volatile Components on an Agilent J&W DB-UI 8270D	504

A

Acids	570
Acrylate Impurities I	627
Acrylate Impurities II	628
Acrylates	628
Agilent's Ultra Inert Test Probe Mixture	532
Alcohol Beverage Standard	568
Alcohols I	602
Alcohols II	605
Alcohols III	606
Aldehydes and Acids	611
Aldehydes and Ketones	612
Alditol Acetates	568
Amines and Alcohols	606
Amines and Nitriles	610
Amines in Water	611
Amphetamines and Precursors – TMS Derivatives	636
Anabolic Steroids	648
Analysis of Acetylenes' Mixture	634
Analysis of Amino Alcohols in Water	606
Analysis of Drugs of Abuse in Urine via GC/MS	642
Analysis of Ethanolamines	607
Analysis of Fragrance and Allergens	561
Analysis of Oxygenates in Mixed C4 Streams	595
Analysis of Oxygenates in a C1 to C5 Hydrocarbon Mix	600
Analysis of Polycyclic Aromatic Hydrocarbons	503
Analysis of Process Gas	601
Analysis of Semivolatiles	510
Analysis of Solvents	626
Analysis of Volatile Organic Compounds in Environmental Waters Using the Agilent 7697A Headspace and 7890B/5977A GC/MS	546
Anesthetics	643
Anilines	629
Anticonvulsants	643
Antiepileptic Drugs	644
Antihistamines	644
Aroclors 1016-1268 (without 1221)	514
Aromatic Solvents	604
Aromatics Analysis – ASTM D16 Analytes	593
Aromatics Analysis – Ethylbenzene Impurities	593
Aromatics I	614
Aromatics II	615
Aromatics in Finished Gasoline – ASTM Method D5769	597
Aspirin and Ibuprofen in Methanol	647
Automated Cleanup of PCB extracts from Waste Oil Using 7696A Sample Prep Workbench	535

B

Bacterial Fatty Acid Methyl Esters	571
Barbiturates	637
Baseline Resolution of Air/ CO , CO_2 , and Methane in a Natural Gas Sample	583
Benzodiazepines I	636
Benzodiazepines II	640
Blood Alcohols I (Static Headspace/Split)	638
Blood Alcohols II (Static Headspace/Split)	638
Blood Pollutants I	649
Blood Pollutants II	650
Bourbon	568
Butter Triglycerides I	574
Butter Triglycerides II	574

C

C ₁ to C ₄ Hydrocarbon Mix	580
C ₁ and C ₂ Halocarbons (Freons)	527
CLP Pesticides	514
Canola Oil Margarine Partially Hydrogenated FAMES AOCS Method 1c-89	574
Chiral Compounds in Essential Oils and Fragrances	560
Chlorinated Isooctane	623
Chlorinated Pesticides, EPA Method 508	521
Citrus Flavored Carbonated Beverage (Soda)	567
CLP Pesticides	514
Cold-pressed Orange Oil	564
Column Performance for USP <467> Standards	651
Common Drug Screen	641
Congeners in DIN Method PCBs	533

D

DB-624UI 1 μ L/L Fermented Beverage Standard Mix	554
DB-624UI Organic Acid Performance	509
DB-Select 624 UI for <467> Megabore Early Eluting Peaks	635
Denatured Fuel Ethanol – ASTM D5501	596
Detailed Hydrocarbon Analysis of Petroleum Naphthas Through N-nonane Using ASTM D5134	601
Determination of Chlorophenols in Water and Soil	502
Diesel Analysis	600
Diesel Fuel	503
Dioxins and Dibenzofurans	504
Direct Comparison for Rapid CLP (Contract Laboratory Program) Pesticide Analysis	513
Drug Screen	640

E

EPA 625 Halogenated Pesticides on "1701" Type Phases	528
EPA Air Analysis Compendium Method TO-14 Standard	549

EPA Air Analysis Method TO-15 (1 ppbv standard)	553
EPA Method 504.1 – 1,2 dibromoethane (EDB), 1,2-dibromo-3-chloropropane (DBCP), and 1,2,3-trichloropropane (123TCP)	507
EPA Method 508.1 – Chlorinated Pesticides and Herbicides	520
EPA Method 525.2	518
EPA Method 551	547
EPA Method 551 – Chlorinated Solvents, Trihalomethanes (THMs), and Disinfection Byproducts (DBPs)	509
EPA Method 552.2	543
EPA Volatiles by GC/MS (Split Injector)	517, 548
Essential Oils	557
Esters I	617
Esters II	618
Esters III	618
Ethers	619
Ethoxyethanol	607
Ethylene	584
Ethylene Glycol Mixture	621
Ethylene Oxide	631
Ethylene Oxide Synthetic Standard	594
European Red List Volatiles	547
Extended Analyte List for EPA Method 8021 (ELCD)	544
Extended Hydrocarbon Analysis I	588
Extended Hydrocarbon Analysis II	589
Extended Temperature Program Resolving Congeners 52 and 138	534

F

FAME Standard	572, 573
FAMES	561
Fast Analysis of Aromatic Solvent	576
Fast Analysis of Permanent Gases and CO_2 using Tandem PLOT columns	552
Fast CLP Pesticides	506
Fast Screening of FAME Isomers in Butter	575
Fast Separation of Silanes	633
Fast VOC Analysis	545
Fentanyl	645
Flavor Mixture	563
Formaldehyde Underivatized	613
Formaldehyde, 50 ppb	550
Formaldehyde-DNPH Derivative	613
Fragrance Allergens	563
Fragrance Reference Standard	558, 559
Free Organic Acids/C ₄ -C ₅ Isomers	608
Free Phenols	528
Free Steroids	648

G

Gasoline Unleaded ASTM D5769	599
Glycols I	619
Glycols II	620
Glycols III	620
Glycols/Diols	622

H			
Hallucinogens	646	Organochlorine Pesticides, EPA Method 8081B.....	508
Halocarbons	631	Organophosphorous Pesticides in Apple Matrix.....	531
Halogenated Hydrocarbons I	603	Organophosphorous Pesticide Residues	
Halogenated Hydrocarbons II	622	in Olive Oil Extract.....	562
Halothane	632	Organotin Compounds I.....	537
Herbicides I	526	Organotin Compounds II.....	537
Herbicides II	526	Over-the-Counter Pain Killers – TMS Derivatives	647
High Resolution Phenol Analysis by GC/MS	542	Oxygenates in Gasoline ASTM D5599 (GC-OFID)	595
High Resolution Separation of Xylene Isomers	632		
High Speed VOC, EPA Method 8260	515		
		P	
I		PAHs	540
Impurities in Ethylbenzene	616	PBDEs.....	516
Impurities in Ethylene.....	584	PBDEs by ECD.....	502
Impurities in Mixed Xylenes	632	PCBs by EPA Method 8082.....	534
Impurities in Propylene.....	585	PFBHA Derivative.....	614
Impurities in Styrene.....	616	PONA Mix as Specified by	
Impurities in p-Xylene – ASTM D3798	594	AFNOR Method #2.....	596
Inorganic Gases	605	Peppermint Oil.....	565
Inorganic Hydride Gases.....	633	Perfume.....	560
		Permanent Gases	582
L		Permanent Gases on a Thick Film	
Lavender Oil Characterization.....	556	Molsieve Column.....	552
Lemon Oil	564	Pesticides and Fire Retardants (US EPA 527)	519
		Pesticides in Sunflower Oil	575
M		Pesticides, EPA 508.1	511
Marijuana ($\Delta 9$ -THC) and Major Metabolites – TMS		Phenols	541
Derivatives.....	649	Phenols According to EPA Method 8040	543
Menthol.....	560	Phenols I.....	604
Mercaptans	592	Phenols II.....	630
		Phenols III.....	630
N		Phenoxy Acid Herbicides – Methyl Derivatives,	
N_2O I	551	EPA 8151A.....	512
N_2O II.....	551	Polybrominated Diphenyl Ethers (PBDEs).....	505
N_2O III.....	551	Polyethylene	599
Narcotics	637	Polyethyleneamines.....	609
Narcotics and Adulterants	647	Polymer Additives.....	633
Natural Gas	583	Primary Amines.....	609
Nitrogen Containing Herbicides		Propylene.....	585
(EPA Method 507)	527	Pyrethrins	536
Nitrogen-based Solvents I.....	626	Pyrolysates of Polystyrene.....	617
Nitrogen-based Solvents II.....	627		
Nitrogen/Phosphorus Containing		R	
Pesticides, EPA Method 507	525	Reference Gas Oil.....	598
Noble Gases	582	Refinery Gas.....	590
n-Paraffin Standard	578	Refinery Gas I.....	577
		Regular Unleaded Gasoline	
O		(California Phase 1) – "Normal" GC Run I.....	598
Organic Acids	570, 607	Regular Unleaded Gasoline	
Organochlorine Pesticides.....	521, 525	(California Phase 1) – "Normal" GC Run II.....	598
Organochlorine Pesticides I		Residual Solvents, DMI Diluent	639
EPA Method 8081A.....	529	Residual Solvents, USP 467	650
Organochlorine Pesticides II		Rosemary Oil.....	567
EPA Method 8081A.....	530		
Organochlorine Pesticides III.....	522	S	
Organochlorine Pesticides IV	523	Sedative Hypnotics.....	646
Organochlorine Pesticides		Selected Oxygenates	582
to EPA 625 via GC/MS.....	529	Semivolatile Compounds, US EPA Method 8270.....	538
Organochlorine Pesticides, DB-5/DB-1701P.....	524	Separation of TMS-derivatized Sugars	
		using VF-1ms.....	569
		Separation of cis-trans FAME Isomers	571
		Simulated Distillation.....	597
		Solvents	625
		Solvents I.....	623
		Solvents II.....	624
		Solvents III.....	624
		Solvents IV.....	625
		Spearmint Oil.....	555
		Spearmint Oil (Western)	565
		Strawberry Syrup	569
		Substituted Anilines.....	629
		Sulfur Compounds in Naphtha.....	593
		Sulfur Compounds in Natural Gas –	
		Synthetic Mixture	592
		Sulfur Compounds in Propylene (1 ppm).....	579, 591
		Sulfur Gas Analysis in Light	
		Hydrocarbon Streams I	590
		Sulfur Gas Analysis in Light	
		Hydrocarbon Streams II	591
		Sulfur Gases.....	634
		Sulfur Impurities in Propylene.....	579
		Sulfur Standards in Toluene.....	578
		Sulfur in Air.....	550
		T	
		Tetrachlorodibenzo-p-furans	533
		Tocopherols	645
		Trace Active Amines, 10 ng on-column	608
		Trace Level Polycyclic Aromatic	
		Hydrocarbon (PAH) Analyses	532
		Trace Oxygenates in Light	
		Hydrocarbon Matrices.....	581
		Trace Sulfur Compounds in Methane (50 ppbv).....	581
		Tricyclic Antipsychotics.....	644
		Triethylene Glycol and Impurities	621
		U	
		US EPA Method 8061 (Phthalate Esters)	539
		Underivatized Drugs of Abuse –	
		Agilent Fast Toxicology Analyzer	639
		Unleaded Gasoline	501, 577
		Urine Drug Screen	642
		V	
		Volatile Amines.....	608
		Y	
		Ylang Ylang Oil	566

Compound Index

A

Acenaphthalene	540
Acenaphthalene-d10	519
Acenaphthene	503, 532, 538, 540
Acenaphthene-d10	518, 538, 542
Acenaphthylene	503, 518, 532, 538
Acephate	531, 562
Acetal (acetaldehyde diethyl acetal)	619
Acetaldehyde	550, 553, 568, 581, 600, 612, 625, 638, 649-650
Acetaminophen	641, 647
Acetic acid	509, 568, 570, 607, 611
Acetone	515, 517, 548, 550, 553-554, 558-559, 568, 570, 581-582, 600, 612, 624-626, 632, 638-639, 649-650
Acetone-d6	550
Acetonitrile	515, 624-627, 635, 639, 649-651
Acetyl aldehyde	554
Acetylcedrene	560
6-Acetylcodeine	637, 647
Acetylene	580, 584-590
Acetylsalicylic acid (aspirin)	647
Acifluorfen	512
Acrolein	515, 612, 626-627
Acrylamide	610
Acrylic acid	509
Acrylonitrile	515, 517, 548, 626-627
Air/CO	577, 583, 601
Alachlor	508, 511, 518, 520, 524-526, 529-530
Aldrin	506, 508, 511, 513-514, 518, 520-525, 529-530, 542
Alfentany	645
Allobarbitol	637, 641
Allyl acrylate	628
Allyl alcohol	515, 581
Allyl butyrate	558-559
Allyl chloride	515, 517, 544, 548
Allyl ether	619
Allyl ethyl ether	619
Alpha isomethyl ionone	563
Alphenal	637
Alprazolam	194, 636, 639-640
Ametryn	518, 525-526
2-Aminoazotoluene	629
2-Amino-1-butanol	606
2-Amino-ethanol	606
Aminoethylethanolamine	611
n-(2-Aminoethyl) piperazine	611
2-Aminonaphthalene	629
1-Amino-4-nitronaphthalene	540
5-Amino-1-pentanol	606
1-Amino-2-propaol	606
Amitriptyline	642
Amobarbital	637, 641
Amphetamine	636, 639, 642
Amyl acetate	554, 569, 617-618
Amyl alcohol	554
n-Amyl alcohol	568
Amyl butyrate	569
Amyl cinnamic alcohol	561
Amyl cinnamyl alcohol	563
Amyl cinnamyl aldehyde	561
n-Amyl mercaptan	592

tert-Amyl mercaptan	592
tert-Amyl methyl ether (TAME)	515, 546, 581, 595, 599, 600, 619
n-Amyl salicylate	558-560
5 α -Androstan-17 α -ol-3-one (stanolone)	648
Androsterone	648
5- β -Androsterone	648
Aniline	538, 542, 610, 629
Anisic alcohol	561
Anisyl alcohol	563
Antazoline	644
Anthracene	503, 518, 532, 538, 540
Anthracene-d10	519
Aprobarbital	637, 641
Arabinitol	568
Arabitol	569
Arachidic acid (eicosanoic acid)	570
Arachidic acid methyl ester	572, 573
Arachidonic acid methyl ester	572, 573
Argon	552, 582
Arsine	633
Atraton	525-526
Atrazine	511, 518-520, 525-527, 542
Azinphos-ethyl	562
Azinphos methyl	531, 562
Azobenzene	538
Azulene	540

B

Balan	527
Barbital	637
BDE	505, 519
Behenic acid methyl ester	572, 573
Benactyzine	641
Bentazone	512
Benthiocarb	519
Benz[a]anthracene	505, 518, 532
Benz[a]anthracene-7,12-dione	540
Benzaldehyde	558-559, 563, 612
Benzaldehyde, 3 methoxy	554
1,2-Benzanthracene	540
Benzene	501, 515, 517, 544-546, 548-549, 553, 576-577, 586-589, 593, 596-597, 599, 601, 604, 614-616, 623-626, 632, 635, 649-651
Benzene-d6	597
Benzene ethanol	558-560
Benzidine	538, 629
Benzo[a]anthracene	503, 538
Benzocaine	643
Benzo[b]fluoranthene	503, 505, 518, 532, 538, 542
Benzo[j]fluoranthene	505
Benzo[k]fluoranthene	503, 505, 518, 532, 538, 540
Benzo[b]fluoroanthene	540
Benzo[c]fluorene	505
2,3-Benzofluorene	540
Benzo[g,h,i]perylene	503, 505, 518, 532, 538, 540
Benzo[a]pyrene	503, 518, 532, 538, 540
Benzo[e]pyrene	540
Benzoic acid	538
Benzonitrile	610, 626-627
Benzophenone	558-559
5,6-Benzoquinoline	540

Benzothiophene	593
Benzphetamine	636, 641
Benzo[a]pyrene	505
Benzyl acetate	558-560, 566, 617-618
Benzyl alcohol	538, 558-559, 561, 563, 602, 605
Benzylamine	609-610
Benzyl benzoate	539, 558-561, 563, 566, 569
Benzyl butyl phthalate	518
Benzyl chloride	515, 544, 553
Benzyl cinnamate	561, 563
Benzyl ether	619
n-Benzylmethylamine	610
Benzyl salicylate	558, 560-561, 563, 566
α -Bergamotene	556
trans- α -Bergamotene	564
α -BHC	506, 508, 511, 513-514, 518, 520-525, 529-530
β -BHC	506, 508, 511, 513-514, 518, 520-525, 529-530
δ -BHC	506, 508, 511, 513-514, 518, 520-525, 529-530
γ -BHC	506, 508, 511, 513-514, 518, 520-525, 529-530
BHEB	633
BHT	561, 567, 633
Bifenthrin	519
Bioallethrin	519
Biphenyl	540
β -Bisabolene	559, 564, 566
α -Bisabolol	556
Borneol	556, 558-560, 567
Borneol acetate	556
Bornyl acetate	557-558
α -Bourbonene	555
β -Bourbonene	555, 565
BPE (PB)	536
Bromacil	518-519, 525-527
Bromazepam	636, 640
Bromoacetic acid	543
Bromoacetone	515
4-Bromoaniline	610, 629
Bromobenzene	515, 517, 544-546, 548, 614-615
2-Bromobiphenyl	529-530
Bromochloroacetic acid	543
Bromochloroacetonitrile	509, 547
Bromochlorodifluoromethane	527
Bromochloromethane	515, 517, 544-546, 548-550, 553, 622
Bromodichloroacetic acid	543
Bromodichloromethane	507, 509, 515, 517, 544-548, 553, 622
2-Bromo-4,6-dinitroaniline	629
Bromoethane (ethyl bromide)	553, 622
Bromofluorobenzene	545
4-Bromofluorobenzene	515, 517, 550, 548-549, 553
Bromoform	507, 509, 515, 517, 544-548, 553, 622
Bromofos	575
Bromomethane	515, 517, 544-546, 548-549, 553
1-Bromo-2-nitrobenzene	508, 529-530
Bromopheniramine	641, 644
4-Bromophenyl-phenylether	538
3-Bromopyridine	610
Bromotrifluoromethane	527
Bucizine	644
Bufotenine	646
Butabarbital	637, 641
Butacaine	643
Butachlor	518, 520, 525-526

1,2-Butadiene	586-587, 634	Butyl ethyl ether	619	Cetearyl octanoate	558-559
1,3-Butadiene	553, 584-591, 634	Butyl heptanoate	618	Chloral hydrate	509, 547
Butalbitol	637	1-Butyl mercaptan	592	Chloramben	512
Butanal	611	n-Butyl mercaptan	592	Chlorcyclizine	641, 644
1,3-Butandiol	619, 622	tert-Butyl mercaptan	581, 592	α -Chlordane	506, 508, 511, 513-514, 518, 520-524, 529-530
1,4-Butandiol	619	sec-Butyl methacrylate	628	γ -Chlordane	506, 508, 511, 513-514, 518, 520-524, 529-530
2,3-Butandiol	619	Butyl methyl ether	619	Chlordiazepoxide	636
Butane	599, 601	tert-Butyl methyl ether (MTBE)	553	Chloroacetic acid	543
iso-Butane	583, 590-591	Butyl propionate	618, 627-628	Chloroacetonitrile	515, 517, 548
n-Butane	577, 580, 583-590, 596, 634	4-tert-Butyltoluene	614-615	2-Chloroaniline	610, 629
n-Butane/cis-2-butene	591	Butyl valerate	618	3-Chloroaniline	610, 629
1,3-Butanediol	620	2-Butyne (dimethylacetylene)	586-587	4-Chloroaniline	538, 610, 629
1,4-Butanediol	620	1-Butyne (ethylacetylene)	586-587	Chlorobenzene	515, 517, 544-546, 548-549, 553, 604, 614-615, 623, 625, 651
2,3-Butanediol	620	2-sec-Butyl-4,6-dinitrophenol (dionseb)	542, 543	Chlorobenzene-d ₅	550, 553, 549
2,3-Butanedione (VDK)	554	Butyraldehyde	600, 612	Chlorobenzilate	508, 511, 518, 520-521, 524, 529-530
2,3-Butanedione (diacetyl)	558-559	Butyric acid	570, 607-608, 611	4-Chlorobenzonitrile	610
1-Butanethiol	578-579, 581, 591	iso-Butyric acid	611	2-Chlorobiphenyl	518
Butanol	626-628	n-Butyric acid	509	1-Chlorobutane	515, 517, 548, 603, 622, 649-650
1-Butanol	515, 554, 595, 600, 602, 605-606, 625, 649-650	Butyric acid methyl ester	572, 573	2-Chlorobutane	649-650
2-Butanol	554, 635, 639			Chlorodibromoacetic acid	543
iso-Butanol	606			1-Chloro-1,1-difluoroethane	527
n-Butanol	568, 624, 650			Chlorodifluoromethane	527
sec-Butanol	568, 595, 602, 605, 624, 650			2-Chloro-4,6-dinitroaniline	629
tert-Butanol	595, 602, 605, 624, 650			Chloroethane	515, 517, 527, 544-546, 548-549, 553
2-Butanone (MEK)	515, 517, 548, 550, 553, 600, 612, 626, 639			2-Chloroethanol	515, 631
2-Butanthiol	581			bis(2-Chloroethoxy) methane	538
Butene	617			bis(2-Chloroethyl) ether	538
Butene-1	584, 586-587			2-Chloroethyl vinyl ether	544
1-Butene	584-585, 588-591, 617, 634			1-Chloro-3-fluorobenzene	517, 544, 548
cis-2-Butene	584-590			Chloroform	507, 509, 515, 517, 544-549, 553, 603, 622, 624, 635, 649-650
trans-2-Butene	584-591			1-Chlorohexane	515, 603, 622
1-Butene/methyl acetylene	590			1-Chloro isooctane	623
3-Buten-1-ol	602, 605			3-Chloro isooctane	623
2-Buten-1-ol (crotyl alcohol)	602, 605			4-Chloro isooctane	623
Butethal	637			bis(2-Chloroisopropyl) ether	538, 544
Butylpentytin	537			Chloromethane	515, 517, 527, 544-546, 548-549, 553
2-Butoxyethanol	602, 605			4-Chloro-2-methylaniline	629
bis(2-n-Butoxyethyl) phthalate	539			4-Chloromethyl 2,2'-dimethyl pentane	623
Butyl acetate	617-618, 627-628			2-Chloro-5-methylphenol	541
n-Butyl acetate	625			4-Chloro-2-methylphenol	541
sec-Butyl acetate	617-618			4-Chloro-3-methylphenol	528, 538, 541-543, 604
n-Butylbenzene	517, 544, 546, 548, 576, 593, 597, 604, 632			2-Chloronaphthalene	538
tert-Butyl acetate	617-618			Chloroneb	508, 511, 518, 520-521, 524, 529-530
Butyl acrylate	627-628			2-Chloro-4-nitroaniline	629
n-Butyl acrylate	628			4-Chloro-2-nitroaniline	629
n-Butyl alcohol	581, 625			1-Chloro-4-nitrobenzene	614-615
s-Butyl alcohol	581			4-Chloro-3-nitrobenzotrifluoride	529-530
sec-Butyl alcohol	624-625, 649-650			Chloropentafluoroethane	527
tert-Butyl alcohol	515, 581, 625, 649-650			2-Chlorophenol	502, 528, 538, 541-543, 604, 630
Butyl aldehyde	554, 581			3-Chlorophenol	502
Butylate	518, 525			4-Chlorophenol	502, 541
Butylbenzene	515, 577, 608, 614-615			4-Chlorophenyl-phenyl ether	538
s-Butylbenzene	576, 593			Chloropicrin	509, 547
sec-Butylbenzene	515, 517, 544-546, 548, 577, 593, 604, 614-615			Chloroprene	544
tert-Butylbenzene	515, 517, 544-546, 548, 576, 593, 601, 604, 614-615			2-Chloropropane	544
Butyl benzyl phthalate	538, 539			2-Chloropropene	594
Butyl caproate	618			3-Chloropropene (allyl chloride)	553, 603, 622
Butyl cellosolve	625			3-Chloropropionitrile	515
Butyl cellosolve acetate	625			Chloropropylate	508, 529-530
Butylene	601				
Butyl ether	619				

C

Cadinene	564
δ -Cadinene	557-559, 566
γ -Cadinene	556-557
Caffeine	636, 640-641, 647
Camphene	556, 558-559, 560, 564, 567
Camphor	556, 558-559, 597
Capric acid methyl ester	561, 572, 573
Caproic acid methyl ester	561, 572, 570, 573, 607
Caprylic acid methyl ester	561, 572, 573
Captafol	508, 529-530
Captan	524, 529-530
d-Carane	557
Carbamazepine	640, 643
Carbazole	538, 540
Carbepoxide 10/11	644
Carbinoxamine	644
Carbon dioxide	549, 552, 577, 583, 601, 605
CO/air	552, 601
Carbon disulfide	517, 546, 548, 550, 553, 578, 579, 591
Carbon monoxide	582
Carbon tetrachloride	507, 509, 515, 517, 544-549, 553, 603, 622, 635, 649-651
Carbonyl sulfide (COS)	550, 577-579, 581, 591, 605, 634
Carbophenothion	524, 562
Carboxin	518, 525-526
tau-Cardinol	556
3-Carene	556
Carfentany	645
Carvacrol	560
cis-Carveol	563, 565
trans-Carveol	555, 563, 565
Carvone	565
l-Carvone	555
Carvone phenylethyl acetate	558
cis-Caryyl acetate	555, 565
Caryophyllene	556
α -Caryophyllene	566
β -Caryophyllene	555, 557-559, 564-566
Caryophyllene oxide	556
Celestolide	558
Cellosolve acetate	625
Cetearyl decanoate	558-559

Chloropyrifos.....	575
4-Chlorostyrene.....	614-615
2-Chloro-1,1,1,2-tetrafluoroethane.....	527
Chlorothalonil.....	508, 511, 518, 520-521, 529-530, 542
2-Chlorotoluene.....	515, 517, 544-546, 548, 604, 614-615
3-Chlorotoluene.....	546, 604, 614-615
4-Chlorotoluene.....	515, 517, 544, 546, 548, 604, 614-615, 623
Chlorotrifluoromethane.....	527
Chlorpheniramine.....	641, 644
Chlorpropham.....	518, 525
Chlorprothixene.....	644
Chlorpyrifos.....	524, 526, 531, 562
4-Chlortestosterone-17-acetate (clostebol).....	648
5- α -Cholestane.....	648
Cholesterol.....	648
Chrysene.....	503, 505, 518, 532, 538, 540
Chrysene-d12.....	518-519, 538, 542
1,8-Cineol.....	555, 558-559, 565
Cineole.....	560
Cinerin I.....	536
Cinerin II.....	536
Cinnamaldehyde.....	561, 563
trans-Cinnamaldehyde.....	560
Cinnamic alcohol.....	563
Cinnamyl acetate.....	563, 566
trans-Cinnamyl acetate.....	566
Cinnamyl alcohol.....	561
Cinnamyl cinnamate.....	558-559
Cinnamyl phenyl acetate.....	558-559
Cinnanzine.....	644
Citral.....	563
Citronellal.....	558, 564
Citronellic acid.....	509
Citronellol.....	557-561, 563, 569
Citronellyl acetate.....	558-559
Citronellyl formate.....	557-558
Citronellyl propionate.....	558
Citronellyl tiglate.....	558-559
cis-Citronellyl tiglate.....	563
trans-Citronellyl tiglate.....	563
Clemizole.....	644
Clobazam.....	640
Clonazam.....	636
Clonazepam.....	194, 636, 639-640
Cocaine.....	639-642, 647
Codeine.....	637, 639, 641-642, 647
Commamyl acetate.....	560
α -Copaene.....	555, 558-559, 566
Coprostane (5- β -cholestane).....	648
Cotinine.....	642
Coumarin.....	560-561, 563
2-Cresol.....	630
3-Cresol.....	630
4-Cresol.....	630
m-Cresol.....	542-543, 630
o-Cresol.....	542-543, 630
p-Cresol.....	542-543, 630
Crotonaldehyde.....	515, 612
Cumene.....	576, 593, 616, 623, 632
Cumic aldehyde.....	556
Cyanazine.....	511, 518, 520, 526
2-Cyanopyridine.....	610
3-Cyanopyridine.....	610

Cyclizine.....	644
Cycloate.....	518, 525
Cyclohexane.....	553, 576-577, 588-589, 593, 623-624, 626, 635, 639, 650-651
Cyclohexanol.....	602, 605-606
Cyclohexanone.....	612
2-Cyclohexyl-4,6-dinitrophenol.....	541
Cyclohexyl methacrylate.....	628
Cyclopenta[c,d]pyrene.....	505
Cyclopentane.....	588-589, 596, 601
Cyclopentanol.....	602, 605
Cyclopentanone.....	612
Cyclopentene.....	588-589
1-Cyclopentene.....	596
Cyclopentylbarbital.....	637
Cyclopropane.....	584-587
Cyheptamide.....	642
o-Cymene.....	556
p-Cymene.....	558, 559, 567
r-Cymene.....	555, 564-565

D

2,4-D.....	512
Dalapon.....	512, 543
2,4-DB.....	512
4,4'-DDD.....	506, 508, 511, 513, 514, 520-521, 525, 529
o,p'-DDD.....	524, 575
p,p'-DDD.....	518, 521-524, 529-530, 575
4,4'-DDE.....	506, 508, 511, 513, 514, 520-521, 525, 529
p,p'-DDE.....	518, 521-524, 529-530, 575
o,p'-DDE.....	524, 575
4,4'-DDT.....	506, 508, 511, 513, 514, 520-521, 525, 528, 529, 542
o,p'-DDT.....	524
p,p'-DDT.....	518, 521-524, 529-530, 575
Decabromobiphenyl ether (decaBDE).....	502
Decachlorobiphenyl.....	506, 508, 513-514, 520-524, 534
Decanal.....	564, 611
1,10-Decandiol.....	619
Decane.....	249, 251-252, 503, 576, 593, 599, 623
n-Decane.....	249, 260, 532, 577, 596-597
1,10-Decanediol.....	619, 620, 622
Decanoic acid.....	570, 611
Decanol.....	567
1-Decanol.....	602, 605
n-Decylamine.....	609-610
DEF.....	518
1-Dehydro-17- α -methyltestosterone (methandrostenolone).....	648
Dehydroisoandrosterone (prasterone).....	648
1-Dehydrotestosterone acetate.....	648
1-Dehydrotestosterone benzoate.....	648
1-Dehydrotestosterone (boldenone).....	648
1-Dehydrotestosterone undecylenate.....	648
Delorazepam.....	640
Demoxepam.....	636, 640
n-Des.....	644
Desalkyl aurazepam.....	636
Desipramine.....	640
Desmethyldiazepam.....	640
n-Desmethyl methsuximide.....	643
Dextromethorphan.....	637
Diacetone alcohol.....	625
Diacetyl.....	568

Diacetylmorphine (heroin).....	194, 637, 639-641, 647
Di-allate A.....	529-530
Di-allate B.....	529-530
2,4-Diaminoanisole.....	629
2,4-Diaminotoluene.....	629
3,4-Diaminotoluene.....	610
Diamyl phthalate.....	539
Diazepam.....	636, 639-642
Diazinon.....	518, 525, 531, 562
Dibenz[a,h]anthracene.....	503, 505, 518, 532, 538
1,2,3,4-Dibenzanthracene.....	540
1,2,5,6-Dibenzanthracene.....	540
Dibenzo[a,e]pyrene.....	505
Dibenzo[a,h]pyrene.....	505
Dibenzo[a,i]pyrene.....	505
Dibenzo[a,l]pyrene.....	505
Dibenzofuran.....	538, 540
Dibenzo-p-dioxin.....	540
Dibenzothiophene.....	540
Dibenzyl ether.....	558-559
Dibenzyl phthalate.....	539
Diborane.....	633
Dibromoacetic acid.....	543
Dibromoacetone.....	509, 547
1,2-Dibromobenzene.....	549
4,4'-Dibromobiphenyl.....	511
Dibromochloromethane.....	507, 509, 515, 517, 544-548, 553
1,2-Dibromo-3-chloropropane (DBCP).....	507-509, 515, 517, 529-530, 544-548, 603, 622
1,2-Dibromoethane (EDB).....	507, 509, 515, 517, 544, 546-549, 553, 603, 622
Dibromofluoromethane.....	515
Dibromomethane.....	515, 517, 544-546, 548, 622
1,2-Dibromomethane.....	545
α , α -Dibromo-m-xylene.....	529-530
2,6-Dibromo-4-nitroaniline.....	629
4,4'-Dibromooctafluorobiphenyl.....	512
2,4-Dibromophenol.....	541
1,2-Dibromopropane.....	547
2,3-Dibromopropionic acid.....	543
1,2-Dibromo-1,1,2,2-tetrafluoroethane.....	527
Dibucaine.....	643
Dibutyl chlorendate.....	529-530
tert-Dibutyl disulfide.....	592
Dibutyl ether.....	627-628
Dibutylpentyltin.....	537
Dicamba.....	512
Dichlobenil.....	526
Dichloro.....	508
Dichloroacetic acid.....	543
Dichloroacetone.....	509, 547
3,4-Dichloroaniline.....	610, 629
1,2-Dichlorobenzene.....	515, 517, 538, 544-546, 548-549, 553, 604, 614-615, 623
1,3-Dichlorobenzene.....	256, 258, 515, 517, 538, 544-546, 548-549, 553, 604, 614-615, 623
1,4-Dichlorobenzene.....	256, 258, 515, 517, 538, 544-546, 548-549, 553, 604, 614-615, 623
1,4-Dichlorobenzene-d4.....	515, 538, 542
3,3'-Dichlorobenzidine.....	538, 542, 629
3,5-Dichlorobenzoic acid.....	512
2,3-Dichlorobiphenyl.....	518

p,p'-Dichlorobiphenyl.....	524	Diethylbenzene isomer.....	593	Dimethyl phthalate.....	518, 538, 539
cis-Dichlorobutene.....	515	Diethyl disulfide.....	578	Dimethyl sulfide.....	550, 578-579, 581, 591-592, 626
cis-1,4-Dichlorobutene.....	544	Diethylene glycol.....	620-621	Dimethyl sulfoxide (DMSO).....	626-627, 639, 649-650
trans-Dichlorobutene.....	515	Diethylene glycol monobutyl ether.....	619, 620	Dimethyl tetrachloroterephthalate	
trans-1,4-Dichloro-2-butene.....	517, 548, 603, 622	Diethylene glycol monoethyl ether.....	619, 620	(DCPA).....	508, 511-512, 518, 520-521, 524, 526, 529-530
Dichlorodifluoromethane.....	515, 517, 527, 544, 546, 548, 549, 553	Diethylenetriamine.....	611	Dimethyltryptamine.....	646
Dichlorodimethyl silane.....	633	Diethyl ether.....	517, 548, 581, 600, 626, 632, 635, 649-650	Di-n-butyl phthalate.....	518, 538, 539
Dichloroethane-d4.....	515	1,2-Diethyl-4-ethylbenzene.....	596	2,4-Dinitroaniline.....	629
1,1-Dichloroethane.....	515, 517, 544-549, 553, 603, 622	1,3-Diethyl-5-ethylbenzene.....	596	2,2'-Dinitrobiphenyl.....	540
1,2-Dichloroethane.....	515, 517, 544, 546-549, 553, 603, 622, 635, 651	Diethyl formamide (DEF).....	518	2,7-Dinitrofluorene.....	540
cis-1,2-Dichloroethane.....	651	Diethyl phthalate.....	518, 538-539, 560	4,6-Dinitro-2-methyl phenol.....	538, 542
trans-1,2-Dichloroethane.....	651	Diethyl sulfide.....	578, 581	1,3-Dinitronaphthalene.....	540
1,1-Dichloroethene.....	515, 517, 544-546, 548-549, 553, 603, 622, 651	Diethyltryptamine.....	646	1,5-Dinitronaphthalene.....	540
cis-1,2-Dichloroethene.....	515, 517, 544, 546, 548-549, 553	1,4-Difluorobenzene.....	550, 517, 548-549, 553	2,4-Dinitrophenol.....	528, 538, 541-543, 604, 630
trans-1,2-Dichloroethene.....	515, 517, 544-546, 548, 553, 622	1,1-Difluoroethane.....	527	2,5-Dinitrophenol.....	541
1,1-Dichloroethylene.....	547, 635	Difolatan.....	526	2,4-Dinitrotoluene.....	518, 538, 542
cis-1,2-Dichloroethylene.....	635	Diglyme.....	619, 620	2,6-Dinitrotoluene.....	518, 538
trans-1,2-Dichloroethylene.....	635	Dihexyl phthalate.....	539	Di-n-octyl phthalate.....	538, 539
1,1-Dichloro-1-fluoroethane.....	527	9,10-Dihydroanthracene.....	540	Dinonyl phthalate.....	539
Dichlorofluoromethane.....	527, 545	trans-Dihydro carveol acetate.....	555	Dinoseb.....	512, 541
Dichloromethane.....	527, 546, 626, 633, 635, 651	Dihydro carveone.....	555, 565	1,3-Dioxalane.....	619
Dichloromethyl silane.....	633	trans-Dihydro carvyl.....	565	1,4-Dioxane.....	515, 553, 576, 593, 619, 624-626, 635, 650-651
2,6-Dichloro-4-nitroaniline.....	629	Dihydrocodeine.....	637	Diphenamid.....	518, 525-526
2,3-Dichlorophenol.....	502, 541	Dihydropentaborane.....	633	Diphenhydramine.....	644
2,4-Dichlorophenol.....	502, 528, 538, 541-543, 604, 630	Diisobutyl phthalate.....	539	Diphenylaniline.....	610
2,5-Dichlorophenol.....	502, 541	Diisopropylamine.....	610	9,10-Diphenylanthracene.....	540
2,6-Dichlorophenol.....	502, 541-543	1,3-Diisopropylbenzene.....	604, 614-615	Diphenyl isophthalate.....	539
3,4-Dichlorophenol.....	502, 541	1,4-Diisopropylbenzene.....	604, 614-615	Diphenyl oxide.....	558-559
3,5-Dichlorophenol.....	502, 541	Diisopropyl ether (DIPE).....	515, 581, 595, 600	Diphenyl phthalate.....	539
Dichloroprop.....	512	Dimenhydrinate.....	641, 644	Diphenyl pyraline.....	644
1,1-Dichloropropane.....	622	Dimethoate.....	519, 531, 562	Diphenyl sulfide.....	578
1,2-Dichloropropane.....	515, 517, 544-546, 548-549, 553, 603, 622	3,3'-Dimethoxybenzidine.....	629	Dipropylene glycol.....	619-620
1,3-Dichloropropane.....	515, 517, 544, 546, 548, 622	1,2-Dimethoxyethane.....	595, 639	Dipropyl ether.....	600
2,2-Dichloropropane.....	515, 517, 544-546, 548, 622	Dimethoxyethane (DME).....	595	Disulfoton.....	518, 525
1,3-Dichloro-2-propanol.....	515	n,n-Dimethylacetamide.....	625, 639	Diuron.....	526
1,1-Dichloro-2-propanone.....	509, 517, 547-548	Dimethylamine.....	608	cis-13,16-Docosadienoic acid methyl ester.....	572, 573
1,1-Dichloropropene.....	515, 517, 544, 548	Dimethylamphetamine.....	636	cis-4,7,10,13,16,19-Docosahexaenoic acid methyl ester.....	572, 573
1,1-Dichloro-1-propene.....	546	2,4-Dimethylaniline.....	610	Docusane.....	503
cis-1,2-Dichloropropene.....	603	2,6-Dimethylaniline.....	610, 629	Dodecahydrotriphenylene.....	540
cis-1,3-Dichloropropene.....	515, 517, 544-546, 548-549, 553, 622	7,12-Dimethylbenz[a]anthracene.....	540	γ-Dodecalactone.....	558-559
trans-1,2-Dichloropropene.....	603	1,2-Dimethylbenzene.....	560	Dodecane.....	249, 251-252, 503, 576-577, 593, 623
trans-1,3-Dichloropropene.....	515, 517, 544-546, 548-549, 553, 622	n,n-Dimethylbenzylamine.....	610	n-Dodecane.....	249, 260, 582, 596-597
Dichlorotetrafluoroethane.....	553	2,2-Dimethylbutane.....	588-589, 596, 599	Dodecanoic acid.....	570
1,2-Dichloro-1,1,2,2-tetrafluoroethane.....	527, 549	2,3-Dimethylbutane.....	588-589, 596, 599	Dodecanol.....	558-559
2,2-Dichloro-1,1,1-trifluoroethane.....	527	Dimethyl disulfide.....	578	Dodecenal.....	564
Dichlorvos.....	518, 525	Dimethyl ether (DME).....	581, 601	n-Dotriacontane.....	597
Dicyclohexylamine.....	609-610	Dimethylformamide (DMF).....	624, 626-627, 650	Doxylamine.....	641, 644
Dicyclohexyl phthalate.....	539	n,n-Dimethylformamide.....	625-626, 639	Droperidol.....	646
Dieldrin.....	506, 508, 511, 513-514, 518, 520-525, 529-530, 575	2,6-Dimethylhept-5-enal (melonal).....	558, 559	Dursban (chlorpyrifos).....	518, 519
Diethanolamine (DEA).....	606, 607	2,2-Dimethylhexane.....	517, 548, 599		
Diethylamine.....	610	1,3-Dimethyl-2-imidazolidinone (DMI).....	639		
n,n-Diethylaniline.....	610	2,6-Dimethylnaphthalene.....	540		
2,6-Diethylaniline.....	610	1,2-Dimethyl-2-nitrobenzene.....	519		
Diethylbenzene.....	593	1,3-Dimethyl-2-nitrobenzene.....	518, 542		
1,2-Diethylbenzene.....	576, 597	2,2-Dimethylpentane.....	596, 599		
1,3-Diethylbenzene.....	576, 596	2,3-Dimethylpentane.....	588-589, 599		
1,4-Diethylbenzene.....	597	2,4-Dimethylpentane.....	588-589, 596, 599		
		Dimethylpentyltin.....	537		
		3,6-Dimethylphenanthrene.....	540		
		2,3-Dimethylphenol.....	541		
		2,4-Dimethylphenol.....	528, 538, 541-543, 604		
		2,5-Dimethylphenol.....	541		
		2,6-Dimethylphenol.....	541		
		3,4-Dimethylphenol.....	541		

E

cis-11,14-Eicosadienoic acid methyl ester.....	572, 573
Eicosane.....	503, 558
n-Eicosane.....	597
cis-5,8,11,14,17-Eicosapentaenoic acid methyl ester.....	572, 573
cis-11,14,17-Eicosatrienoic acid methyl ester.....	572, 573
cis-8,11,14-Eicosatrienoic acid methyl ester.....	572, 573
cis-11-Eicosenoic acid methyl ester.....	572, 573
Elaidic acid methyl ester.....	572, 573

Endosulfan I.....	506, 508, 511, 513-514, 518, 520-525, 529-530	Ethyl cellosolve.....	625	Farnesyl acetate.....	566
Endosulfan II.....	506, 508, 511, 513-514, 518, 520-525, 529-530	Ethyl decanoate.....	558-559	Fenamiphos.....	518-519, 525
α -Endosulfan.....	575	Ethyl dodecanoate.....	558-559	Fenarimol.....	518, 525
β -Endosulfan.....	575	Ethylene.....	577, 580, 584-585, 588-590, 601, 634	Fenchone.....	563
Endosulfan sulfate.....	506, 508, 511, 513-514, 518, 520-525, 529-530, 575	Ethylene/acetylene.....	591	Fenchyl acetate.....	558
Endrin.....	506, 508, 511, 513-514, 518, 520-525, 529-530, 542	Ethylenediamine.....	611	Fenitrothion.....	562, 575
Endrin aldehyde.....	506, 508, 511, 513-514, 518, 520-525, 528-530	Ethylene.....	595	Fentanyl.....	640, 645-646
Endrin ketone.....	506, 508, 513-514, 518, 521-524, 529-530	Ethylene glycol.....	602, 605, 619-622, 649-650	Fenthion.....	562
Ephedrine.....	636	Ethylene glycol monoacetate.....	607	Fenvalerate.....	519
Epiandrosterone (trans-androsterone).....	648	Ethylene glycol monobutyl ether.....	619-620	Florazone.....	558-559
Epichlorohydrin.....	515, 619	Ethylene glycol monoethyl ether.....	619-620	Flunitrazepam.....	194, 636, 639-640
EPTC.....	518, 525-526	Ethylene glycol monoformate.....	607	Fluoranthene.....	503, 532, 538, 540
Eptam.....	527	Ethylene glycol monomethyl ether.....	619-620	Fluorene.....	503, 518, 532, 538, 540, 542
Erucic acid methyl ester.....	572, 573	Ethylene glycol phenyl ether.....	622	Fluorobenzene.....	515, 517, 544, 545, 548, 614-615
Erythritol.....	569	Ethylene oxide.....	594, 605, 607, 631	2-Fluorobiphenyl.....	538
Esfenvalerate.....	519	Ethyl ether.....	515, 619, 624-625, 650	2-Fluorophenol.....	538
Estazolam.....	636	Ethyl formate.....	568, 607, 617-618, 624-625, 649-650	Flurazepam.....	640, 642
17- α -Estradiol.....	648	Ethyl heptanoate.....	558-559	Fluridone.....	518, 525-526
β -Estradiol.....	648	Ethyl hexadecanoate.....	558-559	Fonofos.....	526
Estriol.....	648	Ethyl hexanoate.....	558-559	Formaldehyde.....	550, 553, 613
Estrone.....	648	2-Ethyl hexanoic acid.....	509	Formaldehyde-DNPH.....	613
Ethane.....	580, 583-585, 588-591, 596, 601, 634	2-Ethyl-1-hexanol.....	602, 605	Formaldehyde-PFBHA.....	614
Ethane/ethylene.....	580	Ethyl isovalerate.....	558-559	Formic acid.....	509, 570, 607
Ethanethiol.....	578	Ethyl mercaptan.....	579, 581, 591-592	Frambinone (raspberry ketone).....	558-559
Ethanol.....	515, 554, 568, 581, 595-596, 600, 602, 605, 623-628, 635, 638-639, 649-650	Ethyl methacrylate.....	515, 517, 548, 628	Fucitol.....	568
Ethanolamine.....	608	Ethyl methyl sulfide.....	578	Furan.....	619
Ethchlorvynol.....	646	Ethyl morphine.....	637	Furfural.....	558-559, 612
Ethinamate.....	646	Ethyl nonanoate.....	558-559	Furfuryl alcohol.....	602, 605
Ethoprop.....	518, 525-526	Ethyl octadecanoate.....	558-559		
Ethosuximide.....	643	Ethyl octanoate.....	558-559		
Ethotoin.....	644	Ethyl parathion.....	575		
2-Ethoxyethanol.....	602, 605, 607, 635, 639	Ethyl pentadecanoate.....	559		
bis(2-Ethoxyethyl) phthalate.....	539	Ethyl pentadecanoate nonadecane.....	558		
2-Ethoxyethyl acetate.....	607, 617-618	Ethyl pentanoate.....	558-559		
2-Ethyl-1-hexanol.....	567	α -Ethylphenethyl alcohol.....	602, 605		
bis(2-Ethylhexyl) adipate.....	518	β -Ethylphenethyl alcohol.....	602, 605		
bis(2-Ethylhexyl) phthalate.....	518, 538-539	2-Ethylphenol.....	630		
Ethyl acetate.....	515, 553-554, 558-559, 568-569, 617-618, 624-628, 635, 639, 649-650	4-Ethylphenol.....	630		
Ethyl acetylene.....	634	Ethyl-3-phenyl oxiran carboxylate.....	569		
Ethyl acrylate.....	617-618, 626-628	Ethyl propanoate.....	554		
Ethyl alcohol.....	607	Ethyl propionate.....	558-559, 617-618, 627-628		
2-(Ethylamino)-ethanol.....	606	Ethyl tetradecanoate.....	558-559		
Ethyl and dimethyl thiophenes.....	593	4-Ethyltoluene.....	553		
Ethylbenzene.....	501, 515, 517, 544-546, 548-549, 553, 576-577, 588-589, 593-594, 596-597, 599, 601, 604, 614-616, 623, 625, 632, 639, 649-651	m-Ethyltoluene.....	576, 593		
Ethylbenzene-d10.....	597	p-Ethyltoluene.....	576, 593		
Ethyl benzoate.....	569, 617-618	Ethyl undecanoate.....	558-559		
Ethyl butanoate.....	554	Ethyl valerate.....	618		
2-Ethyl-1-butanol.....	606	Ethyl vinyl ether.....	619		
Ethyl t-butyl ether (ETBE).....	515, 546, 581, 595, 600, 625,	Etridiazole.....	508, 511, 520-521, 524		
Ethyl butyrate.....	558-559, 569, 618	Eucalyptol.....	556-557		
Ethyl caprate.....	554	Eugenol.....	561, 563		
Ethyl caproate.....	618	iso-Eugenol.....	561		
Ethyl caprylate.....	554	Evernyl.....	558		
		Eugenyl acetate.....	558-559		
		Eugenyl methyl ether.....	559		

Heptachlor epoxide	506, 508, 511, 513-514, 518, 520-525, 529-530, 542
2,2',3,3',4,4',6-Heptachlorobiphenyl	518
Heptadec-1-ene	558
Heptadecane	536
n-Heptadecane	597
Heptadecanoic acid methyl ester	572, 573
cis-10-Heptadecenoic acid methyl ester	572, 573
Heptanal	554, 611-612
1,7-Heptandiol	619, 620
Heptane	576, 593, 601, 616, 623, 625, 639, 649-650
n-Heptane	553, 577, 588-589, 596-597, 624, 635, 650
Heptanoic acid	570, 607, 611
n-Heptanoic acid	509
1-Heptanol	532, 602, 605-606
2-Heptanol	602, 605
3-Heptanol	602, 605
2-Heptanone	612
3-Heptanone	612
4-Heptanone	612
cis-4-Hepten-1-ol	602, 605
trans-2-Hepten-1-ol	602, 605
Hercolyn D (tetrahydro & dihydro methyl abietate)	558
Heroin	194, 637, 639-641, 647
Hexabromobenzene (HBB)	524
Hexabromobiphenyl	519
Hexacosane	503
Hexachlorobenzene	508, 518, 520, 524, 529-530, 538
Hexachlorobezene	511, 521
2,2',4,4',5,6'-Hexachlorobiphenyl	518
Hexachlorobutadiene	515, 517, 538, 544-546, 548, 553, 603, 622
Hexachloro-1,3-butadiene	549
β-Hexachlorocyclohexane	575
δ-Hexachlorocyclohexane	575
γ-Hexachlorocyclohexane	575
Hexachlorocyclopentadiene	508, 511, 518, 538, 542, 603, 622
Hexachloroethane	515, 517, 538, 548, 603, 622
Hexachloropentadiene	529-530
Hexadecane	503
n-Hexadecane	249, 597
Hexadecanoic acid	570
2,4-Hexadienal	557
Hexanal	554, 611-612
1,6-Hexandiol	619
Hexane	517, 548, 599, 601, 623-626, 639, 649-650
n-Hexane	553, 577, 583-585, 588-589, 596-597, 624, 635, 650
1,6-Hexanediol	620
Hexanoic acid	611
Hexanol	249, 251-252, 558-559
1-Hexanol	554, 602, 606
2-Hexanol	602, 605
3-Hexanol	602, 605
2-Hexanone	517, 548, 553, 612
3-Hexanone	612
n-Hexatriacontane	597
Hexazinone	518-519, 525-527
Hexchlorocyclopentadiene	520
cis-2-Hexen-1-ol	602, 605
cis-3-Hexen-1-ol	602, 605
cis-2-Hexene	588-589
trans-2-Hexene	588-589
Hexobarbital	637, 641
Hexyl acetate	558-559
Hexyl butyrate	556
Hexyl cinnamaldehyde	563
Hexyl cinnamic aldehyde	561
Hexylene glycol	558-559
Hexyl 2-ethylhexyl phthalate	539
1-Hexyl mercaptan	592
n-Hexyl mercaptan	592
Hexyl methacrylate	628
α-Humulene	564
Hydrocarbon	616
Hydrocodone	639, 641
Hydrogen oxide	577
Hydrogen sulfide	550, 577-579, 581, 583, 591, 592, 601, 605, 634
Hydroxy acetate	607
Hydroxy citronellal	558-559, 561, 563
4-Hydroxy-4-methyl-2-pentanone	602, 605
Hydroxyphenamate	641
Hydroxypropionitrile	515
Hydroxypropyl acrylate	628
Hydroxyzine	644
Ibogaine	640, 646
Ibuprofen	641, 647
Imipramine	642
Indan	596-597
Indeno[1,2,3-c,d]pyrene	503, 505, 518, 532, 538
Inositol	568
Iodobenzene	614-615
Iodoform	603, 622
Iodomethane	517, 544, 548, 603, 622
Irgafos	633
Irganox	633
Isoamyl acetate	554, 569, 617-618
Isoamyl alcohol	554, 568, 649-650
Isoamyl butyrate	569
Isoamyl mercaptan	592
Isoamyl salicylate	558, 560
Isoborneol	558-560, 567
Isobutane	584-589, 599
Isobutanol	515, 568, 581, 595, 602, 605, 626-628
Iso-butene	634
Isobutyl acetate	554, 558-559, 617-618, 624, 627-628
Isobutyl acrylate	627-628
Isobutyl alcohol	554
Isobutyl aldehyde	554, 581
Isobutylbenzene	577, 604, 614-615
Isobutylene	584-587
Isobutyl mercaptan	592
Isobutyl propionate	627-628
Isobutyraldehyde	582, 600, 612, 625
Isobutyric acid	509, 570, 607-608
Isocaproic acid	570, 607
Isodrin	508, 524, 529-530
Isoeugenol	563
trans-Isoeugenol	563
Isomenthol	557
Isomenthone	555, 557-559
d-Isomethone	565
Isonox	633
Isocetane	577, 588-589, 596-597, 623-624, 635
Isopentane	577, 584-589, 596
Isopentanoic acid	509
Isophorone	518, 538, 542
Isoprene	586-589
Isopropanol	554, 581, 602, 605, 624-625, 638, 650
Isopropyl acetate	617-618, 624, 639
Isopropyl acrylate	627-628
Isopropyl alcohol	515, 553, 624-625, 649-650
Isopropyl amine	649-650
Isopropyl benzene	515, 517, 544-546, 548, 577, 588, 589, 593, 604, 614-616
Isopropyl ether	619
Isopropyl mercaptan	592
Isopropyl myristate	558-559
2-Isopropylphenol	630
Isopropyltoluene	515
4-Isopropyltoluene	545-546
p-Isopropyltoluene	517, 544, 548
Isovaleraldehyde	581, 600
Isovaleric acid	570, 607
J	
Jasmodin I	536
Jasmodin II	536
cis-Jasmone	555, 565
K	
Kelthane	508, 529-530
Kepone	519
Ketamine	646, 647
Krypton	582
L	
L-9-Carboxy-11-nor-δ9-THC	649
L-11-Hydroxy-δ9-THC	649
Lactose	569
Lactulose	569
Lauric acid methyl ester	572, 573
Lavandulol	556
Lavandulyl acetate	556
Levulinic acid	509
Lidocaine	640, 643, 647
Lignoceric acid methyl ester	572, 573
Lilial	561, 563
Limonene	555, 558-560, 563-564, 565, 567
δ-Limonene	556-557, 564
S-(-)-Limonene	567
Linalool	555, 558-561, 563-567
β-Linalool	556
Linalool acetate	556
cis-Linalool oxide	558-559
trans-Linalool oxide	558-559
Linalyl acetate	558-560
Linoleic acid methyl ester	572, 573
Linolelaidic acid methyl ester	572, 573
Linolenic acid methyl ester	572, 573
γ-Linolenic acid methyl ester	572, 573
Lofentanyll	645
Lorazepam	636, 639-640
Lormetazepam	636
Lyrall	561
Lyrall 1	563
Lyrall 2	563
Lysergic acid diethylamide (LSD)	640, 646

M

β-Maaliene.....	557	3-Methyl-2-buten-1-ol.....	602, 605	Methyl 3-hydroxydodecanoate.....	571
Malathion.....	519, 531, 562, 575	2-Methyl-1-butene.....	596	Methyl 2-hydroxyhexadecanoate.....	571
Manitol.....	569	2-Methyl-2-butene.....	584, 588-589, 596	Methyl 2-hydroxytetradecanoate.....	571
Mannitol.....	568	3-Methyl-1-butene.....	586-589	Methyl 3-hydroxytetradecanoate.....	571
Meclizine.....	644	trans-2-Methyl-2-butenic acid.....	608	Methyl isobutyl ketone.....	624-625
Medazepam.....	636, 640	2-Methylbutyl acetate.....	617-618	Methyl laurate.....	571
Menthofuran.....	557	2-Methylbutyl alcohol.....	558-559	Methyl mercaptan.....	550, 579, 581, 591-592
Menthol (MeOH).....	557, 560, 565, 577, 595	3-Methylbutyl alcohol.....	558-559	Methyl methacrylate.....	515, 517, 548, 627-628
Menthone.....	557-559	Methyl t-butyl ether (MTBE).....	515, 595, 639, 649-650	Methyl 14-methylhexadecanoate.....	571
Menthyl acetate.....	565	Methyl butyrate.....	558-559	Methyl 12-methyltetradecanoate.....	571
Meperidine.....	639, 642	2-Methyl butyric acid.....	608	Methyl monoethanolamine (MMEA).....	606
Mephobarbital.....	637	Methyl chavicol.....	561	Methyl myristate.....	571
Mepivacaine.....	643	3-Methylcholanthrene.....	540	1-Methylnaphthalene.....	540, 597
Meprobamate.....	641, 646	2-Methyl-4-chlorophenoxyacetic acid (MCPA).....	512	2-Methylnaphthalene.....	538, 540, 596-597
Merphos.....	518, 525	Methyl chlorophenoxypropionic acid (MCPP).....	512	2-Methyl-5-nitroaniline.....	629
Mescaline.....	646	5-Methylchrysene.....	505	Methyl nonadecanoate.....	571
Mesitylene.....	623	Methyl-cresol.....	558	Methyl octadecanoate.....	618
Mesterolone.....	648	Methyl-p-cresol.....	559	Methyl octine carbonate.....	561, 563
Methacrolein.....	612, 626-627	Methyl-r-cresol.....	566	Methyl oleate.....	571
Methacrylic acid.....	509	Methylcyclohexane.....	588-589, 599, 635, 639, 651	Methyl palmitate.....	571
Methacrylonitrile.....	515, 517, 548, 626-627	Methylcyclopentane.....	577, 588-589, 596, 599	Methyl palmitoleate.....	571
Methadone.....	639-640, 642	1-Methyl-1-cyclopentene.....	596	Methyl paraoxon.....	518
Methamidophos.....	531, 562	Methyl decanoate.....	618	Methyl parathion.....	575
Methamphetamine.....	636, 639, 642	Methyl-2,4-dichlorophenylacetate.....	512	Methyl pentadecanoate.....	571
Methane.....	552, 577, 580, 582-585, 588-591, 601, 626, 634-635	Methyl diethanolamine (MDEA).....	606	2-Methylpentane.....	588-589, 596, 599
Methanethiol.....	578	4-Methyl-2,5-dimethoxyamphetamine (STP).....	636	3-Methylpentane.....	501, 588-589, 596, 599, 623
Methanol.....	554, 568, 581-582, 595-596, 600-602, 605-606, 613, 623-624, 626-628, 635, 638-639, 649-651	2-Methyl-4,6-dinitrophenol.....	528, 541-543, 604	2-Methyl-2 pentanol.....	581
Methapyriline.....	641, 644	Methyl dodecanoate.....	618	3-Methyl-3-pentanol.....	606
Methaqualone.....	640, 646	Methyl eicosenoate.....	618	4-Methyl-2-pentanol.....	602, 605-606
Methidathion.....	531, 562	Methyl elaidate.....	571	2-Methyl-3-pentanone.....	612
Methofuran.....	565	4,4'-Methylenedianiline.....	629	4-Methyl-2-pentanone (MIBK).....	517, 548, 553, 612
Methohexital.....	637	3,4-Methyl enedioxyamphetamine.....	642	2-Methyl-1-pentene.....	588-589, 596
Methone.....	565	3,4-Methyl enedioxyethylamphetamine.....	642	4-Methyl-1-pentene.....	588-589
Methoprene.....	536	3,4-Methyl enedioxyethylamphetamine.....	642	bis(4-Methyl-2-pentyl) phthalate.....	539
Methoxychlor.....	506, 508, 511, 513-514, 518, 520-524, 529-530	Methyl ethyl ketone (2-butanone).....	649, 650	α-Methylphenyl alcohol.....	602
2-Methoxyethanol.....	602, 605, 626	Methyl isobutyl ketone (MIBK).....	639	2-Methylphenol.....	538, 541
bis(4-methoxyethyl) phthalate.....	539	α-Methyl ionone.....	560	3-Methylphenol.....	541
2-Methoxy-5-methylaniline.....	629	Methyl cis-9,10-methylene octadecanoate.....	571	4-Methylphenol.....	538, 541
1-Methoxy-2-propanol.....	639	Methyl cis-9,10-methyl hexadecanoate.....	571	5-Methyl-5-phenylhydantoin.....	643
Methsuximide.....	643	Methylene chloride.....	515, 517, 544-545, 548-549, 553, 603, 622, 624-625, 639, 649-650	2-Methyl-1-propanethiol.....	578-579, 581, 591
Methyl acetate.....	617-618, 627-628, 635	Methylenedioxyamphetamine (MDA).....	636, 639	2-Methyl-2-propanethiol.....	578-579, 591
Methyl acetylene.....	591	Methylenedioxyethylamphetamine (MDE).....	636, 639	2-Methyl-1-propanol (isobutanol).....	600
Methyl acrylate.....	515, 517, 548, 627-628	Methylenedioxyethylamphetamine (MDMA).....	636, 639	2-Methyl-2-propanol (tert-butanol).....	600, 639
Methyl alcohol.....	625	Methyl ephedrine.....	636	Methyl propionate.....	617-618, 627-628
Methylamine.....	608	1-Methyl-2-ethylbenzene.....	597	1-Methyl-1-propnaethiol.....	579, 591
4-Methylaminorex.....	646	1-Methyl-3-ethylbenzene.....	597	1-Methyl-3-propylbenzene.....	596
r-Methylansiole.....	566	1-Methyl-4-ethylbenzene.....	597	1-Methyl-1-propyl mercaptan.....	592
2-Methylanthracene.....	540	Methyl ethyl ketone (MEK).....	581, 624-625, 650	2-Methyl-1propyl mercaptan.....	592
9-Methylanthracene.....	540	Methyl ethyl sulfide.....	581	2-Methyl-2-propyl mercaptan.....	592
Methyl arachidate.....	571	Methyl eugenol.....	561	1-Methyl-2-pyrrolidine.....	610
Methyl benzoate.....	558-559, 566, 617-618	2-Methylfluoranthene.....	540	1-Methyl-2-pyrrolidinone.....	626-627
5-Methylbenzo[b]thiophene.....	578	Methyl formate.....	617-618	n-Methylpyrrolidone.....	639
3-Methylbenzothiophene.....	578	Methyl γ ionone.....	561	Methyl stearate.....	571
2-Methyl butanal.....	611	3-0-Methylglucose 1.....	569	α-Methylstyrene.....	576, 593, 614-616
2-Methylbutane.....	599	3-0-Methylglucose 2.....	569	4-Methylstyrene.....	614-615
2-Methyl-1-butanol.....	568, 602, 605, 611	Methyl heptadecanoate.....	571	Methyl sulfide.....	634
2-Methyl-2-butanol.....	602, 605	2-Methylheptane.....	577, 588-589, 596	Methyl tert-butyl ether (MTBE).....	517, 546, 548, 581-582, 597, 600, 619, 623
3-Methyl-1-butanol.....	568, 602, 605, 625, 635	3-Methylheptane.....	588-589, 596	17α-Methyltestosterone.....	648
3-Methyl-2-butanone.....	612, 635	4-Methylheptane.....	577	Methyl tetradecanoate.....	618
2-Methyl-3-buten-2-ol.....	602, 605	Methyl heptene carbonate.....	561, 563	2-Methylthiophene.....	578
		Methyl hexadecanoate.....	618	3-Methylthiophene.....	578
		2-Methylhexane.....	588-589, 599	Methyl thiophenes.....	593
		3-Methylhexane.....	588-589, 599	Methyl tridecanoate.....	571
		Methyl 2-hydroxydecanoate.....	571	Methyltripentyltin.....	537
		Methyl 2-hydroxydodecanoate.....	571		

Methyl undecanoate.....	571
Metolachlor.....	511, 518, 520, 525-526
Metribuzin.....	511, 518, 520, 525-526
Mevinphos.....	518, 525, 531, 542
MGK-264.....	518, 525
Endo-MGK.....	264, 536
Exo-MGK.....	264, 536
Mirex.....	508, 519, 524, 529-530
Molinate.....	518, 525-526
6-Monoacetylmorphine.....	637
Mono ethylene glycol.....	606
Mono-ethanolamine (MEA).....	607
Monuron.....	526
Morphine.....	637, 641, 647
Musk T (ethylene brassylate).....	558-559
Musk ketone.....	560
Musk xylene.....	560
Myrcene.....	555, 558-560, 564-565
β -Myrcene.....	556-557
Myristic acid methyl ester.....	572, 573
Myristoleic acid methyl ester.....	572, 573

N

Naled.....	531
Naphthalene.....	249, 251-252, 501, 503, 515, 517, 532, 538, 540, 542, 544-546, 548, 577, 596-597, 623
Naphthalene-d8.....	538, 597
1-Naphthol.....	630
Napropamide.....	518, 525-526
Nefopam.....	643
Neomenthol.....	560, 565
Neon.....	552, 582
Neral.....	558-559, 563-564
Nerol.....	557
Nerol acetate.....	556
cis-Nerolidol.....	558
trans-Nerolidol.....	558
Nervonic acid methyl ester.....	572, 573
Neryl acetate.....	558-559, 564
Nicotinamine.....	636
Nicotine.....	639-641, 647
Nitrazepam.....	194, 636, 639
2-Nitroaniline.....	538, 610, 629
3-Nitroaniline.....	538, 610, 629
4-Nitroaniline.....	538, 610, 629
Nitrobenzene.....	515, 517, 538, 548, 614-615
Nitrobenzene-d5.....	538
2-Nitrobiphenyl.....	540
3-Nitrobiphenyl.....	540
4-Nitrobiphenyl.....	540
Nitrofen.....	508
Nitrogen.....	552, 582, 605, 632
Nitromethane.....	635
1-Nitronaphthalene.....	540
2-Nitronaphthalene.....	540
Nitrophen.....	519
2-Nitrophenol.....	528, 538, 541-543, 604, 630
3-Nitrophenol.....	541
4-Nitrophenol.....	512, 528, 538, 541-543, 604
2-Nitropropane.....	515, 517, 548
n-Nitrosodimethylamine.....	538, 542
n-Nitrosodiphenylamine.....	538
n-Nitroso-di-n-propylamine.....	538
2-Nitrotoluene.....	604, 614-615

3-Nitrotoluene.....	604, 614-615
4-Nitrotoluene.....	604, 614-615
trans-Nonachlor.....	508, 518, 529-530
Nonadecane.....	559
Nonadec-1-ene.....	558-559
Nonanal.....	564
γ -Nonalactone.....	558-559
1,9-Nonandiol.....	619
Nonane.....	576, 593, 599, 601, 623
n-Nonane.....	532, 577, 588-589, 596-597
1,9-Nonanediol.....	620, 622
Nonanol.....	567
1-Nonanol.....	602, 605
Nonyl aldehyde.....	612
Nonylamine.....	608
n-Nonylamine.....	609-610
Nootkatone.....	564
Norcodeine.....	637
Nordazepam.....	636
Norethandrolone.....	648
Norflurazon.....	518-519, 525-526
Normorphine.....	637
19-Nortestosterone (nandrolone).....	648
19-Nortestosterone-17-decanoate.....	648
19-Nortestosterone-17-propionate.....	648

O

cis-Ocimene.....	555, 565
trans-Ocimene.....	555
β -cis-Ocimene.....	556
β -trans-Ocimene.....	556
Octacosane.....	503
n-Octacosane.....	597
2,2',3,3',4,5',6,6'-Octachlorobiphenyl.....	518
1,2,3,4,6,7,8,9-Octachlorodibenzodioxin.....	504
1,2,3,4,6,7,8,9-Octachlorodibenzofurans.....	504
Octadecane.....	503, 536
n-Octadecane.....	597
Octadecanoic acid.....	570
Octanal.....	554, 564, 611
1,8-Octandiol.....	619
Octane.....	576, 593, 599, 601, 623
isooctane.....	501
n-Octane.....	532, 577, 588-589, 596-597
1,8-Octanediol.....	620
Octanoic acid.....	570, 611
Octanol.....	558-559, 564, 567
1-Octanol.....	602, 605
3-Octanol.....	555, 565
3-Octanone.....	532, 556
1-Octene.....	532
1-Octen-3-ol.....	556, 565
Octen-1-ol acetate.....	556
Octyl acetate.....	558-559
3-Octyl acetate.....	555, 565
Octyl aldehyde.....	612
n-Octylamine.....	609-610
n-Octyl mercaptan.....	592
Oleic acid methyl ester.....	572, 573
Omethoate.....	562
Ordram.....	527
Oxadiazon.....	527
Oxazepam.....	636, 639-640
Oxazopam.....	636
Oxychlorane.....	519

Oxycodone.....	194, 639
Oxydemeton-methyl.....	531
4,4'-Oxydianiline.....	629
Oxygen.....	552, 582
Oxymetholone.....	648
Oxymorphone.....	641

P

Paarlan.....	527
Palmitic acid methyl ester.....	572, 573
Palmitoleic acid methyl ester.....	572, 573
Paraldehyde.....	515
Parathion.....	519, 562
Parathion-methyl.....	562
Pebulate.....	525-526
Pentaborane.....	633
2,2',3,4,4'-Pentabromodiphenyl ether.....	502
2,2',4,4',5-Pentabromodiphenyl ether.....	502
2,2',4,4',6-Pentabromodiphenyl ether.....	502
Pentachlorobenzene.....	614-615
2,2',3',4,6-Pentachlorobiphenyl.....	518
1,2,3,7,8-Pentachlorodibenzofuran.....	504
2,3,4,7,8-Pentachlorodibenzofuran.....	504
Pentachloroethane.....	515, 517, 548, 603, 622
Pentachloronitrobenzene.....	508, 511, 524, 529-530
Pentachlorophenol.....	512, 518, 528, 538, 541-543, 604
n-Pentadecane.....	249, 582, 597
Pentadecanoic acid methyl ester.....	572, 573
cis-10-Pentadecenoic acid methyl ester.....	572, 573
cis-1,3-Pentadiene.....	586-587
trans-1,3-Pentadiene.....	586-587
Pentafluorobenzene.....	515, 517, 548
Pentafluoroethane.....	527
Pentamethylbenzene.....	597
Pentamethyl disiloxane.....	633
Pentanal.....	611
1,5-Pentandiol.....	619
Pentane.....	599, 601, 603, 617, 625-626
iso-Pentane.....	583, 590-591, 601
n-Pentane.....	577, 583-591, 596-597
1,2-Pentandiol.....	532
1,5-Pentandiol.....	620
2,3-Pentanedione (VDK).....	554
2,3-Pentanedione (acetyl propionyl).....	558-559
n-Pentanoic acid.....	509
1-Pentanol.....	554, 602, 605-606, 625, 639
2-Pentanol.....	602, 605
3-Pentanol.....	554, 602, 605
2-Pentanone.....	515, 612, 624
3-Pentanone.....	612
Pentene-1.....	584, 586-587
1-Pentene.....	584-585, 588-589, 596, 625
cis-2-Pentene.....	584, 586-589
trans-2-Pentene.....	584, 586-589
4-Pentenoic acid.....	608
trans-2-Pentenoic acid.....	608
trans-3-Pentenoic acid.....	608
1-Penten-3-ol.....	602, 605
2-Penten-1-ol.....	602, 605
Pentobarbital.....	637, 641
Pentyl ether.....	619
1-Pentyl mercaptan.....	592
cis-Permethrin.....	508, 511, 518, 520-521, 524, 529-530

trans-Permethrin.....508, 511, 518, 520-521,
524, 529-530
Perthane508, 529-530
Perylene.....540
Perylene-d12.....518-519, 538, 542
 α -Phellandrene564
Phenacetin636, 641
Phenanthrene503, 518, 532, 538, 540
Phenanthrene-d10.....518, 538, 542
Phencyclidine (PCP).....639, 640, 642, 646
Phenethylalcohol567
Pheniramine.....644
Phenmetrazine641
Phenobarbital637, 640-641, 643-644
Phenol249, 251-252, 502, 528,
538, 541-543, 604, 630
Phenol-d5538
2-Phenoxyethanol.....602, 605
Phensuximide643
Phentermine636, 639
Phenyl acetaldehyde.....561, 563
Phenylacetone.....636
Phenylacetylene.....576, 593
9-Phenylanthracene540
Phenyl ephedrine.....636
1-Phenyl ethyl acetate.....554
Phenylethylamine610
Phenyl ethyl cinnamate558-559
Phenylethylmalonamide.....643
Phenyl ethyl tiglate558-559
2-Phenyl-naphthalene540
Phenyltoin.....644
Phenyltoloxamine644
Phenytoloin.....643
Phorate526
Phosmet.....531
Phosphine633
Phylethylmalonamide (PEMA)644
Picloram512
4-Picoline532
 α -Pinene555-559, 564-567
 β -Pinene555, 557-560, 564-565, 567
Piperazine611
Piperitone565
Pirimiphos-methyl.....562
Polyethylene glycol621
Polymethoxyflavones.....564
Prazepam636, 640
Prebane (terbutryn).....518
Primidone643-644
Procaine.....641, 643, 647
Progesterone648
Promazine.....644
Promethazine.....644
Prometon518, 525-526
Prometryn.....518-519, 525-526
Promofos575
Pronamide518, 525-526
Propachlor.....508, 511, 518, 520-521, 524,
526-527, 529-530
Propadiene584-587, 634
Propanal.....553
2-Propanamine625
1,3-Propanediol.....619
1,2-Propanediol (propylene glycol).....619
Propane.....580, 583-591, 596, 599, 601, 617

1,2-Propanediol (propylene glycol).....620
1,3-Propanediol620
1-Propanethiol578-579, 591
2-Propanethiol578-579, 581, 591
Propanil526
Propanoic acid.....611
Propanol550
1-Propanol515, 554, 595, 600, 602, 605-606,
625-626, 635, 638-639, 649-650
2-Propanol595, 626, 635, 639
n-Propanol568, 581, 624, 650
Propargyl alcohol515
Propazine.....518-519, 525-527
Propene.....553
2-Propen-1-ol (allyl alcohol).....602, 605
Propiolactone.....515
Propiomazine.....646
Propionaldehyde581, 600, 612
Propionic acid.....509, 568, 570, 607
1-Propionic acid532
Propionitrile.....515, 517, 548, 610, 626-627
Proposyphene.....642
Propyl acetate.....554, 617-618, 635
n-Propyl acetate.....625
Propylbenzene515, 576-577, 593, 614-616, 623
n-Propylbenzene517, 532, 544-546, 548, 588-589,
593, 596-597, 604, 616
Propyl benzoate617-618
Propyl butyrate618
Propyl caproate.....618
Propylene.....577, 579-580, 584-591, 601, 617
Propylene glycol.....602, 605
Propyl ether.....581, 619
n-Propyl ether.....625
Propyl formate.....617-618
1-Propyl mercaptan.....592
2-Propyl mercaptan.....592
n-Propyl mercaptan.....592
2-Propyl pentanoic acid.....509
Propyl propionate.....617-618
Propyl valerate.....618
Propyne580, 584-587, 588-589
2-Propyn-1-ol (propargyl alcohol).....602, 605
Prowl.....527
Psilocin.....646
Pulegone.....565
Pyrazole.....610
Pyrene.....503, 518, 532, 538, 540
Pyrene-d10518
Pyrethrin I536
Pyrethrin II.....536
Pyridazine610
Pyridine.....515, 538, 610, 624-627, 635, 650
Pyrilamine.....644
Pyrimidine610
Pirimiphos ethyl575
Pirimiphos methyl.....575
Pyrrithyldione.....646
Pyrrbutamine641

Q

Quinine.....647

R

Rhamnitol568-569
Rhamnose 1569
Rhamnose 2569
Ribitol.....568-569
Ribose 2569
Ro-Neet.....527
Rosatol (rosetone).....559
Rosatol (rosetone) geranyl butyrate558
Rose oxide, cis-rose.....558-559
Rose oxide, trans-rose558-559

S

Sabinene.....555, 558, 564-565
trans-Sabinene hydrate.....555, 565
Salicylamide643
 α -Santolene556
Secbumeton526
Secobarbital637, 641
Selenide633
Sencor.....527
Simazine.....511, 518, 520, 525-527, 542
Simetryn.....518, 525-526
 α -Sinensal.....564
 β -Sinensal564
Sorbitol569
Stanozolol.....648
Stearic acid methyl ester.....572, 573
Stigmasterol.....648
Stirifos525
Strawberry aldehyde.....569
Strychnine639, 647
Styrene515, 517, 544, 546, 548-549, 553, 576,
588-589, 593, 604, 614-617, 623, 632
Styrene dimer617
Styrene oxide.....614, 615
Styrene trimer.....617
Sucrose569
Sufentanyl645
Sulfallate526
Sulfur dioxide.....605, 634
Sulpiride646
Sutan527

T

Talbutal637, 646
Tebuthiuron518, 525
Temazepam194, 636, 639-640
Tenocyclidine (TCP).....646
Terbacil518, 525-527
Terbufos518, 525-526, 542
Terbufos sulfone519
Terbuthylazine.....526
Terbutryn.....525-526
Terphenyl-d14538
Terpinene565
 α -Terpinene555, 564-565
 γ -Terpinene555, 557-559, 564-565, 567
Terpinen-4-ol.....555-556, 558-559, 564-565, 567
 α -Terpineol556, 558-560, 564-565, 567
Terpinolene555, 559, 564-565
 γ -Terpinolene565
Terrazole.....518, 529-530

Testosterone	648	o-Tolualdehyde	612	Triflupromazine	644
2,2',4,4'-Tetrabromodiphenyl ether	502	p-Tolualdehyde	612	Trifluralin	508, 511, 518, 520-521, 524, 529-530, 542
2,3',4,4'-Tetrabromodiphenyl ether	502	Toluene	249, 251-252, 501, 515, 517, 544-546, 548-549, 553, 576, 577, 586-589, 593-594, 596-597, 599, 601, 604, 614-617, 623-625, 632, 635, 639, 649-651	Triglyme (triethylene glycol dimethyl ether)	619, 620
2,3',4',6'-Tetrabromodiphenyl ether	502			Trimethylamine	608
Tetraborane	633			2,4,5-Trimethylaniline	629
Tetrabutyltin	537			1,2,3-Trimethylbenzene	
Tetracaine	643	Toluene-d8	515, 549	(hemimellitene)	546, 577, 596-597, 599, 601, 614-615
1,2,4,5-Tetrachlorobenzene	614-615	o-Toluidine	629	1,2,4-Trimethylbenzene	
2,2',4,4'-Tetrachlorobiphenyl	518	Tonalid	558-560	(pseudocumene)	501, 515, 517, 544-546, 548-549, 553, 577, 588-589, 596-597, 614-615, 623
2,3,7,8-Tetrachlorodibenzodioxin (TCDD)	504	Trazodone	639	1,3,5-Trimethylbenzene	
2,3,7,8-Tetrachlorodibenzofuran (TCDF)	504	Treflan	527	(mesitylene)	515, 517, 544-546, 548-549, 553, 588-589, 596-597, 614-615
Tetrachloroethane	507, 509	2,2',4-TriBDE (BDE-17)	502	Trimethylchloro silane	633
1,1,1,2-Tetrachloroethane	515, 517, 544-546, 548, 603, 622	2,4,4'-TriBDE (BDE-28)	502	3,5,5-Trimethylhexanol	558-559
1,1,2,2-Tetrachloroethane	515, 517, 544, 546, 548-549, 553, 603, 622	Triacetanamine (TAA)	595	2,2,4-Trimethylpentane (isooctane)	553
Tetrachloroethene	515, 517, 544-549, 553, 622	Triadimefon	518, 525	2,3,3-Trimethylpentane	577, 596
Tetrachloroethylene	547	Triazolam	636, 640	2,3,4-Trimethylpentane	588-589
Tetrachloro-m-xylene (TCMS)	506, 508, 513-514, 520, 524, 529-530, 534	Tribromoacetic acid	543	2,3,5-Trimethylphenol	630
2,4,5,6-Tetrachloro-m-xylene	521-523	2,4,6-Tribromophenol	538	Trimethyl phosphate	532
2,3,4,5-Tetrachlorophenol	502, 541	Tributylpentyltin	537	Trimipramine	640
2,3,4,6-Tetrachlorophenol	502, 541	Trichloroacetic acid	543	2,4,6-Trinitrophenol	604
2,3,5,6-Tetrachlorophenol	502, 541-543	Trichloroacetonitrile	509, 547	1,3,5-Trioxane impurity	554
Tetrachlorovinphos (stirifos)	518	2,4,5-Trichloroaniline	629	Tripelennamine	644
n-Tetracontane	597	1,2,3-Trichlorobenzene	515, 517, 544-546, 548, 614-615	Triphenyl phosphate	518, 519, 562
Tetracosane	503	1,2,4-Trichlorobenzene	515, 517, 538, 544-546, 548-549, 553, 614-615	Triprolidine	644
n-Tetracosane	597	1,3,5-Trichlorobenzene	546, 614-615		
Tetradecane	503, 596	2,4,5-Trichlorobiphenyl	518		
n-Tetradecane	249, 582, 597	Trichloroethane	507, 509, 624		
Tetradecanoic acid	570	1,1-Trichloroethane	649-650		
Tetraethylenepentamine	611	1,1,1-Trichloroethane	507, 509, 515, 517, 544-549, 553, 603, 622, 635, 651		
1,1,1,2-Tetrafluoroethane	527	1,1,2-Trichloroethane	507, 509, 515, 517, 544-546, 548-549, 553, 603, 622		
1,1,2,2-Tetrafluoroethane	527	Trichloroethene	515, 517, 544-549, 553, 622		
Tetrahydrocannabinol	639	Trichloroethylene	547, 624, 626, 635, 650		
Δ9-tetrahydrocannabinol (THC)	649	Trichlorofluoromethane	515, 517, 527, 544-546, 548, 549, 553		
1,2,3,4-Tetrahydrofluoranthene	540	1,1,2-trichloro-1,2,2-trifluoroethane	549		
Tetrahydrofuran (THF)	517, 548, 553, 619, 624, 626, 639, 650, 651	Trichloromethane	625-626		
Tetrahydropyran	619	Trichloromethyl silane	633		
1,2,3,5-Tetramethylbenzene	596-597	2,3,4-Trichlorophenol	502, 541-543		
1,2,4,5-Tetramethylbenzene	577, 596-597	2,3,5-Trichlorophenol	502, 541-543		
Tetrapentyltin	537	2,3,6-Trichlorophenol	502, 541-543		
Thenylidiamine	641	2,4,5-Trichlorophenol	502, 538, 541-543		
Thiamine pyrophosphate (TPP)	531	2,4,6-Trichlorophenol	502, 528, 538, 541-543, 630		
Thiamylal	637	3,4,5-Trichlorophenol	502, 541		
Thiazopyr	519	2,4,5-Trichlorophenoxyacetic acid	512		
4,4'-Thiodianiline	629	1,2,3-Trichloropropane	507, 509, 515, 517, 543-544, 546, 548, 603, 622		
Thiopental	637	1,1,1-Trichloropropanone	547		
Thiophene	578-579, 581, 591, 592, 593	1,1,1-Trichloro-2-propanone	509		
Thioridazine	644	1,1,2-Trichloro-1,2,2-trifluoroethane	527, 553		
Thonzylamine	644	1,1,2-Trichlorotrifluoroethane	603, 622		
Threhalose	569	Tricosanoic acid methyl ester	572, 573		
Threitol	569	Tricyclazole	518, 525-526		
Thujone	563	Tridecane	576-577, 593		
α-Thujone	558, 564	n-Tridecane	249, 582		
β-Thujone	564	Tridecanoic acid methyl ester	572, 573		
Tillam	518, 527	Triethanolamine (TEA)	607, 608		
Tinuvin P	633	Triethylamine	610, 626-627		
α-Tocopherol	645	Triethylene glycol	621		
β-Tocopherol	645	Triethylenetetramine	611		
δ-Tocopherol	645	1,1,1-Trifluoroethane	527		
γ-Tocopherol	645	Trifluoromethane	527		
Tolban	527	Trifluorotoluene	544		
o-Tolidine	629				
m-Tolualdehyde	612				

U

Undecane	576, 593, 599, 623
n-Undecane	249, 260, 577, 596-597
Undecanoic acid methyl ester	572, 573

V

Valencene	564
Valeraldehyde (pentanal)	581, 600, 612
Valeric acid (pentanoic acid)	570, 607, 608, 611
iso-Valeric acid	611
Valerolactone	608
Vanillin	563
Verapamil	639
Vernam	527
Vernolate	518, 525-526
Vertenex	558-559
Vinclozolin	519, 575
Vinyl acetate	553, 617-618
Vinyl chloride	515, 517, 544-546, 548-549, 553
Vinyl-trimethyl silane	633
Viridiflorol	555, 565

W

Water	583, 595, 607
-------	---------------

X

Xenon	582
m-Xylene	515, 517, 544-546, 548-549, 553, 576-577, 588-589, 593-594, 596-597, 599, 601, 604, 614-616, 623, 625, 632
m,p-Xylene	501, 639, 649-651

o-Xylene	501, 515, 517, 544-546, 548-549, 553, 576-577, 588-589, 593-594, 596-597, 599, 601, 604, 614-616, 623, 625, 632, 639, 649-651
p-Xylene	515, 517, 544-546, 548-549, 553, 576-577, 588-589, 593-594, 596-597, 599, 601, 604, 614-616, 623, 625, 632
2,3-Xylenol	630
2,4-Xylenol	630
2,5-Xylenol	630
2,6-Xylenol	630
3,4-Xylenol	630
3,5-Xylenol	630
Xylose 1	569
Xylose 2	569

Ordering Information

Easy Ordering Terms and Conditions

Discounts and Delivery

Agilent Technologies specializes in fast delivery. In the US, if you call before 2 PM EST, we will ship your order that day. You may also request overnight express delivery before 6 PM EST and you will have your order the next day. Volume discounts on a variety of individual products are offered when the entire quantity is shipped to one address at one time.

A shipping and handling fee will be added to your order unless the purchase is over \$2000 US for orders place online or over \$4000 for orders place via phone. Special shipping (i.e., overnight in the US) is available in most regions at an additional cost.

Agilent is required to collect all state and local sales taxes unless the buyer's tax-exempt certificate is on file with Agilent Technologies. Please be prepared to provide a copy if it is not on file, when placing your order.

Please check with your Agilent Customer Service Representative, local Authorized Distributor, or the Agilent website for current prices, special offers, promotions and discounts when placing your order.

Satisfaction Guaranteed

If you are not satisfied with your Agilent product within the first 60 days, you may return your purchase in its original condition for a full refund or credit. A return policy statement is included in every Agilent shipment and posted under Product Information on the website. In the US and Canada, please call for a Return Authorization form and return instructions at **1-800-227-9770**. If your Agilent product was purchased from a distributor, please contact the distributor.

Shipping Damages

If items are damages in transit, please follow the instructions below:

- If a shipment is visibly damaged upon arrival, do not accept it until the person making the delivery has endorsed the bill of lading with statement for the extent of the damage.
- If any damage is found after unpacking, retain all cartons and inner packaging and immediately request an inspection from the carrier.
- Notify the Agilent Customer Contact Center at **1-800-227-9770** about the damaged shipment so that we can make the appropriate sales adjustment and/or provide you with return instructions (Sales order number, product number and quantity damaged will be needed).

Easy Ways To Order

- Phone: **1-800-227-9770** (option 1, 1) in the US and Canada — Mon-Fri, 8AM to 8PM EST
- Fax: **1-302-633-8901** in the US
- Email: **cag_sales-na@agilent.com** in the US and Canada
- Online: **www.agilent.com/chem** in the US and Canada

Payment Options

- In the US, Visa, MasterCard, Discover and American Express are accepted with a minimum order of \$20 (not applicable in all countries).
- Email ePay@agilent.com to make an electronic payment using the ACH/EFT (Automated Clearing House/Electronic Funds Transfer) method.
- Establish a charge account through your Agilent Customer Service Representative or Your Local Agilent sales office. An account number will be assigned to you for charging your purchases. Payment terms are net 30 days from the invoice date. All orders are subject to credit approval.

We will be happy to supply a price quote via, phone, email or fax if you need it in writing.

Warranties

All Agilent Technologies products in this catalog are designed and manufactured to stringent standards under the Agilent quality system registered to ISO 9001. At Agilent, we back every product with a 90-day warranty and a money-back guarantee. If Agilent receives notice of defects during the warranty period. Agilent shall, at its option, either repair or replace products which prove to be defective. If Agilent is unable, within a reasonable time, to repair or replace any product to a condition as warranted, the buyer shall be entitled to a refund of the purchase price upon return of the product to Agilent. The warranty period for each product begins on the day of shipment.

This warranty shall not apply to any defect, failure, or damage caused by improper use or improper or inadequate maintenance or care. This warranty is exclusive and no other warranty, whether written or oral, is expressed or implied. Agilent specifically disclaims the implied warranties of merchantability and fitness for particular purposes. The remedies provided herein are the buyer's sole and exclusive remedies. In no event shall Agilent be liable for direct, indirect, special, incidental, or consequential damages (including loss of profits) whether based on contract, tort, or any other legal theory.

Agilent Technologies Order Form

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.



Agilent Technologies

Order Date		Purchase Order or Credit Card Number & Expiration		Taxable Y or N?	Yes	No
Name:				If No, please provide Certificate #		
Title:						
Phone:				Fax:		
Company:				Email:		
		Shipping Address		Billing Address		
Company:						
Street:						
Room/Bldg/Dept:						
City:						
State/Province/Country:						
Zip/Postal Code:						
Deliver to:						
Part Number	Description			Quantity	Price	Total Cost
Special Instructions:					Subtotal:	
					Tax:	
					Total:	
For Assistance: Please call the Agilent Technologies Customer Contact Center at 1-800-227-9770 (U.S. and Canada) Email Node: cag_sales-na@agilent.com Fax Number: 302-633-8901 U.S. and Canada Note: 1. All pricing, tax, discount, and availability information is subject to verification by Agilent Technologies. 2. Shipping and handling is free for orders over \$2000 US.						

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.

WWW.AGILENT.COM/CHEM/GC

Agilent Technologies Order Form

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.



Agilent Technologies

Order Date		Purchase Order or Credit Card Number & Expiration		Taxable Y or N?	Yes	No
Name:				If No, please provide Certificate #		
Title:						
Phone:				Fax:		
Company:				Email:		
		Shipping Address		Billing Address		
Company:						
Street:						
Room/Bldg/Dept:						
City:						
State/Province/Country:						
Zip/Postal Code:						
Deliver to:						
Part Number	Description		Quantity	Price	Total Cost	
Special Instructions:				Subtotal:		
				Tax:		
				Total:		
For Assistance: Please call the Agilent Technologies Customer Contact Center at 1-800-227-9770 (U.S. and Canada) Email Node: cag_sales-na@agilent.com Fax Number: 302-633-8901 U.S. and Canada Note: 1. All pricing, tax, discount, and availability information is subject to verification by Agilent Technologies. 2. Shipping and handling is free for orders over \$2000 US.						

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.

WWW.AGILENT.COM/CHEM/GC

Agilent Technologies Order Form

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.



Agilent Technologies

Order Date		Purchase Order or Credit Card Number & Expiration		Taxable Y or N?	Yes	No
Name:				If No, please provide Certificate #		
Title:						
Phone:				Fax:		
Company:				Email:		
		Shipping Address		Billing Address		
Company:						
Street:						
Room/Bldg/Dept:						
City:						
State/Province/Country:						
Zip/Postal Code:						
Deliver to:						
Part Number	Description		Quantity	Price	Total Cost	
Special Instructions:				Subtotal:		
				Tax:		
				Total:		
For Assistance: Please call the Agilent Technologies Customer Contact Center at 1-800-227-9770 (U.S. and Canada) Email Node: cag_sales-na@agilent.com Fax Number: 302-633-8901 U.S. and Canada Note: 1. All pricing, tax, discount, and availability information is subject to verification by Agilent Technologies. 2. Shipping and handling is free for orders over \$2000 US.						

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.

WWW.AGILENT.COM/CHEM/GC

Agilent Technologies Order Form

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.



Agilent Technologies

Order Date		Purchase Order or Credit Card Number & Expiration		Taxable Y or N?	Yes	No
Name:				If No, please provide Certificate #		
Title:						
Phone:				Fax:		
Company:				Email:		
		Shipping Address		Billing Address		
Company:						
Street:						
Room/Bldg/Dept:						
City:						
State/Province/Country:						
Zip/Postal Code:						
Deliver to:						
Part Number	Description			Quantity	Price	Total Cost
Special Instructions:					Subtotal:	
					Tax:	
					Total:	
For Assistance: Please call the Agilent Technologies Customer Contact Center at 1-800-227-9770 (U.S. and Canada) Email Node: cag_sales-na@agilent.com Fax Number: 302-633-8901 U.S. and Canada Note: 1. All pricing, tax, discount, and availability information is subject to verification by Agilent Technologies. 2. Shipping and handling is free for orders over \$2000 US.						

Outside the U.S. and Canada, please contact your local Agilent office or Authorized Distributor when ordering.

WWW.AGILENT.COM/CHEM/GC















ENSURING AN INERT GC FLOW PATH HAS NEVER BEEN MORE CRITICAL

For labs that need to perform trace level analysis on very active compounds, **Agilent Inert Flow Path solutions** ensure a reliably inert flow path for higher sensitivity, accuracy, and reproducibility. These analyses can now be done with the utmost confidence because each lot of Agilent Inert Flow Path components are qualified using chromatographic QC tests to provide reliable, consistent inertness that spans the entire flow path.

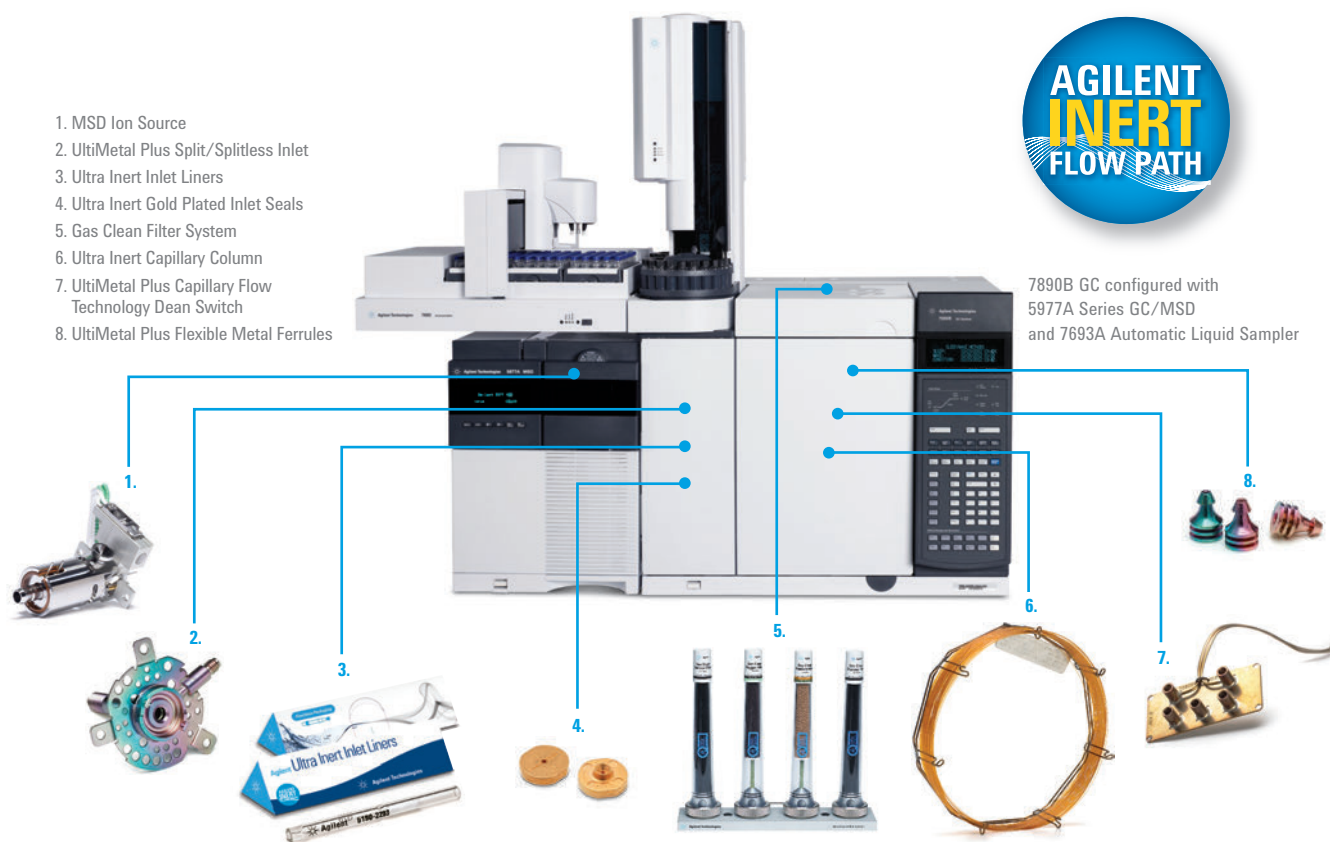

Learn more about creating an inert flow path, and order the poster at agilent.com/chem/inert

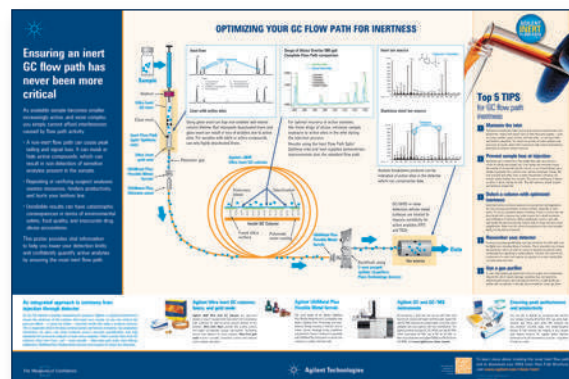

20
15|16

For more information

Buy online:

www.agilent.com/chem/store

Contact us:

www.agilent.com/chem/contactus

This information is subject to change without notice.

© Agilent Technologies, Inc. 2014

Published in December 2014

5991-5213EN



Agilent Technologies

